

# **Service Manual**

## **Agilent 7890 Series Gas Chromatograph**



**Agilent Technologies**

# Notices

© Agilent Technologies, Inc. 2013

No part of this manual may be reproduced in any form or by any means (including electronic storage and retrieval or translation into a foreign language) without prior agreement and written consent from Agilent Technologies, Inc. as governed by United States and international copyright laws.

## Manual Part Number

G3430-90006

## Edition

Ninth edition, August 2013

Eighth edition, May 2013

Seventh edition, November 2011

Printed in USA

Agilent Technologies, Inc.

2850 Centerville Road

Wilmington, DE 19808-1610 USA

## Warranty

The material contained in this document is provided "as is," and is subject to being changed, without notice, in future editions. Further, to the maximum extent permitted by applicable law, Agilent disclaims all warranties, either express or implied, with regard to this manual and any information contained herein, including but not limited to the implied warranties of merchantability and fitness for a particular purpose. Agilent shall not be liable for errors or for incidental or consequential damages in connection with the furnishing, use, or performance of this document or of any information contained herein. Should Agilent and the user have a separate written agreement with warranty terms covering the material in this document that conflict with these terms, the warranty terms in the separate agreement shall control.

## Safety Notices

### **WARNING**

**A warning calls attention to a condition or possible situation that could cause injury to the user.**

### **CAUTION**

A caution calls attention to a condition or possible situation that could damage or destroy the product or the user's work.

# Contents

## 1 Before Servicing the Instrument

Safety	20
Hand tools	20
Specialized tools for GC support	24
Electronic tools for GC support	26
Pneumatic tools for GC support	28
For Checking Inlet Pressure	30

## 2 Inlets

Inlet overview	34
Split/Splitless Inlet	35
Replacing the split/splitless inlet	35
Replacing the split/splitless inlet top weldment assembly	38
Replacing the split/splitless inlet split vent line	39
Replacing the split/splitless inlet split vent filter canister and valve assembly	40
Reinstalling the split/splitless inlet	42
Multimode Inlet	43
Replacing the multimode inlet	44
Replacing the multimode inlet top weldment assembly	46
Replacing the multimode inlet split vent line	47
Replacing the multimode inlet split vent filter canister and valve assembly	48
Replacing a Multimode CO2 cryo valve	50
Replacing the Multimode N2 cryo valve	52
Removing the MMI cryo/PC board bracket	54
Replacing the MMI heater/cryo PC board	56
Cool On-Column Inlet	58
Replacing the COC inlet	58
Replacing the COC heater/sensor assembly	60
Replacing an inlet cryo cooling tube	62
Purged Packed Inlet	65
Replacing the Purged Packed inlet	65
Replacing the Purged Packed inlet heater	70
Programmed Temperature Vaporization Inlet (PTV)	73
Replacing the PTV inlet	74
Replacing the PTV inlet cryo jacket	77

Replacing the PTV CO2 or LN2 inlet cryo assembly	79
Replacing a PTV CO2 cryo valve	80
Replacing the PTV N2 cryo valve	81
Removing the PTV cryo/PC board bracket	83
Replacing the PTV heater/cryo PC board	85
Volatiles Interface	87
Replacing or cleaning the volatiles interface	87
Replacing the volatiles interface heater/sensor assembly	90
Calibrating the interface	92
Solvent Vapor Exit Accessory	93
Replacing the SVE valve/fitting assembly	94
Replacing the SVE bleed restrictor column	96
Replacing the SVE tri-column assembly	97
Replacing the SVE pre-column assembly	98

### 3 Detectors

Detector overview	100
About the detector signal boards	100
Accessing detector signal boards	103
Replacing detector signal boards	104
Detector signal board, bad checksum error	107
Flame Ionization Detector (FID)	113
Selecting an FID jet	113
Replacing the FID interconnect assembly or spring	115
To replace the interconnect spring only	117
To replace the interconnect (and spring)	117
Replacing the heater	118
Replacing the FID signal board	122
Replacing the FID electrometer	123
Replacing the FID glow plug (ignitor)	125
Replacing the entire FID	127
Ultrasonic cleaning the FID	130
Nitrogen Phosphorus Detector (NPD)	132
Selecting an NPD jet	132
Removing the NPD electrometer	134
Removing the NPD electrometer	134
Replacing the entire NPD	136
Replacing an NPD signal board	138
Replacing the NPD heater	140
Replacing the NPD electrometer	142



Cleaning the NPD jet and collector	144
Cleaning the NPD	146
Thermal Conductivity Detector (TCD)	148
Replacing the TCD detector cell	148
Replacing the TCD switching valve	150
Micro-cell Electron Capture Detector (uECD)	152
Replacing the uECD heater/sensor assembly	152
Replacing the uECD makeup gas adapter	154
Frequency test	156
Leak test	158
Performing a radioactivity leak test (wipe test)	159
Flame Photometric Detector (FPD) 7890A	160
Preparing for maintenance	160
Disassembling the FPD	161
Rebuilding the FPD optics assembly	171
Reassembling the FPD	178
Replacing the FPD	186
Bake out and run checkout test	188
Replacing the heater/sensor assemblies	189
Replacing the FPD signal board	194

## 4 EPC modules

EPC overview	196
Repairing EPC modules	197
Accessing EPC modules	199
Replacing the EPC modules	200
Replacing EPC module proportional valves	203
Inlet modules	205
Inlet modules	206
Proportional valves	207
Detector modules	208
Detector modules	209
Proportional valves	210
Pneumatics Control Module (PCM)	211
Calibrating the PCM interface	212
Changing the PCM pressure channel frit	213
Correcting PCM leaks	215
Potential PCM leak areas	216
Aux EPCs	217

Frits for auxiliary pressure control channels	218
Procedure for replacing frit	220
Updating PID Constants	221
Updating PID constants	222
No positive feedback on reconnect	226
Updating PID constants available for use by Instrument Utilities	227

## 5 Performance Verification

About Chromatographic Checkout	229
To Prepare for Chromatographic Checkout	230
To Check FID Performance	232
Typical values	237
To Check NPD Performance	238
Typical values	243
To Check TCD Performance	244
Typical values	249
To Check uECD Performance	250
Typical values	255
To Check FPD Performance (Sample 5188-5953)	256
Preparation	256
Phosphorus performance	257
Typical values	261
Sulfur performance	262
Typical values	263
To Check FPD Performance (Sample 5188-5245, Japan)	264
Preparation	264
Phosphorus performance	265
Typical values	269
Sulfur performance	270
Typical values	271

## 6 Configuration

Configuration overview	273
Inlet example	274
Detector example	275
PCM example	276
AUX example	277
Heater assignments	278
GC modules and the communications buss	280

To configure an MSD transfer line	283
To configure a valve box	286
MSD transfer line heater error	286
Initial steps	286
If that doesn't work	287
Replacing a GC module	288
Removing a GC module	289
Changing the GC configuration	290
Agilent data system methods	290
Methods stored in GC	290
Configuration locks	291
Installing new devices	292
Configuring time	294

## 7 External Cabling

Connectors on the back of the GC	296
Remote Start/Stop	297
APG Remote Control	297
Electrical specifications	297
Suggested Drive Circuits	298
The Remote Connector	298
Signal Descriptions	299
Timing Diagram	299
Cable pinouts, remote start/stop, general use	300
Cable pinouts, GC to 3395B/3396C Integrator	301
Cable pinouts, GC to 35900C, D, E/MSD/Sampler, 2 meters	302
Cable pinouts, GC to Headspace Sampler	303
Cable pinouts, APG Remote Y-cable	304
External event	305
Automatic sampler for GC	306
Cable pinouts, GC to G2613A Injector	306
Cable pinouts, GC to G4513A and 7650 Injector	307
Cable pinouts, GC to G2614A Tray	308
Cable pinouts, GC to G4514A Tray	309
Analog signal outputs	310
Analog cable: GC to 3395A/B or 3396B/C Integrators and 35900 C/D/E	
Analog to Digital Interface instrument	310
Analog cable: general use	311

Analog signal inputs	312
General purpose AIB cable	312
XCD Chemiluminescence detector AIB cable	312

## 8 Internal Cabling

Internal cabling overview	314
MSD transfer line heater cable	315
Third detector heated zone cable	316
Communication harness	317
Power connection table	318
Data communication table	319
Motor Harness & AC Control Assembly	320
Keyboard and display connection board (original style only)	322
AC power inlet assemblies	323
Keyboard and display connection board (original style only)	324
AC power inlet assemblies	325
Inlet/Detector harness	326
Auxiliary zone/Valve box harness	328
Third detector EPC communication cable	330
Keyboard/Display harness	331
Keyboard interconnect board	332
Heated zone extension cable	333
Ignitor cable, FID	334
Ignitor cable, FPD	335
RS-232 cable, ALS controller	336
PTV thermocouple cable	337
NPD power cable	339

## 9 Mainframe

Mainframe overview	341
Covers and Fans	342
Removing and replacing the covers	343
Detector top cover	343
The pneumatics top cover	343
Replacing the pneumatics top cover buttons	343
The electronics top cover	343
The side covers	343

The rear covers	344
The side cover on the third detector assembly	344
Replacing the inlet cooling fan	345
Replacing the EPC cooling fan	346
Replacing the oven bezel	348
Replacing the oven door	349
Replacing the oven door button	350
Replacing the keypad assembly and display board	351
Replacing the keypad assembly (original style)	351
Replacing only the display board (original style)	354
Replacing the keypad assembly (new style)	356
Replacing only the display board (new style)	357
Replacing the original keypad assembly with the new assembly	359
Oven temperature control	362
Oven temperature troubleshooting	363
Oven does not heat up	363
Oven not cooling or never cools	364
Testing resistance of the heater coil	365
Cryo valve installation/replacement	366
Installing a new cryo valve	366
Replacing an existing cryo valve	368
Replacing the oven sensor	371
Replacing the oven door sensor	372
Replacing the oven shroud, oven fan, and oven fan motor	375
Replacing the oven flapper assembly	381
Replacing the flapper or flapper motor	383
Replacing Components Inside the Electronics Panel	385
Accessing the analog and power board	385
Replacing the ALS board	386
Replacing the fuse on the ALS controller board	389
Troubleshooting	390
Removing the valve bracket	391
Replacing the logic board	392
Check GC configuration	395
Replacing the analog and power board	396
Replacing Components Inside the Lower Rear Metal Cover	400
Replacing AC board fuse sets	401
Replacing the AC board	402

Replacing the oven triac	404
Replacing the power transformer	406

## 10 Valves

Valves	409
Introduction	410
Valco W-series minivalves	411
Valve bodies	411
Valve rotors	412
Preload assembly (not shown)	412
Valco W-series minivalve	413
Valve bodies	413
Valve rotors	414
General purpose valves	415
Gas sample valves	415
Gas sample loops	415
Adjustable restrictor valves	415
Liquid sample valves	416
Valve Box	417
Installing the valve box	417
Removing the valve box assembly	422
Actuators	423
Installing the actuators	423
Installing the actuator on a heated valve box	426
Installing the actuator on a side mount bracket	429
Valve actuator alignment	430
Valve actuator drivers	431
Solenoid differences	431
Which solenoid valve is installed?	431
Replacing an original solenoid valve	431
Installing a new valve stack assembly	432
Original style	432
Assemble the valve driver block	432
Install the bracket and cabling (original style)	433
Installing the valve driver block	434
New style	436
Replacing a solenoid valve (new style)	436
Installing the valve driver (solenoid) bracket, cabling, and elbow fitting (new style)	441
Typical Valve Configurations	442
Option numbers	443
Custom Plumbing (diagram required), Option 700 or 730	443

Gas Sampling Option, Option 701 or 731	443
Column Isolation, Option 702 or 732	444
Two Stream Selection (Requires Gas Sampling), Option 703 or 733	444
Backflush to Detector, Option 704 or 734	445
Backflush a Precolumn to Vent, Option 705 or 735	445
Column Selection (Unused Column Isolated), Option 706 or 736	446
Sequence Reverse, Option 707 or 737	446
Sequence Reverse with Backflush of Column 1, Option 708 or 738	447
Custom Plumbing (Diagram Required), Option 800 or 830	447
Gas Sampling with Backflush of Precolumn to Vent, Option 801 or 831	448
Gas Sampling with Backflush to Detector, Option 802 or 832	448
Gas Sampling of Alternate Streams, Option 803 or 833	449
Gas Sampling with Sequence Reverse, Option 804 or 834	449
Gas Sampling with Sequence Reverse and Backflush of Column 1, Option 805 or 835	450
Column Selection with Backflush to Vent, Option 806 or 836	450
Liquid Sampling, Options 850, 852, 854, and 856	451
Troubleshooting	452
Chromatographic symptoms	452
Loss of sensitivity or excessive drift	453
Loss of peaks in specific areas of the chromatogram	453
Extraneous peaks	453
Peak broadening and tailing	454
Baseline shifts	454
Baseline upsets	454
Variation in peak area and retention time	455
Pressure check	455

## 11 Capillary flow technology (CFT) devices

CFT overview	458
Ultimate union	460
Tools and kits	460
Metal ferrules and nuts	461
Assembling nut, ferrule and swaging nut	461
Connecting and disconnecting columns	464
Mounting the union	464

## 12 Electrical

Power options	466
Regular oven	466
Fast oven	467
Converting the power option	469
Procedure	469
Remote start/stop connection	470
Remote control	471
Measuring inlet and detector heater and PRT resistance	472
Sensor resistance	474
Temperature sensor resistance	475
Oven heater coil resistance	476
Setting the instrument power configuration	477
Line voltage configuration plug	478
Ceramic and glass fuses	479
AC power board schematic	480
Oven shroud	481

## 13 Firmware

Firmware overview	483
To Update GC Firmware	484
Troubleshooting problems	486
Firmware history	487
7890 A.01.10.2	488
7890 A.01.10.1	489
7890 A.01.10	490
7890 A.01.09.2	490
7890 A.01.09.1	491
7890 A.01.09	491
7890 A.01.08	492
7890 A.01.07	492
7890A.01.04	492
7890A.01.03	494
7890A.01.02.581	494
Confirm firmware update	495
Problems and solutions	496
Possible solutions to unusual behaviors	496



Clearing battery-backed PF-RAM (7890A only)	497
Partial clear	497
Total clear	497
Set the default method	500
Default method listing	501

## 14 Firmware Instructions for Channel Partner Devices

Firmware overview for channel partners	508
Tuning a small heated zone	509
Thermal Aux Configuration	510
Heater and cryo connection locations	510
Connections assigned to a thermal zone	511
Auxiliary Type Selection	512
User Configurable Heater	513
Custom heater parameters	513
Programmable temperature zones	514
Cryo control model	515
Configuring the thermal Aux zone (Firmware A.01.07)	516
Configuring the thermal Aux zone (Firmware A.01.08 and A.01.09)	519
Valve Box Configuration	521
Heater connection locations	522
Custom heater parameters	523
Custom Heater Conf.	523
View Custom Heater Conf.	523
The PID and Mass fields	523
The Power field	523
Configuring the Valve Box Zone	524
Inlet configuration	525
Hardware configuration	526
EPC module installation	527
Heater and Cryo installation	528
The Inlet configuration display	529
Gas type	529
Cryo Parameters	529
Cryo type	529
Cryo	529
Use Cryo temperature	529
Cryo timeout	529
Cryo fault	529
Custom heater parameters	530
Configured hardware	531

Configuring a custom inlet (Firmware A.01.05)	532
Configuring a custom inlet (Firmware A.01.07)	534
Configuration procedure	537
Configuring a custom inlet (Firmware A.01.08 and A.01.09)	538
Cryo focus valve configuration	540
Operation	540
Assigned connections	540
Cryo trap configuration	541
Detector configuration	543
FID Signal 1 Detector Hardware Connections	544
AIB Signal 1 Detector Hardware Connections	546
Aux Detector Hardware Connections	548
EPC module installation	548
Detector Board Installation	549
Heater installation	549
Configuring a custom detector	550

## 15 Service Mode Diagnostics

Service mode overview	553
Detector diagnostics	554
Signal board diagnostics	555
Multiplexed ADC	556
Pneumatics	557
Power diagnostics	558

## 16 Accessories and options

## 17 Printed circuit boards

Summary of PCB Part Numbers	568
AC power board	569
Analog and power board	570
Logic and communications board (7890A)	572
Logic and communications board (7890B)	573
FID signal board	574
TCD signal board	575
uECD signal board	576
FPD signal board	577

NPD signal board	578
FID electrometer board	579
uECD electrometer board	580
NPD electrometer board	581
Analog input board	582


## 18 Illustrated Parts Breakdown

Inlets	584
Split/Splitless Inlet (SSL)	585
Split/Splitless Inlet, Top	586
SSL Split Vent Trap	587
Split Splitless Inlet Body	588
Split Splitless Inlet Warmer	590
Multimode Inlet (MMI)	591
Multimode Inlet Body	594
MMI Top Insert Assembly	596
MMI Split Vent Trap	598
7890A MMI Temperature controller	599
7890B MMI Temperature controller	600
7890A MMI Cryo assemblies	601
7890B MMI Cryo assemblies	602
Purged Packed Inlet (PP)	603
Purged Packed Inlet Upper Body	604
Purged Packed Inlet Lower Parts	606
Cool on-column inlet (COC)	608
Cool on-column inlet body	609
Programmed Temperature Vaporization Inlet (PTV)	611
PTV Inlet Body	613
PTV Septumless head	615
PTV Split vent	616
PTV Temperature controller	617
PTV CO2 cryo assembly	618
PTV LN2 cryo assembly, 7890A	619
PTV LN2 cryo assembly, 7890B	620
Volatiles Interface (VI)	622
Volatiles Interface (VI)	623
Volatiles Interface Trap	625
Detectors	627
Flame Ionization Detector (FID)	628
FID detector	629
FID Base Assembly	631

PM kits for the FID	633
Flame Photometric Detector (FPD)	634
FPD inert transfer line parts	636
FPD ignitor and heat shield assembly	638
PMT and bracket assemblies	640
FPD lens assembly	643
FPD Covers	644
FPD <sup>+</sup> (High Temperature FPD)	646
FPD <sup>+</sup> Assemblies	648
Microcell Electron Capture Detector (uECD)	654
Microcell Electron Capture Detector (uECD) parts	655
uECD Side Mount (3rd Detector) Parts	657
Nitrogen Phosphorus Detector (NPD)	659
Nitrogen Phosphorus Detector (NPD) parts	660
Thermal Conductivity Detector (TCD)	662
Thermal Conductivity Detector (TCD) parts	665
EPC Module Parts	669
Lookup table	669
Inlet modules	670
Detector modules	671
PCM modules	672
Aux EPCs	673
Nickel Catalyst Accessory	674
Oven	676
Oven assembly	677
Oven flapper assembly	680
CO2 cryogenic cooling	682
Liquid nitrogen (LN <sub>2</sub> ) cryogenic cooling	684
Valves	685
Valve box assembly	686
Valve driver assembly	688
Valve driver assembly (original style)	689
Valve driver assembly (new style)	691
Valve actuator assembly (1 of 2)	694
Valve actuator assembly (2 of 2)	696
Miscellaneous Valve Parts	698

Electronics and Fans	700
Electronics carrier	701
Keyboard Assemblies	705
7890A Keyboard and Display (original style)	705
7890A Keyboard and Display (new style)	707
7890B Keyboard and Display	709
AC circuit board components	711
AC power	713
Power cords	715
Main transformer	716
Analog input board (G1556A accessory)	718
Chassis fans	719
Covers	720
Plastic covers	721
Back covers	723
Auto-injector mounting and parking posts	724
7890A Side Carrier	726
7890B Side Carrier	727
Mainframe	728
Tools	729
Shipping Crates	729
Hand tools	730
Specialized tools for GC support	734
Electronic tools for GC support	736
Pneumatic tools for GC support	738
For Checking Inlet Pressure	740



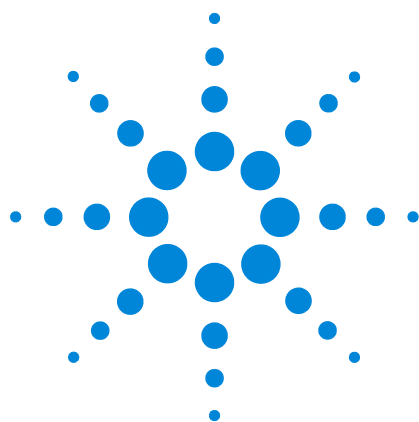

## 1 Before Servicing the Instrument

Safety	20
Hand tools	20
Specialized tools for GC support	24
Electronic tools for GC support	26
Pneumatic tools for GC support	28



## Safety

Before servicing the GC, observe the safety precautions in the [7890 Series GC Safety Manual](#).

## Hand tools

We assume you have a well-equipped toolbox, but some of the special tools listed here will be helpful.

**Table 1** 7890 Series GC tool bag, 5182-3456

Description	Part no.
Wrench, 1/4-inch to 5/16-inch	8710-0510
Wrench, open-end, 7/16-inch and 9/16-inch	8710-0803
Wrench, open-end, 1/2-inch and 7/16-inch	8710-0806
Wrench, open-end, 7/16-inch and 3/8-inch	8710-0972
Driver, nut, 7-mm	8710-1217
Cutter, tube	8710-1709
Screwdriver, Torx T-20	8710-1615
Screwdriver, Torx T-10	8710-1623
Tool bag	9222-3253
Tool, ferrule, .25–.32 mm, 2 x 3 x 3/4-inch	RFT-2500
Tool, ferrule, .53 mm	RFT-5300
Lens, magnifying	G2855-40001
Septum tool, knurled handle	450-1000





### Drivers

- T-20 Torx Driver — 5182-3465
- T-20 Torx Key (for close-quarters work) — 8710-1807
- T-10 Torx Driver — 5182-3466

- T-10 Torx Key (for close-quarters work) – 8710-2140
- #1 Pozidriv screwdriver – 8710-0899
- #2 Pozidriv screwdriver – 8710-0900
- Flat-blade screwdriver – 8730-0008
- 1/4-inch Nut Driver – 8710-1561 (FID Jets)
- 7-mm Nut Driver – 8710-1217

### **Wrenches**

- Wrench, 1/4-inch to 5/16-inch – 8710-0510
- Wrench, open-end, 7/16-inch and 9/16-inch – 8710-0803
- Wrench, open-end, 1/2-inch and 7/16-inch – 8710-0806
- Wrench, open-end, 7/16-inch and 3/8-inch – 8710-0972
- Wrench, adjustable, 12-inch – 8710-1712
- Wrench, open-end, 3/16-inch x 1/4-inch – 8710-2618
- Wrench, open-end, 9/16-inch x 5/8-inch – 8720-0010
- Wrench, 7-mm x 8-mm – 1340407010
- Wrench, 6-mm x 7-mm – 1340407011
- Wrench, 8-mm x 10-mm – 1340407012
- Wrench, 6-mm x 6-mm – 8710-2156
- Wrench, 5-mm x 5-mm – 8710-2157

### **Cutters/Crimpers/Strippers**

- Diagonal Cutters
- Crimper/Wire Strippers
- Tubing Cutters
  - Restek #20193 <http://www.restekcorp.com>
  - Agilent 8710-1709
  - Agilent precision tubing cutter for 1/16-in. SS tubing 5190-1442
  - Alternate for Europe:

HICHROM Ltd - Part No: HI-196

Description: Tubing Cutters - Price: £17.00

1 The Markham Centre Station Road – Theale, Berkshire

RG7 4PE, UK Telephone: +44 (0) 118 930 3660

[mailto: salesw@hichrom.co.uk](mailto:salesw@hichrom.co.uk)

### **Pliers**

- Channel lock pliers
- Slip-joint pliers
- Needle nose pliers
- Agilent needle nose pliers – 8710-0004
- Agilent Truarc pliers – 8710-0018

## Specialized tools for GC support



- FID Flow Tool – 19301-60660
- FID Cleaning Kit – 9301-0985, (.010-inch Guitar String)
- Injection Port Cleaning Kit – 480-0003

- 12 Piece File Kit – RSF-1200
- (Ferrule Removal tool, not shown – 440-1000)
- Fiberglass Tape – 0460-0186
- PTFE Tape – 0460-0016
- White Gloves – 8650-0030




- NPT Adapter – 0100-0118
- 5/16-inch by #20 Thread Chaser (Die), to clean 1/8-inch Swagelok threads – Other Supplier
- On/Off Valve, 1/8-inch Ball Valve – 0100-2144
- Exacto Knife
- Metric (8710-0641) and English Allen Wrench – Hex Keys
- Pin Vise and small drill bit
- Tweezers – 8710-0007
- Micro Probes (Sharp Object) – RMP-5005 – Qty 5
- Inspection Mirror – 707-0027
- 1.25-inch socket and ratchet for removing FID brass retaining nut
- Column cutter wafer, 4/pk – 5181-8836
- 6890 inlet wrench – 19251-00100
- 7890 inlet wrench (for SSL and MMI) – G3452-20512



## Electronic tools for GC support



- Digital Multimeter – (Agilent U1231A True RMS shown)
- Power Outlet Test Tool (Radio Shack/Sears/Electrical Supply)
- Static Strap – 9300-0969 (small) or 9300-0970 (large), 5 ft cord/clip – 9300-0980
- Wrist strap, disposable, 4-LG 1-W – 9300-1408
- Various jumpers and clip leads – electronic supply
- ESD mat – 9300-1484
- Useful Cables – not shown:
  - Crossover LAN Cable – 5183-4649
  - 9 Pin RS-232 Null Modem Cable – G1530-60600
  - USB/RS-232 adapter cable – 8121-1013
- Sensor pin extractor – 8710-1542

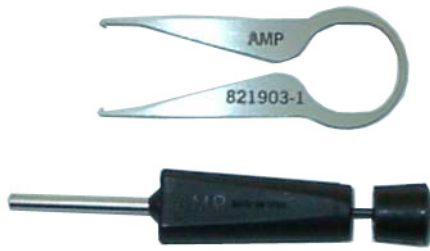

## Pneumatic tools for GC support







- Electronic Leak Detector
  - G3388B (or G3388A)
  - 5182-9646, 120 VAC
  - 5182-9648, 220 VAC (Agilent)
- Rotameter, Porter Model 65
 

Tube 1: 10-70 cc/min He, Tube 2: 85-850 cc/min Air

[http://www.porterinstrument.com/indust\\_21.php](http://www.porterinstrument.com/indust_21.php)
- Electronic Flow Meter
  - Preferred flow meter: Precision gas flow meter, flow range 5-500 mL/min – 5067-0223



- ADM 1000 – 220-1170
- ADM 2000 – 220-1171-U (Mass Flow Version)
- Electronic "Mass" Flow Meter – Flow Tracker
  - Model 1000 Flow Only – 5183-4779
  - Model 2000 with Leak Detector – 5183-4780
  - Universal Power Adapter – 5183-4781
- White Silicone Tubing:
  - 4 meters – 701-0016
  - 7 cm – 220-1179
- PTFE tubing, 1/8-inch, 1 m long – G1290-80050

## For Checking Inlet Pressure

- Quality Analog Pressure Gauge (0-100 or 0-60 PSIG) adapted with a Headspace Probe (301-013-HSP) – need NPT adapter – 0100-0118 and TCD Ferrules Set – 5180-4103, 5182-9673, 5182-3477
- Digital Pressure Gauge (0-100 PSIG) – Omega – Model HHP-201 Adapted with 1/8-inch NPT adapter 0100-0111, SS Capillary 5021-1831, Reducing Ferrule – 0100-1342 and 530  $\mu$ M Column Ferrule – 5080-8773



NPD Flow adapter - G1534-60640



G1530-20610 - FID/NPD Jet Plug



"No-Hole" Column Ferrules – 5181-7458



FID Flow adapter - 19301-60660



5060-9055 uECD/TCD Detector Plug




Plastic 1/8-inch Swagelok Cap – P/N 0100-2414 Qty 3

Replacement Septum Purge Fitting for EPC Module – P/N G2630-61720 (Not included in kit)

## 1 Before Servicing the Instrument



Blanking plug, finger tight, 5020-8294. (Can be used with any ferrule to duplicate using a no-hole ferrule.)

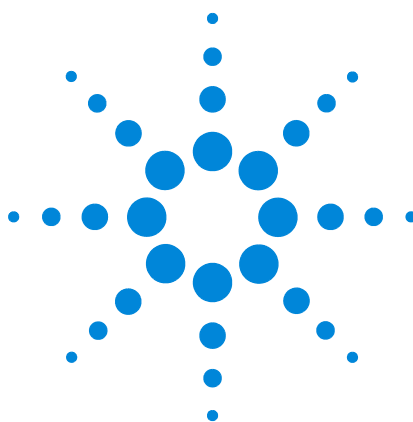

## 2 Inlets

Split/Splitless Inlet	35
Multimode Inlet	43
Cool On-Column Inlet	58
Purged Packed Inlet	65
Programmed Temperature Vaporization Inlet (PTV)	73
Volatiles Interface	87
Solvent Vapor Exit Accessory	93



## Inlet overview

Inlets are a means of introducing a sample into the carrier gas stream and then into the analyzing column. The term "Inlets" includes:

- Injection ports, for use with a syringe
- Sampling valves, both liquid and gas, for high sample size reproducibility
- Vapor inlets, in which the sample is vaporized in an external sampler, then swept into the GC. Headspace and thermal desorbers are examples.
- Solid samplers
- and other devices.

## Split/Splitless Inlet

This is the general-purpose inlet for use with capillary columns. It vaporizes the sample in a heated liner, then delivers all (splitless mode) or a specified fraction (split mode) of the vapor to the column.

### Replacing the split/splitless inlet

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

#### **WARNING**

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

---

#### **CAUTION**

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.


---

#### **WARNING**

**Inlets, detectors, and the oven are insulated with fibrous materials which may cause irritation to skin, eyes, and/or mucous membranes. Always wear gloves when working with the insulation. Additionally, if the insulation is flaky/crumbly, wear protective eyewear and a suitable breathing mask and/or respirator.**


---

- 3 Inside the oven, remove interferences as needed: column(s), the inlet insulating cup, and so forth.
- 4 Remove the split/splitless top weldment assembly. See [“Replacing the split/splitless inlet top weldment assembly”](#).
- 5 At the top of the GC, remove screws from the injection port top cover and remove the cover to expose the inlet assembly.

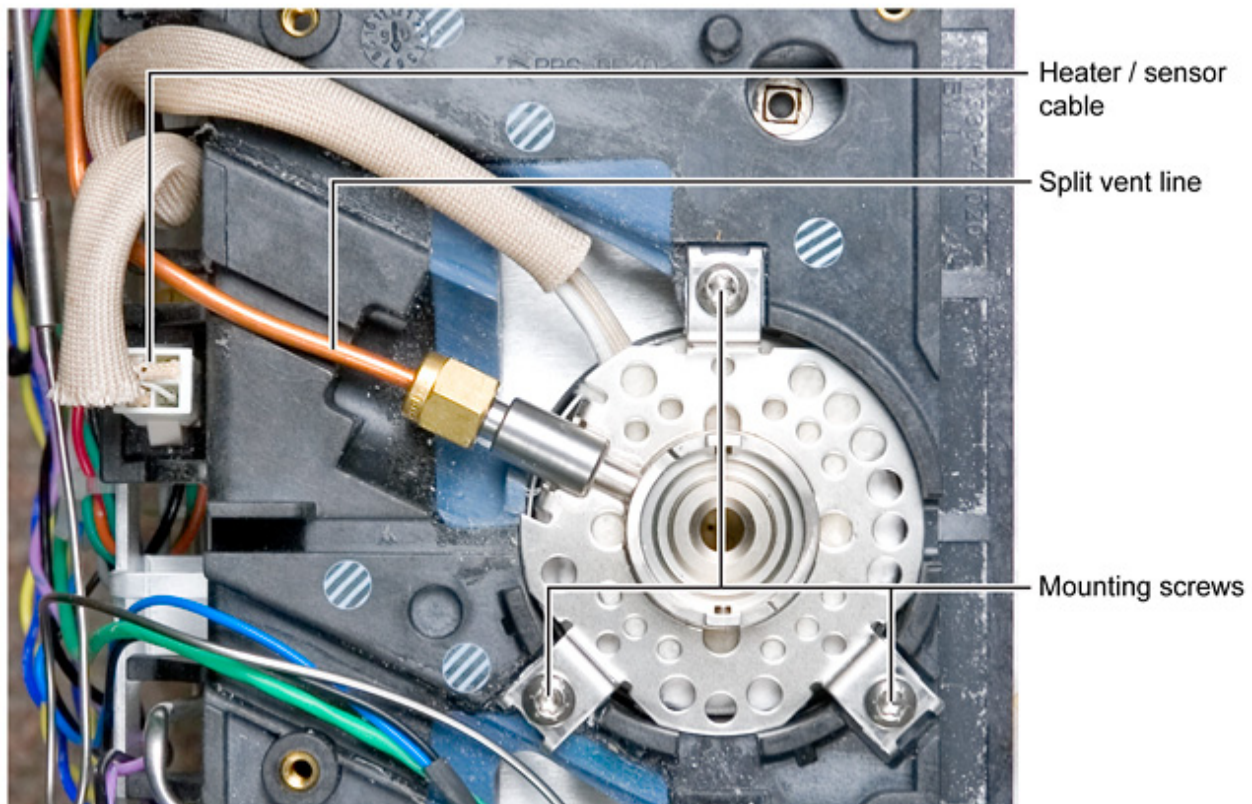

- 6 Disconnect the split vent line.
- 7 Disconnect the sensor/heater cable and work it back through any interfering wiring, tubing, and/or GC frame members.
- 8 Remove 3 screws to release the inlet body assembly with its insulation and cup and lift the assembly from the GC.



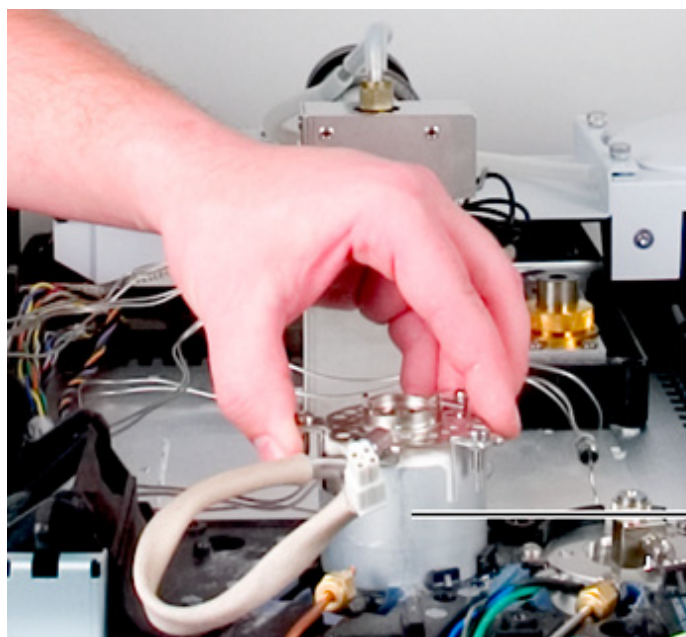

Split / splitless inlet assembly

- 9 Replacing the inlet assembly is the reverse of these steps. Pay attention to rotational orientations of the inlet body, the preformed insulation and its cup, and positions of the split vent connection and the heater/sensor cable as you fit the assembly into the GC.


## Replacing the split/splitless inlet top weldment assembly

- 1 Turn off all gas flows at their sources.
- 2 Turn off the GC main power switch.

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

- 3 Rotate the locking tab to release the top weldment assembly from the split/splitless inlet.

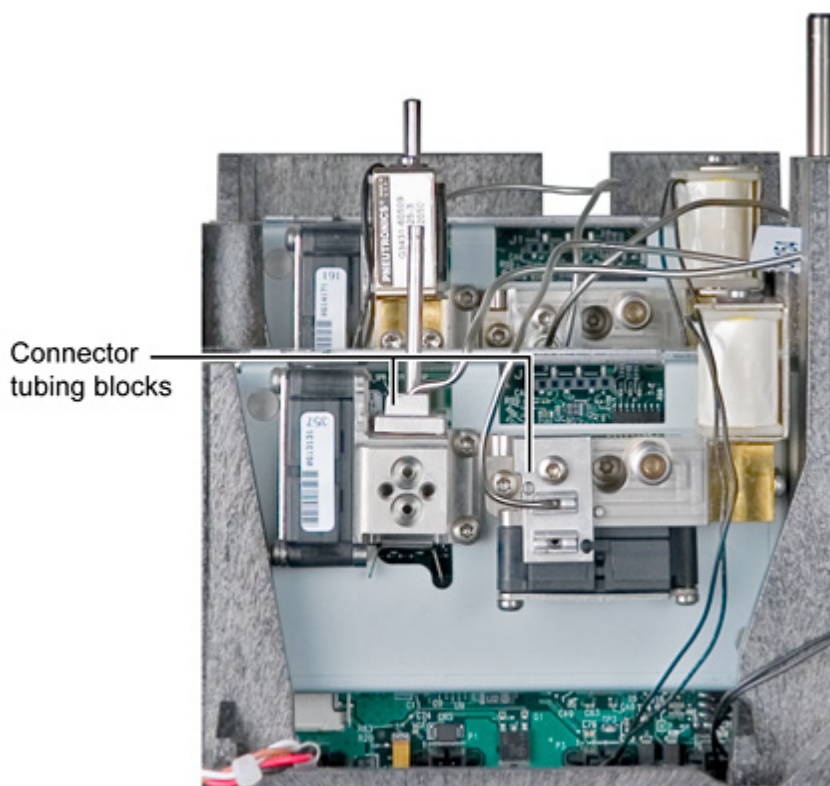

- 4 Disconnect tubing connector blocks at the EPC module locations (T-10 screwdriver). Take care to make sure sealing O-rings are not dislodged or damaged.
- 5 Guide the top weldment assembly carefully through any interfering wiring and/or tubing to remove it from the GC.
- 6 Reassembly is the reverse of these steps.

## Replacing the split/splitless inlet split vent line


- 1 Turn off all gas flows at their sources.
- 2 Turn off the GC main power switch.

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 3 Remove the top EPC module cover and the left side panel on the GC.

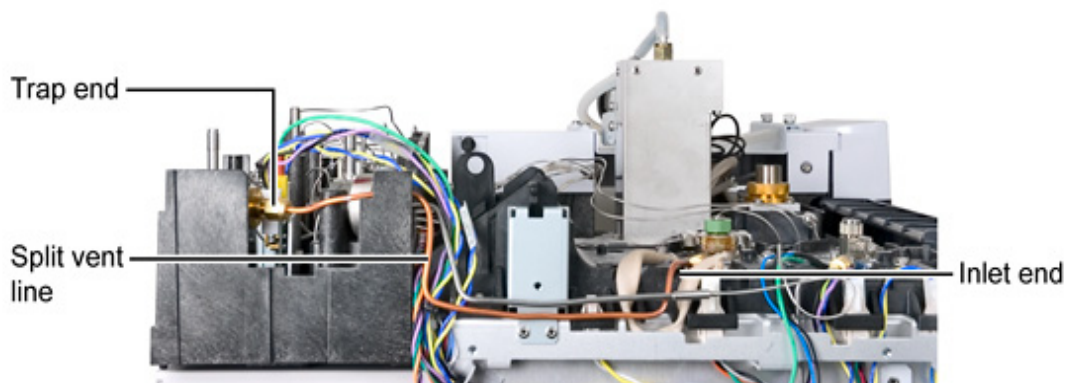

- 4 Remove the inlet cover.
- 5 Disconnect the split vent line from the inlet and from the trap canister. At the trap end, use two wrenches against each other to prevent the trap body from rotating.

### CAUTION

When removing/attaching the split vent line at its trap canister, always use one wrench to support the canister fitting and one to tighten the nut. Failure to do this could break the seal within the canister.

---

- 6 Remove the split vent line for cleaning or replacement.
- 7 Replacement is the reverse of these steps.

## Replacing the split/splitless inlet split vent filter canister and valve assembly

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.


### WARNING

Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

- 3 Remove the top EPC module cover.

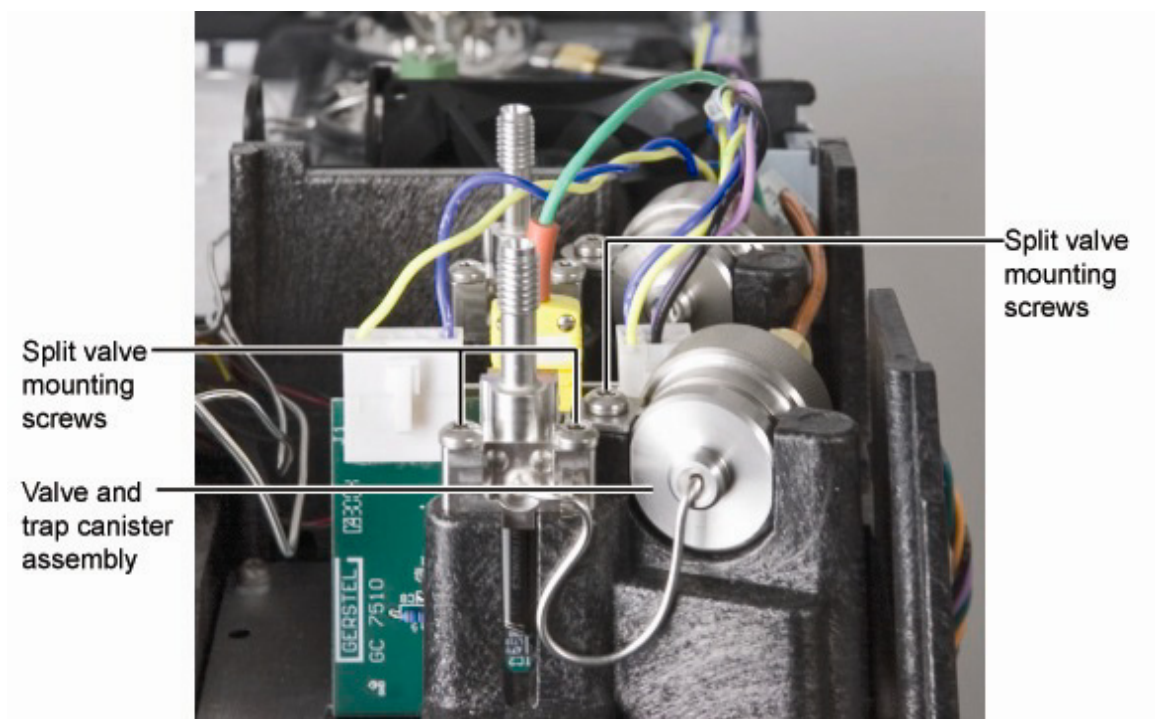

- 4 Disconnect the cable to the split valve from the inlet EPC module. See “Replacing the split/splitless inlet split vent line”.

- 5 Disconnect the split vent line. See [“Replacing the EPC modules”](#).
- 6 Loosen the screw securing the canister retainer enough to rotate it aside. Loosen 2 screws to release the split valve assembly.
- 7 The canister and split valve are removed as a unit. Note that the split valve is now easily accessible for replacement if needed.
- 8 Reassembly is the reverse of these steps.

## Reinstalling the split/splitless inlet

- 1 Make sure the heater/sensor assembly is installed and the inlet insulation sleeve is in place.
- 2 Install a column nut and blank ferrule on the bottom of the inlet to prevent insulation contamination, and place the inlet into the inlet carrier.

**NOTE**

Make sure the insulation is properly seated around the inlet and that the heater/sensor wiring harness insulation sleeve is tucked under the top inlet plate.

- 3 Retighten the three screws (Torx T-20) to secure the top inlet weldment plate to the inlet carrier.
- 4 Reconnect the split vent flow line.
- 5 Reinstall the top insert assembly (with septum and carrier lines attached). Line up the tab on the bottom of the top insert assembly with the hole on the insert weldment assembly so that the locking tab is on the left. Push down to connect. Slide the locking tab to the back.
- 6 Route the septum purge and carrier gas lines along the oven side of the inlet carrier and into the chassis slots leading to the EPC module.
- 7 Seat the heater/sensor leads into the channel on the inlet carrier.
- 8 Reconnect the heater/sensor assembly into the provided connector (front or back) on the left side of the GC.
- 9 Reinstall the insulated thermal cup and insulation in the GC oven.
- 10 Reinstall the inlet cover.

## Multimode Inlet

The multimode inlet (MMI) is a programmable temperature vaporization inlet based off of the standard split-splitless inlet. It has six modes of operation: split, splitless, pulsed split, pulsed splitless, solvent vent, and direct. The inlet has the same internal configuration and dimensions of the split/splitless inlet, and uses the same liners, consumables, ferrules, O-rings, and septa as the split/splitless inlet.

The differences include:

- Cooling jacket for CO<sub>2</sub>, LN<sub>2</sub>, and air cooling
- Free-spinning nut for column connection
- No gold seal

Like the PTV, the MMI requires a thermocouple PCB and additional harness connections.

**7890A and 7890A+:** Due to space limitations (you can only have one thermocouple PCB in the GC), only one MMI or PTV is allowed per GC.

**7890B:** Two MMIs can be installed.

## Replacing the multimode inlet

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

### WARNING

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

---

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---


### WARNING

**Inlets, detectors, and the oven are insulated with fibrous materials which may cause irritation to skin, eyes, and/or mucous membranes. Always wear gloves when working with the insulation. Additionally, if the insulation is flaky/crumbly, wear protective eyewear and a suitable breathing mask and/or respirator.**

---

- 3 Inside the oven, remove interferences as needed: column(s), the inlet insulating cup, and so forth.
- 4 Remove the multimode top weldment assembly. See [“Replacing the multimode inlet split vent line”](#).
- 5 At the top of the GC, remove screws from the injection port top cover and remove the cover to expose the inlet assembly.



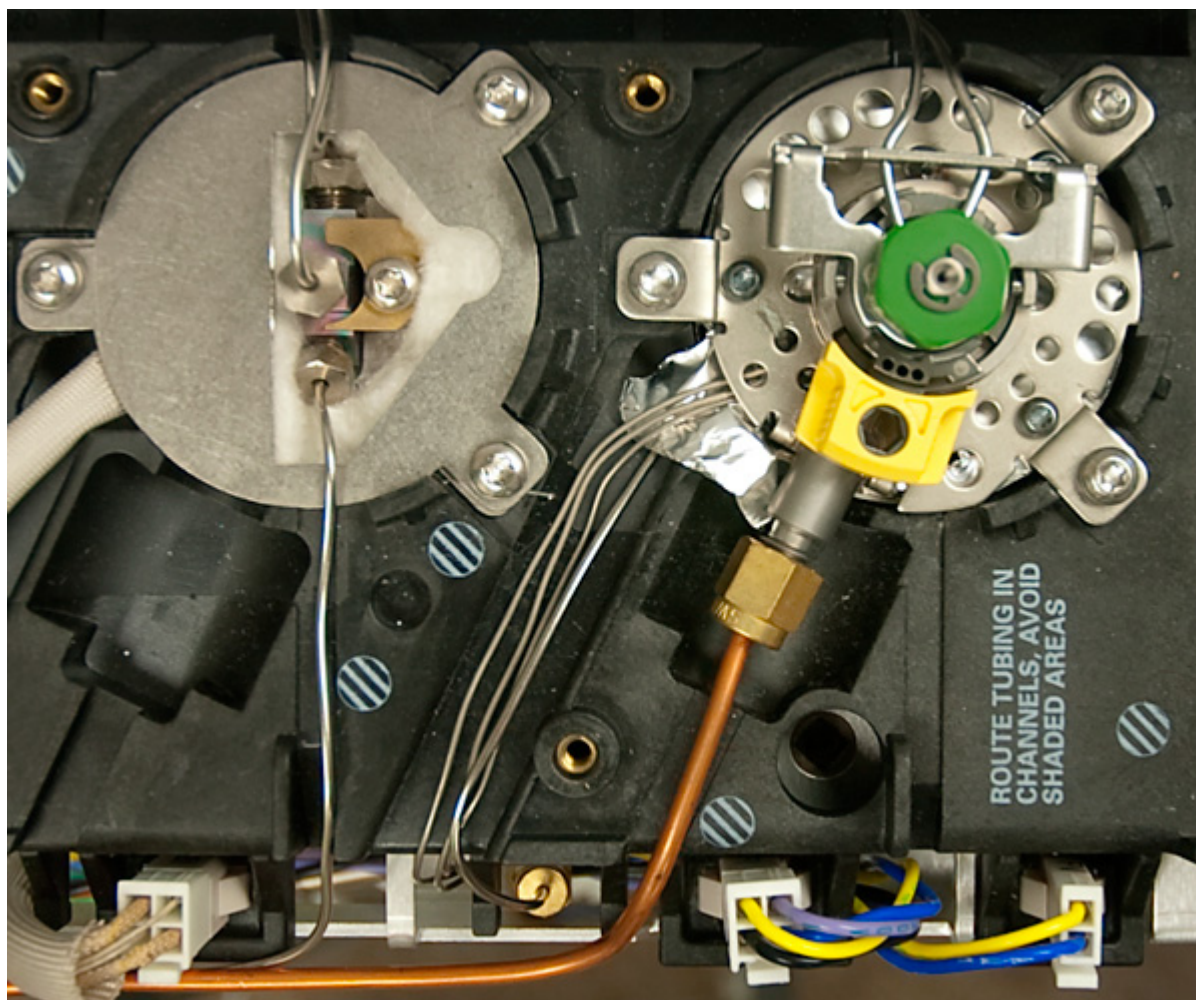

- 6 Disconnect the split vent line.
- 7 Disconnect the sensor/heater cable and work it back through any interfering wiring, tubing, and/or GC frame members.
- 8 Disconnect the plug labeled PCB from the wiring harness.
- 9 Disconnect the cryo line (use two wrenches against each other to prevent rotation).
- 10 Remove 3 screws to release the inlet body assembly with its insulation and cup and lift the assembly from the GC.
- 11 Replacing the inlet assembly is the reverse of these steps. Pay attention to rotational orientations of the inlet body, the preformed insulation and its cup, and positions of the split vent connection and the heater/sensor cable as you fit the assembly into the GC.


## Replacing the multimode inlet top weldment assembly

- 1 Turn off all gas flows at their sources.
- 2 Turn off the GC main power switch.

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

- 3 Rotate the locking tab to release the top weldment assembly from the multimode inlet.

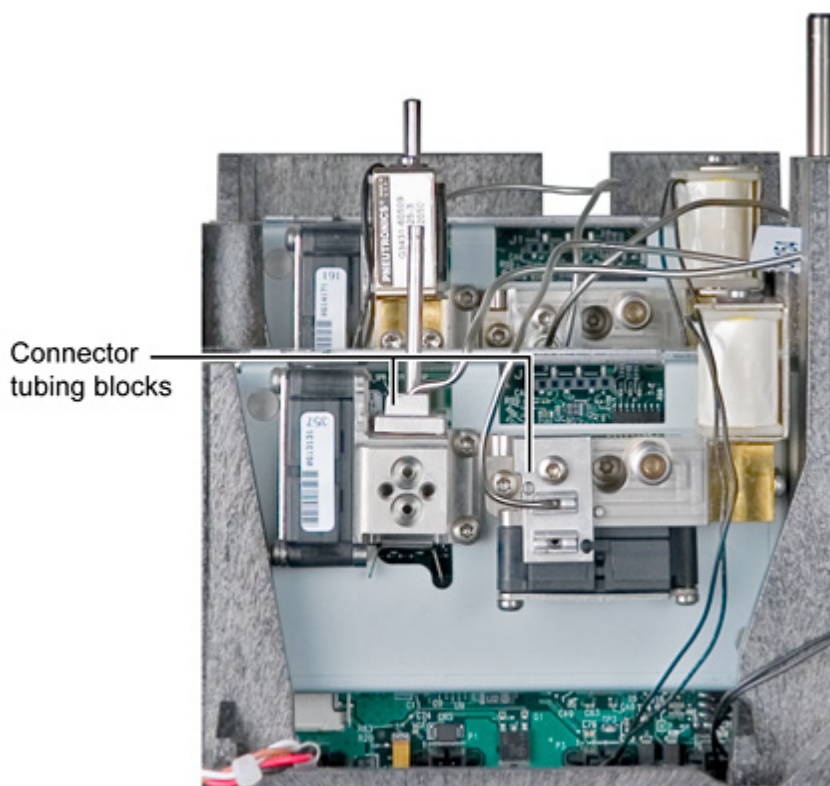

- 4 Disconnect tubing connector blocks at the EPC module locations (T-10 screwdriver). Take care to make sure sealing O-rings are not dislodged or damaged.
- 5 Guide the top weldment assembly carefully through any interfering wiring and/or tubing to remove it from the GC.
- 6 Reassembly is the reverse of these steps.

## Replacing the multimode inlet split vent line


- 1 Turn off all gas flows at their sources.
- 2 Turn off the GC main power switch.

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 3 Remove the top EPC module cover and the left side panel on the GC.

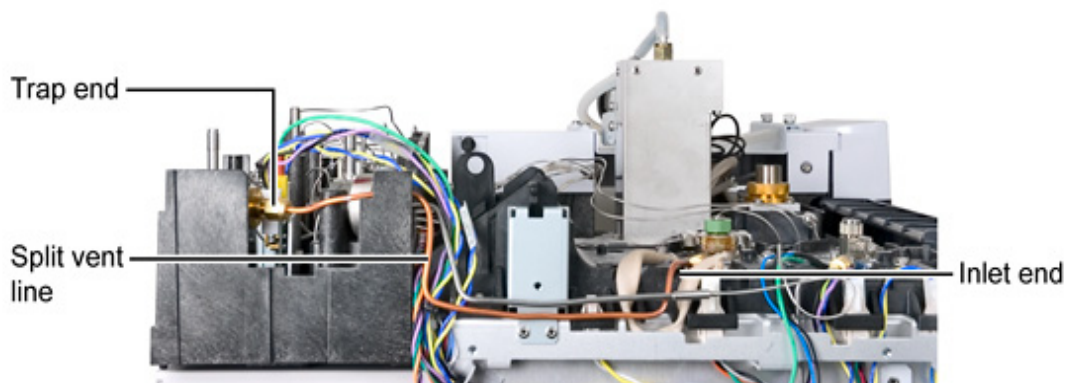

- 4 Remove the inlet cover.
- 5 Disconnect the split vent line from the inlet and from the trap canister. At the trap end, use two wrenches against each other to prevent the trap body from rotating.

### CAUTION

When removing/attaching the split vent line at its trap canister, always use one wrench to support the canister fitting and one to tighten the nut. Failure to do this could break the seal within the canister.

---

- 6 Remove the split vent line for cleaning or replacement.
- 7 Replacement is the reverse of these steps.

## Replacing the multimode inlet split vent filter canister and valve assembly

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.


### WARNING

Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

- 3 Remove the top EPC module cover.

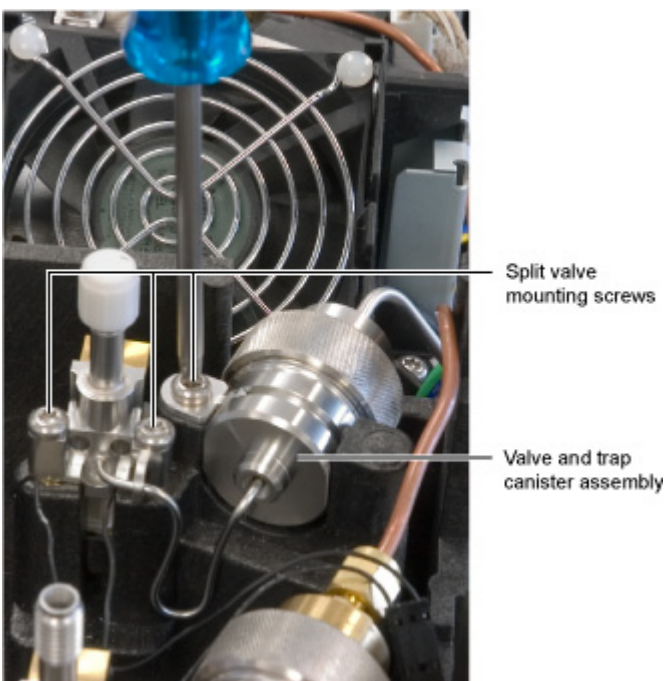

- 4 Disconnect the cable to the split valve from the inlet EPC module. See [“Replacing the multimode inlet split vent line”](#).
- 5 Disconnect the split vent line. See [“Replacing the EPC modules”](#).

- 6 Loosen the screw securing the canister retainer enough to rotate it aside. Loosen 2 screws to release the split valve assembly.
- 7 The canister and split valve are removed as a unit. Note that the split valve is now easily accessible for replacement if needed.
- 8 Reassembly is the reverse of these steps.

## Replacing a Multimode CO2 cryo valve

### **WARNING**

**High pressure or extremely cold gas. Before replacing the cryo valve, turn off the cryo line at the source, and depressurize the cryo line to the GC.**

---

If changing the GC configuration, see [“Changing the GC configuration”](#) for important information regarding GC methods. Then proceed with the steps below.

- 1 Cool down the oven and all heated zones.
- 2 Turn off all gas flows at their sources.
- 3 Turn off the GC main power switch and disconnect its power cord.

### **WARNING**

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**


---

### **CAUTION**

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---



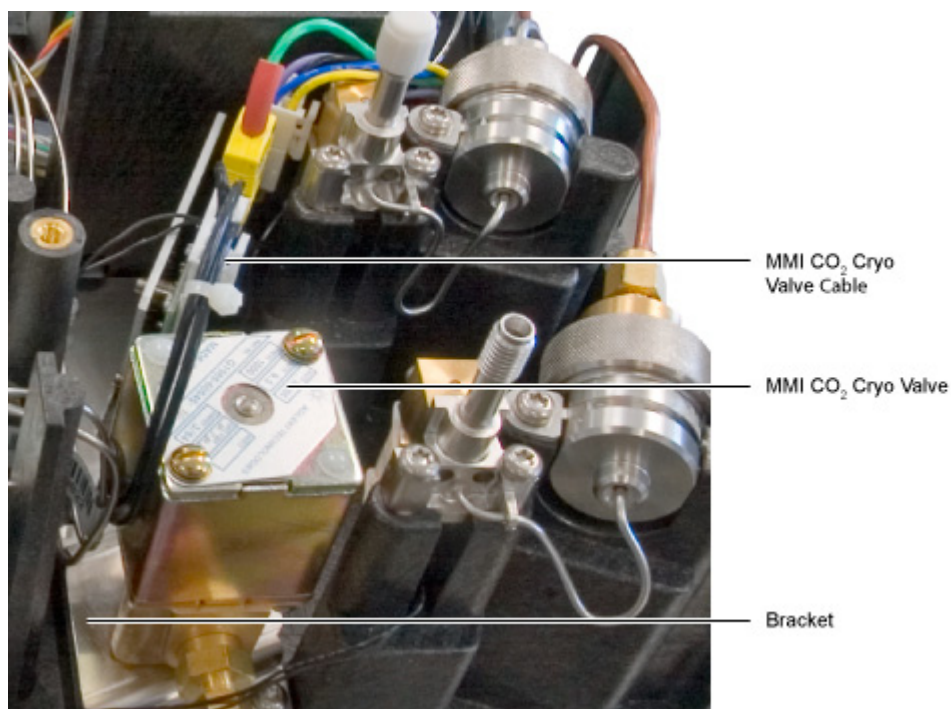

- 4 Remove the valve via its two mounting screws accessible under the removed bracket assembly. See [“Removing the MMI cryo/PC board bracket”](#) for details.
- 5 When reinstalling, make sure the valve is oriented correctly on the bracket with its cryo feed input connection towards the rear of the GC and that its cable is reconnected to the nearby PC board.

## Replacing the Multimode N2 cryo valve

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.


### WARNING

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

- 3 Remove the inlet cover and GC left side panel.
- 4 Remove the N2 cryo feed tube between the cryo valve and the multimode inlet.
- 5 Remove two screws from the multimode N2 cryo feed cover plate. Remove the cover plate.

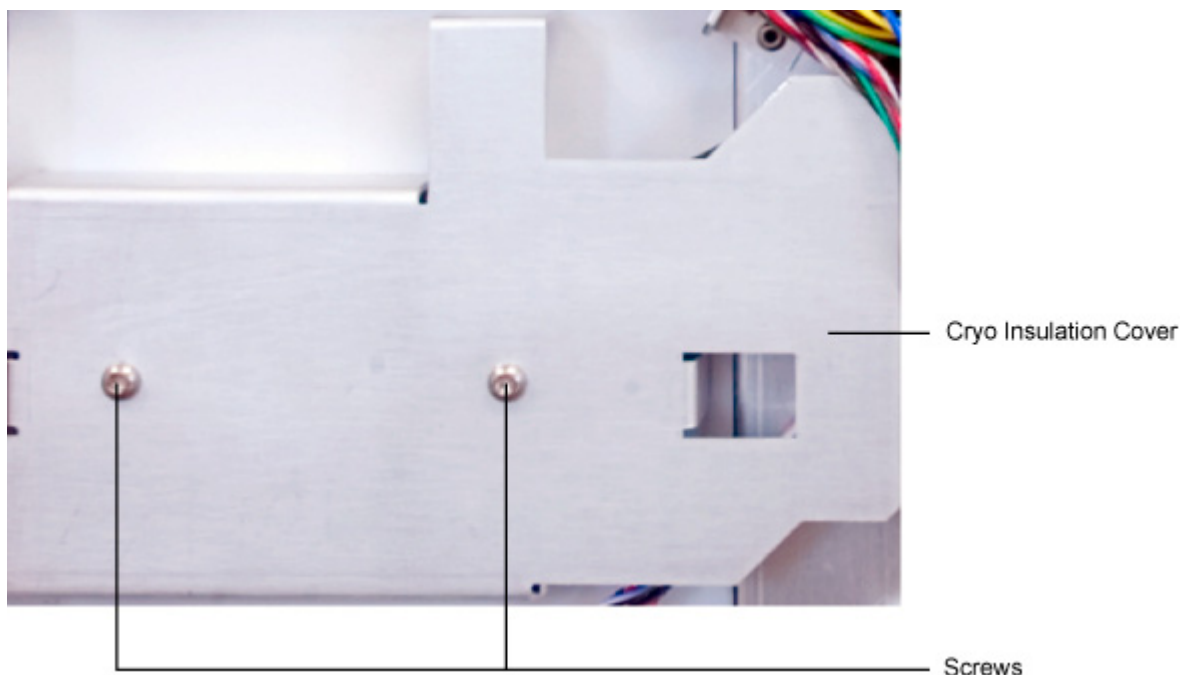


- 6 Disconnect the multimode N2 cryo valve cable. If you intend to replace this cable with a new one, use the existing cable as



a model while you install the new cable, then remove the old cable.

- 7 Remove the valve bracket and cryo valve by removing three screws.

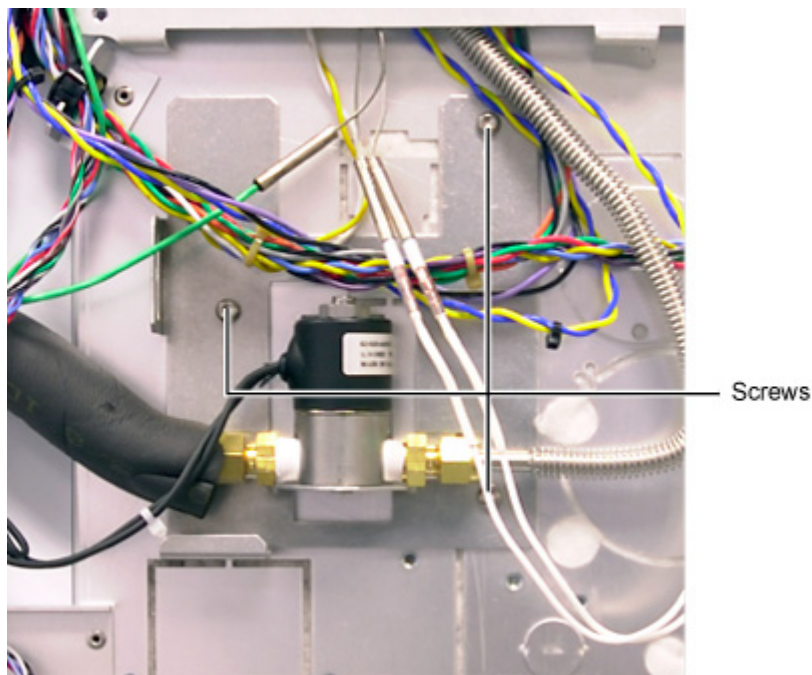

- 8 Replace the N2 cryo valve assembly by reversing these steps.
- 9 Reinstall the N2 cryo tube between the valve and the inlet.

## Removing the MMI cryo/PC board bracket

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.


### WARNING

Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

- 3 Remove the top EPC module cover.

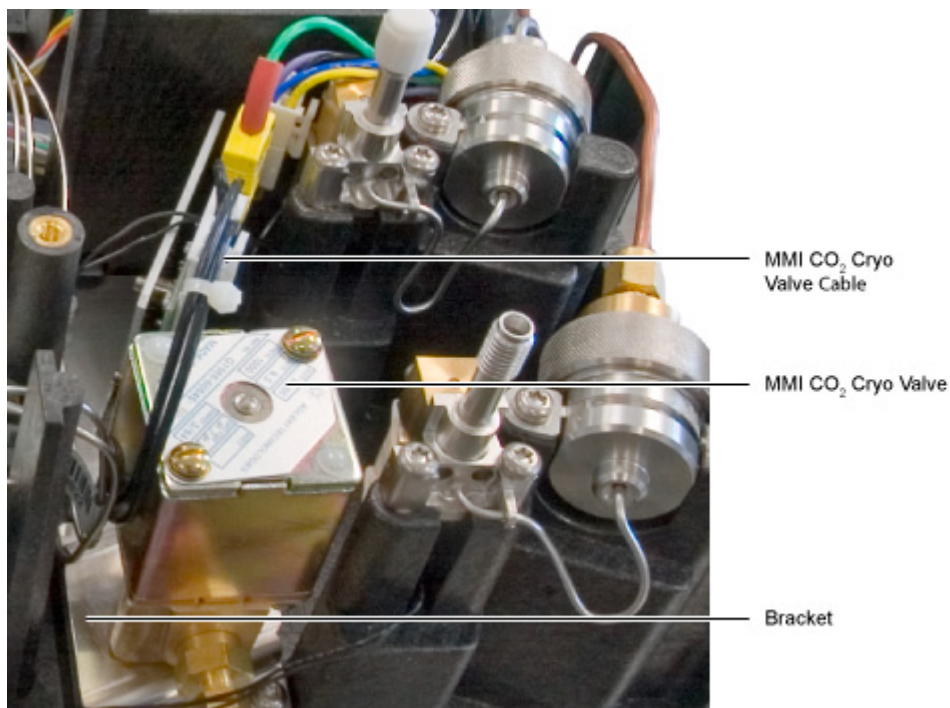

- 4 Disconnect cables to the multimode heater/cryo PC board. See [“Replacing the MMI heater/cryo PC board”](#).
- 5 If a multimode CO<sub>2</sub> valve is present, disconnect its cryo feed lines (CO<sub>2</sub> in / CO<sub>2</sub> out).
- 6 Remove two screws securing the bracket to the GC.

- 7 Lift the bracket assembly out of the GC working it around any interfering wiring and/or tubing.

**NOTE**

The multimode CO<sub>2</sub> cryo valve, if present, is also mounted on this bracket.

---

- 8 Reassembly is the reverse of these steps.

## Replacing the MMI heater/cryo PC board

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

### WARNING

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**


---

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 3 Remove the top EPC module cover.

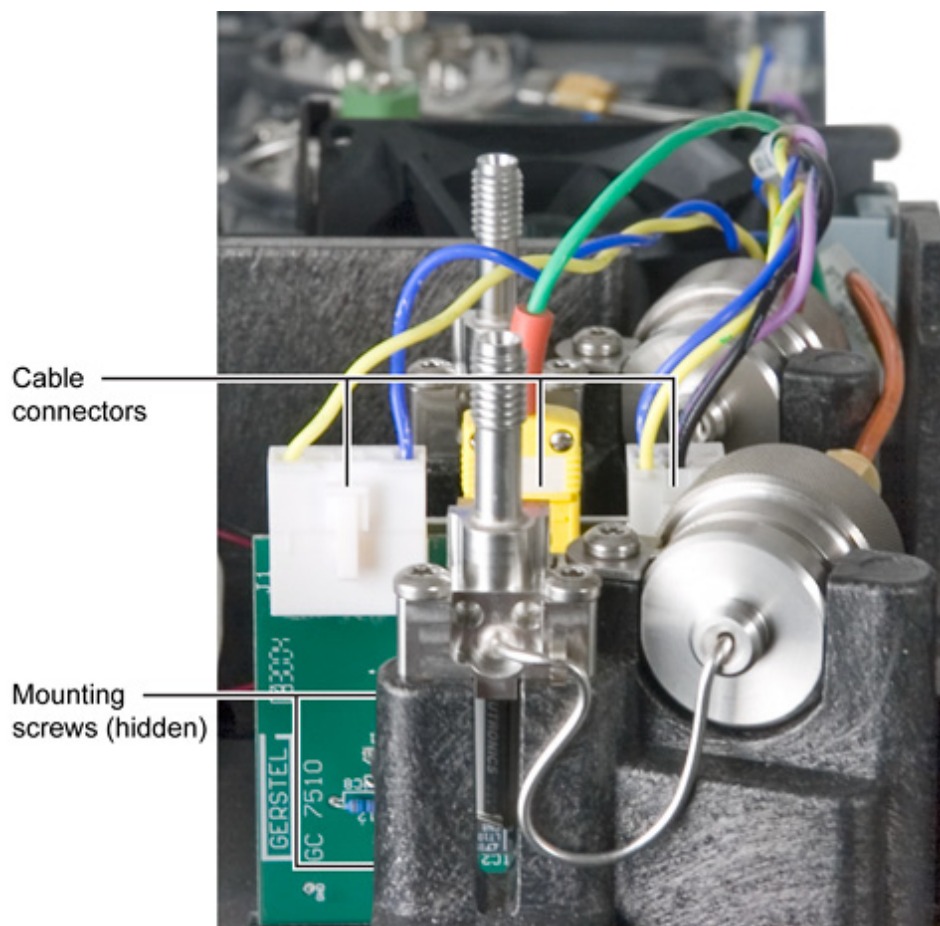

- 4 To remove the heater/cryo PC board, first remove its three cables: from left-to-right, respectively, the LN2 (liquid nitrogen) cryo valve cable (if present), the thermocouple cable, and the heater/sensor cable.
- 5 Next, remove two T-20 screws securing the board onto its bracket.
- 6 Lift the board straight up and out of the GC.
- 7 Installing the new board is the reverse of these steps.

## Cool On-Column Inlet

Cool On-Column (COC) inlets allow the injection syringe to deposit the liquid sample directly into the capillary column. This is accomplished by the use of an insert which aligns the syringe with the capillary column and the syringe needle.

### Replacing the COC inlet

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

#### **WARNING**

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**


---

#### **CAUTION**

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 3 From inside the oven, remove the column from the inlet fitting.
- 4 From the top of the GC, remove the inlet carrier cover (or the tray/tray bracket, if installed) and the left side cover.
- 5 Unclip the heater/sensor leads from the connector to the left of the inlet carrier.
- 6 Disconnect the inlet plumbing and remove from the chassis slot leading to the EPC flow module.
- 7 Disconnect the inlet plumbing pneumatics blocks from the EPC flow module (Torx T-10).
- 8 Use a Torx T-20 screwdriver to loosen the three captive screws that attach the inlet weldment to the top of the inlet carrier.
- 9 Slide the inlet up out of the inlet carrier. If necessary, you can also slide the insulation sleeve off the bottom of the inlet.
- 10 Reinstallation is the reverse of these steps. Use new O-rings when reconnecting the pneumatic tubing to the EPC flow module.

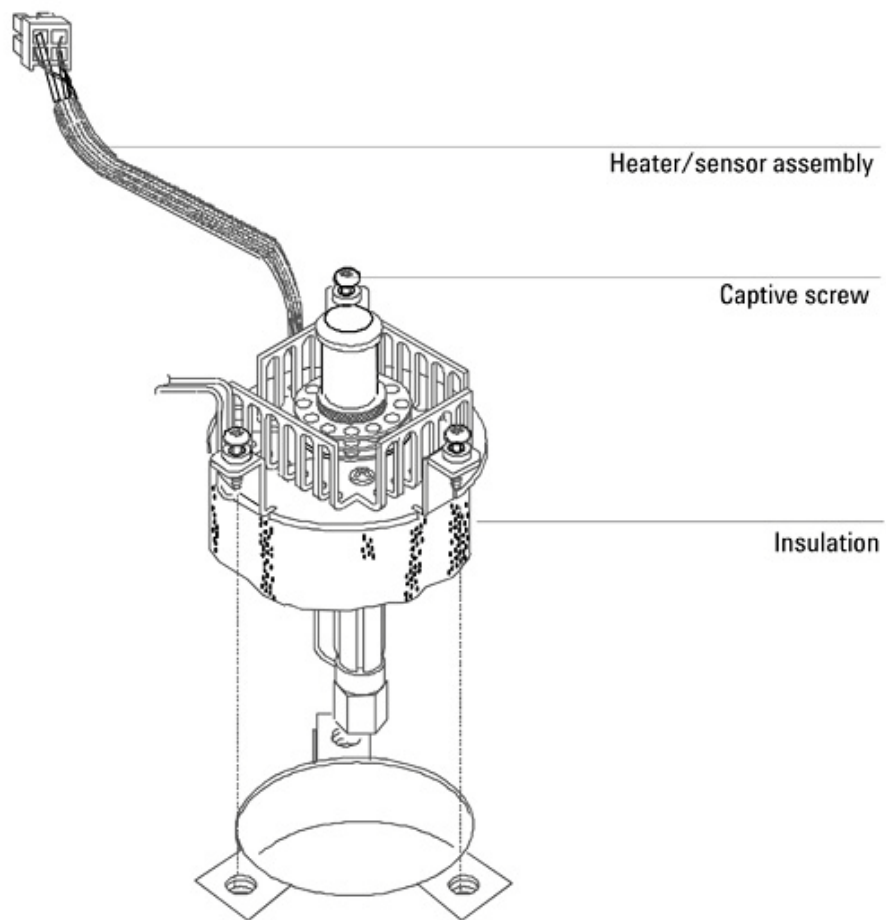

## Replacing the COC heater/sensor assembly

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

### WARNING

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

---

### WARNING

**Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.**

---

- 3 If necessary, remove the septum nut, cooling tower and/or needle guide to provide access to the two screws in the top of the cooling fin.


### NOTE

If desired, you can remove the entire inlet for better access.

---

- 4 Remove the screws securing the cooling fin to the inlet weldment and remove the fin.
- 5 Lift the heater/sensor leads out of the weldment channel and lift the assembly out of the inlet.
- 6 Install the new heater/sensor assembly and reassemble the inlet. You may need to use tweezers to seat the cable back in the channel and fully seat the heater/sensor in the weldment.



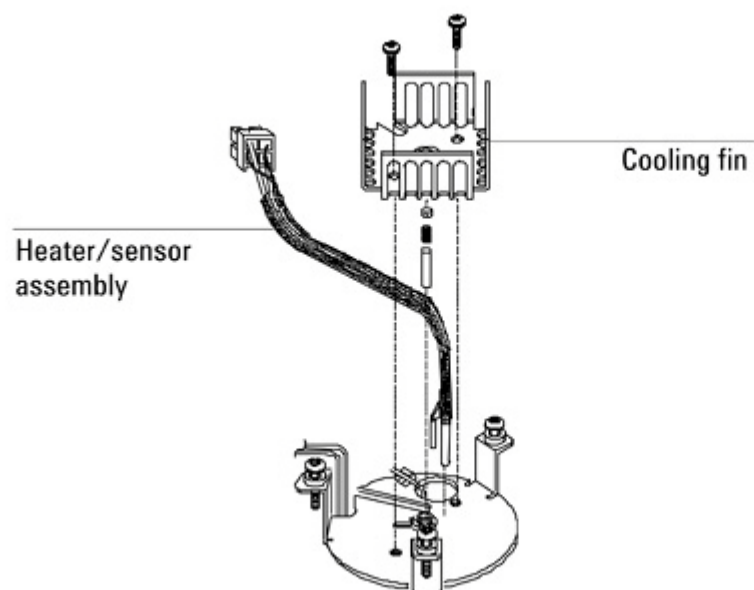

## Replacing an inlet cryo cooling tube

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

### WARNING

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

---

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

### NOTE

Unless otherwise noted, the installation process is the same for both LC02 and LN2 oven cryo valves, and for front or back inlet locations.

---

### Preparation of the inlet

- 1 Remove the septum retaining nut carefully: there is a small septum and very small coil spring held in place by the nut.

### NOTE

These parts, especially the coil spring, are easily lost. For their safety, it is recommended that the spring, septum, and septum retainer nut be kept together in a safe place for later reassembly.

---

- 2 Next, remove the heat sink by removing its two screws.
- 3 **If** the inlet is secured into the inlet chassis, release its three mounting screws and pull it straight upwards enough to expose the full length of the inlet body.
- 4 **If**, as you raise the inlet body, there is a large insulation pad either on the inlet body or in the inlet's mounting hole, it must be carefully removed to preserve its integrity.
- 5 Carefully inspect both the inlet mounting flange and the existing end of the inlet cryo cooling tube to be replaced:
  - Note the large, round hole in the flange - this is the pass through for the nozzle end of the inlet cryo cooling tube.

- Also note the smaller, square hole with a raised metal tab adjacent to the round hole - this is a location and alignment guide. The double right-angle bent portion of the inlet cryo cooling tube must be routed over the square hole and against the tab when installation is complete.
- 6 Remove the old tube from the inlet and entirely from the GC. Use two wrenches against each other at the T-fitting to avoid twisting tube.
  - 7 Install the new inlet cryo cooling tube into the opening in the mounting flange, as noted in step 5.
  - 8 Carefully work the insulation pad onto both the inlet body and inlet cryo cooling nozzle such that you end with both the inlet cryo cooling tube and the tube on the inlet body in the slit provided in the insulation pad (if two precut slits are present, use one for each tube). Work the pad upwards until it contacts the inlet mounting flange.

**NOTE**

In the next step, maintain control over the inlet cryo cooling tube as you insert and secure the inlet: the tube must not be allowed to slip deeper into the inlet cavity than as described in the inspection detail item in steps 3-8: basically, the tube's double right-angle bend must remain in its described position at all times.

- 9 Lower the inlet into its prepared location with its three screws aligned with the holes in the inlet chassis and its heater/sensor cable sitting in the trough in the inlet chassis. While doing so, mark where to bend the inlet cryo cooling tube into its channel across the inlet chassis.
- 10 Remove the inlet and bend the cryo tube at the marked location to pass between the inlet chassis and the GC frame.
- 11 Reinsert the inlet while guiding the inlet cryo cooling tube through the channel between the inlet chassis and the GC frame. Again, maintain control over the inlet cryo cooling tube inside the inlet, as noted above.
- 12 Start each mounting screw one at a time to insure the inlet is aligned and screws are properly threaded. Then tighten each screw in turn evenly until snug and the inlet secured.
- 13 From the left side of the GC, locate the inlet's heater/sensor cable and its corresponding connector plug. Tuck the heater/sensor cable underneath the routing tabs at the side of the GC and connect it to the heater connector.

- 14** In the following order, replace the heat sink, coil spring, septum, and septum retaining nut on the inlet. Route the inlet cryo cooling tube to the tee attached to the cryo valve.
- 15** Tighten the swage fitting where the inlet cryo cooling tube connects to the T-fitting at the oven cryo valve. Use two wrenches against each other to avoid twisting tubing.

This completes replacement of a inlet cryo cooling tube for a cool on-column capillary inlet.

## Purged Packed Inlet

The purged packed column inlet controls column flow by means of a forward pressure/flow, electronic proportional control valve. A forward pressure regulator controls flow out the septum purge vent. The inlet can be used for packed or capillary columns. Based on the column you have configured, the inlet can be operated in either a mass flow controlled mode or a pressure controlled mode.

### Replacing the Purged Packed inlet

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

---

**WARNING**

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

---

**CAUTION**

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 3 Turn off all gas flows at their sources.

---


**WARNING**

**Inlets, detectors, and the oven are insulated with fibrous materials which may cause irritation to skin, eyes, and/or mucous membranes. Always wear gloves when working with the insulation. Additionally, if the insulation is flaky/crumbly, wear protective eyewear and a suitable breathing mask and/or respirator.**


---

- 4 Remove the insulating cup and insulation from inside the oven.
- 5 Disconnect the column from the inlet and cap the column to minimize contamination.
- 6 Remove the inlet cover.

- 7 Note and record plumbing and wiring layouts and connectivity.
- 8 Remove the assembly stepwise. You will need Torx T-10 and T-20 drivers.
- 9 Disconnect the heater/sensor cable.
- 10 Disconnect the 1/16-inch tubing pair from the EPC module.
- 11 Remove the three top screws retaining the inlet to the GC.
- 12 Gently lift the assembly up and out.
- 13 Insert replacement inlet.
- 14 Reassembly is the reverse of these steps.
- 15 Use new O-rings during EPC reconnection.
- 16 Refer to the figures and photos in this section for guidance.

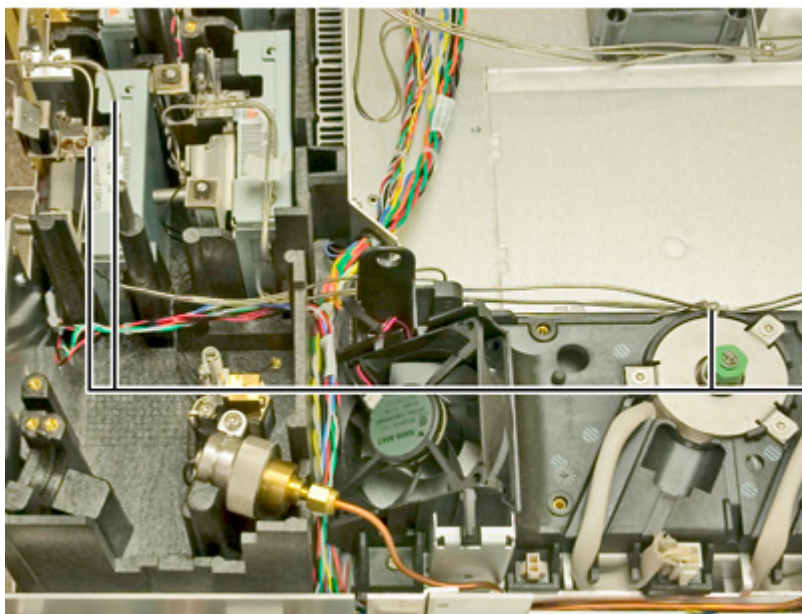

Gas tubing  
to/from PP inlet

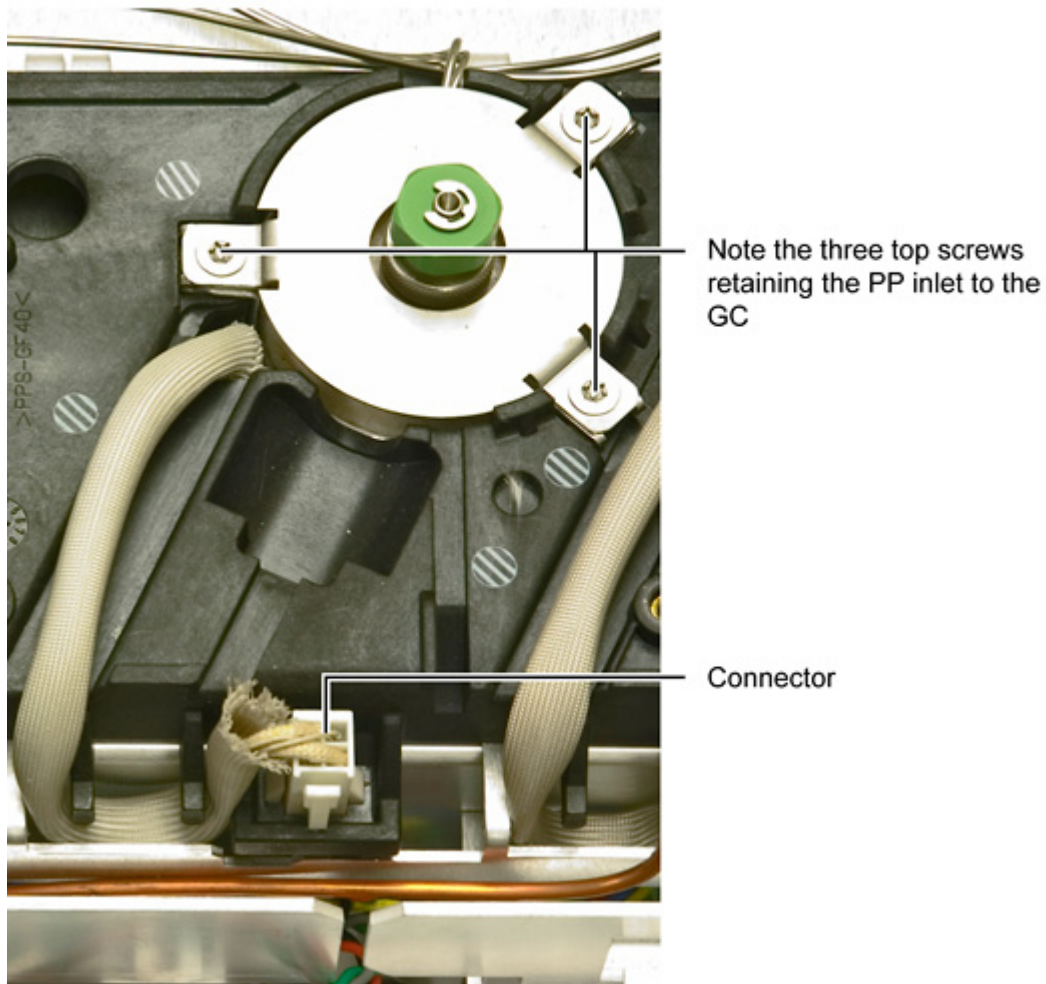

If the EPC module is disconnected, the inlet can be removed and replaced.

## 2 Inlets





This housing is designed to accommodate a variety of Agilent inlets.

Note fibrous insulation within housing. See earlier warning.





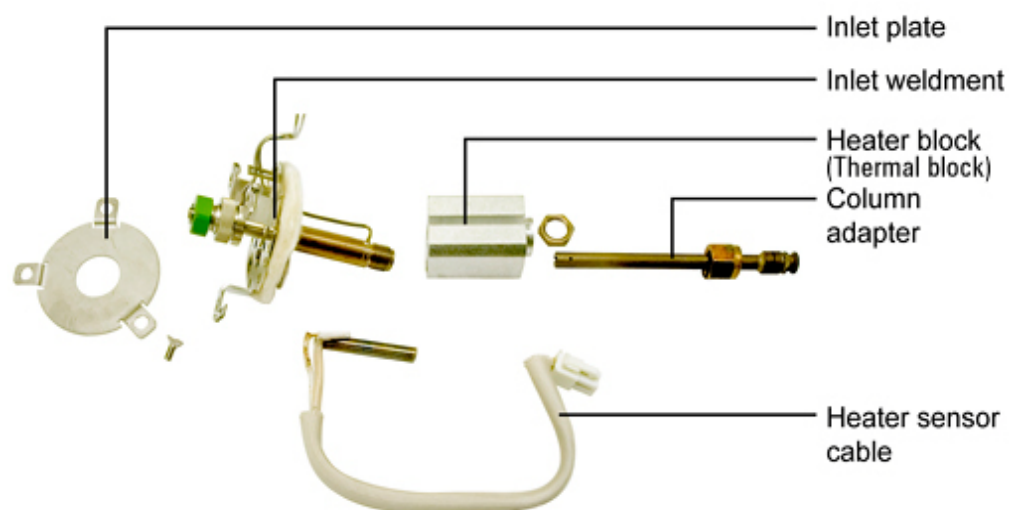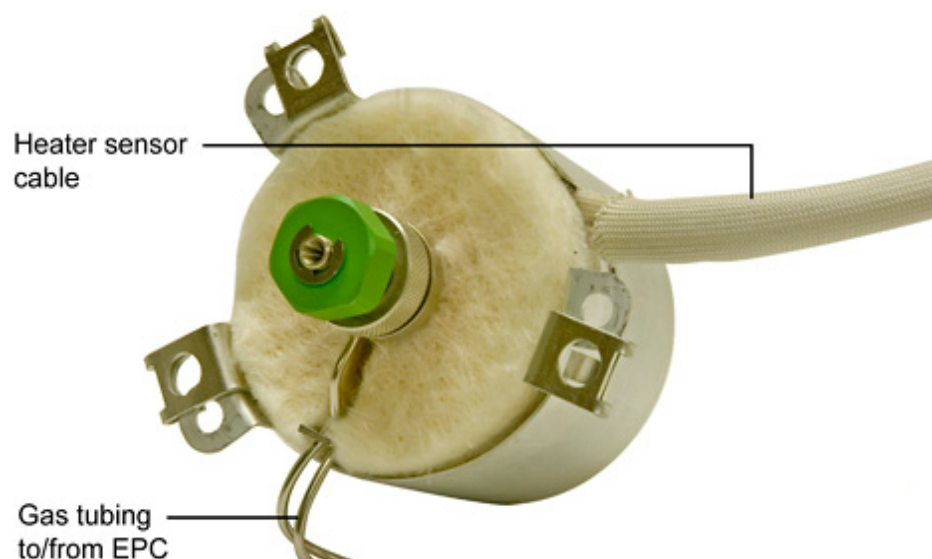

## Replacing the Purged Packed inlet heater

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

### **WARNING**

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

---

### **CAUTION**

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.


---

### **WARNING**

**Inlets, detectors, and the oven are insulated with fibrous materials which may cause irritation to skin, eyes, and/or mucous membranes. Always wear gloves when working with the insulation. Additionally, if the insulation is flaky/crumbly, wear protective eyewear and a suitable breathing mask and/or respirator.**

---

Heater replacement requires removal of the inlet. Refer to [“Replacing the Purged Packed inlet”](#) for removal instructions.

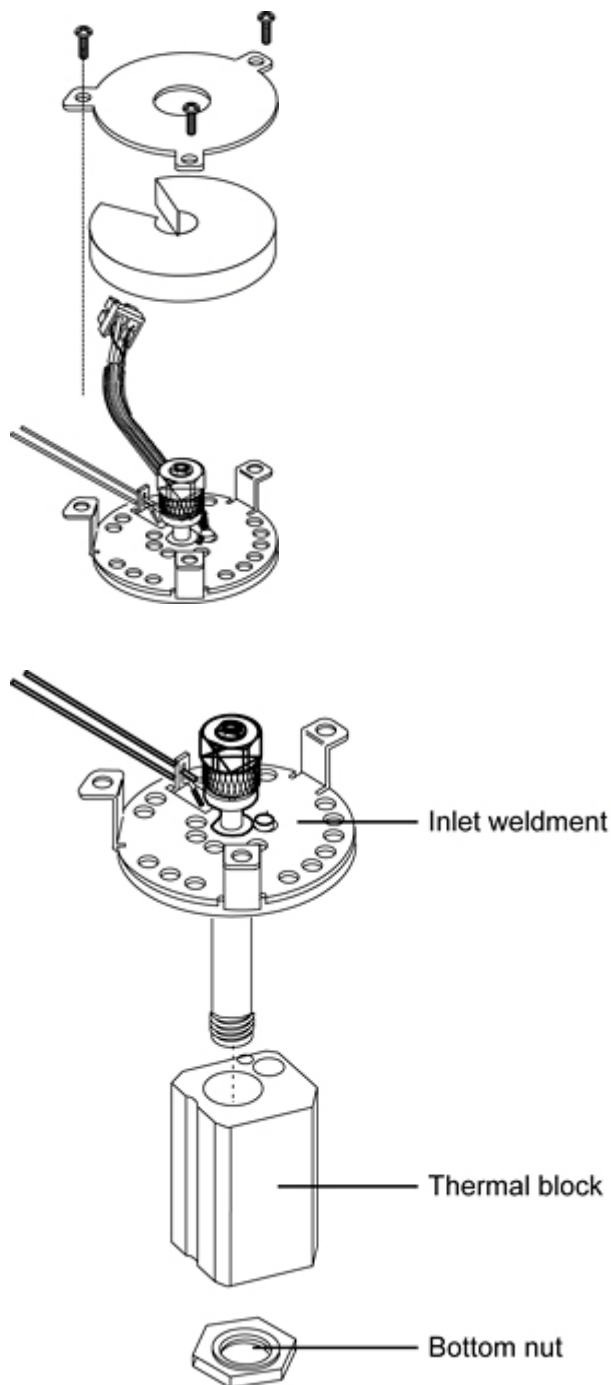

With the inlet removed:

- 1 Remove the tubing nut and column adapter.
- 2 Remove the bottom nut holding the thermal block.
- 3 Slide the thermal block off the inlet weldment.
- 4 Carefully remove the heater and sensor from the thermal (heater) block.

- 5 Replace heater/sensor cable assembly with a new unit.
- 6 Reassembly is the reverse of these steps.

## Programmed Temperature Vaporization Inlet (PTV)

In addition to the PTV consumables (inlet adapters, columns, PTFE ferrules, and septa), the replacement parts in the PTV inlet assembly are:

- The PTV inlet
- PTV weldment and front trap assemblies
- The filter
- The head assembly (septum or septumless)
- The cryo shroud on the inlet body (CO<sub>2</sub> or liquid nitrogen)
- The PTV thermocouple board
- The O-rings and restrictors in the gas fitting assembly.

## Replacing the PTV inlet

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

### WARNING

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

---

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---


The PTV inlet system is removed and replaced as a single unit consisting of the inlet and split vent trap assemblies connected via its split vent line.

### WARNING

**Inlets, detectors, and the oven are insulated with fibrous materials which may cause irritation to skin, eyes, and/or mucous membranes. Always wear gloves when working with the insulation. Additionally, if the insulation is flaky/crumbly, wear protective eyewear and a suitable breathing mask and/or respirator.**

---

- 1 Remove the top EPC module cover.
- 2 Remove the inlet cover and GC left side panel.
- 3 At the inlet, if present, disconnect the cryo line (use two wrenches against each other to prevent rotation).

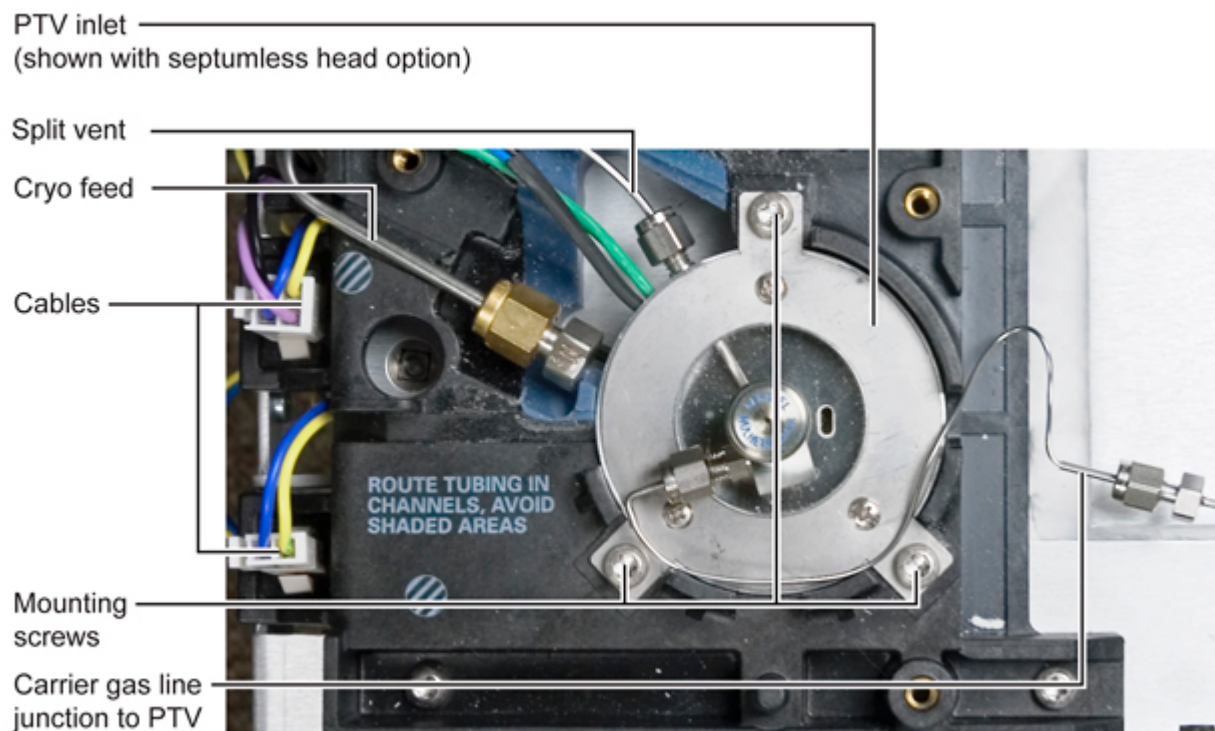**CAUTION**

When removing/attaching the cryo line use two wrenches: one wrench to support one fitting and the other to loosen/tighten the nut. Failure to do this could damage the fitting/tubing.

- 4 At the inlet, also disconnect the carrier gas line. If equipped with septum purge, disconnect the purge line from the inlet.
- 5 At the PTV EPC module, disconnect the split vent valve cable from the lower right corner of the PC board underneath the module exposed by lifting the module.
- 6 Disconnect inlet thermocouple and heater/sensor cables at the PTV heater/cryo PC board and work them back through other wiring and tubing to the inlet. Disconnect any additional cables at the inlet. Note that, if fitted, the PTV LN2 cryo feed tube from the LN2 cryo valve to the inlet may need to be entirely removed from the GC to provide room for cable connectors.
- 7 Release the inlet assembly by removing three screws and lifting the entire assembly from the GC.

- 8 In reassembly, pay attention to the following:
  - Maintain correct orientations of insulation and the inlet body.
  - Ensure that the inlet's insulation aligns properly with the tube for the cryo attachment, if present.
- 9 Reconnect PTV inlet cables.
- 10 As necessary, reconnect the split vent line, cryo line, carrier gas line, and septum purge (if present). Use new O-rings as needed.



## Replacing the PTV inlet cryo jacket

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

### WARNING

Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.

---

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.


---

### WARNING

Inlets, detectors, and the oven are insulated with fibrous materials which may cause irritation to skin, eyes, and/or mucous membranes. Always wear gloves when working with the insulation. Additionally, if the insulation is flaky/crumbly, wear protective eyewear and a suitable breathing mask and/or respirator.

---

- 3 Remove the PTV inlet assembly. See [“Replacing the PTV inlet”](#).

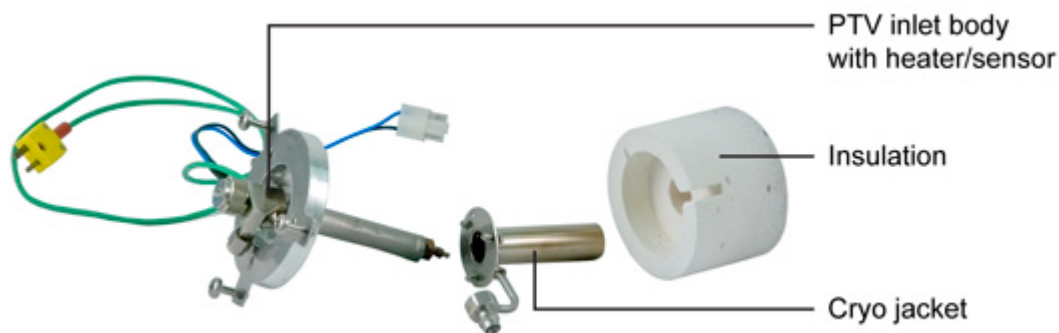

- 4 Remove three screws from the underside of the inlet body assembly securing the cryo jacket to the inlet body.

---

**NOTE**

The CO2 cryo jacket (marked "CO2" on its inlet fitting) is different from the LN2 jacket (marked "N2" on its inlet fitting).

---

- 5 Reassembly is the reverse of these steps.

---


**NOTE**

The inlet body with its heater and sensor are not individually replaceable parts.

---

## Replacing the PTV CO2 or LN2 inlet cryo assembly

- 1 Once the PTV with cooling assembly is removed, the cryo assembly on the inlet can be replaced.
- 2 Remove the PTV inlet as described under [“Replacing the PTV inlet”](#).
- 3 Remove the three screws on the cryo assembly from the PTV inlet, and remove the cryo assembly.
- 4 Install the new cryo assembly over the inlet body. Be sure to align the assembly so that the fitting is closest to the wires.
- 5 Reassembly is the reverse of these steps.
- 6 Check for leaks.

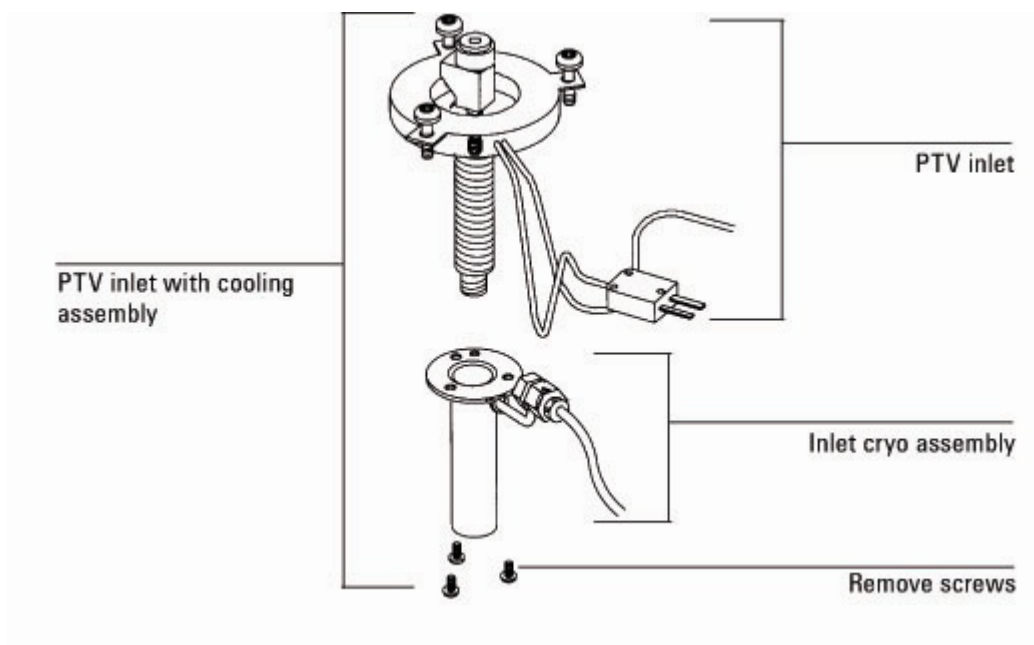

## Replacing a PTV CO<sub>2</sub> cryo valve

If changing the GC configuration, see [“Changing the GC configuration”](#) for important information regarding GC methods. Then proceed with the steps below.


- 1 Cool down the oven and all heated zones.
- 2 Turn off all gas flows at their sources.
- 3 Turn off the GC main power switch and disconnect its power cord.

### WARNING

Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

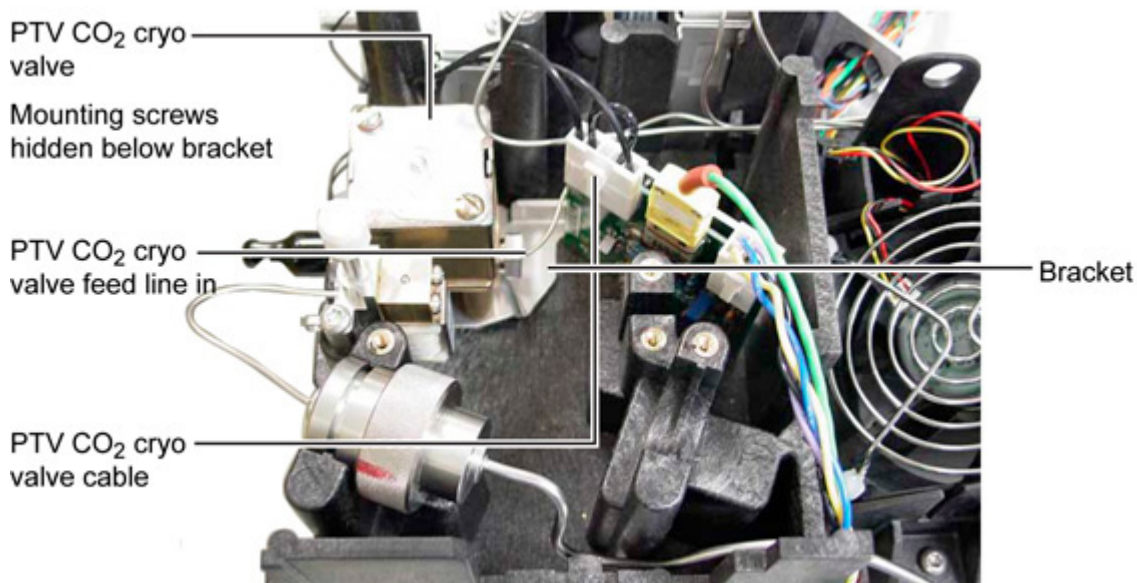

- 4 Remove the valve via its two mounting screws accessible under the removed bracket assembly. See [“Removing the PTV cryo/PC board bracket”](#) for details.
- 5 In reinstallation, make sure the valve is oriented correctly on the bracket with its cryo feed input connection towards the rear of the GC and that its cable is reconnected to the nearby PC board.

## Replacing the PTV N<sub>2</sub> cryo valve

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.


### WARNING

Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

- 3 Remove the inlet cover and GC left side panel.

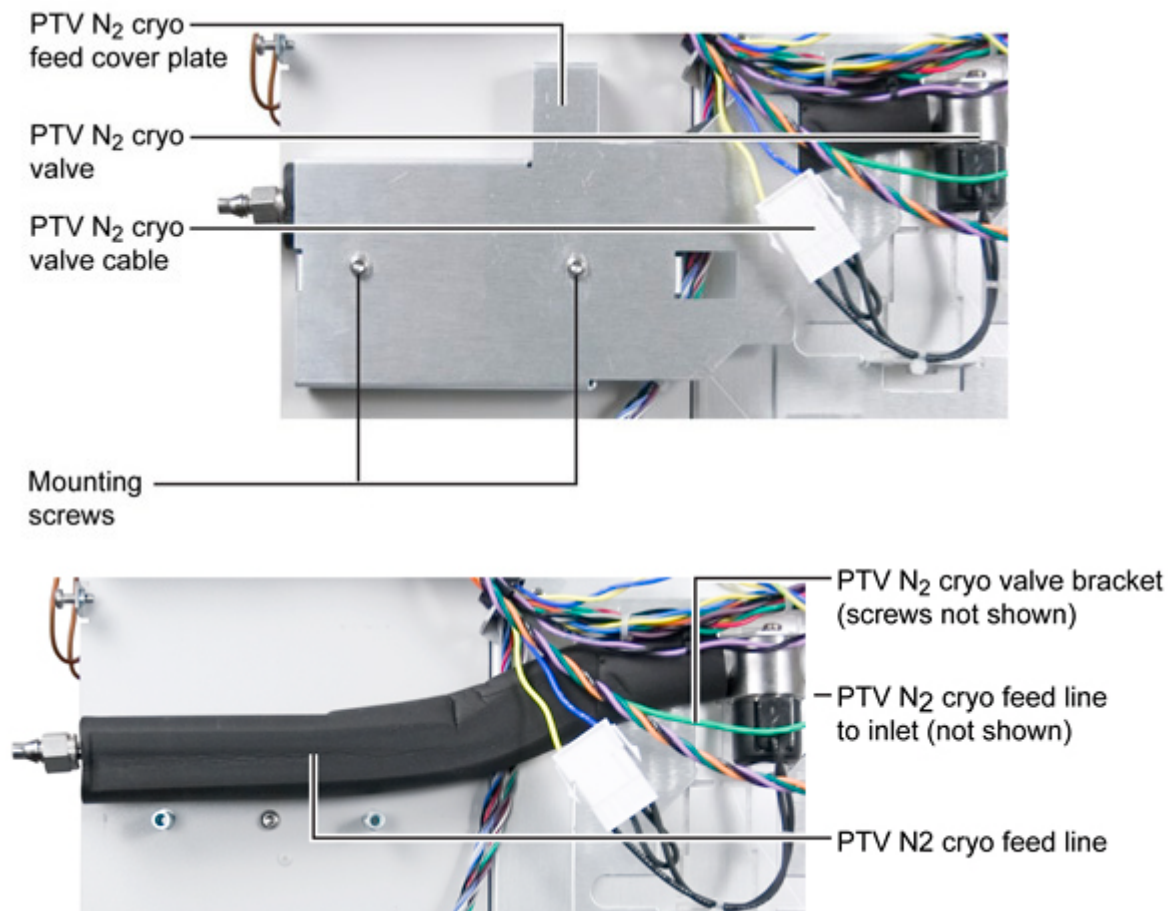

- 4 Remove the N2 cryo feed tube between the cryo valve and the PTV inlet.
- 5 Remove two screws from the PTV N2 cryo feed cover plate. Remove the cover plate.
- 6 Disconnect the PTV N2 cryo valve cable. If you intend to replace this cable with a new one, use the existing cable as a model while you install the new cable, then remove the old cable.
- 7 Remove the valve bracket and cryo valve by three screws.
- 8 Replace the N2 cryo valve assembly by reversing these steps.
- 9 Reinstall the N2 cryo tube between the valve and the inlet.

## Removing the PTV cryo/PC board bracket

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.


### WARNING

Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

- 3 Remove the top EPC module cover.

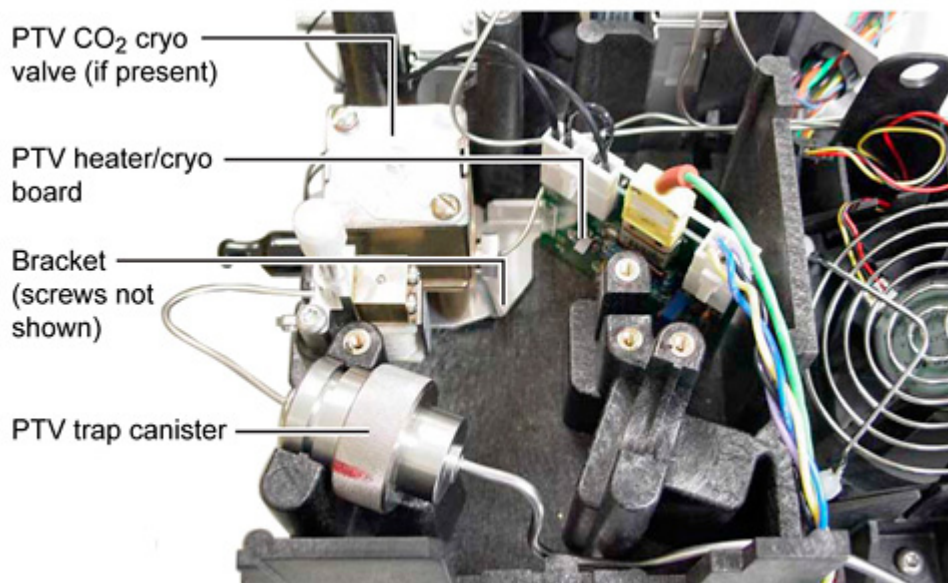

- 4 Disconnect cables to the PTV heater/cryo PC board. See [“Replacing the PTV heater/cryo PC board”](#).
- 5 If a PTV CO<sub>2</sub> valve is present, disconnect its cryo feed lines (CO<sub>2</sub> in / CO<sub>2</sub> out).
- 6 Remove two screws securing the bracket to the GC.
- 7 Lift the bracket assembly out of the GC working it around any interfering wiring and/or tubing.

**NOTE**

The PTV CO2 cryo valve, if present, is also mounted on this bracket.

---

**8** Reassembly is the reverse of these steps.



## Replacing the PTV heater/cryo PC board

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.


### WARNING

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

- 3 Remove the top EPC module cover.

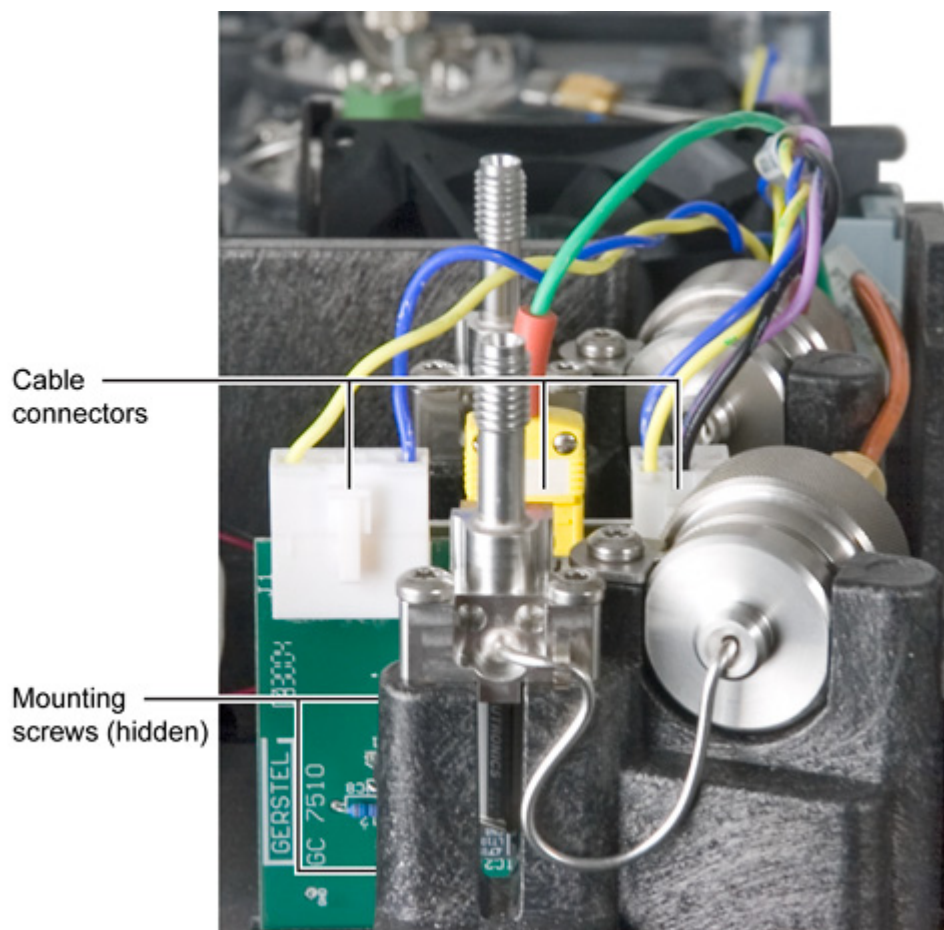

- 4 To remove the heater/cryo PC board, first remove its three cables: from left-to-right, respectively, the LN2 (liquid nitrogen) cryo valve cable (if present), the thermocouple cable, and the heater/sensor cable.
- 5 Next, remove two T-20 screws securing the board onto its bracket.
- 6 Lift the board straight up and out of the GC.
- 7 Installing the new board is the reverse of these steps.

## Volatiles Interface

The Volatiles Interface (VI) uses a controlled stream of gas to transfer a vaporized sample from an external sampler (headspace, purge and trap, thermal desorber, etc.) to a capillary column.

### Replacing or cleaning the volatiles interface

- 1 Turn off all gas flows at their sources.
- 2 Turn off the GC main power switch.

#### CAUTION


Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.


- 3 Disconnect the transfer line. Loosen the nut with a 1/4-inch wrench and remove the line. Remove the clamping plate from the interface by loosening the captive screw with a T-20 Torx screwdriver. Put the plate in a safe place.
- 4 Lift the interface out of the heater block.
- 5 If a column is installed, remove it.
- 6 Remove the split and pressure sensing lines by loosening the hex nuts with the wrench.
- 7 Clean or replace the interface. If cleaning the interface, sonicate it twice and then rinse.
- 8 Reinstall the split line and pressure sensing lines and finger tighten the hex nuts. Tighten the hex nuts an additional 1/4 turn with the wrench.
- 9 Reinstall the column in the interface.
- 10 Place the interface in the heater block. Replace the clamping plate and tighten the screw until snug. Do not overtighten.
- 11 Reinstall the transfer line. Finger tighten the nut and then tighten an additional 1/4 turn with the wrench.

#### NOTE

If the transfer line is from a TMR-8900 Purge and Trap, install the transfer line support nut assembly up and inside the metal sleeve of the heated line assembly to prevent damage to the fused silica line.

- 12 After the column is installed at both the interface and the detector, establish a flow of carrier gas through the interface and maintain it for 10 to 15 minutes. Check for leaks. Heat the interface to operating temperatures and retighten the fittings, if necessary.

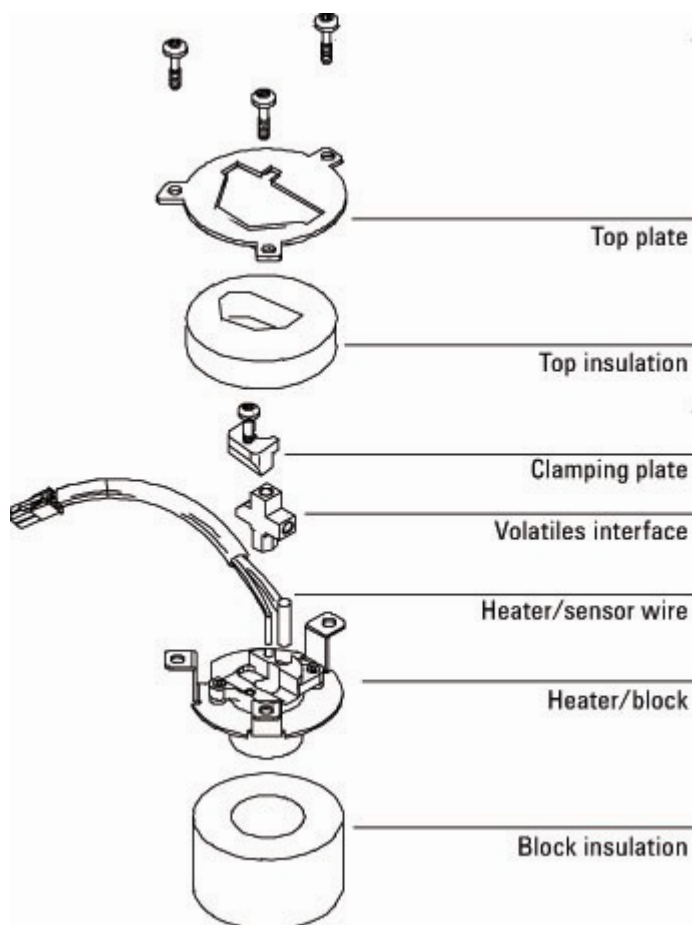

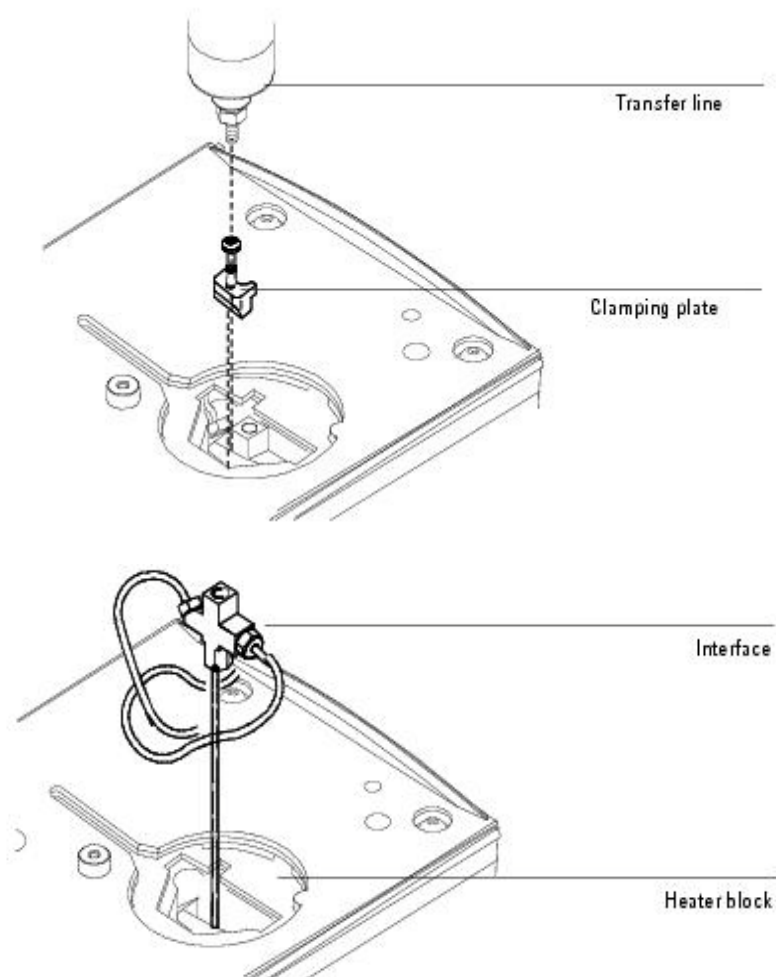

## Replacing the volatiles interface heater/sensor assembly

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

**WARNING**

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**


---

**CAUTION**

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 3 Disconnect the transfer line. Loosen the nut with a 1/4-inch wrench and remove the line.
- 4 Remove the heater/sensor wire connector from the GC connection.
- 5 Remove the three screws in the top plate of the volatiles assembly which mount it in the GC using a T-20 Torx screwdriver. Loosen each screw a little at a time.
- 6 Remove the top plate and the top insulation from the GC.
- 7 Remove the heater/sensor assembly and replace.
- 8 Reinstall the top insulation and the top plate. Align the volatiles interface with the mounting holes.
- 9 Reinstall the three Torx screws. Tighten each screw once with the T-20 Torx screwdriver until the interface is properly aligned. Tighten each screw again until snug.
- 10 Reinstall the heater/sensor wire connector in the GC. Route the wire as shown.
- 11 Reinstall the transfer line using a 1/4-inch wrench and check for leaks.

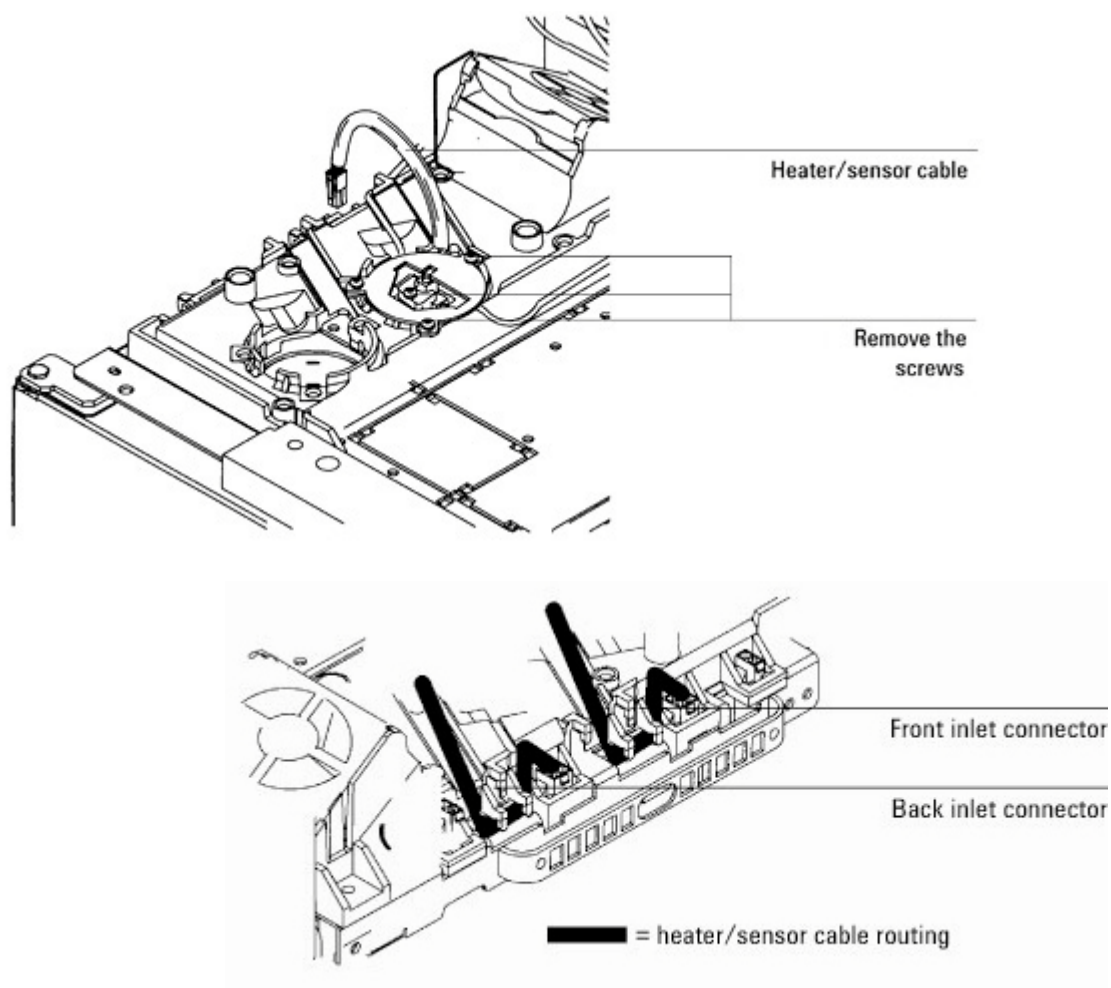

**Figure 1** Cable routing

## Calibrating the interface

The interface's flow module contains a pressure sensor that must be zeroed after it is installed on your GC. This calibration procedure ensures an accurate interface pressure display.

Do not connect the carrier gas to the flow module until you have zeroed the interface's pressure sensor.

- 1 Plug in the GC and turn it on, if you haven't already done so.
- 2 Wait 15 minutes. This allows the GC to reach thermal equilibrium.
- 3 Zero the interface's pressure sensor:
  - a Press **[Options]**, scroll to **Calibration** and press **[Enter]**.
  - b Scroll to the module to be zeroed and press **[Enter]**.
  - c Scroll to a zero line and press **[Info]**. The GC will remind you of the conditions necessary for zeroing that specific sensor.

**Flow sensors.** Verify that the gas is connected and flowing (turned on).


**Pressure sensors.** Disconnect the gas supply line at the back of the GC. Turning it off is not adequate; the valve may leak.

- d Press **[On/Yes]** to zero or **[Clear]** to cancel.
- 4 Turn off the GC.
- 5 Plumb the carrier gas to the flow module. If you need help with this step, see the GC site preparation and installation manuals.
- 6 Turn on the GC.
- 7 Configure the column and carrier gas.



## Solvent Vapor Exit Accessory

This accessory vents the carrier gas and solvent vapor to waste for a controlled time after injection. This greatly reduces the size of the solvent peak.

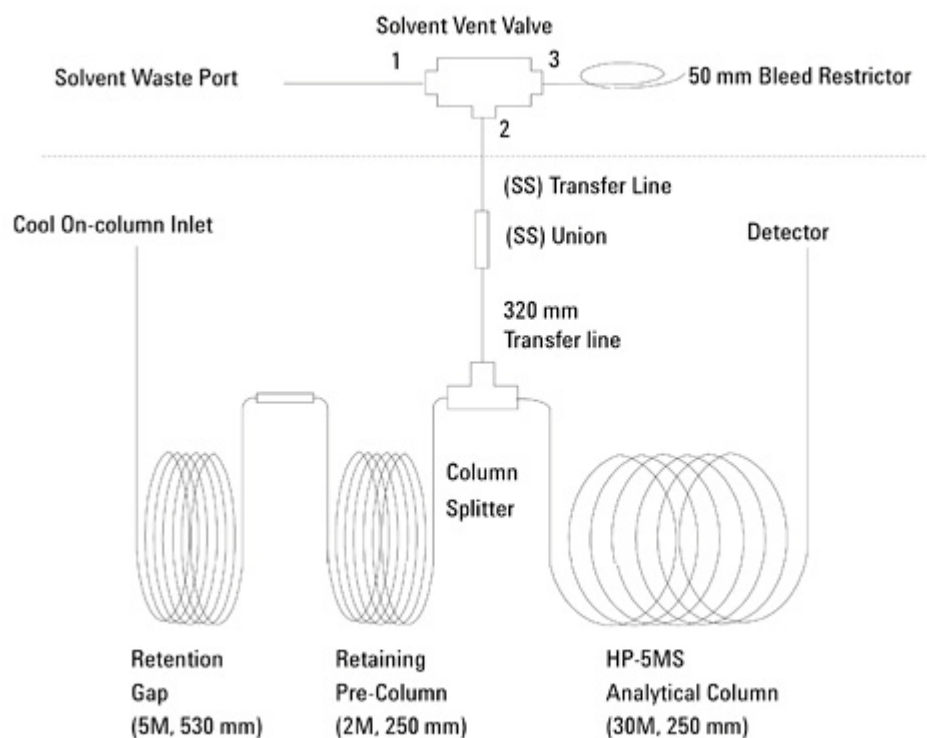

## Replacing the SVE valve/fitting assembly

**WARNING**

Hazardous voltages are present in the mainframe when the GC power cord is plugged in. Avoid a potentially dangerous shock hazard by unplugging the power cord before removing the side panels.

---

**CAUTION**

Prevent electrostatic voltages from damaging the GC by using an ESD wriststrap.


---

- 1 Turn off the GC and unplug the power cord. Allow time for all heated zones to cool and then turn off supply gases at their sources.
- 2 Remove the top cover, the pneumatics cover, the electronics carrier cover, and the right side cover.
- 3 Remove the bleed restrictor column and the solvent vent waste line from the valve fitting assembly.
- 4 Inside the oven, disconnect the transfer line from the stainless union on the 1/16-inch stainless steel tubing. Carefully remove the tubing from the oven through the cutout in the top of the oven.
- 5 Disconnect the valve driver cable from the valve jumper cable or valve driver assembly, as applicable.
- 6 Remove the mounting screw in the valve/fitting assembly and remove the assembly from the GC.
- 7 Install the new valve/fitting assembly and mount it to the GC oven with the screw. Cover the open end of the union to avoid contamination and route the 1/16-inch stainless steel tubing and union into the oven. Repack the insulation around the cabling. Connect the valve driver cable to the valve jumper cable or valve driver assembly, as applicable.
- 8 Use a new graphite/Vespel ferrule to reconnect the transfer line to the union.
- 9 Examine the bleed restrictor column. If the column is damaged, replace it with a new 0.5 m length of 0.050 mm column, installing a new fitting and ferrule. Be sure to trim 5 to 10 mm off the end of the new column after installing the new fitting and ferrule.
- 10 Reattach the solvent waste vent line.

**WARNING**

Because a significant amount of solvent is vented through the SVE valve assembly, it is important that the bleed restrictor and the solvent vent are connected to an appropriate laboratory ventilation system.

11 Check for leaks.

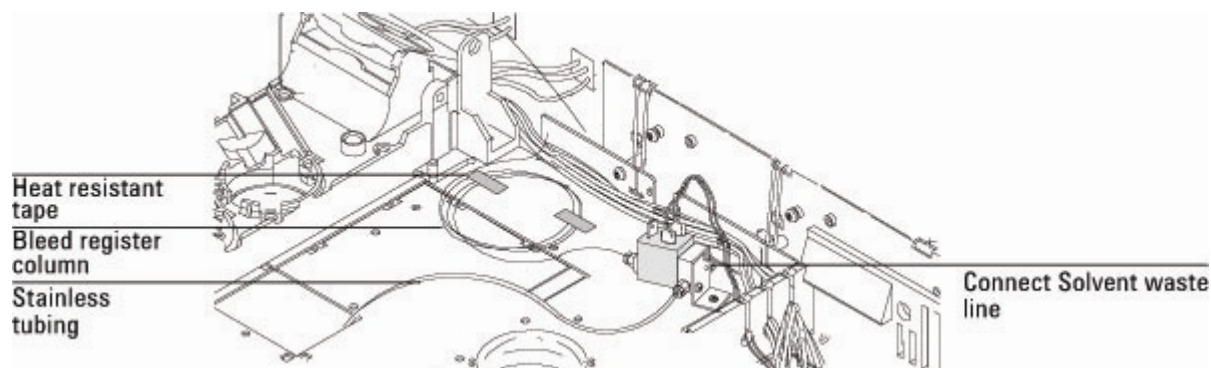

## **Replacing the SVE bleed restrictor column**

- 1** Turn off the GC and unplug the power cord. Allow time for all heated zones to cool and then turn off supply gases at their sources.
- 2** Remove the top cover and the right side cover. If necessary, also remove the pneumatics cover and the electronics carrier cover.
- 3** Using a wrench, remove the old bleed restrictor column.
- 4** Cut a 0.5 m section off of the new 50  $\mu$ m bleed column for use as the restrictor.
- 5** Attach a male fitting and ferrule to the restrictor column, then trim 5 to 10 mm from that end of the column.
- 6** Connect the 0.5 m x 50  $\mu$ m bleed restrictor to the SVE valve.
- 7** Make a loop (or loops) in the excess column protruding from the valve and move it to an unobstructive position.
  - Coil it in a loose coil (approximately 6 inch diameter).
  - Secure the coil with heat resistant tape or equivalent.

## Replacing the SVE tri-column assembly

- 1 Turn off all gas flows at their sources.
- 2 Turn off the GC main power switch.

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 3 Disconnect the Tri-column assembly from the COC inlet, detector and (SS) union and remove the column assembly from the GC oven.
- 4 Install the new SVE Tri-column assembly onto a column hanger and hang the assembly inside the GC oven. Position the column assembly so that the end of the retention gap (530  $\mu\text{m}$ ) is located under the COC inlet and the end of the analytical column (250  $\mu\text{m}$ ) is under the detector.
- 5 Connect the retention gap to the COC inlet using a graphite ferrule and a column nut.
- 6 If using a MSD, connect the analytical column (HP-5MS) to the MSD using a column nut (part no. 05988-20066) and a graphite/Vespel ferrule (part no. 5062-3508). If using any other type of detector, connect the analytical column (HP-5MS) to the GC detector using a column nut (part no. 5181-8830) and a graphite/Vespel ferrule (part no. 5062-9527).
- 7 Remove the stainless ferrules from the union. Use a graphite/Vespel ferrule to connect the 0.5 m  $\times$  320  $\mu\text{m}$  transfer line from the quartz Y-splitter of the Tri-column assembly to the unused end of the stainless union located inside the GC oven.
- 8 Check all of the connections for leaks using an electronic leak detector.

## Replacing the SVE pre-column assembly

- 1 Turn off all gas flows at their sources.
- 2 Turn off the GC main power switch.

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---


- 3 Disconnect the Tri-column assembly from the COC inlet, the detector and the stainless steel union and remove the column assembly from the GC oven.
- 4 Using a column cutter, cut the transfer column and the analytical columns as close to the quartz Y-splitter as possible. Be sure to cut the columns straight.
- 5 Using the column cutter, trim the ends of the tri-column assembly ensuring that the cuts are straight. Then use a lint free wipe and methanol to clean any dirt and finger prints from approximately three to four centimeters from each of the five column ends.
- 6 Being careful not to touch the ends of the columns, insert them into the quartz splitter as far as possible, making sure that the column ends are making contact with the inside walls of the splitter. Sufficient pressure should be applied to give a good seal. Too much pressure, on the other hand, can damage the polyimide layer or even the column and result in leakage. Test to see that the column has been installed correctly by trying to pull the column out of the splitter using medium force. If the column comes out easily, trim the end and insert again.

### NOTE

A concentric circle inside the splitter should result if the column is installed correctly.

---

- 7 After connecting the column ends into the quartz splitter, reconnect the tri-column assembly inside the GC oven and increase the column head pressure incrementally to the desired pressure. An incremental increase is better than immediately applying the total column head pressure to avoid a "shock" that can loosen the connection.

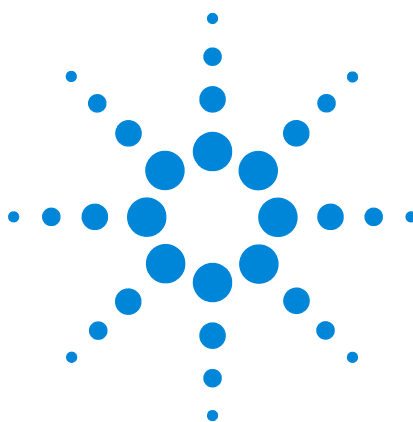

## 3 Detectors

Detector overview	100
About the detector signal boards	100
Accessing detector signal boards	103
Replacing detector signal boards	104
Detector signal board, bad checksum error	107
Flame Ionization Detector (FID)	113
Nitrogen Phosphorus Detector (NPD)	132
Thermal Conductivity Detector (TCD)	148
Micro-cell Electron Capture Detector (uECD)	152
Flame Photometric Detector (FPD) 7890A	160



## Detector overview

A detector monitors the gas stream exiting from the analyzing column. Its electrical output changes when the composition of the gas does. This section deals with the most widely used detectors.

## About the detector signal boards

The GC can control up to 4 detector signal boards and simultaneously process their digital signals. Except for a TCD, all detector boards are mounted in the electronics carrier on the right-hand side of the GC. A TCD signal board can be mounted in the optional 3rd detector carrier.

Sometimes other problems seem like a signal board failure. Replace a signal board only after ruling out other possible causes as follows:


- 1 Verify the signal board is receiving 24 VDC power, the green LED is ON, and the board appears configured in the GC display.
- 2 Check diagnostics for faults caused by components connected to the signal board, such as, ignitors, electrometers, filaments, and PMTs.
- 3 Replace the signal board only after eliminating the possible causes listed in steps 1 and 2.

Repairs should be verified by the following:

- Typical detector signal baseline output and noise.

There are 6 types of detector signal boards. Their part numbers and possible locations are listed in the table below.



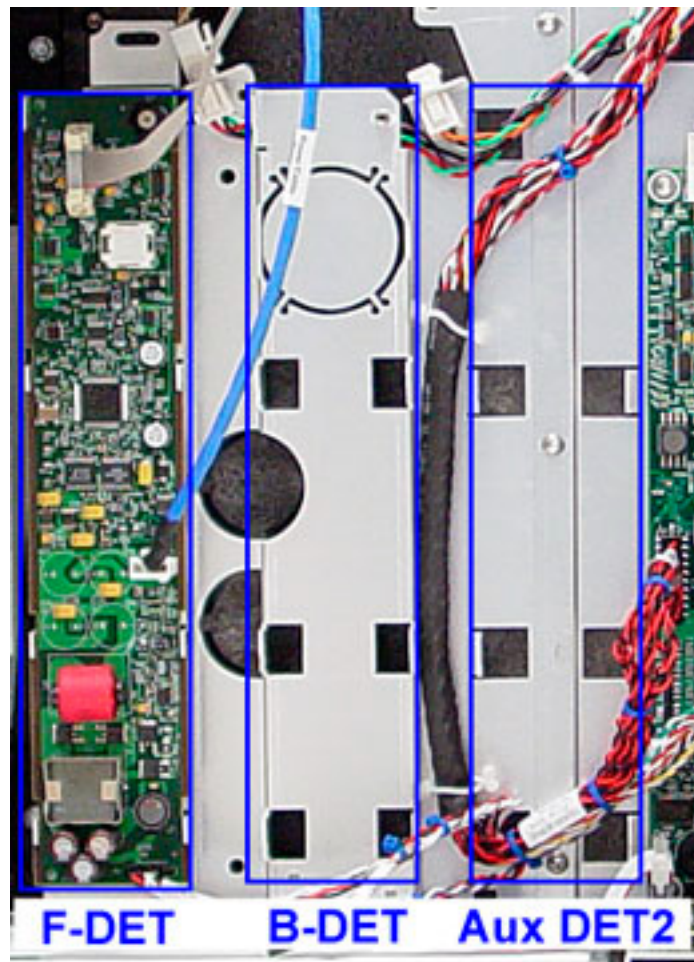

**Table 2** Detector signal boards

Description	Part number	Slot
AIB	G3456-60010	Aux DET1
uECD	G3433-61020	F-DET or B-DET
FID	G3431-60020	F-DET or B-DET
FPD	G3435-60010	F-DET or B-DET or Aux DET2
NPD	G3434-60020	F-DET or B-DET
TCD	G3432-60010	F-DET or B-DET or Aux DET1

**Table 3** Detector signal electrical connectors

Description	P1	P2	P3	P4	Other
AIB	Communication buss			J1 Signal input	P5
uECD	Communication buss	Electrometer			
FID	Communication buss	Ignitor		Electrometer	
FPD	Communication buss	Ignitor		Signal input	P7 S11 High voltage source for PMT
NPD	Communication buss	Bead current		Electrometer	
TCD	Communication buss		Switching valve	J4 Filaments	

## Accessing detector signal boards

The front and back detector signal boards are located under the right-hand side panel.

- 1 Loosen the T-20 captured screws until they do not engage the electronics carrier.
- 2 Slide the panel toward the back of the GC and lift off. Place the tip of the driver in one of the top ventilation slots and push back.

The 3rd detector signal board is located in the 3rd detector carrier on the left side of the GC.

- 1 Remove the 2 T-20 screws from the side panel.
- 2 Lift the panel (both side and top) off of the carrier.
- 3 Remove the screw from the detector signal board cover. Remove the cover.

## Replacing detector signal boards

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

---

**WARNING**

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**


---

**CAUTION**


Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 3 Remove the right-hand side panel for the electronics carrier.

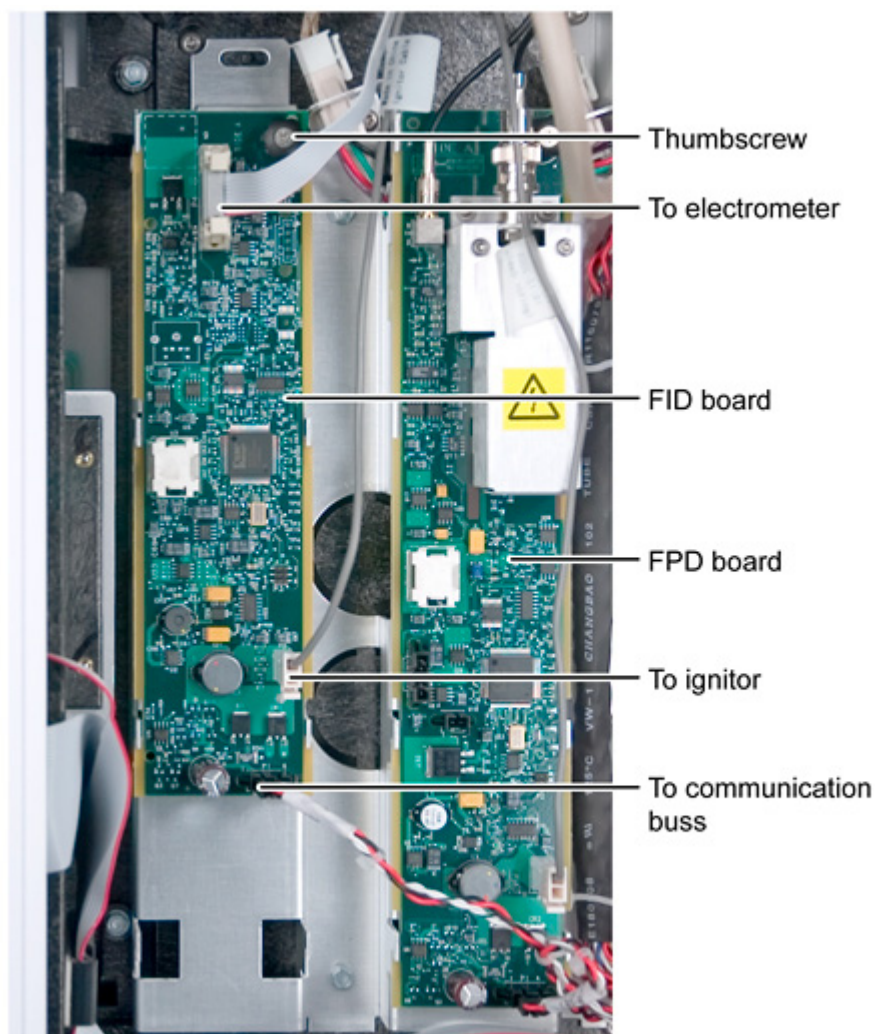

- 4 Disconnect the electrical cables to the board.
- 5 Turn the thumb screw counter-clockwise and lift the screw to disengage it from the sheet metal.

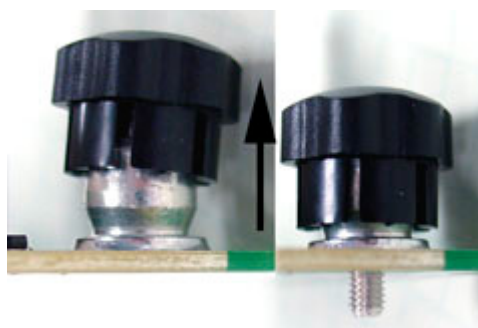

- 7 When replacing the signal board, short boards, such as the FID and uECD, slide into 3 slots. Long boards, such as NPD and FPD, slide into 4 slots.
- 8 Push in on the thumbscrew to engage the threads. Turn the screw clockwise until it is tight.

**CAUTION**

The grounding point for the signal board is the thumbscrew. If this screw is not secured, the detector signal will be noisy and may cause damage to the board.

---

- 9 Connect the electrical cables to the board.
- 10 Turn ON the GC and reconfigure the detector signal board. See [“Replacing a GC module”](#) for details.

## Detector signal board, bad checksum error

It is very possible that you will get error messages at power on. The messages can be a fault or exception. They are often due to a checksum error with the detector signal boards. Use the following instructions to add a serial number and manufacturing date to the signal board eeprom.

Detector signal boards are programmed and labelled with a unique ID number along with a manufacturing date. The following is an example of an FPD signal board. The serial number is **STI190698769**.



In some rare occasions, the board ID number or manufacturing date is not programmed into the board's EEPROM. In these cases, the GC will post an exception or a "bad checksum error." In this state, the detector is not usable.

Before you start, you can confirm the problem by trying to view the information. Navigate to **Service Mode > Diagnostics > Front Detector > Signal board**. At the bottom of the table, you will see the manufacturing date and ID number.


You can perform the following procedure to insert the ID number and/or date into the EEPROM of the board. You will need:

- PC with a telnet program, such as, **HyperTerminal**.
- LAN cable or cross-over cable.
- LAN connection to the GC and IP address for this GC.
- The ID number of the board.


### Procedure

- 1 Power off the GC, remove the detector signal board, and record the serial number.
- 2 Reinstall the detector signal board and power ON the GC.

- 3 If necessary, configure the GC and the PC for a LAN connection.
- 4 Unlock the keyboard.
- 5 Launch **HyperTerminal** or other telnet program.
- 6 Name the connection and click **OK**.

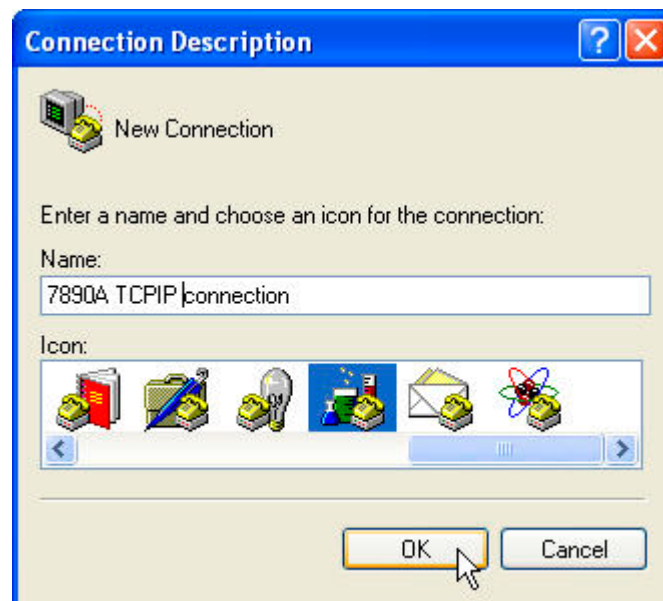


- 7 Select TCP/IP (Winsock) for the connection.
- 8 Enter the IP address of the instrument and a Port number of **9101**.

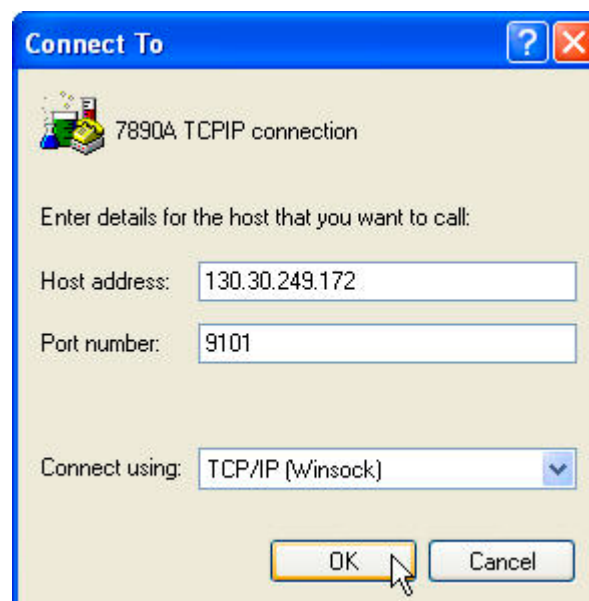


- 9 Under the Connection properties select ASCII Setup. Check the items shown below.



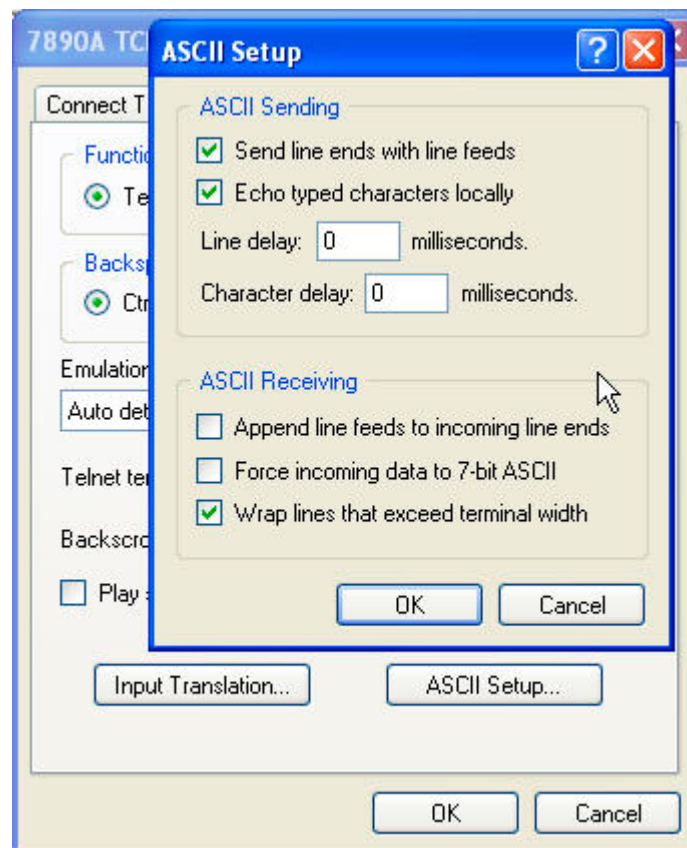

- 10 The GC will automatically reply with this prompt. Enter **go** and press Enter.

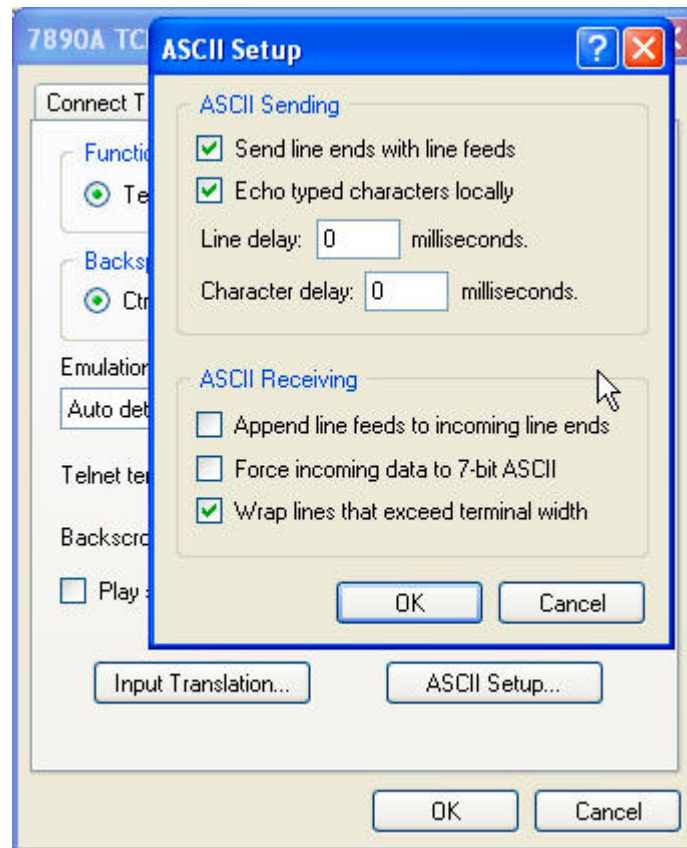

- 11 Type **go** and press Enter.
- 12 For an ID number, enter the following string and press enter. The memory location is 1.


```
dvcLES<space>< 9 or 8 or 7 or 6
><space>1<space>(string)<serial number>
```

- 9 for FRONT DET
- 8 for BACK DET
- 7 for AUX DET2
- 6 for AUX DET1 or 3rd detector

- 13 For a manufacturing date, enter the following string and press enter. The memory location is 7.

```
dvcLES<space><9 or 8 or 7 or
6><space>7<space>(int)1134699051
```

If the process is successful, the GC will respond with a 0.

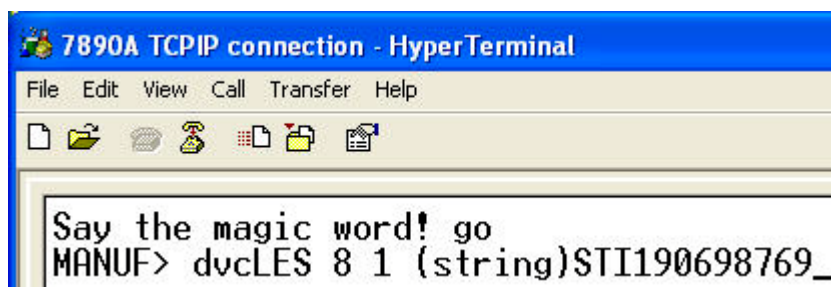

- 14 After writing the information to the eeprom, confirm the ID number and date by entering these strings and pressing enter.


dvcRES<space>< 9 or 8 or 6 ><space>1<space>(string) and

dvcRES<space><9 or 8 or 7 or 6><space>7<space>(int)

```
MANUF> dvcRES 8 1 (string)
dvcRES 0 1 STI190698769
MANUF>
```

The next steps reconfigure the detector signal board.

- 15 Press [Config] [Front Det] or [Back Det].

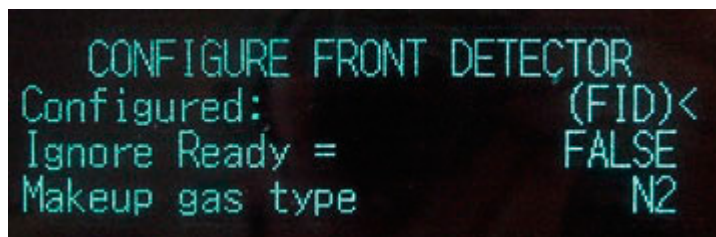

- 16 Press [Mode/Type] and press [Enter].




- 17 The GC displays a message prompting to power cycle. Do not power cycle.

- 18 Press [Mode/Type] again and press [Enter].



- 19 Power cycle the GC. Your detector is reconfigured.

If you do not perform steps 14 through 18, the GC will display the following messages at power ON.

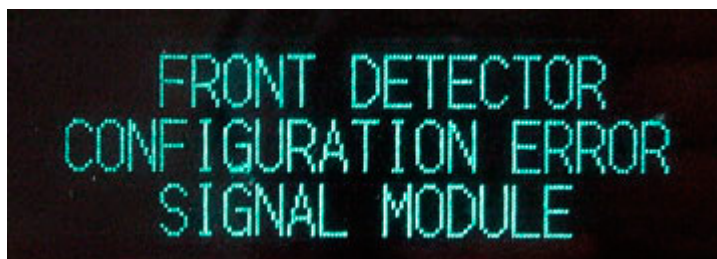

FRONT DETECTOR  
CONFIGURATION ERROR  
SIGNAL MODULE

and




FRONT DETECTOR (FID) BROKEN  
F-DET (OK) 68 watts 0.0<  
Signal: (OK) FID  
Epc: (OK) FID

## Flame Ionization Detector (FID)

The FID is the most widely used of the GC detectors. It responds to almost all organic compounds (there are a few exceptions), has good sensitivity and a wide linear range, and is easy to use.

### Selecting an FID jet

Open the oven door and locate the column connection fitting at the base of the detector. It will look like either a capillary optimized fitting or an adaptable fitting.

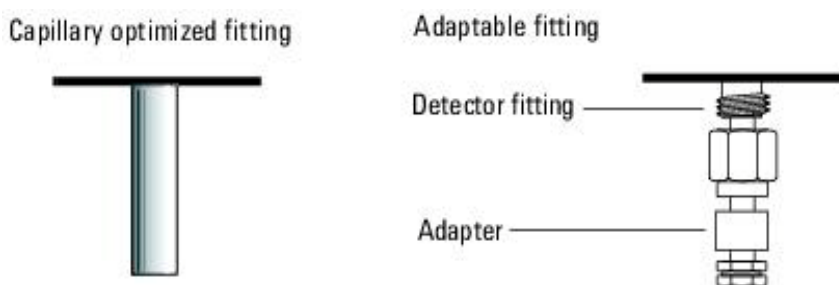

- If you have an application that tends to clog the jet, select a jet with a wider tip id.
- When using packed columns in high column-bleed applications, the jet tends to clog with silicon dioxide.
- In simulated distillation applications, the high-boiling hydrocarbons tend to clog the jet.

For capillary optimized fittings, select a jet from the table below, [Jets for capillary optimized fittings](#).

**Table 4** Jets for capillary optimized fittings

Figure 2 ID	Jet type	Part number	Jet tip id	Length
1	Capillary	G1531-80560	0.29 mm (0.011 inch)	48 mm
2	High-temperature (use with simulated distillation)	G1531-80620	0.47 mm (0.018 inch)	48 mm




**Figure 2** Capillary optimized jets

For adaptable fittings, select a jet from the table below, [Jets for capillary adaptable fittings](#).

**Table 5** Jets for capillary adaptable fittings

Capillary adaptable jets ID	Jet type	Part number	Jet tip id	Length
1	Capillary	19244-80560	0.29 mm (0.011 inch)	61.5 mm
2	Capillary, high- temperature (use with simulated distillation)	19244-80620	0.47 mm (0.018 inch)	61.5 mm
3	Packed	18710-20119	0.46 mm (0.018 inch)	63.6 mm
4	Packed, wide-bore (use with high-bleed applications)	18789-80070	0.76 mm (0.030 inch)	63.6 mm

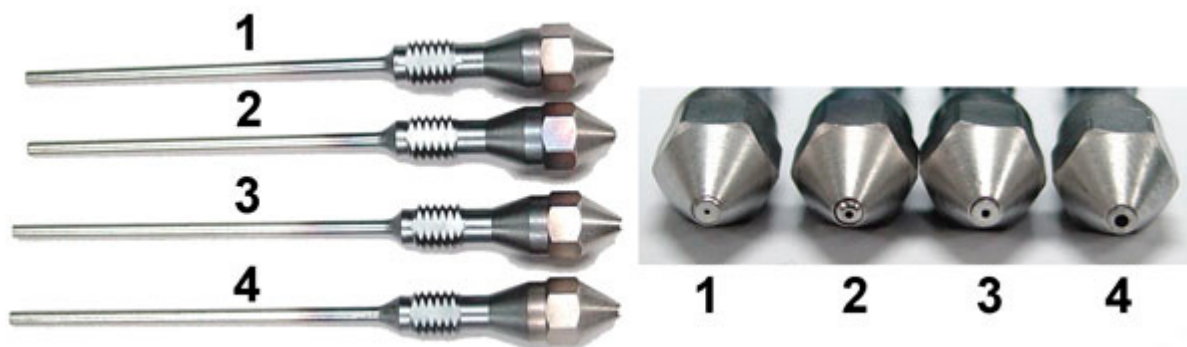

**Figure 3** Capillary adaptable jets

## Replacing the FID interconnect assembly or spring

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

### WARNING

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

---

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 3 Gather the following:
  - New FID collector assembly. (See [“PM kits for the FID”](#).)
  - T-20 Torx screwdriver
  - 1/4-inch nut driver
  - Tweezers
  - Lint-free gloves
- 4 Shut off the detector and the detector gases and let the detector cool.


### CAUTION

To avoid contaminating the FID, wear clean, lint-free gloves when handling the collector assembly.


---



- 5 Disconnect the ignitor cable assembly.

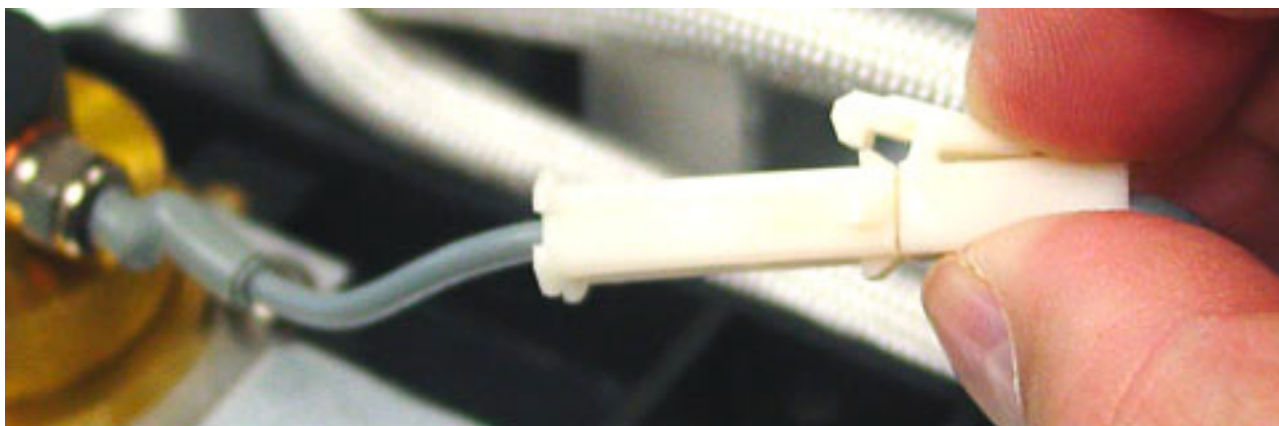

- 6 Remove the three screws holding the collector assembly to the mounting pallet.

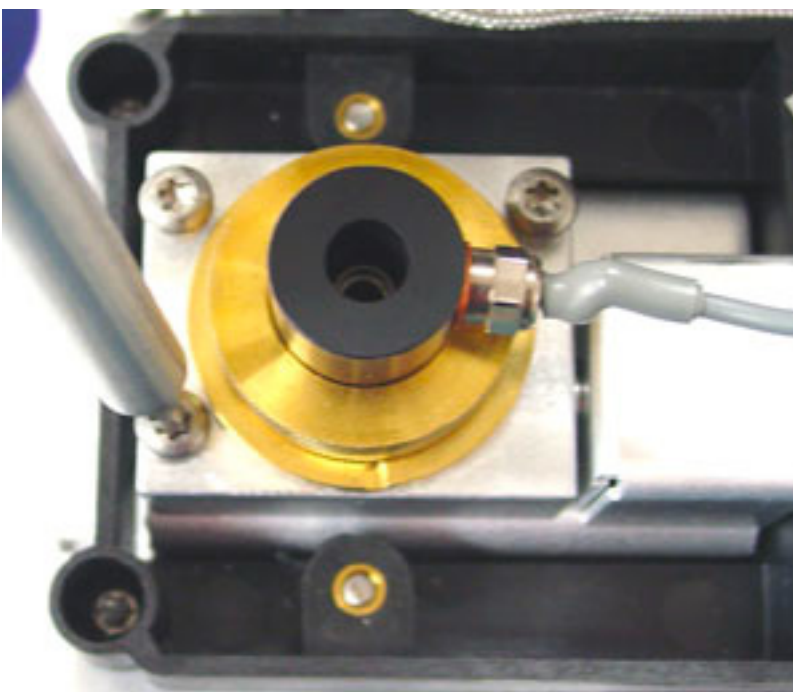

**CAUTION**

The next step exposes the interconnect spring. Be careful not to touch or disfigure the spring while working on the FID. Any dirt or bending will reduce the sensitivity of your detector.

---



- 7 Lift and remove the assembly from the pallet.



#### **To replace the interconnect spring only**

- 1 Pull the spring off the end of the interconnect. Push a new spring on.
- 2 Reassemble the detector.

#### **To replace the interconnect (and spring)**

- 1 Remove the screws at each end of the electrometer body. DO NOT loosen the screw in the center of the electrometer cover.
- 2 The interconnect is held by a small spring clip. Remove the screw holding this clip.
- 3 Disconnect the ribbon cable that connects the electrometer to the signal board. Remove the electrometer.
- 4 There is a hex section where the interconnect enters the electrometer. Use a wrench to loosen and remove the interconnect.
- 5 Insert a new interconnect and tighten it firmly (do not overtighten!).
- 6 Reassemble the detector.

## Replacing the heater

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

### **WARNING**

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

---

### **CAUTION**

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 3 Turn off all gas flows at their sources.


### **WARNING**

**Inlets, detectors, and the oven are insulated with fibrous materials which may cause irritation to skin, eyes, and/or mucous membranes. Always wear gloves when working with the insulation. Additionally, if the insulation is flaky/crumbly, wear protective eyewear and a suitable breathing mask and/or respirator.**

---


- 4 Disconnect the column from the bottom of the detector.
- 5 Remove the right side electronics panel.
- 6 Access the FID module by unscrewing the single T-20 screw holding the GC cover over the detector(s).
- 7 Raise it out of the way or remove it temporarily.

- 8 Disconnect the heater/sensor cable from the heater/sensor connector.

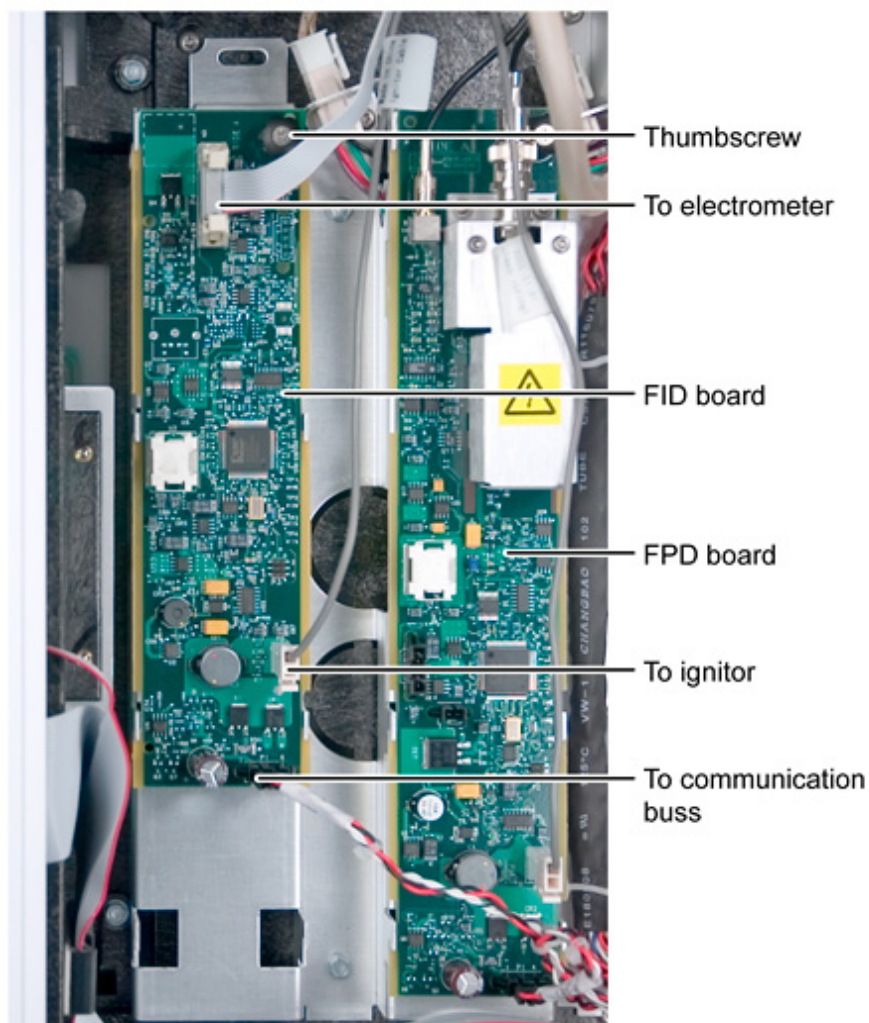


- 9 Disconnect the ignitor cable from the FID signal board to the ignitor castle.
- 10 Disconnect the electrometer cable from the FID signal board.

- 11** Unscrew the three captive T-10 screws retaining the FID flow blocks to the EPC module.

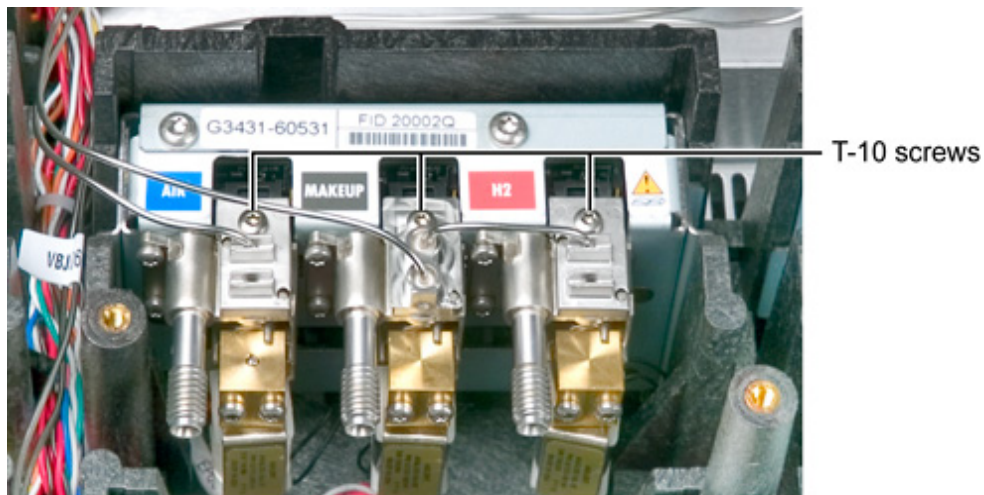


- 12** Remove the four T-20 screws retaining the FID assembly to the GC.

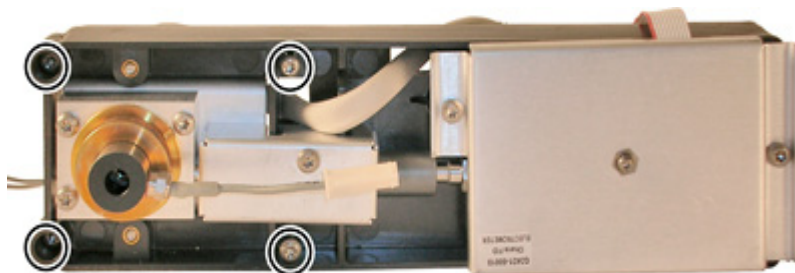

- 13** Lift the entire FID assembly up and off.  
**14** Lay the assembly down on a clean flat surface.  
**15** Examine the underside of the detector assembly.

**16** Remove the insulated cover over the heater block.

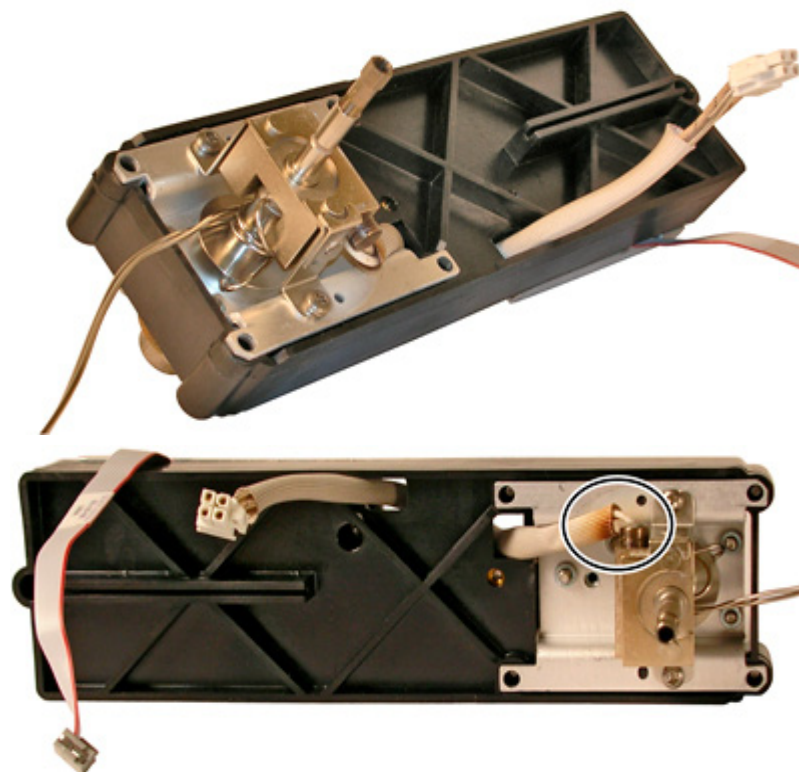

**17** Loosen the heater/sensor cable sufficiently to allow for the lifting and removal of the heater from the heater block.

**18** Carefully remove the defective heater/sensor cable and replace it with a new heater/sensor cable.

**19** Reassembly is the reverse of these steps.

## Replacing the FID signal board

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

**WARNING**

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

---

**CAUTION**

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 3 Turn off all gas flows at their sources.

**WARNING**

**Inlets, detectors, and the oven are insulated with fibrous materials which may cause irritation to skin, eyes, and/or mucous membranes. Always wear gloves when working with the insulation. Additionally, if the insulation is flaky/crumbly, wear protective eyewear and a suitable breathing mask and/or respirator.**

---

- 4 Remove the right side electronics panel from the GC.
- 5 Lift up the GC detector cover to access the FID.
- 6 Disconnect all electrical connections between the FID and the signal board.
- 7 Rotate the FID board thumbscrew counterclockwise, lift the board up and then downward to remove it from the GC.
- 8 Replace the defective board with a new FID board.
- 9 Restore the FID electrical connections.
- 10 Reassembly is the reverse of these steps.
- 11 Reconfigure the GC using the front panel keypad, supplying the requested information.

## Replacing the FID electrometer

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

### WARNING

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

---

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 3 Turn off all gas flows at their sources.

### WARNING


**Inlets, detectors, and the oven are insulated with fibrous materials which may cause irritation to skin, eyes, and/or mucous membranes. Always wear gloves when working with the insulation. Additionally, if the insulation is flaky/crumbly, wear protective eyewear and a suitable breathing mask and/or respirator.**

---

- 4 Remove the right side electronics cover.
- 5 Disconnect the electrometer ribbon cable from the FID signal board.
- 6 Remove the Torx T-20 screw and clamp over the interconnect tube.
- 7 Remove the two Torx T-20 screws from each end of the electrometer.
- 8 Lift the electrometer up and away.



### 3 Detectors

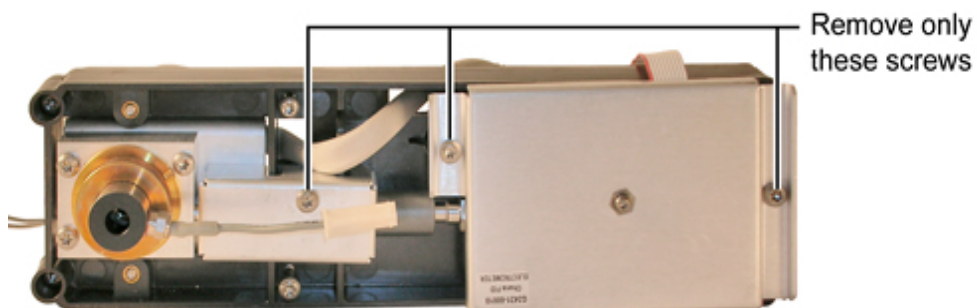

- 9 Replace the defective electrometer with a new one, making sure that there is firm contact between the interconnect spring and the collector.
- 10 Reassembly is the reverse of these steps.
- 11 After reassembly, turn the GC power on and reconfigure the modified detector using the front panel keypad, supplying the requested information.



## Replacing the FID glow plug (ignitor)

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

### **WARNING**

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

---

### **CAUTION**

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.


---

### **WARNING**


**Inlets, detectors, and the oven are insulated with fibrous materials which may cause irritation to skin, eyes, and/or mucous membranes. Always wear gloves when working with the insulation. Additionally, if the insulation is flaky/crumbly, wear protective eyewear and a suitable breathing mask and/or respirator.**

---

- 3 Remove the right side electronics panel and identify the FID signal board.
- 4 Disconnect the ignitor cable from the FID signal board.

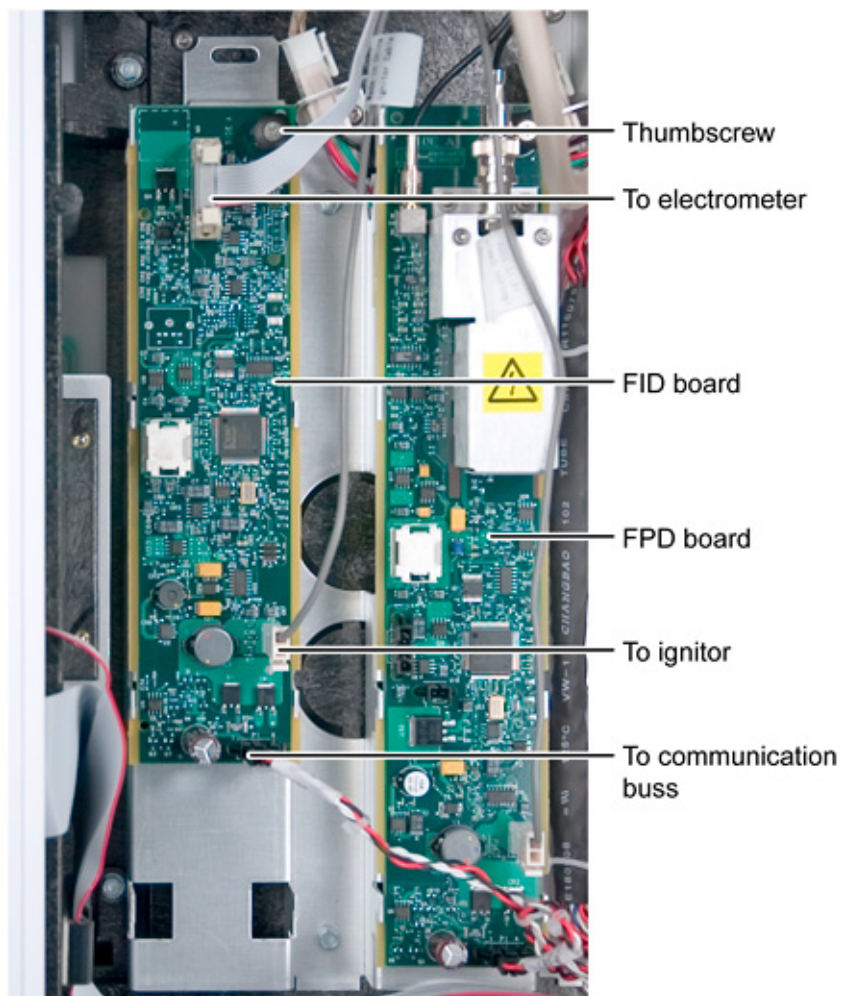

- 5 Lift up or remove the GC detector cover to access the FID.
- 6 Unclip the ignitor cable from the cable extension to the FID board.
- 7 From the top, unscrew the defective ignitor cable from the ignitor castle.
- 8 Replace the defective cable with a new ignitor cable.
- 9 Reassembly is the reverse of these steps.

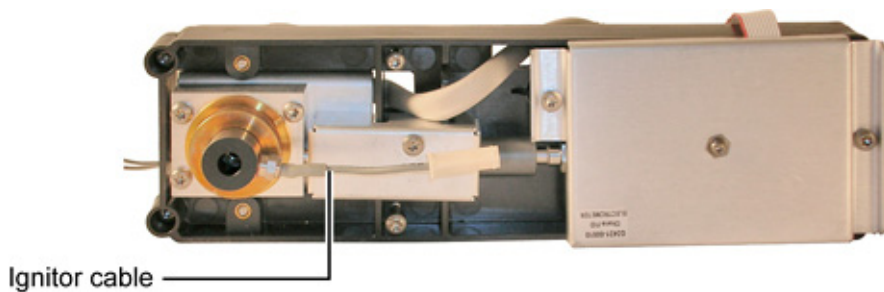

## Replacing the entire FID

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

### WARNING

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

---

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 3 Turn off all gas flows at their sources.


### WARNING

**Inlets, detectors, and the oven are insulated with fibrous materials which may cause irritation to skin, eyes, and/or mucous membranes. Always wear gloves when working with the insulation. Additionally, if the insulation is flaky/crumbly, wear protective eyewear and a suitable breathing mask and/or respirator.**


---

- 4 Disconnect the column from the bottom of the detector.
- 5 Remove the right side electronics panel.
- 6 Access the FID module by unscrewing the single T-20 screw holding the GC cover over the detector(s).
- 7 Raise the cover out of the way or remove it temporarily.
- 8 Disconnect the heater/sensor cable to the FID.
- 9 Disconnect the ignitor cable from the FID signal board.
- 10 Disconnect the electrometer cable from the FID signal board.


### 3 Detectors

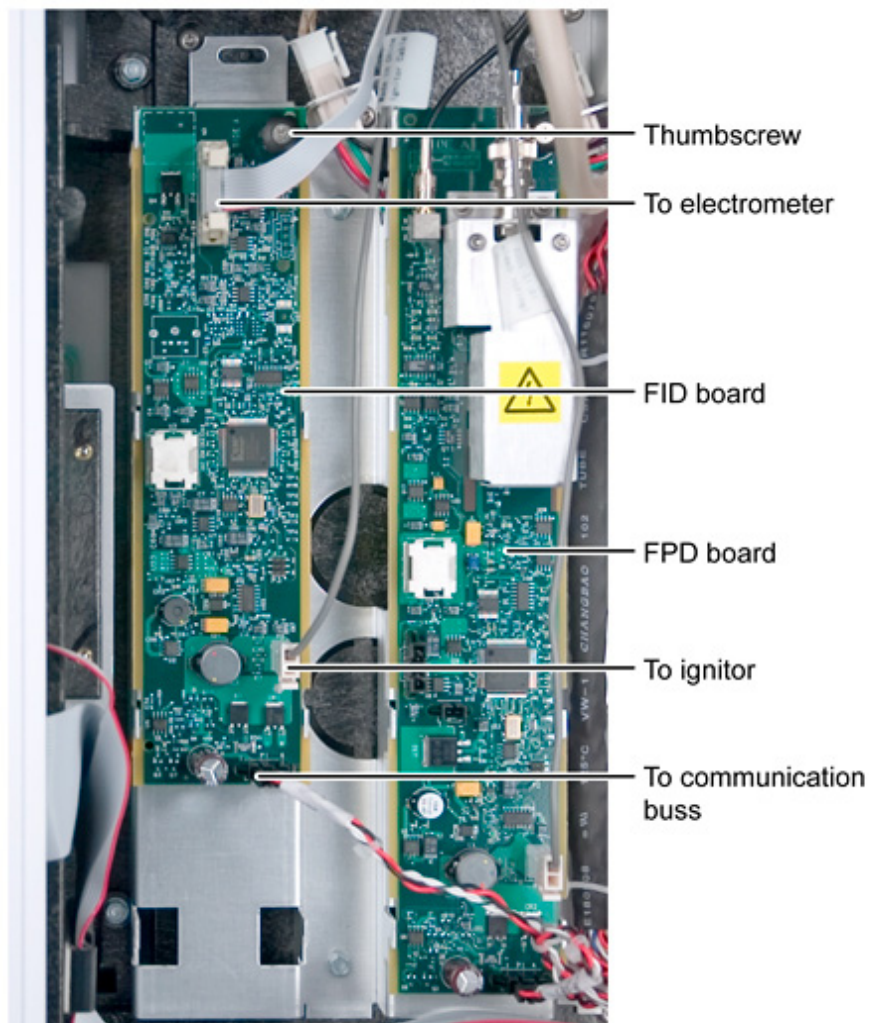

**11** Unscrew the three captive T-10 screws retaining the tubing connection blocks to the EPC module.

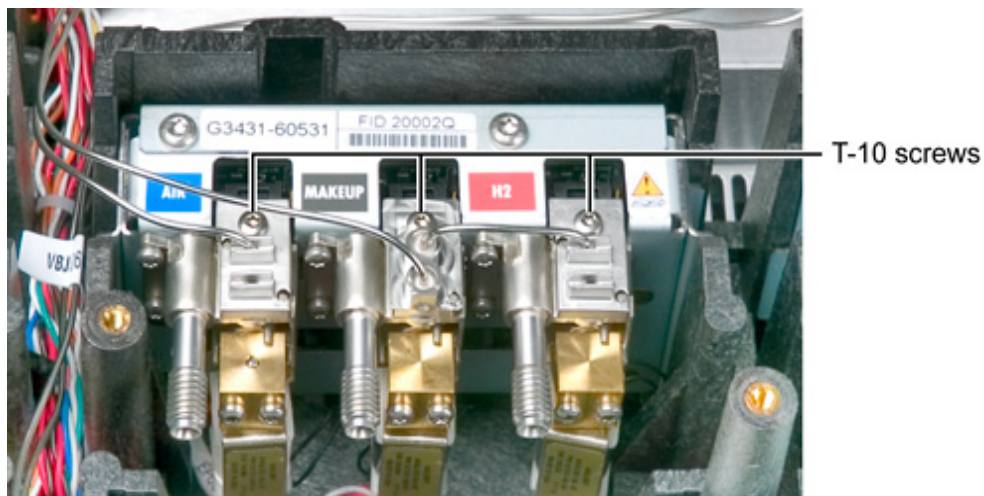

- 12** Remove the four T-20 screws retaining the FID assembly to the mounting pallet.

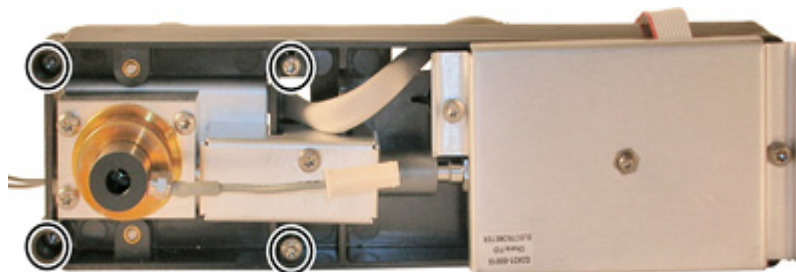

- 13** Lift the entire FID assembly up and off.
- 14** At this point the entire unit can be replaced with a new or repaired unit.
- 15** Reassembly is the reverse of these steps.
- 16** After completing the reassembly turn the power on and reconfigure the GC from the front panel keypad.

## Ultrasonic cleaning the FID

The collector requires occasional cleaning to remove deposits (usually white silica from column bleed, or black, carbonaceous soot). Deposits reduce sensitivity and cause chromatographic noise and spikes. The cleaning procedure presented here suggests you use an ultrasonic bath to clean the collector and other parts of the detector. However, if your collector is not too dirty, it may be sufficient to scrub it with a nylon brush and then use a burst of compressed air or nitrogen to blow stray particles away.

**WARNING**

**This procedure summarizes the general steps for cleaning the parts. You need to follow the standard safety practices of your laboratory for handling chemicals. For example, wear the appropriate safety eye glasses, lab coat, and gloves.**

---

**CAUTION**

Scratches on the FID jet affects its performance. If you choose to clean the jet, be sure not to scratch or bend the jet.

---

**CAUTION**

The FID castle is coated with a layer of PTFE. Ultrasonic cleaning of this part for more than 5 to 10 seconds will damage this coating.

---

It is often more convenient to replace dirty jets with new ones than to clean them, especially jets that have been badly contaminated.

**1** Gather the following:

- Small ultrasonic cleaning bath
- Aqueous detergent
- GC-grade methanol in a PTFE wash bottle. Wash bottles made of other materials usually contain plasticizer contaminants
- Dry, filtered, compressed air or nitrogen
- Clean cloth

**2** If you are cleaning the jet, run the cleaning wire through the tip of the jet. Run it back and forth a few times until it moves smoothly. Be careful not to scratch the jet.





- 3 Place the parts in your glassware, cover them with the aqueous detergent, and place them in the ultrasonic bath.
- 4 If you are cleaning the castle, sonicate for only 5 to 10 seconds. Push the wire brush through the collector.

For other parts, sonicate for 1 to 5 minutes.

- 5 Remove the parts from the bath.
- 6 Push the wire brush through the collector.




- 7 Rinse the parts in tap water or distilled water.
- 8 Rinse all surfaces of the parts in methanol; inside and outside surfaces. To insure good rinsing use either a PTFE wash bottle or a beaker.
- 9 From this point on, handle the parts only with forceps (or tweezers). Remove the jet from the bath and rinse it thoroughly with hot tap water and then with a small amount of methanol.
- 10 Blow the jet dry with a burst of compressed air or nitrogen and then place the jet on a clean cloth to air dry.

## Nitrogen Phosphorus Detector (NPD)

The NPD is a variation on the FID, in that the sample is burned and the resulting ions are collected. However, the hydrogen/air ratio is set to suppress carbon ionization and an alkali salt enhances nitrogen and phosphorus ionization.

### Selecting an NPD jet

Open the oven door and locate the column connection fitting at the base of the detector. It will look like either a capillary optimized fitting or an adaptable fitting.

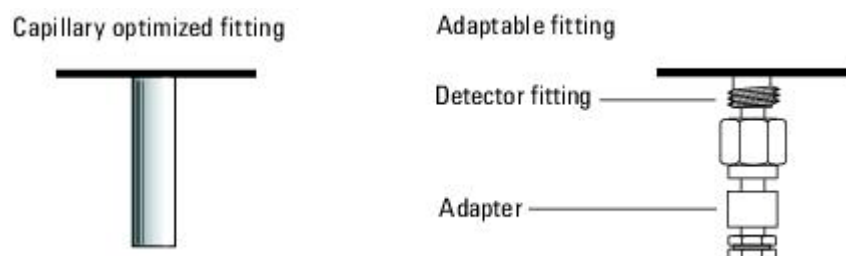


- If you have an application that tends to clog the jet, select a jet with a wider tip id.
- When using packed columns in high column-bleed applications, the jet tends to clog with silicon dioxide.

For capillary optimized fittings, select one of the following from the table below, [Jets for capillary optimized fittings](#).

**Table 6** Jets for capillary optimized fittings

Capillary optimized NPD jets ID	Jet type	Part number	Jet tip id	Length
1	Capillary with extended jet (recommended)	G1534-80580	0.29 mm (0.011 inch)	51.5 mm
2	Capillary	G1531-80560	0.29 mm (0.011 inch)	43 mm
3	High-temperature	G1531-80620	0.47 mm (0.018 inch)	43 mm



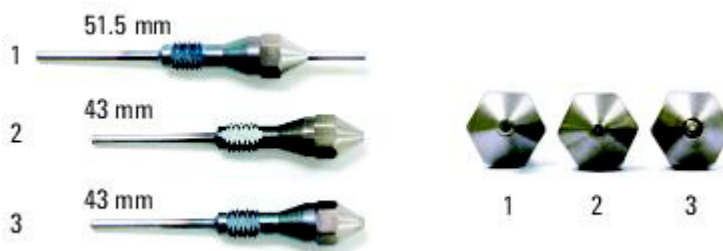

**Figure 4** Capillary optimized NPD jets

For adaptable fittings, select one of the following from the table below, [Jets for adaptable fittings](#).

**Table 7** Jets for adaptable fittings

Adaptable NPD jets ID	Jet type	Part number	Jet tip id	Length
1	Capillary with extended jet (recommended)	G1534-80590	0.29 mm (0.11 inch)	70.5 mm
2	Capillary	19244-80560	0.29 mm (0.011 inch)	61.5 mm
3	Capillary, high- temperature	19244-80620	0.47 mm (0.018 inch)	61.5 mm
4	Packed	18710-20119	0.46 mm (0.018 inch)	63.6 mm



**Figure 5** Adaptable NPD jets

## Removing the NPD electrometer

### Removing the NPD electrometer

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

---

**NOTE**

Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.


---

**NOTE**

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 3 Remove both the electronics top cover and the right side cover.
- 4 Lift up the hinged detector tower cover and remove the Torx T-20 screw and the clamp on the electrical interconnect.

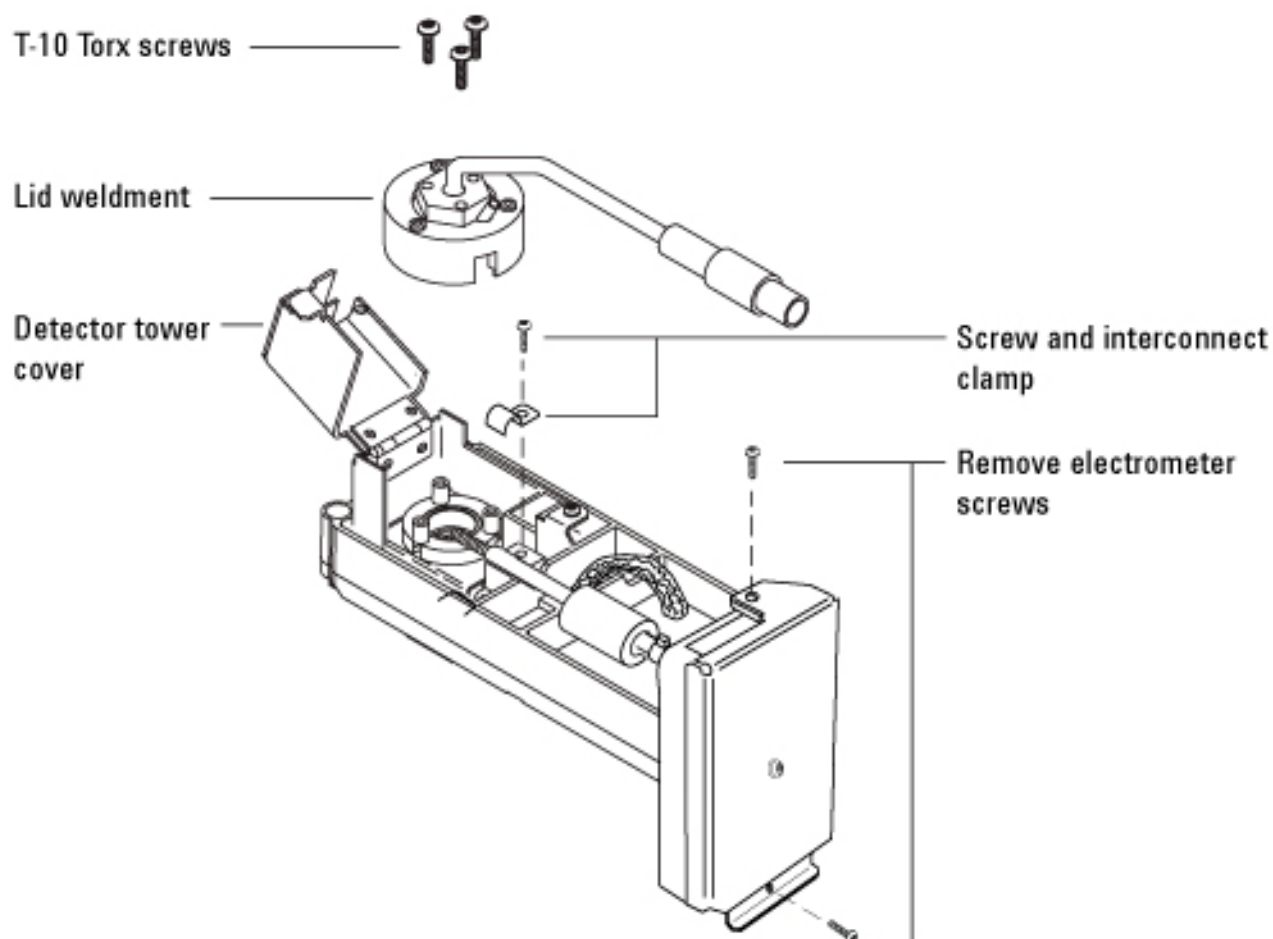

- 5 Loosen the three Torx T-10 screws on the lid weldment and remove the lid.
- 6 Remove one Torx T-20 screw from each end of the electrometer. (You do not need to remove the screw on the top of the electrometer that holds the cover on.)
- 7 Unlock and detach the electrometer's ribbon cable from the detector's interface board and lift the electrometer from the detector pallet.
- 8 Reassembly is the reverse of these steps.

## Replacing the entire NPD

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

### WARNING

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

---

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

### WARNING

**Inlets, detectors, and the oven are insulated with fibrous materials which may cause irritation to skin, eyes, and/or mucous membranes. Always wear gloves when working with the insulation. Additionally, if the insulation is flaky/crumbly, wear protective eyewear and a suitable breathing mask and/or respirator.**

---

- 3 Remove the right side electronics panel.
- 4 Disconnect the column from the bottom of the detector.
- 5 To access the NPD module unscrew the single T-20 screw holding the GC cover over the detector(s).
- 6 Raise it out of the way or remove it temporarily.
- 7 Disconnect the power/sensor cable to the NPD.
- 8 Disconnect the electrometer cable from the NPD signal board.
- 9 Remove the three T-10 screws fastening the NPD pneumatic tubing assembly to the EPC module.
- 10 Remove the four T-20 screws retaining the NPD pallet.
- 11 Lift the entire NPD unit up and off.
- 12 At this point the entire unit can be replaced.
- 13 Reassembly is the reverse of these steps.

- 14 Restore power and reconfigure the GC using the keypad, supplying the requested information.
- 15 Note that the replacement of some NPD components does not require the complete removal of this detector.

## Replacing an NPD signal board

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

**WARNING**

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

---

**CAUTION**

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 3 Turn off all gas flows at their sources.

**WARNING**

**Inlets, detectors, and the oven are insulated with fibrous materials which may cause irritation to skin, eyes, and/or mucous membranes. Always wear gloves when working with the insulation. Additionally, if the insulation is flaky/crumbly, wear protective eyewear and a suitable breathing mask and/or respirator.**


---

**CAUTION**

Wear lint-free gloves to minimize source contamination when servicing the bead assembly.

---

- 4 Disconnect three electrical connections (to electrometer, to heater, and the 4-wire communication buss).
- 5 Turn thumb screw counterclockwise, lift screw to disengage, slide board down, and then lift out. (Note that the grounding point for the board is at the thumbscrew.)
- 6 Insert the replacement board.
- 7 Reassembly is the reverse of these steps.
- 8 Reconfigure the modified GC from the front panel keypad, supplying requested information.

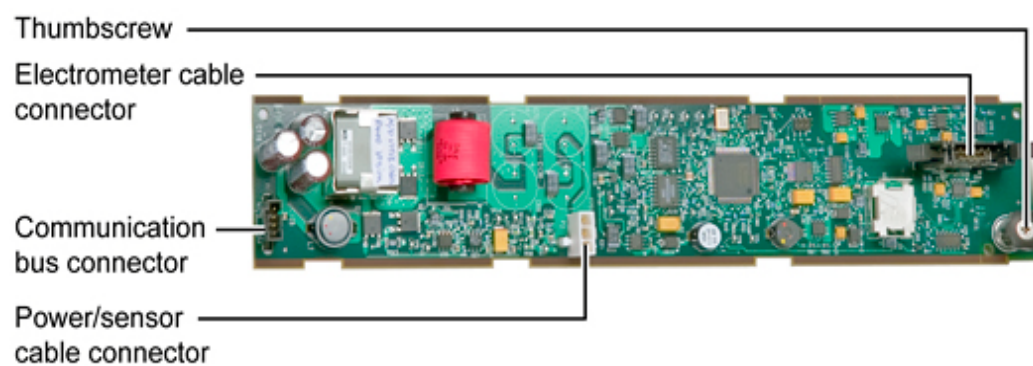

## Replacing the NPD heater

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

**WARNING**

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

---

**CAUTION**

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 3 Turn off the oven and all heated zones and let them cool.
- 4 Turn off all gas flows at their sources.
- 5 Turn off the GC main power switch and disconnect its power cord.


**WARNING**

**Inlets, detectors, and the oven are insulated with fibrous materials which may cause irritation to skin, eyes, and/or mucous membranes. Always wear gloves when working with the insulation. Additionally, if the insulation is flaky/crumbly, wear protective eyewear and a suitable breathing mask and/or respirator.**

---

- 6 Disconnect the column from the bottom of the detector.
- 7 Remove the right side mainframe panel.
- 8 Access the NPD module by unscrewing the single T-20 screw holding the GC cover over the detector(s).
- 9 Raise it out of the way or remove it temporarily.
- 10 Disconnect the power/sensor cable to the NPD detector.
- 11 Disconnect the electrometer cable from the NPD signal board.
- 12 Remove the three T-10 screws fastening the NPD pneumatic tubing connectors to the EPC module.



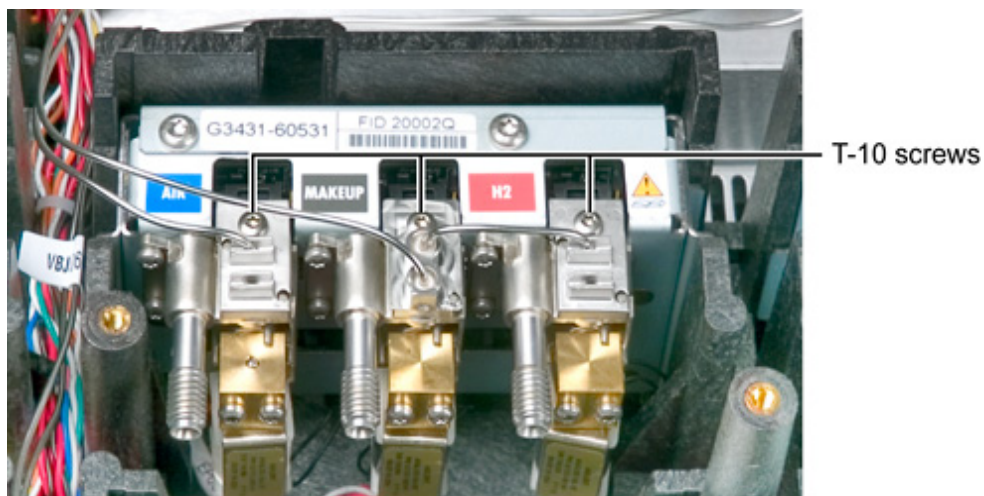

- 13 Remove the four T-20 screws retaining the NPD assembly to the GC.
- 14 Remove the NPD detector assembly from the GC.
- 15 Examine the underside of the detector assembly.
- 16 Cover or plug the sample inlet tube while servicing this unit to limit contamination.
- 17 Remove the insulated cover over the heater/detector block.
- 18 Loosen the heater cable sufficiently to allow for the lifting off and removal of the heater from the detector block.
- 19 Remove and replace the entire heater/cable assembly.
- 20 Reassembly is the reverse of these steps.

## Replacing the NPD electrometer

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

### WARNING

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

---

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 3 Turn off all gas flows at their sources.

### WARNING

**Inlets, detectors, and the oven are insulated with fibrous materials which may cause irritation to skin, eyes, and/or mucous membranes. Always wear gloves when working with the insulation. Additionally, if the insulation is flaky/crumbly, wear protective eyewear and a suitable breathing mask and/or respirator.**

---



### CAUTION

Wear lint-free gloves to minimize source contamination when servicing the bead assembly.

---

- 4 Remove both the electronics top cover and the right side cover.
- 5 Lift up the hinged NPD cover and remove the T-20 screw and J-clamp on the interconnect assembly.
- 6 Loosen the three Torx T-10 screws on the lid weldment and remove the lid.
- 7 Remove one Torx T-20 screw from each end of the electrometer. Do not remove the screw on the top of the electrometer that holds the cover on.
- 8 Disconnect the power/sensor cable from the bead assembly and from the NPD signal board.

- 9 Unlock and detach the electrometer's ribbon cable from the NPD signal board and lift the electrometer up from the pallet.
- 10 Reassembly is the reverse of these steps.
- 11 Reconfigure the modified GC from the front panel keypad, supplying requested information.

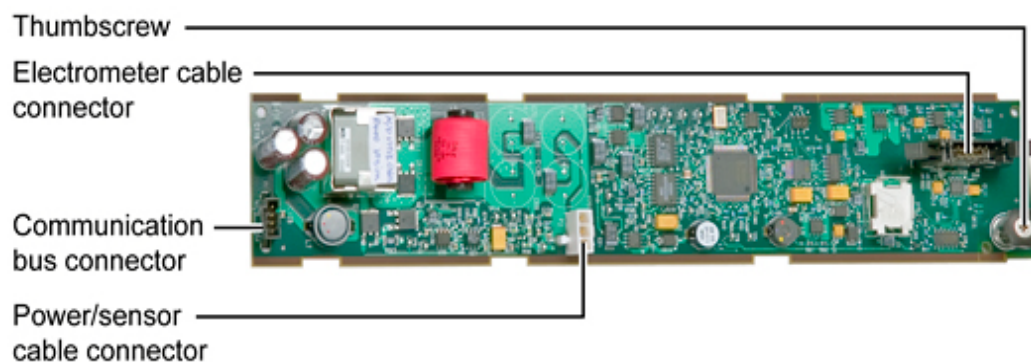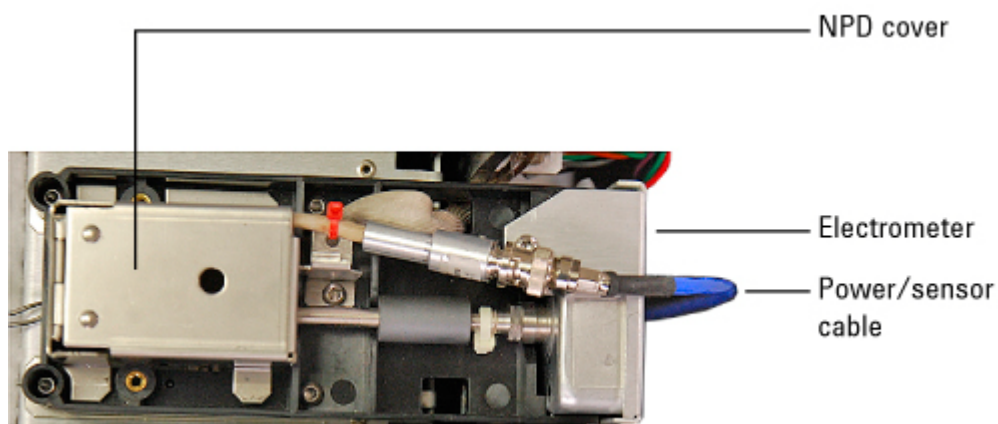

## Cleaning the NPD jet and collector

The collector requires occasional cleaning to remove deposits (usually white silica from column bleed, or black, carbonaceous soot). Deposits reduce sensitivity and cause chromatographic noise and spikes. The cleaning procedure presented here suggests you use an ultrasonic bath to clean the collector and other parts of the detector. However, if your collector is not too dirty, it may be sufficient to scrub it with a nylon brush and then use a burst of compressed air or nitrogen to blow stray particles away.

### WARNING

**This procedure summarizes the general steps for cleaning the parts. You need to follow the standard safety practices of your laboratory for handling chemicals. For example, wear the appropriate safety eye glasses, lab coat, and gloves.**

---

### CAUTION

Scratches on the NPD jet affects its performance. If you choose to clean the jet, be sure not to scratch or bend the jet.

---

It is often more convenient to replace dirty jets with new ones than to clean them, especially jets that have been badly contaminated.

**1** Gather the following:

- Small ultrasonic cleaning bath
- Aqueous detergent
- GC-grade methanol in a PTFE wash bottle. Wash bottles made of other materials usually contain plasticizer contaminates.
- Dry, filtered, compressed air or nitrogen
- Clean cloth

**2** If you are cleaning the jet, run the cleaning wire through the tip of the jet. Run it back and forth a few times until it moves smoothly. Be careful not to scratch the jet.



- 3 Place the parts in your glassware, cover them with the aqueous detergent, and place them in the ultrasonic bath for 1 to 5 minutes.
- 4 Remove the parts from the bath and push the wire brush through the collector.
- 5 Rinse the parts in tap water or distilled water.
- 6 Rinse all surfaces of the parts in methanol; inside and outside surfaces. To insure good rinsing use either a PTFE wash bottle or a beaker.
- 7 From this point on, handle the parts only with forceps (or tweezers). Remove the jet from the bath and rinse it thoroughly with hot tap water and then with a small amount of methanol.
- 8 Blow the jet dry with a burst of compressed air or nitrogen and then place the jet on a clean cloth to air dry.

## Cleaning the NPD

- 1 Turn off all gas flows at their sources.
- 2 Turn off the GC main power switch.

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

### WARNING

**Inlets, detectors, and the oven are insulated with fibrous materials which may cause irritation to skin, eyes, and/or mucous membranes. Always wear gloves when working with the insulation. Additionally, if the insulation is flaky/crumbly, wear protective eyewear and a suitable breathing mask and/or respirator.**

---

### CAUTION

Wear lint-free gloves to minimize source contamination when servicing the bead assembly.

---

- 3 Remove the right side mainframe panel.
- 4 To access the NPD module unscrew the single T-20 screw holding the GC cover over the detector(s).
- 5 Raise it out of the way or remove it temporarily.
- 6 Disconnect the power/sensor cable to the NPD.
- 7 Disconnect the electrometer cable from the NPD signal board.
- 8 Remove the three T-20 screws fastening the NPD pneumatic tubing assembly to the EPC module.
- 9 Remove the four T-20 screws retaining the NPD assembly to the GC.
- 10 Lift entire unit up and off.
- 11 Place unit on a clean surface.
- 12 Lift up the hinged metal cover over the detector assembly.
- 13 Remove the three T-10 screws retaining the bead assembly to the housing and then remove the bead assembly to prevent accidental damage.
- 14 Do not touch the ceramic bead.

- 15 Apply protective cap over the bead assembly, if available.
- 16 Loosen the three captive T-20 screws retaining the lid.
- 17 Lift the lid up and off.
- 18 Using a pair of tweezers, carefully remove the three O-rings, two alumina insulators, and the collector assembly.
- 19 At this point one can choose to clean only the parts removed or continue the disassembly to remove, clean, or replace the jet.
- 20 Use a 1/4-inch deep-socket nut driver to unscrew the jet.
- 21 Use slight side pressure on the nut driver to drag the jet up out of the detector weldment and then use tweezers to grasp and remove the jet.
- 22 Reassembly is the reverse of these steps.

## Thermal Conductivity Detector (TCD)

The TCD is sometimes called the "universal" detector because it responds to anything that is not the carrier gas.

### Replacing the TCD detector cell

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

#### **WARNING**

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

---

#### **CAUTION**

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 3 Turn off any detector gases at their supply.
- 4 Shut off the detector and the detector gases and let the detector cool.
- 5 Remove the detector cover, the electronics cover, the right side cover, and the rear top panel.
- 6 Disconnect the TCD filament leads which run from the detector to the detector interface card. Disconnect the wires from the detector interface card using a small flat blade screwdriver to push down on the connector tabs while you pull out the wires.
- 7 Remove the cover and the insulation.
- 8 Use a T-20 Torx screwdriver to remove the two screws securing the detector cell to the aluminum detector carrier bracket and lift the cell from the bracket.
- 9 Slide the heater/sensor assembly from the detector cell.
- 10 Slide the PRT out of the detector cell.



11 Install the new cell and reassemble the detector. Reassembly is the reverse of these steps with the following additional considerations:


- Before replacing the insulation, place the cap that came with your detector over the detector vent to prevent plugging the vent with insulation. Remove the cap once the insulation is in place.

12 When replacing the detector cover, make sure that the hole in the top of the cover is positioned over the detector vent and that the filaments, plumbing and heater/sensor leads including the  $\Delta$ PRT leads are positioned under the appropriate cut-outs on the sides of the cover.

- Make sure that only the filament wire is inserted into the connector, not the insulation sleeve. Then, check the filaments by tugging slightly on them.
- Check the  $\Delta$ PRT sensor lead connections to the detector card by tugging slightly on them.
- Replace the heater/sensor and  $\Delta$ PRT in the sensor holes as shown below.

**CAUTION**

Make sure the  $\Delta$ PRT leads are properly installed. If the  $\Delta$ PRT leads are not properly installed in the detector card, filament burnout can occur.

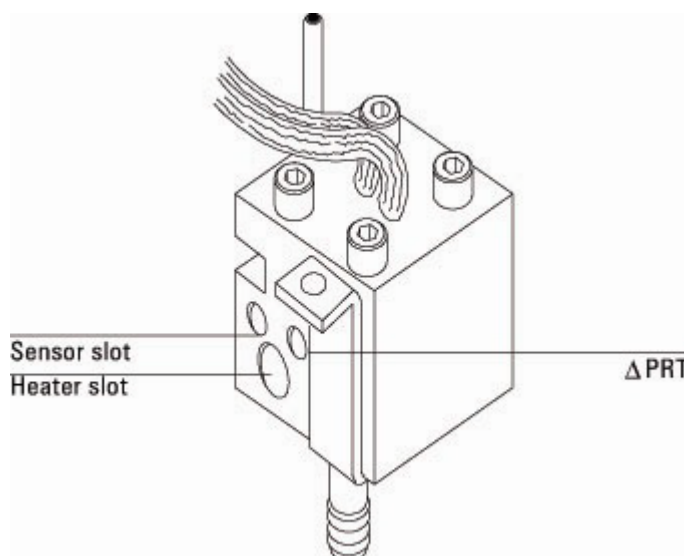


## Replacing the TCD switching valve

- 1 Turn off all gas flows at their sources.
- 2 Turn off the GC main power switch.

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

- 3 Remove the EPC module cover and the module retainer. With a side-mounted TCD, remove the left side cover.
- 4 It is possible, but quite difficult, to replace the switching valve without further disassembly. We recommend raising the module for better accessibility. With a side-mounted TCD, remove the module from its bracket.

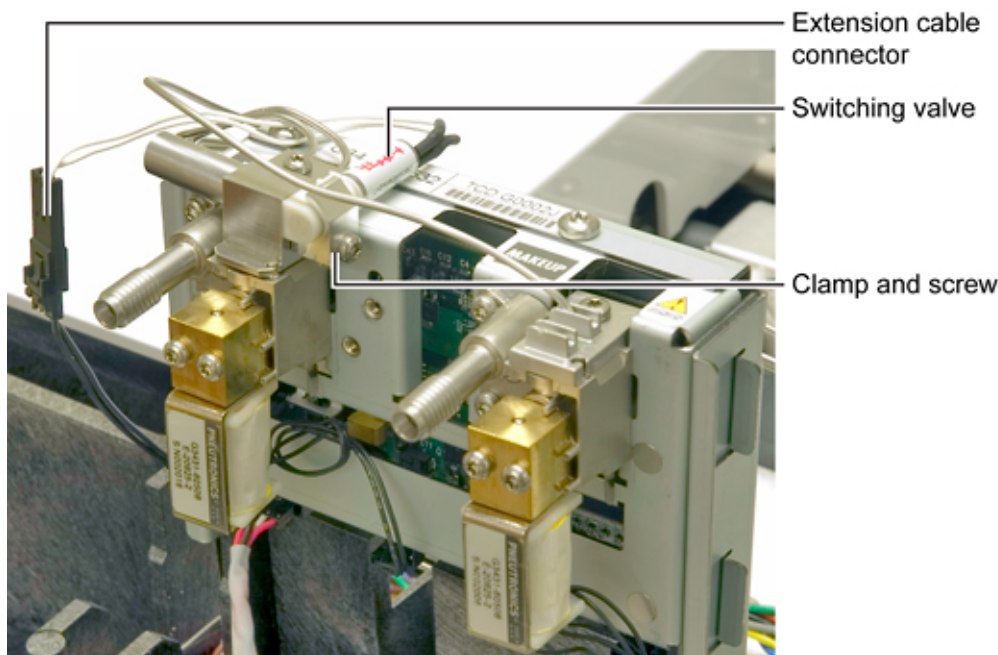

- 5 Trace the wires from the valve to the extension cable connector. Disconnect them.
- 6 Loosen the screw holding the clamp and slide the old valve out.
- 7 Examine the new valve. The plastic part has a flat side with 3 small O-rings. These must be placed over the 3 holes in the piece to which the screw connects.
- 8 Slide the new valve into position (flat side against the metal) and place the clamp over it.

- 9 Align the flat end of the valve with the flat metal surface next to it. Adjust the valve position until it lies flat against the adjacent metal.
- 10 Tighten the clamp screw. There should be no gap between the plastic and the metal. If there is, loosen the clamp screw and repeat the adjustment.
- 11 Reconnect the wires and remount the module.
- 12 Restore the covers.

## Micro-cell Electron Capture Detector (uECD)

This detector simply ignores most compounds, but responds with enormous sensitivity to electron-accepting species such as the halogens.


### Replacing the uECD heater/sensor assembly


After removing the uECD detector from the GC, you can further disassemble it to replace the heater/sensor assembly.

#### **WARNING**

**The ECD micro-cell contains radioactive  $^{63}\text{Ni}$ . To reduce the risk of exposure, wear disposable gloves while handling the ECD micro-cell. When you are finished, dispose of the gloves and wash your hands with soap and water.**

- 1 Remove the detector.
- 2 Loosen the locking tab screw on top of the detector, slide the locking tab back, and pivot it out of the way.
- 3 Lift the thermal cover up and carefully slide it off the anode assembly.
- 4 Remove the two screws holding the upper heated block onto the assembly. Lift the block over the anode lead and remove.
- 5 Slide the heater and sensor out of the lower heated block.

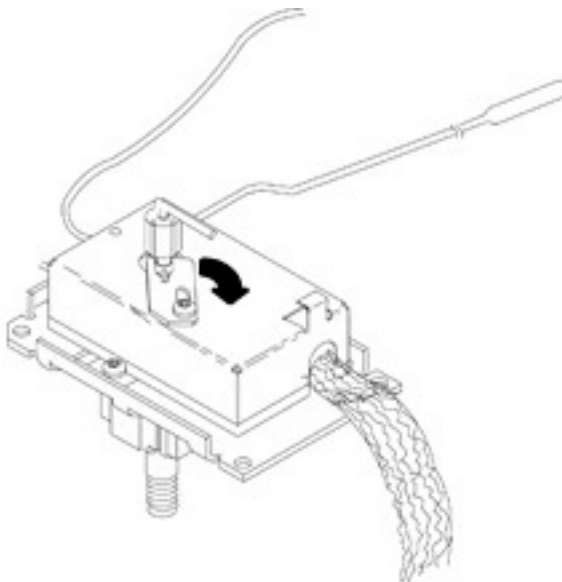


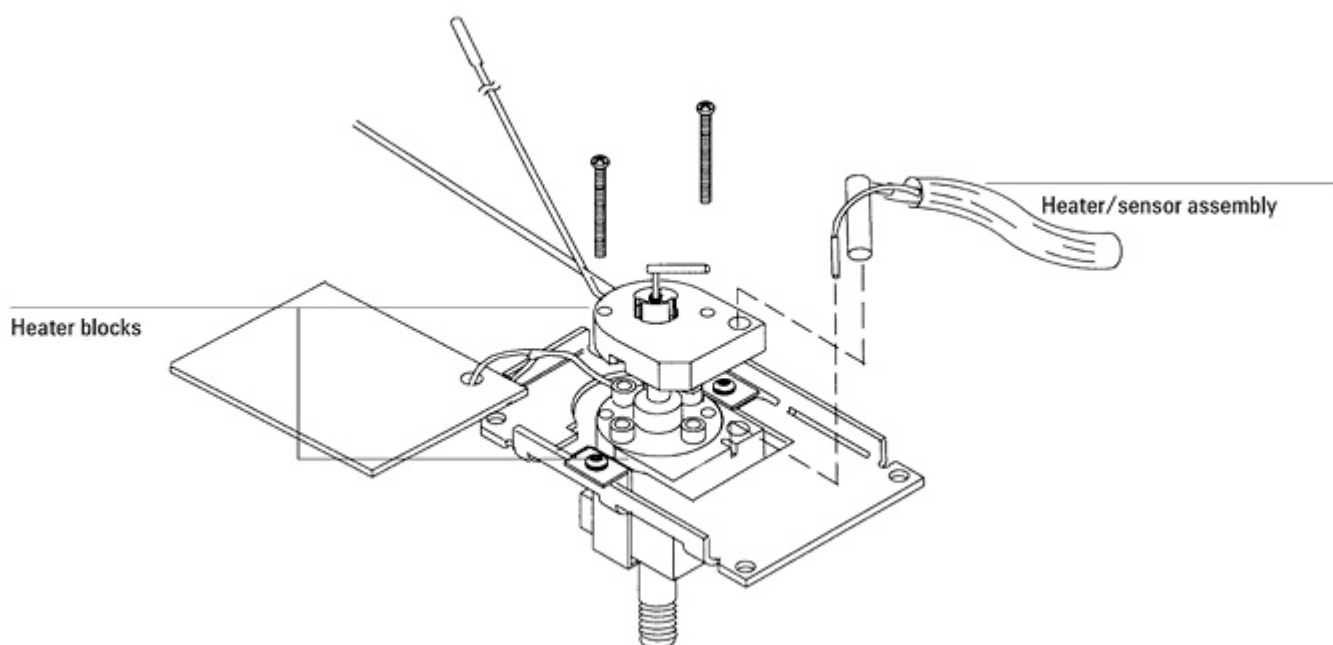

## Replacing the uECD makeup gas adapter

After removing the uECD detector from the GC, you can further disassemble it to replace the makeup gas adapter. The makeup gas adapter consists of a line from the detector pneumatics manifold that carries makeup gas to a weldment that screws into the bottom of the uECD detector. From there, the makeup gas sweeps past the end of the column and carries the column effluent into the uECD cell.

- 1 Remove the detector.
- 2 Remove the Torx T-20 screw holding the pneumatics block(s) to the detector manifold.
- 3 Slide the makeup gas adapter up and out of the GC.
- 4 When re-installing the makeup gas adapter, ensure the following:
  - Approximately 6 inches of the makeup gas line resides in the oven after installation.
  - The makeup gas line is bent into a coil or loop (inside the oven) that loops around the bottom of the detector weldment and makeup gas adapter.
  - The end of the column does not protrude from the top of the makeup gas adapter. For most columns (outer diameter > 0.20 mm), insert the column as far as it will go into the gigabore liner. If using a microbore column (id < 0.20 mm) that passes completely through the liner, position the column so that the total length from the back of the capillary nut to the end of the column is about 69 to 71 mm.

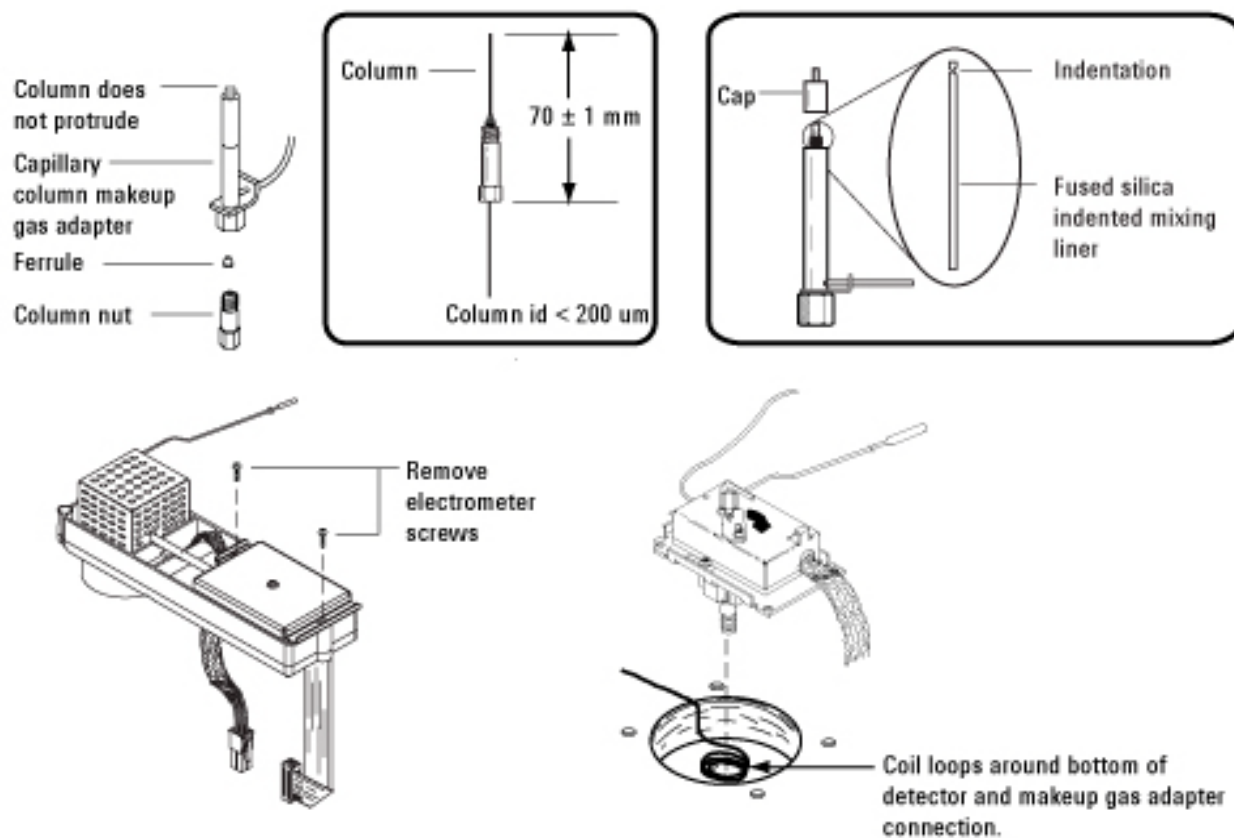

## Frequency test

Perform this test to make sure that the base frequency for the uECD during a blank run indicates a relatively contaminant-free system. See also the related troubleshooting manual topic on uECD Life Cycle.

### NOTE

It may take 24 hours for the uECD baseline to completely stabilize, especially if you are starting with a cold system and want to assure high-sensitivity operation. Therefore, for the most accurate results, run the detector at normal operating conditions for as long as possible (at least 2 hours and up to 24 hours) before running the frequency test. If you will be injecting into an unused inlet, you must use low-bleed septa. Make sure to condition new septa before use in an inlet for several hours with 1 to 5 mL/min carrier flow.

- 1 Make sure you are using normal operating conditions and that *at least* 2 hours have elapsed since the last run.
- 2 Turn on the uECD and the corresponding signal.
- 3 Check the displayed "Output":
  - $< 25$  = uECD frequency is acceptable
  - $\geq 1000$  = Contaminants in system

### NOTE

Each display count equals a frequency of 1 Hertz (e.g., a display reading of 100 = 100 Hz.).

- 4 If the uECD frequency indicates contamination ( $\geq 1000$ ) check for the following:
  - Contaminated carrier gas trap(s) and or supply—replace carrier gas supply tank and any traps on the carrier supply line.
  - Insufficient column conditioning—fully condition the column.
  - Contaminated detector—bake out the detector.
  - Column, inlet and/or septum bleed—clean the inlet/replace the septum with a conditioned, low bleed septum.



- Leaks—perform leak tests on both the inlet and detector systems.
- Anode current leakage—make sure the anode contacts are clean. Make sure the anode nut is tight.

## Leak test

**NOTE**

Once you have determined that the flow system components upstream from the detector (gas supply tubing, inlet, column fittings) are leak free, perform the following uECD detector leak test.

---

- 1 With the GC on and operating normally, set the oven, detector, and inlet temperatures to ambient.
- 2 Turn off the uECD and then turn off the inlet pressure.
- 3 Turn off the anode and makeup gas flows.
- 4 Cap the uECD exhaust vent with a vent plug (part no. 5060-9055).
- 5 Set carrier gas pressure at the inlet corresponding to the uECD to 15 psi (103 kPa).
- 6 Wait until the system reaches the setpoint pressure and then turn off the pressure and monitor the actual pressure value for at least 10 minutes.
- 7 Check for pressure drop:
  - If the pressure stays stable or drops only 0.5 psi, you can consider the uECD leakfree.
  - If the pressure drops more than 0.5 psi, you have a leak.

If you are sure none of the upstream flow system components are leaking, check for leaks at the column fitting and plugged inlet. If you find leaks, tighten the fittings and repeat the leak test.

**NOTE**

If you can find no other leaks, the uECD itself is probably leaking. The uECD cannot be disassembled without special license from the Nuclear Regulatory Commission or Agreement State Licensing Agency (USA only). Return the leaking uECD to Agilent for disposal.

---

## **Performing a radioactivity leak test (wipe test)**

Micro-cell ECDs must be tested for radioactive leakage at least every 6 months. Records of tests and results must be maintained for possible inspection by the Nuclear Regulatory Commission and/or responsible state agency. More frequent tests may be conducted when necessary.

The procedure used is the wipe test. A Wipe Test Kit (part no. 18713-60050) is supplied with each new uECD. Refer to the information card supplied in the Wipe Test Kit for instructions on performing the wipe test.

## Flame Photometric Detector (FPD) 7890A

The sample is burned in a hydrogen-rich flame. Excited species rise into a cooler zone above the flame, decay, and give off characteristic radiation. This is filtered and measured by a high-gain photomultiplier.

The intense yellow carbon radiation is blocked by a shield around the flame. Filters select either sulfur or phosphorus radiation.

### Preparing for maintenance

- 1 Turn off all gas flows at their sources.
- 2 Turn off the GC main power switch.

#### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

#### CAUTION

To prevent damage to the column or columns, remove the columns from the GC.

---


Additional tasks in preparation include:

- If you are replacing the transfer line, turn off the source gases to the FPD detector.
- If you are not replacing the transfer line, plug base of transfer line with a 1/8-inch Swagelok nut to keep it clean.


## Disassembling the FPD

Our objective is to disassemble the detector and prevent it from getting dirty. We recommend using lint free gloves during most of these steps.


- 1 Disconnect the spring and remove the photomultiplier tube (PMT). Place it in a safe place away from the light. (For a dual FPD, repeat for the second PMT.)

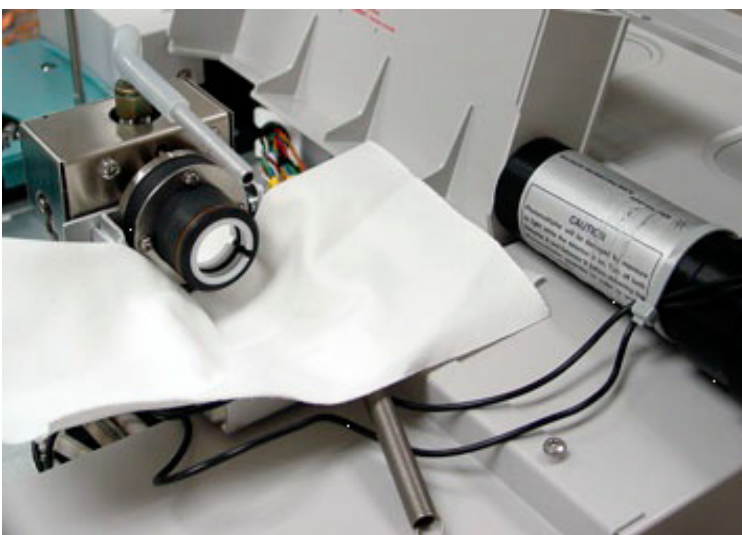

- 2 Remove and set aside the filter in a lint free cloth. For sulfur, remove the plastic spacer and then the filter. Try tapping the side of the detector or using the edge of a cotton swab. The phosphorus filter sits closer to the end of the housing.

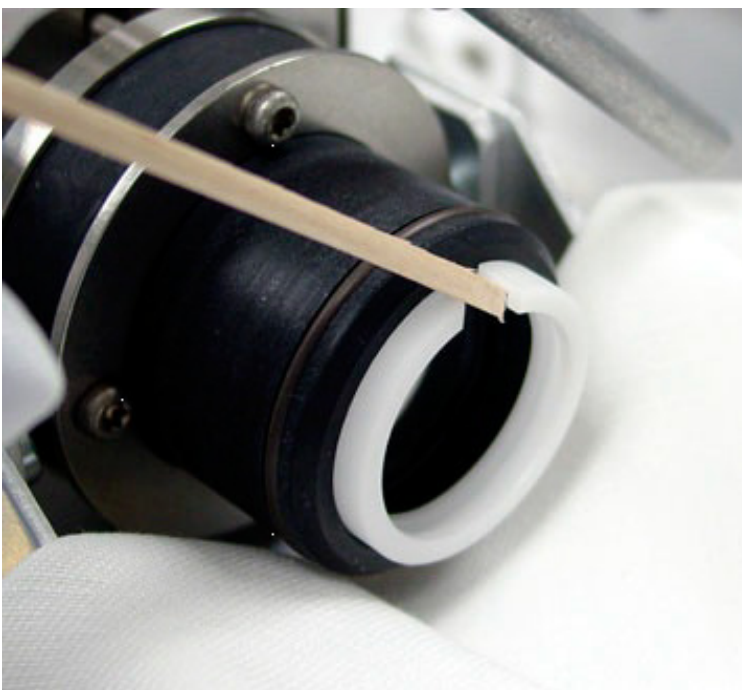


- 3 Remove the vent or exhaust tube with a 9/16-inch wrench.

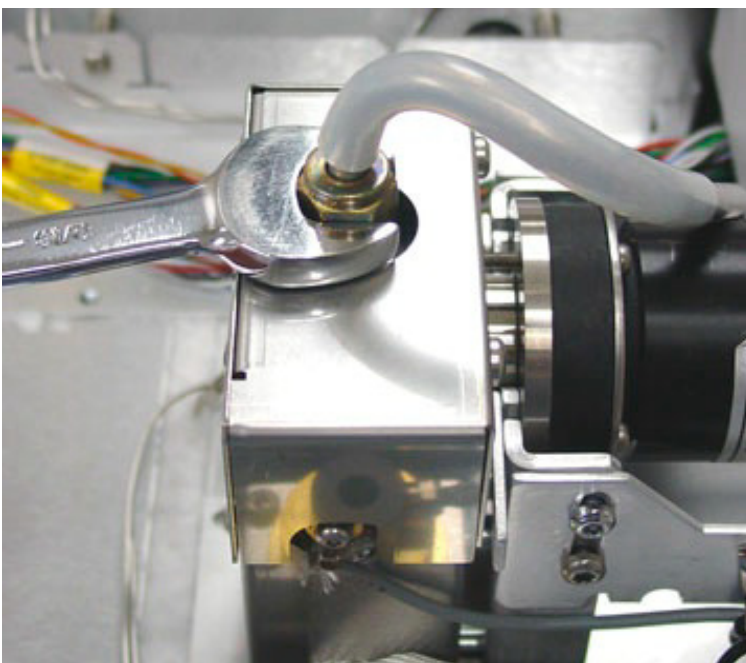

- 4 Remove the 4 screws and detector cover with a Torx T-20 driver.






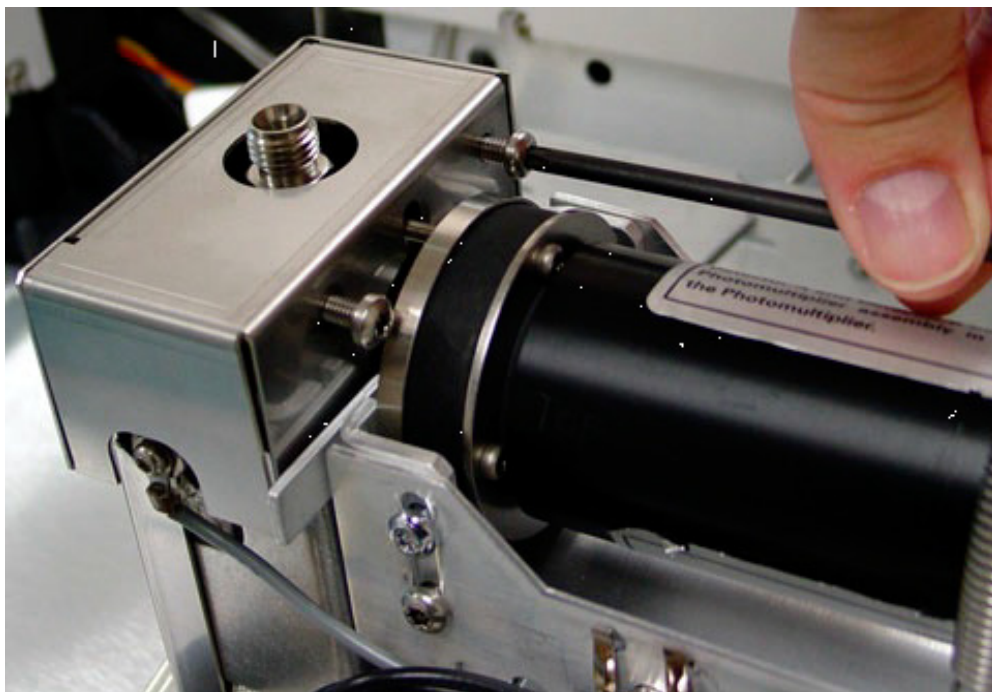

- 5 Loosen one screw and remove the ignitor collar with a Torx T-10 driver. The collar of your ignitor may have 1 or 2 screws.
- 6 Remove the heater/sensor assembly from the emission block.





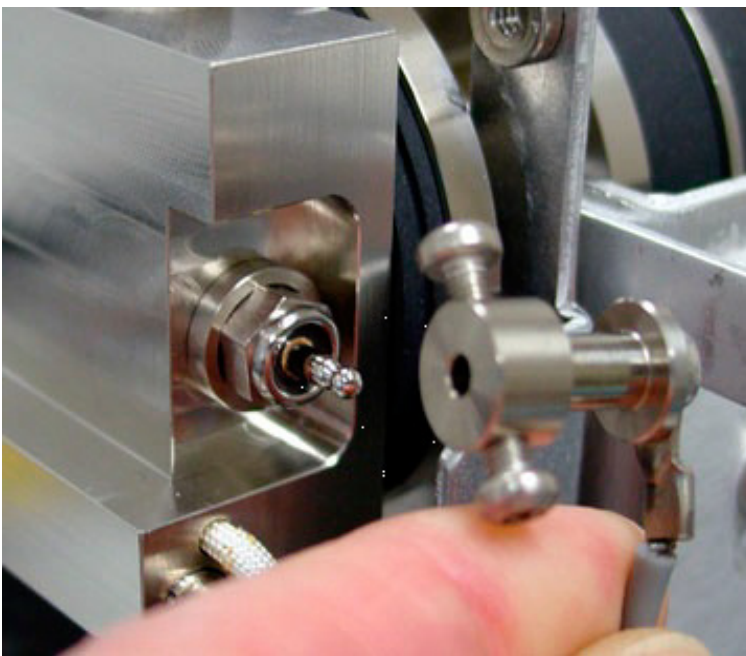

- 7** The heater/sensor assembly for the dual FPD attaches with a screw.



- 8** Remove one screw and retainer for the heater/sensor assembly in the transfer line with a Torx T-10 Key.




- 9** Remove the heater/sensor assembly.


- 10 Move the heater/sensor assemblies and ignitor cable to the front of the oven top.



- 11 Remove two 7-mm mounting nuts from the base of the transfer line.

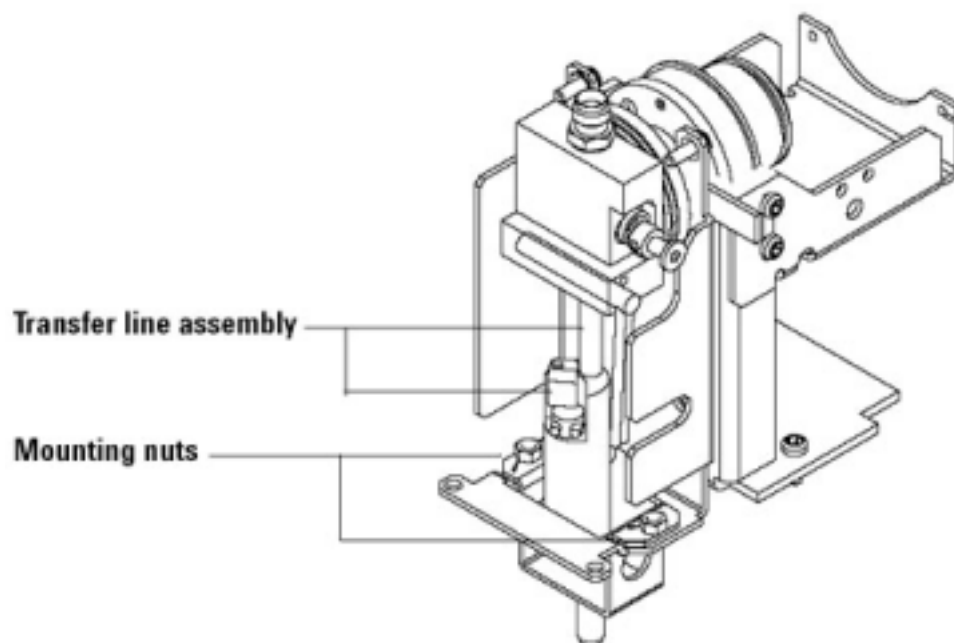

- 12** Loosen 3 screws that hold the optics assembly to the detector bracket with a Torx T-10 driver.

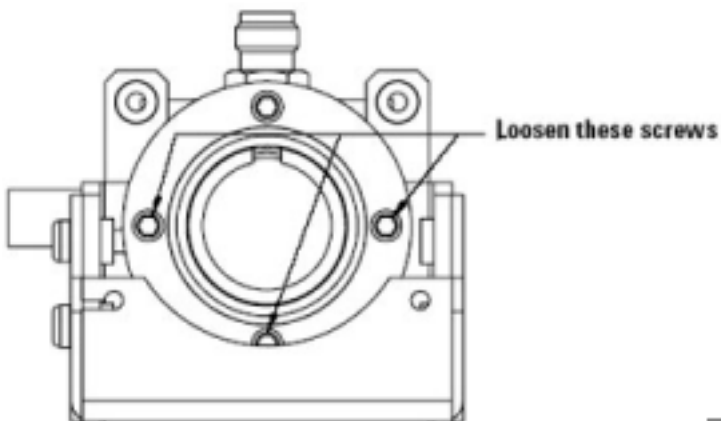


- 13** The objective of this step is to separate the optics from the transfer line. Do not flex the tubing where it is brazed to the transfer line weldment. Grasp the transfer line with your left hand and the optics assembly with your right hand. Lift them just high enough so that the 3 screws that you loosened in step 11 are above the bracket.





- 14** Twist and lift the optics while holding the transfer line stationary. Carefully separate the assemblies. The O-ring on the transfer jet is compressed against the inside of the emissions chamber. This is the resistance that you feel.
- 15** Place the optic assembly on a lint free cloth.

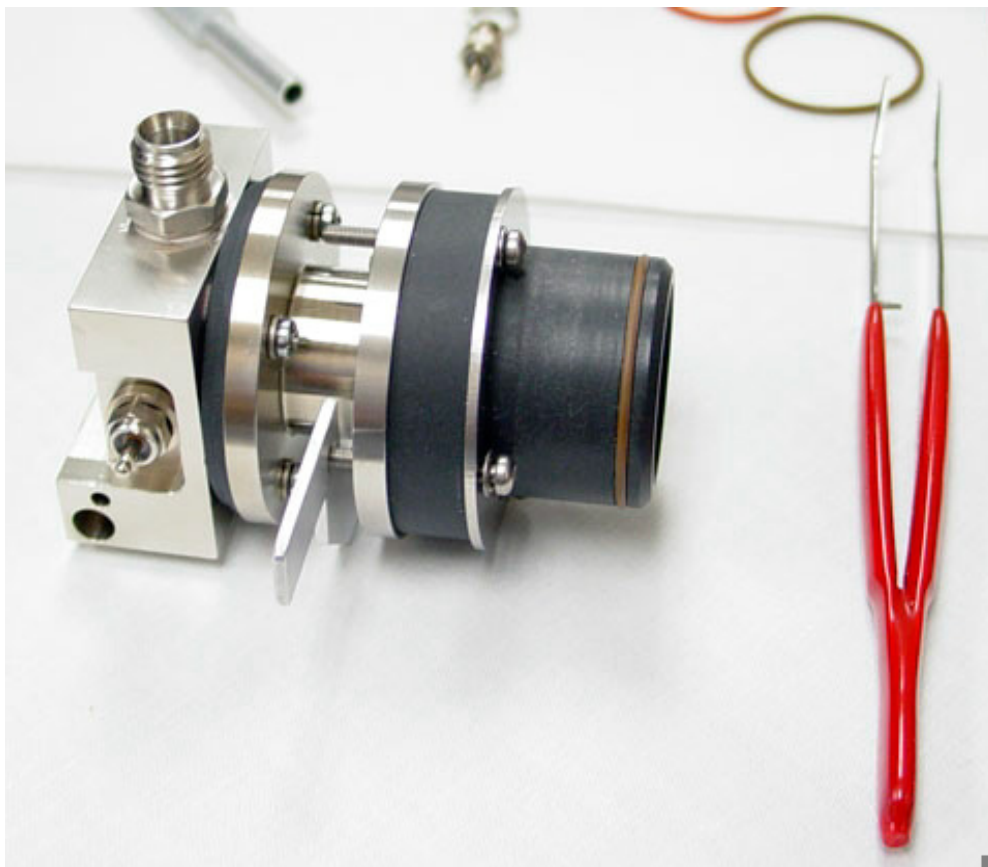

- 16** If you are not replacing the transfer line, skip to the section [“Rebuilding the FPD optics assembly”](#). If you are replacing the transfer line, continue with the next section.




## Rebuilding the FPD optics assembly

Use lint-free gloves when handling the optics assembly and O-ring seal.

- 1 Completely loosen the 4 T-10 screws from the clamp and coupling.




- 2 Set the Y-shaped clamp to one side.
- 3 Carefully separate the lens housing from the stainless steel coupling. The focusing lens and O-ring usually stick to the coupling, but could stay with the lens housing. Try to keep the housing, flange ring, screws and washers in place.

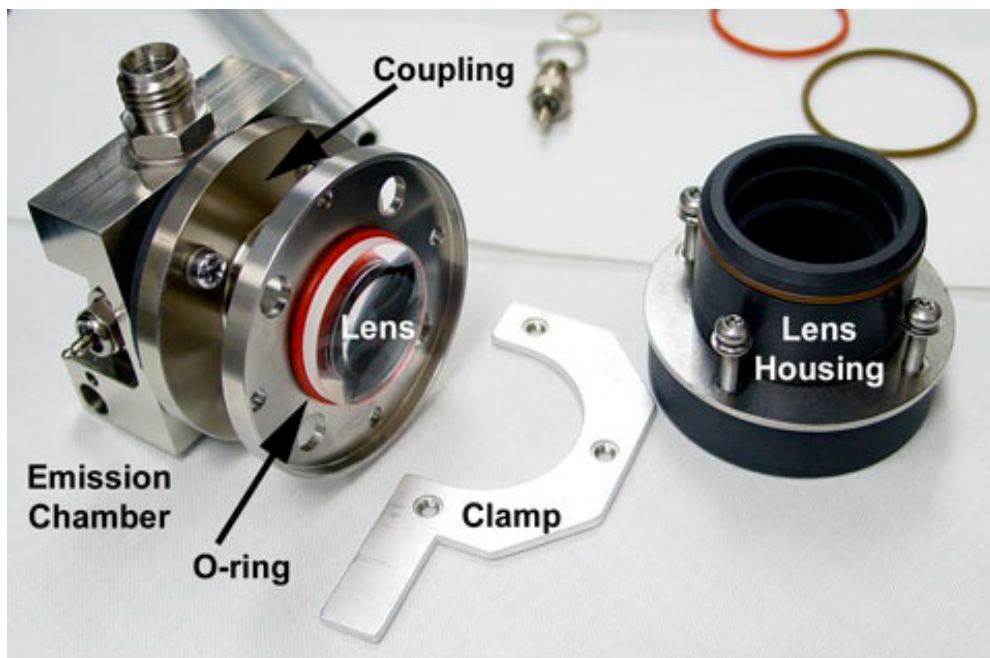

- 4 Remove the brown O-ring on the lens housing with a small pick. Roll the new seal over the edge of the housing until it sits in the groove.
- 5 Remove the coupling and heat shield disk from the old emission chamber with a No 1 Pozidrive or Philips screwdriver. Try to keep the screws, washers, and parts together.






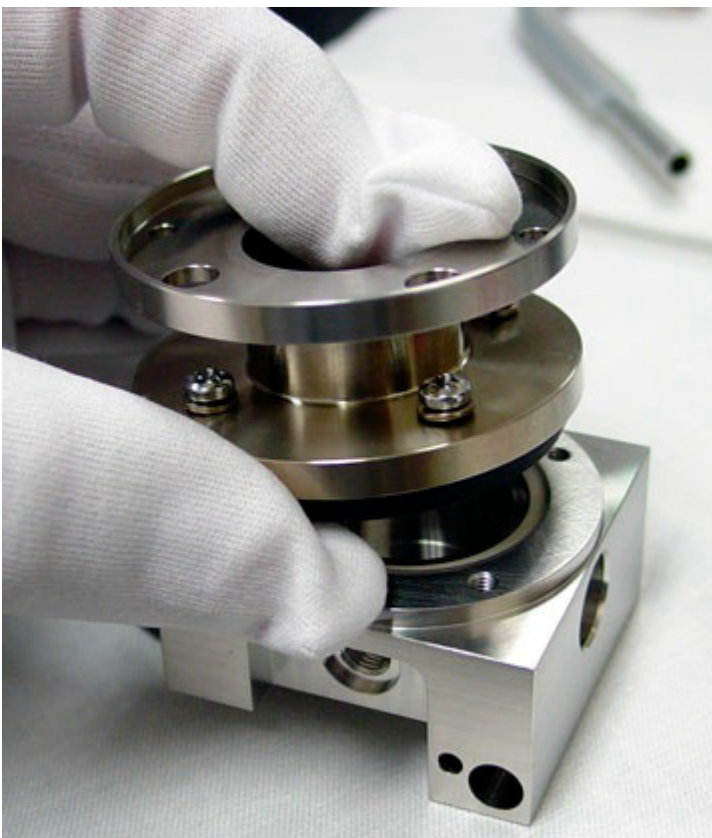



- 6 Assemble the ignitor parts. Slip the stainless steel spacer over the glow plug, followed by the O-ring. Do not use the copper washer.
- 7 Screw the glow plug into the new emission chamber and tighten with a 5/16-inch wrench.
- 8 Assemble the new emission chamber, seal, and window.



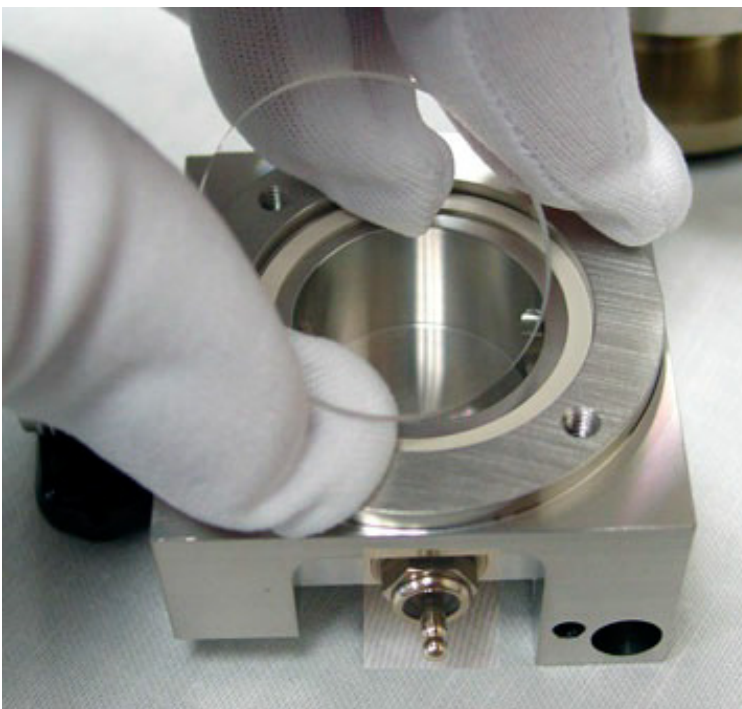



- 9** Position the heat shield and coupling onto the emission chamber so the threaded hole is aligned with the threaded fitting for the exhaust tube.
- 10** Reconnect the parts using the 4 screws and washers. Tighten opposite screws, similar to a wheel on an automobile, to insure a good seal.

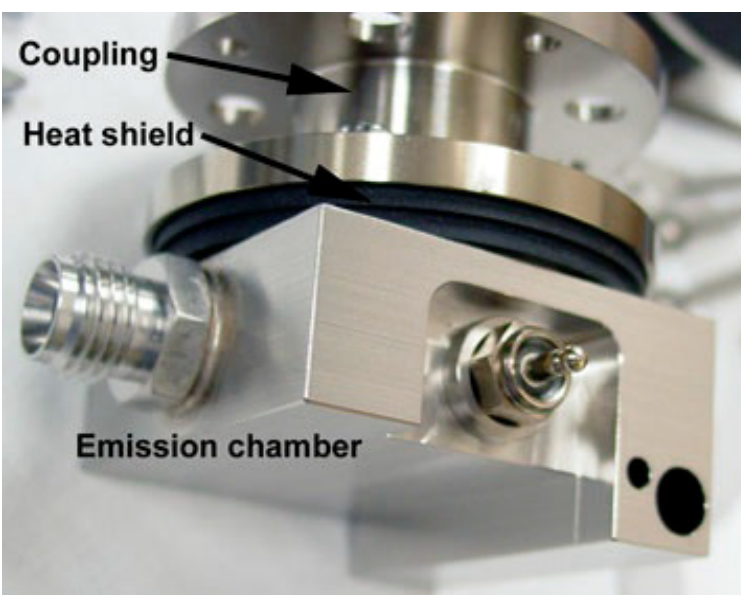

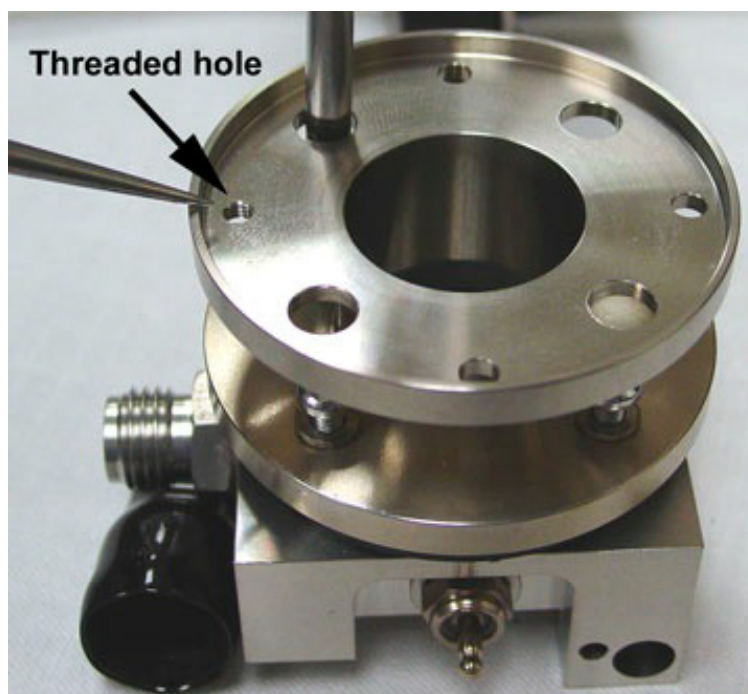

- 11** Turn the lens housing so that it faces the coupling. Insert the lens and rust colored O-ring.

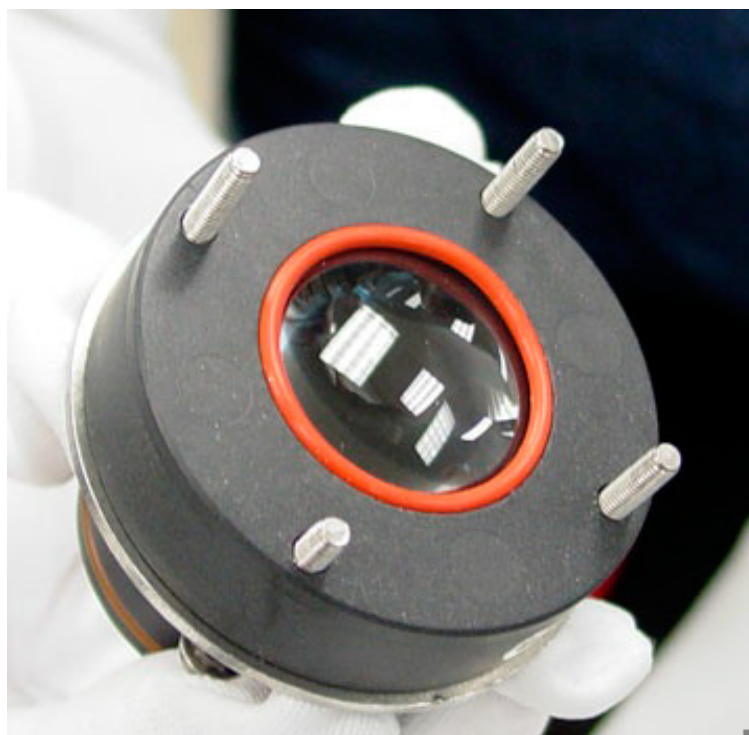


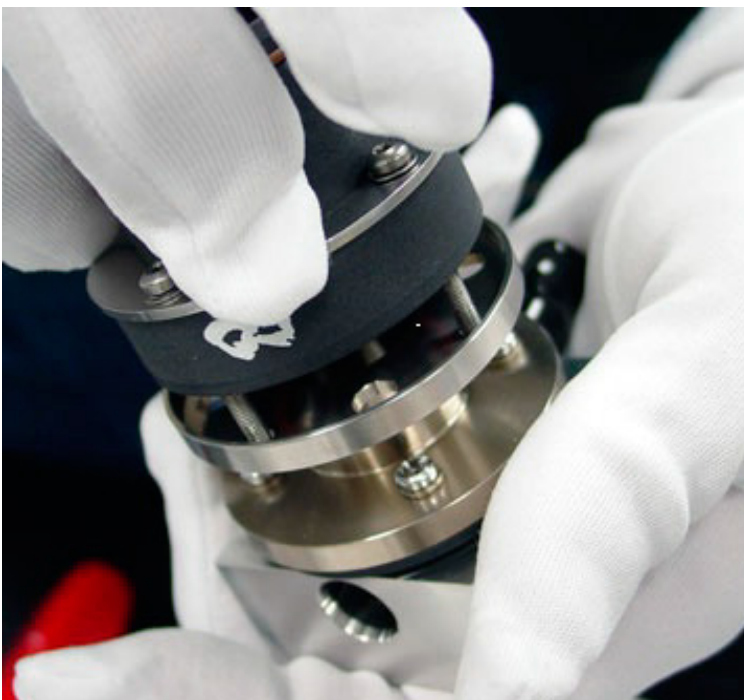

- 12** Place the housing onto the coupling and fasten the top screw just tight enough to keep the parts together.



- 13** Hold the y-shaped clamp and start the remaining screws into the clamp. Do not tighten until [step 7](#) of “Reassembling the FPD”.

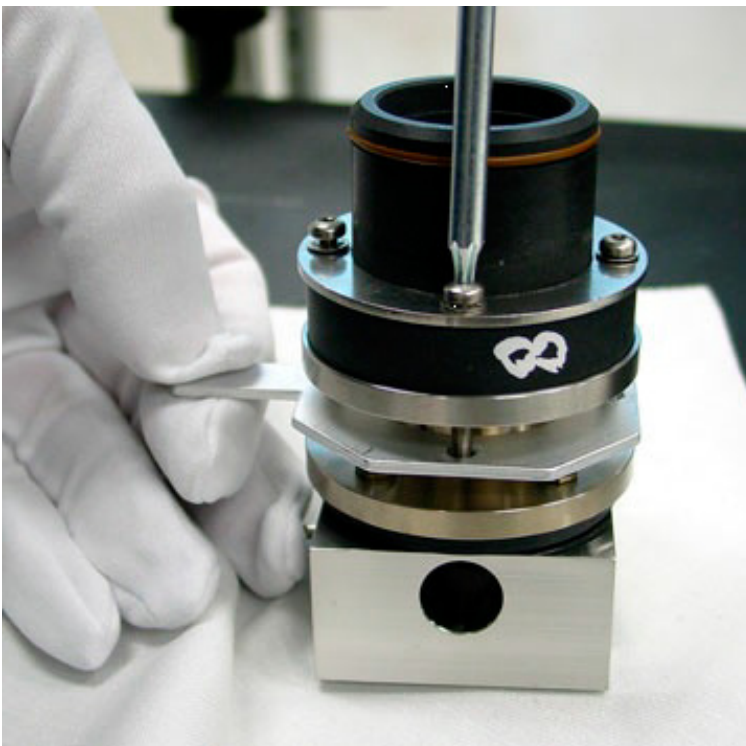



## Reassembling the FPD

This part of the procedure reassembles the optics assembly with the transfer line, reconnects both to the detector bracket, replaces the heater/sensor and ignitor wires, replaces the covers and PMT.

- 1 Replace or install the O-ring on the transfer line jet. Remove the O-ring with a small pick. Place the new O-ring over the jet. Roll it into the groove below the jet.




- 2 Reposition the detector optics above the bracket.
- 3 The objective of this step is to reattach the optics to the transfer line and insure a good seal between the O-ring and the emissions chamber. Grasp the transfer line with your left


hand and the optics assembly with your right hand. Push them together while twisting back and forth.

- 4 Before lowering the assemblies into the bracket, make sure there is no gap between the transfer line and the emission chamber.



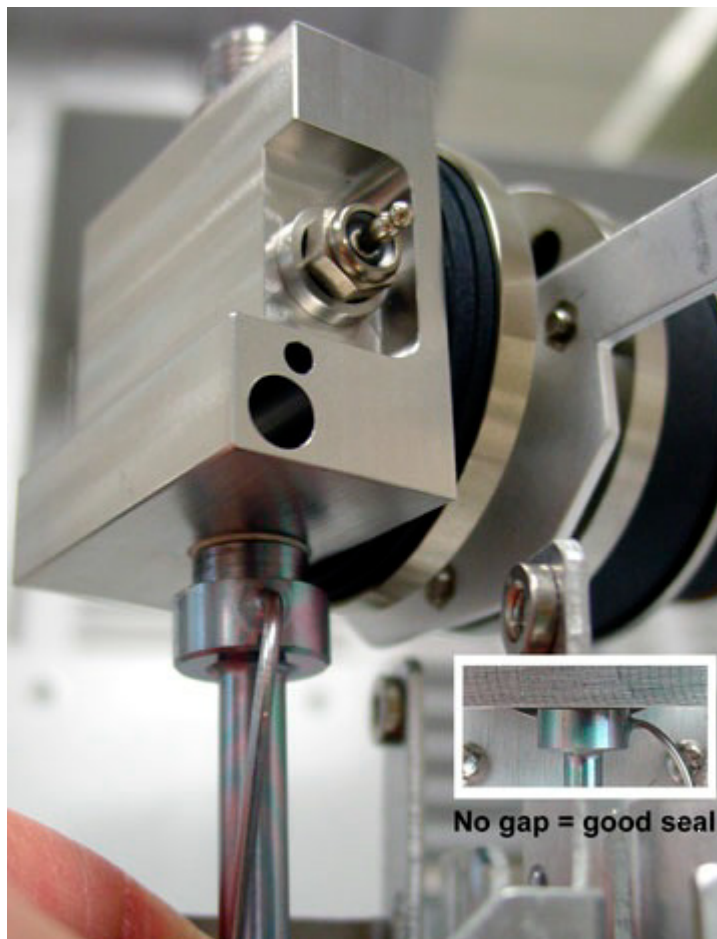

- 5 Lower the assemblies into the bracket. Line up the holes in the transfer line base with the threaded posts in the bracket. Line up the 3 screws and the clamp with the notches in the detector bracket.

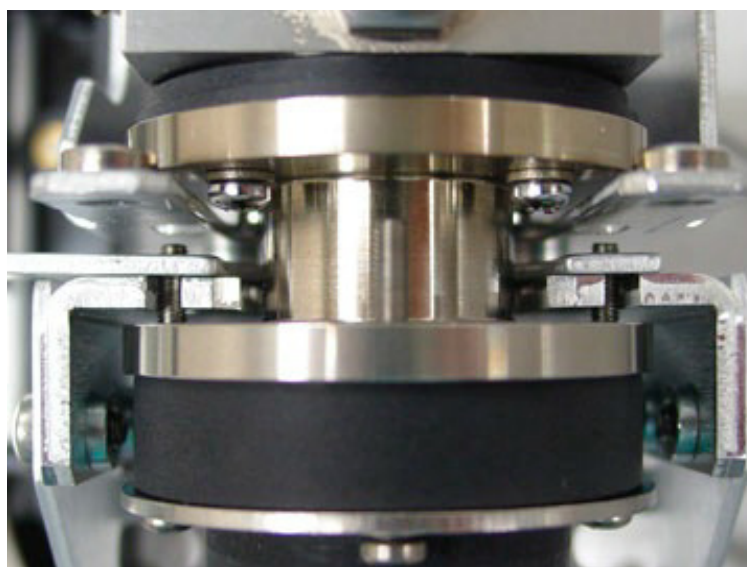




- 6** Reattach the transfer line to the bracket with the 7-mm nuts.
- 7** Tighten the clamp against the detector bracket. Tighten the 3 bottom screws on the optics assembly.
- 8** Insert the heater/sensor assembly into the transfer line. Make sure the sensor is at the bottom of the hole.
- 9** Reinstall the retainer and screw.









- 10 Insert the heater/sensor assembly into the emission chamber.
- 11 Reconnect the ignitor wire to the glow plug.

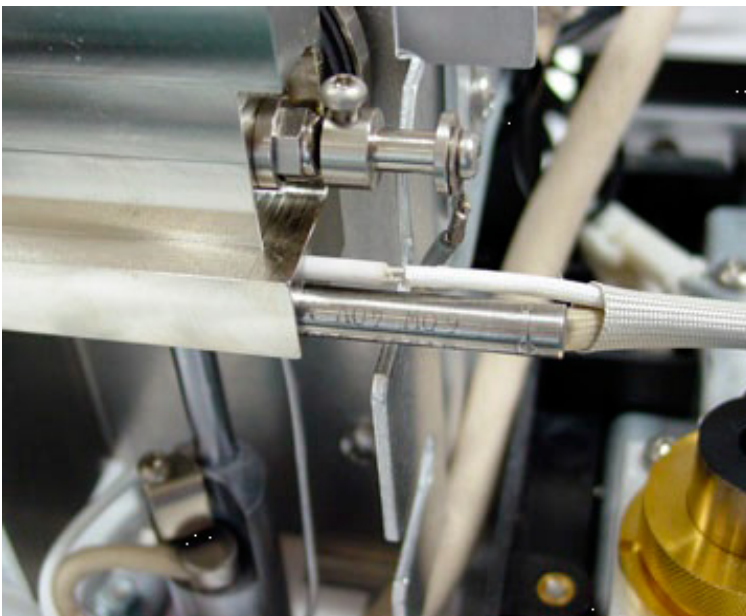


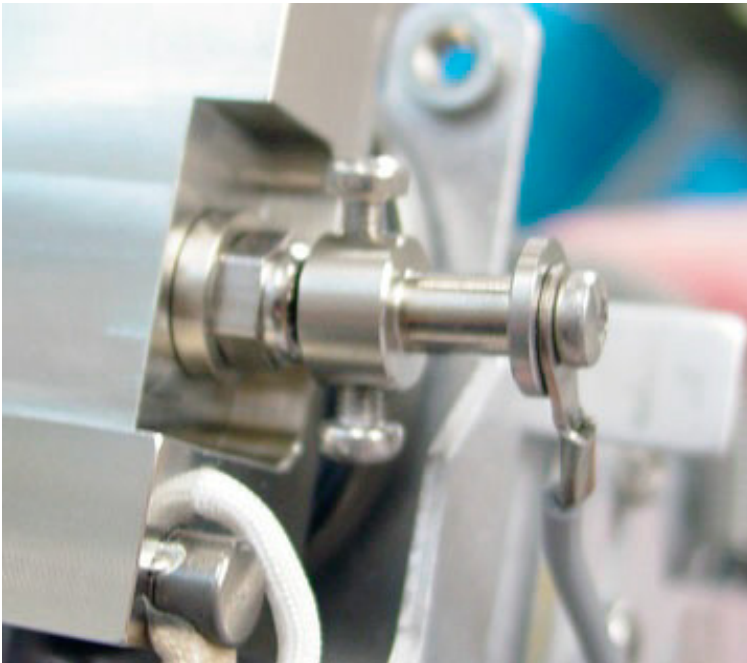


- 12** Replace the cover. First, start the 2 screws on the right-hand side of the cover. Second, start and tighten the screws at the base on the left-hand side. Third, tighten the screws on the right-hand side.

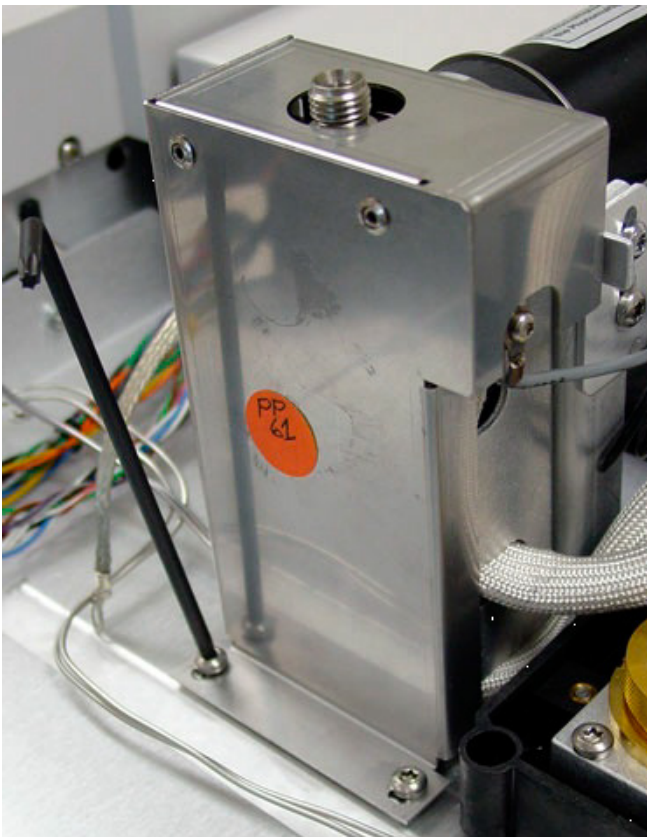



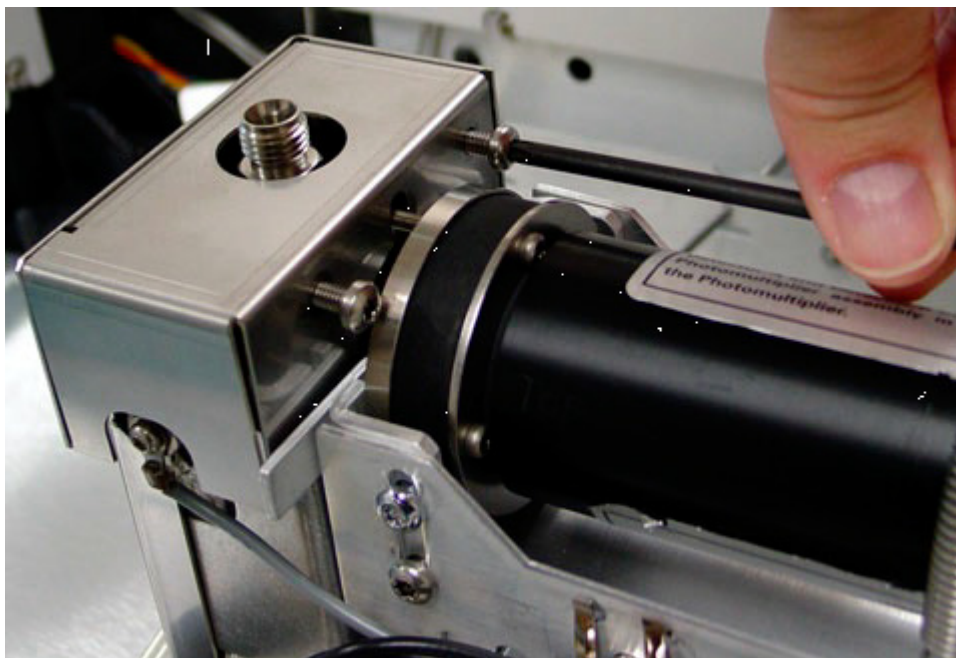


- 13** Replace the filter. If you are doing the Agilent checkout, run the phosphorus filter before the sulfur filter with spacer.
- 14** Replace the photomultiplier tube (PMT). Connect the spring.
- 15** Install the new vent tube. Make sure it is tight to prevent light leaks.

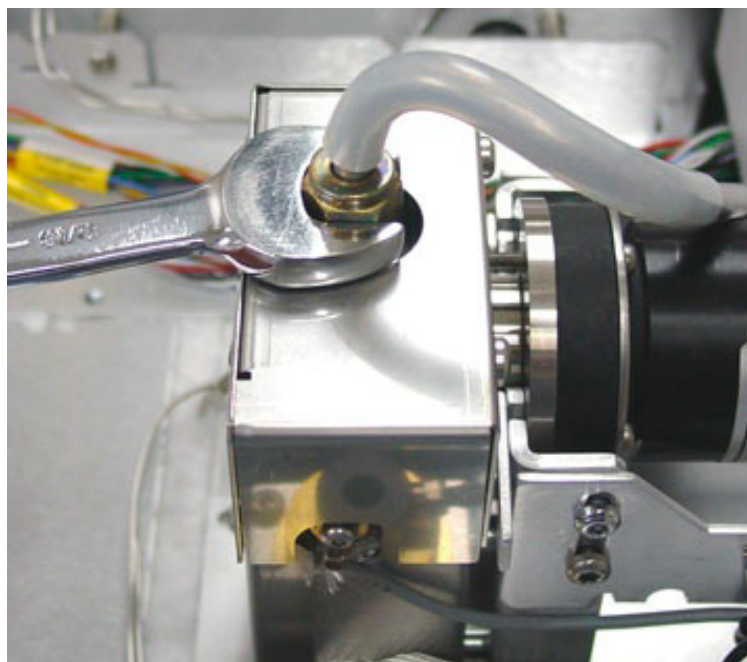


## Replacing the FPD

- 1 Before this can be done you must first disconnect some cables originating from the FPD board in the electronics panel.
- 2 From the side panel side, disconnect heater sensors to the Aux heater and to the back detector, the ignitor cable, the high voltage connector, and the PMT signal cable.

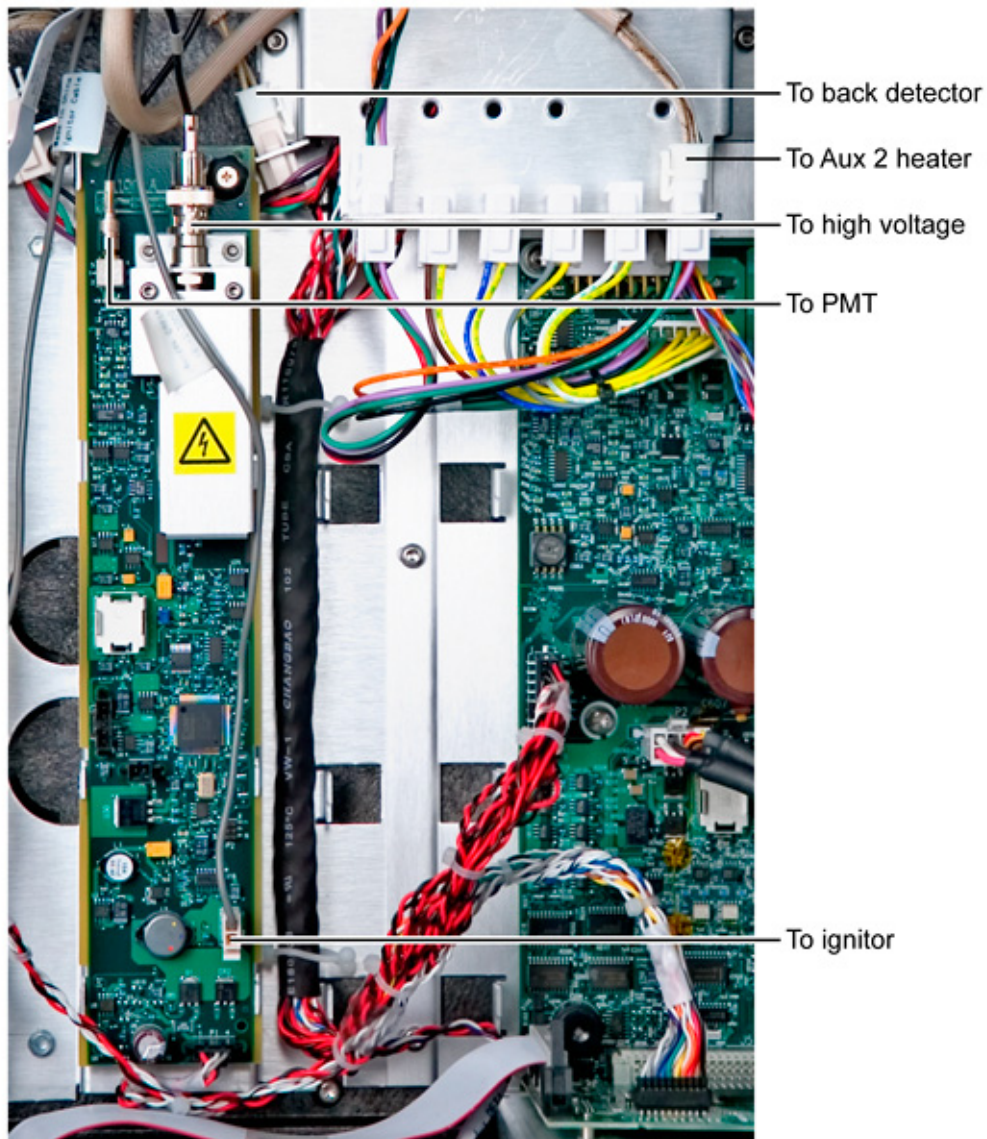

- 3 Remove the top four T-20 screws retaining the PMT assembly.

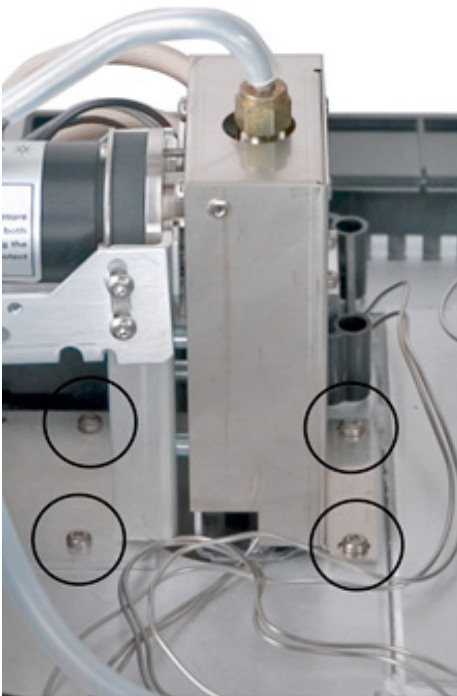

- 4 Lift up the entire assembly and lay it down on the top GC surface.
- 5 Loosen the one screw holding the adjacent cover and pivot the cover up and off.
- 6 The FPD is now ready for replacement, if required.

If not, reassembly is the reverse of these steps.

### **Bake out and run checkout test**

This part of the procedure bakes out the detector and restores conditions. The detector output will level off in about 1 hour after you restore the conditions. The detector output will continue to drift down slowly for about 24 hours.

- 1** Reinstall the column or columns.
- 2** Reinstall the drain tube to the FPD exhaust tube.
- 3** Restore the gases to the GC.
- 4** Restore the power to the GC.
- 5** Restore the conditions to the inlets and detectors, but turn off the flame of the FPD to prevent condensation.
- 6** After the FPD is at temperature for about 10 minutes, light the flame.
- 7** Bake out the detector. Set the oven and detector temperatures to 250°C for 15 minutes. Insure there is adequate gas flows through the column.
- 8** After the bakeout, restore the oven and FPD conditions.
- 9** Allow the output to level off.
- 10** Run a checkout and compare the results.



## Replacing the heater/sensor assemblies

Replacing the heater/sensor assemblies takes about 30 minutes. It requires you to turn off the GC, remove the vent, and the detector cover. (If using a dual-wavelength FPD--not available on all GC models--also remove the left PMT assembly.)

### CAUTION

When turning the GC off, turn off the flame first to prevent condensation from dripping into the jet and column.

You may wish to replace or check the FPD heater/sensor assemblies for the following reasons:

- One or both of the heaters or sensors are defective.
- The actual temperature reading on the display of the heaters is cycling more than 1 °C.

### Materials needed

- G1535-60610 Heater/Sensor assembly with short lead for the emissions block assembly.
- G1535-60620 Heater/Sensor assembly with long lead for the transfer line.


### Procedure

- 1 Turn off the flame.
- 2 Turn off the GC.
- 3 Remove the vent assembly and cover.
- 4 Put on an ESD wrist strap. Remove the right-side cover to access the GC electronics. Disconnect the heater/sensor leads from the auxiliary heater board and the connectors above the main board.
- 5 Carefully pull the cables up onto the top of the GC.

### Transfer line

- 1 Use a Torx T-10 driver to remove the screw and retainer clip holding the lower heater/sensor assembly. Remove the heater and sensor from the transfer line.
- 2 Remove the protective cap from the temperature sensor of the heater/sensor assembly with the short cable.
- 3 Insert the heater and sensor into the transfer line. Make sure the sensor is seated at the bottom of the hole. If not, the AUX temperature will wander above and below the detector setpoint.

- 4** Position the retainer clip over the heater/sensor assembly and install the screw.

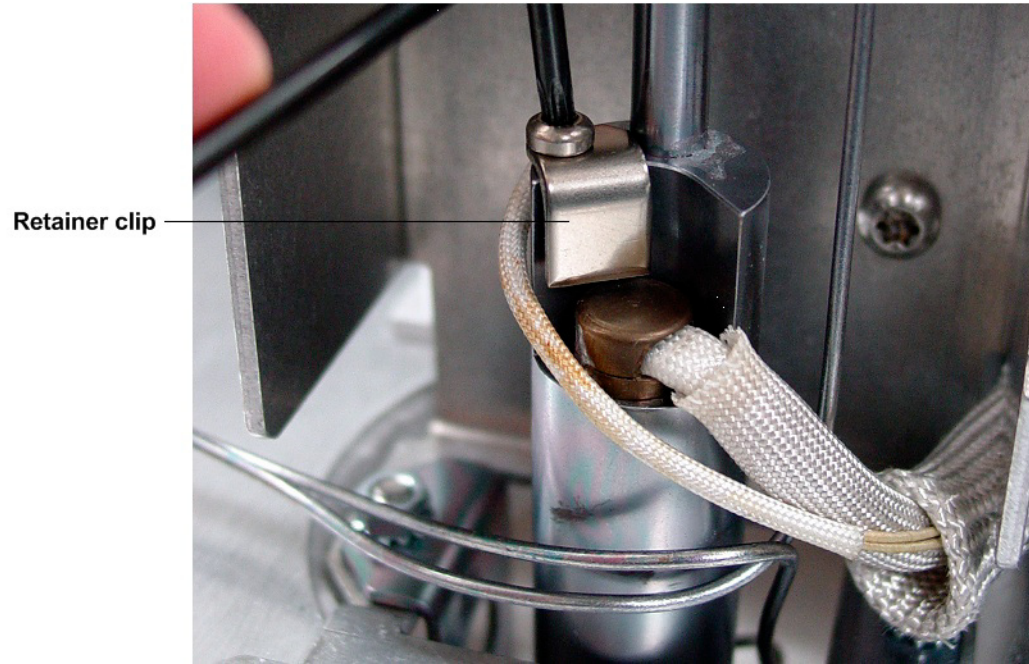



### **Emissions block assembly**


- 1** Remove the upper heater and temperature sensor from the emissions block assembly.
- 2** Remove the protective cap from the sensor of the heater/sensor assembly with the long cable.
- 3** Install the upper heater and sensor in the emissions block assembly.



#### Closing up

- 1 Route the heater/sensor cables out of the bracket as shown.
- 2 Replace the cover, vent assembly, and, on the dual-wavelength detector, the left PMT assembly. Refer to Accessing Heaters and Ignitor.
- 3 Carefully thread the heater/sensor cables into the electronics compartment.
- 4 Put on an ESD wrist strap. Connect the short cable as shown in Back-top panel cutouts. Connect the long cable to the auxiliary heater board as shown in Long cable connectors and leads.
- 5 Replace the right-side cover to the GC electronics compartment.
- 6 Turn on the GC. Confirm that the flame is off.
- 7 Restore the operating conditions.
- 8 Wait 20 minutes for the detector to heat up, then ignite the flame.

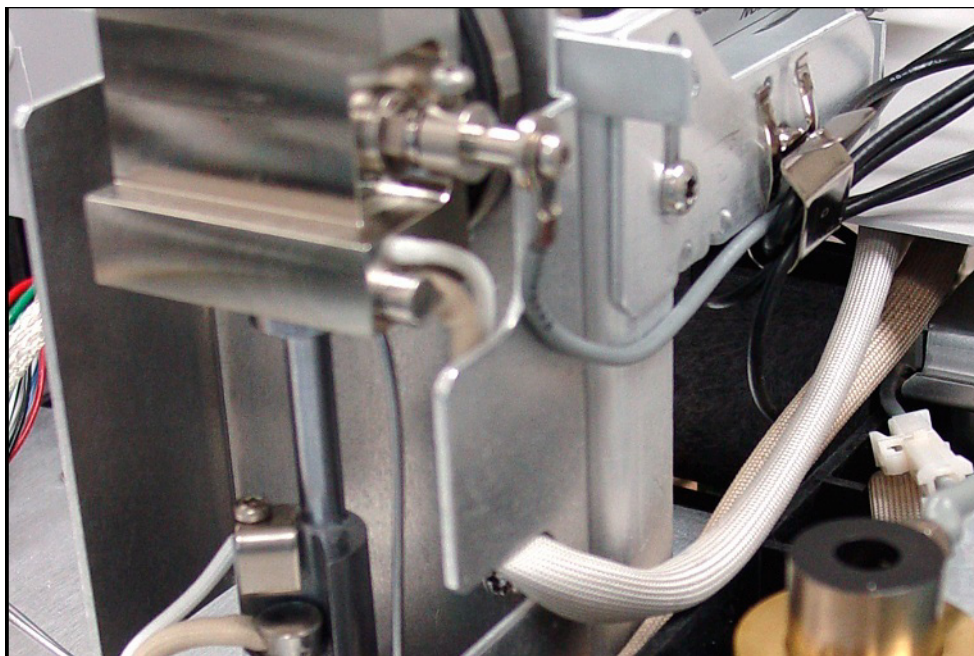

**Table 8** Long cable connectors and leads

Detector type	Location	Connector
Single-wavelength	Front	A1
Single-wavelength	Back	A2
Dual-wavelength	Back	A2

## Replacing the FPD signal board

If changing the GC configuration, see “[Changing the GC configuration](#)” for important information regarding GC methods. Then proceed with the steps below.

- 1 Cool down the oven and all heated zones.
- 2 Turn off all gas flows at their sources.
- 3 Turn off the GC main power switch and disconnect its power cord.

### WARNING

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**


---

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 4 Disconnect all attached electrical connections to the FPD board.
- 5 Turn thumb screw counterclockwise, lift screw to disengage, slide board down, and then lift out. (Note that the grounding point for the board is at the thumbscrew).
- 6 Replace board if required. (Note that the FPD board is longer than the FID board and is attached to the full length of the GC board receptacle.)
- 7 Replace connections in reverse order of disassembly.
- 8 Reconfigure the modified GC from the front panel keypad supplying requested information.

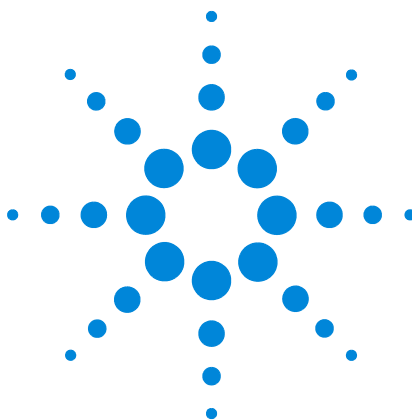

## 4 EPC modules

EPC overview	196
Repairing EPC modules	197
Accessing EPC modules	199
Replacing the EPC modules	200
Replacing EPC module proportional valves	203
Inlet modules	205
Detector modules	208
Pneumatics Control Module (PCM)	211
Calibrating the PCM interface	212
Changing the PCM pressure channel frit	213
Correcting PCM leaks	215
Aux EPCs	217
Frits for auxiliary pressure control channels	218
Updating PID Constants	221



## EPC overview

An Electronic Pneumatic Controller (EPC) module senses gas pressure or mass flow and controls that pressure or flow.

There are four kinds of EPC modules:

- Inlet flow modules
- Detector flow modules
- Pneumatic control modules (PCM)
- Aux EPCs

Inlet and detector modules are each designed for specific inlets and detectors.


A PCM has two control channels. One is intended for carrier gas. The other may function as either a forward- or back-pressure regulator, depending on the external plumbing.

An Aux EPC provides three independent forward-pressure regulation channels.



## Repairing EPC modules

The 7890 Series GC can control up to 6 EPC modules. Each has a particular purpose and is limited to a particular location in the EPC carrier. A TCD or AUX EPC module can be located in the optional 3rd detector carrier.

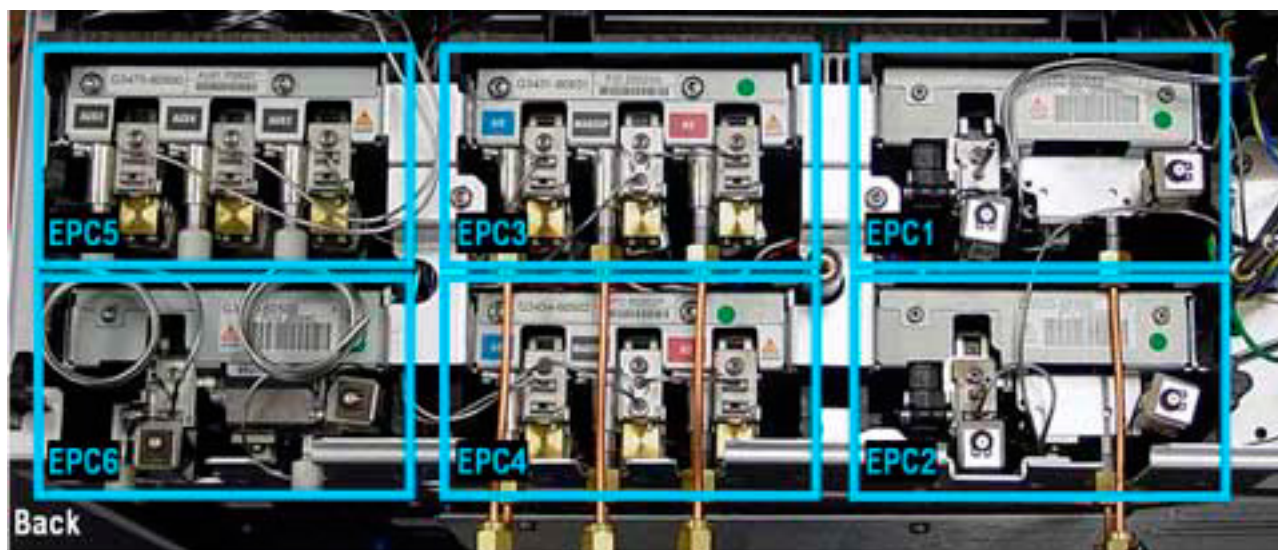

Repairs on EPC modules should be in this order:

- 1 Verify the EPC module is receiving 24 VDC power, the green LED is ON, and the EPC appears configured in the GC display.
- 2 Compare the Communication buss connection with an adjacent module that has a good power and signal connection.
- 3 Replace proportional valve or valves for flow or pressure problems.
- 4 Replace the module only after 1 and 2 (when dealing with power issues).

Repairs should be verified by the following:

- Zero the pressure channels before performing any verification tests.
- Inlets: With the inlet capped, perform pressure decay test. With column installed, compare calculated flow reading from the display with measured flow for column, septum purge, and split vent.

- Detectors: With the detector capped, compare calculated flow readings from the display with measured flow for each gas.

Failure modes for a proportional valve include:

- Sticking solenoid plunger causing variable flows.
- Electronic open causing no flow. Valve fail in a closed position.

## Accessing EPC modules

- 1 Turn off all gas flows at their sources.
- 2 Turn off the GC main power switch.

**CAUTION**

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 3 Remove the top rear cover by pressing in the two black side tabs and lifting the cover up and off.
- 4 Disconnect the gas plumbing to the module that is being replaced. If the module is towards the back (slots 2, 4, or 6), the gas plumbing on the front module may need to be removed.

**NOTE**

In some cases, it may be necessary to remove the upper rear cover plate. Do this by first loosening the three retaining T-20 screws along the bottom and then lifting the cover plate up and off.

---

## Replacing the EPC modules

- 1 Access the EPC module. See [“Accessing EPC modules”](#).

### **WARNING**

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**


---

### **CAUTION**


Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 2 EPC modules are secured by one or two clamps on the top edge. Remove the one or two retaining clamps holding the EPC module in place.
- 3 Remove the 1/8 inch Swagelok gas supply tubing connections (1, 2, or 3) from the EPC module.
- 4 Lift the EPC module up and remove the communication buss connector. The SSL, PTV, and VI EPC modules also require removal of an electrical connector to their split vent valves. This illustration (a split/splitless module) shows these features.

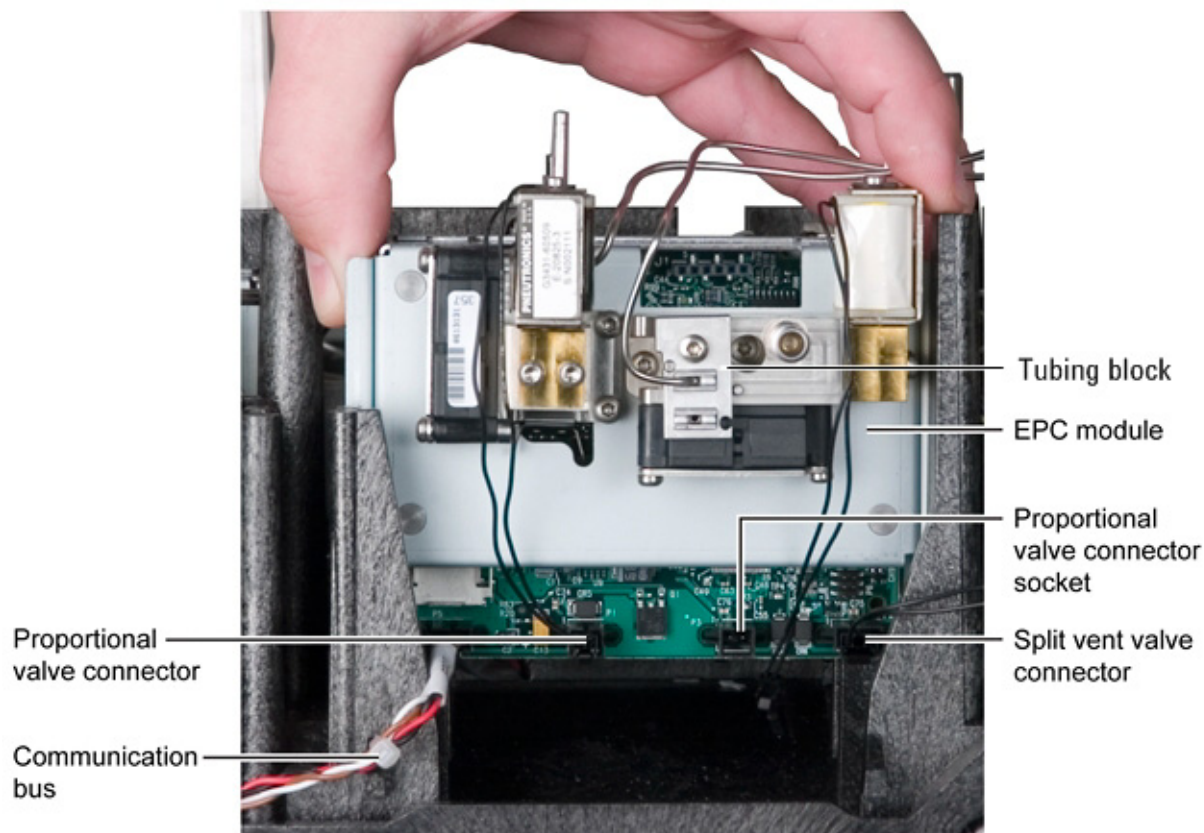

- 5 Remove the 1, 2, or 3 tubing block using a T-10 driver. This is the EPC module for an FID detector.

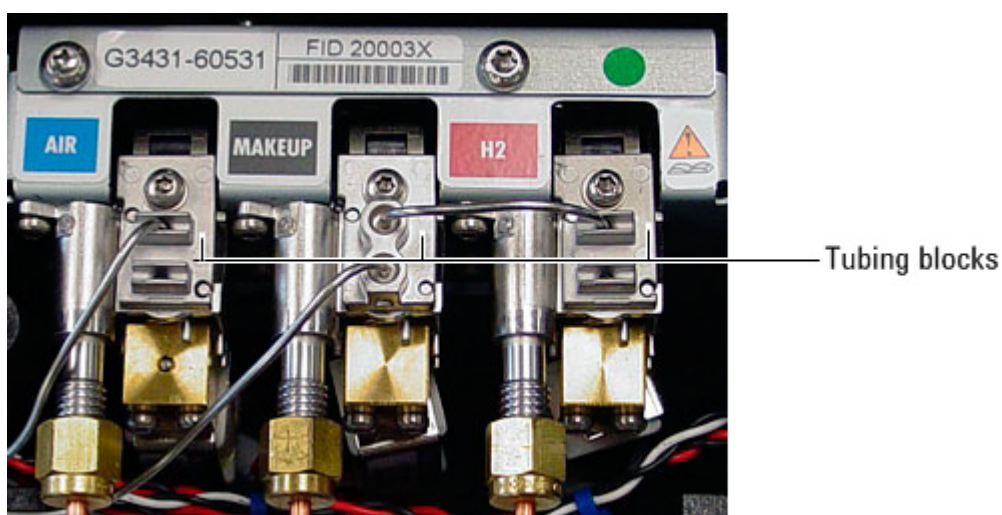


- 6 Remove the EPC module by lifting straight up.
- 7 Install the new EPC module. Remember to plug in the Communication buss. If a SSL, PTV, or VI inlet module, plug in the split vent valve.

- 8 Remove the O-rings for the pneumatic tubing blocks. Clean the sealing surface under the O-rings with a lint-free cloth. Install new O-rings.
- 9 Install the new pneumatic tubing blocks and tighten the screws.
- 10 Turn ON the GC and configure the EPC module. See [“Replacing a GC module”](#) for details.
- 11 Zero the flow and/or pressure sensors (**Options / Calibration**).
- 12 Verify the flow or pressure control.


## Replacing EPC module proportional valves

This is a generalized procedure for replacing all G3440A and G3440B EPC proportional valves. Select the valve and valve part number from “[EPC Module Parts](#)”. Valves come in a kit with spare O-rings and screws.

- 1 To protect the column and inlet liner, place the instrument into Inlet or GC maintenance mode. Let instrument go ready.
- 2 Turn off instrument.
- 3 If you need to lift the EPC module to access the wire connectors, remove 1 or 2 retaining brackets.
- 4 Lift the EPC module and disconnect the 2-wire connector from the printed circuit board. The figure shows an Aux EPC module, which has three proportional valves.

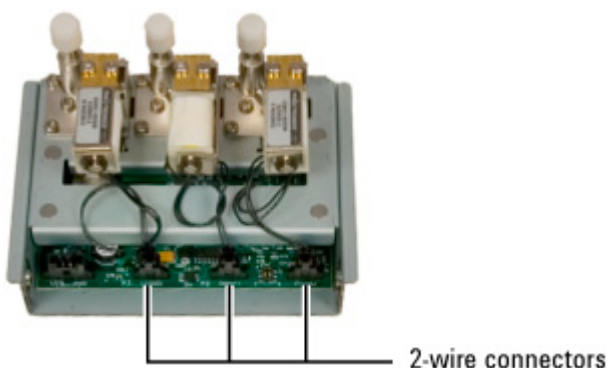


- 5 Using a T-10 driver, remove the two retaining screws holding the proportional valve to the EPC module. Save and reuse the screws if you did not get them in the kit.

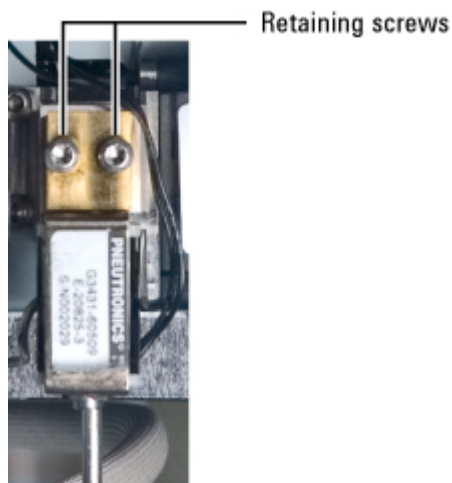


- 6 Remove the valve and the O-rings. Clean the seat of the O-ring on the flow block with a lint free cloth.

- 7 Insert 2 new O-rings.

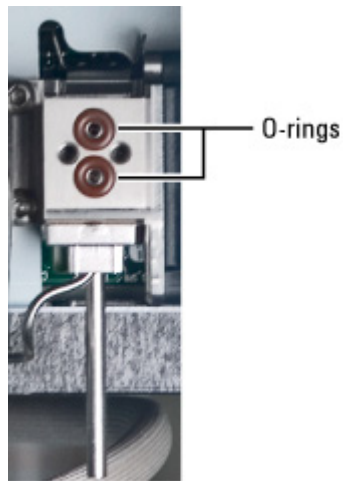

- 8 Except for the split vent valve, tie a loose knot in the proportional valve wire to shorten the wire.

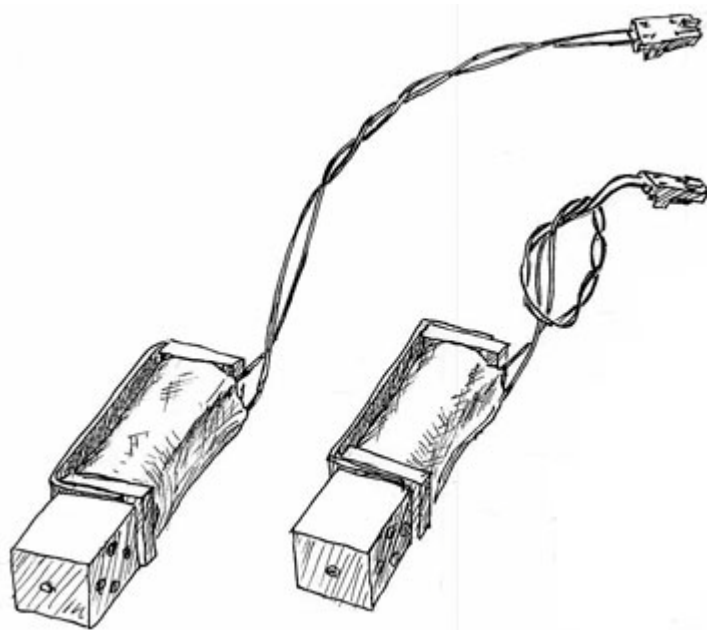


- 9 Attach the new proportional valve, the electrical connector, the EPC retaining clamps, and body panels.



## Inlet modules

These modules are essentially carrier gas regulators. The SSL, MMI, VI, and PTV modules have a split vent valve that sits next to the module..

This is a split/splitless inlet module:

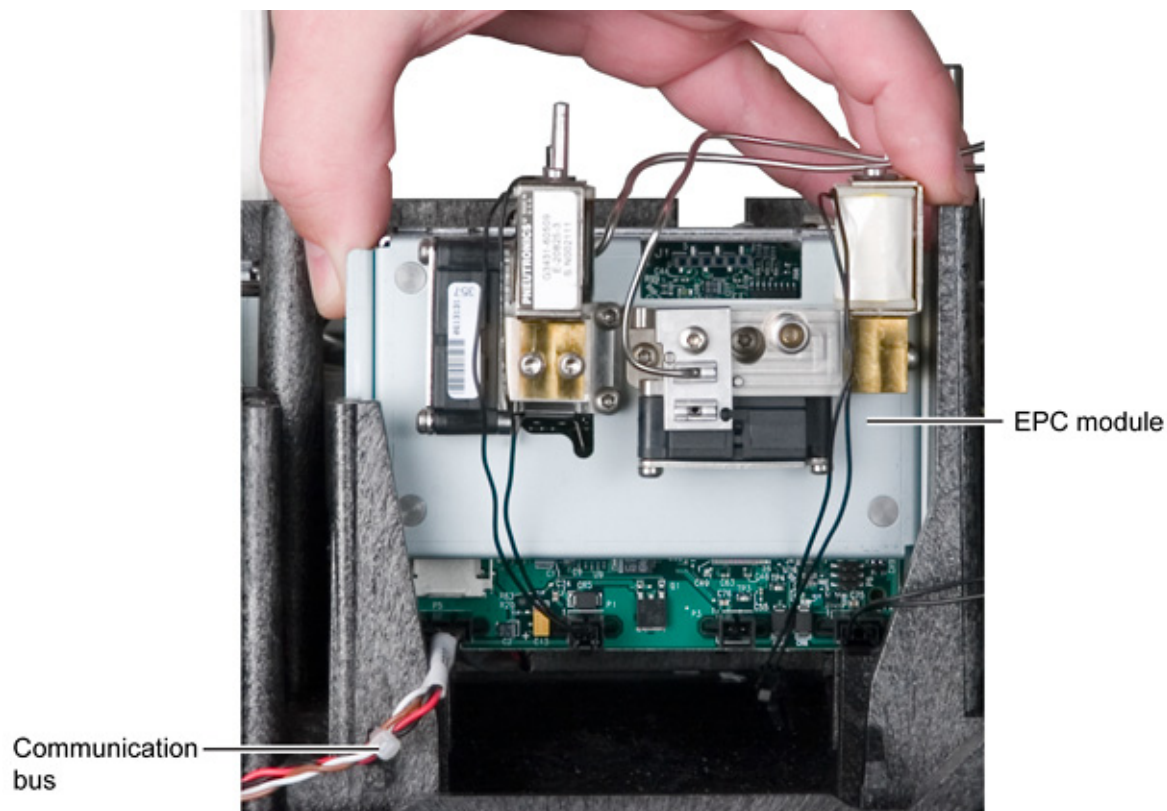

It consists of

- One or more proportional valves
- Gas tubing connections
- A printed circuit board
- A connection to the communications buss in the lower left corner

**Inlet modules**

Description	Part number	O-rings	Slot
COC	G3454-60554	3	1 or 2
PP	G3451-60551	3	1 or 2
PTV	G3500-60500	3	1 or 2
S/SL 100 PSI	G3452-60552	3	1 or 2
S/SL 150 PSI	G3452-60510	3	1 or 2
VI	G3504-60501	4	1 or 2
MMI	G3510-60500	3	1 or 2

The type of proportional valve varies with the module.


## Proportional valves

Description	Carrier	Septum purge	Split vent
COC	G3430-67013	G3430-67014	
PP	G3430-67013	G3430-67014	
PTV	G3430-67013	G3430-67014	G3430-67015
S/SL 100 PSI	G3430-67013	G3430-67014	G3430-67015
S/SL 150 PSI	G3430-67013	G3430-67014	G3430-67015
VI	G3430-67013	G3430-67014	G3430-67015
MMI	G3430-67013	G3430-67014	G3430-67015

Valves require 2 O-rings each (5180-4181, 12/pkg).

## Detector modules

Detector modules often supply multiple detector gases and a makeup gas. They always control flow by controlling pressure through a fixed restrictor. This is a typical FID module. FID, NPD and FPD modules control Air, Hydrogen, and column makeup. The Hydrogen and makeup are tied together and mixed. The Air has a separate connection.

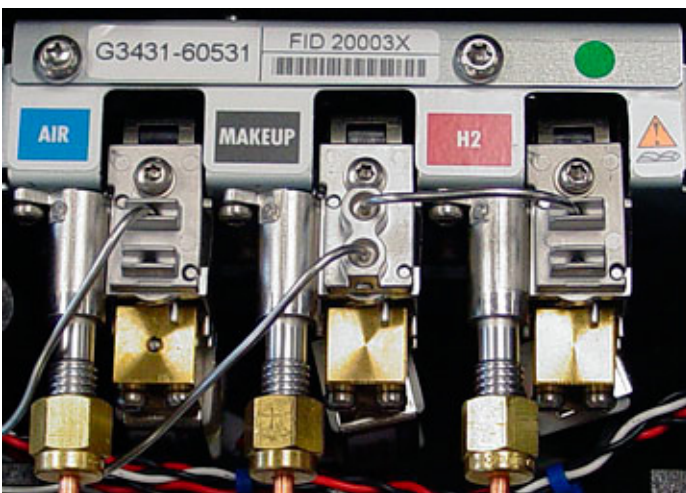

The TCD module supplies both carrier and reference gases and has a switching valve in the reference gas supply.

## Detector modules

Description	Part number	Slot
uECD	G3433-60533	3 or 4
FID	G3431-60531	3 or 4
FPD	G3435-60535	3 or 4
NPD	G3434-60502	3 or 4
TCD	G3432-60532	3 or 4 or AUX Det 2

## Proportional valves

Description	Air	Makeup	Hydrogen	Reference
uECD		G3430-67013		
FID	G3430-67013	G3430-67014	G3430-67014	
FPD	G3430-67013	G3430-67013	G3430-67013	
NPD	G3430-67013	G3430-67014	G3430-67014	
TCD		G3430-67013		G3430-67013

Valves require 2 O-rings each (5180-4181, 12/pkg).

## Pneumatics Control Module (PCM)

This module contains two channels.

- The PCM carrier channel has both pressure and flow sensors. It is functionally identical to the Purged Packed inlet flow module.
- The PCM auxiliary channel has only a pressure sensor. When plumbed normally (gas in at the threaded connector, out via the tubing), it provides forward pressure regulation. If the plumbing is reversed, it operates as a back-pressure regulator.

**Table 9** PCM module

Description	Part number	O-rings	Slot
PCM	G3476-60501	3	1, 2, 5, or 6

**Table 10** Proportional valves

Description	Carrier	AUX
PCM	G3430-67013	G3430-67016

## Calibrating the PCM interface

The interface's flow module contains a pressure sensor that must be zeroed after it is installed on the GC. Calibration ensures an accurate interface pressure display.

Do not connect the carrier gas to the flow module until you have zeroed the interface's pressure sensor. For more information on zeroing the pressure sensor, refer to the 7890A Advanced User Guide or 7890B Operation Manual.

- 1 If the gas supply is connected to the GC, turn off the supply at the source, then disconnect the supply line from the PCM inlet fitting.
- 2 Turn on the GC and wait 15 minutes to allow it to reach thermal equilibrium.
- 3 When the GC has reached thermal equilibrium, press **[Options]**, scroll to **Calibration** and press **[Enter]**.
- 4 Scroll to the module to be zeroed and press **[Enter]**.
- 5 Scroll to a zero line and press **[Info]**. The GC will remind you of the conditions necessary for zeroing that specific sensor.

**Flow sensors.** Verify that the gas is connected and flowing (turned on).

**Pressure sensors.** Disconnect the gas supply line at the back of the GC. Turning it off is not adequate; the valve may leak.

- 6 Press **[On/Yes]** to zero or **[Clear]** to cancel.
- 7 Turn off the GC.
- 8 Plumb the carrier gas to the flow module.
- 9 Turn on the GC.


If you were calibrating the flow sensor after replacing the PCM, check for leaks.




## Changing the PCM pressure channel frit

To change a frit:

- 1 Gather the following:
  - O-rings, package of 12, 5180-4181
  - Restrictor
- 2 Turn off the gas supply to the channel.
- 3 Remove the screw holding the tubing block with the output tubing.
- 4 Remove the tubing block. Remove the frit and O-ring. Be careful to avoid scratching the metal surfaces.
- 5 Remove the other O-ring as well. Replace it with a new O-ring.

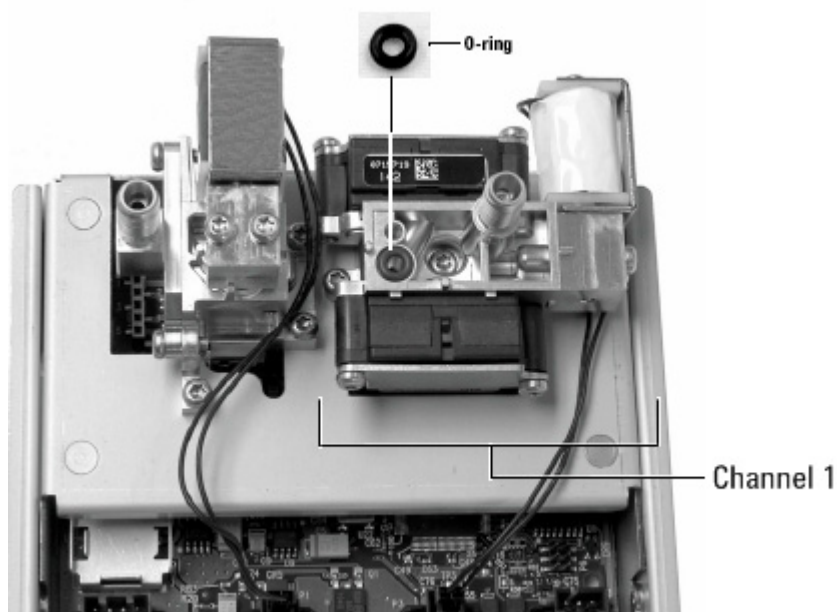

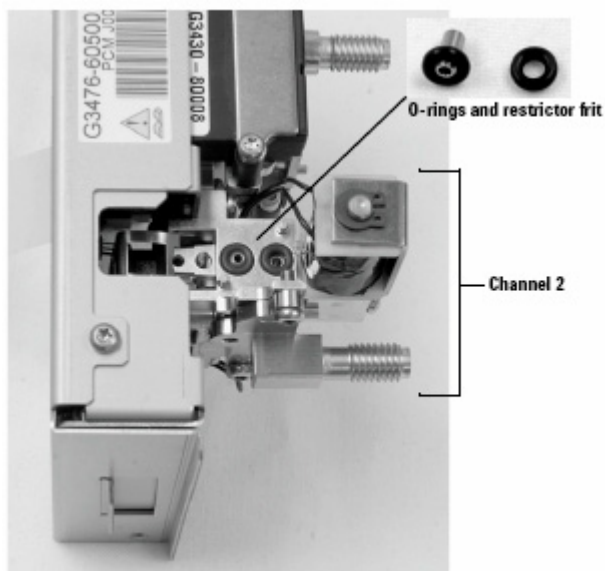

- 6 Place a new O-ring on the new frit and press it down into the block.
- 7 Place the block on the module and tighten the screw firmly.
- 8 Restore the gas supply.

## Correcting PCM leaks

- 1 Use the electronic leak detector to check all areas of the PCM that are potential sources of a leak.
- 2 Tighten any connections which are leaking and retest.

**NOTE**

If the pressure drop is now 0.05 psi/min or less, you can consider the interface system leak-free.

---

- 3 If the interface still leaks, continue to check for and correct leaks.

**Potential PCM leak areas**

Check the following areas when checking for leaks in the PCM.

**Flow channel**

- Two O-rings behind the block where the pneumatic lines enter the PCM
- Two O-rings for the valve in the PCM

**Pressure channel**

- Flatness of tubing block fitting
- Two O-rings behind the block where the pneumatic lines enter the PCM
- Two O-rings for the valve in the PCM

**Forward pressure valve leaks**

Occasionally an increase in pressure, rather than a decrease, may be observed. This is usually due to slight leakage into the module across the forward pressure control proportional valve. Although slight leaks of this nature do not create chromatographic problems, they may obscure other small leaks that do cause problems by allowing air into the system. The valves can leak at about 0.2 mL/min and be within specification.

To check for internal valve leakage (when leak testing the PCM only):

- 1 Remove the supply pressure at the carrier inlet fitting, and quickly cap the fitting with a solid 1/8-inch Vespel plug and a Swagelok nut.
- 2 Check the actual pressure on the display and monitor it for 5 minutes. Pressure loss should not be greater than 0.5 psi.

## Aux EPCs

Up to three of these modules can be installed, for a possible total of 9 pressure-regulated channels. Each channel has a block connector covering a frit with O-ring for the gas output via the block connector.

The module is shipped with a brown (low resistance) frit in each channel.

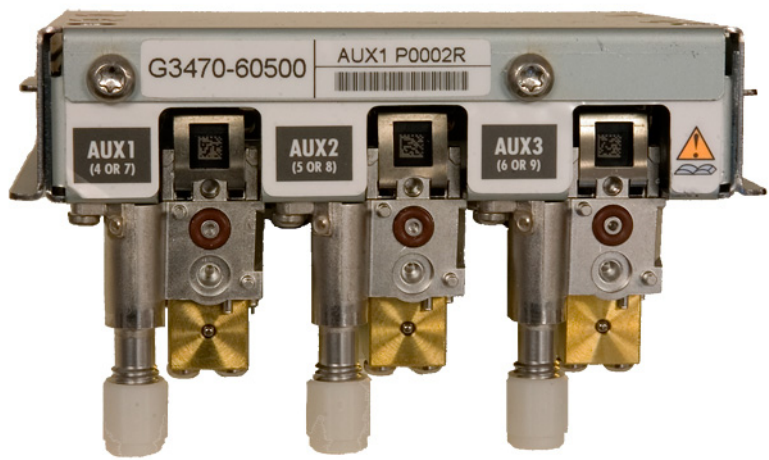


Table 11 Module

Description	Part number	Slot
AUX	G3470-60501	5 or 6 or AUX

Proportional valves

Table 12 Proportional valves

Description	Air, Makeup, Hydrogen
AUX EPC	G3430-67013


Table 13 Tubing weldment

Description	Part number
Pneumatic tubing with block connector and captive screw	G3430-60550

## Frits for auxiliary pressure control channels

The auxiliary pressure control channels in the Aux EPC module are used to provide a flow to various devices including a purged splitter or a QuickSwap or a headspace sampler. The PCM Aux channel can be used to provide a flow in forward pressure control mode.

A replaceable frit is used in these auxiliary channels to optimize control for the range of flows. The image below shows one of the three channels.

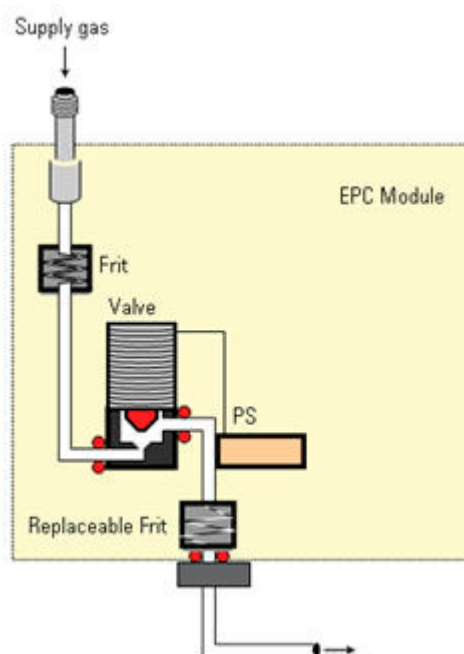**NOTE**

Different frits may require different PID constants. See [“Updating PID Constants”](#).

This table lists the 4 types of frits, their flow characteristics, and common uses.

Kit	Component	Marking	Flow characteristic	Resistance	Often used with
No	G3430-80061	1 ring Brown	400 +/- 30 SCCM AIR @ 40 PSIG	Low	FID Air, QuickSwap, Splitter, Deans Switch
Yes	G3430-80062	2 rings Red	30 +/- 1.5 SCCM H2 @ 15 PSIG	Medium	FID Hydrogen
Yes	G3430-80063	3 rings Blue	3.33 +/- 0.3 SCCM H2 @ 15 PSIG	High	NPD Hydrogen
Yes	G3430-20011	None	No restriction	Zero	Headspace vial pressurization

Notes for this table

- The Air frit, G3430-80061, ships in each channel in the AUX EPC module. No frit ships in the PCM Aux channel.
- Yes indicates this part is included in the AUX module ship kit and the PCM ship kit.
- Always use a new O-ring. 5180-4181, O-rings 12/pk

### Procedure for replacing frit

- 1 Select the appropriate frit.
- 2 Gather a new O-ring, a T-10 Torx driver, and a tweezer.
- 3 Loosen the captured screw for the channel that you want to change.
- 4 Lift up the pneumatic tubing block.
- 5 Replace the O-ring and frit.
- 6 Replace the pneumatic block and tight the screw.





## Updating PID Constants

Pneumatic PID constants are the coefficients used by the EPC modules to control flows and pressures. Certain gases or special applications (such as QuickSwap, Headspace vial pressurization, or splitter and backflush applications) require different PIDs than provided at the factory. For the Agilent 7890 Series GC, you must use G4600AA Agilent Instrument Utilities software to download new PIDs to the EPC modules.

The table below summarizes custom PID values required for selected applications. Note that if updating an Aux EPC module, you will need to change the frit for the channel used.

Application	Module	AUX Frit	Select Available PID Values
QuickSwap	AUX EPC	G3430-80061 brown or 1 ring	<b>Quickswap</b>
Purged Splitter and Deans Switch when using Backflush	AUX EPC	G3430-20011 no color or no rings	<b>Quickswap</b>
Purged Splitter and Deans Switch	AUX EPC	G3430-80061 brown or 1 ring	<b>Standard</b>
Headspace vial pressurization	AUX EPC	G3430-20011 no color or no rings	<b>AUX_EPC_Headspace</b>
Headspace sampling loop	PCM in backpressure control	NA	<b>PCM_Headspace</b>

### NOTE

The driver for the 7890Series GC will calculate backflush conditions. The calculation of **Hold Times** assumes that a zero restrictor is installed in the AUX EPC channel. If you use an FID Air restrictor the actual hold time will be less than what the software displays.

## Updating PID constants

Follow these steps to update the PID constants for any EPC module (inlet, detector, AUX, or PCM). Use these same directions to apply special PIDs for use with QuickSwap, Capillary Flow Technology devices, and Headspace samplers, or to restore factory defaults.

The PID Constants screen displays (for the selected instrument):

- **Current**—lists devices (and the associated Module for each) installed on the selected GC
- **Version**—version number of the currently-installed PID values for the selected module
- **Available**—PID files available for upload to the selected module
- **Selected Version**—version number of the Available PID values selected for update
- **Status**—indicates (after an update) whether or not the PID values were uploaded to the module successfully

To set your PID values:

- 1 Close any online data system session for this instrument.
- 2 Open Instrument Utilities.
- 3 Select **Firmware Update** from the software explorer pane.
- 4 Select the PID Constants tab.

- 5 From the **Instruments** list, select the GC with the module(s) you wish to update, and click **Connect**. The list of modules for this GC will display. The list may not appear immediately.

**Tools**


PID Constants

This feature requires a live connection to the instrument which may already be in use

Disconnect

Module	Current	Version	Available	Selected Version	Status	Set Values
EPC1	<INLET_SS>	500	SS_100_Inlet	514		Update
EPC2	No device present					Update
EPC3	No device present					Update
EPC4	No device present					Update
EPC5	No device present					Update
EPC6 - Aux 1	<AUX_APC>	513	AUX_EPC	513		Update
EPC6 - Aux 2	<AUX_APC>	513	AUX_EPC	513		Update
EPC6 - Aux 3	<AUX_APC>	513	AUX_EPC	513		Update
Aux Detector	No device present					Update

- 6 Locate the **Module** you wish to update. If updating an EPC AUX or PCM module, determine the specific channel to be updated.
- If you are unsure of the channel to be updated, remove the pneumatics cover and trace the wiring from the device to the EPC module. Use the figure below to determine the module's EPC location.

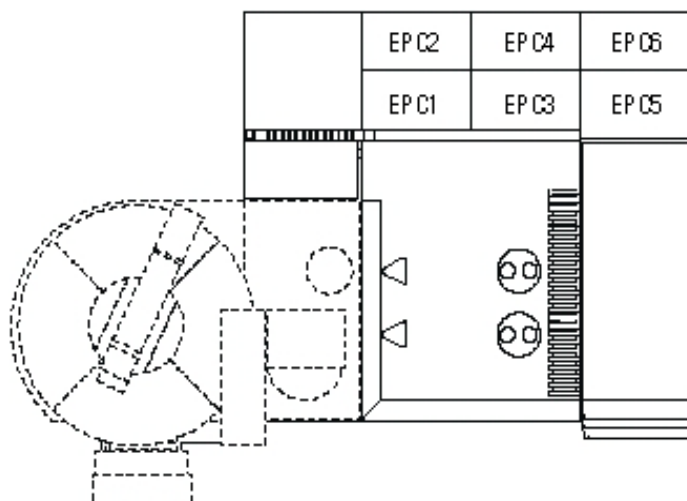

## 4 EPC modules

- 7 Select the row of the EPC module (and channel) that you need to update. This is most often an AUX\_EPC module.

**Tools**

PID Constants

This feature requires a live connection to the instrument which may already be in use

Disconnect

Module	Current	Version	Available	Selected Version	Status	Set Values
EPC1	<INLET_SS>	500	SS_100_Inlet	514		Update
EPC2	No device present					Update
EPC3	No device present					Update
EPC4	No device present					Update
EPC5	No device present					Update
EPC6 - Aux 1	<AUX_APC>	513	AUX_EPC	513		Update
EPC6 - Aux 2	<AUX_APC>	513	AUX_EPC	513		Update
EPC6 - Aux 3	<AUX_APC>	513	AUX_EPC	513		Update
Aux Detector	No device present					Update

- 8 From the **Available** selection box, select the appropriate set of values to download. Note the Selected Version number that appears. In the example below, we are updating Aux 1 with QuickSwap PIDs.

**Tools**

Vapor Volume Calculator Pressure Flow Calculator Method Translator PID Constants

This feature requires a live connection to the instrument which may already be in use

Disconnect

Module	Current	Version	Available	Selected Version	Status	Set Values
EPC1	<INLET_SS>	500	SS_100_Inlet	514		Update
EPC2	No device present					Update
EPC3	No device present					Update
EPC4	No device present					Update
EPC5	No device present					Update
EPC6 - Aux 1	<AUX_APC>	513	QuickSwap	514		Update
EPC6 - Aux 2	<AUX_APC>	513	AUX_EPC	513		Update
EPC6 - Aux 3	<AUX_APC>	513	AUX_EPC	513		Update
Aux Detector	No device present					Update

- 9 When finished specifying your desired values, click **Update** in the row(s) of the module(s) you wish to update. This will download the PIDs to the appropriate EPC modules.

- 10 If the download is successful, the status will read **Values Set** (see below). If the file is not accepted by Instrument Utilities, the status will read **Invalid File**. Note that these Status messages are the only feedback on the success of the PID update.

The following screen capture illustrates updating PIDs for QuickSwap.

**Tools**

PID Constants

This feature requires a live connection to the instrument which may already be in use

Disconnect

Module	Current	Version	Available	Selected Version	Status	Set Values
EPC1	<INLET_SS>	500	SS_100_Inlet	514		Update
EPC2	No device present					Update
EPC3	No device present					Update
EPC4	No device present					Update
EPC5	No device present					Update
EPC6 - Aux 1	<AUX_APC>	513	QuickSwap	514	Values Set	Update
EPC6 - Aux 2	<AUX_APC>	513	AUX_EPC	513		Update
EPC6 - Aux 3	<AUX_APC>	513	AUX_EPC	513		Update
Aux Detector	No device present					Update

- 11 The software writes a log entry for the instrument.
- 12 Power cycle the GC to activate the new PIDs.
- 13 If using Instrument Utilities, close the software to remove your CE key from the customer's system.

## No positive feedback on reconnect

### NOTE

The ChemStation driver for the 7890A will calculate backflush conditions. The calculation of **Hold Times** assumes that a zero restrictor is installed in the AUX EPC channel. If you use an FID Air restrictor the actual hold time will be less than what the software displays. The only feedback on the success of the update is the message Values Set and an entry in the software log file.

If you reconnect, you will not see any confirmation that you updated a PID. For example, if you update one or more AUX EPC channels to QuickSwap PID values, the screen does not tell you that the channel has QuickSwap values.

However, the **Version** number may have changed for all channels. This number is read from the EPC module. There is only one version number in each EPC module. If you update one channel of an AUX EPC module, the version number will be the latest update. Because there are not a version number for each channel, it will look like you updated all channels. This is not correct. The current version does not provide you positive feedback that you updated one channel. Here is an example screen after the update of one channel to QuickSwap.

**Tools**

PID Constants

This feature requires a live connection to the instrument which may already be in use

Disconnect

Module	Current	Version	Available	Selected Version	Status	Set Values
EPC1	<INLET_SS>	500	SS_100_Inlet	514		Update
EPC2	No device present					Update
EPC3	No device present					Update
EPC4	No device present					Update
EPC5	No device present					Update
EPC6 - Aux 1	<AUX_APC>	514	AUX_EPC	513		Update
EPC6 - Aux 2	<AUX_APC>	514	AUX_EPC	513		Update
EPC6 - Aux 3	<AUX_APC>	514	AUX_EPC	513		Update
Aux Detector	No device present					Update

The **Selected Version** is read from the eeprom file on your PC.


## Updating PID constants available for use by Instrument Utilities

If at any time the PID constant values are updated, follow the procedure below to update the files within the Instrument Utilities directory. Perform this procedure **before** accessing the PID Constants interface.

- 1 Obtain the latest EEPROM files.

Browse to where the Instrument Utilities is installed (default is **C:\Program Files\Agilent\Lab Monitor**).

- 1 In this Instrument Utilities directory, browse to **firmware\pneumatic-eprom**.
- 2 If the latest EEPROM file dates are newer than the files in this directory:
  - a Delete the old .xml files.
  - b Copy the new .xml files into the **firmware\pneumatic-eprom** directory. The directory will look similar to the one below.

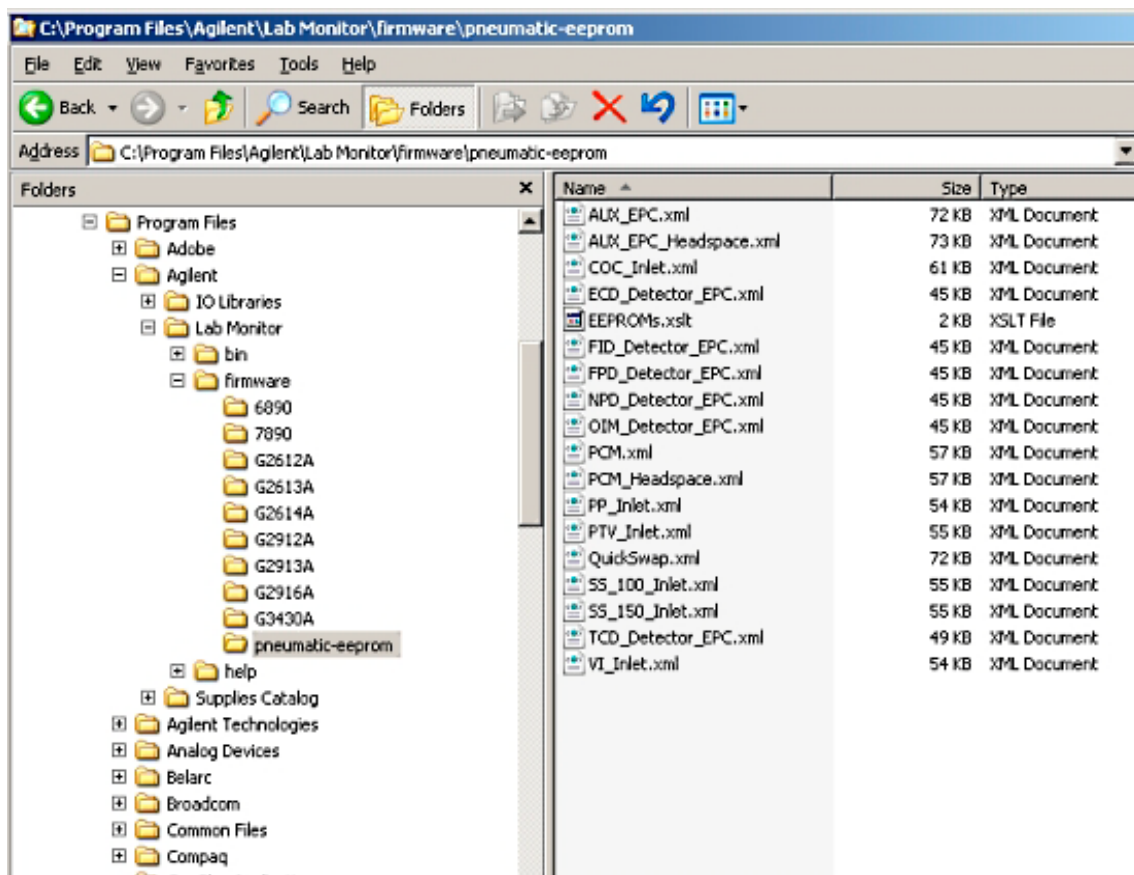

- c See "Updating PID Constants" to download these values to the EPC modules.



## 5 Performance Verification

About Chromatographic Checkout	229
To Prepare for Chromatographic Checkout	230
To Check FID Performance	232
To Check NPD Performance	238
To Check TCD Performance	244
To Check uECD Performance	250
To Check FPD Performance (Sample 5188-5953)	256
To Check FPD Performance (Sample 5188-5245, Japan)	264





## About Chromatographic Checkout

The tests described in this section provide basic confirmation that the GC and detector can perform comparably to factory condition. However, as detectors and the other parts of the GC age, detector performance can change. The results presented here represent typical outputs for typical operating conditions and are not specifications.

The tests assume the following:

- Use of an automatic liquid sampler. If not available, use a suitable manual syringe instead of the syringe listed.
- Use of a 10- $\mu$ L syringe in most cases. However, a 5- $\mu$ L syringe is an acceptable substitute for the 1- $\mu$ L injections described here.
- Use of the septa and other hardware (liners, jets, adapters, and so forth) described. If you substitute other hardware, performance can vary.

## To Prepare for Chromatographic Checkout

Because of the differences in chromatographic performance associated with different consumables, Agilent strongly recommends using the parts listed here for all checkout tests. Agilent also recommends installing new consumable parts whenever the quality of the installed ones is not known. For example, installing a new liner and septum ensures that they will not contribute any contamination to the results.

- 1 Check the indicators/dates on any gas supply traps. Replace/recondition expended traps.
- 2 Install new consumable parts for the inlet and prepare the correct injector syringe (and needle, as needed).

**Table 14** Recommended parts for checkout

Recommended part for checkout	Part number
<b>Split splitless inlet</b>	
Syringe, 10- $\mu$ L	5181-1267
O-ring	5188-5365
Septum	5183-4757
Liner	5062-3587 or 5181-3316
<b>Multimode inlet</b>	
Syringe, 10- $\mu$ L	5181-1267
O-ring	5188-6405
Septum	5183-4757
Liner	5188-6568
<b>Packed column inlet</b>	
Syringe, 10- $\mu$ L	5181-1267
O-ring	5080-8898
Septum	5183-4757
<b>Cool on-column inlet</b>	
Septum	5183-4758
Septum nut	19245-80521
Syringe, 5- $\mu$ L on-column	5182-0836
0.32 mm needle for 5- $\mu$ L syringe	5182-0831

Recommended part for checkout	Part number
7683B Needle support assembly for .250/.320 mm injections	G2913-60977
Insert, fused silica, 0.32-mm id	19245-20525
<b>PTV inlet</b>	
Syringe, 10- $\mu$ L -- for septum head	5181-1267
Syringe, 10- $\mu$ L, 23/42/HP -- for septumless head	5181-8809
Inlet adapter, Graphpak-2M	5182-9761
Silver seal for Graphpak-2M	5182-9763
Glass liner, multibaffle	5183-2037
PTFE ferrule (septumless head)	5182-9748
Microseal replacement (if installed)	5182-3444
Ferrule, Graphpak-3D	5182-9749

## To Check FID Performance

- 1 Gather the following:
  - Evaluation column, HP-5 30 m – 0.32 mm – 0.25  $\mu$ m (19091J-413)
  - FID performance evaluation (checkout) sample (5188-5372)
  - Chromatographic-grade isooctane (wash solvent)
  - 4-mL solvent and waste bottles or equivalent for autoinjector
  - 2-mL sample vials or equivalent for sample
  - Inlet and injector hardware
- 2 Verify the following:
  - Capillary column jet installed. If not, select and install a capillary column jet.
  - Capillary column adapter installed (adaptable FID only). If not, install it.
  - Chromatographic-grade gases plumbed and configured: helium as carrier gas, nitrogen, hydrogen, and air.
  - Empty waste vials loaded in sample turret.
  - 4-mL solvent vial with diffusion cap filled with isooctane and inserted in Solvent A injector position.
- 3 Replace consumable parts (liner, septum, traps, syringe, and so forth) as needed for the checkout.
- 4 Install the evaluation column. (See the procedure for the SS, MMI, PP, COC, or PTV in the Maintenance manual.)
  - Bake out the evaluation column for at least 30 min at 180 °C. (See the procedure for the SS, MMI, PP, COC, or PTV in the Maintenance manual.)
  - Be sure to configure the column.
- 5 Check the FID baseline output. The output should be between 5 pA and 20 pA and relatively stable. (If using a gas generator or ultra pure gas, the signal may stabilize below 5 pA.) If the output is outside this range or unstable, resolve this problem before continuing.
- 6 If the output is too low:
  - Check that the electrometer is on.
  - Check that the flame is lit.
- 7 Create or load a method with these parameter values.

**Table 15** FID checkout conditions

<b>Column and sample</b>	
Type	HP-5, 30 m x 0.32 mm x 0.25 µm (19091J-413)
Sample	FID checkout 5188-5372
Column flow	6.5 mL/min
Column mode	Constant flow
<b>Split/splitless inlet</b>	
Temperature	250 °C
Mode	Splitless
Purge flow	40 mL/min
Purge time	0.5 min
Septum purge	3 mL/min
Gas saver	Off
<b>Multimode inlet</b>	
Mode	Splitless
Inlet temperature	75 °C
Initial time	0.1 min
Rate 1	720 °C/min
Final temp 1	250°C
Final time 1	5.0 min
Purge time	1.0 min
Purge flow	40 mL/min
Septum purge	3 mL/min
<b>Purged packed inlet</b>	
Temperature	250 °C
Septum purge	3 mL/min
<b>Cool on-column inlet</b>	
Temperature	Oven Track
Septum purge	15 mL/min
<b>PTV inlet</b>	
Mode	Splitless
Inlet temperature	75 °C

<b>Column and sample</b>	
Initial time	0.1 min
Rate 1	720 °C/min
Final temp 1	350 °C
Final time 1	2 min
Rate 2	100 °C/min
Final temp 2	250 °C
Final time 2	0 min
Purge time	0.5 min
Purge flow	40 mL/min
Septum purge	3 mL/min
<b>Detector</b>	
Temperature	300 °C
H2 flow	30 mL/min
Air flow	400 mL/min
Makeup flow (N2)	25 mL/min
Lit offset	Typically 2 pA
<b>Oven</b>	
Initial temp	75 °C
Initial time	0.5 min
Rate 1	20 °C/min
Final temp	190 °C
Final time	0 min
<b>ALS settings (if installed)</b>	
Sample washes	2
Sample pumps	6
Sample wash volume	8
Injection volume	1 µL
Syringe size	10 µL
Solvent A pre washes	2
Solvent A pre washes	2
Solvent A wash volume	8

Column and sample	
Solvent B pre washes	0
Solvent B post washes	0
Solvent B wash volume	0
Injection mode (7693A)	Normal
Airgap volume (7693A)	0.20
Viscosity delay	0 s
Inject dispense speed (7693A)	6000
Plunger speed	Fast (Slow for cool on-column inlet)
PreInjection dwell	0
PostInjection dwell	0
Manual injection	
Injection volume	1 µL
Data system	
Path	Setting
Method and Run Control View: Instrument>Edit Parameters>Options	Keep instrument keyboard locked after method is loaded? = No
Method and Run Control View: Instrument>Edit Parameters>Signals	Assign Signal 1 to the detector Choose Save Data All Data Rate = 5 Hz


- 8 If using a data system, prepare the system to perform one run using the loaded checkout method. Make sure that the data system will output a chromatogram.

- 9 Start the run.

If performing an injection using an autosampler, start the run using the data system or press **[Start]** on the GC.

If performing a manual injection (with or without a data system):

- a Press **[Prep Run]** to prepare the inlet for splitless injection.
  - b When the GC becomes ready, inject 1 µL of the checkout sample and press **[Start]** on the GC.
- 10 The following chromatogram shows typical results for a new detector with new consumable parts installed and nitrogen makeup gas.

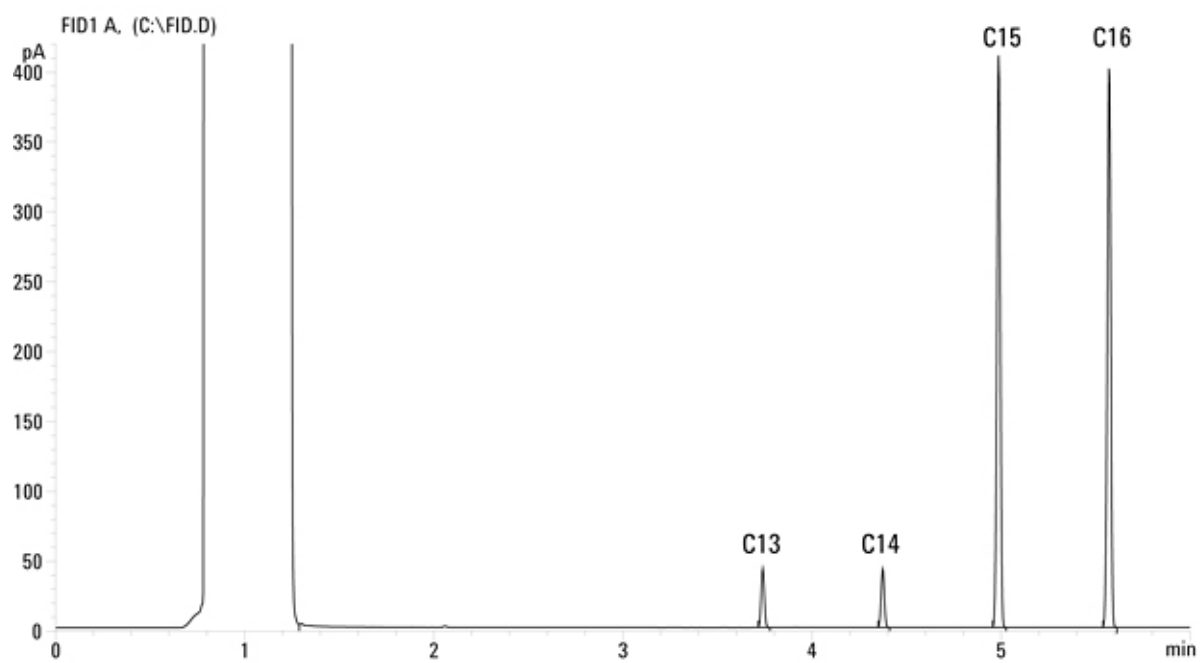



## Typical values

The following are typical values for the detector results. The results for a detector depend on site preparation related factors. These are not pass/fail criteria.

Components	Typical area counts
C13 and C14	46 pA-sec (area counts)
C15 and C16	460 pA-sec

FID ASTM noise  $\leq .0382$  pA

Repeatability (6 runs, first deleted to avoid first run anomalies):  
RSD% of RT = 0.01%;  $<0.0008$  m of area  $<1\%$ .

Split: C13/C16 and C14/C16 sensitivity ratios:  $>0.95$  and  $<1.05$

## To Check NPD Performance

- 1 Gather the following:
  - Evaluation column, HP-5 30 m – 0.32 mm – 0.25  $\mu$ m (19091J-413)
  - NPD performance evaluation (checkout) sample (18789-60060)
  - 4-mL solvent and waste bottles or equivalent for autoinjector
  - Chromatographic-grade isooctane
  - 2-mL sample vials or equivalent for sample
  - Inlet and injector hardware
- 2 Verify the following:
  - Capillary column jet installed. If not, select and install a capillary column jet.
  - Capillary column adapter installed. If not, install it.
  - Chromatographic-grade gases plumbed and configured: helium as carrier gas, nitrogen, hydrogen, and air.
  - Empty waste vials loaded in sample turret.
  - 4-mL vial with diffusion cap filled with isooctane and inserted in Solvent A injector position.
- 3 If present, remove any protective caps from the inlet manifold vents.
- 4 Replace consumable parts (liner, septum, traps, syringe, and so forth) as needed for the checkout.
- 5 Install the evaluation column. (See the procedure for the SS, MMI, PP, COC, or PTV in the Maintenance manual.)
  - Bake out the evaluation column for at least 30 min at 180 °C. (See the procedure for the SS, MMI, PP, COC, or PTV in the Maintenance manual.)
  - Be sure to configure the column.
- 6 Create or load a method with these parameter values.

**Table 16** NPD checkout conditions

<b>Column and sample</b>	
Type	HP-5, 30 m x 0.32 mm x 0.25 µm (19091J-413)
Sample	NPD checkout 18789-60060
Column mode	Constant flow
Column pressure	6.5 mL/min (helium)
Split/splitless inlet	
Temperature	200 °C
Mode	Splitless
Purge flow	60 mL/min
Purge time	0.75 min
Septum purge	3 mL/min
Multimode inlet	
Mode	Splitless
Inlet temperature	60 °C
Initial time	0.75 min
Rate 1	720 °C/min
Final temp 1	350 °C
Final time 1	2 min
Purge time	1.0 min
Purge flow	60 mL/min
Septum purge	3 mL/min
Packed column inlet	
Temperature	200 °C
Septum purge	3 mL/min
Cool on-column inlet	
Temperature	Oven track
Septum purge	15 mL/min
PTV inlet	
Mode	Splitless
Inlet temperature	60 °C
Initial time	0.1 min

Column and sample	
Rate 1	720 °C/min
Final temp 1	350 °C
Final time 1	2 min
Rate 2	100 °C/min
Final temp 2	250 °C
Final time 2	0 min
Purge time	0.75 min
Purge flow	60 mL/min
Septum purge	3 mL/min
Detector	
Temperature	300 °C
H2 flow	3 mL/min
Air flow	60 mL/min
Makeup flow (N2)	Makeup + column = 10 mL/min
Output	30 display units (30 pA)
Oven	
Initial temp	60 °C
Initial time	0 min
Rate 1	20 °C/min
Final temp	200 °C
Final time	3 min
ALS settings (if installed)	
Sample washes	2
Sample pumps	6
Injection volume	1 µL
Syringe size	10 µL
PreInj Solvent A Washes	2
PreInj Solvent B Washes	0
PostInj Solvent A Washes	2
PostInj Solvent B Washes	0
Viscosity delay	0

Column and sample	
Plunger speed	Fast (Slow for cool on-column inlet)
PreInjection dwell	0
PostInjection dwell	0
Manual injection	
Injection volume	1 µL
Data system	
Path	Setting
Method and Run Control View: Instrument>Edit Parameters>Options	Keep instrument keyboard locked after method is loaded? = No
Method and Run Control View: Instrument>Edit Parameters>Signals	Assign Signal 1 to the detector Choose Save Data All Data Rate = 5 Hz


- 7 If using a data system, prepare the system to perform one run using the loaded checkout method. Make sure that the data system will output a chromatogram.

- 8 Start the run.

If performing an injection using an autosampler, start the run using the data system or press **[Start]** on the GC.

If performing a manual injection (with or without a data system):

- a Press **[Prep Run]** to prepare the inlet for splitless injection.
  - b When the GC becomes ready, inject 1 µL of the checkout sample and press **[Start]** on the GC.
- 9 The following chromatogram shows typical results for a new detector with new consumable parts installed.

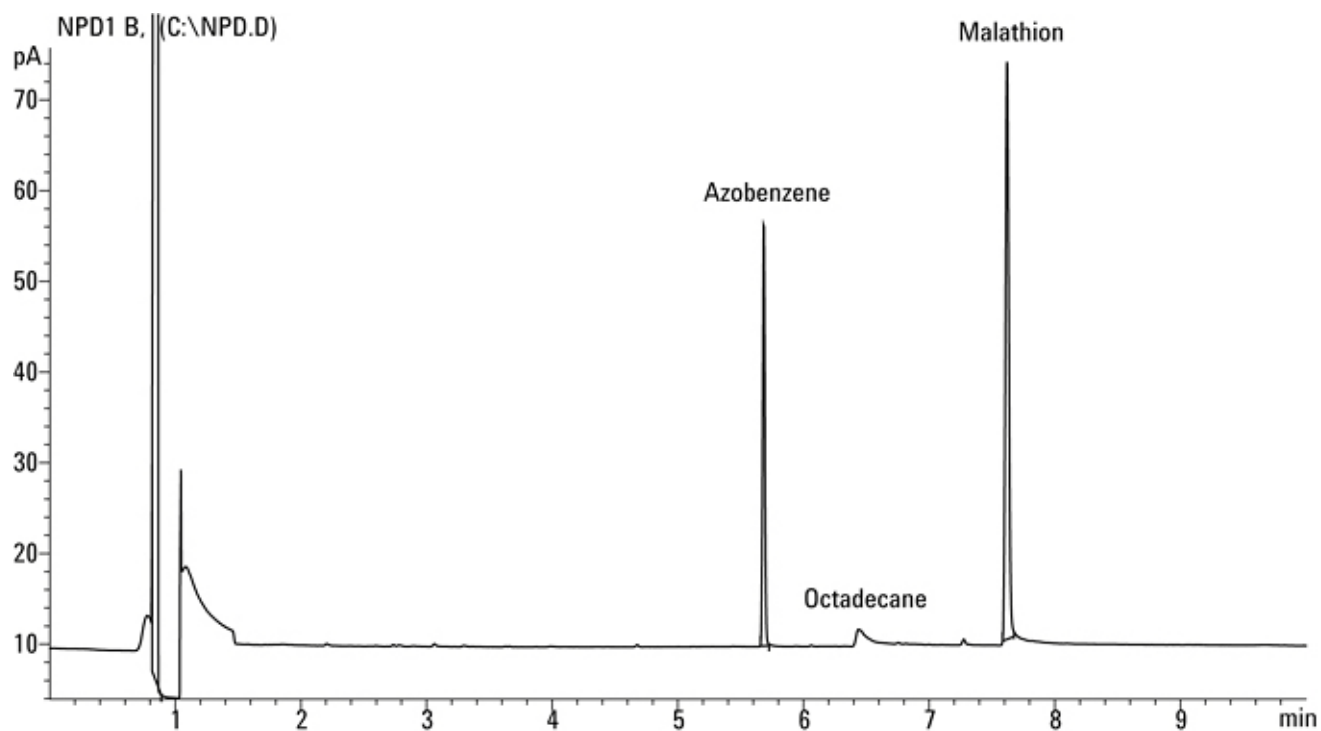

## Typical values

The following are typical values for the detector results. The results for a detector depend on site preparation related factors. These are not pass/fail criteria.

### Agilent ChemStation

		SS/MMI/COC/PTV/PP
Sensitivity	Azobenzene area counts	38
	Malathion area counts	71

### Agilent Chemstation—ASTM noise

Measured with GC isothermal @ 100°C

10 minute blank run, noise range(s) >1 minute

Performance + noise report

NPD ASTM noise ? .0765 pA

### 3396B or 3396C or 3397A integrator sample conditions with analog input [GC Analog output, Range 0]

>>Set DATE and TIME<<

Zero	10	Pk wd	0.04
Att 2^	7	Thresh	7
Cht sp	1	[Stop time 12.0]	
Ar Rej	1000		

		SS/MMI/COC/PTV/PP
Sensitivity	Azobenzene area counts	306,000
	Malathion area counts	575,000

## To Check TCD Performance

- 1 Gather the following:
  - Evaluation column, HP-5 30 m x 0.32 mm x 0.25  $\mu$ m (19091J-413)
  - FID/TCD performance evaluation (checkout) sample (18710-60170)
  - 4-mL solvent and waste bottles or equivalent for autoinjector
  - Chromatographic-grade hexane
  - 2-mL sample vials or equivalent for sample
  - Chromatographic-grade helium as carrier, makeup, and reference gas
  - Inlet and injector hardware
- 2 Verify the following:
  - Chromatographic-grade gases plumbed and configured: helium as carrier gas, nitrogen, hydrogen, and air.
  - Empty waste vials loaded in sample turret.
  - 4-mL solvent vial with diffusion cap filled with isooctane and inserted in Solvent A injector position.
- 3 Replace consumable parts (liner, septum, traps, syringe, and so forth) as needed for the checkout.
- 4 Install the evaluation column. (See the procedure for the SS, MMI, PP, COC, or PTV in the Maintenance manual.)
  - Bake out the evaluation column for at least 30 min at 180 °C. (See the procedure for the SS, MMI, PP, COC, or PTV in the Maintenance manual.)
  - Be sure to configure the column.
- 5 Create or load a method with these parameter values.



**Table 17** TCD Checkout Conditions

<b>Column and sample</b>	
Type	HP-5, 30 m — 0.32 mm — 0.25 µm (19091J-413)
Sample	FID/TCD checkout 18710-60170
Column flow	6.5 mL/min
Column mode	Constant flow
<b>Split/splitless inlet</b>	
Temperature	250 °C
Mode	Splitless
Purge flow	60 mL/min
Purge time	0.75 min
Septum purge	3 mL/min
<b>Multimode inlet</b>	
Mode	Splitless
Inlet temperature	40 °C
Initial time	0.1 min
Rate 1	720 °C/min
Final temp 1	350 °C
Final time 1	2 min
Purge time	1.0 min
Purge flow	40 mL/min
Septum purge	3 mL/min
<b>Purged packed inlet</b>	
Temperature	250 °C
Septum purge	3 mL/min
<b>Cool on-column inlet</b>	
Temperature	Oven track
Septum purge	15 mL/min
<b>PTV inlet</b>	
Mode	Splitless
Inlet temperature	40 °C
Initial time	0.1 min

<b>Column and sample</b>	
Rate 1	720 °C/min
Final temp 1	350 °C
Final time 1	2 min
Rate 2	100 °C/min
Final temp 2	250 °C
Final time 2	0 min
Purge time	0.5 min
Purge flow	40 mL/min
Septum purge	3 mL/min
<b>Detector</b>	
Temperature	300 °C
Reference flow (He)	20 mL/min
Makeup flow (He)	2 mL/min
Baseline output	< 30 display counts on Agilent ChemStation (< 750 µV)
<b>Oven</b>	
Initial temp	40 °C
Initial time	0 min
Rate 1	25 °C/min
Final temp	90 °C
Final time	0 min
Rate 2	15 °C/min
Final temp	170 °C
Final time	2 min
<b>ALS settings (if installed)</b>	
Sample washes	2
Sample pumps	6
Injection volume	1 µL
Syringe size	10 µL
PreInj Solvent A Washes	2
PreInj Solvent B Washes	0


Column and sample	
PostInj Solvent A Washes	2
PostInj Solvent B Washes	0
Viscosity delay	0
Plunger speed	Fast (Slow for cool on-column inlet)
PreInjection dwell	0
PostInjection dwell	0
Manual injection	
Injection volume	1 $\mu$ L
Data system	
Path	Setting
Method and Run Control View: Instrument>Edit Parameters>Options	Keep instrument keyboard locked after method is loaded? = No
Method and Run Control View: Instrument>Edit Parameters>Signals	Assign Signal 1 to the detector Choose Save Data All Data Rate = 5 Hz

- 6 Display the signal output. A stable output at any value between 12.5 and 750  $\mu$ V (inclusive) is acceptable.
  - If the baseline output is < 0.5 display units (< 12.5  $\mu$ V), verify that the detector filament is on. If the offset is still <0.5 display units (< 12.5  $\mu$ V), your detector requires service.
  - If baseline output is > 30 display units (> 750  $\mu$ V), there may be chemical contamination contributing to the signal. Bakeout the TCD. If repeated cleanings do not give an acceptable signal, check gas purity. Use higher purity gases and/or install traps.
- 7 If using a data system, prepare the system to perform one run using the loaded checkout method. Make sure that the data system will output a chromatogram.
- 8 Start the run.

If performing an injection using an autosampler, start the run using the data system or press **[Start]** on the GC.

If performing a manual injection (with or without a data system):

- a** Press [**Prep Run**] to prepare the inlet for splitless injection.
  - b** When the GC becomes ready, inject 1  $\mu\text{L}$  of the checkout sample and press [**Start**] on the GC.
- 9** The following chromatogram shows typical results for a new detector with new consumable parts installed.

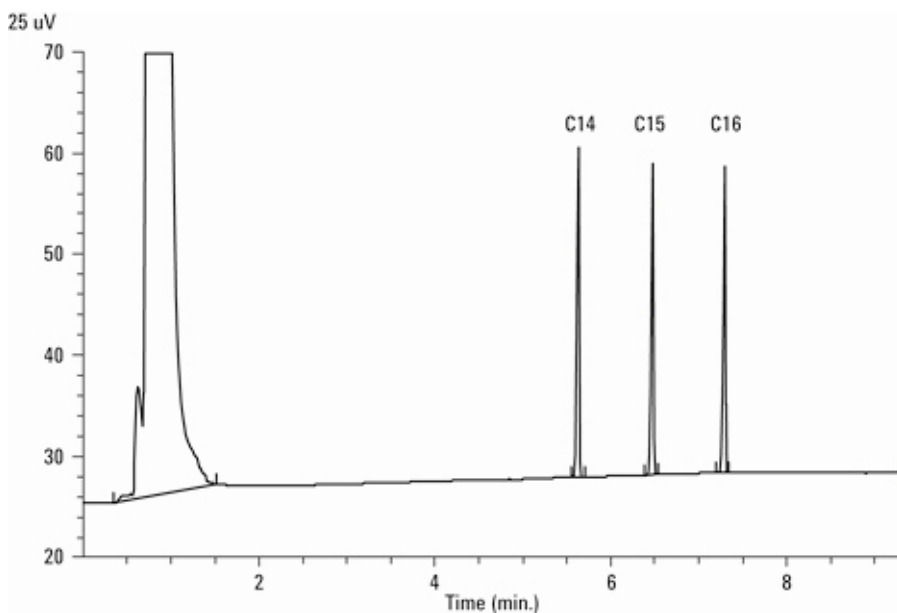

## Typical values

The following are typical values for the detector results. The results for a detector depend on site preparation related factors. These are not pass/fail criteria.

### Agilent ChemStation

		SS/MMI/COC/PTV	PP
Sensitivity	C <sub>14</sub> , C <sub>15</sub> , C <sub>16</sub> area counts	73	65
Discrimination	C <sub>14</sub> /C <sub>16</sub> area ratio	1.00 ± 0.10	1.00 ± 0.10

### Agilent ChemStation—ASTM noise

Measured with GC isothermal at 100 °C

10 minute blank run, noise range(s) > 1 minute

Performance + noise report

TCD ASTM noise -0.05733 display units (25 µV/display unit)

Detector signal set to 5 Hertz

### 3396B or 3396C or 3397A integrator sample conditions with analog input (GC Analog output, Range 0)

>>Set DATE and TIME<<			
Zero	10	Pk wd	0.04
Att^2	7	Thresh	5
Cht sp	1	Stop time	8.75
Ar Rej	1000		

		SS/MMI/COC/PTV	PP
Sensitivity	C <sub>14</sub> , C <sub>15</sub> , C <sub>16</sub> area counts	585,000	515,000
Discrimination	C <sub>14</sub> /C <sub>16</sub> area ratio	1.00±0.10	1.00±0.10

## To Check uECD Performance

- 1 Gather the following:
  - Evaluation column, HP-5 30 m x 0.32 mm x 0.25  $\mu$ m (19091J-413)
  - uECD performance evaluation (checkout) sample (18713-60040, Japan: 5183-0379)
  - 4-mL solvent and waste bottles or equivalent for autoinjector
  - Chromatographic-grade isooctane
  - 2-mL sample vials or equivalent for sample
  - Inlet and injector hardware
- 2 Verify the following:
  - Clean fused silica indented mixing liner installed. If not, install it.
  - Chromatographic-grade gases plumbed and configured: helium.
  - Empty waste vials loaded in sample turret.
  - 4-mL vial with diffusion cap filled with hexane and inserted in Solvent A injector position.
- 3 Replace consumable parts (liner, septum, traps, syringe, and so forth) as needed for the checkout.
- 4 Install the evaluation column. (See the procedure for the SS, MMI, PP, COC, or PTV in the Maintenance manual.)
  - Bake out the evaluation column for at least 30 minutes at 180 °C. (See the procedure for the SS, MMI, PP, COC, or PTV in the Maintenance manual.)
  - Be sure to configure the column.
- 5 Display the signal output to determine baseline output. A stable baseline output at any value between 0.5 and 1000 Hz (ChemStation display units) (inclusive) is acceptable.
  - If the baseline output is < 0.5 Hz, verify that the electrometer is on. If the offset is still < 0.5 Hz, your detector requires service.
  - If the baseline output is > 1000 Hz, there may be chemical contamination contributing to the signal. Bakeout the uECD. If repeated cleanings do not give an acceptable signal, check gas purity. Use higher purity gases and/or install traps.
- 6 Create or load a method with the parameter values listed in the table below.

**Table 18**  $\mu$ ECD checkout conditions

<b>Column and sample</b>	
Type	HP-5, 30 m — 0.32 mm — 0.25 $\mu$ m (19091J-413)
Sample	$\mu$ ECD checkout (18713-60040 or Japan: 5183-0379)
Column mode	Constant flow
Column pressure	6.5 mL/min (helium)
<b>Split/splitless inlet</b>	
Temperature	200 °C
Mode	Splitless
Purge flow	60 mL/min
Purge time	0.75 min
Septum purge	3 mL/min
<b>Multimode inlet</b>	
Mode	Splitless
Inlet temperature	80 °C
Initial time	0.1 min
Rate 1	720 °C/min
Final temp 1	250 °C
Final time 1	5 min
Purge time	1.0 min
Purge flow	60 mL/min
Septum purge	3 mL/min
<b>Purged packed inlet</b>	
Temperature	200 °C
Septum purge	3 mL/min
<b>Cool on-column inlet</b>	
Temperature	Oven track
Septum purge	15 mL/min
<b>PTV inlet</b>	
Mode	Splitless
Inlet temperature	80 °C

<b>Column and sample</b>	
Initial time	0.1 min
Rate 1	720 °C/min
Final temp 1	350 °C
Final time 1	2 min
Rate 2	100 °C/min
Final temp 2	250 °C
Final time 2	0 min
Purge time	0.75 min
Purge flow	60 mL/min
Septum purge	3 mL/min
<b>Detector</b>	
Temperature	300 °C
Makeup flow (N2)	30 mL/min (constant + makeup)
Baseline output	Should be < 1000 display counts in Agilent ChemStation (< 1000 Hz)
<b>Oven</b>	
Initial temp	80 °C
Initial time	0 min
Rate 1	15 °C/min
Final temp	180 °C
Final time	10 min
<b>ALS settings (if installed)</b>	
Sample washes	2
Sample pumps	6
Injection volume	1 µL
Syringe size	10 µL
PreInj Solvent A Washes	2
PreInj Solvent B Washes	0
PostInj Solvent A Washes	2
PostInj Solvent B Washes	0
Viscosity delay	0
Plunger speed	Fast (Slow for cool on-column inlet)



Column and sample	
PreInjection dwell	0
PostInjection dwell	0
Manual injection	
Injection volume	1 µL
Data system	
Path	Setting
Method and Run Control View: Instrument>Edit Parameters>Options	Keep instrument keyboard locked after method is loaded? = No
Method and Run Control View: Instrument>Edit Parameters>Signals	Assign Signal 1 to the detector Choose Save Data All Data Rate = 5 Hz


- 7 If using a data system, prepare the system to perform one run using the loaded checkout method. Make sure that the data system will output a chromatogram.

- 8 Start the run.

If performing an injection using an autosampler, start the run using the data system or press **[Start]** on the GC.

If performing a manual injection (with or without a data system):

- a Press **[Prep Run]** to prepare the inlet for splitless injection.
  - b When the GC becomes ready, inject 1 µL of the checkout sample and press **[Start]** on the GC.
- 9 The following chromatogram shows typical results for a new detector with new consumable parts installed.

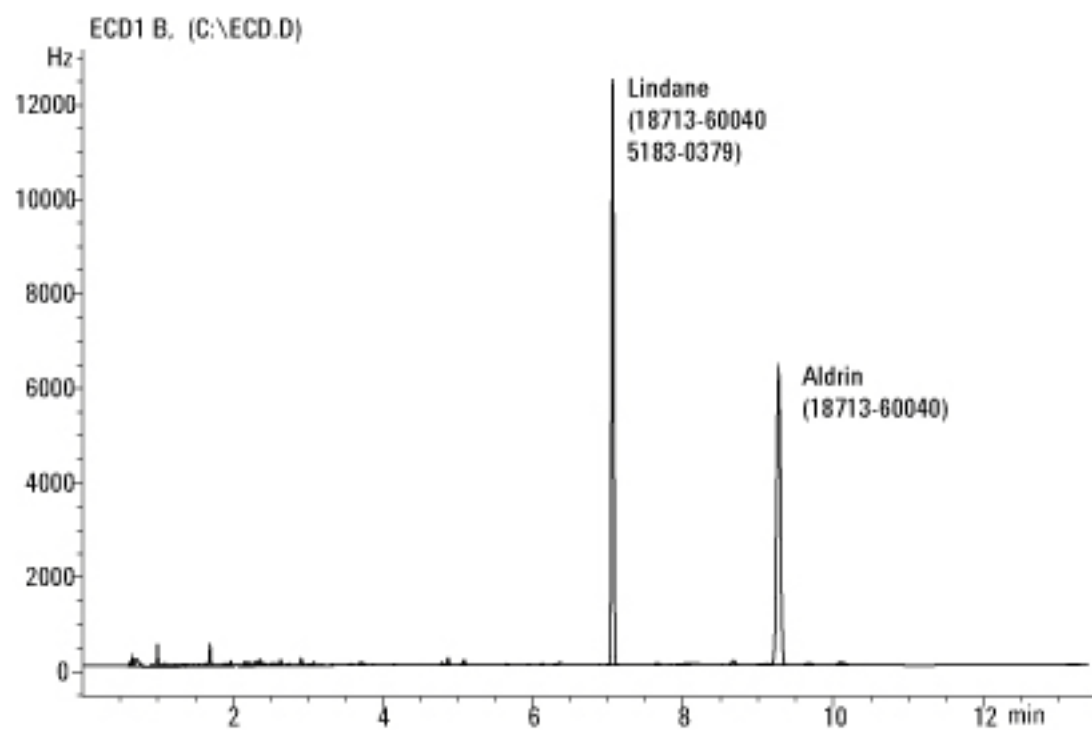

## Typical values

The following are typical values for the detector results. The results for a detector depend on site preparation related factors. These are not pass/fail criteria.

### Agilent ChemStation Noise ([N3])

Measured by ChemStation, ASTM [N3]	< 3 Hz
------------------------------------	--------

### Signal

Area, Lindane peak	7500 × [N3]
--------------------	-------------

**3396B or 3396C or 3397A integrator sample conditions with analog input [GC Analog output, Range 0]**

>>Set DATE and TIME<<

### GC Analog output

### Noise ([N1])

Measured at [GC Range 0],

[3396 Attn. 3, Chart speed 1, zero 50]	<37 mm for 1 min. measurement
--	-------------------------------

### Signal

GC Range 6

### 3396

Zero	10	Pk wd	0.04
Att 2^	5	Thresh	5
Cht sp	1	Stop time	11.0
Ar Rej	100000		

Area, Lindane peak	47,950 × [N1]
--------------------	---------------

## To Check FPD Performance (Sample 5188-5953)

If sample FPD checkout sample 5188-5953 is available, use it to test FPD performance.

To check FPD performance, first check the phosphorus performance, then the sulfur performance.

### Preparation

- 1 Gather the following:
  - Evaluation column, HP-5 30 m x 0.32 mm x 0.25  $\mu$ m (19091J-413)
  - FPD performance evaluation (checkout) sample (5188-5953)
  - Phosphorus filter
  - Sulfur filter and filter spacer
  - 4-mL solvent and waste bottles or equivalent for autoinjector
  - 2-mL sample vials or equivalent for sample
  - Chromatographic-grade isooctane for syringe wash solvent
  - Inlet and injector hardware
- 2 Verify the following:
  - Capillary column adapter installed. If not, install it.
  - Chromatographic-grade gases plumbed and configured: helium as carrier gas, nitrogen, hydrogen, and air.
  - Empty waste vials loaded in sample turret.
  - 4-mL vial with diffusion cap filled with isooctane and inserted in Solvent A injector position.
- 3 Verify that the Lit Offset is set appropriately. Typically, it should be about 2.0 pA for the checkout method.
- 4 Replace consumable parts (liner, septum, traps, syringe, and so forth) as needed for the checkout.
- 5 Install the evaluation column. (See the procedure for the SS, MMI, PP, COC, or PTV in the Maintenance manual.)
  - Set the oven, inlet, and detector to 250 °C and bake out for at 15 min. (See the procedure for the SS, MMI, PP, COC, or PTV in the Maintenance manual.)

Be sure to configure the column.

## Phosphorus performance

- 1 If it is not already installed, install the phosphorus filter.
- 2 Create or load a method with the parameter values listed in the table below.

**Table 19** FPD Checkout Conditions (P)

Column and sample	
Type	HP-5, 30 m x 0.32 mm x 0.25 µm (19091J-413)
Sample	FPD checkout (5188-5953)
Column mode	Constant pressure
Column pressure	25 psi
Split/splitless inlet	
Temperature	200 °C Split/splitless
Mode	Splitless
Purge flow	60 mL/min
Purge time	0.75 min
Septum purge	3 mL/min
Multimode	
Mode	Splitless
Inlet temperature	75 °C
Initial time	0.1 min
Rate 1	720 °C/min
Final temp 1	250 °C
Final time 1	5.0 min
Purge time	1.0 min
Purge flow	60 mL/min
Septum purge	3 mL/min
Purged packed inlet	
Temperature	200 °C
Septum purge	3 mL/min
Cool on-column inlet	
Temperature	Oven track (cool on-column)

Column and sample	
Septum purge	15 mL/min
PTV inlet	
Mode	Splitless
Inlet temperature	75 °C
Initial time	0.1 min
Rate 1	720 °C/min
Final temp 1	350 °C
Final time 1	2 min
Rate 2	100 °C/min
Final temp 2	250 °C
Final time 2	0 min
Purge time	0.75 min
Purge flow	60 mL/min
Septum purge	3 mL/min
Detector	
Temperature	200 °C (On)
Hydrogen flow	75 mL/min (On)
Air (Oxidizer) flow	100 mL/min (On)
Mode	Constant makeup flow OFF
Makeup flow	60 mL/min (On)
Makeup gas type	Nitrogen
Flame	On
Lit offset	Typically 2 pA
Electrometer	On
Oven	
Initial temp	70 °C
Initial time	0 min
Rate 1	25 °C/min
Final temp 1	150 °C
Final time 1	0 min
Rate 2	5 °C/min

Column and sample	
Final temp 2	190 °C
Final time 2	4 min
ALS settings (if installed)	
Sample washes	2
Sample pumps	6
Injection volume	1 µL
Syringe size	10 µL
PreInj Solvent A Washes	2
PreInj Solvent B Washes	0
PostInj Solvent A Washes	2
PostInj Solvent B Washes	0
Viscosity delay	0
Plunger speed	Fast (Slow for cool on-column inlet)
PreInjection dwell	0
PostInjection dwell	0
Manual injection	
Injection volume	1 µL

Path	Setting
Method and Run Control View: Instrument>Edit Parameters>Options	Keep instrument keyboard locked after method is loaded? = No
Method and Run Control View: Instrument>Edit Parameters>Signals	Assign Signal 1 to the detector Choose Save Data All Data Rate = 5 Hz

- 3 Ignite the FPD flame, if not lit.

- 4** Display the signal output and monitor. This output typically runs between 40 and 55 but can be as high as 70. Wait for the output to level off. This takes approximately 1 hour.

If the baseline output is too high:

- Check column installation. If installed high, the stationary phase burns out and increases measured output.
- Check for leaks.
- Bake out the detector and column at 250 °C.
- Wrong flows set for installed filter.

If the baseline output is zero, verify the electrometer is on and the flame is lit.


- 5** If using a data system, prepare the system to perform one run using the loaded checkout method. Make sure that the data system will output a chromatogram.
- 6** Start the run.

If performing an injection using an autosampler, start the run using the data system or press **[Start]** on the GC.

If performing a manual injection (with or without a data system):

- a** Press **[Run]** to prepare the inlet for splitless injection.
  - b** When the GC Prep becomes ready, inject 1 µL of the checkout sample and press **[Start]** on the GC.
- 7** The following chromatogram shows typical results for a new detector with new consumable parts installed.



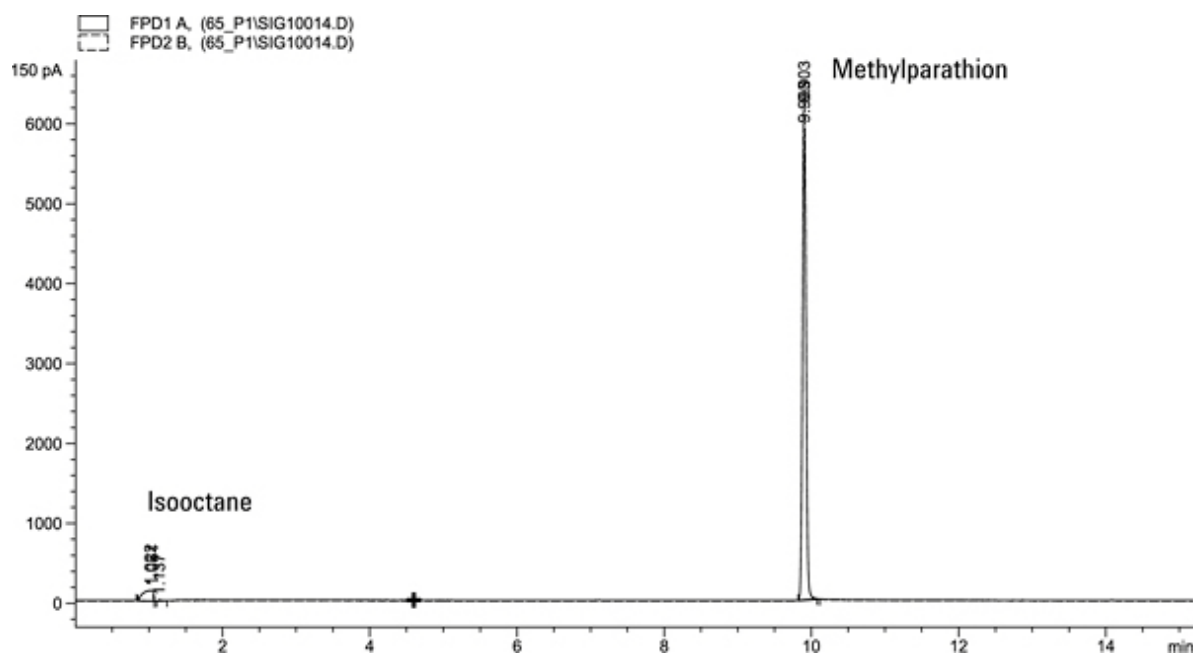

### Typical values

The following are typical values for the detector results. The results for a detector depend on site preparation related factors. These are not pass/fail criteria.

FPD P filter	Typical range after 24 hours	Limits at installation
MDL (pg/sec)	0.06 to 0.08	$\leq 0.10$
Peak area	19000 to 32000	$\geq 19000$
Signal height	5000 to 11000	--
Noise	1.6 to 3.0	$\leq 4$
Half-width (min)	0.05 to 0.07	--
Output	34 to 80	$\leq 80$

## Sulfur performance

- 1 Install the sulfur filter and filter spacer.
- 2 Make the following method parameter changes.

Sulfur method parameters (S)

Parameter	Value ( mL/min)
H2 flow	50
Air flow	60

- 3 Display the signal output and monitor. This output typically runs between 50 and 60 but can be as high as 70. Wait for the output to level off. This takes approximately 1 hour.

If the baseline output is too high:

- Check column installation. If installed high, the stationary phase burns out and increases measured output.
- Check for leaks.
- Bake out the detector and column at 250 °C.
- Wrong flows set for installed filter.


If the baseline output is zero, verify the electrometer is on and the flame is lit.

- 4 If using a data system, prepare the system to perform one run using the loaded checkout method. Make sure that the data system will output a chromatogram.
- 5 Start the run.

If performing an injection using an autosampler, start the run using the data system or press **[Start]** on the GC.

If performing a manual injection (with or without a data system):

- a Press **[Prep Run]** to prepare the inlet for splitless injection.
  - b When the GC becomes ready, inject 1 µL of the checkout sample and press **[Start]** on the GC.
- 6 The following chromatogram shows typical results for a new detector with new consumable parts installed.

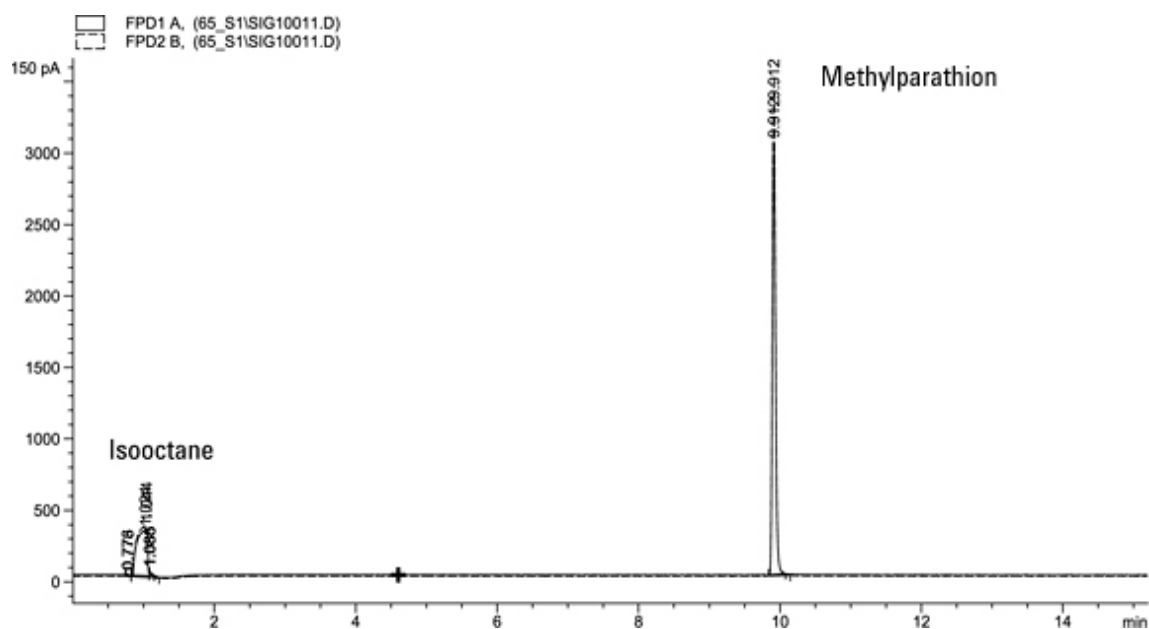

### Typical values

The following are typical values for the detector results. The results for a detector depend on site preparation related factors. These are not pass/fail criteria.

FPD S filter	Typical range after 24 hours	Limits at installation
MDL (pg/sec)	3.8 to 5	$\leq 6$
Peak area	8000 to 19000	$\geq 8000$
Signal height	2500 to 6000	--
Noise	2 to 4	$\leq 5$
Half-width (min)	0.06 to 0.08	--
Output	34 to 65	$\leq 70$

## To Check FPD Performance (Sample 5188-5245, Japan)

To verify FPD performance, first check the phosphorus performance, then the sulfur performance.

### Preparation

- 1 Gather the following:
  - Evaluation column, DB5 15 m x 0.32 mm x 1.0  $\mu$ m (123-5513)
  - FPD performance evaluation (checkout) sample (5188-5245)
  - Phosphorus filter
  - Sulfur filter and filter spacer
  - 4-mL solvent and waste bottles or equivalent for autoinjector
  - 2-mL sample vials or equivalent for sample.
  - Chromatographic-grade isooctane for syringe wash solvent.
  - Inlet and injector hardware
- 2 Verify the following:
  - Capillary column adapter installed. If not, install it.
  - Chromatographic-grade gases plumbed and configured: helium as carrier gas, nitrogen, hydrogen, and air.
  - Empty waste vials loaded in sample turret.
  - 4-mL vial with diffusion cap filled with isooctane and inserted in Solvent A injector position.
- 3 Verify the lit offset is set appropriately. Typically, it should be about 2.0 pA for the checkout method.
- 4 Replace consumable parts (liner, septum, traps, syringe, and so forth) as needed for the checkout.
- 5 Install the evaluation column. (See the procedure for the SS, MMI, PP, COC, or PTV in the Maintenance manual.)
  - Set the oven, inlet, and detector to 250 °C and bake out for at 15 min. (See the procedure for the SS, MMI, PP, COC, or PTV in the Maintenance manual.)
  - Configure the column.

## Phosphorus performance

- 1 If it is not already installed, install the phosphorus filter.
- 2 Create or load a method with these parameter values.

**Table 20** FPD Phosphorus Checkout Conditions

Column and sample	
Type	DB-5MS, 15 m x 0.320 µm x 1.0 µm (123-5513)
Sample	FPD checkout (5188-5245)
Column mode	Constant flow
Column pressure	7.5 mL/min
<b>Split/splitless inlet</b>	
Temperature	250 °C Split/splitless
Mode	Splitless
Total purge flow	69.5 mL/min
Purge flow	60 mL/min
Purge time	0.75 min
Septum purge	3 mL/min
<b>Multimode</b>	
Mode	Splitless
Inlet temperature	75 °C
Initial time	0.1 min
Rate 1	720 °C/min
Final temp 1	250 °C
Final time 1	5.0 min
Purge time	1.0 min
Purge flow	60 mL/min
Septum purge	3 mL/min
<b>Purged packed inlet</b>	
Temperature	250 °C
Septum purge	3 mL/min

<b>Column and sample</b>	
<b>Cool on-column inlet</b>	
Temperature	Oven track
Septum purge	15 mL/min
<b>PTV inlet</b>	
Mode	Splitless
Inlet temperature	80 °C
Initial time	0.1 min
Rate 1	720 °C/min
Final temp 1	350 °C
Final time 1	2 min
Rate 2	100 °C/min
Final temp 2	250 °C
Final time 2	0 min
Purge time	0.75 min
Purge flow	60 mL/min
Septum purge	3 mL/min
<b>Detector</b>	
Temperature	200 °C (On)
Hydrogen flow	75.0 mL/min (On)
Air (oxidizer) flow	100.0 mL/min (On)
Mode	Constant makeup flow Off
Makeup flow	60.0 mL/min (On)
Makeup gas type	Nitrogen
Flame	On
Lit offset	Typically 2 pA
Electrometer	On
<b>Oven</b>	
Initial temp	70 °C
Initial time	0 min
Rate 1	10 °C/min
Final temp	105 °C

Column and sample	
Final time	0 min
Rate 2	20 °C/min
Final temp 2	190 °C
Final time 2	7.25 min for sulfur 12.25 min for phosphorus
ALS settings (if installed)	
Sample washes	2
Sample pumps	6
Injection volume	1 µL
Syringe size	10 µL
PreInj Solvent A Washes	2
PreInj Solvent B Washes	0
PostInj Solvent A Washes	2
PostInj Solvent B Washes	0
Viscosity delay	0
Plunger speed	Fast (Slow for cool on-column inlet)
PreInjection dwell	0
PostInjection dwell	0
Manual injection	
Injection volume	1 µL
Data System	

Path	Setting
Method and Run Control View: Instrument>Edit Parameters>Options	Keep instrument keyboard locked after method is loaded? = No
Method and Run Control View: Instrument>Edit Parameters>Signals	Assign Signal 1 to the detector Choose Save Data All Data Rate = 5 Hz

- 3 Ignite the FPD flame, if not lit.

- 4 Display the signal output and monitor. This output typically runs between 40 and 55 but can be as high as 70. Wait for the output to level off. This takes approximately 1 hour.

If the baseline output is too high:

- Check column installation. If installed high, the stationary phase burns out and increases measured output.
- Check for leaks.
- Bake out the detector and column at 250 °C.
- Wrong flows set for installed filter.
- If the baseline output is zero, verify the electrometer is on and the flame is lit.

- 5 If using a data system, prepare the system to perform one run using the loaded checkout method. Make sure that the data system will output a chromatogram.

If performing an injection using an autosampler, start the run using the data system or press [Start] on the GC.

If performing a manual injection (with or without a data system):

- a Press [**Prep Run**] to prepare the inlet for splitless injection.
- b When the GC becomes ready, inject 1 µL of the checkout sample and press [Start] on the GC.

- 6 Start the run.


If performing an injection using an autosampler, start the run using the data system or press [**Start**] on the GC.

If performing a manual injection (with or without a data system):

- a Press [**Prep Run**] to prepare the inlet for splitless injection.
- b When the GC becomes ready, inject 1 µL of the checkout sample and press [**Start**] on the GC.

- 7 The following chromatogram shows typical results for a new detector with new consumable parts installed.



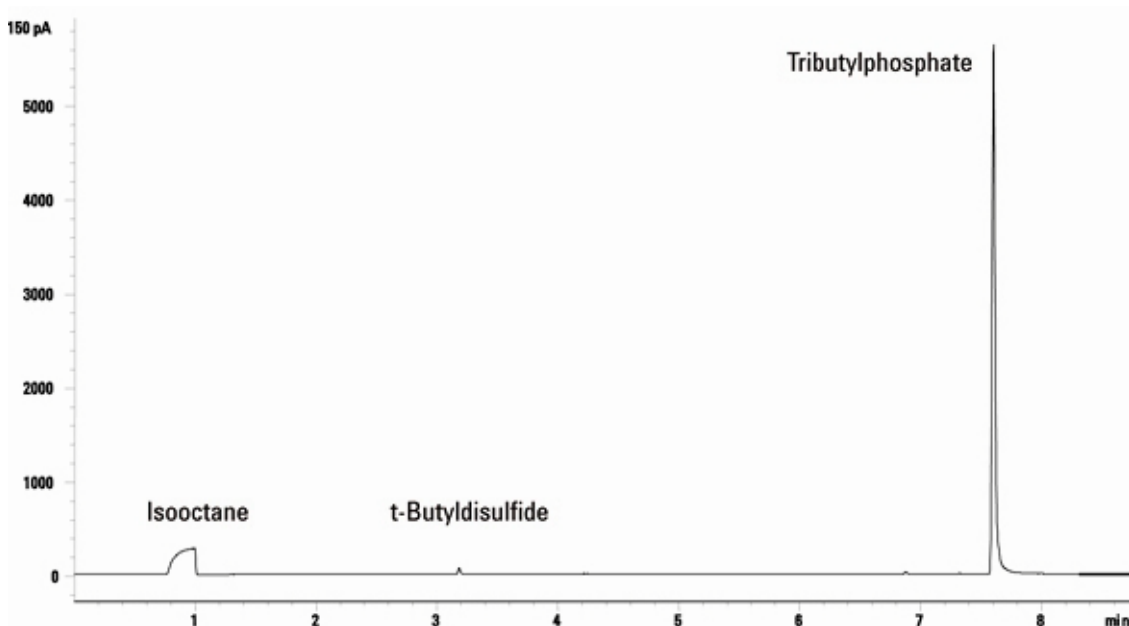

## Typical values

The following are typical values for the detector results. The results for a detector depend on site preparation related factors. These are not pass/fail criteria.

FPD P filter	Typical range after 24 hours	Limits at installation
MDL (pg/sec)	0.06 to 0.08	$\leq 0.10$
Peak area	19000 to 32000	$\geq 19000$
Signal height	5000 to 11000	--
Noise	1.6 to 3.0	$\leq 4$
Half-width (min)	0.05 to 0.07	--
Output	34 to 80	$\leq 80$

## Sulfur performance

- 1 Install the sulfur filter.
- 2 Make the following method parameter changes.

Sulfur method parameters

Parameter	Value ( mL/min)
H2 flow	50
Air flow	60

- 3 Ignite the FPD flame, if not lit.
- 4 Display the signal output and monitor. This output typically runs between 50 and 60 but can be as high as 70. Wait for the output to level off. This takes approximately 2 hours.

If the baseline output is too high:

- Check column installation. If installed high, the stationary phase burns out and increases measured output.
- Check for leaks.
- Bake out the detector and column at 250 °C.
- Wrong flows set for installed filter.

If the baseline output is zero, verify the electrometer is on and the flame is lit.


- 5 If using a data system, prepare the data system to perform one run using the loaded checkout method. Make sure the data system will output a chromatogram.
- 6 Start the run.

If performing an injection using an autosampler, start the run using the system or press [**Start**] on the GC.

If performing a manual injection (with or without a data system):

- Press [**Prep Run**] to prepare the inlet for splitless injection.
- When the GC becomes ready, inject 1 µL of the checkout sample and press [**Start**] on the GC.


- 7 The following chromatogram shows typical results for a new detector with new consumable parts installed.

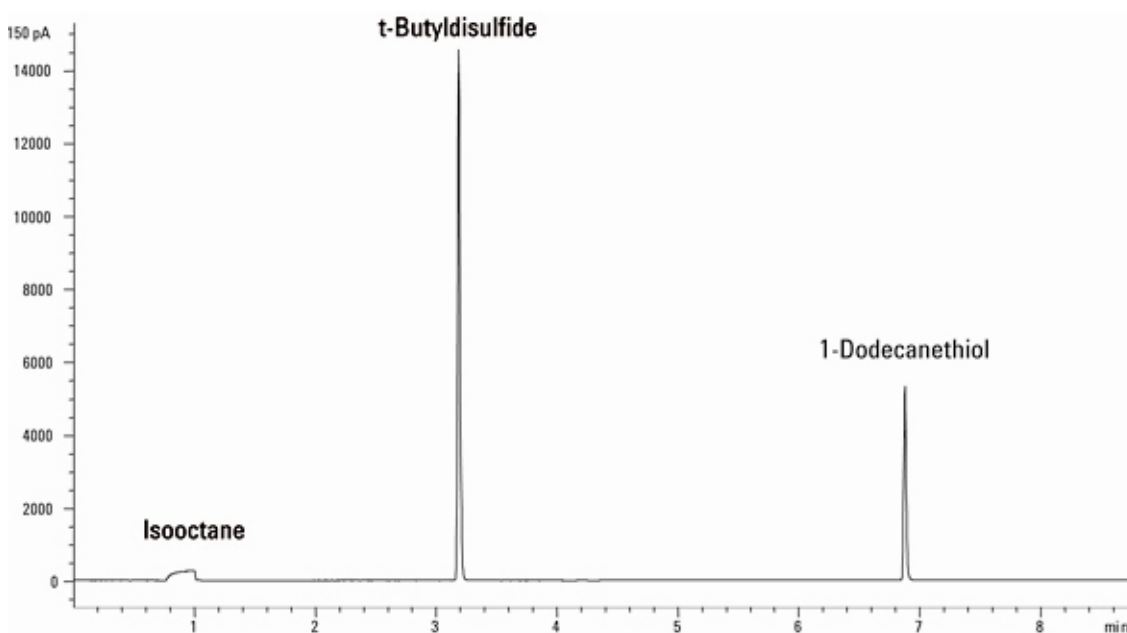

### Typical values

The following are typical values for the detector results. The results for a detector depend on site preparation related factors. These are not pass/fail criteria.

FPD S filter	Typical range after 24 hours	Limits at installation
MDL (pg/s)	3.8 to 5	$\leq 6$
Peak area	8000 to 19000	$\geq 8000$
Signal height	2500 to 6000	--
Noise	2 to 4	$\leq 5$
Half-width (min)	0.06 to 0.08	--
Output	34 to 65	$\leq 70$

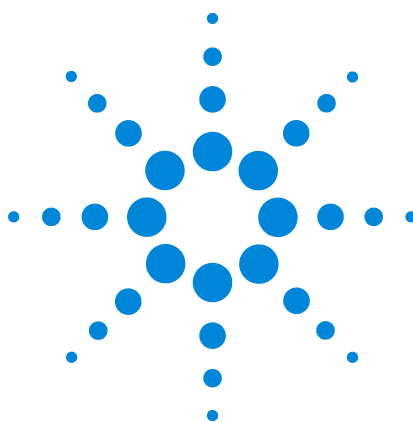

## 6 Configuration

Configuration overview	273
GC modules and the communications buss	280
To configure an MSD transfer line	283
To configure a valve box	286
Replacing a GC module	288
Removing a GC module	289
Changing the GC configuration	290
Configuration locks	291
Installing new devices	292
Configuring time	



## Configuration overview

There are 2 kinds of configuration information: hard and soft.

- Hard configuration includes what hardware options are installed on the GC: the type of inlets and detectors, where they are located, and what auxiliary heaters are connected.
- Soft configuration includes application specific information: carrier gas type, make up gas type, inlet control modes, column dimensions.

The GC modules that make up an instrument (i.e., EPC modules, detector boards, and heaters) are initially configured during manufacturing. This information is stored in the signal board memory.

- An inlet is configured as an EPC module and a heater.
- A detector is configured as an EPC module, a detector signal board, and a heater.
- A pneumatic control module (PCM) is configured as an EPC module and can be associated with a heater.
- An Aux EPC module (AUX) is configured as 3 channels of auxiliary pressure control and can be associated with a heater.
- A valve box is configured with one or two heaters.
- An AUX thermal zone is configured with a heater. For example, an MSD transfer line or a nickel-catalyst accessory.

There are many custom configuration options available. Most are explained in G3430-90024, Firmware A.01.05 Instructions For Channel Partner Devices.

You can configure the GC modules by pressing [Config] and then pressing the key for the GC module, for example, [Config] [Front Det]. Under the configuration menu, you will find these features:

- Install or remove the module.
- If you wish to ignore a ready or not ready status.
- The type of gas that the EPC module is controlling.
- The baseline output to use for reigniting the flame.

## Inlet example

The configuration control table for a front SSL Inlet includes 5 items.

CONFIGURE FRONT INLET	
Configured:	Split/Splitless
Ignore Ready =	FALSE
Makeup gas type	He
[ EPC1 ] = (INLET)	(SS)
FINLET (OK) 68 watts	250.0

This list interprets each line of the control table.

- **Configured:** The type of module configured as a front inlet is an split/splitless capillary inlet.
- **Ignore Ready:** If TRUE, the GC will ignore a Not Ready status. If the inlet actuals do not match their setpoints, the GC will go ready. The ignored items are displayed under Status.
- This line indicates the EPC module. In this case, it is plugged into the 4-wire connector labeled EPC1 and is an inlet EPC module for an split/splitless inlet.
- This line describes the heater and sensor. In this case, there is a 64 watt heater cartridge associated with the front inlet and the sensor is reading 250.0°C.

## Detector example

The configuration control table for a front FID detector includes 7 items.

CONFIGURE FRONT DETECTOR	
Configured:	(FID)
Ignore Ready =	FALSE
Makeup gas type	N2
Lit offset	2.0
[ F-DET ] = (SIGNAL)	(FID)
[ EPC3 ] = (DET_EPC)	(FID)
F-DET (OK) 64 watts	300.0

This list interprets each line of the control table.

- Configured: The type of module configured as a front detector is an FID.
- Ignore Ready: If TRUE, the GC will ignore a Not Ready status. If the detector actuals do not match their setpoints, the GC will go ready. The ignored items are displayed under Status.
- Lit offset: The detector will try to reignite if the output drops to this level.
- This line indicates the signal board. In this case, it is plugged into the 4-wire communication buss connector labelled F-DET and is an FID signal board.
- This line indicates the EPC module. In this case, it is plugged into the 4-wire connector labeled EPC3 and is a detector EPC module for an FID.
- This line describes the heater and sensor. In this case, there is a 64 watt heater cartridge associated with the front detector and the sensor is reading 300.0°C.

## PCM example

The configuration control table for an pneumatic control Module (PCM) module includes 6 items. For custom hardware, the configuration could include a heater.

```
CONFIGURE PCM C
Configured:                PCM
Ignore Ready =             FALSE
Gas type                   N2
Aux Gas type               N2
Aux Mode                   Forward Pressure
[ EPC6 ] = (AUX_EPC)       (PCM)
```

This list interprets each line of the control table.

- Configured: The type of module configured as a pneumatic control module or PCM.
- Ignore Ready: If TRUE, the GC will ignore a Not Ready status. If the PCM actuals do not match their setpoints, the GC will go ready. The ignored items are displayed under Status.
- The next two lines indicate what gas type is associated with each channel. The Gas type refers to channel one. The options for an both channels are: Helium, Hydrogen, Nitrogen, Argon methane 5%, and Air.
- This line indicates the Aux Mode for the pressure control channel. This can be Forward Pressure or Back Pressure control.
- This line indicates the EPC module. In this case, it is plugged into the 4-wire connector labelled EPC6 and is an AUX EPC module for a PCM.
- There are no heaters configured in this example. However, you can associate a PCM with a front or back inlet heater.

```
CONFIGURE PCM C
Install EPC6               Htr = NONE
Install EPC6               Htr = FINLET
Install EPC6               Htr = BINLET
```



## AUX example

The configuration control table for an Aux EPC module includes 6 items. For custom hardware, the configuration could include a heater.


CONFIGURE APC 1,2,3	
Configured:	Aux Pressure
Ignore Ready =	FALSE
Chan 1 Gas type	N2
Chan 2 Gas type	N2
Chan 3 Gas type	N2
[ EPC6 ] = (AUX_EPC)	(AUX)

This list interprets each line of the control table.

- Configured: The type of module configured as Aux EPC channels 1, 2 and 3 is Aux Pressure module.
- Ignore Ready: If TRUE, the GC will ignore a Not Ready status. If the AUX actuals do not match their setpoints, the GC will go ready. The ignored items are displayed under Status.
- The next three lines indicate what gas type is associated with each channel. The options for an AUX channel are: Helium, Hydrogen, Nitrogen, Argon methane 5%, and Air
- This line indicates the EPC module. In this case, it is plugged into the 4-wire connector labeled EPC6 and is an AUX EPC module labeled AUX.
- There are no heaters associated with an AUX module in this position.

## Heater assignments

There are 4 standard connection for heaters and sensors. With a valve bracket and harness, you add 2 more for a total of 6. The following picture and table describe the harness labels and typical uses.

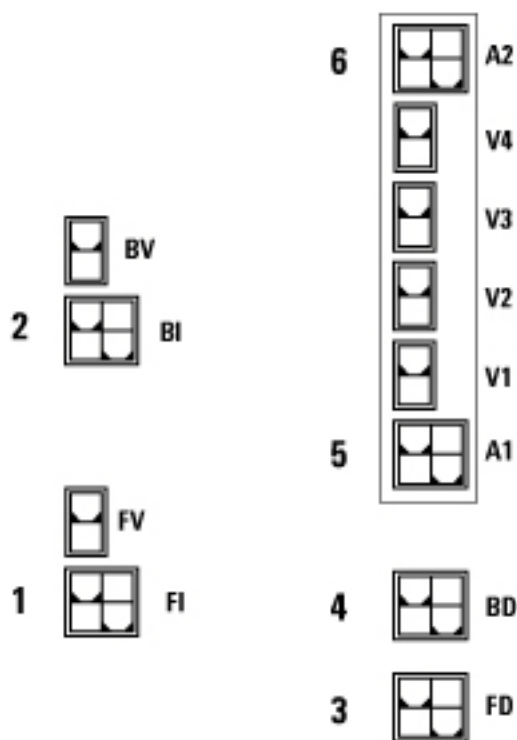

Number	Label	Typical uses
1	FI	Front inlet
2	BI	Back inlet
3	FD	Front detector
4	BD	Back detector
5	A1	Valve box, 3rd detector, Front FPD
6	A2	MSD transfer line, Valve box, Back FPD

For our specials group and for our channel partners, there are many options for configuring the heaters. These options are based on connecting a heater and configuring it to a type of module to enable a certain type of behavior. See the

configuration topics in Firmware Instructions for Channel Partner Devices.

- As an Inlet or detector, a small zone is isothermal.
- As an Thermal Aux, a small zone can be programmed.

This table describes this flexibility available in firmware.


Controlling	Number	Heater type	Configuration	Setpoints
Front inlet	1	Inlet	[Config][Front Inlet]	[Front Inlet]
Back inlet	2	Inlet Nickel catalyst	[Config][Back Inlet] [Config][Thermal Aux 1]	[Back Inlet] [Aux Temp]
Front detector	3	Detector	[Config][Front Det]	[Front Det]
Back detector	4	Detector Custom heater	[Config][Back Det] [Config][Thermal Aux 2]	[Back Det] [Aux Temp]
Thermal Aux 1	5	3rd Detector Front FPD	[Config][Aux Det] [Config][Back Det]	[Aux Det] [Back Det]
Thermal Aux 2	6	MSD transfer line Back FPD	[Config][Aux Temp] [Config][Back Det]	[Aux Temp] [Back Det]
Valve box	5 6 5 and 6	Valve box Valve box Valve box	[Config][Thermal Aux 2] [Config][Back Det] [Config][Valve Box]	[Aux Temp] [Back Det] [Valve Box]

Notes:


- 1 You can only configure 6 heaters. Each heater can only be configured to one component type.
- 2 Front FPD uses heaters 3 and 5.
- 3 Back FPD uses heaters 4 and 6.
- 4 FPDs can be configured for one or two heater versions.

## GC modules and the communications buss


This GC provides a flexible architecture. Along with the flexibility, the configuration of the GC can be more complex. The flexible architecture is represented in this schematic. Above each arrow is a label that represents an address on a communications buss.

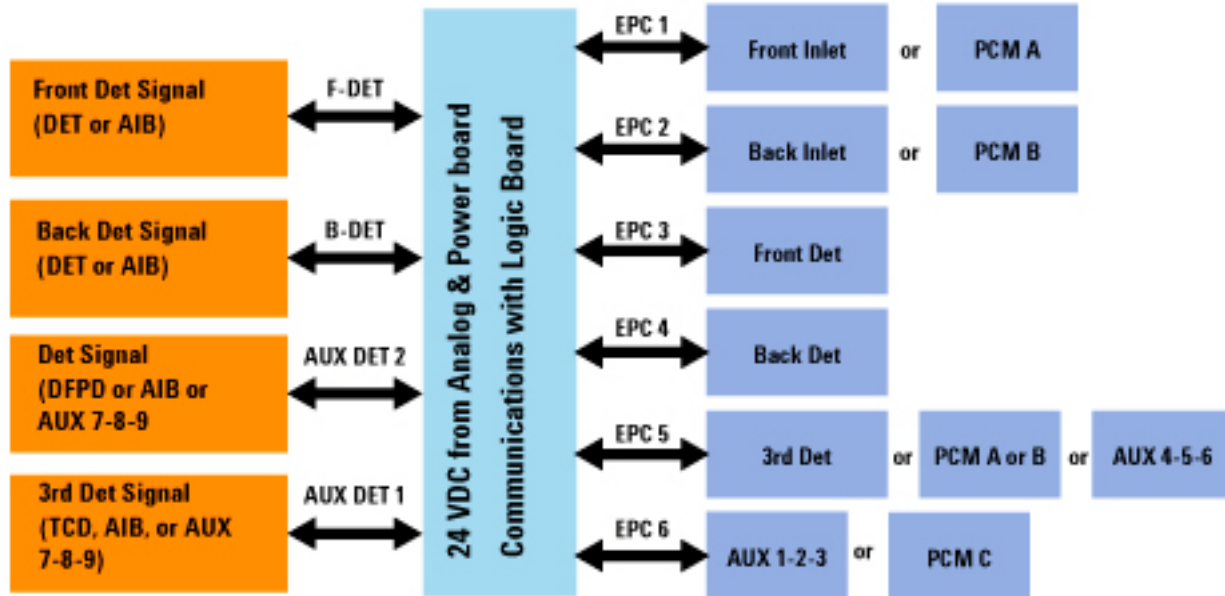

The 4-wire communications buss is the nervous system of the GC. It connects to up to 6 EPC modules and up to 4 detector signal boards. On each of the 10 connectors, there is a label to help identify what can be plugged in.

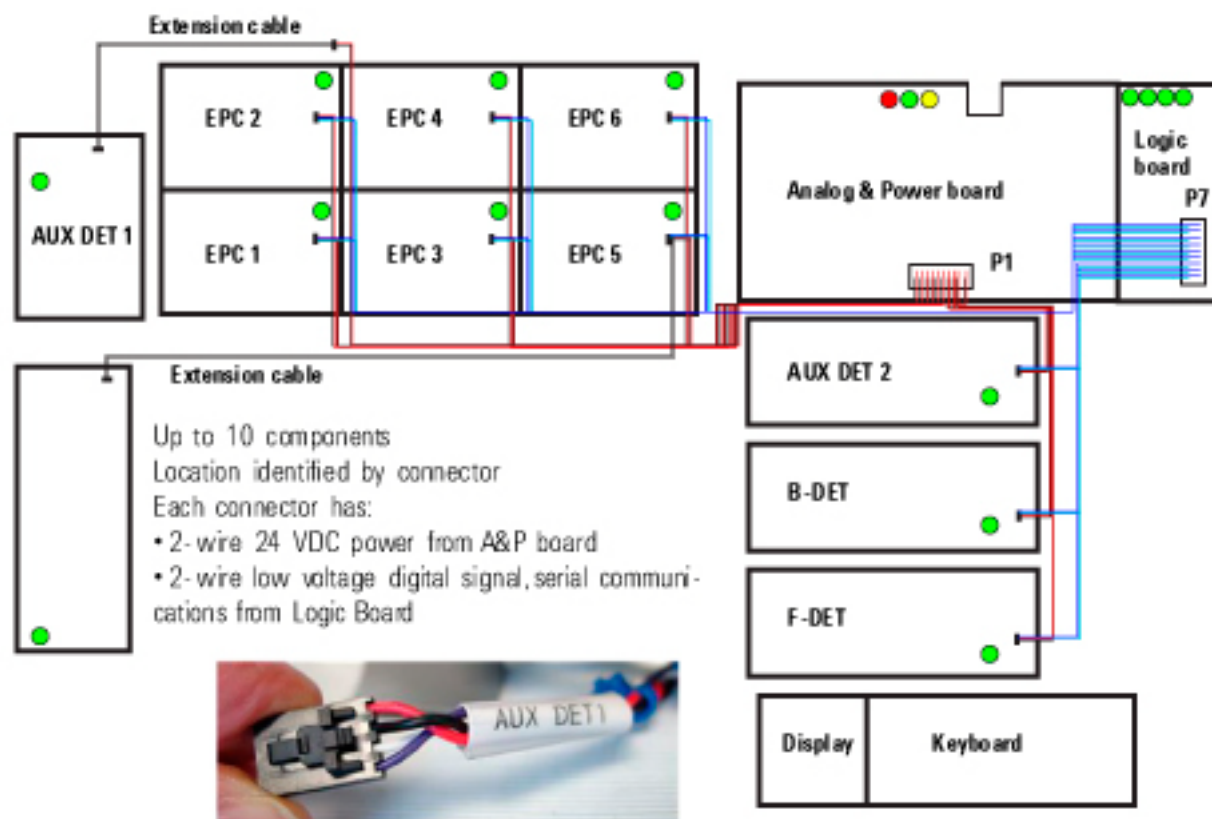


The red and black wires on the buss provide power to the module. The black wire is ground. The red wire is switched, unregulated +24 VDC. All of these wires plug into the Analog and Power Board. The white and colored wires provide serial communications between the signal board and the module.

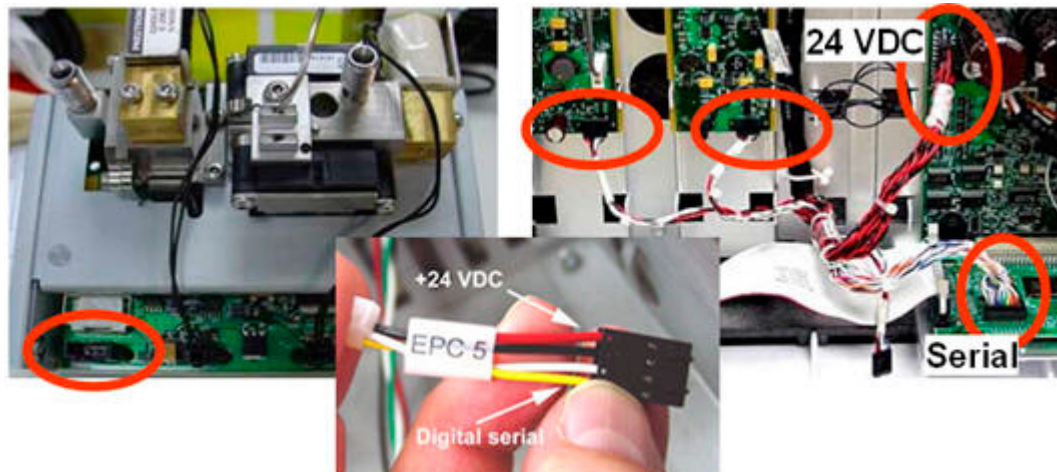

There are 2 extension cables for the GC modules of the 3<sup>rd</sup> detector:


- For the GC modules, there is an extension cable for both the EPC module and the signal board; G3432-60505.
- For the heater, there is an extension heater cable to connect to the AUX 2 heater; G3432-60506.

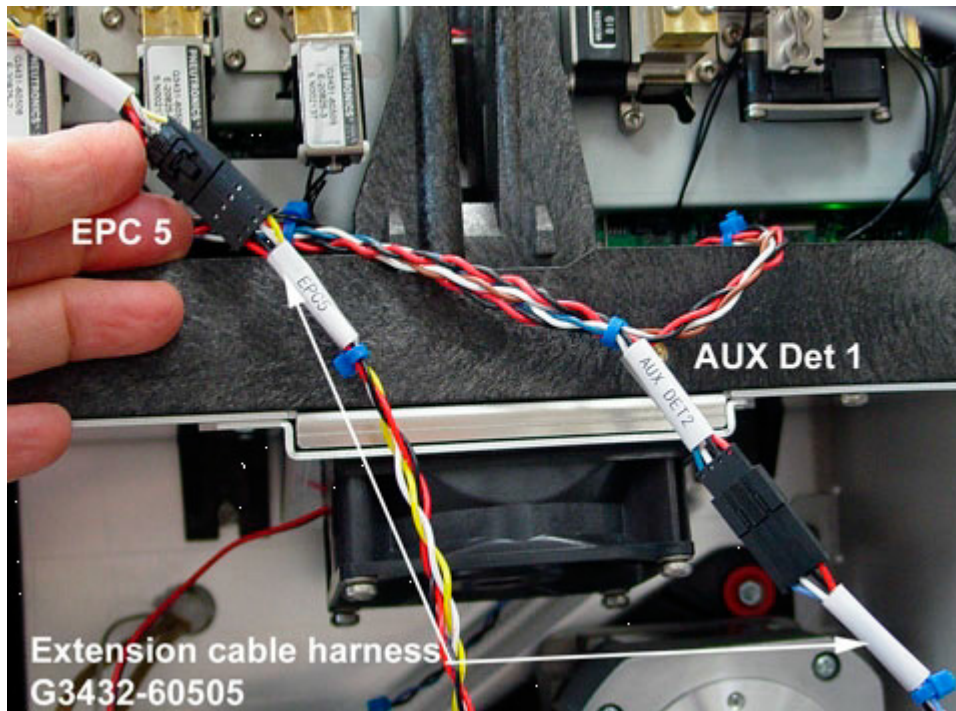


## To configure an MSD transfer line

The factory does NOT configure any MSD transfer line. Thus you may find, on installation, that the MSD heater will not heat, that you cannot set a temperature, and that the transfer line is not seen by the software. The transfer line heater may be configured as a valve box heater or other device.

These instructions explain how to go from the configuration shipped to you which may look like this (the valve box case)

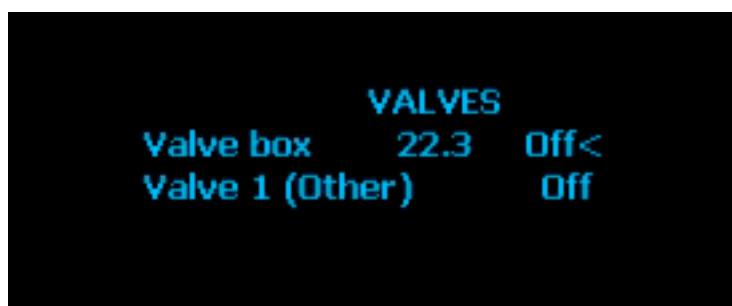


to the configuration that you need for the transfer line.

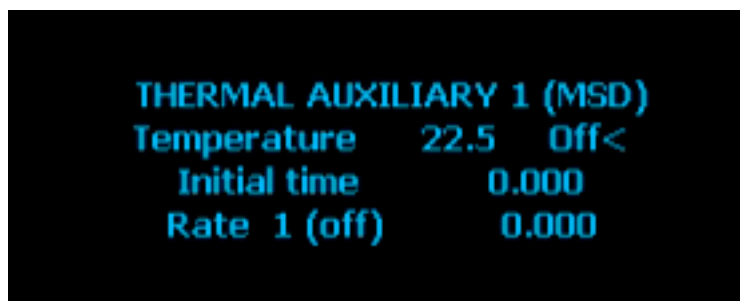

To set up the MSD transfer line requires the user to un-configure where the aux heater has been configured.


- 1 Press Config Aux Temp #. Scroll to the Valve Box and press Enter. CONFIGURE VALVE BOX will be displayed.

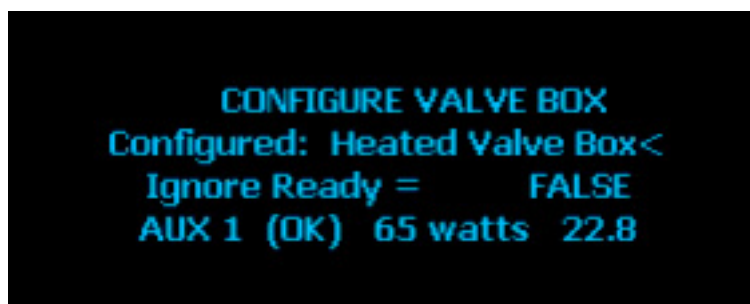


- 2 Press Mode/Type. This will display the Remove module message.



- 3 Press Enter. This displays the following message.

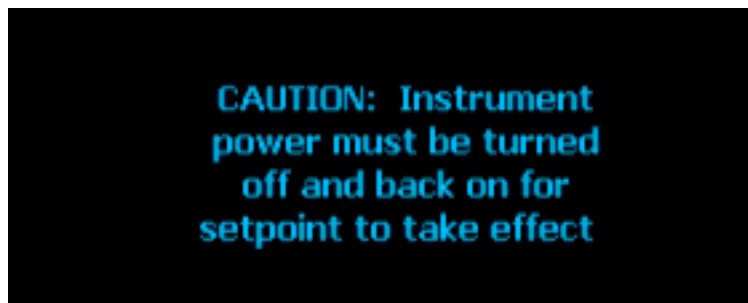

- 4 While a power cycle is required for the change to take effect, it does not need to be done until the MSD transfer line configuration is complete. Press Clear to return to CONFIGURE VALVE BOX. Return to Configure Aux Temp and select the Thermal Aux zone corresponding to the MSD transfer line connection.
  - Thermal Aux 1 if the heater cable is connected to connector 5.
  - Thermal Aux 2 if the heater cable is connected to connector 6.
- 5 Using Thermal Aux 1 as the example, CONFIGURE AUX HEATER 1 is displayed after selecting Thermal Aux 1 and pressing Enter.

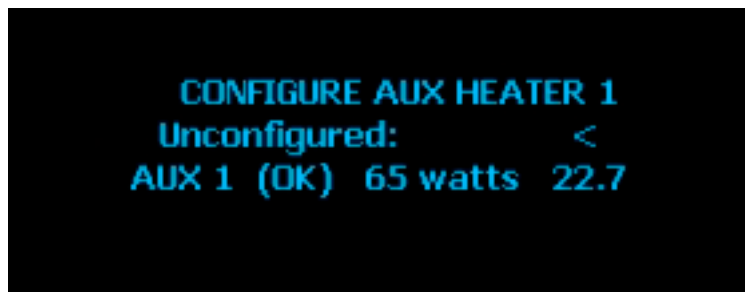




- 6 Press Mode/Type. The following is displayed.



```

CONFIGURE AUX HEATER 1
Install Heater    AUX 1 <
  
```


- 7 Press Enter. The Caution message appears again. Press Clear. This displays the following.

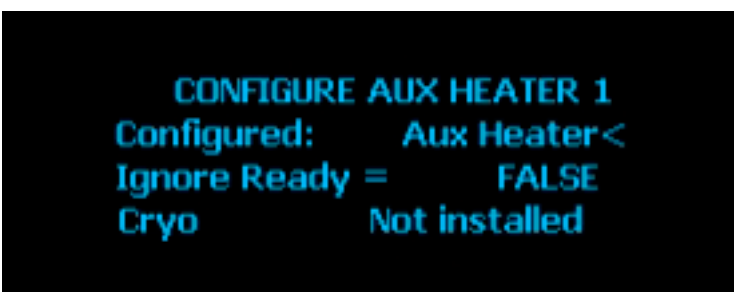

```

CONFIGURE AUX HEATER 1
Configured:      Aux Heater<
Ignore Ready =   FALSE
Cryo             Not installed
  
```


- 8 Scroll to Auxiliary type and press Mode/ Type. Scroll to MSD transfer line.

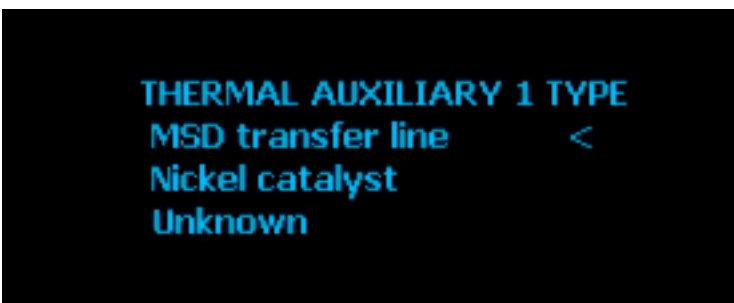

```

THERMAL AUXILIARY 1 TYPE
MSD transfer line    <
Nickel catalyst
Unknown
  
```

- 9 Press Enter to select. This shows the Caution message again. Now power cycle the instrument.
- 10 Press Thermal Aux # and the appropriate number if necessary. This shows the thermal aux that was configured.

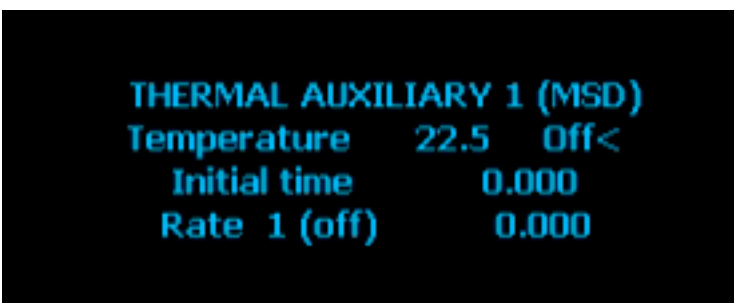

```

THERMAL AUXILIARY 1 (MSD)
Temperature  22.5  Off<
Initial time    0.000
Rate 1 (off)    0.000
  
```

## To configure a valve box

The common approach is to configure the valve box heater to Valve Box. Alternates are to use Thermal Aux 2 or Thermal Aux 1. A valve box can house one or two heaters.

### MSD transfer line heater error

If you are updating a GC with an MSD transfer line heater, you may observe a heater fault or error associated with Thermal Aux 2.

Objective: Configure the Thermal Aux 2 heater, plugged into the Aux 2 position, as an MSD transfer line heater. This may be difficult because the heater may be configured as something else. You may have to unconfigure it first.

Assumptions: The heater cable is plugged into the 5975A MSD and into the Aux 2 position of the GC. Firmware is version A.01.04 or higher.

### Initial steps

- 1 Press **Options**. Scroll to Keyboard & Display and press **Enter**.
- 2 Turn Hard configuration lock off.
- 3 Press **Config**. Scroll to Thermal Aux 2 and press **Enter**.
- 4 You should see Configured Aux Heater on line 2 and Auxiliary type: MSD transfer line on line 5.
- 5 If line 2 is correct, scroll to line 5. If the type is wrong, press **Mode/Type** and select MSD transfer line. Confirm the change by pressing **Aux Temp**.
- 6 If line 2 indicates that nothing is configured, press **Mode/Type** to see if anything is available for configuration. If nothing is, go to the next steps. If Thermal Aux 2 is available, press **Enter** to select it, then power cycle the GC.

## If that doesn't work

Examine the other heaters that could be configured as Thermal Aux 2. These are Valve Box or Thermal Aux 1.

- 1 Press **Config**. Scroll to Valve Box and examine line 2.
- 2 If it is configured, press **Mode/Type** and then **Enter** to Remove the module.
- 3 Navigate back to Thermal Aux 2 to see if you can configure it as Aux 2.
- 4 Continue in this way until you find where Aux 2 is configured and unconfigure it.

## Replacing a GC module

When replacing an EPC module or a detector signal board, you will have to update the configuration of the GC. This is because the configuration information includes the unique ID# and manufacturing date for these parts.

To configure the new GC module

- 1 Replace the module and power on the GC. The GC will list the Exceptions to the configuration and indicate that these modules are broken or missing.
- 2 Unlock the “[Configuration locks](#)”.
- 3 Press [**Config**] and select the GC module that you replaced.
- 4 Press [**Mode/Type**]. The GC displays **Remove module**. Press [**Enter**].
- 5 A Caution message tells you to reboot the GC.
- 6 Press [**Mode/Type**]. The GC displays **Install module**. Press [**Enter**].

### CAUTION

Be sure that you select the correct module configuration. Incorrect selections can cause operation problems. For example, for an FPD select the 2 heater version.

- 7 A Caution message tells you to reboot the GC.
- 8 Press [**Options**]. Scroll to Communications and press [**Enter**].
- 9 Scroll to Reboot the GC. Press [**Yes/On**] twice.
- 10 Navigate to the configuration table for the GC module.
- 11 Verify that the configuration is correct. Modify any additional configuration items, such as, Gas type or Lit offset.

## Removing a GC module

It is easy to unconfigure a module and make it appear to be uninstalled without disconnecting any cables. This is useful if a module fails and the customer still wants to use the other modules. For example, if the back inlet fails, the customer can unconfigure the back inlet and continue to use the front inlet and detectors.

- 1 Press **[Options]** and select Keyboard & Display.
- 2 Press **[Off / No]** to turn the Configuration Lock off.
- 3 Press **[Config]** and the key for the GC module to remove.
- 4 From the line labeled Configure, press **[Mode/Type]**.
- 5 Press **[Enter]**. The GC will display a message requesting you to power cycle the instrument.
- 6 Power cycle the instrument.
- 7 Verify the removal by pressing the key for the GC module. The GC should display an uninstalled message.

## Changing the GC configuration

The GC configuration is changed when you do any of the following:

- Add or replace a detector signal board.
- Add or replace an EPC module.
- Replace the GC logic board.
- Add or remove a tray or auto-injector.
- Add or remove a heater.
- Change any soft configuration items, such as, gas types, heater types, and Lit offset.

Changing the GC configuration impacts data system methods and methods stored in the GC.

### Agilent data system methods

During the first connection and when loading any of the methods created with the old configuration, the data system will note the change in configuration and force you to resolve the differences. If you do not resolve differences in both the configuration and method parameter screens the method will not load. Run method or run sequence will not start the run.

- Use the GC Configuration Connection screen to get the new configuration. Review and update any configuration setpoints.
- Use the Edit parameters screen to review and update any method setpoints.

### Methods stored in GC

Methods stored in the GC that do not match the new configuration must be updated. Make the necessary additions or modifications. Store this method over the old method.

When you load a method, the GC compares the method against the current hardware it detects (inlet types, detector types, PCM, heaters, valves, and so forth). A method which does not match is not loaded. Default setpoints for this configuration are loaded. Before adding new hardware or changing the installed hardware type, verify that any methods stored in the GC have been documented.

## Configuration locks

When trying to change the GC configuration—as by adding, removing, or changing an Inlet module—you may encounter this or a similar message:

### **Configuration is Locked**

#### **Go to keyboard options to unlock**

To unlock the configuration

- 1** Press [**Options**].
- 2** Select **Keyboard and Display**.
- 3** Select **Hard Configuration lock**, then press [**No/Off**].

When you power cycle the GC, the configuration lock is automatically set to On.

When you use the Reboot GC feature under Options | Communications, the configuration lock is not changed.

## Installing new devices

When a new device is added to the GC, it must be described so that the GC knows what it is and how to use it.

This procedure, using a PCM module as an example, illustrates the process.

- 1 Replace all instrument covers in the reverse order in which they were removed.
- 2 Plug in the GC and turn on the power.
- 3 Unlock the GC keypad.
  - a Press [**Options**].
  - b Select **Keyboard & Display**.
  - c Select Hard Configuration lock, then press [**No/Off**].
- 4 Press [**Configure**].
- 5 Scroll to select the entry for the new PCM module (**PCM A**, **PCM B**, or **PCM C**) and press [**Enter**].
- 6 On the GC keypad, press [**Mode/Type**].
- 7 The GC presents available configuration options.
  - If using the PCM for something other than an inlet (for example, to provide carrier gas to a splitter), scroll to **Install EPC# Htr = None** (where # indicated the slot chosen).
  - If using the PCM for inlet carrier gas supply, refer to the inlet manufacturer's documentation for the correct configuration choice.
- 8 Press [**Enter**]. A caution message appears instructing you to reboot.
- 9 Reboot the GC.
  - a Press [**Options**].
  - b Scroll to **Communications** and press [**Enter**].
  - c Scroll to **Reboot the GC?** and press [**Yes/On**] twice to reboot the GC and have the changes take effect. This also resets the Hard Configuration lock.
- 10 Zero the pressure sensors.




- 11 Configure the gas types and control mode.
  - a Press [**Config**], then scroll to the correct **PCM** entry (**PCM A**, **PCM B**, or **PCM C**).
  - b Scroll to **Gas type** and use the [**Mode/Type**] key to set the gas type connected to the PCM's Channel **1**.
  - c Scroll to **Aux Gas type** and use the [**Mode/Type**] key to set the gas type connected to the PCM's Channel **2**.
  - d Scroll to **Aux Mode** and use the [**Mode/Type**] key to set the pressure control mode (forward or backward) for the PCM's Channel **2**.

## Configuring time

To localize the date/time stamp in the Run and Event logs, use this procedure.

- 1** Press [**Config**][**Time**].
- 2** Enter the offset from Greenwich Mean Time using 24-hour format. Press [**Enter**].
- 3** Enter the local time (24 hour format). Press [**Enter**].
- 4** Enter the date, in ddmmyy format. Press [**Enter**].

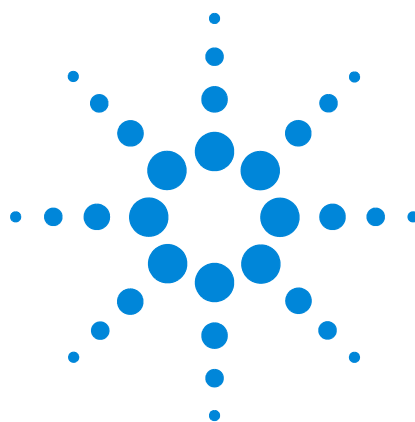

## 7 External Cabling

Connectors on the back of the GC 296

Remote Start/Stop 297

External event 305


Automatic sampler for GC 306

Analog signal outputs 310

Analog signal inputs 312



## Connectors on the back of the GC

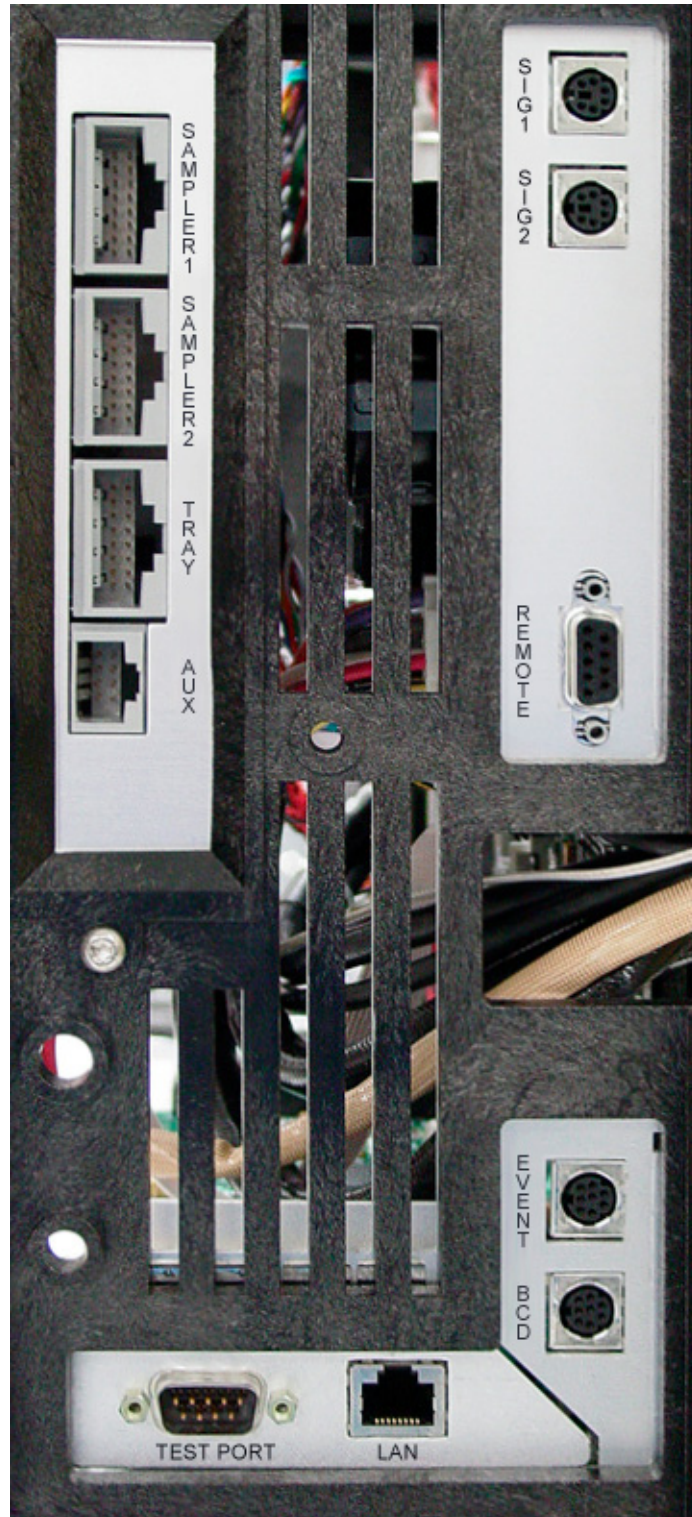

## Remote Start/Stop

These cables are used to link multiple instruments together, so that they start and stop together.

### APG Remote Control

Remote control allows easy connection between single instruments or systems to ensure coordinated analysis with simple coupling requirements. For example, you might have an integrator, automatic sampler, and a gas chromatograph connected with Remote cables.

You can synchronize a maximum of ten instruments using Remote cables. Control of analysis is maintained by instrument readiness - READY for next analysis, followed by START of run and optional STOP of run triggered on the respective lines. In addition, PREPARE may be issued to initiate pre-run activities. All devices connected to APG remote are connected in parallel. For example, if one device is NOT READY then all devices on the bus are also NOT READY.

### Electrical specifications

The APG signals are a modified open collector type. The signal levels are generally TTL levels (low voltage is logic zero, high voltage is logic one) but the open circuit voltage will be between 2.5 to 3.7 Volts. The typical voltage is 3 Volts. A voltage over 2.2 volts will be interpreted as a high logic state while a voltage below 0.4 volts will be interpreted as a low logic state. These levels provide some margin over the specifications of the devices used.

The pull-up resistance, connected to the open-circuit voltage, is in the range of about 1K ohms to 1.5K ohms. For a logic-low state, for a single device on the bus, the minimum current you must be able to sink is 3.3 milliamps. Since devices are connected in parallel, when you have multiple devices this minimum current must be multiplied by the number of devices attached on the bus. The maximum voltage for a low-input state = 0.4V.

The bus is passively pulled high. Leakage current out of a port must be less than 0.2 milliamps to keep the voltage from being pulled lower than 2.2 volts. Higher leakage current may cause the state to be interpreted as a low.

Over-voltage protection - APG Remote connections are clamped by a zener diode to 5.6 Volts. Exceeding this voltage will damage the circuit (main board).

Suggested Drive Circuits

A signal on the APG bus may be driven by another APG device or by one of the following circuits.


A relay, with one side connected to ground, when closed will set a logic-low state.

An NPN transistor, with the emitter connected to ground and the collector connected to the signal line will set a logic-low state if proper base current is supplied.

An open-collector logic gate will perform this same function.

A low-side drive IC will also work, but Darlington-type drivers should be avoided as they will not meet the low-side voltage requirement of less than 0.4V.

The Remote Connector

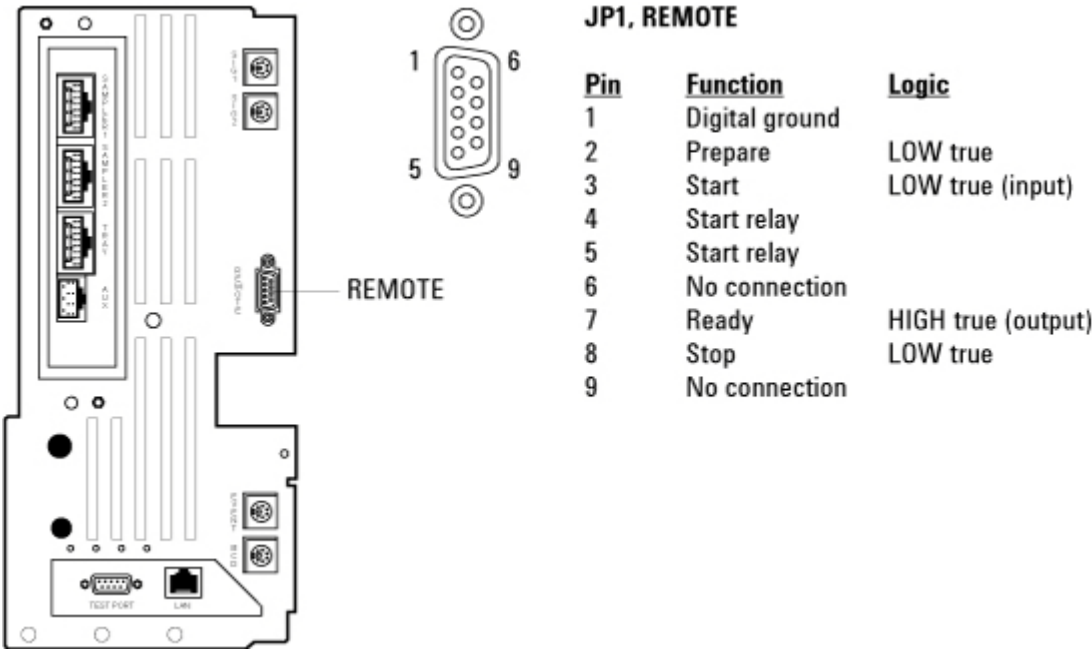

## Signal Descriptions


**Prepare (Low True)**—Request to prepare for analysis. Receiver is any module performing pre-analysis activities. For example, shorting pin 2 to ground will put the GC into "Prep-Run" state. This is useful for Splitless Mode to prepare the inlet for injection or when using "Gas Saver Mode." This function is not needed by Agilent autosampler systems.

**Ready (High True)**— If The Ready Line is high ( $> 2.2$  VDC) then the system is ready for next analysis. Receiver is any sequence controller.

**Start (Low True)**—Request to start run/timetable. Receiver is any module performing runtime-controlled activities.

**Stop (Low True)**—Request to reach system ready state as soon as possible (for example, stop run, abort or finish, and stop injection). Receiver is any module performing runtime-controlled activities. Normally this line is not connected, if the GC oven program is used to control the method "Stop" time.


## Timing Diagram

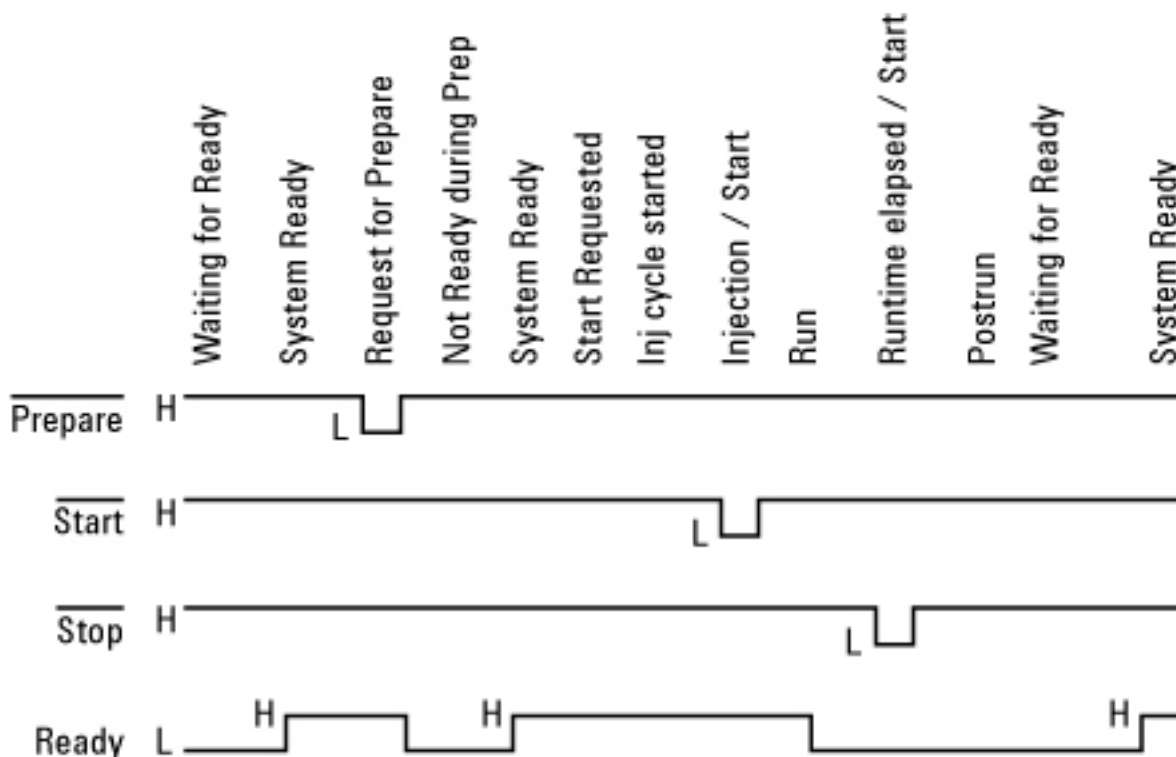

## Cable pinouts, remote start/stop, general use

Part number 35900-60670

Part no. 35900-60670

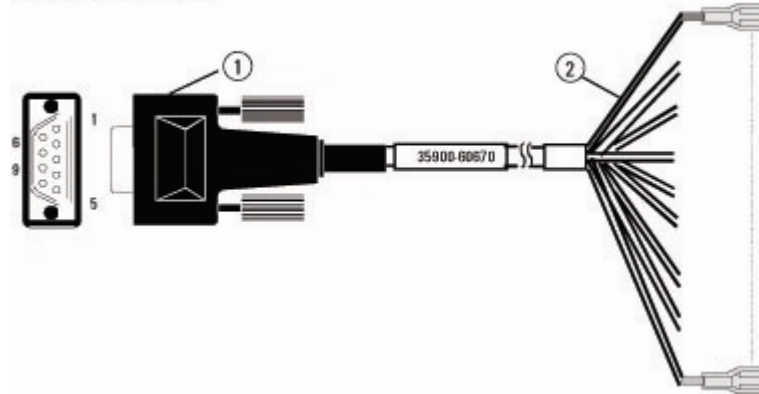


Connector 1 9 pin (male)	Signal name	Connector 2 (spade lugs)
1	GND	Black
2	Prepare (low true)	White
3	Start (low true)	Red
4	Start relay	Green
5	Start relay	Brown
6	No connection	Blue
7	Ready (high true input)	Orange
8	Stop (low true)	Yellow
9	No connection	Violet



## Cable pinouts, GC to 3395B/3396C Integrator

Part number 03396-61010

Part no. 03396-61010

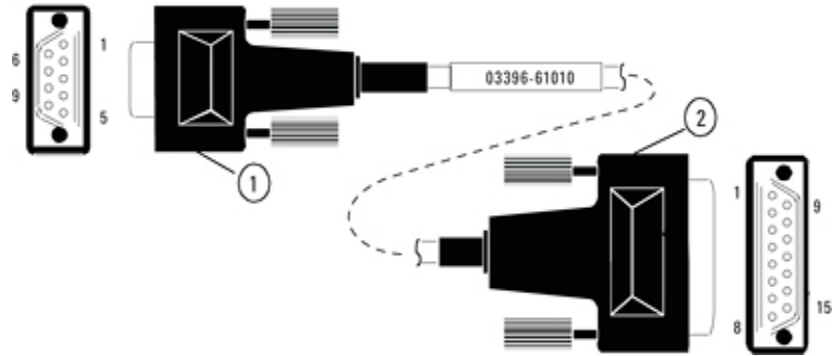


Connector 1 9 pin (male)	Signal name	Connector 2 15 pin (male)
1	GND	9 — Ground
2	Prepare	No connection
3	Start	3 — Start in
4	Start relay	No connection
5	Start relay	No connection
6	No connection	No connection
7	Ready	14 — Ready out
8	Stop	4 — STOP2 In
9	No connection	No connection

## Cable pinouts, GC to 35900C, D, E/MSD/Sampler, 2 meters

Part number G1530-60930

Part no. G1530-60930

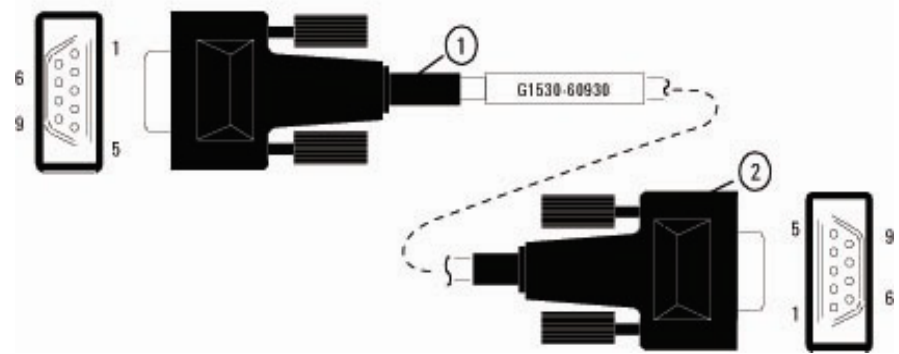


Connector 1 (male)	Signal name	Connector 2 (male)
1	GND	1
2	Prepare	2
3	Start	3
4	Start relay	4
5	Start relay	5
6	No connection	6
7	Ready	7
8	Stop	8
9	No connection	9

## Cable pinouts, GC to Headspace Sampler

Part number G1290-60575

Part no. G1290-60575

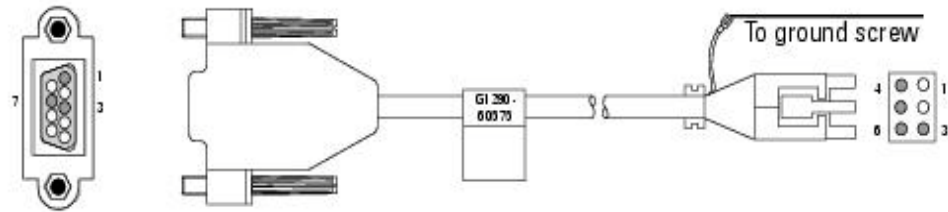


Connector 1 3-pin male	Signal name	Connector 4-pin male
1	GND	4, 5
3	Start	6
7	Ready	3

## Cable pinouts, APG Remote Y-cable

Part number 31530-61020

Part no. G1530-61200

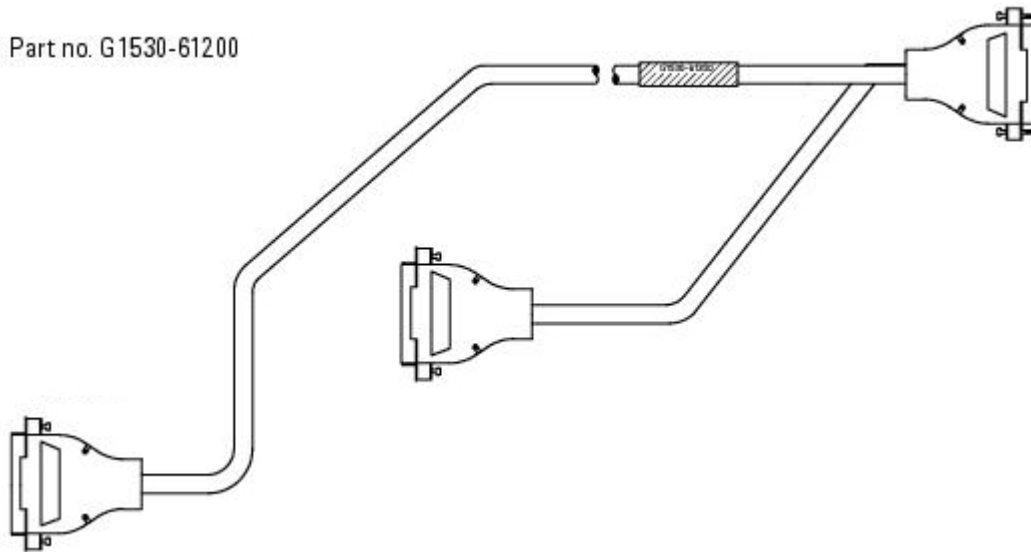

All connectors 9 pin (male)	Signal name
1	GND
2	Prepare
3	Start
4	Start relay
5	Start relay
6	No connection
7	Ready
8	Stop
9	No connection


## External event

Two passive relay contact closures and two 24-volt control outputs are available for controlling external devices. Devices connected to the passive contact closures must be connected to their own power source.

### CAUTION

The external event connector is similar to the BCD connector. Plugging an external event cable into the BCD connector can damage the A&P board.

Part number G1530-60590

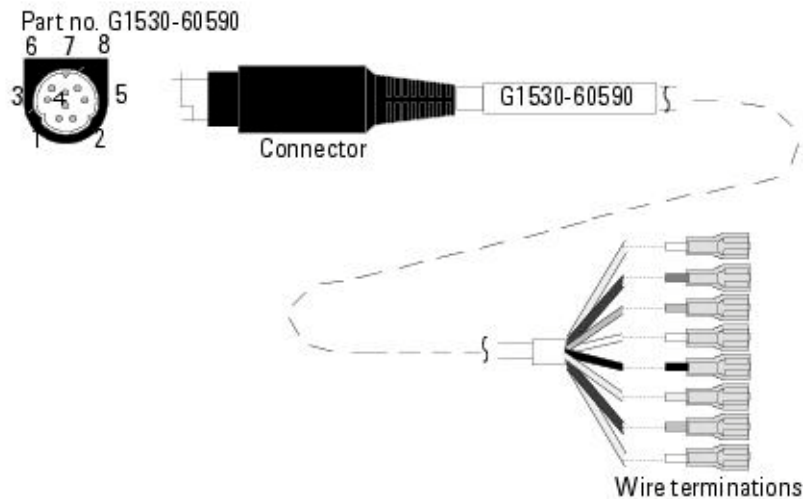


Connector	Signal name	Maximum rating	Wire terminations	Corresponds to valve #
24 volt control output				
1	24 volt output 1	150 mA	Yellow	5
2	24 volt output 2	150 mA	Black	6
3	Ground		Red	
4	Ground		White	
Relay contact closures (Normally open)				
5	Contact closure 1	48V ac/dc, 250 mA	Orange	7
6	Contact closure 1		Green	7
7	Contact closure 2	48 V ac/dc, 250 mA	Brown or Violet	8
8	Contact closure 2		Blue	8

## Automatic sampler for GC

There are three connectors, Sampler 1, Sampler 2, and Tray, for the 7683 Autosampler system.

### Cable pinouts, GC to G2613A Injector


Part no. G2613-60590

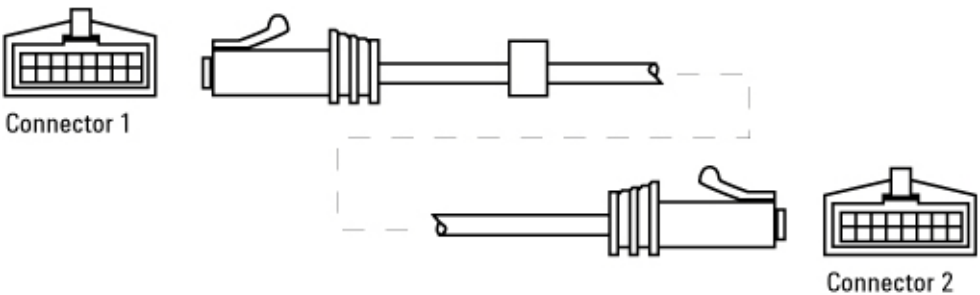

Connector 1 —16 pin (female)		Connector 2 — 16 pin (female)	
A1	RXD	A1	TXD
A2	CTS	A2	RTS
A3	DSR	A3	DTR
A4	NCTL Reset	A4	NCTL Reset
A5	GND	A5	GND
A6	VAC 1	A6	VAC 1
A7	GND	A7	GND
A8	VAC 2	A8	VAC 2
B1	TXD	B1	RXD
B2	RTS	B2	CTS
B3	DTR	B3	DSR
B4	GND	B4	GND
B5	GND	B5	GND
B6	VAC 1	B6	VAC 1
B7	GND	B7	GND
B8	VAC 2	B8	VAC 2

## Cable pinouts, GC to G4513A and 7650 Injector


Part no. G4514-60610

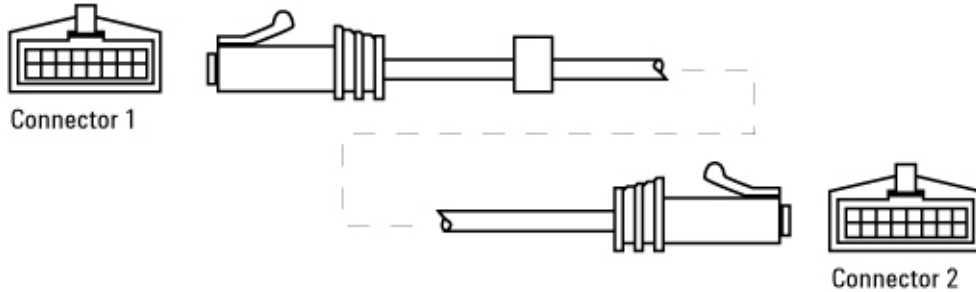

Connector 1 —16 pin (female)		Connector 2 — 16 pin (female)	
A1	RXD	A1	TXD
A2	CTS	A2	RTS
A3	DSR	A3	DTR
A4	NCTL Reset	A4	NCTL Reset
A5	GND	A5	GND
A6	VAC 1	A6	VAC 1
A7	GND	A7	GND
A8	VAC 2	A8	VAC 2
B1	TXD	B1	RXD
B2	RTS	B2	CTS
B3	DTR	B3	DSR
B4	GND	B4	GND
B5	GND	B5	GND
B6	VAC 1	B6	VAC 1
B7	GND	B7	GND
B8	VAC 2	B8	VAC 2

Cable pinouts, GC to G2614A Tray

Part no. G2614-60610

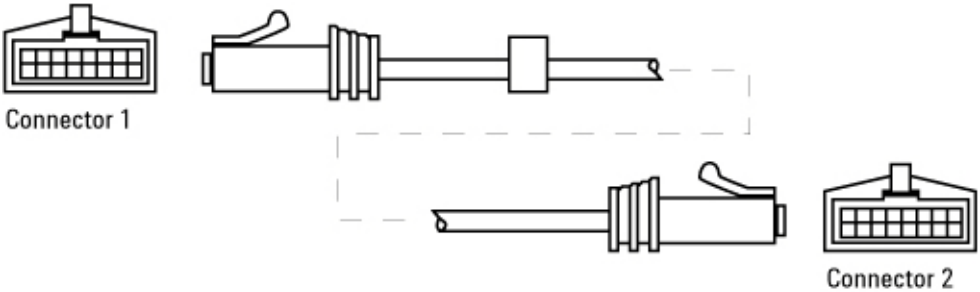


Connector 1 —16 pin (female)		Connector 2 — 16 pin (female)	
A1	RXD	A1	TXD
A2	CTS	A2	RTS
A3	DSR	A3	DTR
A4	M Reset	A4	M Reset
A5	GND	A5	GND
A6	VAC 1	A6	VAC 1
A7	GND	A7	GND
A8	VAC 2	A8	VAC 2
B1	TXD	B1	RXD
B2	RTS	B2	CTS
B3	DTR	B3	DSR
B4	GND	B4	GND
B5	GND	B5	GND
B6	VAC 1	B6	VAC 1
B7	GND	B7	GND
B8	VAC 2	B8	VAC 2



## Cable pinouts, GC to G4514A Tray

Part no. G4514-60610

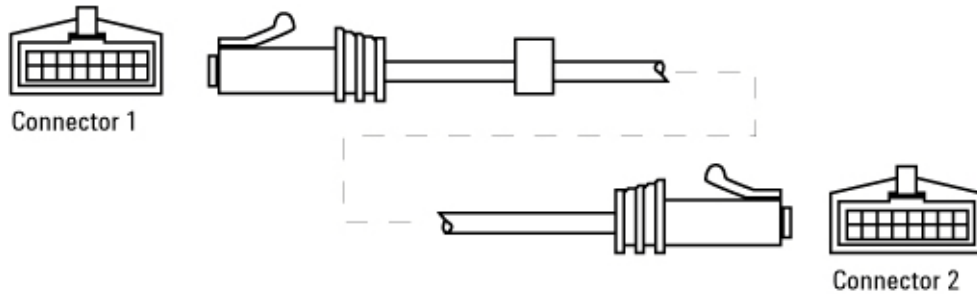

Connector 1 —16 pin (female)		Connector 2 — 16 pin (female)	
A1	RXD	A1	TXD
A2	CTS	A2	RTS
A3	DSR	A3	DTR
A4	M Reset	A4	M Reset
A5	GND	A5	GND
A6	VAC 1	A6	VAC 1
A7	GND	A7	GND
A8	VAC 2	A8	VAC 2
B1	TXD	B1	RXD
B2	RTS	B2	CTS
B3	DTR	B3	DSR
B4	GND	B4	GND
B5	GND	B5	GND
B6	VAC 1	B6	VAC 1
B7	GND	B7	GND
B8	VAC 2	B8	VAC 2


# Analog signal outputs

There are two channels of analog output available on the back panel, labeled *Sig 1* and *Sig 2*. Two cables are available – one for 3395/6 series integrators and one for general use.

## Analog cable: GC to 3395A/B or 3396B/C Integrators and 35900 C/D/E Analog to Digital Interface instrument

Part number G1530-60570

Part no. G1530-60570

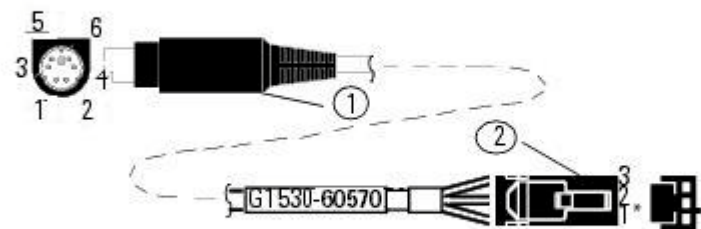


\*1 next to triangle etched on connector

Connector 1	Signal Name	Color	Connector 2
4	1 V	Black	3
2	Common	White	2
Shell	Ground	Orange	1

## Analog cable: general use

Part number G1530-60560

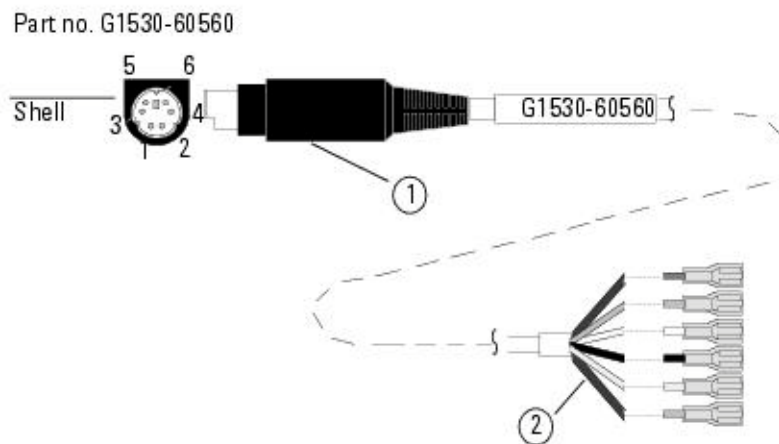


Connector 1	Signal name	Connector 2 — quick disconnects
1	not used	Brown
2	0 to 1 V, 0 to 10 V(–)	White
3	not used	Red
4	1 V (+)	Black
6	10 V (+)	Blue
Shell	Ground	Orange

## Analog signal inputs

The data from external detectors, such as the chemiluminescence detectors, can be imported into the GC using the appropriate signal input cable and the analog input board.

### General purpose AIB cable

#### G1530-60560 Signal cable, general purpose, 6 pins

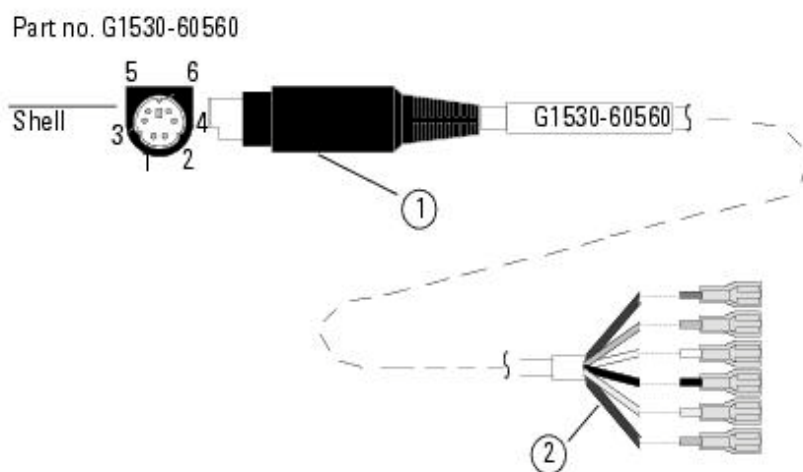


Connector 1	Signal name	Connector 2 — quick disconnects
1	not used	Brown
2	0 to 1 V, 0 to 10 V(–)	White
3	not used	Red
4	1 V (+)	Black
6	10 V (+)	Blue
Shell	Ground	Orange

### XCD Chemiluminescence detector AIB cable

#### G6600-60800 Recorder Cable



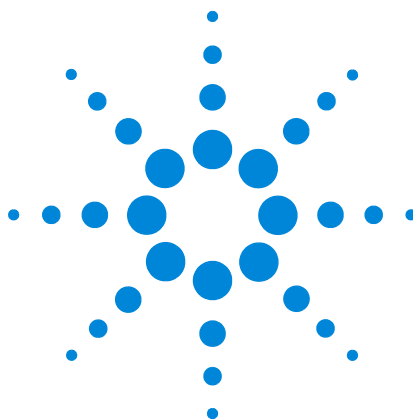

## 8 Internal Cabling

Internal cabling overview	314
MSD transfer line heater cable	315
Third detector heated zone cable	316
Communication harness	317
Motor Harness & AC Control Assembly	320
Keyboard and display connection board (original style only)	324
AC power inlet assemblies	325
Inlet/Detector harness	326
Auxiliary zone/Valve box harness	328
Third detector EPC communication cable	330
Keyboard/Display harness	331
Keyboard interconnect board	332
Heated zone extension cable	333
Ignitor cable, FID	334
Ignitor cable, FPD	335
RS-232 cable, ALS controller	336
PTV thermocouple cable	337
NPD power cable	339

This section covers wiring harnesses in the gas chromatograph. The following wiring harnesses are described:

- Communication harness
- Motor harness and AC control assembly
- Third detector EPC communication cable
- Inlet/Detector harness
- Aux zone / valve box cable
- Third detector heated zone cable
- Keyboard/Display harness
- Interconnect board, keyboard
- Ignitor cable, FID
- Ignitor cable, FPD
- RS-232 cable, ALS controller
- PTV thermocouple cable
- NPD power cable
- MSD transfer line heater cable



## **Internal cabling overview**

The various modules that make up the GC are powered and interconnected by a set of cables and wiring harnesses. These modules include inlets, detectors, temperature controllers, printed circuit boards, and power supplies.

## MSD transfer line heater cable

**Part number G1530-60790**


This cable provides heater power and temperature sensing for an MSD transfer line heater.



### Third detector heated zone cable

**Part number G3432-60506**

This is a 4-conductor cable that provides heater power and sensor connections for a third detector mounted on the left side of the GC.

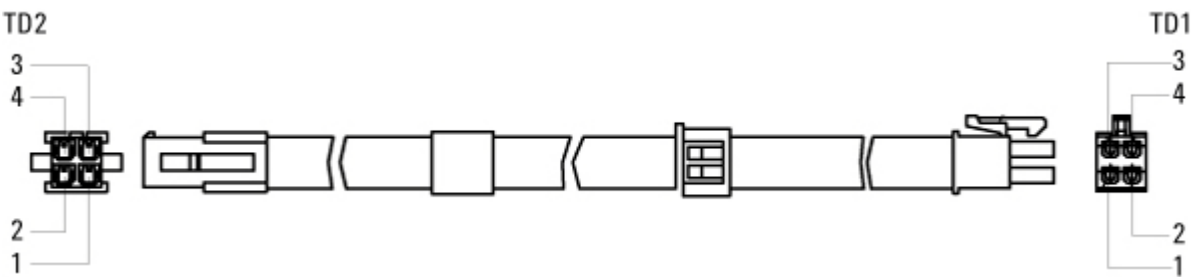

**Table 21** Connection table


TD1 connector pin	Wire color	Twisted pair?	TD2 connector pin
1	Orange	Yes	1
3	Black		3
2	Grey	Yes	2
4	Green		4



## Communication harness

### Part number G3430-60513

This harness provides 24 VDC from the Analog and Power board to all EPC modules. It also carries communications between the EPC modules and the signal board. The harness connector labeled "power" connects to P1 on the Analog and Power board. The harness connector labeled "Comm Bus" connects to P7 on the logic board.

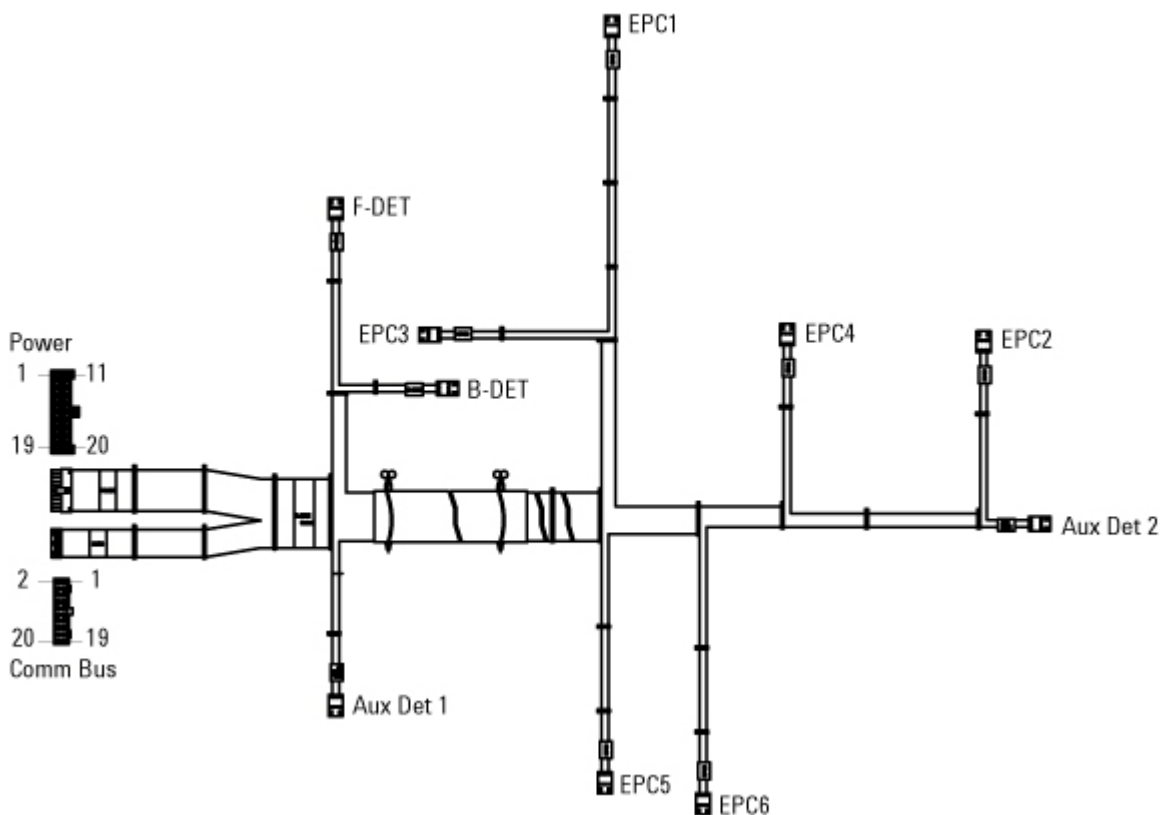

## Power connection table

Black wires, that are not part of twisted pairs, are GROUND for pin 2 of the associated EPC.

Red wires, that are not part of twisted pairs, supply + 24 VDC to pin 1 of the associated EPC module.

Power connector pin	Wire color	Signal	EPC (pin)
1	Red	+ 24 VDC	EPC 1 (1)
2	Red	+ 24 VDC	EPC 2 (1)
3	Red	+ 24 VDC	EPC 3 (1)
4	Red	+ 24 VDC	EPC 4 (1)
5	Red	+ 24 VDC	EPC 5 (1)
6	Red	+ 24 VDC	EPC 6 (1)
7	Red	+ 24 VDC	AUX DET 2 (1)
8	Red	+ 24 VDC	AUX DET 1 (1)
9	Red	+ 24 VDC	B DET (1)
10	Red	+ 24 VDC	F DET (1)
11	Black	GROUND	EPC 1 (2)
12	Black	GROUND	EPC 2 (2)
13	Black	GROUND	EPC 3 (2)
14	Black	GROUND	EPC 4 (2)
15	Black	GROUND	EPC 5 (2)
16	Black	GROUND	EPC 6 (2)
17	Black	GROUND	AUX DET 2 (2)
18	Black	GROUND	AUX DET 1 (2)
19	Black	GROUND	B DET (2)
20	Black	GROUND	F DET (2)

## Data communication table

Data are transmitted to and from the EPC modules on twisted pairs of wires.


Comm Bus connector pin	Wire color	Twisted pair?	EPC (pin)
1	White	Yes	EPC 1 (3)
2	Black		EPC 1 (4)
3	White	Yes	EPC 2 (3)
4	Brown		EPC 2 (4)
5	White	Yes	EPC 3 (3)
6	Red		EPC 3 (4)
7	White	Yes	EPC 4 (3)
8	Orange		EPC 4 (4)
9	White	Yes	EPC 5 (3)
10	Yellow		EPC 5 (4)
11	White	Yes	EPC 6 (3)
12	Green		EPC 6 (4)
13	White	Yes	AUX DET 2 (3)
14	Blue		AUX DET 2 (4)
15	White	Yes	AUX DET 1 (3)
16	Violet		AUX DET 1 (4)
17	White	Yes	B DET (3)
18	Grey		B DET (4)
19	White	Yes	F DET (3)
20	White/Black		F DET (4)

## Motor Harness & AC Control Assembly

### Part number G3430-60512

This harness provides 24 VDC power from the Analog and Power board at P18 for the following functions:

- Flapper motor (6 pin, 2 x 3 connector)
- Inlet fan motor (3 pin)
- Oven sensor (2 pin, shortest)
- Pneumatics fan motor (2 pin)
- Cryogenic cooling valve for the oven (5 pin)
- Shut-off solenoid on the AC board (5 pin, 3 used)

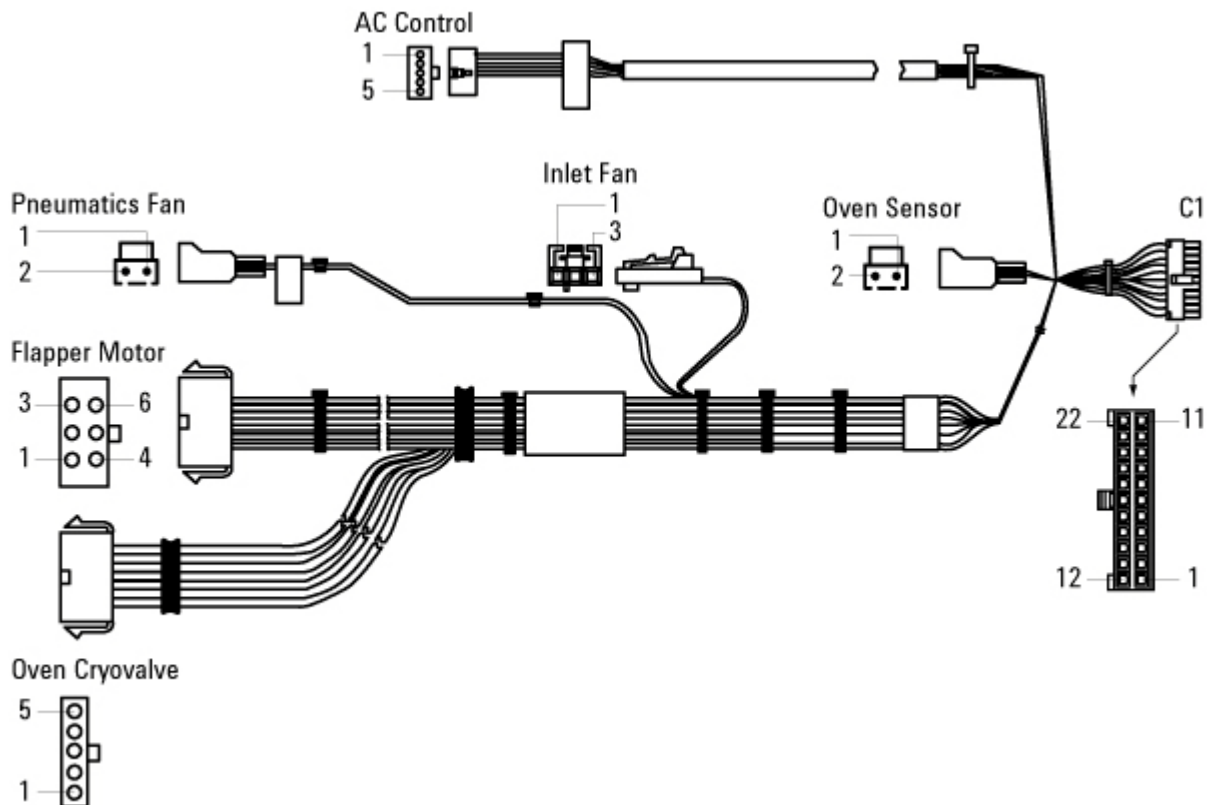


**Table 22** Connection table

<b>C1 connector pin</b>	<b>Wire color</b>	<b>Destination</b>
1	Black	Flapper motor 2
2	White/Blue	Flapper motor 4
3	Red	Flapper motor 6
4	Orange	Inlet fan 2
5	White/Green	Inlet fan 3
6	Blue	Pneumatics fan 2
7	Green	Oven cryovalve 3
8	Brown	Oven cryovalve 4
9	White	AC control 3
10	no connection	
11	Violet	Oven sensor 2
12	White	Flapper motor 1
13	White/Black	Flapper motor 3
14	Blue	Flapper motor 5
15	Grey	Oven cryovalve 5
16	Yellow	Inlet fan 1
17	Brown	Pneumatics fan 1
18	Violet	Oven cryovalve 2
19	White/Red	Oven cryovalve 5
20	Black	AC control 2
21	Red	AC control 1
22	Black	Oven sensor 1

### Keyboard and display connection board (original style only)

This board combines the oven door sensor (P2) with the keyboard and display ribbon cable (P1) and connects them to the keyboard and display harness (P3).

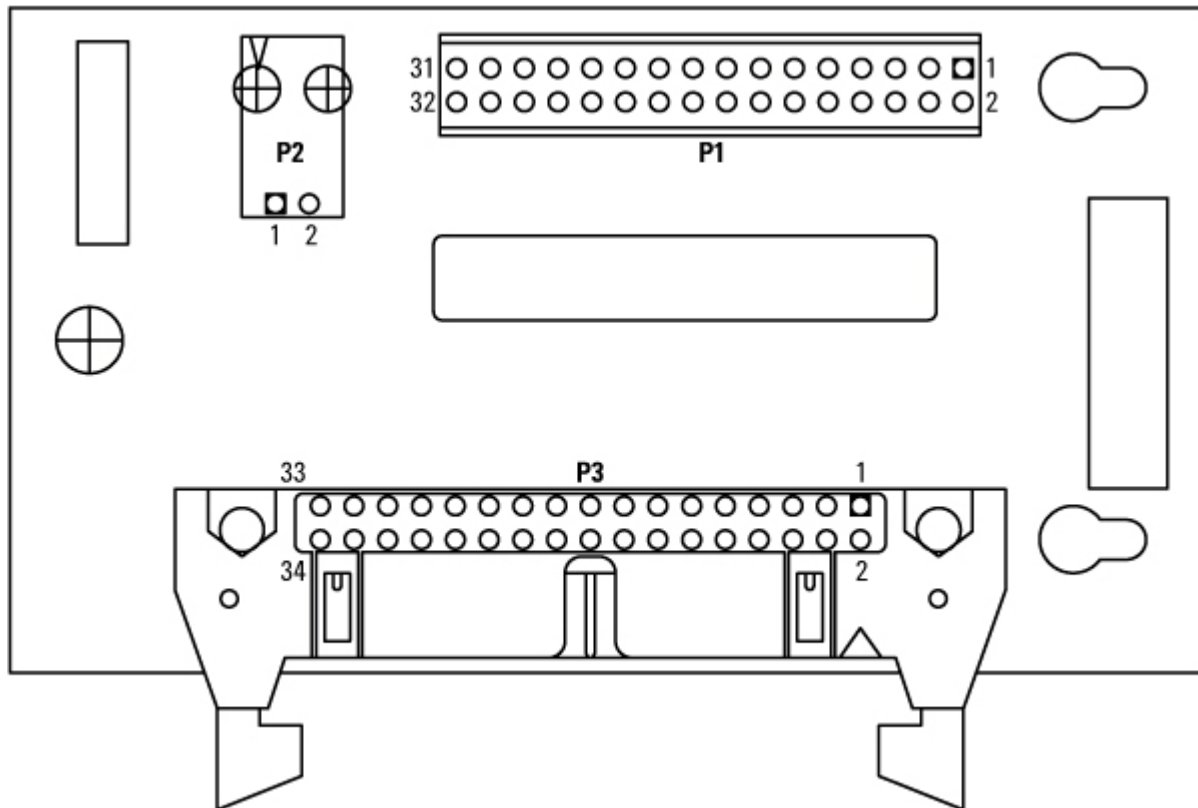

P1 pins are connected to the corresponding pins on P3 (P2-1 to P3-1, etc.).


P2 pins are connected to P3 pins 33 and 34.

AC power inlet assemblies

These assemblies provide the connection between the external power cord and the internal AC power board.


C20 assembly (for C19 plug)

Part number G1530-61550

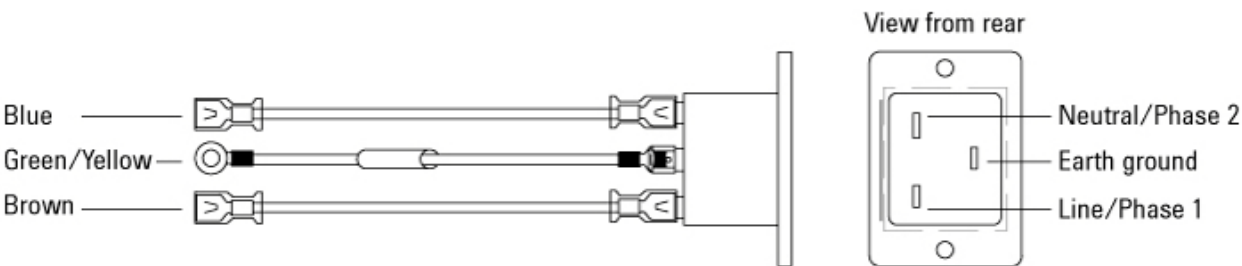


C14 assembly (for C13 plug)

Part number G1530-61560

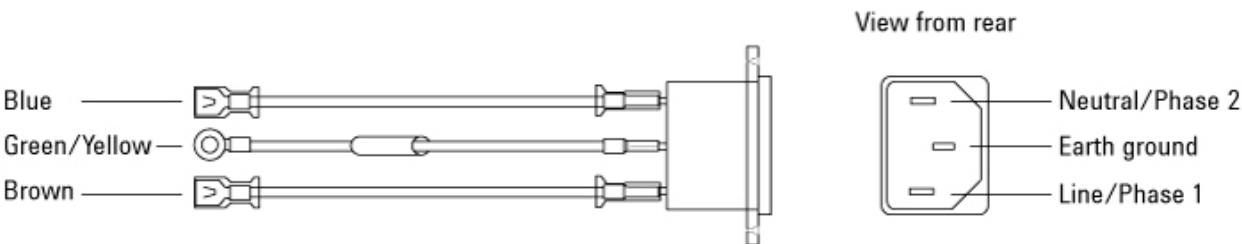

## Keyboard and display connection board (original style only)

This board combines the oven door sensor (P2) with the keyboard and display ribbon cable (P1) and connects them to the keyboard and display harness (P3).

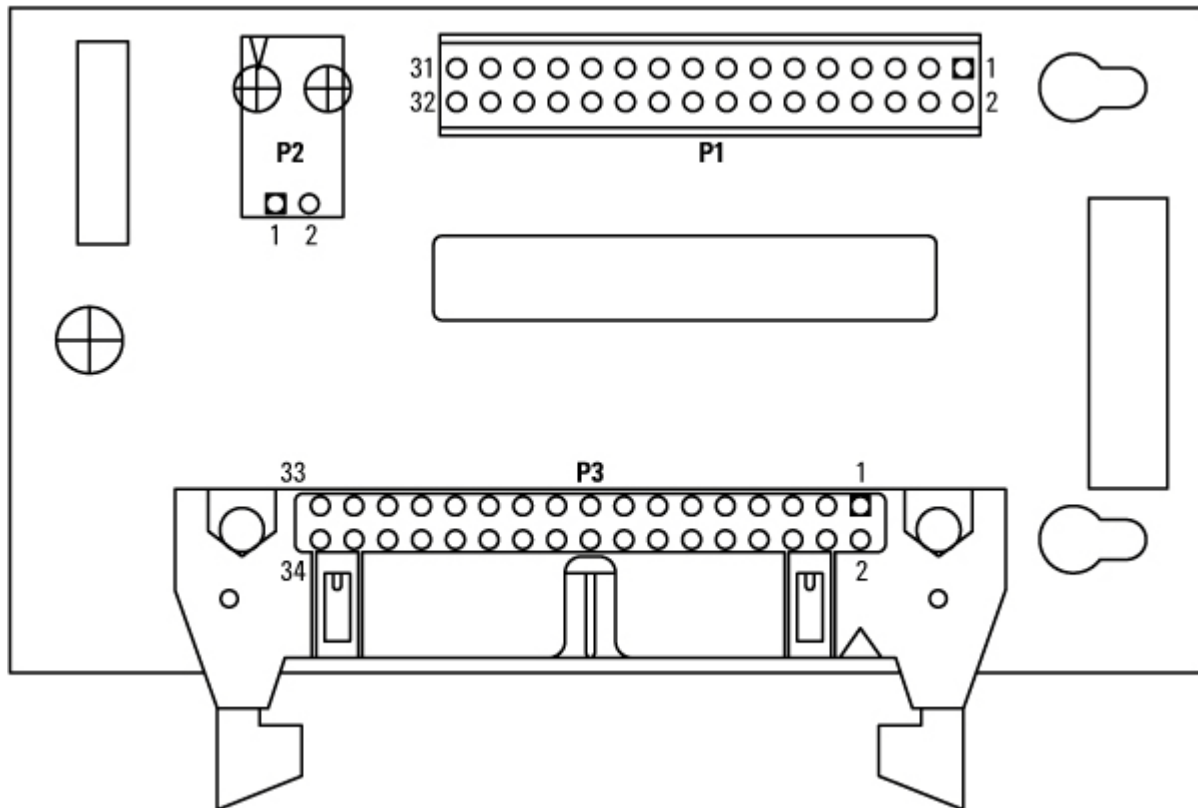

P1 pins are connected to the corresponding pins on P3 (P2-1 to P3-1, etc.).

P2 pins are connected to P3 pins 33 and 34.




# AC power inlet assemblies

These assemblies provide the connection between the external power cord and the internal AC power board.


## C20 assembly (for C19 plug)

Part number G1530-61550

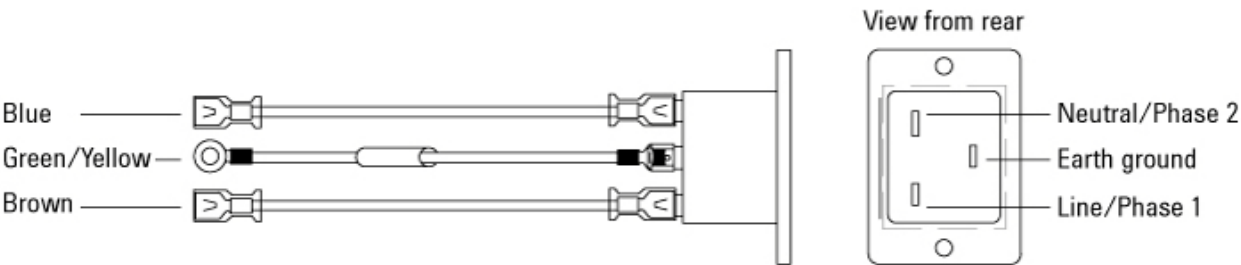

## C14 assembly (for C13 plug)


Part number G1530-61560

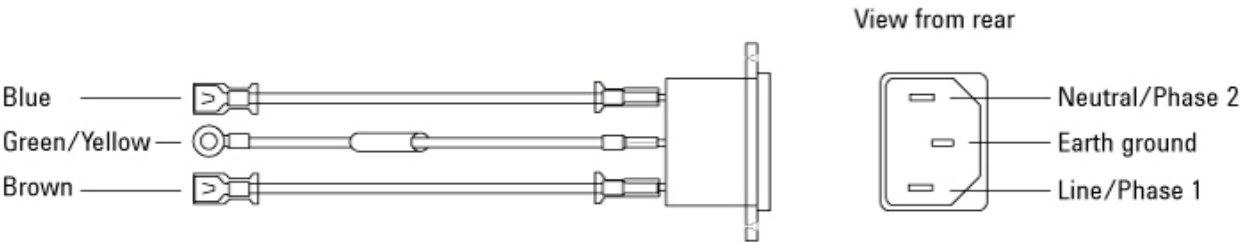

## Inlet/Detector harness

### Part number G3430-60535

This harness connects the inlet and detector EPCs to the logic board at P7.

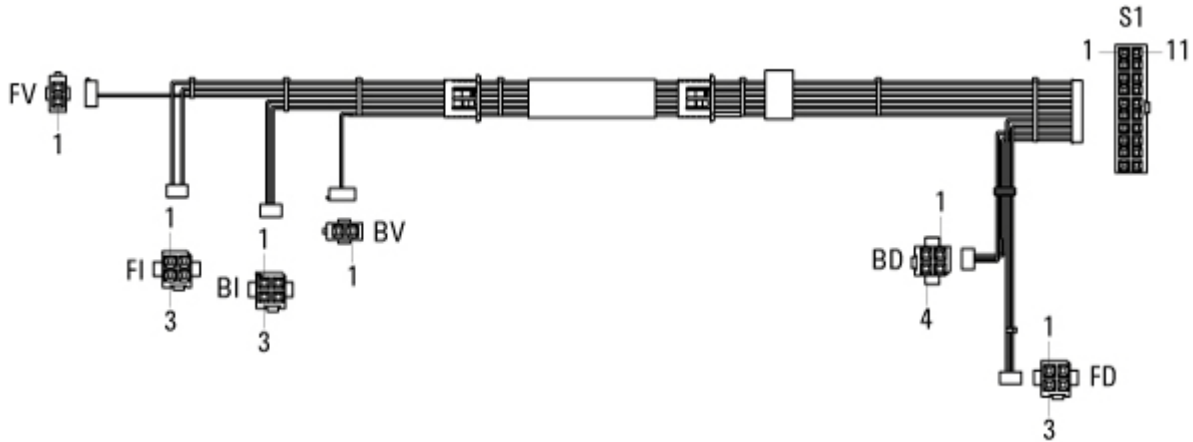

The FI and BI connectors are for inlets; FD and BD are for detectors. FV and BV are two-conductor connectors mounted on the top left edge of the GC near the inlets. They supply 24 VDC to various accessories.

**Table 23** Connection table

S1 pin	Wire color	Signal name	Destination (pin)
1	Black	ADC-GND	F1 (3)
2	Black	ADC-GND	BI (3)
3	Black	ADC-GND	FD (3)
4	Black	ADC-GND	BD (3)
5	Green	C GROUND	FI (4)
6	Green	C GROUND	BI (4)
7	Green	C GROUND	FD (4)
8	Green	C GROUND	BD (4)
9	Yellow	+ 24 V	FV (1)
10	Yellow	+ 24 V	BV (1)
11	Violet	INJA-SENSE	FI (2)


S1 pin	Wire color	Signal name	Destination (pin)
12	Grey	INJB-SENSE	BI (2)
13	Violet	DETA-SENSE	FD (2)
14	Grey	DETB-SENSE	BD (2)
15	Red	INJ-FRONT	FI (1)
16	Orange	INJ-BACK	BI (1)
17	Red	DET-FRONT	FD (1)
18	Orange	DET-BACK	BD (1)
19	Blue	ICRYOSPLT-FRONT	FV (2)
20	White	ICRYOSPLT-BACK	BV (2)

## Auxiliary zone/Valve box harness

### Part number G1530-60660

The Auxiliary Zone/Valve Box harness runs from connector P22 on the A&P board to the valve driver bracket screwed on to the right side of the GC above the A&P board. It powers the valve actuator drivers and two auxiliary heater/sensors.

The A1 and A2 connectors have a tab which helps ID which wire is the heater, sense, or ground.

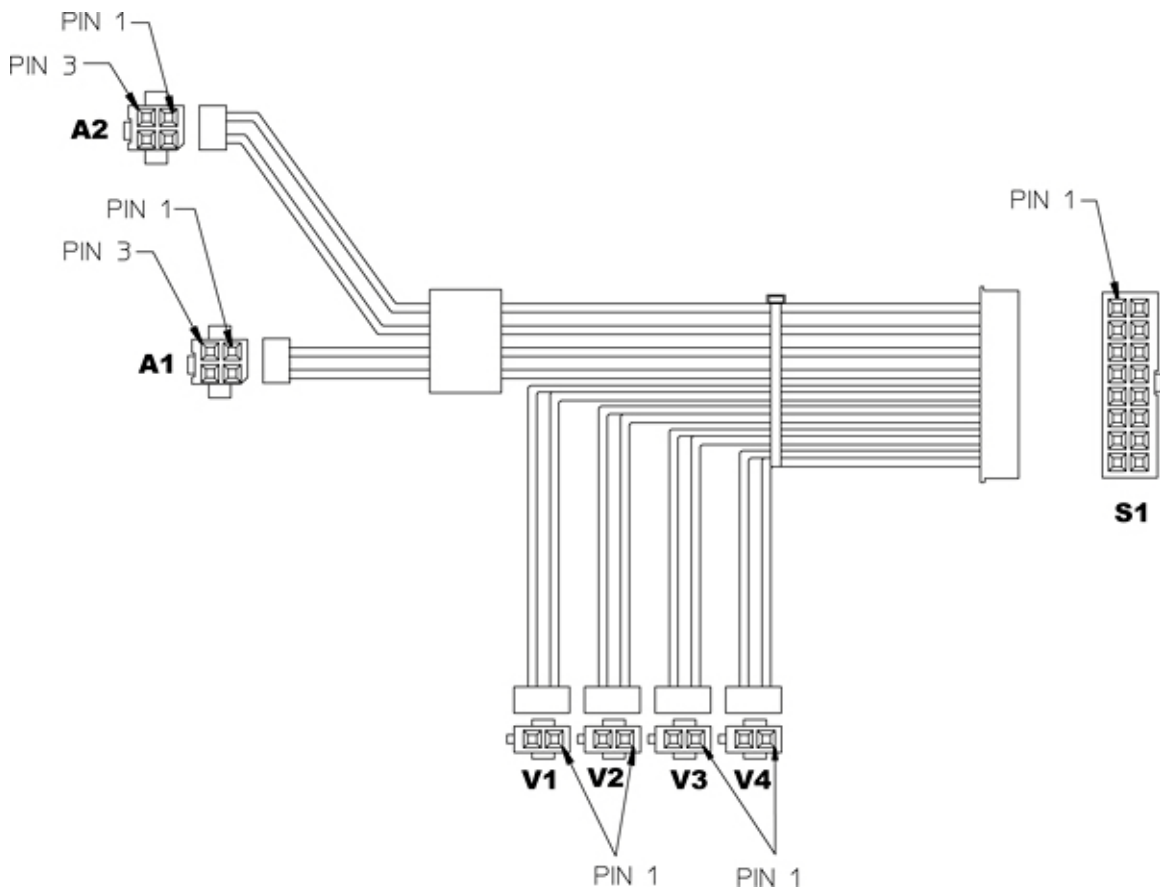

**Table 24** Connection table


Wire color	Pin # on S1	Signal name	Destination (and Pin #)
Violet	1	Aux 1 Sense	A1 (2)
	2	Aux 2 Sense	A2 (2)
Green	3	Heater GND	A1 (4)
	4	Heater GND	A2 (4)
Yellow	5	+24 Volts	V1 (1)
	6	+24 Volts	V2 (1)
	7	+24 Volts	V3 (1)
	8	+24 Volts	V4 (1)
Black	9	Sensor GND	A1 (3)
	10	Sensor GND	A2 (3)
Red	11	Aux 1 Heater	A1 (1)
Orange	12	Aux 2 Heater	A2 (1)
Brown	13	Valve #1	V1 (2)
Blue	14	Valve #2	V2 (2)
Gray	15	Valve #3	V3 (2)
White	16	Valve #4	V4 (2)

## Third detector EPC communication cable

### Part number G3432-60505

This cable extends the communications harness to a third detector mounted on the left side of the GC. It connects to the EPC 6 connector and extends across the top of the GC to the side mount.

When removing and replacing the third detector cable, note carefully how it enters the side detector carrier. You may need to loosen the screws holding the carrier to the GC body.

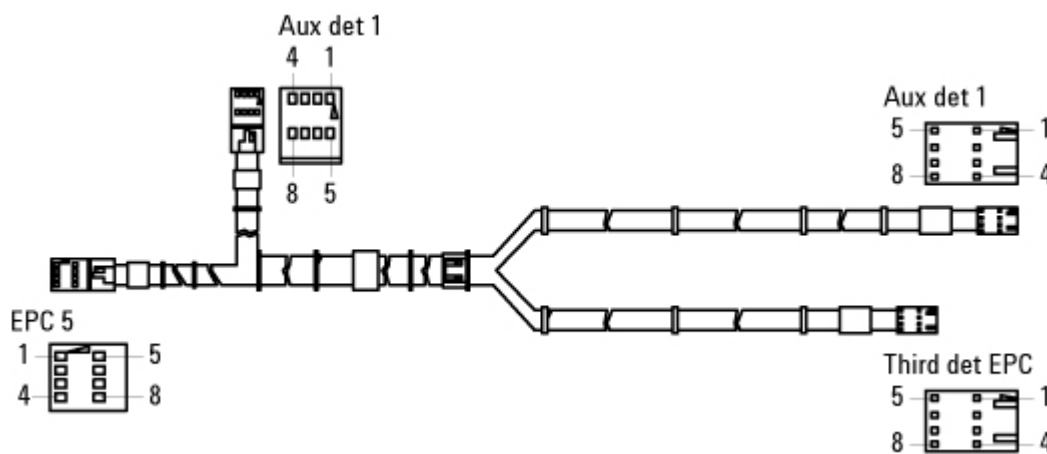

**Figure 6** EPC AreaSide Carrier Area


**Table 25** Connection table

Connector, EPC area (pin)	Wire color	Twisted pair?	Connector, side carrier area (pin)
Aux Det 1 (1)	Red	Yes	Aux Det 1 (1)
Aux Det 1 (2)	Black		Aux Det 1 (2)
Aux Det 1 (3)	White	Yes	Aux Det 1 (3)
Aux Det 1 (4)	Blue		Aux Det 1 (4)
EPC 5 (1)	Red	Yes	Third Det EPC (1)
EPC 5 (2)	Black		Third Det EPC (2)
EPC 5 (3)	White	Yes	Third Det EPC (3)
EPC 5 (4)	Yellow		Third Det EPC (4)

## Keyboard/Display harness

**Part number G3430-60514**


This harness connects the keyboard and display to the logic board at P8.

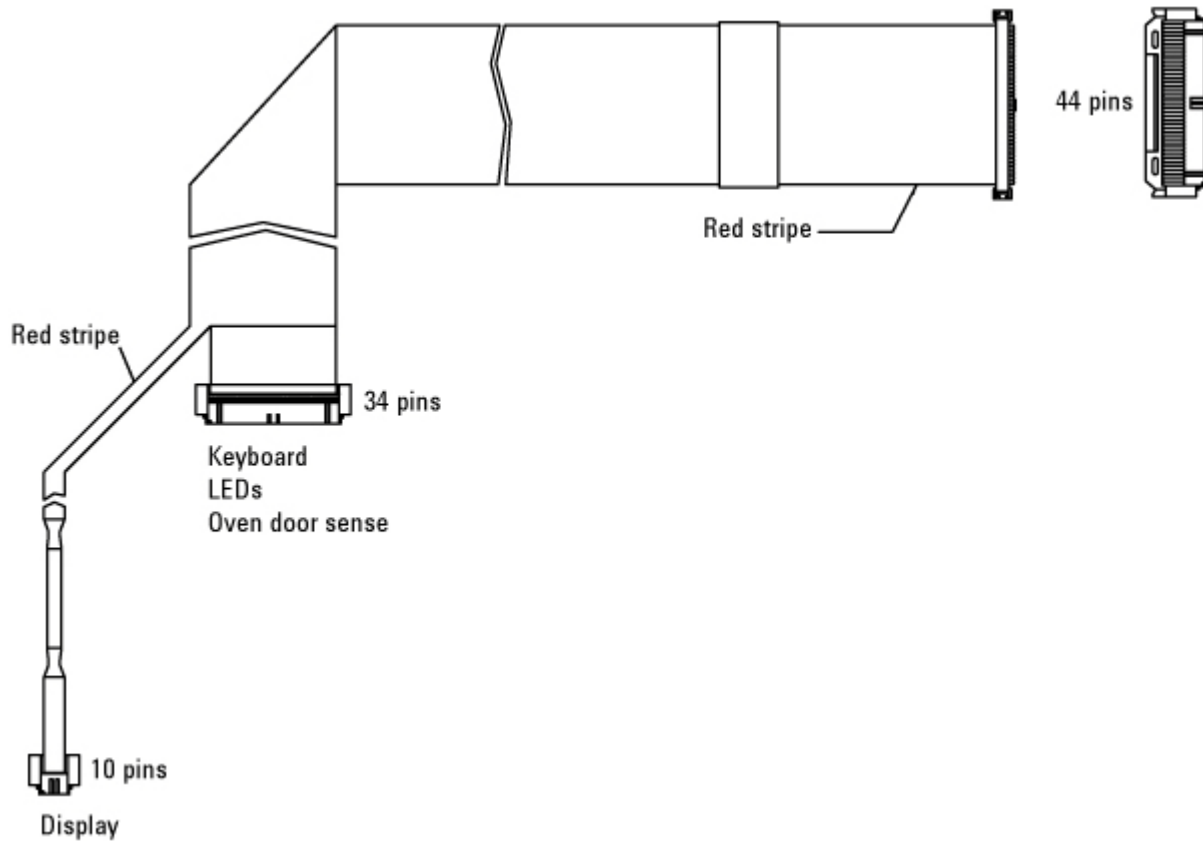

## Keyboard interconnect board

**Part number G3430-60003**

This board, mounted behind the keyboard/display module, connects the keyboard to the keyboard/display harness.

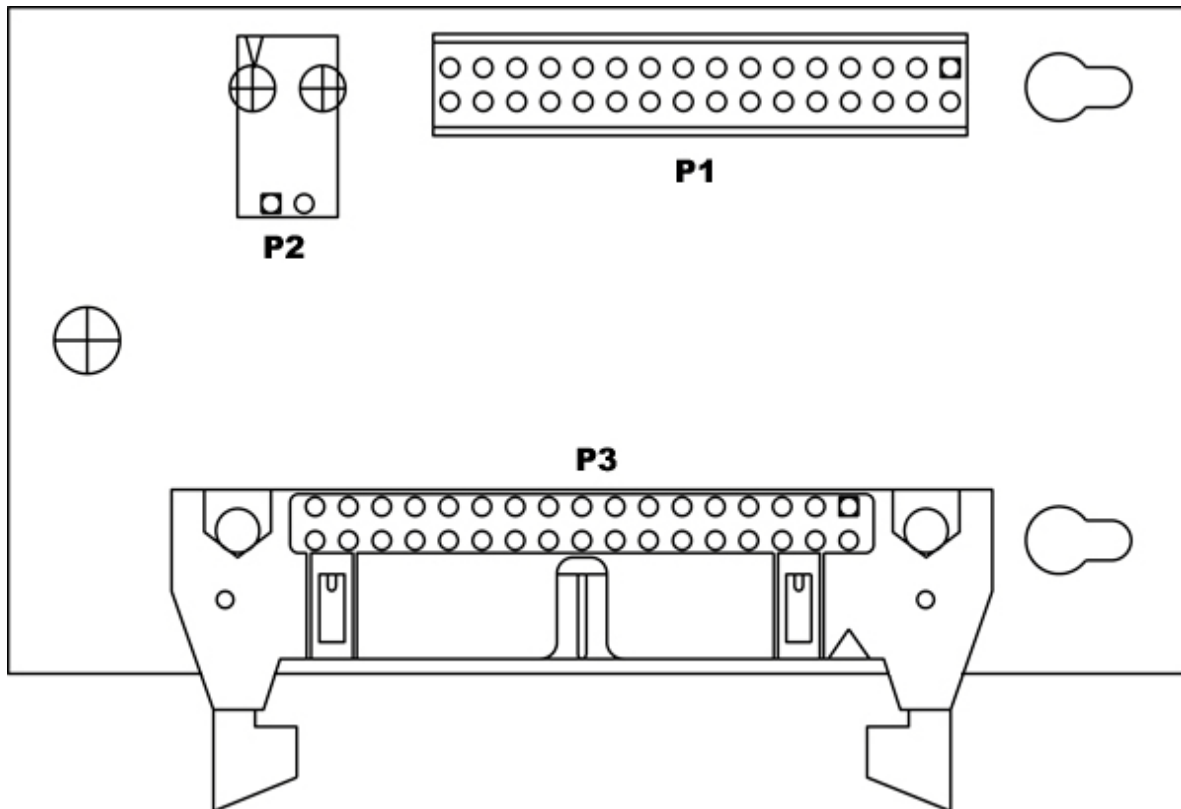

- Flex circuit from keypad attaches to P1.
- Oven door sensor attaches to P2.
- Ribbon cable from logic board attaches to P3.



Heated zone extension cable

Part number G1530-60790




Table 26 Connection table

Connector C1 pin	Wire color	Connector C2 pin
1	Orange	1
2	Violet	2
3	Black	3
4	Green	4

## Ignitor cable, FID

**Part number G3431-60680**

The ignitor cable connects the glow plug ignitor to its power source.

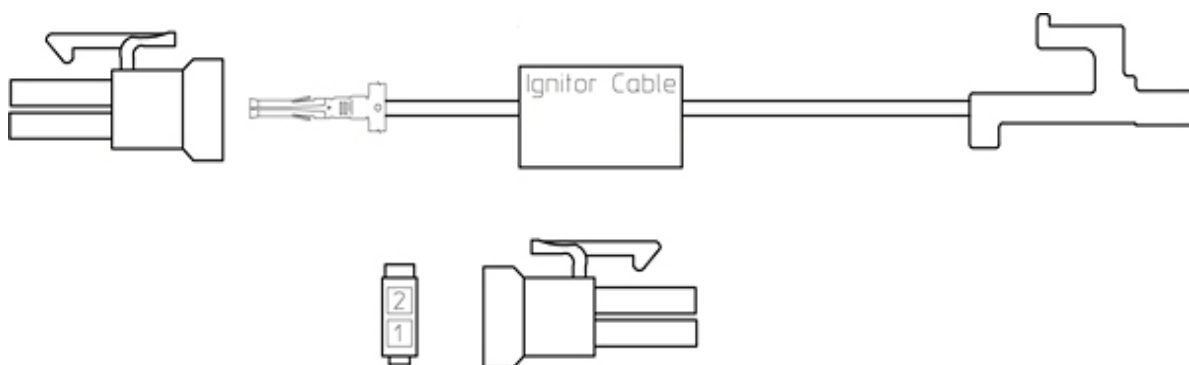

## Ignitor cable, FPD

**Part number G3435-60600**

This cable connects the glow plug ignitor to its power source.



RS-232 cable, ALS controller

Part number G2612-60510

The ALS controller is driven by the analog and power board.  
This cable makes the connection at J8.





J8 pin	Signal
1	+5 V
2	no connection
3	+5 V
4	DSRA
5	Ground
6	+5 V
7	Ground
8	RSTA
9	Ground
10	Ground
11	Ground
12	Ground
13	nStart (ALS Start paralleled with APG Start to FPGA as nAPG_START)
14	Spare.VDSP
15	RxDA
16	Spare.VDSN
17	TxDA
18	DTRA
19	CTSA
20	RTSA

## PTV thermocouple cable

### Part number G2617-60505

The Thermocouple Cable connects the PTV thermocouple board to the PTV inlet heater/sensor cable, and to the valve and power connectors on the inlet/detector wiring harness.

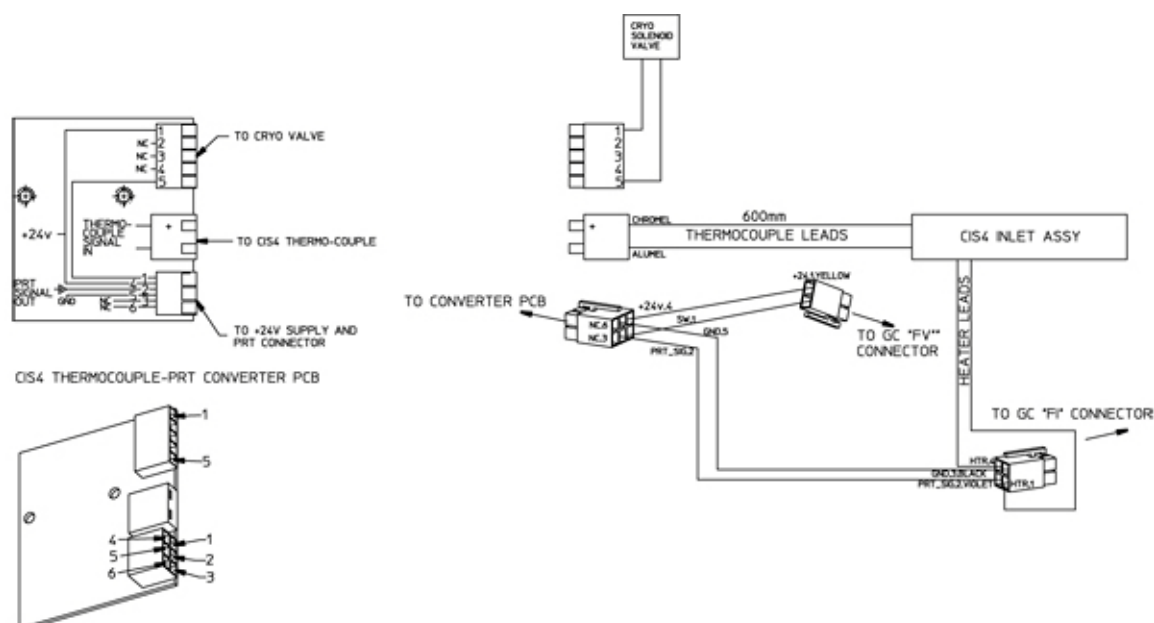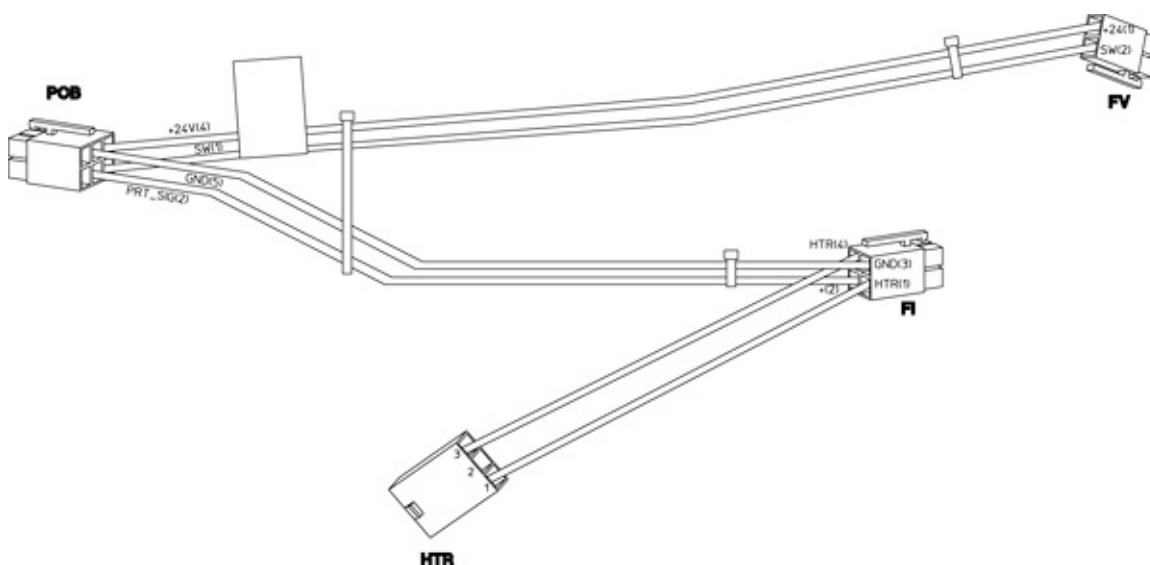

**Table 27** Connection tables

PCB pin	Wire color	Signal	Destination (pin)
1	Blue	ICRYO_SPLT_FRO NT	FV (2)
2	Violet	INJA_SENSE	FI (2)
3	No connection		
4	Yellow	+ 24 VDC	FV (1)
5	Black	ACC-GND	FI (3)
6	No connection		

HTR pin	Wire color	Destination (pin)
1	Yellow	FI (1)
2	No connection	
3	Blue	FI (4)

# NPD power cable

Part number G3434-60600

This cable provides power to the heated bead in the NPD.


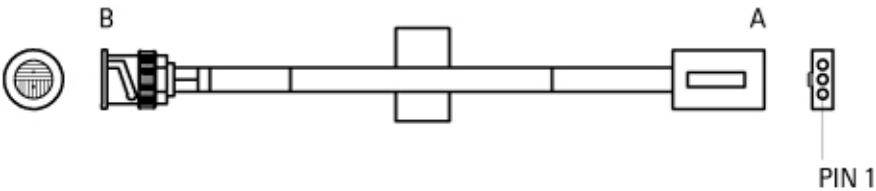


Table 28 Connection table

Connector A pin	Wire color	Connector B pin
1	Braided shield	Shell
2	Blue	Female
3	Clear	Male

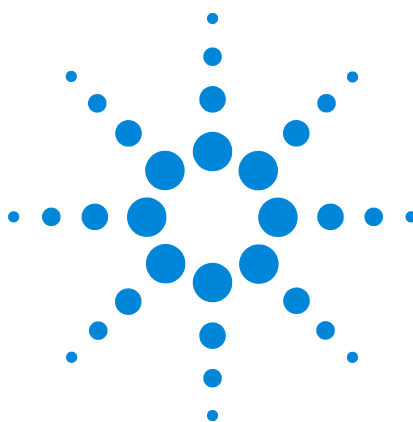

## 9 Mainframe

Mainframe overview	341
Covers and Fans	342
Replacing the oven bezel	348
Replacing the oven door	349
Replacing the oven door button	350
Replacing the keypad assembly and display board	351
Oven temperature control	362
Oven temperature troubleshooting	363
Testing resistance of the heater coil	365
Cryo valve installation/replacement	366
Replacing the oven sensor	371
Replacing the oven door sensor	372
Replacing the oven shroud, oven fan, and oven fan motor	375
Replacing the oven flapper assembly	381
Replacing the flapper or flapper motor	383
Replacing Components Inside the Electronics Panel	385
Replacing Components Inside the Lower Rear Metal Cover	400





## Mainframe overview

The GC mainframe consists of the chassis (baseplate) plus the following:

- A set of metal and plastic covers
- The oven, oven door, and related parts
- Electronic components inside the right side cover. These are the "mind" of the GC.
- Power handling components inside the lower back cover
- Wiring cables and harnesses

Inlets, detectors, valves, and other selected components are mounted on or in the mainframe.

## **Covers and Fans**

Metal and plastic covers protect the GC components and guard the user from exposure to heat and electric shock. The covers also play a part in controlling the flow of cooling air through the instrument. The GC should not be operated unless all covers are in place.

Additional cooling is needed in two areas, the inlets and the EPC pneumatic controllers. The inlet fan is located at the rear of the inlet cover. The pneumatics fan is located on the back of the GC under the center of the pneumatics slots.

## Removing and replacing the covers

### Detector top cover

This cover protects the detectors, valve box, and valve assembly. To remove it, raise it to a vertical position, lift the right side, and disengage the pin on the lower left side.

#### CAUTION

Do not force the cover, either when installing it or closing it. This could break the plastic parts.

To replace the cover, make certain that the slot in the brass bushing (lower right corner) is vertical and that the bushing is fully seated. Installation is then the reverse of removal.

### The pneumatics top cover

The large plastic cover over the pneumatics area on the top, rear of the GC is held by a black button on each end. Press both buttons and lift the cover off.

### Replacing the pneumatics top cover buttons

- 1 Remove the cover from the GC.
- 2 Turn the cover over and remove the two screws holding the button to be replaced.
- 3 Assemble the plastic button and the metal mount. Position the assembly over the two nuts and install the screws.

### The electronics top cover

- 1 Raise or remove the top cover.
- 2 This exposes a screw near the left front side of the electronics cover. Loosen this screw.
- 3 Raise the cover to the vertical position.
- 4 To remove the cover completely after raising it to the vertical position, Tilt it to the right to disengage the left end, then slide it to the left and remove it.

### The side covers

The left and right side covers are each held by 2 captive screws at the bottom and a hook at the rear. Unfasten the screws, slide the panel to the rear, and remove it.

### **The rear covers**

The upper and lower rear covers are held by several screws each, but all are visible from behind the GC. Note that some of the screws are in slots rather than holes; these screws do not need to be completely removed.

### **The side cover on the third detector assembly**

Remove 2 screws at the bottom of the panel, then lift and remove the combined side and top panel.


## Replacing the inlet cooling fan

- 1 Cool down the oven and all heated zones to below 70 °C to avoid the creation of active sites.
- 2 Turn off the GC main power switch.

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

- 3 Remove the top rear cover by pressing in on the side retaining clips and lifting off.

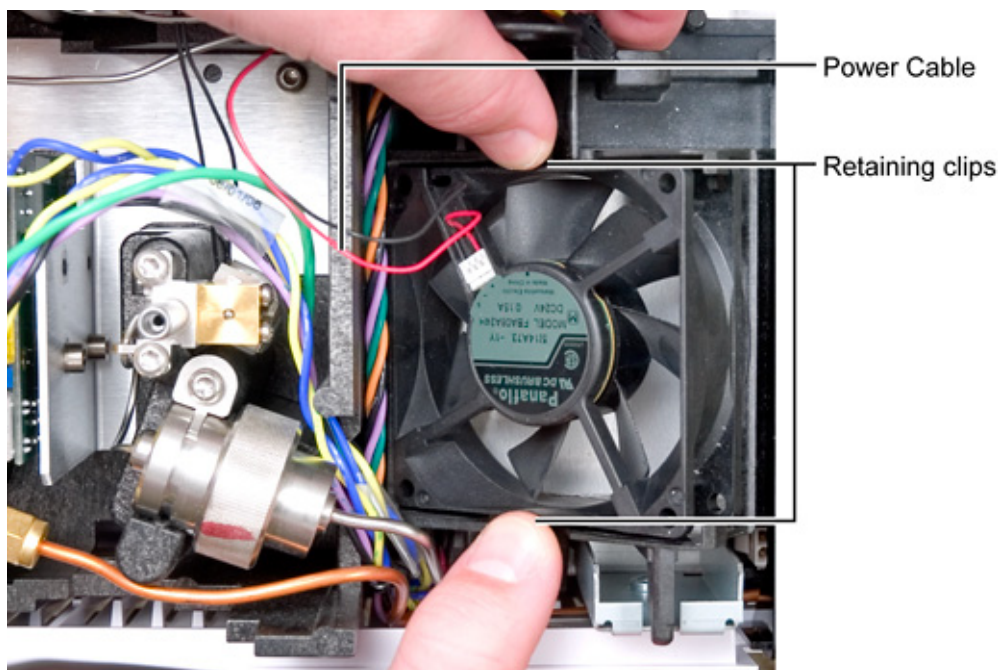

- 4 Disconnect its cable from the GC wiring harness.
- 5 Remove the fan from its clips.
- 6 In reassembly, make sure the label side of the fan faces upwards (flow must be in the upwards direction). Also, the cable should exit nearest its connector on the GC wiring harness.

## Replacing the EPC cooling fan

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

### WARNING

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**


---

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 3 Remove the lower rear cover on the GC.
- 4 The fan is located below the EPC module compartment; begin by disconnecting the fan's cable.

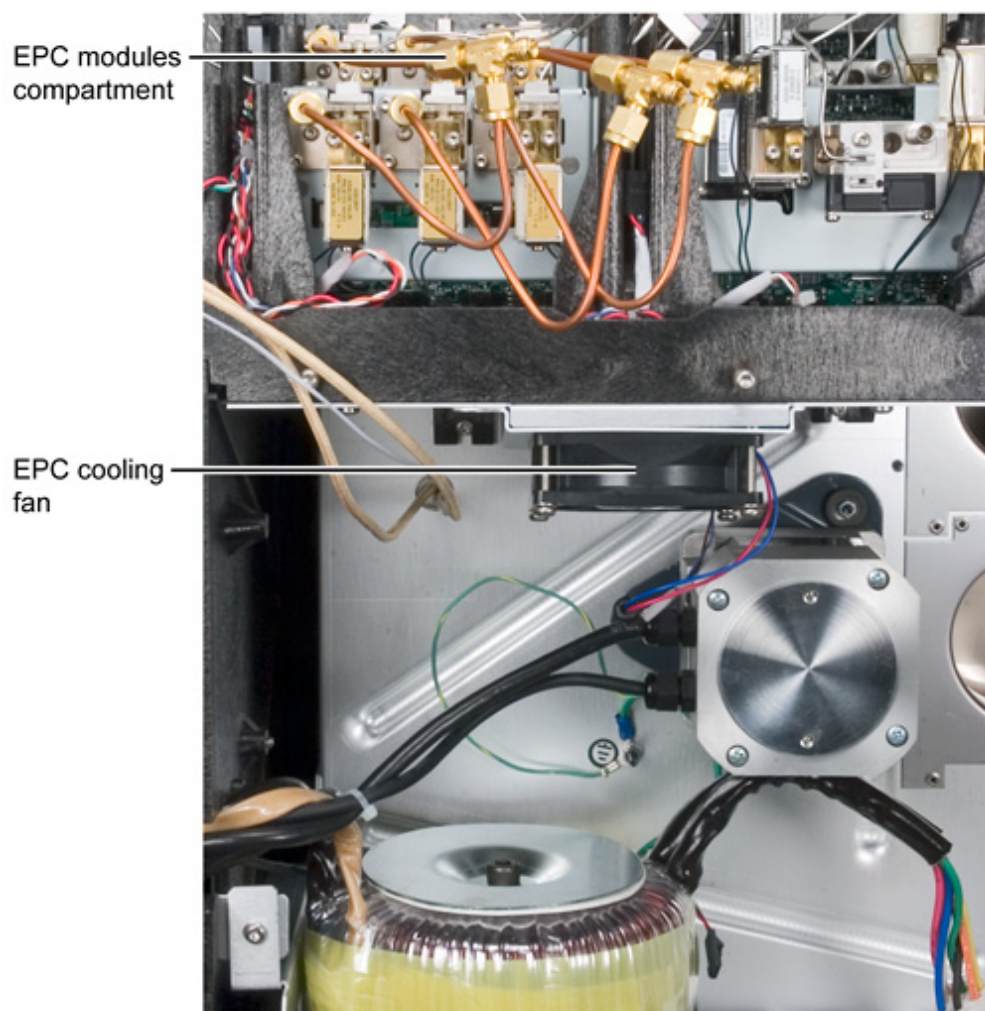

- 5 Remove 4 screws.
- 6 In replacement, be careful to not overtighten the screws such that the fan's plastic frame is cracked. Also, in replacing the screws, it may be easiest to start with the one at the right forward corner due to interference caused by the nearby oven fan motor.

**NOTE**

In reassembly, make sure fan orientation is such that air flow is directed **upwards** through the EPC modules.

## Replacing the oven bezel

**WARNING**

**Be careful! The oven and/or bezel may be hot enough to cause burns. If it is hot, wear gloves to protect your hands.**

---

The bezel consists of a plastic cover over a metal plate. It is held by two screws that are reached through holes in the top of the bezel.

- 1 Unscrew the two mounting screws.
- 2 Lift the bezel and metal plate off the GC.
- 3 When reinstalling, note that the lip on the bezel extends forward over the top of the oven door. Be very certain that the metal plate is between the plastic part and the top of the GC.



## Replacing the oven door

**WARNING**

**Be careful! The oven and/or inlet may be hot enough to cause burns. If the inlet is hot, wear gloves to protect your hands.**

---

- 1 If an ALS is installed, remove it.
- 2 Remove the GC pneumatics cover and the top back panel.
- 3 Remove the six T-20 screws retaining the inlet cover, lift off and remove the cover.
- 4 Open the oven door. This exposes the top of the shaft that attaches the door to the GC.
- 5 The hinge shaft threads completely through its bracket, so that the threads cannot engage and loosen during use. To remove the shaft, use a flat-head screwdriver to lift the shaft from the bottom until it contacts the threads. Hold in place while using a T-20 driver to loosen the shaft. Turn the T-20 driver a few times to engage the threads on the bottom of the shaft with the top threaded plate of the bottom door hinge.
- 6 While supporting the weight of the door, use pliers to pull the door shaft up and out of the door.
- 7 Remove the door.
- 8 Installation of the new door is the reverse of these steps.

## Replacing the oven door button

**WARNING**

The oven door button mechanism is under the mainframe. To reach it, you must either tip the GC on its side or slide it part way off the bench so that the right front corner is accessible. Be very certain that the GC is stable and cannot fall.

---

- 1 Move the GC so that the bottom of the mainframe near the button is accessible.
- 2 Open the oven door and operate the button several times. Note that the operating lever from the button is behind the vertical lever that moves the latch,
- 3 To remove the mechanism, remove the two screws that hold it and pull it down out of the mainframe.
- 4 Replacement is the reverse of these steps. Be certain that the button lever is *behind* the latch lever.

## Replacing the keypad assembly and display board

### Replacing the keypad assembly (original style)

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

#### **WARNING**

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**


---

#### **CAUTION**

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.


---

- 3 Remove three screws with ¼-inch nut heads from the detector electronics panel side. Two of the screws are near the upper left and lower left corners of the front detector signal board; the other is on the bottom of the GC in a slotted hole.
- 4 Disconnect the large keypad ribbon cable from the inner connector board.
- 5 Remove the T-20 screw from the bracket near the top of the keypad assembly to release it from the GC body.
- 6 From the front, slide the keypad panel forward.
- 7 At this point there are two servicing options: replacement of the entire keypad assembly or just the display. [See next topic for the replacement of the display only.]

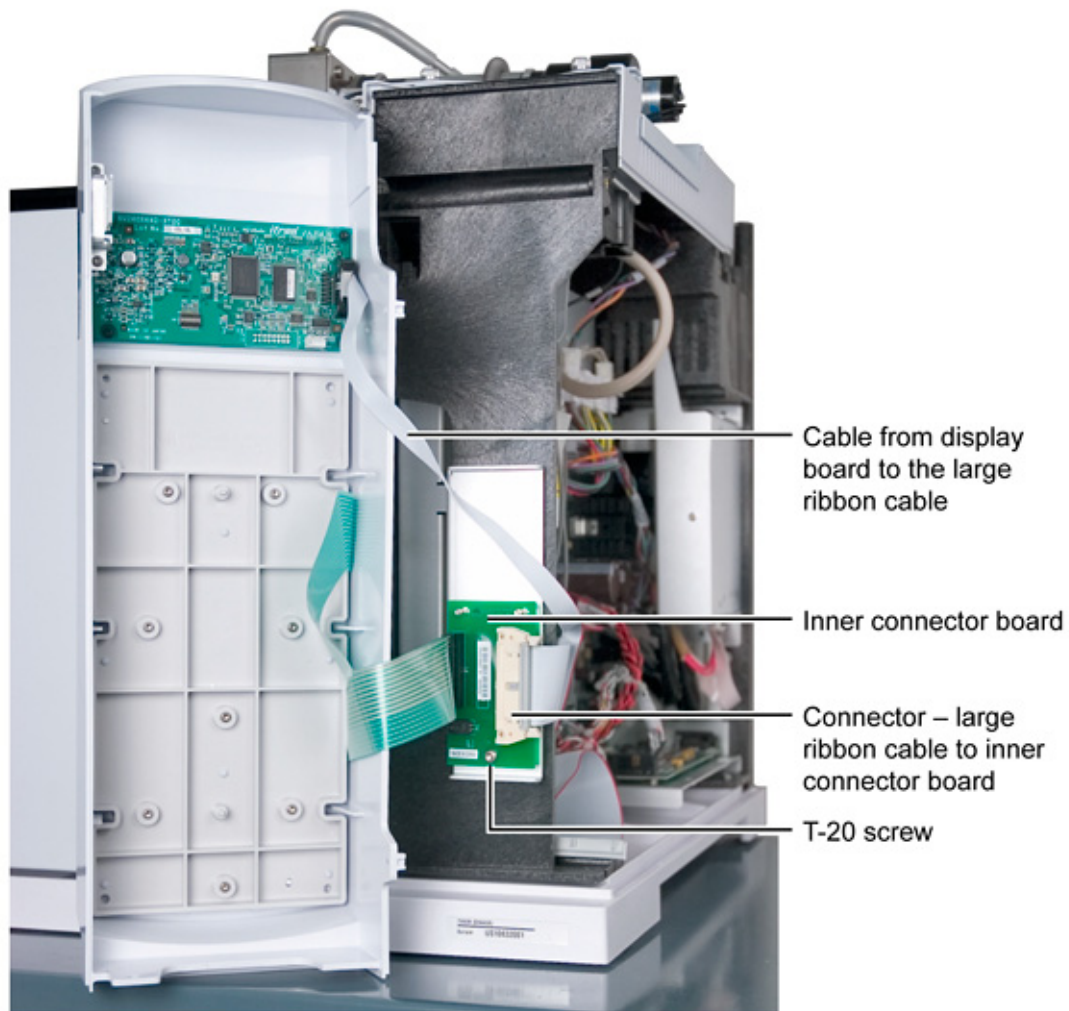

- 8 Disconnect the ribbon cable from the display board.
- 9 Disconnect the large ribbon cable from the inner connector board.
- 10 Disconnect the cable (at the lower left of the inner connector board) from the oven door safety switch.
- 11 Remove the inner connector board from the GC.

**12** Replace the entire keypad assembly with a new one.

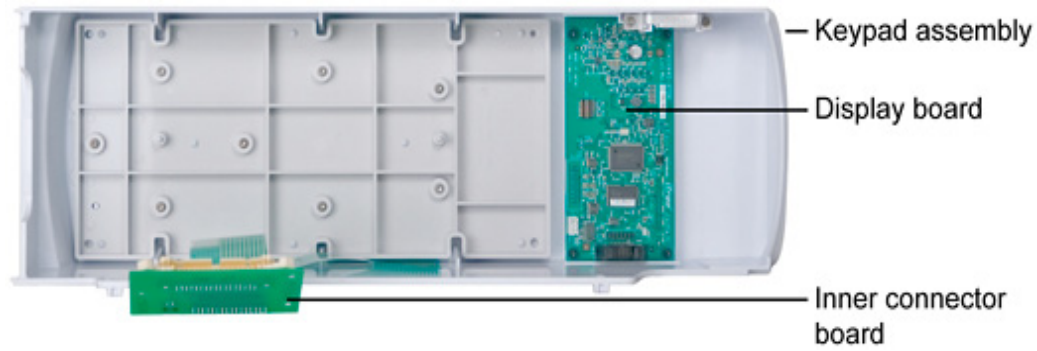

**13** Reassembly is the reverse of these steps.

## Replacing only the display board (original style)

If changing the GC configuration, see “[Changing the GC configuration](#)” for important information regarding GC methods. Then proceed with the steps below.

- 1 Cool down the oven and all heated zones.
- 2 Turn off all gas flows at their sources.
- 3 Turn off the GC main power switch and disconnect its power cord.

### WARNING

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**


---

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 4 Perform steps 4 through 8 of the “[Replacing the keypad assembly \(original style\)](#)” procedure.
- 5 Disconnect the connector cable from the display board to the inner connection board.
- 6 Remove the four screws holding the display assembly to the keypad panel.
- 7 Replace the display panel, making sure it is oriented as shown in the figure, with the connector receptacle on the right.

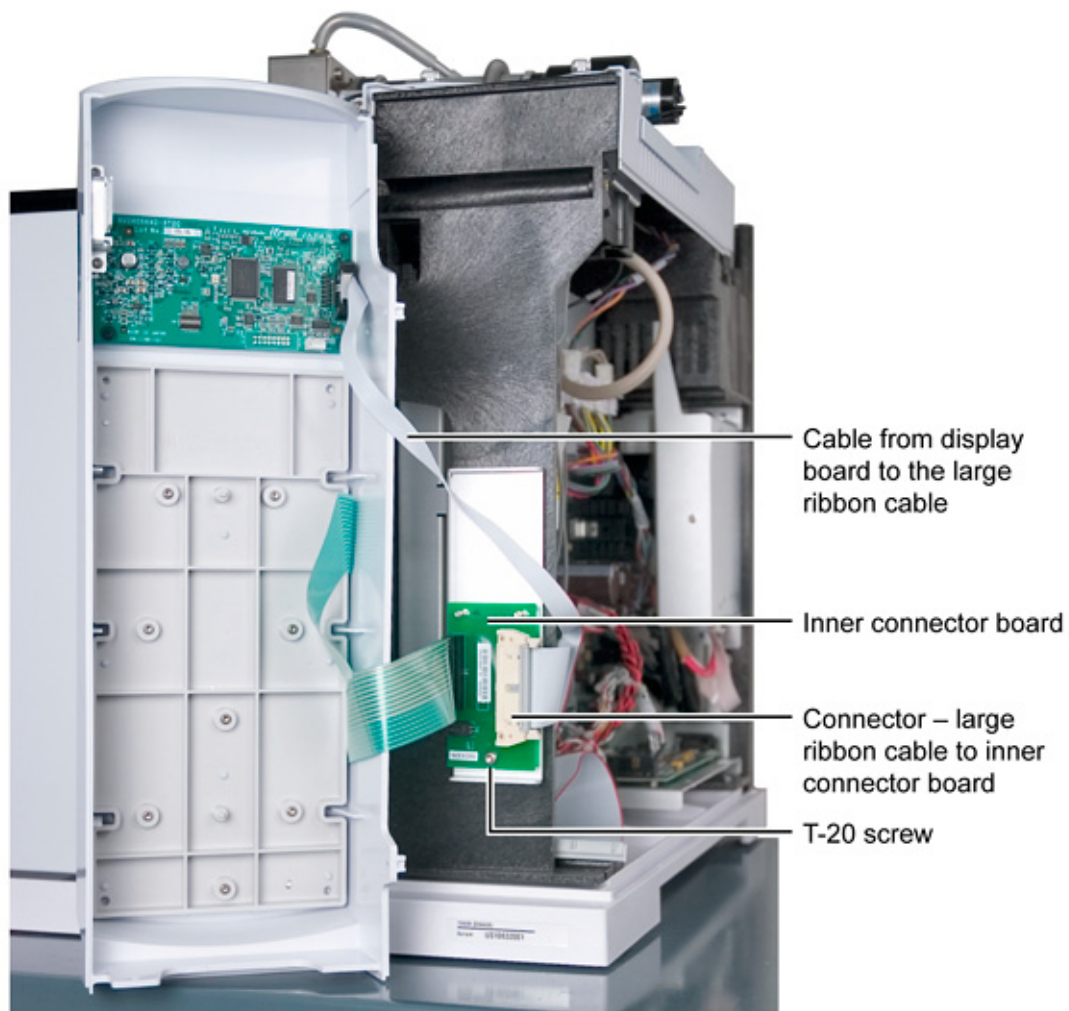

- 8 Reassemble all screws and cables in the reverse order of disassembly.

## Replacing the keypad assembly (new style)

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

### WARNING

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

---

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 3 Open the oven door to facilitate the removal of the front keypad assembly.
- 4 Raise the detector top cover.
- 5 Remove three screws with ¼-inch nut heads from the detector electronics panel side. Two of the screws are near the upper left and lower left corners of the front detector signal board; the other is on the bottom of the GC in a slotted hole.
- 6 Remove the T-20 screw from the bracket near the top of the keypad assembly to release it from the GC body.
- 7 From the front, slide the keypad panel forward.
- 8 At this point there are two servicing options: replacement of the entire keypad assembly or just the display. [See next topic for the replacement of the display only.]
- 9 Reassembly is the reverse of these steps.



## Replacing only the display board (new style)

If changing the GC configuration, see [“Changing the GC configuration”](#) for important information regarding GC methods. Then proceed with the steps below.

- 1 Cool down the oven and all heated zones.
- 2 Turn off all gas flows at their sources.
- 3 Turn off the GC main power switch and disconnect its power cord.

### WARNING

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**


---

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 4 Perform steps 4 through 8 of the [“Replacing the keypad assembly \(new style\)”](#) procedure.
- 5 Disconnect the connector cable from the display board.
- 6 Remove the four screws holding the display assembly to the keypad panel.
- 7 Replace the display panel, making sure it is oriented as shown in the figure, with the connector receptacle on the right.

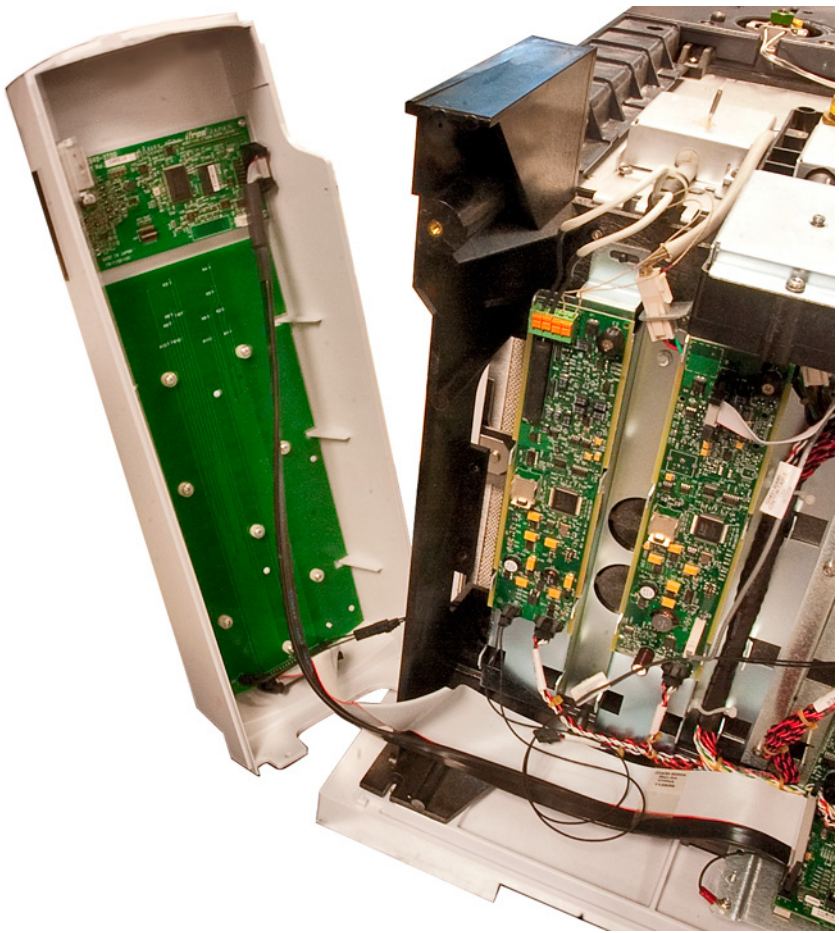

- 8 Reassemble all screws and cables in the reverse order of disassembly.

## Replacing the original keypad assembly with the new assembly

### NOTE

When replacing the original assembly with the new assembly, the only new parts needed are the new Keypad Assembly (G3430-60555) and the Door Sensor Interconnect Cable (G3430-60561). All other original parts are reusable.

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

### WARNING


**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

- 3 Remove three screws with ¼-inch nut heads from the detector electronics panel side. Two of the screws are near the upper left and lower left corners of the front detector signal board; the other is on the bottom of the GC in a slotted hole.
- 4 Disconnect the large keypad ribbon cable from the inner connector board.
- 5 Remove the T-20 screw from the bracket near the top of the keypad assembly to release it from the GC body.
- 6 From the front, slide the keypad panel forward.
- 7 Disconnect the ribbon cable from the display board.
- 8 Disconnect the large ribbon cable from the inner connector board.
- 9 Disconnect the cable (at the lower left of the inner connector board) from the oven door safety switch.
- 10 Remove the inner connector board from the GC.
- 11 Remove the inner connector board bracket from the GC.

- 12** Connect the door sensor interconnect cable to the door sensor cable (to extend the cable length).
- 13** If the new keyboard assembly does not include a display, remove the display from the old keyboard assembly and install it onto the new one.
- 14** Connect the small ribbon cable to the display board.
- 15** If using the new ribbon cable, connect the grounding nut to the GC chassis.
- 16** Connect the large ribbon cable from the the main board to the new keypad assembly ribbon cable.
- 17** Set the new keypad panel in place.
- 18** From the front, slide the keypad panel backward.
- 19** Install the T-20 screw in the bracket near the top of the keypad assembly to fasten it to the GC body.
- 20** Install three screws with 1/4-inch nut heads on the detector electronics panel side. Two of the screws are located near the upper left and lower left corners of the front detector signal board; the other is on the bottom of the GC in a slotted hole.
- 21** Lower the detector top cover, and close the oven door.

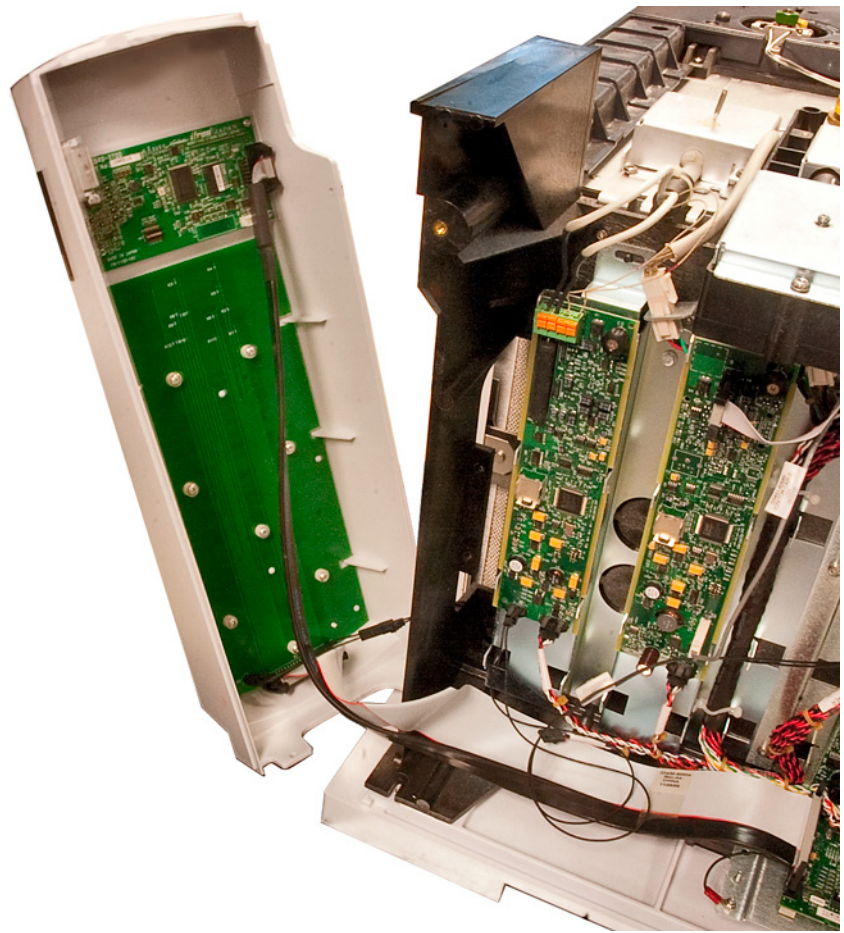

## Oven temperature control

The oven temperature is controlled by a combination of the shroud heater, the flapper assembly, and (optionally) a cryogenic cooling assembly.

### CAUTION

Temperature control and safety require a tight seal between the oven door and the oven body. DO NOT place anything (such as thermocouple leads) between these two parts. The small gap created by a wire allows hot air to escape, which can damage or melt nearby parts.

---

## Oven temperature troubleshooting

See the [7890 GC Troubleshooting manual](#) for more detailed procedures and details. The table below lists summary information for reference.

**WARNING**

**The exhaust coming from the back of the GC is very hot. Keep hands and face away from the exhaust vent.**

---

### Oven does not heat up

- 1 Press [**Status**] to check for errors to report to Agilent.
- 2 Power cycle the GC.
- 3 Check oven flapper operation.
  - a Increase oven temperature by at least 20 degrees.
  - b Verify that the oven flaps in the back of the GC are **closed**. Listen to verify that the fan is operating.
    - If the fan is on but the flap is open, the flap motor or the analog & power board is not functioning properly.
    - If the fan is off but the flap is closed, the fan motor or the analog & power board is not functioning properly.
    - If the fan is off and the flap is open, the analog & power board is probably bad.

## Oven not cooling or never cools

Check oven flapper operation.

- 1 Decrease oven temperature by at least 20 degrees.
- 2 Verify that the oven flaps in the back of the GC are **open**.  
Listen to verify that the fan is operating.
  - If the fan is on but the flap is closed, the flap motor or the analog & power board is not functioning properly.
  - If the fan is off but the flap is open, the fan motor or the analog & power board is not functioning properly.
  - If the fan is off and the flap is closed, the analog & power board is probably bad.

If using cryo cooling:

- Check for sufficient cryo coolant.
- Check if operating limits have been exceeded.
- Bad cryo valve
- Clogged cryo cooling nozzle
- Cryo cooling tank does not have proper "dip tube." See the [GC, GC/MS, and ALS Site Preparation Guide](#).
- Bad A&P board



## Testing resistance of the heater coil

If you believe that your heater coil is cracked or otherwise damaged and has caused an open circuit, you can check it by measuring its resistance.

### WARNING

**Before proceeding, turn off the main power switch and unplug the power cord.**

### CAUTION

Make sure you are properly grounded with an ESD strap before continuing.

- 1 Turn the instrument power off.
- 2 Disconnect the oven heater leads (P3, P4) from the AC power board.
- 3 Use an ohmmeter to measure resistance at the **connectors**.

Acceptable resistance ranges (in ohms) are given below. Acceptable resistances range from the nominal value for a new, cold heater to +5% from the nominal value.

### NOTE

Resistance goes up approximately 3% after heating the coil.

Nominal cold heater resistances, ohms	Standard oven	Fast-ramp oven
120 V	9.07–9.52	n/a
200 V	n/a	17.78–18.7
220 V	n/a	21.51–22.6
230 V	33.06–34.71	23.51–24.7
240 V	n/a	25.60–26.9
n/a = not available		

## Cryo valve installation/replacement

A cryogenic valve allows liquid nitrogen or CO<sub>2</sub> to be dispersed into the double-walled plenum of the GC oven where the fan blows the vapors into the oven itself. Use the procedures that follow to install a new cryogenic valve or to replace an existing cryogenic valve or nozzle.

### Installing a new cryo valve

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

#### **WARNING**

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

---

#### **CAUTION**

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 3 Remove the left side cover on the instrument and remove the plug for the cryo valve.

If a side mount carrier is installed on the left side of the instrument, remove its left cover by removing the two bottom screws, sliding the panel towards the back of the instrument and lifting it off.

- 4 Remove the knockout on the left side of the instrument. Use a screwdriver to pry it out.
- 5 Insert the cryogenic valve probe through the insulation into the oven plenum.
- 6 Screw the valve to the side of the GC oven using the two Torx T-20 screws provided.
- 7 If you are installing COC cryo blast for the front and/or rear inlet(s), connect it at this time. See [“Cool On-Column Inlet”](#) for details.
- 8 Plumb the valve to the liquid nitrogen or carbon dioxide source.

- 9 Plug the cryo valve's wiring harness into the five pin connector to the left of the valve.
- 10 Route the cryo tubing through the cutout in the GC side cover and reinstall the cover.

If a side mount carrier is installed, route the cryo tubing out of the slot in the back of the carrier.

## Replacing an existing cryo valve

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

**WARNING**

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

---

**CAUTION**

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 3 Remove the left side cover on the instrument.


If a side mount carrier is installed, remove its left cover by removing the two bottom screws, sliding the panel towards the back of the instrument and lifting it off.

- 4 Shut off the cryo fluid supply and crack the valve at the supply to release any residual pressure.
- 5 Use a 9/16-inch wrench to disconnect the cryogenic fluid supply tube at the cryo valve.
- 6 Disconnect the cryo valve's wiring harness from the connector to the left of the valve.
- 7 Remove the two Torx T-20 screws holding the old cryo valve and pull it straight out from the instrument.
- 8 To remove the valve from the bracket, disconnect any cryo blast plumbing from the valve at the Swagelok tee. Remove the two Torx T-20 screws on the bottom of the bracket.
- 9 To replace the cryogenic nozzle, proceed as follows:
  - a Use a 9/16-inch wrench to unscrew the old nozzle. Discard the nozzle.
  - b Wrap the threads of the new nozzle with PTFE tape, being careful not to cover the first two threads of the nozzle.
  - c Screw on the new nozzle and tighten firmly with a 9/16-inch wrench.
- 10 Insert the cryogenic valve probe through the insulation into the oven plenum.


- 11 Screw the valve assembly back on to the side of the GC oven using the two Torx T-20 screws.
- 12 If you are installing COC cryo blast for the front and/or rear inlet(s), connect it at this time. See [“Cool On-Column Inlet”](#) for details.
- 13 Plumb the valve to the liquid nitrogen or carbon dioxide source.
- 14 Plug the cryo valve's wiring harness into the connector to the left of the valve.
- 15 Route the cryo tubing through the cutout in the GC's side cover and reinstall the cover.

If a side mount carrier is installed, route the cryo tubing out of the slot in the back of the carrier.

**Removing the valve from the valve bracket (N2 valve shown)**

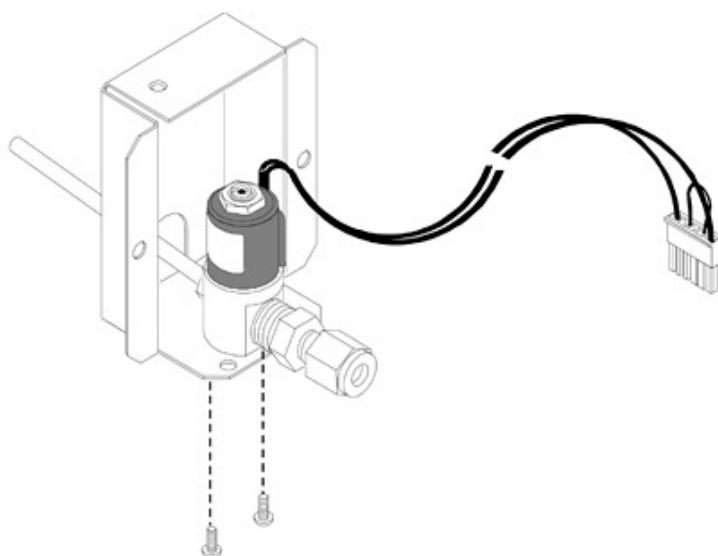

**Cryo blast attachments**

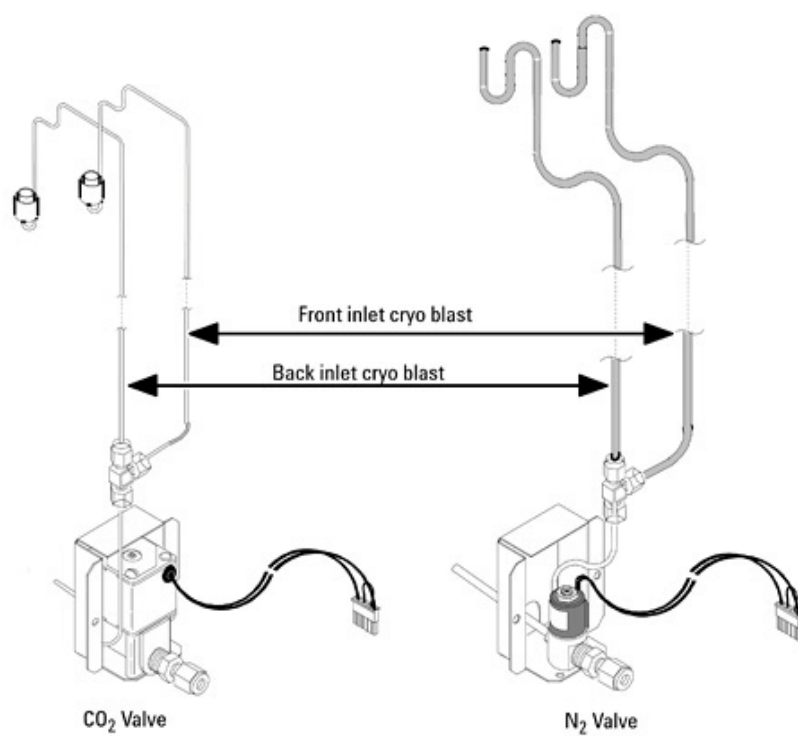

## Replacing the oven sensor

The oven temperature sensor can be replaced after removing the oven shroud.

**WARNING**

**Before proceeding, turn off the oven and let it and any heated zones cool. Turn off the main power switch and unplug the power cord.**

---

- 1 Remove the oven shroud as described in Replacing the oven shroud assembly.
- 2 Loosen the two Torx T-20 screws securing the sensor retainer to the back of the shroud.
- 3 Slide the old sensor out of the retainer. Thread a new sensor through the opening in the rear of the oven.
- 4 Slide the new sensor under the retainer. One of the grill holes on the front of the shroud is stamped. Make sure that the end of the sensor is positioned behind the stamped hole before tightening the retainer.

## Replacing the oven door sensor

The oven door sensor installs above the oven door button, and connects to a cable that runs to the GC keyboard and display.

**WARNING**

**Before proceeding, turn off the oven and let it and any heated zones cool. Turn off the main power switch and unplug the power cord.**

---


To remove the old sensor:

- 1 Prepare the GC. To replace the over door sensor, you will need to tilt the GC slightly. You may need to vent and disconnect any MS devices or otherwise move instruments that block access or otherwise prevent you from lifting the right side of the GC enough to reach underneath it.
- 2 Remove the display and keyboard assembly. See [“Replacing the keypad assembly and display board”](#). (As part of this step, you will disconnect the door sensor cable from the display.)
- 3 Pull the oven door sensor cable back through the electronics carrier.
- 4 Remove the 1/4-inch socket head screw that secures the oven door sensor to the GC.
- 5 Pull the oven door sensor cable through and remove the sensor from the GC.




To install the new sensor:

- 1 Feed the oven door sensor cable down through the hole behind where the sensor mounts.

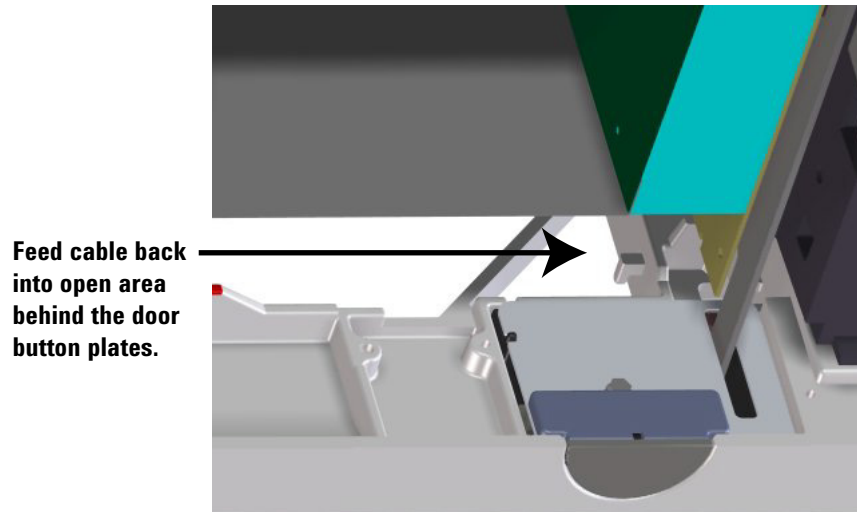


- 2 With one hand, tilt the GC (raise the right side) enough to be able to reach under it and grab the sensor cable. (Prop up if needed.)
- 3 While holding the GC up, pull the cable down and through the opening, then feed it back up and into the space behind the oven door latch plate. Gently let the GC down.

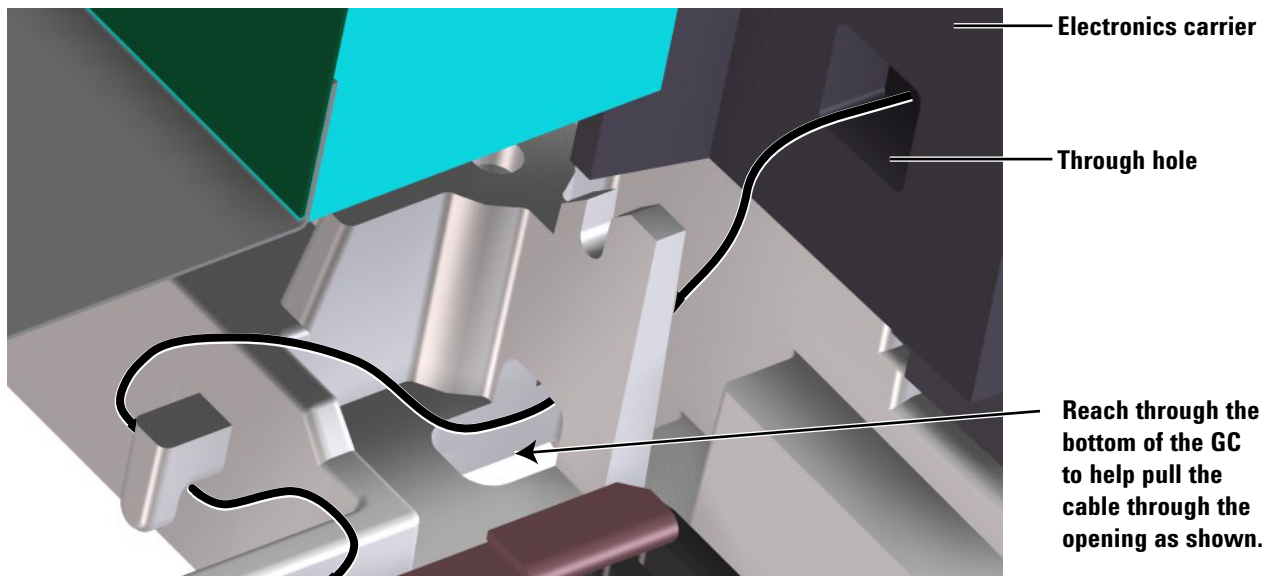


- 4 Use tweezers or similar to pull the cable up and then feed the cable through the opening in the electronics carrier.
- 5 Install the oven door sensor using the 1/4-inch socket head screw. Keep the sensor face square (flat) to the oven door.

The face of the sensor should be parallel to the outer edge of the GC frame.

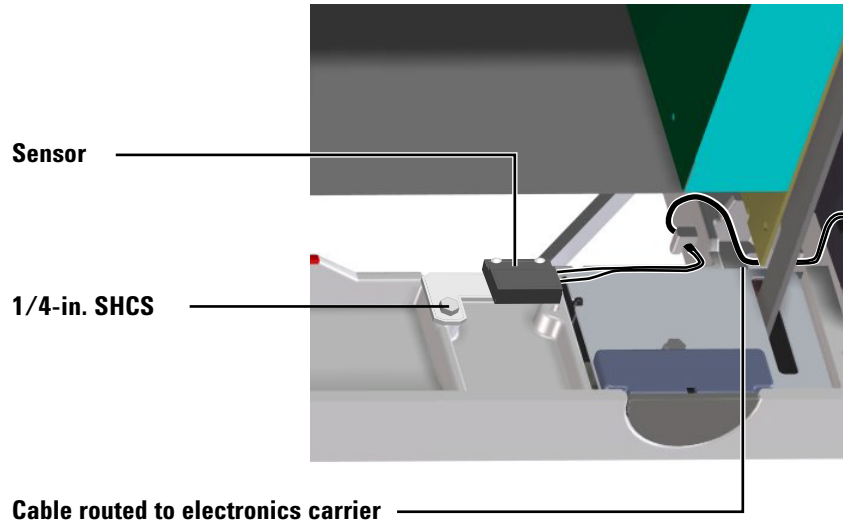

- 6 With one hand, tilt the GC (raise the right side) and use your other hand to set the cable wires in the retaining clip on the GC frame. Make sure the cable wires route as shown in the figures. Proper routing prevents the GC frame and oven door button from clipping the cable.

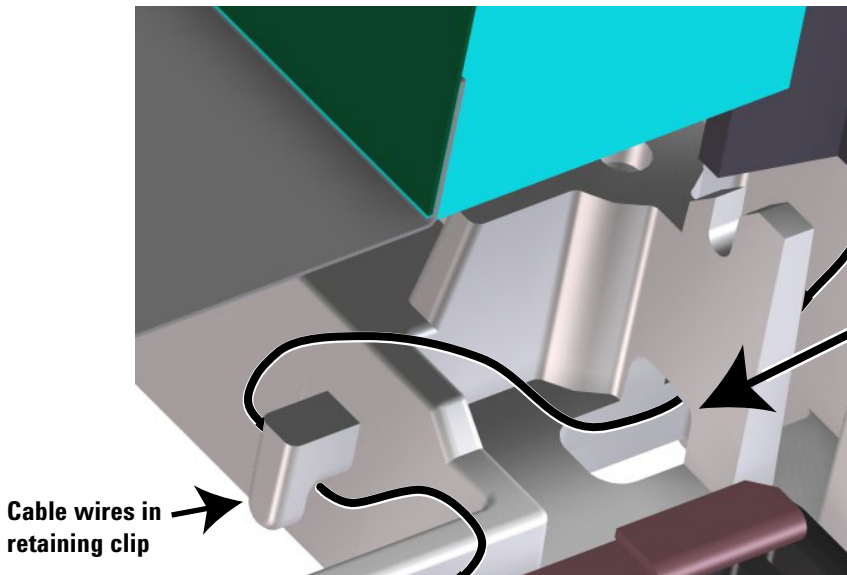

- 7 Reinstall the keyboard and display assembly. Check oven door operation.
- 8 Turn on the GC and verify that the door sensor functions properly.
- 9 Return the system to operating conditions.

## Replacing the oven shroud, oven fan, and oven fan motor

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

**WARNING**

Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.


---

**CAUTION**

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.


---

- 3 At the front of the GC, remove all obstructing items from inside the oven (column(s), column hanger(s), inlet/detector insulation cups, and so forth).

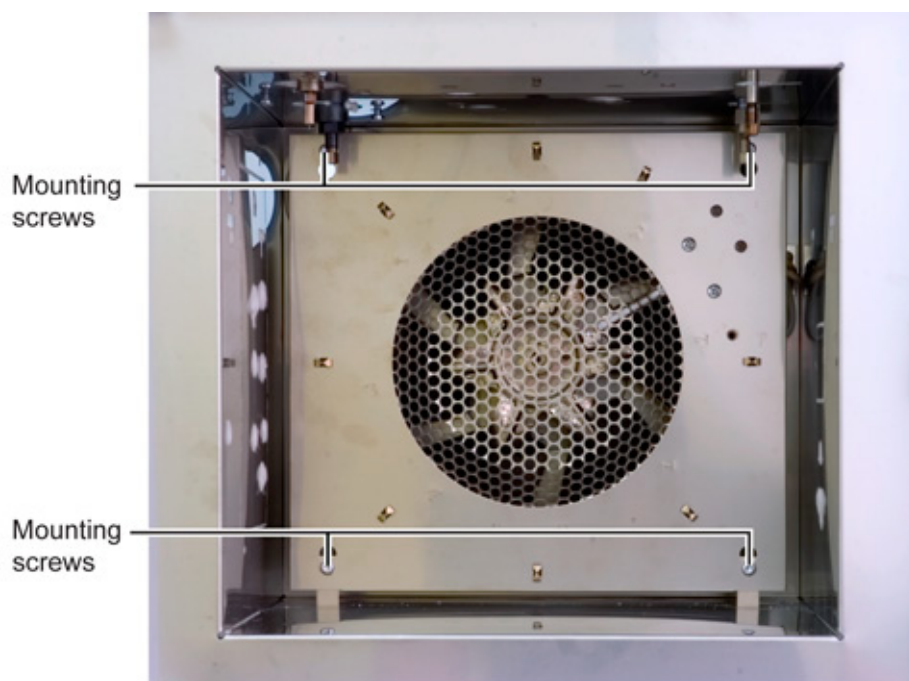

- 4 At the rear of the GC, disconnect the oven sensor cable at the Analog and Power board and feed the cable back to the rear

of the GC. Also disconnect the two oven heater leads at the AC board.

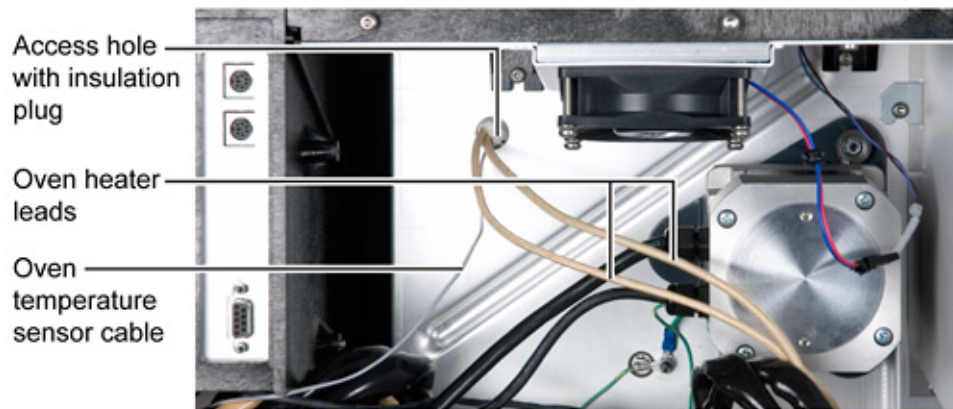

- 5 From the oven, remove the two lower screws securing the oven shroud and loosen the upper two. These are special stainless steel screws: the ones removed must not be confused with any other screws.
- 6 Lift the shroud over the upper screws and work it out of the oven.


**NOTE**

Connectors for the oven heater wires and for the oven sensor must be carefully worked through the access hole at the rear of the oven. There is also a 'plug' of insulation which may fall from the access hole during this process: it must be replaced in reassembly from the rear of the GC.


- 7 Inspect the oven heater and its supporting circular quartz glass rod: they must be undamaged in any way.

**NOTE**

Be very careful to not disturb the oven temperature sensor: it is fragile and easily damaged. Also, its physical location is critical to proper oven temperature control: do not alter its position in any way.

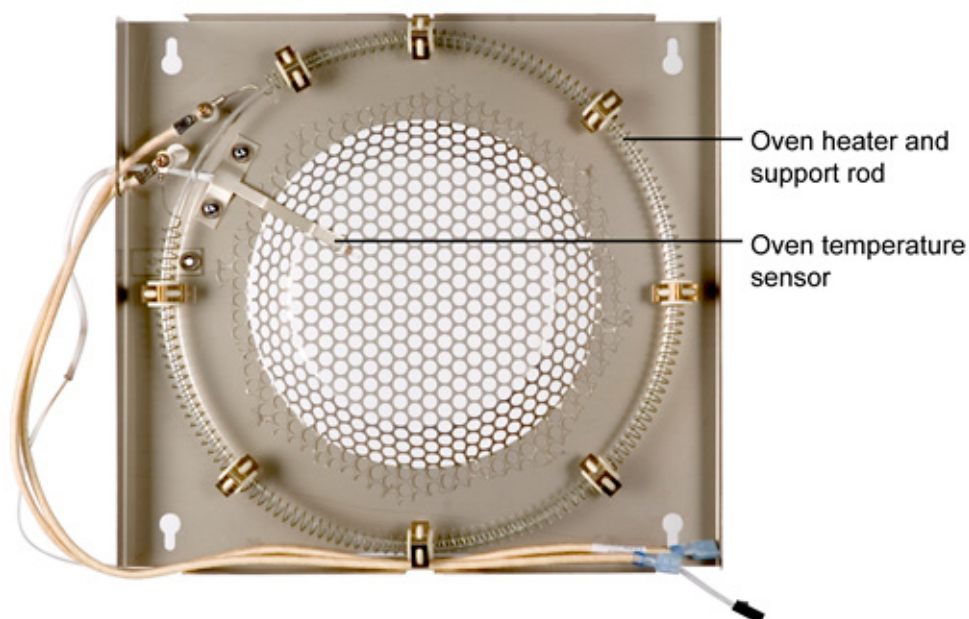


- 8** In replacement, insert oven heater leads first through the access hole, followed by the oven sensor cable.

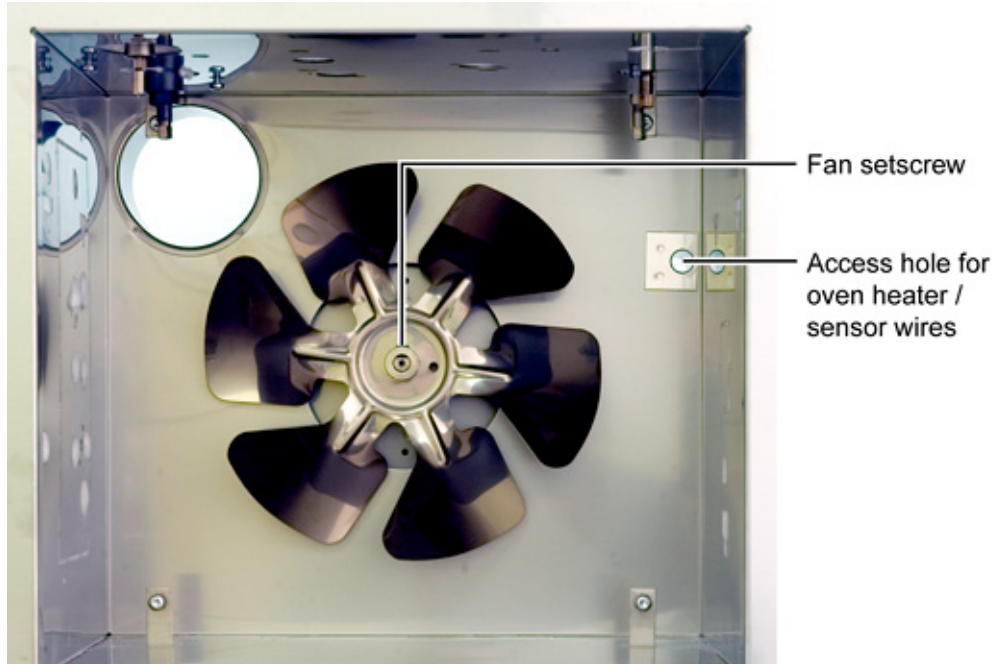

- 9** Before removing the oven fan, note its hub's physical location on the fan motor shaft so as to return it to the same position in reassembly.
- 10** Being careful to not deform its blades, remove the fan by loosening the 5/32-inch Allen setscrew on its hub and then pulling the fan from the motor shaft.



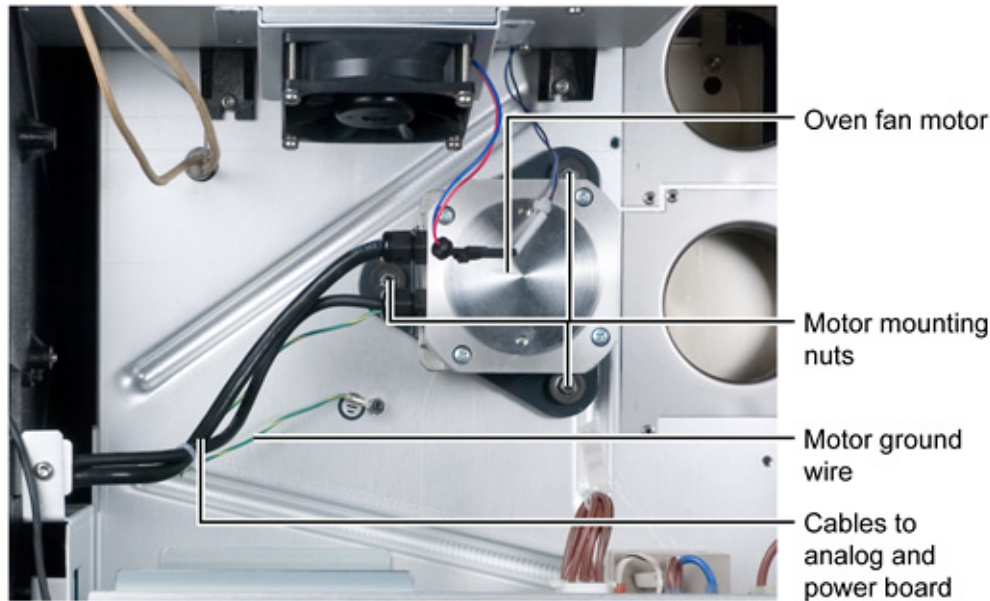

- 11 At the rear of the GC, disconnect the fan motor ground lead. Then disconnect the cable to the Analog and Power board and work the cable back to the rear of the GC.

### WARNING

**The oven fan motor is heavy! Make sure you have a firm grip on the motor before removing / replacing it. Always handle the motor by the body. Do not carry the motor by the shaft.**

- 12 Remove three 7-mm nuts with lockwashers and cone washers to release the motor. Remove the top nut last to maintain mechanical support for the motor. Manually support the motor as the top nut is removed.
- 13 Reassembly is the reverse of these steps with the following considerations:
- Each motor mount must have both an inner and outer cone washer
  - Orient the motor such that its power cable exits towards the electronic side of the GC
  - Reinstall the top nut first to support weight of the motor
  - After motor reassembly into the GC, check that its shaft extends perpendicularly from the back oven wall. If not perpendicular, there may be mounting part(s) missing or misplaced, and/or mounting nut(s) untightened.

- In fan replacement, make sure the fan's hub is returned to its original location, or slightly rearward from being exactly flush with the end of the motor shaft. Also make sure the setscrew is tightened against the flat side of the motor shaft.
- For the oven heater/sensor shroud, if not already done, inspect the oven heater and its supporting circular quartz glass rod: they must be undamaged in any way
- After the oven heater/sensor shroud is reinstalled, remember to return the plug of insulation in the access hole at the back of the oven.




## Replacing the oven flapper assembly


**CAUTION**

Components can be damaged by electrostatic discharge (ESD). Be sure to wear a grounded ESD strap while performing this procedure.


- 1 Remove rear GC cover.
- 2 Remove left side GC cover.
- 3 Disconnect flapper motor cable and pull it back to the motor.

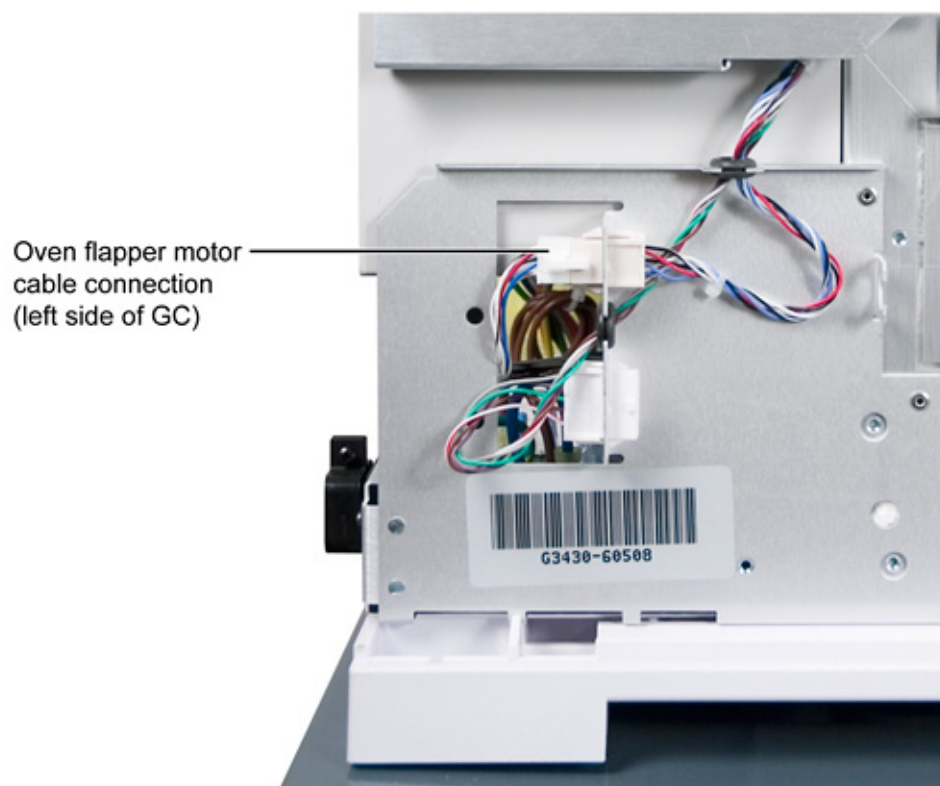

- 4 If present, remove the PTV N2 cryo feed line cover by removing two screws. If necessary, also remove the bracket behind the feed line cover.

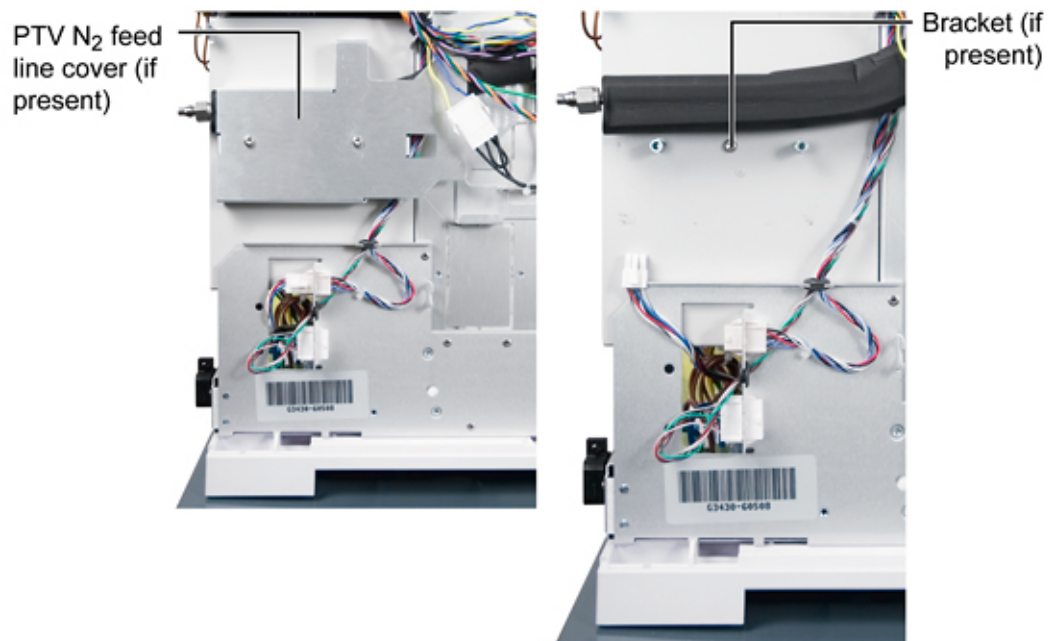

- 5 Loosen three captured T-20 screws to remove the flapper assembly. To support the weight of the assembly, the top-most screw should be removed last.
- 6 Pull the assembly straight out.

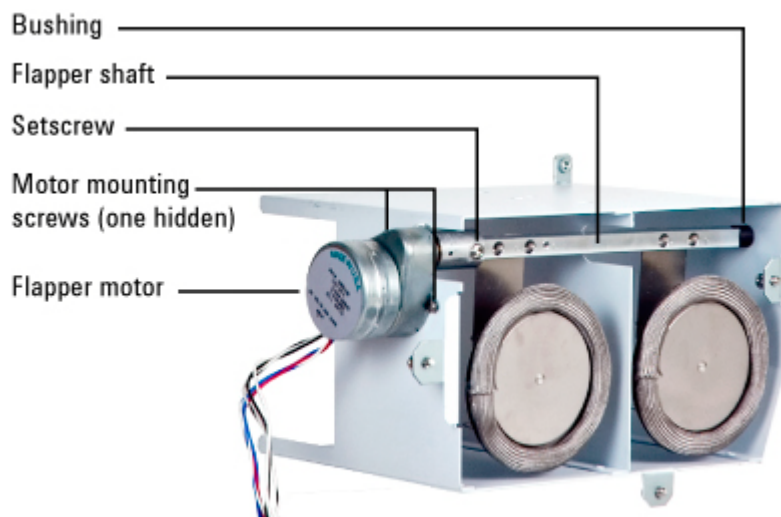


- 7 In replacing the assembly, steps are the reverse of removal steps.

Install the top-most screw first to support the weight of the assembly


## Replacing the flapper or flapper motor

This disassembly procedure enables you to replace either the flapper or flapper motor.

- 1 Remove the flapper assembly. See “[Replacing the oven flapper assembly](#)”.
- 2 Remove the T-20 set screw on the coupling between motor and flapper shafts.

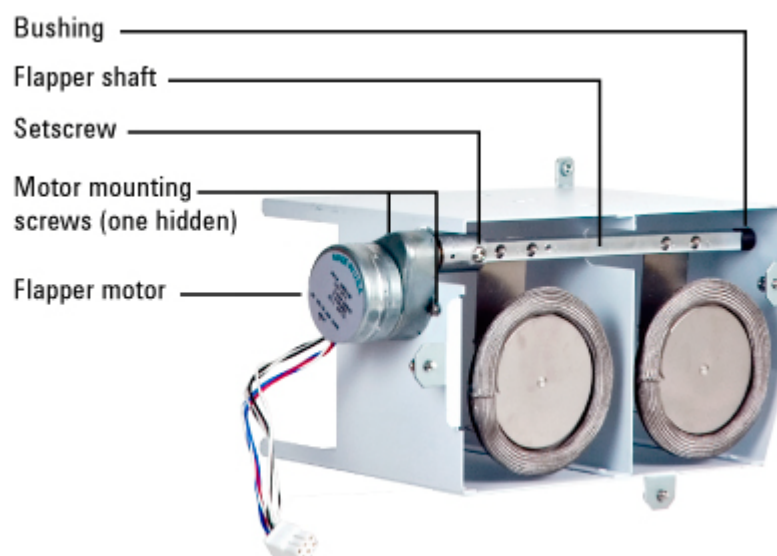


- 3 Remove three T-10 screws securing the motor to the frame.

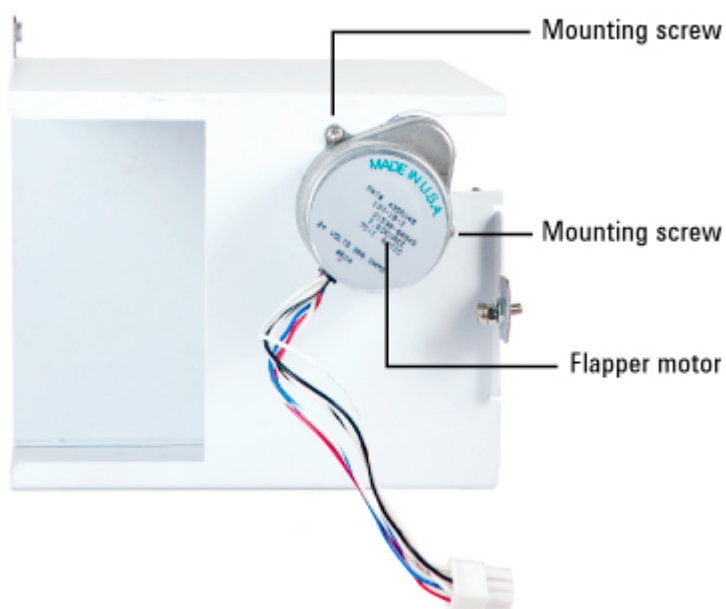

- 4 Slide the motor off the flapper shaft and replace it with a new motor.

If you are replacing the flapper, slide the flapper out of the duct assembly.

Be careful to not lose the bushing at the upper end of the flapper shaft.

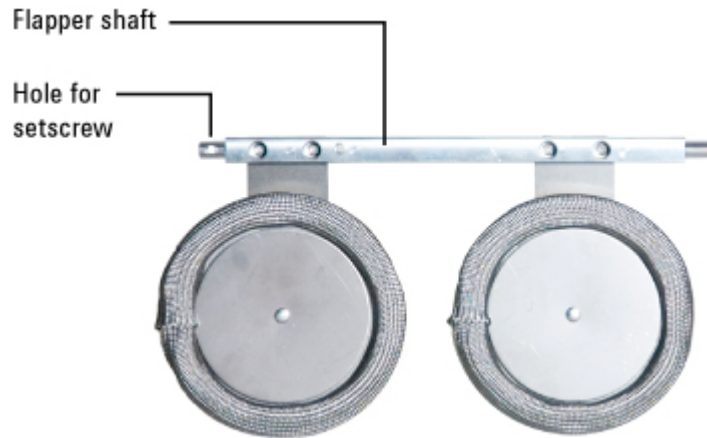

In reassembly, do the following:

- Insure the upper part of the flapper shaft is inserted into its bushing at the top of the duct before securing motor mounting screws.
- Line up the hole in the flapper shaft with the set screw on the coupler mounted on the motor shaft.

Verify the flapper operation by:


- Observing the flapper behavior at turn on. The flapper self-adjusts by closing.
- Set the oven temperature to 70C. The flapper should close.
- Set the oven temperature to 45C. The flapper should go to full open.

## Replacing Components Inside the Electronics Panel

### Accessing the analog and power board

To access the analog and power board in the rear of the electronics carrier, you must first remove:

- The ALS controller board
- The valve bracket
- The logic board

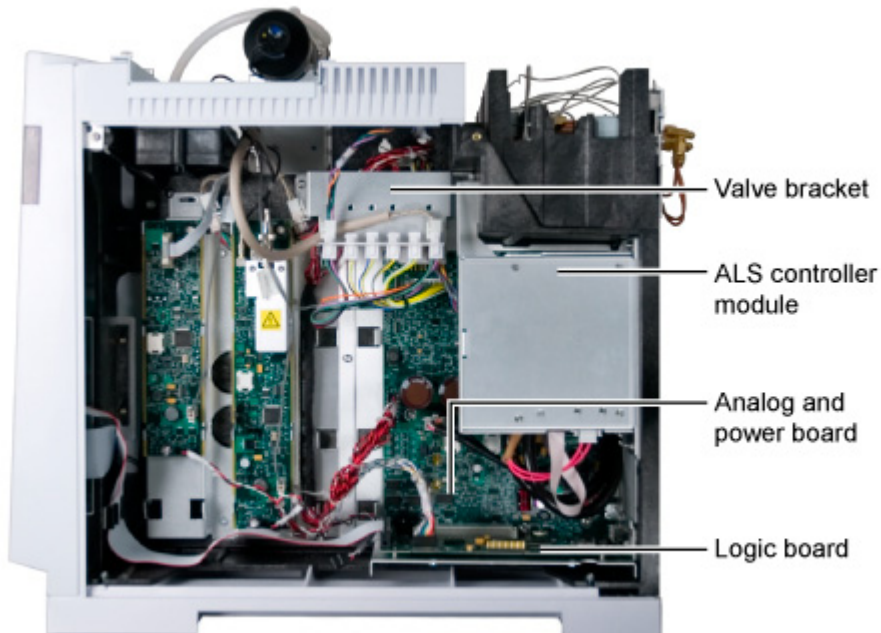

## Replacing the ALS board

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

### WARNING

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**


---

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---


- 3 Remove the right side panel.
- 4 Remove 2 T-20 screws holding the assembly to the back of the electronics carrier.

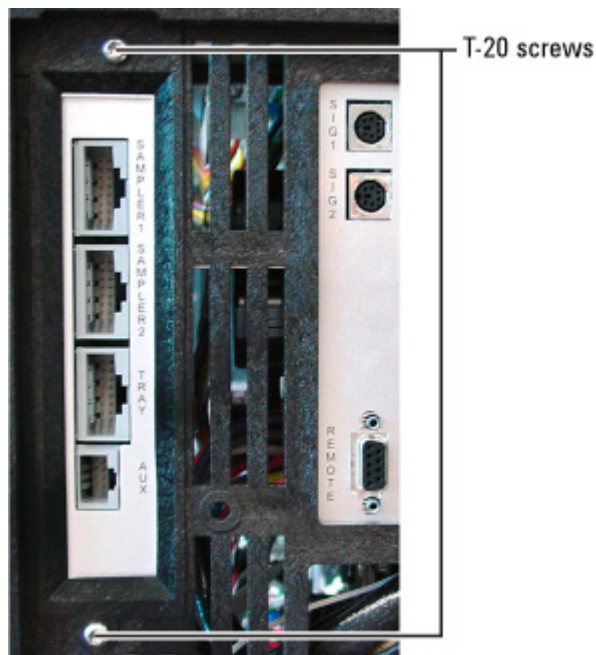


- 5 Slide out the ALS assembly.
- 6 Pinch and pull the 2-wire power connection at the bottom of the assembly and set the assembly aside.



- 7 Carefully pinch and disconnect the flex cable from J8 on the Analog & Power board.

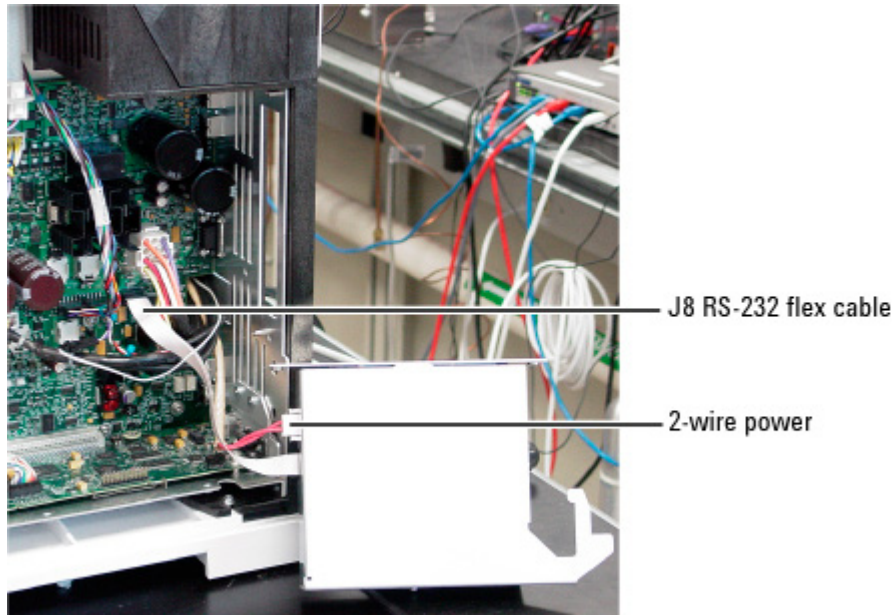

- 8 Pull the ALS module free from the back panel.
- 9 Connect the flex cable of the new ALS module to the Analog & Power board.
- 10 Slide the ALS module hook into the 4th slot.

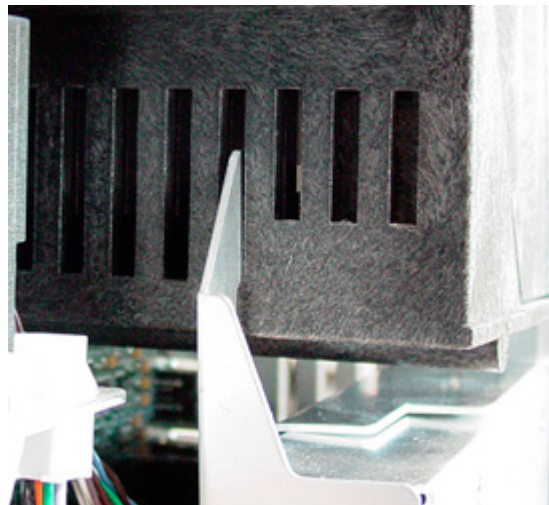

- 11 Connect the 2-wire power cable.
- 12 Secure the assembly with 2 screws.
- 13 Plug in the ALS modules and power on the GC.

- 14 Press **[Service Mode]** and select **Diagnostics, ALS Status**. Verify that the ALS controller model number and firmware revision appear on the display.
- 15 Update the ALS controller board firmware if necessary.

For more information, see the “[Electronics carrier](#)” IPB.



## Replacing the fuse on the ALS controller board

A small fuse on the ALS controller board protects it from over-currents originating in the ALS components (Injector, for example), or shorts in their cables. The fuse is associated with the 42 VAC circuit. You can easily replace this fuse.

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

### WARNING

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**


---

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 3 Gather ESD protection, T-10 driver, and tweezer.
- 4 Remove the ALS controller board. See [“Replacing the ALS board”](#).
- 5 Unscrew the T-10 screw and remove the cover.
- 6 Using the tweezer, carefully replace the fuse.

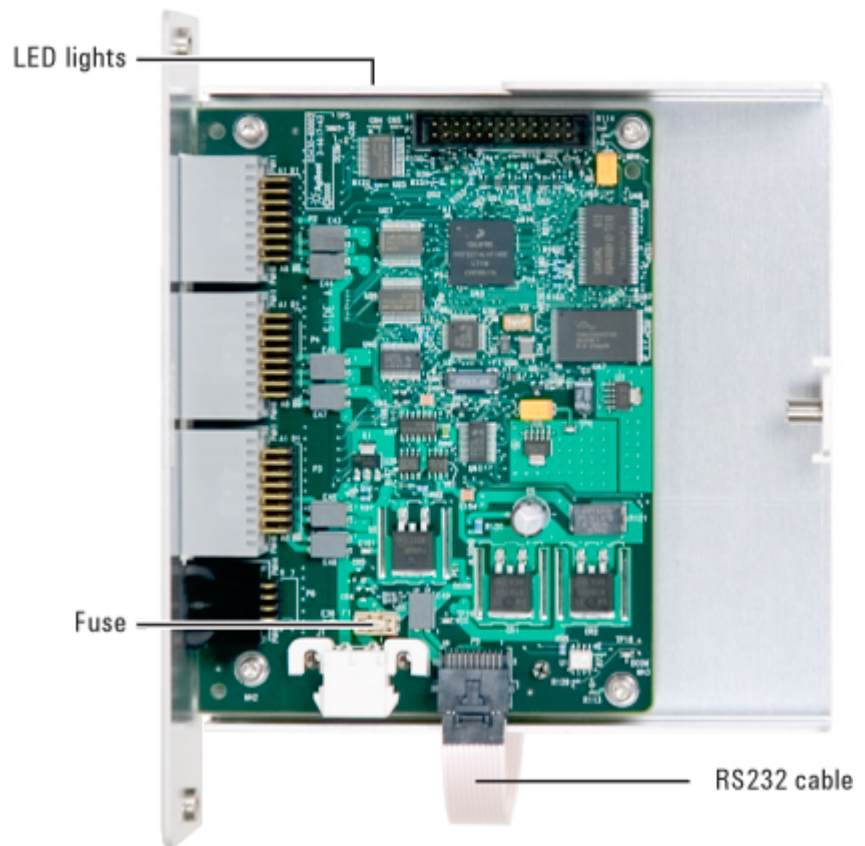

- 7 Reinstall the cover. Insure that the connectors are visible.
- 8 Reinstall the ALS controller board.
- 9 Before plugging in any ALS components (for example injectors), power on the GC. Check Service Mode | Diagnostics | ALS status to make sure the model number and firmware revision of the board are displayed.
- 10 If appropriate, update the firmware to the ALS controller.

### Troubleshooting

- If the ALS Controller is not installed, check the cable connections. Consider replacing the assembly.
- If the firmware revision is blank, the ALS controller board may not be programmed correctly.

## Removing the valve bracket

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.


### WARNING

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

- 3 Remove the GC right side cover.
- 4 Remove the ALS board. See [“Replacing the ALS board”](#).
- 5 Remove all cables connected to the upper side of the valve bracket (the illustration shows 2).
- 6 Loosen the screws at the top left and top right corners of the bracket. Remove the bracket.
- 7 Disconnect the wiring harness from the Analog and Power board at P22.

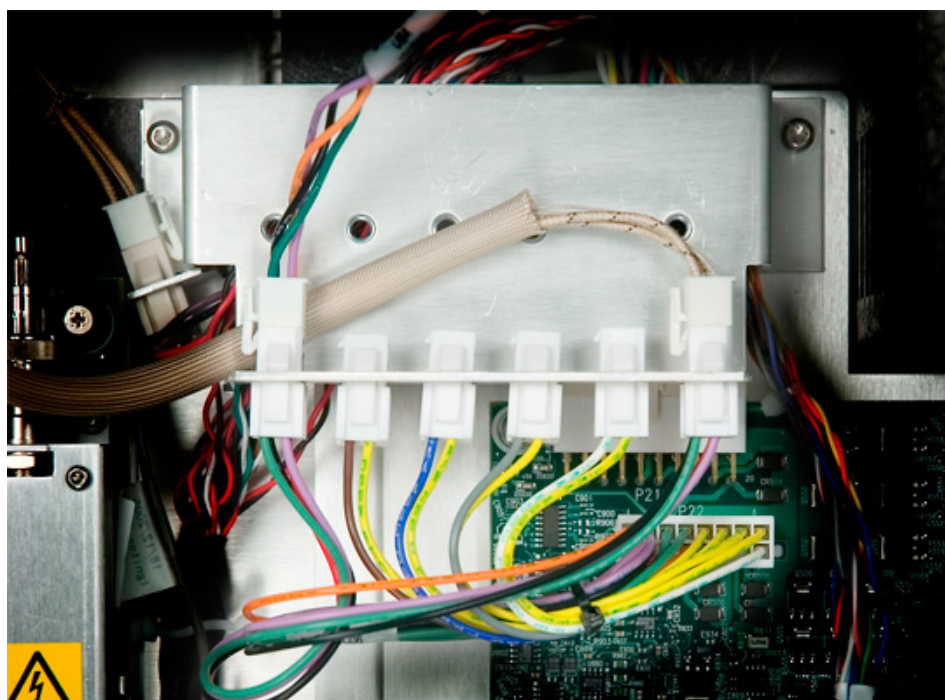

## Replacing the logic board

The logic board contains information specific to the GC and its configuration, including methods, sequences, serial number, date of manufacture, the Gateway, Subnet Mask, and IP address of the LAN, logs, oven type, AUX heater configuration(s), and so forth. These items must be reentered after the new logic board is installed.

If the GC is functional, use the keypad to view and record the following information:

- Configuration-required
  - Gas types
  - Cryogenic cooling, if present
  - Heater assignments: valve box, MSD transfer line, others
  - Valve types
  - If the configuration includes non-Agilent components or specials, you may need to re-enter custom heater or pneumatic PIDs. View and record this information.
- 1 Cool down the oven.
  - 2 Turn off the GC main power switch and disconnect its power cord.

---

**WARNING**

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

---


**CAUTION**

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---


- 3 Remove the right side electronics panel cover.

- 4 Disconnect the ribbon cable to the keyboard/display.

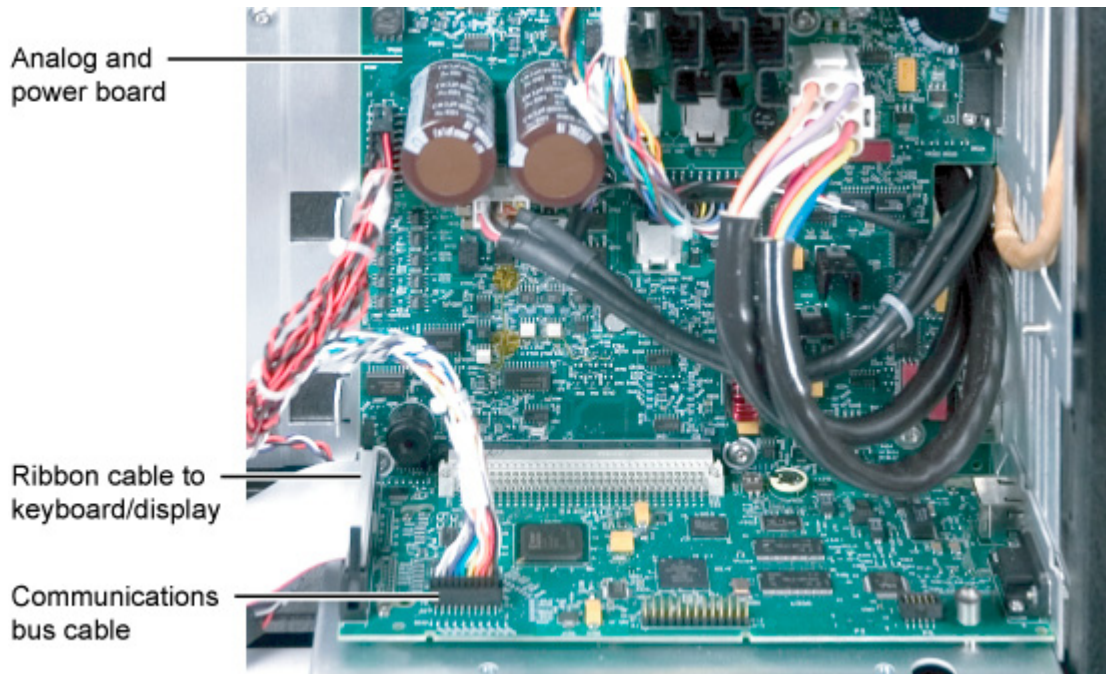


- 5 Disconnect the communications bus cable and cable to the keypad and display. If present, also remove the external LAN cable.
- 6 For the G3430-61010, G3430-61020, or G3430-610101 logic and communication board, loosen the thumbscrew in the bottom right corner of the board. If you have a G3430-61000 logic and communication board remove the two screws that secure the logic board to the rear panel.

G3430-610101 (G3430-61010 and G3430-61020 are similar):

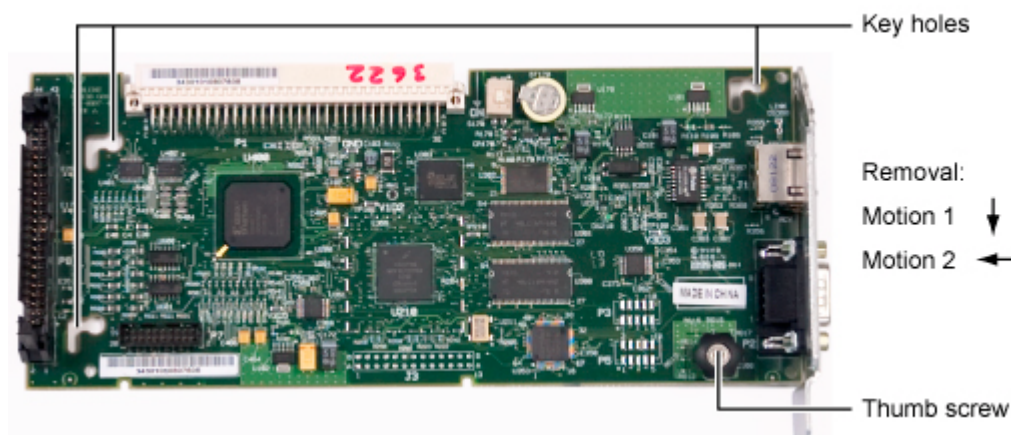




G3430-60100:

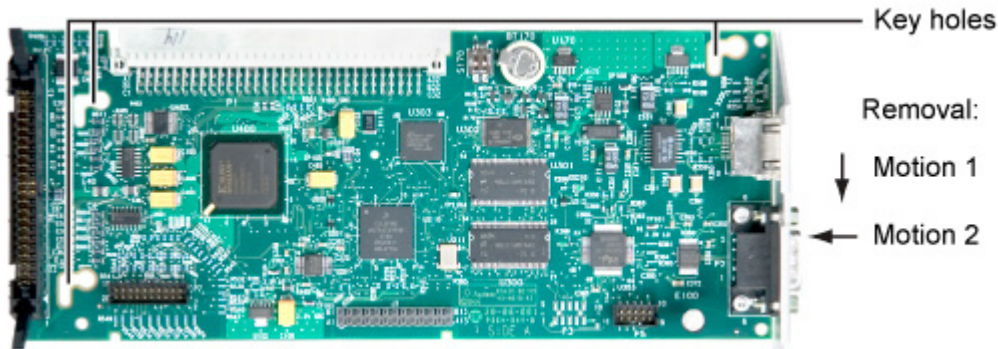


- 7 Slide the logic board toward you while using a rocking motion to unplug the board from its analog and power board receptacle.
- 8 Once free of the analog and power board, slide the logic board to the left and lift it out.
- 9 For 7890B, main board G3430-61010 does not use a battery backup. Skip to the next step.

For 7890A, verify the battery backup is ON. With respect to the new logic boards (G3430-60120 and G3430-60101), verify battery backup switch is ON (change if needed). For the G3430-60100 board, switch the two battery-connect switches found next to the battery: BOTH switches must be in positions away from the analog and power board (see above).

G3430-60101:

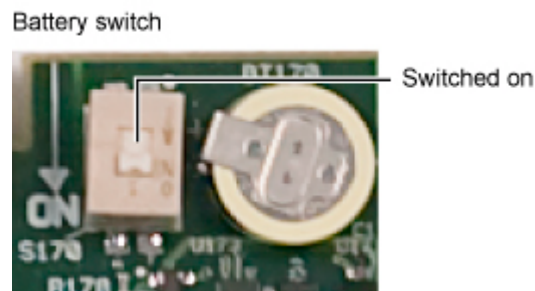

G3430-60120 and G3430-60100:

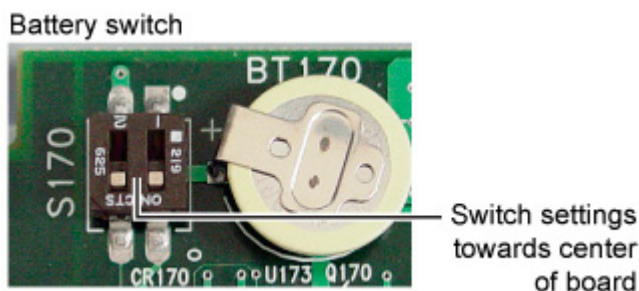

#### NOTE

If switches are incorrectly-set (battery disconnected), critical internal GC information will be lost any time AC power is not present.

- 10 Reinstallation is the reversal of removal steps.
- 11 Restore power to the GC.
- 12 The zero offsets for the flow and pressure sensors are stored on the logic board rather than in the EPC modules. You must now zero all sensors. Navigate to Service Mode | Diagnostics | Electronics | Pneumatics to do this.
- 13 You must now restore information stored on the logic board specific to the GC and its configuration.
- 14 Press Service mode. Scroll to Diagnostics and press [**Enter**].
- 15 Scroll to Instrument Status. Press [**Enter**].
- 16 Press [.] [.] [**Mode/Type**], select the country, and press [**Enter**].
- 17 Press [.] [.] and enter the 8-digit serial number; press [**Enter**].
- 18 Scroll to Mfr date. Press [.] [.] and enter 6 digits in ddmmyy format. Press [**Enter**].
- 19 You have changed the configuration of the GC. See [“Changing the GC configuration”](#) for important information regarding GC methods.

#### Check GC configuration

- Gas types
- Cryogenic cooling
- Heater assignments: valve box, MSD transfer line, others
- Valve types
- If the configuration includes non-Agilent components or specials, you may need to re-enter custom heater or pneumatic PIDs.

## Replacing the analog and power board

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

### WARNING

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

---


### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

---

- 3 Remove the right side electronics panel cover.  
If present, remove the ALS controller module. See [“Replacing the ALS board”](#).
- 4 Remove the logic board. See [“Replacing the logic board”](#).
- 5 If present, disconnect valve box heater/sensor(s) from the valve bracket and remove the bracket. See [“Removing the valve bracket”](#).
- 6 Disconnect cables at all board receptacles. G3430-60150 Analog and power board shown below (G3430-60151 is similar).




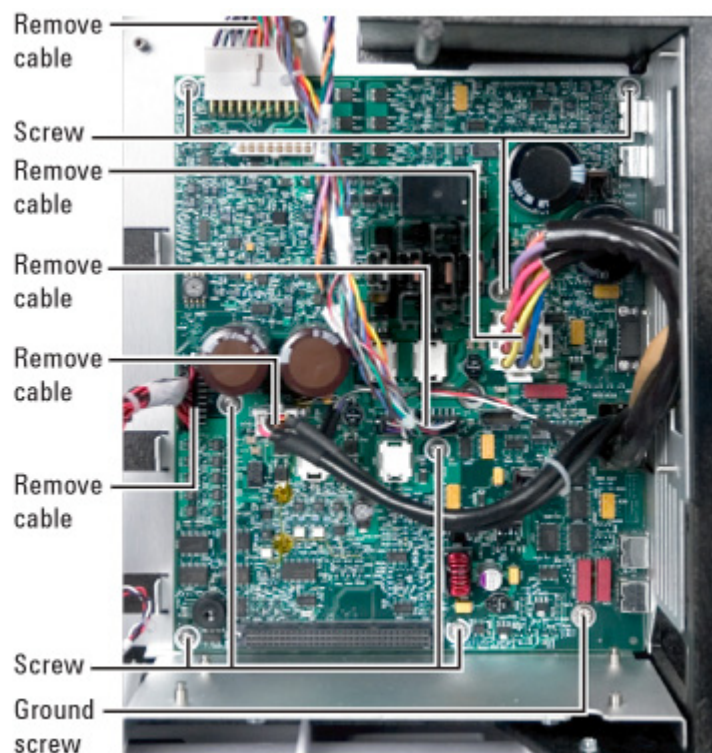

**7** Remove 8 T-20 screws:


- One is a long grounding screw found at the lower right corner of the board
- One screw is difficult to locate next to and above the transformer receptacle near the right center of the board
- In all cases, locate the specific screw first before you attempt to remove it: a misplaced screwdriver may damage nearby board components.

**8** Slide the board slightly outward and to its left to release it from the GC. G3430-60151 Analog and power board shown below (G3430-60150 is similar).

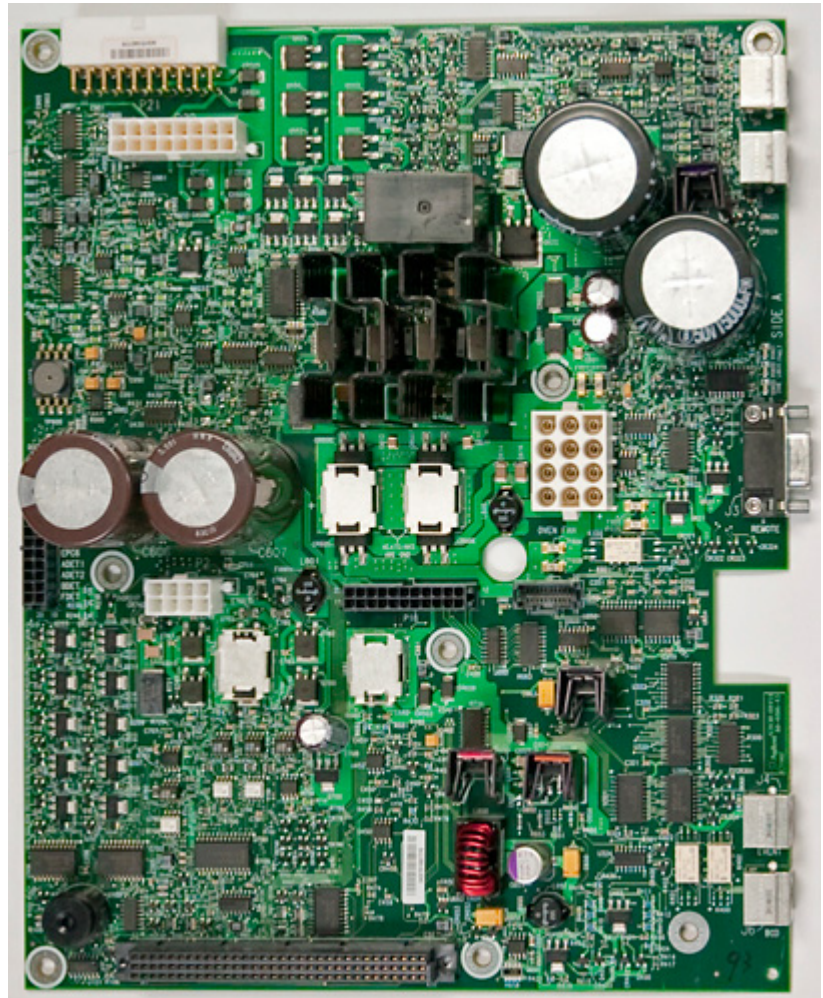

9 Installation of the new board is essentially the reverse of removal steps with the following considerations:

- Make sure the long grounding screw is returned to its proper location
- Be careful in returning all screws to their proper locations in that you avoid accidental contact with / damage to nearby board components

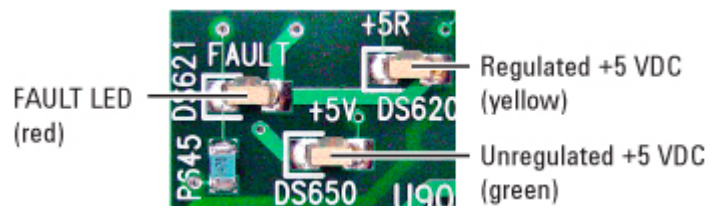

10 For the G3430-60150 Analog and power board there are three LEDs at the upper left corner of the board: red, yellow, and

green. The LEDs for the G3430-60151 Analog and power board are located in the middle on the right next to the **REMOTE** connector. When power is restored, yellow and green LEDs should be lit (on) indicating, respectively, that regulated and unregulated +5 VDC supplies are functional. The red LED is lit only under FAULT conditions.


**NOTE**

Note: The G3430-60151 analog and power board requires GC firmware version A.01.09 or greater. If you have replaced a G3430-60150 board with the newer G3430-60151 board, and the GC firmware is < version A.01.09, the GC will boot into manufacturing test mode (MMON). To use the GC, connect to it using Agilent Instrument Utility or Lab Advisor (see Using Lab Advisor Tools: firmware update) and update the firmware.

---

## Replacing Components Inside the Lower Rear Metal Cover

Item	Description	Qty
1	GC chassis	
2	Pneumatics fan	1
3	Screw	3
4	Flapper assembly	1
5	Oven fan motor	1
6	Transformer	1
7	AC power board	1

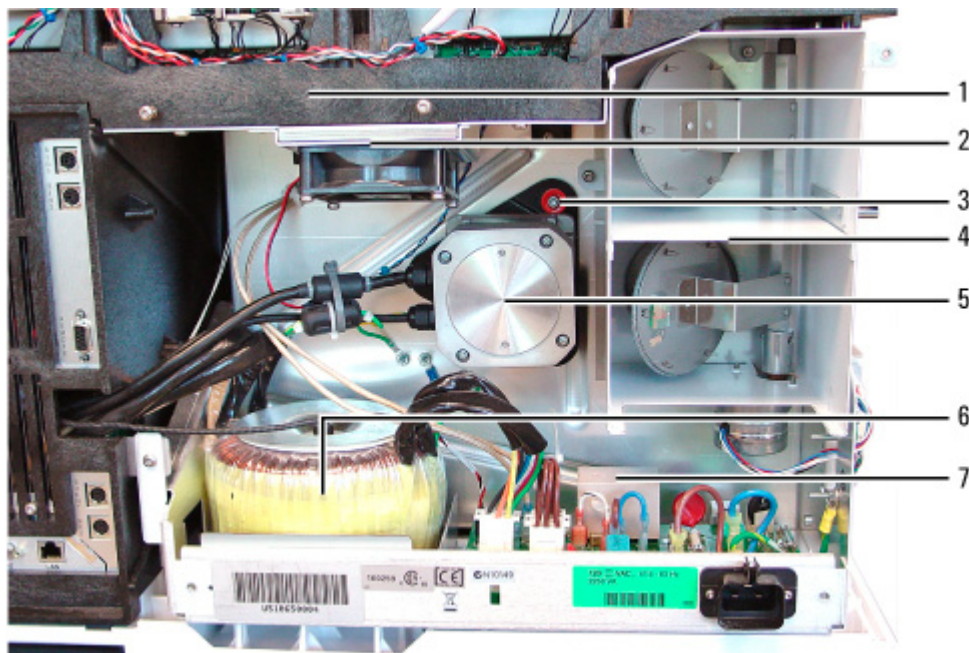



## Replacing AC board fuse sets

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.


### WARNING

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

- 3 Remove the lower rear metal cover from the GC.
- 4 Replace the necessary fuse set. Use of a fuse puller is recommended for ease of removal.

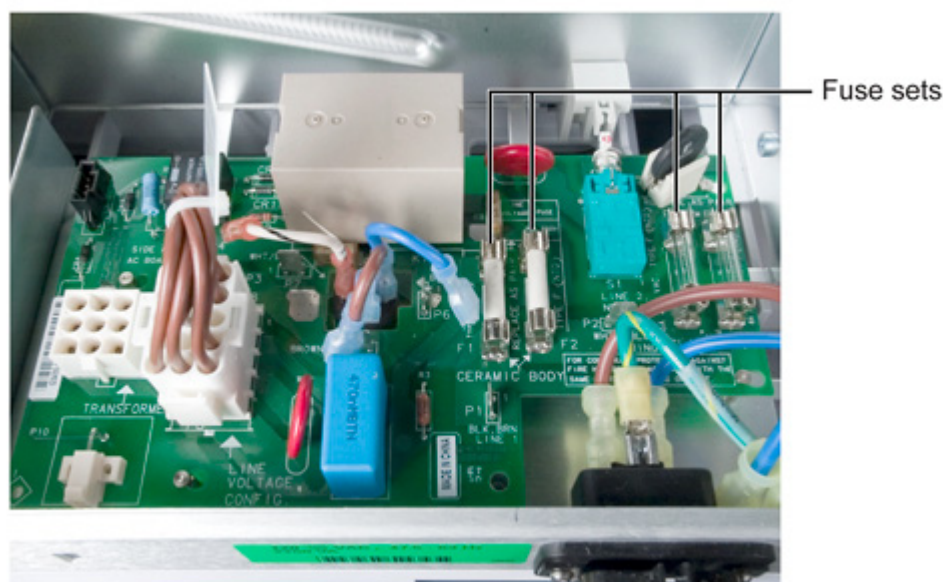

### NOTE

The two sets of fuses on the AC board are each always replaced as a set.

## Replacing the AC board

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.


### WARNING

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

- 3 Remove the rear lower metal cover on the GC.
- 4 The AC board is located below the oven flapper assembly. Disconnect cables (4), 2 wires to AC power receptacle, and 3 wires to the oven triac located below the board. The configuration plug on the original board must be removed and moved to the new replacement board.

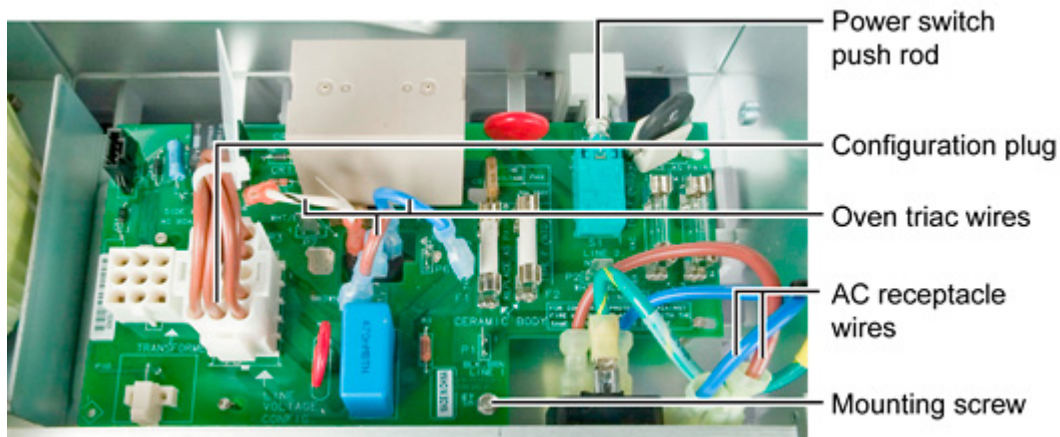

- 5 At the front right corner of the board, disconnect the on/off power switch's connecting push rod by pushing the rod forwards towards the front of the GC until it releases from the switch on the board.
- 6 Rotate the mounting thumbscrew counterclockwise until an audible "click" sound is heard. Then slide the board slightly towards the center of the GC and lift it off its standoffs.

**7** In reassembly, pay attention to the following:

- Do not forget to reconnect the power switch push rod.
- Make sure the new board is located on its standoffs. Then slide it slightly away from the center of the GC such that the mounting thumbscrew is properly positioned.
- Engage the thumbscrew threads by pressing it downwards gently while rotating clockwise.

## Replacing the oven triac

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.



### WARNING

Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

- 3 Remove the AC board. See [“Replacing the AC board”](#).

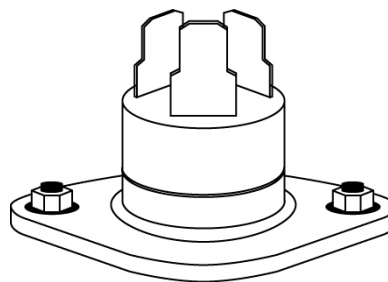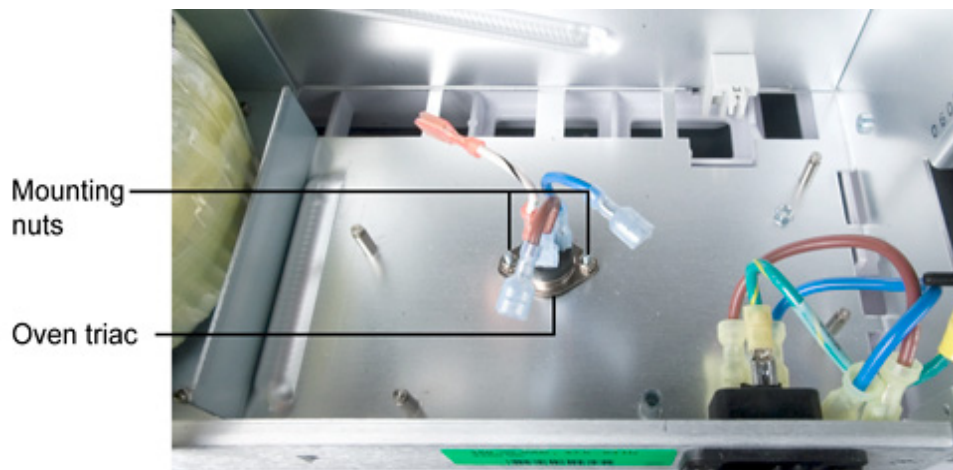



**CAUTION**

Carefully note BOTH physical orientation of the triac AND which of its three leads connects to each of its terminals. A mistake here in reassembly may destroy the new triac and/or result in a lack of proper oven control.

---

- 4 Remove three wires from the triac.
- 5 Remove two 5.5-mm hex nuts with captured lockwashers securing the triac.
- 6 Lift the triac straight up to remove it. Also, if necessary, remove the thin piece of plastic film found under the triac.
- 7 In reinstallation, remember to first install the piece of plastic film supplied with the new triac.

## Replacing the power transformer

- 1 Cool down the oven.
- 2 Turn off the GC main power switch and disconnect its power cord.

### WARNING

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

---

### CAUTION

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.


---

### WARNING

**The transformer is very heavy! Make sure you have a firm grip on the transformer before removing / replacing it.**

---

- 3 Remove the lower rear metal cover on the GC.

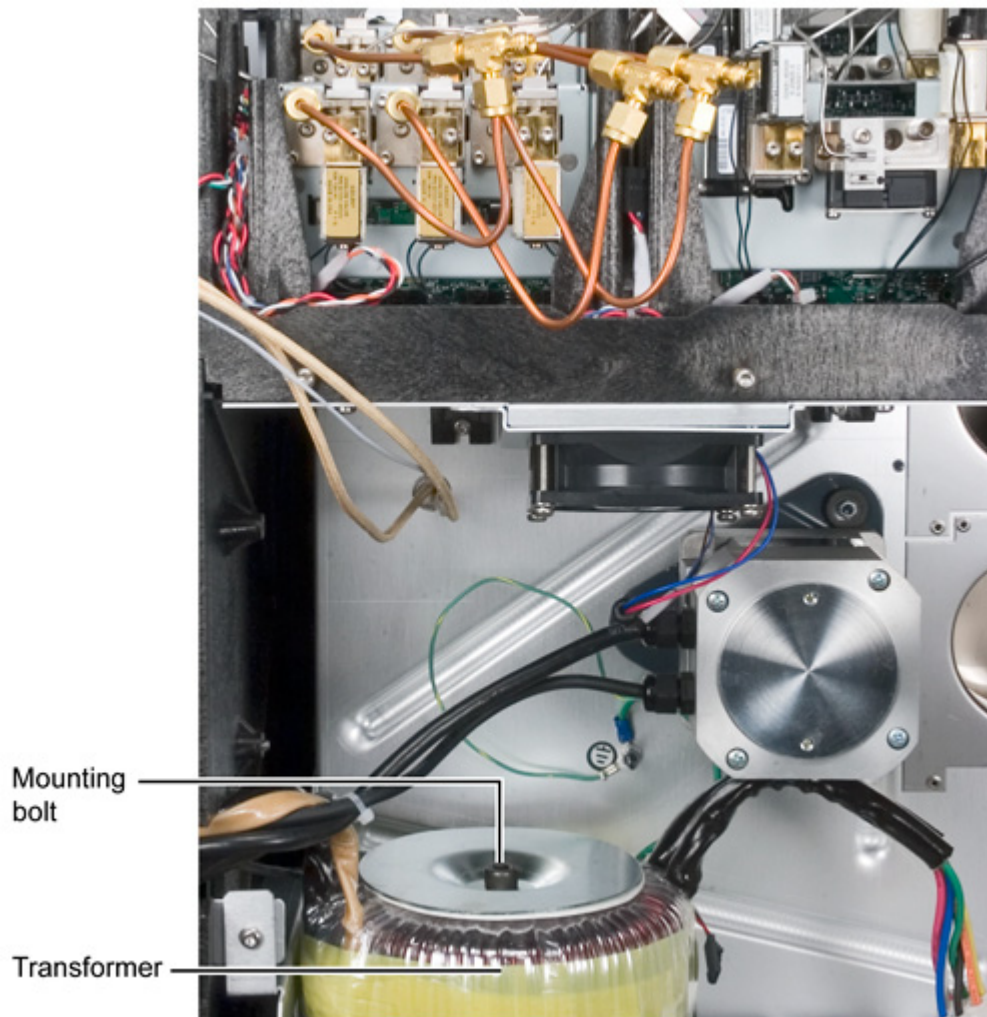

- 4 Disconnect the cable from the transformer to the AC board. Also disconnect the ground wire to the GC frame. Finally, disconnect the cable harness to the Analog and Power Board and feed the harness and connector back to the rear of the GC. Feed the two smaller cable connectors through first to maximize room to pass through the largest of the three cable connectors.
- 5 Note the orientation of the old transformer with respect to the various cables, then remove the Allen bolt to release the transformer. Carefully lift it from the GC.
- 6 In replacement, make sure the lower pad is in place, then orient the new transformer onto the pad.
- 7 Replace the upper pad and associated hardware. In restoring cables, feed the large cable and connector through to the electronics side of the GC first.



## 10 Valves

Valves 409

Valve Box 417

Actuators 423

Valve actuator drivers 431

Typical Valve Configurations 442

Troubleshooting 452



## Valves

Valves are used to modify the carrier gas flow during an analysis. Among their uses are:

- Highly reproducible injections, both liquid and gas
- Changing the order of columns
- Backflushing to eliminate high-boilers
- Column selection

and many others. Some examples are included in this section.

In this GC, valves are rotated by air-driven actuators. This keeps electrical activity away from where hot sample might leak. The air is controlled by solenoid valves located behind the right side cover.

## Introduction

This document does not provide instruction for first time installation of any of the options discussed. The add-on sheets, which accompany the various options, exist for just this purpose, and should be referenced when performing a first time installation.

The valves described in this manual are W-series minivalves, manufactured by VALCO Instruments Co, Houston, Texas. The valve body is made of Nitronic-60 nickel steel with 1/16-inch fittings.

Proper instrument operation will prolong the life of the valve system. Read all the accompanying information and avoid the following operational abuses:

- Exceeding the specified temperature and pressure ranges
- Plugging a valve with column packing or sample precipitation
- Scoring valve surfaces with column packing or particulates in liquid or gas sample
- Contaminating the system with samples (non-eluting materials) or poor quality support gases

Because valves operate best at a constant temperature, the GC places valves in their own heated compartment.


### WARNING

**To reduce the fire hazard when sampling flammable gases or liquids under pressure, operators should routinely make pressure-leak tests of the plumbing, fitting and valves. Both valve positions should be checked. Depending on the nature and pressure of the sample stream, periodic pressure leak test and visual inspection should be made since wear or use could cause leaks to develop. Leaks may occur inside the valve box and be concealed from the operators view.**

---

## Valco W-series minivalves

Valves consist of a driver, valve body, rotor, and preload assembly.

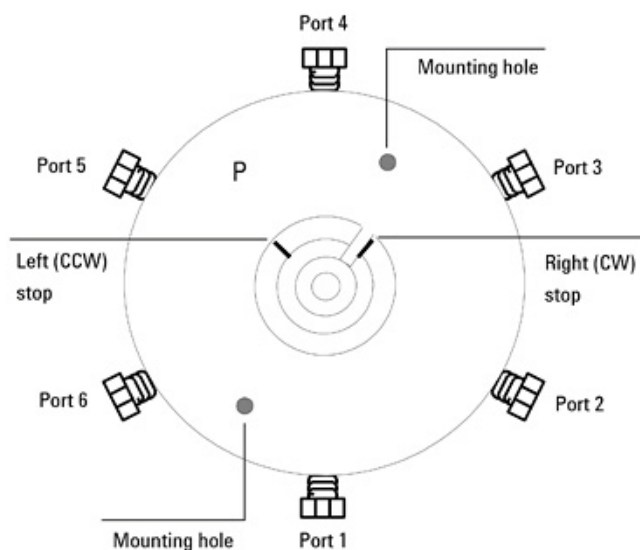

### Valve bodies

Body parts are made from Nitronics 60 nickel steel. If required, the valve may also be produced from Hastelloy C. External tubing (plumbing) is connected to the valve body ports by ferrules and fittings provided with the instrument.

The left (CCW, counterclockwise) and right (CW, clockwise) stops on general purpose valve bodies limit rotor rotation so the correct flow path results when the index pin is close to or against either stop of the index lip.

### CAUTION

Intermediate positions of the rotor may result in an interrupted flow path which could cause damage to the valve or other components in the chromatograph.

### Valve rotors

Rotor type can be identified by color:

- An **off-white** rotor is made of a PTFE composite and may be used from room temperature to 200°C.
- A **black** rotor is made of polyimide and may be used from 100 to 350°C.

#### CAUTION

The life of a valve is shortened if used outside its specified temperature range. Do not mix rotor types in the same system.

---

The rotor seats on a highly polished conical surface. When properly seated, the polished surface prevents leakage around the rotor and between non-selected ports. The finish precludes adsorption of most GC samples.

The rotor assembly is a one-piece part with an integral molded and machined conical hub and the parts necessary for proper seating. The sample contacts only the PTFE composite (low temperature) or polyimide (high temperature) as well as the stainless steel of the valve.


Grooves in the rotor surface form the paths between specific ports. The index pin prevents rotation beyond either stop of the index lip. Valve ports are connected by the grooves only when the index pin is close to or against either stop. Intermediate positions result in flow shutoff through the valve and possible damage if left in this position.

### Preload assembly (not shown)

This assembly consists of a body, a spring, and an adjusting nut. It holds the rotor in the valve body.



### Valco W-series minivalve

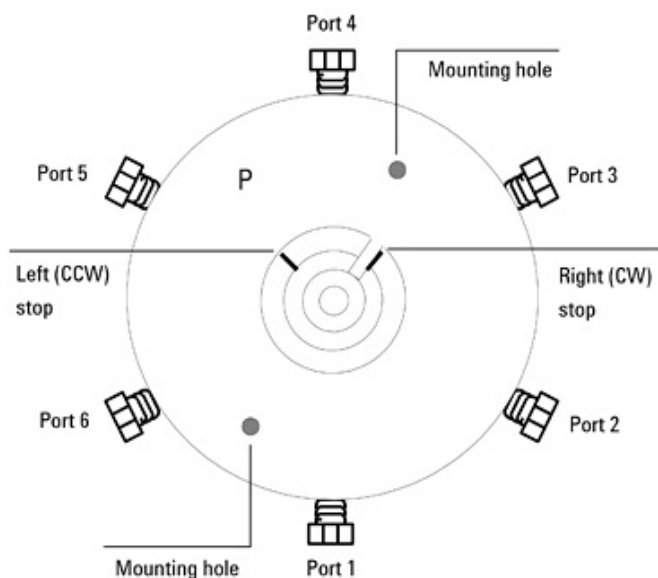

### Valve bodies

Body parts are made from Nitronics 60 nickel steel. If required, the valve may also be produced from Hastelloy C. External tubing (plumbing) is connected to the valve body ports by ferrules and fittings provided with the instrument.

The left (CCW, counterclockwise) and right (CW, clockwise) stops on general purpose valve bodies limit rotor rotation so the correct flow path results when the index pin is close to or against either stop of the index lip.

### CAUTION

Intermediate positions of the rotor may result in an interrupted flow path which could cause damage to the valve or other components in the chromatograph.

### Valve rotors

Rotor type can be identified by color:

- An **off-white** rotor is made of a PTFE composite and may be used from room temperature to 200°C.
- A **black** rotor is made of polyimide and may be used from 100 to 350°C.

#### CAUTION

The life of a valve is shortened if used outside its specified temperature range Do not mix rotor types in the same system.

---

The rotor seats on a highly polished conical surface. When properly seated, the polished surface prevents leakage around the rotor and between non-selected ports. The finish precludes adsorption of most GC samples.

The rotor assembly is a one-piece part with an integral molded and machined conical hub and the parts necessary for proper seating. The sample contacts only the PTFE composite (low temperature) or polyimide (high temperature) as well as the stainless steel of the valve.

Grooves in the rotor surface form the paths between specific ports. The index pin prevents rotation beyond either stop of the index lip. Valve ports are connected by the grooves only when the index pin is close to or against either stop. Intermediate positions result in flow shutoff through the valve and possible damage if left in this position.


## General purpose valves

The 6- and 10-port general purpose valves are suitable for column switching, isolation, backflushing, and other uses as well as gas sampling.

## Gas sample valves

The standard gas sample valves have 1/16-inch zero dead volume fittings and an internal port diameter of 0.016-inch.

### 6-port valve (actuator side view) showing flow path grooves

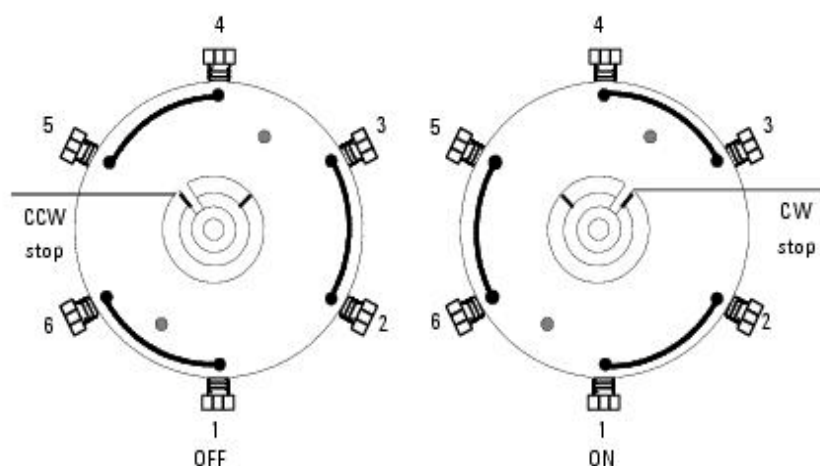

## Gas sample loops

A 0.25 mL sample loop is included with all gas sampling valve systems. 10 mL and 5 mL loops occupy one valve position, limiting the number of valves that can be housed in a valve compartment.

## Adjustable restrictor valves

Adjustable restrictors are used to balance flow resistance between the two valve positions. They are available with ambient to 225°C (part no. 0101-0633) or ambient to 350°C (part no. 0101-0948) operating ranges.

## Liquid sample valves

Agilent Technologies offers 4-port LSVs with 0.2, 0.5, or 1  $\mu\text{L}$  internal loops.

These valves are designed for liquefied gases under pressure such as ethane, propane, butane, LNG, etc. They are not intended for nonvolatile liquids (at room conditions) where a concealed leak could allow an accumulation or pool of liquid to form that may present a significant fire hazard.

All liquid sample valves have 1/16-inch fittings. The 0.5  $\mu\text{L}$  size is available in low- (1000 psig limit) and high-pressure (5000 psig) versions; the other sizes are low-pressure (1000 psig) only.

An adjustable restrictor may be used on the sample outlet line to maintain internal sample pressure to keep a compressed gas liquefied.

### Temperature ranges

1/16 inch PTFE rotor valves	0 to 200 °C
Adjustable restrictor valves	Ambient to 225 °C

### CAUTION

The life of an LSV is shortened if used outside its specified pressure and temperature ranges. Highly dangerous leaks can occur if the valve box temperature ever exceeds the specified temperature limits.

### Adjustable restrictor valve

The adjustable restrictor supplied with a liquid sample valve is designed for temperatures up to 225°C; it is NOT compatible with high-temperature valves.

## Valve Box

### Installing the valve box

The valve box sits on top of the GC and contains the valves and plumbing, heated zones and sensors, and the insulation.

- 1 Use a pair of diagonal cutters to remove the metal cutout on the top of the GC.

**WARNING****Be careful of sharp edges!**

Cut toward the outside edge of the tabs so that the valve box will fit correctly.

- 2 Screw the aluminum standoffs into the valve box bottom plate using one Torx T-20 screw for each standoff.
- 3 Screw the valve box bottom plate onto the top of the GC using three Torx T-20 screws.
- 4 There are six holes in the plate for up to four valves. The outside holes on each side of the plate are for valves and the inside holes are for plumbing into the GC. Punch out the pre-perforated insulation from the holes you plan to use.
- 5 Use two Torx T-20 screws to mount the heater block(s) on the valve box bottom plate.
- 6 Install the heater/sensor into the heater block(s) and secure the leads to the valve box bottom plate with the U-clamp and two Torx T-10 screws.

**NOTE**

Make sure that the sensor is seated all the way in the heater block.

- 7 Insert the valve(s) into the appropriate hole(s) in the heater block and secure each valve with two long Torx T-10 screws.
- 8 Plumb the valves in the appropriate configuration. A variety of valve configurations are diagrammed in the “[Typical Valve Configurations](#)” section in this chapter.
- 9 Orient all of the valves in the OFF (CCW) position.


- 10** Place the valve box insulation in the valve box top. Make sure the perforations in the insulation line up properly with the valve box top.
- 11** Insert the insulation retainer plate in the bottom of the valve box top and secure it to the threaded studs using the two 5.5-mm nuts.

Tighten the nuts until they are flush with the top of the stud. Make sure the holes in the retainer plate line up with the perforations in the insulation.


- 12** While wearing a pair of protective gloves and holding the valve box top over a waste receptacle, punch or cut out the insulation from the appropriate holes.
- 13** Use two Torx T-20 screws to secure the valve box top over the installed valves, making sure the heater/sensor leads are routed under the appropriate cutout(s).

Make sure the valve box top is oriented correctly with holes punched out over the installed valves.


- 14** Install the actuators as described in the “[Actuators](#)” section of this chapter.

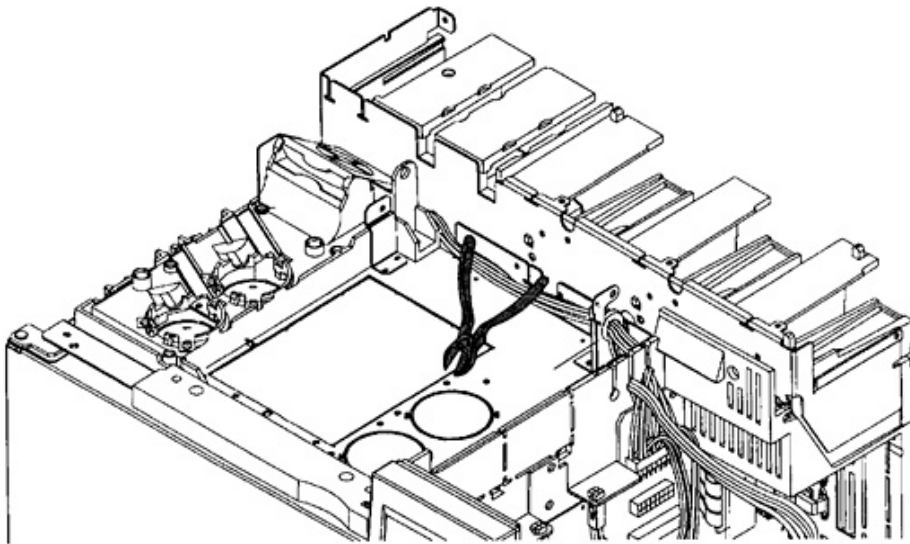


**Figure 7** Removing the sheet metal cutout. 6890 shown.

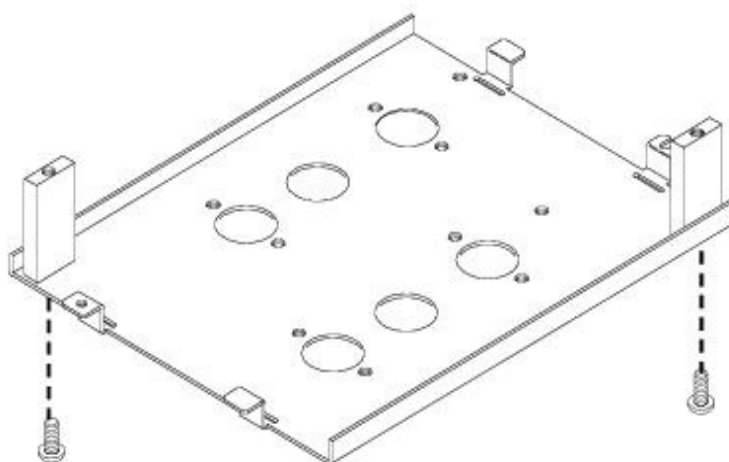


**Figure 8** Attaching the standoffs to the bottom plate

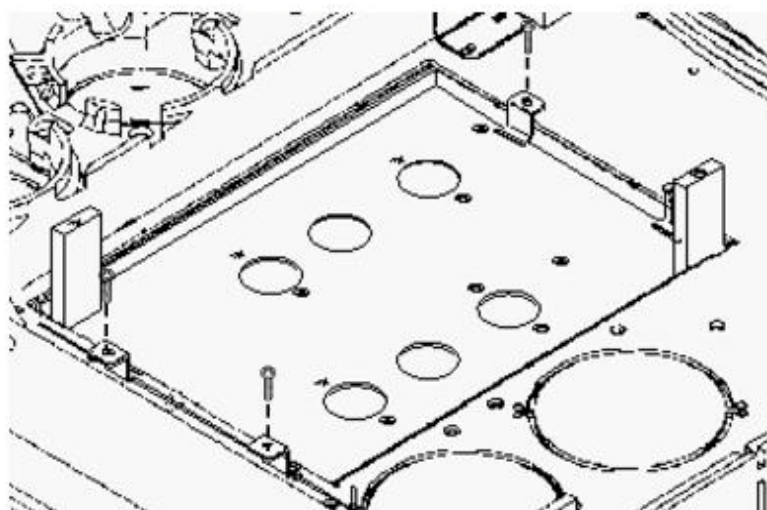

**Figure 9** Attaching the bottom plate to the GC

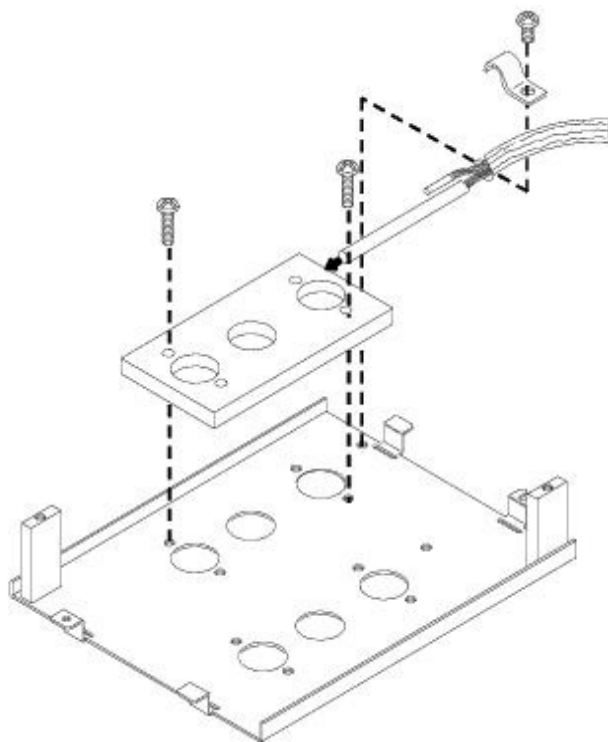


**Figure 10** Attaching heater block to the bottom plate

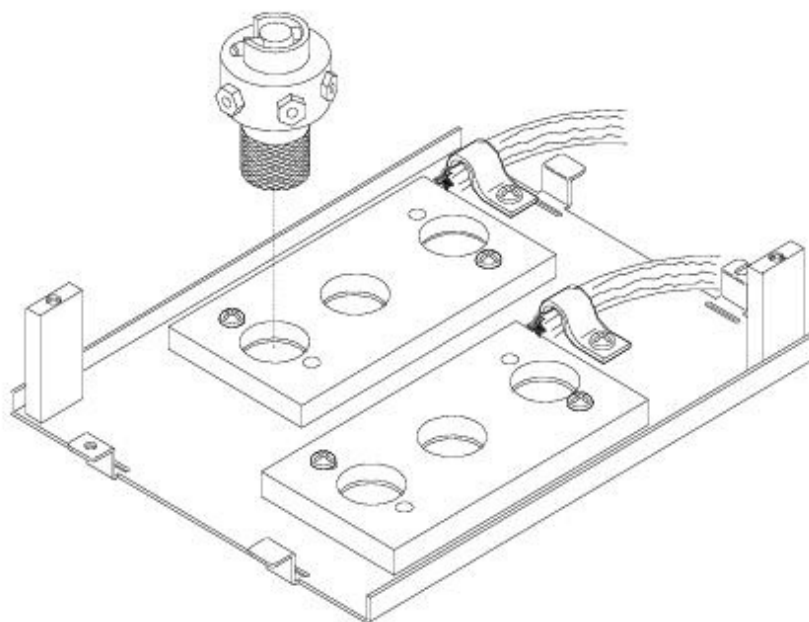


**Figure 11** Installing valves



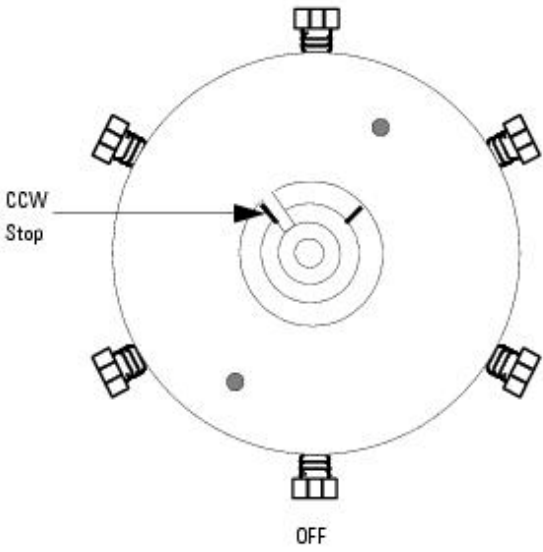

**Figure 12** Valve in the OFF position (top view)

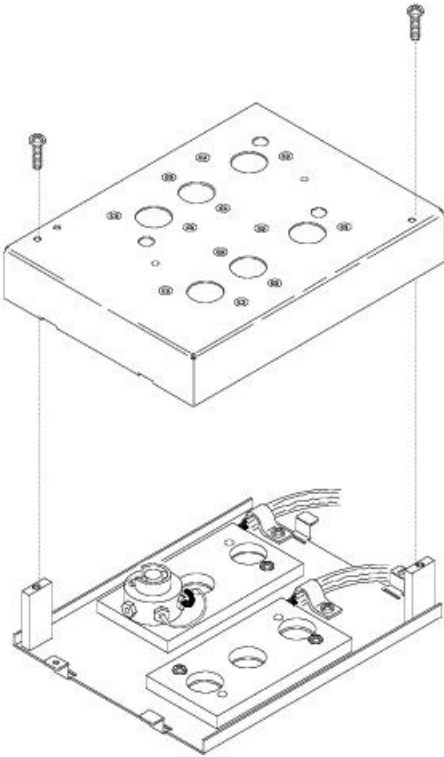

**Figure 13** Installing the valve box top assembly

## Removing the valve box assembly

- 1 Place the main power switch in the off position.

### WARNING

**Hazardous voltages are present in the instrument when the power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before working on the instrument.**

---

- 2 Unplug the line power cord from its receptacle.
- 3 Allow some time for the oven and heated zones to cool.
- 4 When the oven has cooled, turn off all gas supplies.
- 5 Switch the solenoid valve off so the actuator is in its fully extended position (piston rod extended).
- 6 If variable restrictors are present, remove their mounting hardware in the following order: two Torx T-20 screws, hex nut, and mounting bracket for each restrictor valve.
- 7 Remove the two Torx T-20 screws securing the valve box top assembly to the standoffs. Lift the valve box top assembly straight off the valve box. Be careful not to move the valve rotor index pin from its "at rest" position.

### NOTE

If valve/actuator alignment is to be made, see ["Valve actuator alignment"](#) in this section.

---

- 8 To reassemble: Align the two mounting holes in the valve box top assembly with the standoffs in the valve box. Lower the box top assembly until it rests on the standoffs.
- 9 Secure the valve box top assembly with two Torx T-20 mounting screws. Tighten these screws firmly. Reinstall hardware for variable restrictors if present.
- 10 Exercise the valve(s) on and off a few times to verify operation.


## Actuators

Valve rotors are driven by air-operated actuators. The air is controlled by solenoid valves inside the right side panel.


### Installing the actuators

The actuators use pneumatic pressure (40 to 70 psi) to switch the valves between their two positions.


- 1 After installing the valves and valve box as described in the “Valve Box” section, you can install the valve actuators.
- 2 Set each actuator to the appropriate degree of rotation. Move the grenade style pin to the hole on the actuator marked with the correct degree of rotation, as shown below:
  - Four port valves—Place the pin in the 90° hole
  - Six port valves—Place the pin in the 60° hole
  - Ten port valves—Same as six port valves, but with the tubular 36° actuator limiter on the pin.
- 3 Mount an actuator over each valve installed using two Torx T-20 screws.
- 4 Engage the actuator drive shaft coupler with the valve.
  - a Loosen the hex nut on the actuator near the drive shaft.
  - b Slide the shaft down.
  - c Insert a flat bladed screw driver in the slot on the top of the actuator and turn the shaft back and forth until you feel the coupler engage the valve.
  - d Tighten the set screw.
- 5 Install the valve actuator drivers.

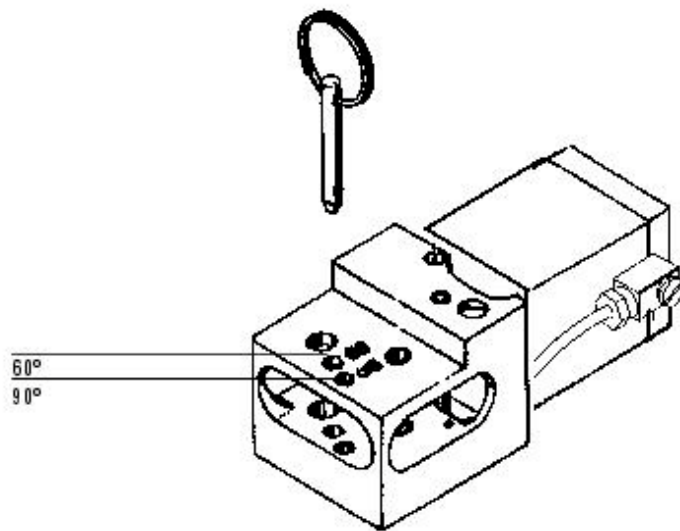


**Mounting the actuator on the valve box**

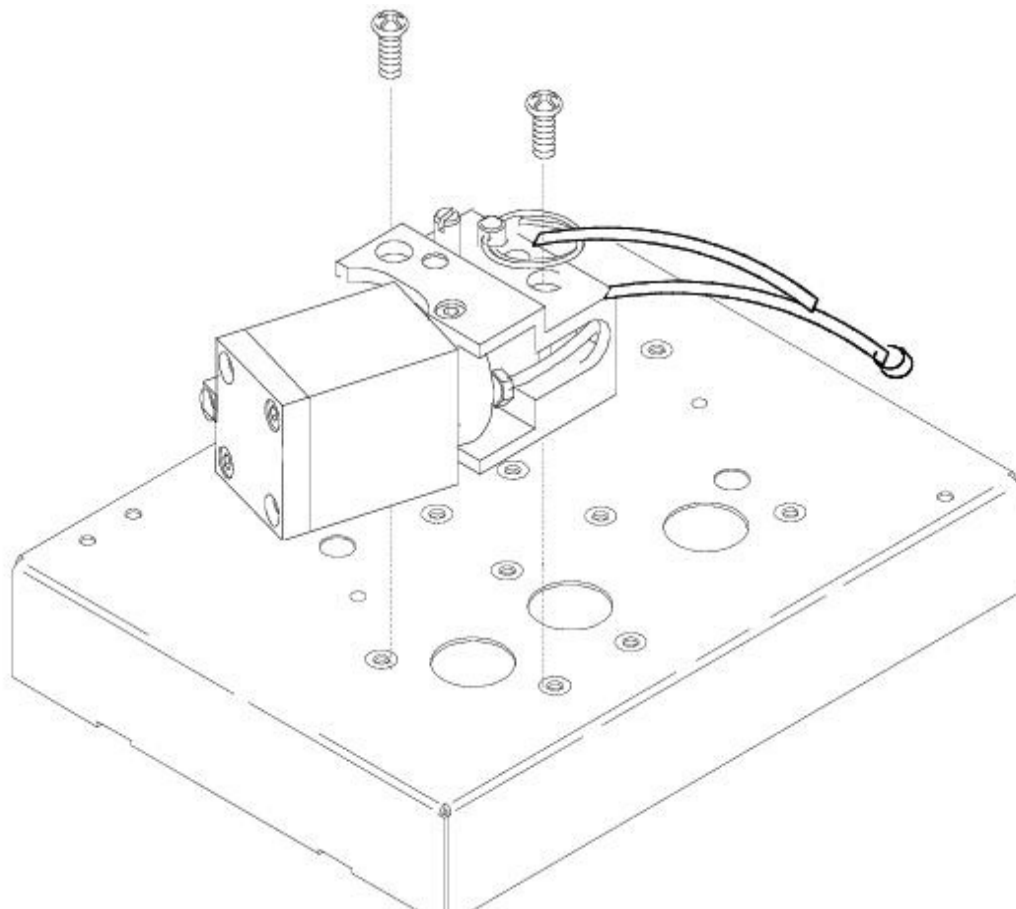

**Engaging the actuator drive shaft with the valve**

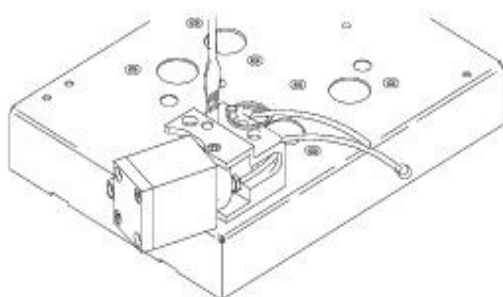

### Setting the actuator's degree of rotation

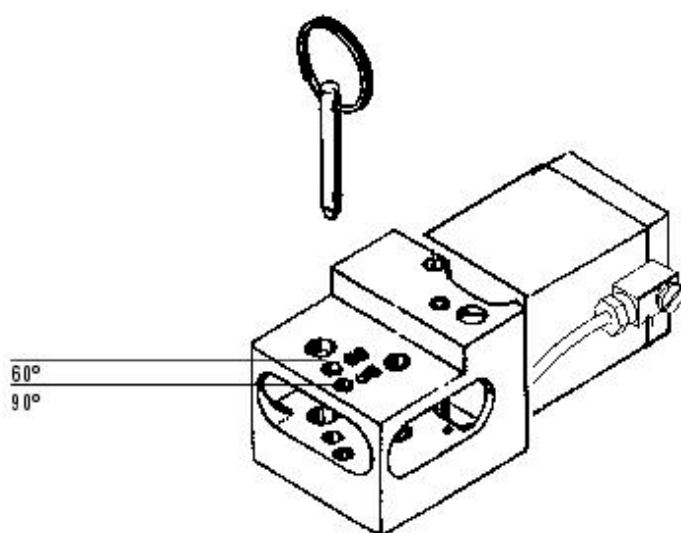


### Installing the actuator on a heated valve box

#### WARNING

The valve box is insulated with a fibrous material that can cause irritation to the skin, eyes, and mucous membranes. Always wear gloves when working with the insulation. If the insulation is flaky/crumbly, wear protective eyewear and a respirator.

Before continuing, install the heated zone kit (if purchased) as described in the instructions included with that accessory. Also install any new rotary valves into the valve box. Leave the GC covers off and the power cord disconnected.


- 1 Hold the upper valve box over a trash can, then remove the pre-cut insulation from the appropriate hole in the upper valve box. Remove any excess or loose pieces from around the hole.

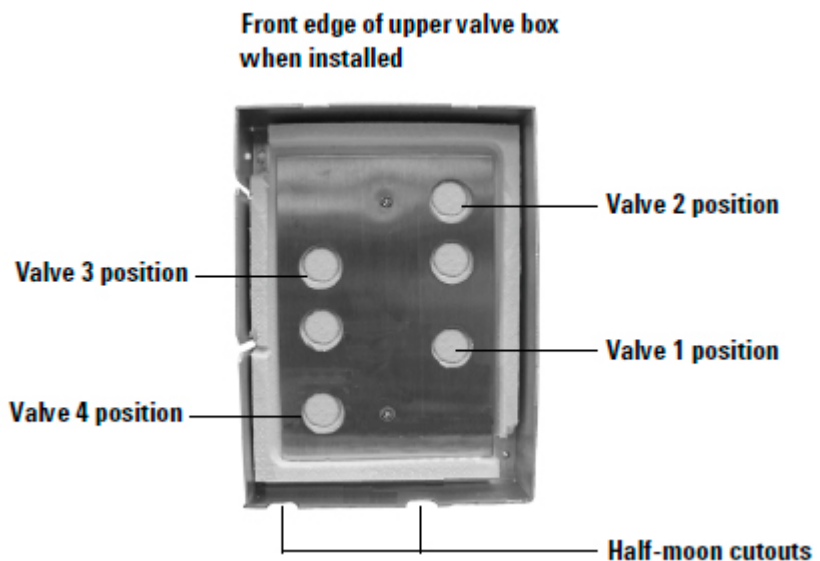


- 2 Verify that all valve rotors are in the full counterclockwise position (valve Off).
- 3 Locate the two half-moon cutouts at the bottom back of the upper valve box. Place the upper valve box on top of the lower valve box assembly, routing the heater/sensor wires through the cutouts. Secure with two T- 20 mounting screws.
- 4 Install the new valve actuator.
  - a Insert the actuator through the valve box cover.
  - b Push the coupling/shaft assembly downward with a flathead screwdriver until the slot on the coupling engages the rotor index pin.

If the coupling and valve do not engage, check that both are fully counterclockwise and try again. If necessary, turn the shaft slightly to engage the coupling.


- c Secure the actuator to the valve box with two T-20 Torx screws provided in the kit.

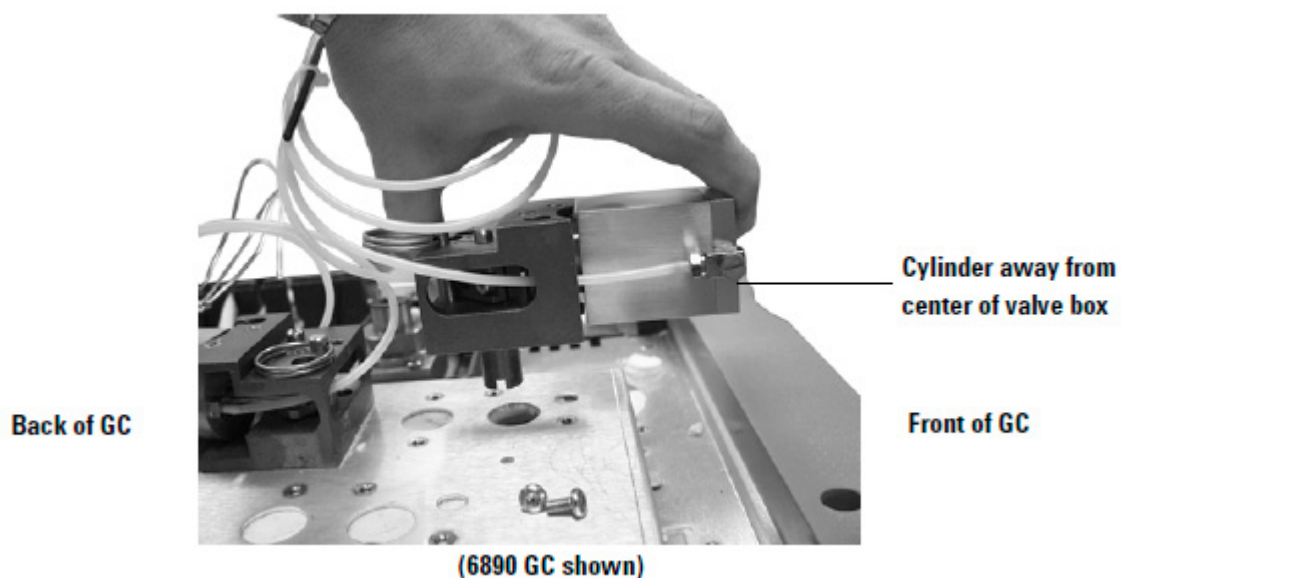

- 5 For each actuator that mates with a newly installed valve:
  - a Loosen the adjustment set screw.

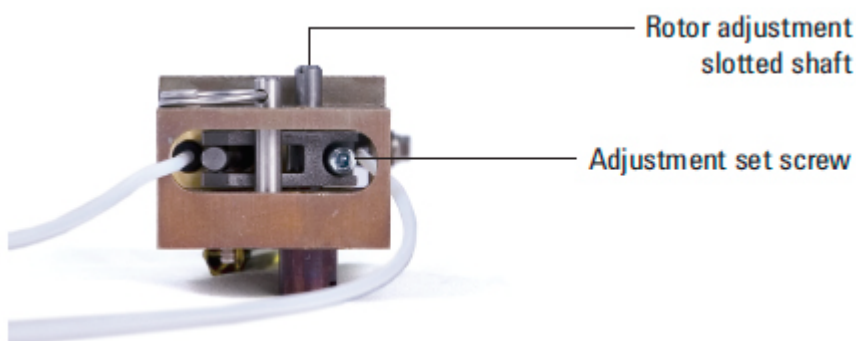

- b Using a flathead screwdriver, turn the rotor adjustment shaft counterclockwise until it stops, then back it off a small amount ( $< 1$  mm) to set one end of the rotor's motion.

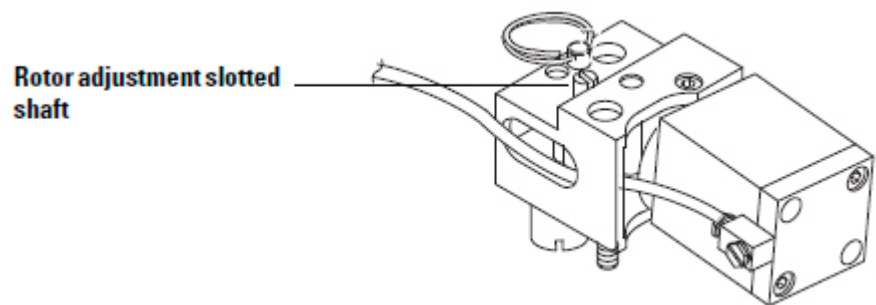


- c Tighten the adjustment set screw.



### Installing the actuator on a side mount bracket

- 1 Install the side mount bracket on the left side of the GC using four machine screws.
- 2 Attach the clamp ring to the back of the bracket with two screws.
- 3 Align the actuator assembly with the upper or lower standoffs. The cylinder must be down. Insert the actuator drive through the hole in the bracket, and secure with two T-20 Torx screws.

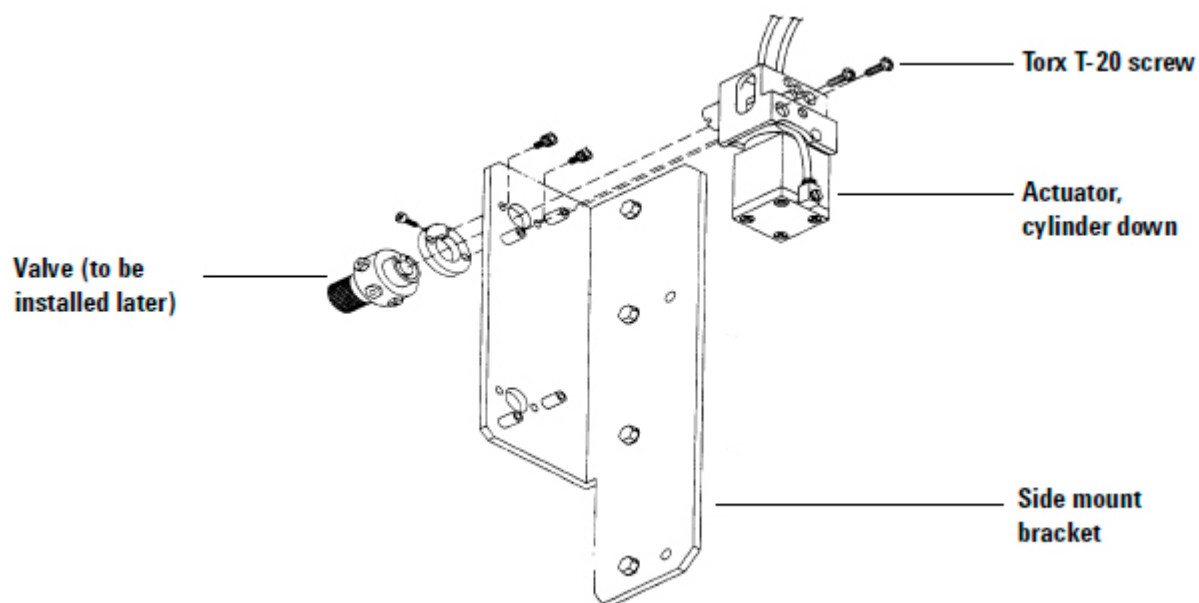

## Valve actuator alignment

- 1 Remove the valve box top assembly. See steps 1 through 3 of “Removing the valve box assembly” for the procedure.
- 2 Loosen the actuator link arm lock screw at each actuator with a 3 mm hex key wrench so that the coupling/shaft assembly is free to rotate. Push the coupling shaft fully into the actuator.
- 3 Turn the valve rotor index pin of each valve counterclockwise until it is 0.010 inch (0.25 mm) from the counterclockwise (left-hand) valve stop.
- 4 Reinstall the valve box top assembly.
- 5 Gently rotate and push the coupling/shaft assembly with a blade-type screwdriver until the slot on the coupling fully engages the valve rotor index pin. Repeat this procedure for each valve installed.

**CAUTION**

Use care in performing the following operation so as not to accidentally turn the valve rotor away from its preset (step 3) position.

---

- 6 Make sure that all solenoid valves are turned "off" by the appropriate valve controller. Turn on the air supply to the solenoid valve(s). The piston rod of each actuator will move all the way out to the extended (OFF) position. Very firmly tighten the link arm lock screw for each actuator.
- 7 Install the hardware for any variable restrictors present.

## Valve actuator drivers

### Solenoid differences


Beginning in April 2010, a new type of solenoid valve became available to drive the GC valve actuators. The new solenoid valve, made by Airtek, is not compatible with the original solenoid valve. They function in the same way, but mount differently and require different valve brackets.

#### Which solenoid valve is installed?


Any valve installed on a GC before April 2010 used the original style solenoid. Upgrading to the new solenoid requires a replacement of the entire driver stack assembly (valve bracket, elbow, solenoids, and mounting hardware).

You can easily distinguish between the solenoid valve types, as shown below. The new style solenoid valve has a silver metal fitting instead of a red pressure vent button.

Original style

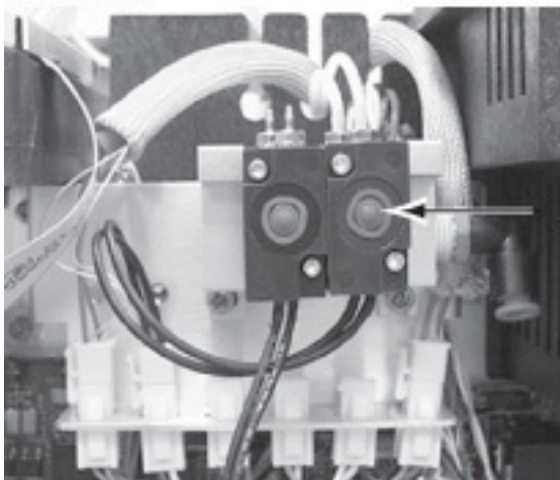

New style

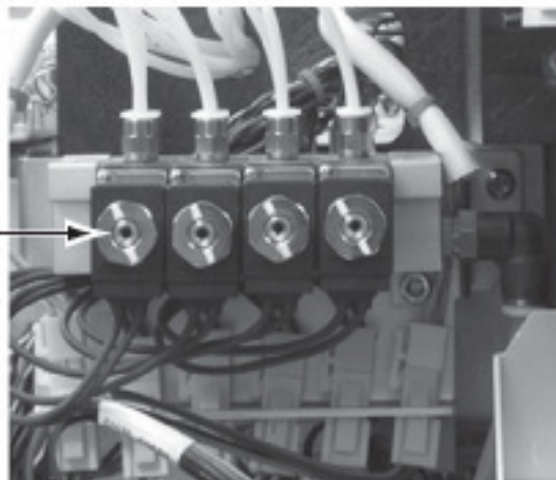

#### Replacing an original solenoid valve

If replacing an original style solenoid valve, first try to order the original replacement part. If that part is not available, then you will need to replace the entire solenoid stack assembly. See the chapter on “[Valves](#)” in the IPB for assembled stack part numbers.

### Installing a new valve stack assembly


If installing a new valve stack assembly, order the new style parts.

## Original style

### Assemble the valve driver block

The valve driver block accommodates up to four valve drivers. A valve driver must be installed for each valve/actuator installed.


- 1 Install two mounting posts on the intake/exhaust endplate (two large threaded holes). Install an O-ring in the supply/exhaust ports on the inside of the plate.
- 2 Slide a valve driver over the mounting posts in the orientation shown. Install two O-rings in the valve driver supply/exhaust ports as shown.

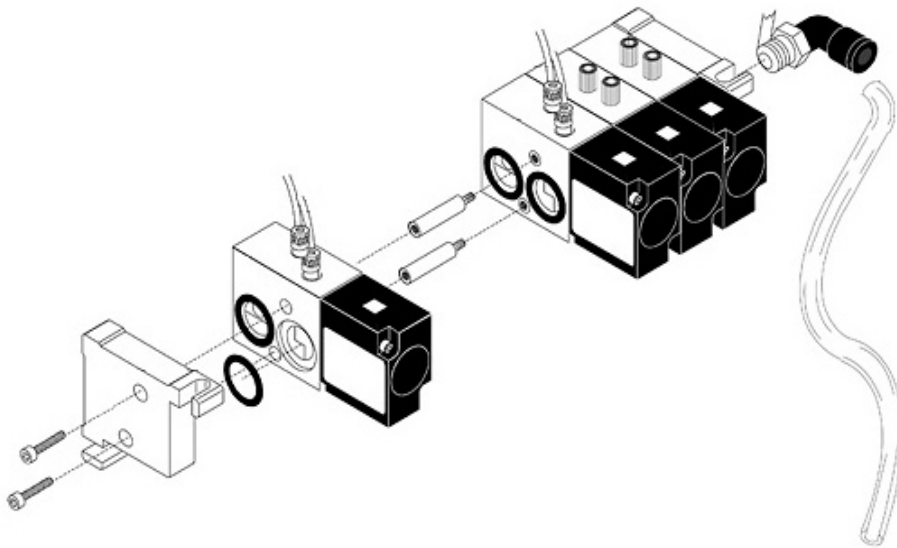

- 3 For each additional valve installed, install two more mounting posts and a valve driver with O-rings in the same manner as the first.
- 4 When all the drivers have been installed, screw on the other end plate with the two hex screws as shown.

### Install the bracket and cabling (original style)

- 1 Install the valve actuators before installing the valve actuator drivers.
- 2 Screw the valve driver bracket into the right side of the GC using the two captive screws.
- 3 Plug the connectors on the valve driver cable harness up through the slots on the valve driver bracket.
- 4 Plug the larger  $2 \times 2$  heater sensor connectors (P1, P2) into the outside slots and the smaller  $1 \times 2$  valve driver connectors (P3 to P6) into the four middle slots.
- 5 Plug in the heater/sensor lead(s) from the valve heater blocks on top of the GC. Thread the heater/sensor lead(s) to the right side of the instrument, through one of the keyhole wiring slots and into the P1 or P2 connector on the actuator bracket.

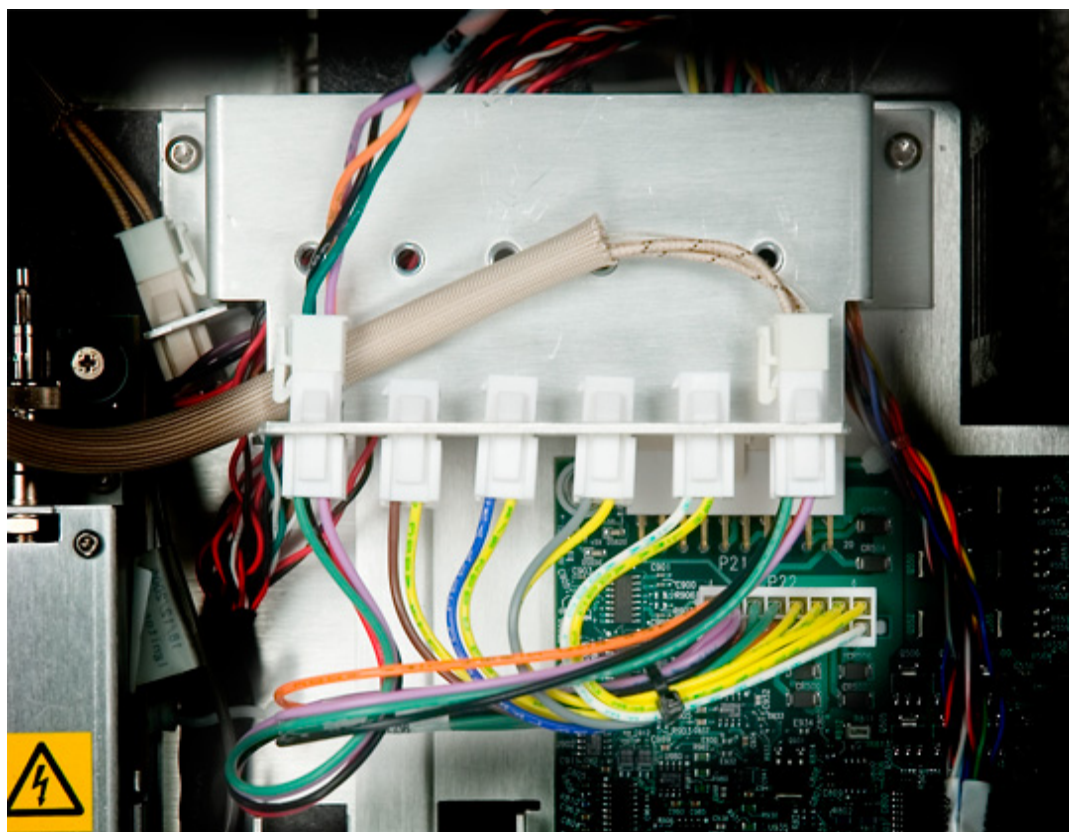

### Installing the valve driver block

- 1 Slide the valve driver block down into the driver bracket until the drivers plug into the connectors.

**NOTE**


**To remove drivers from the driver block**, use a hex wrench to unscrew the two hex screws on the left side of the block. Remove the driver, collapse the block to the width of the remaining drivers and reinstall the hex screws.

---

- 2 Plug the other end of the valve driver cable harness into the P22 connector on the analog and power board.
- 3 Wrap the threaded ends of the 90° elbow fitting in PTFE pipe tape. Screw the fitting into the air supply intake on the side of the valve driver block facing the rear of the GC.

The supply intake is the outside threaded hole, the one farthest from the analog and power board.

- 4 Run a length of 1/4-inch tubing out the hole on the rear of the GC in the lower left corner (when facing the rear of the instrument). Connect this tubing to your air supply.
- 5 Unscrew the screw plugs (if present) on the top of each driver you are using. Replace the screw plugs with tubing connectors.
- 6 Plumb the tubing from each installed actuator to the tubing connectors on the corresponding driver. The tubing from the side of the actuator goes to the connector farthest from the board. Grip the tubing with a piece of sandpaper and push it onto the tubing connector.

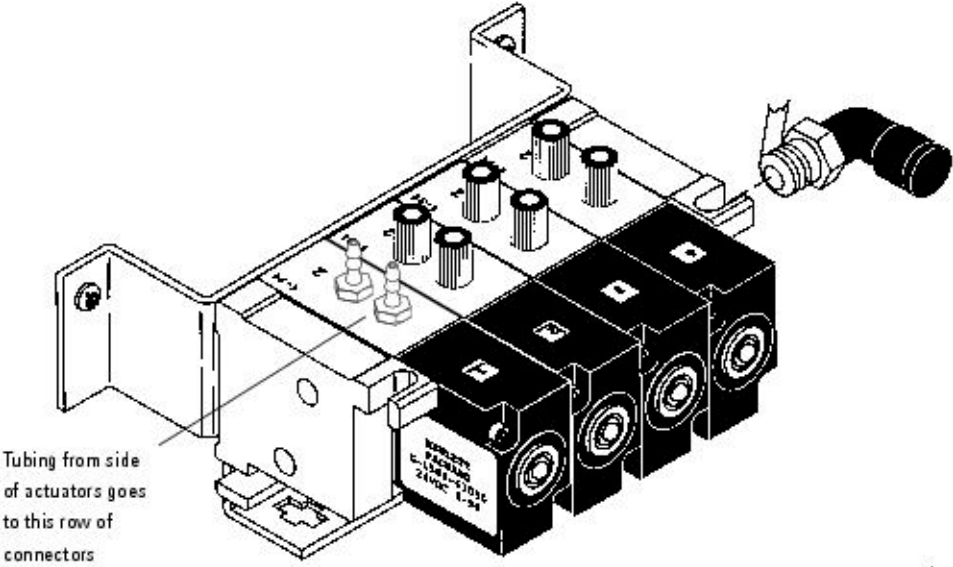

## New style

### Replacing a solenoid valve (new style)

To replace a new style solenoid valve:

- 1 Turn off the valve actuator air supply.

#### **WARNING**

The detector, inlet, oven top, and valve box may be hot enough to cause burns.

---

- 2 Remove the GC right side cover.


#### **WARNING**

Hazardous voltages are present when the GC electronics are exposed.

---


- 3 Press and hold the switch on the bottom of one of the solenoid valves until no residual air pressure remains.

Press to vent  
pressure

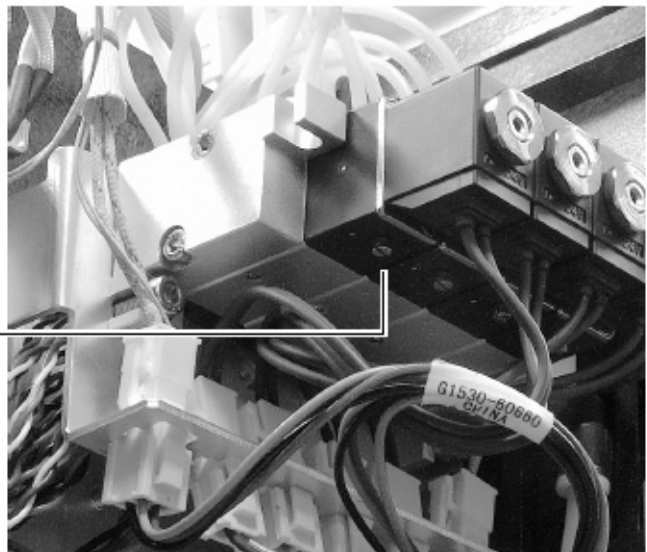


- 4 Turn off the GC and unplug the power cord. If installed, allow the heated valve box to cool.
- 5 Disconnect the valve box heater wires to the A1 connector.




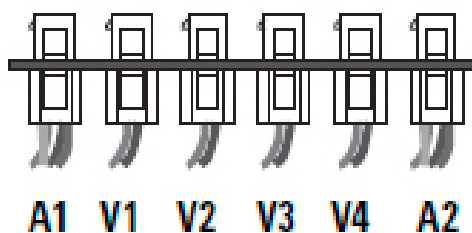


- 6 Disconnect the tubing from the broken or damaged solenoid valve.
- 7 Remove the two screws that hold the end plates and existing valves together.

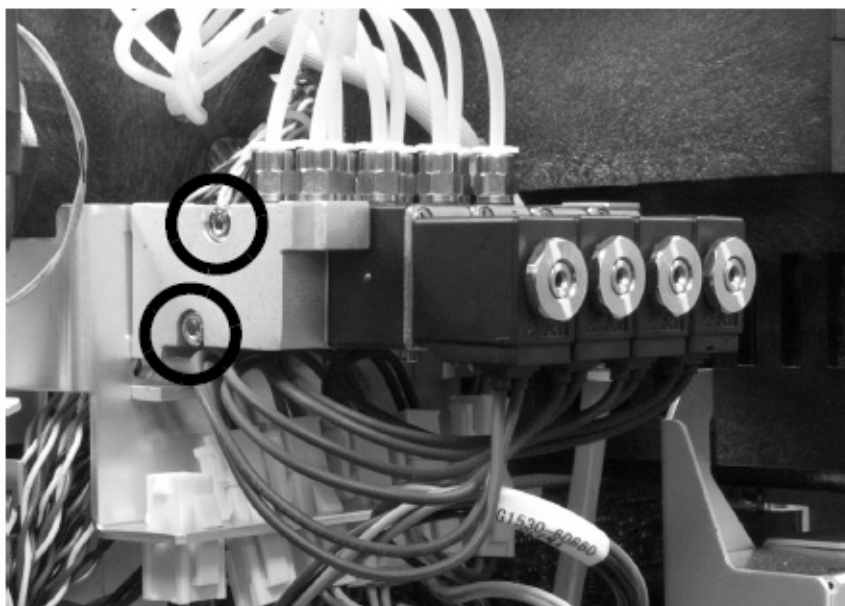

- 8 Use a 4-mm hex key to remove the screw that holds the solenoid stack left end plate in place.

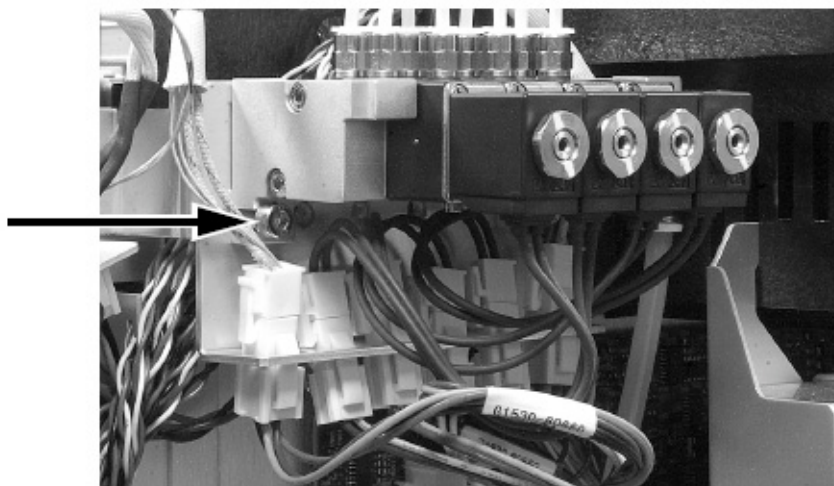


- 9 Remove the broken or damaged solenoid valve wires from the wire connector on the valve bracket.
- 10 Replace the solenoid valve(s).
- 11 Plug the new solenoid valve wires into the appropriate valve bracket connectors. See the diagram below.

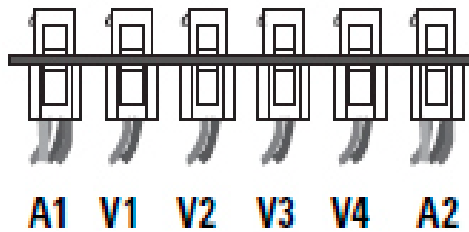


- 12 Reassemble the solenoid stack including the new solenoid valve and left end plate. While holding the valves in one hand, align them as best as possible, then slide a screw through the left end plate and solenoids.
  - Use 25 mm screws for 1 solenoid
  - Use 50 mm screws for 2 solenoids
  - Use 70 mm screws for 3 solenoids
  - Use 85 mm screws for 4 solenoids

Repeat for the second screw. Both screws should extend a little bit out of the last solenoid.
- 13 Line up the stack screws with the holes in the right end plate and assemble. Before tightening, align the solenoid valves so their tops are even.
- 14 Use the screw removed in step 8 to install the left end plate to the valve bracket. Be careful not to trap any wires behind the plate.
- 15 Reconnect the valve box heater wires to the A1 connector.
- 16 Connect the tubing from the actuator to the new solenoid valve by pushing it into the press fitting.
  - The tubing from the side of the actuator connects to the outside fitting on the solenoid (the fitting furthest from the electronics board, labeled 2).
  - Trim the tubing as needed to keep installation neat.

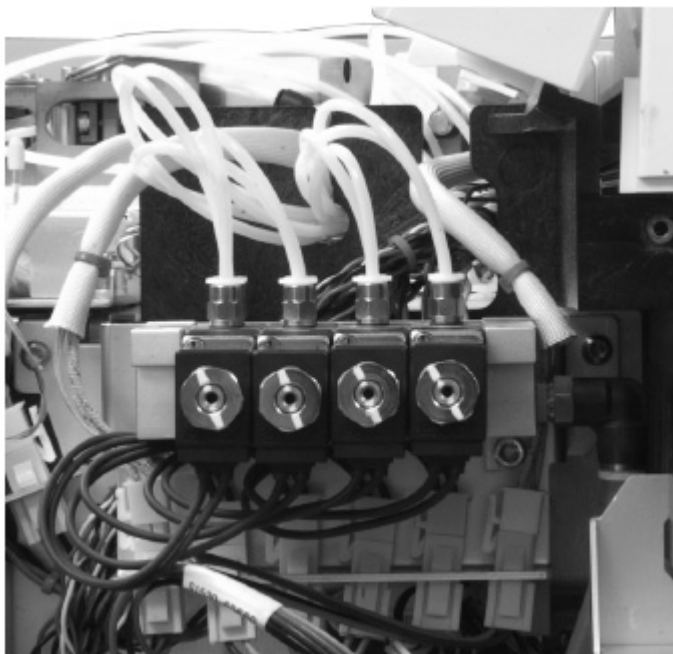

- 17** If installing new tubing, cut the tie-wrap around the tubing leading from the actuator and route the tubing as shown below.

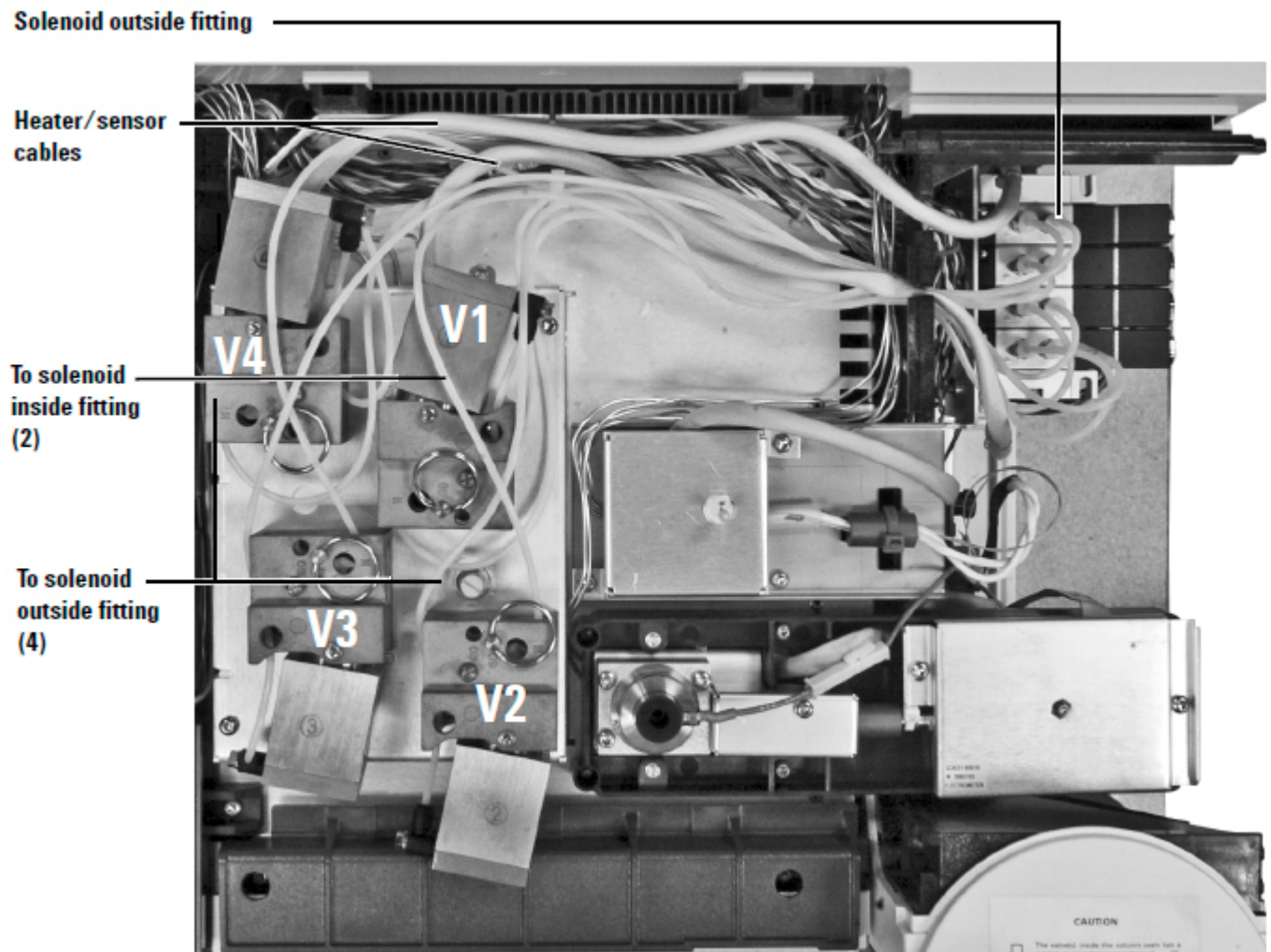

18 Turn on the valve actuator air.

### WARNING

Hazardous voltages are present when the GC electronics are exposed.

- 19 Connect the GC power cord and turn on the GC.
- 20 Press the switch under each solenoid to activate it and pressurize the system.
- 21 The solenoid valve replacement is complete. If no other accessories or parts require installation, reinstall all GC covers.

### Installing the valve driver (solenoid) bracket, cabling, and elbow fitting (new style)


- 1 Install the valve actuators before installing the valve actuator drivers (solenoids).
- 2 Screw the valve driver bracket into the right side of the GC using the two captive screws.
- 3 Plug the valve driver bracket cable harness into P22 on the GC board.
- 4 Plug the connectors on the valve driver cable harness up through the slots on the valve driver bracket.

Plug the larger  $2 \times 2$  heater sensor connectors A1 and A2 into the outside slots and the smaller  $1 \times 2$  valve driver connectors V1 to V4 into the four middle slots.


- 5 Wrap the threaded ends of the  $90^\circ$  elbow fitting in PTFE pipe tape. Screw the fitting into the air supply intake on the side of the valve driver block facing the rear of the GC.

The supply intake is the outside threaded hole, the one farthest from the analog and power board.

- 6 Run a length of 1/4-inch tubing out the hole on the rear of the GC in the lower left corner (when facing the rear of the instrument). Connect this tubing to your air supply.

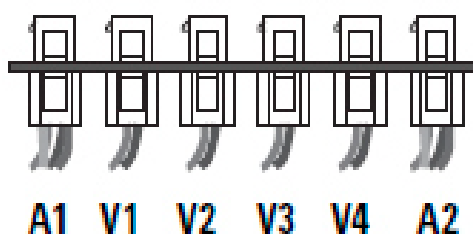

## Typical Valve Configurations

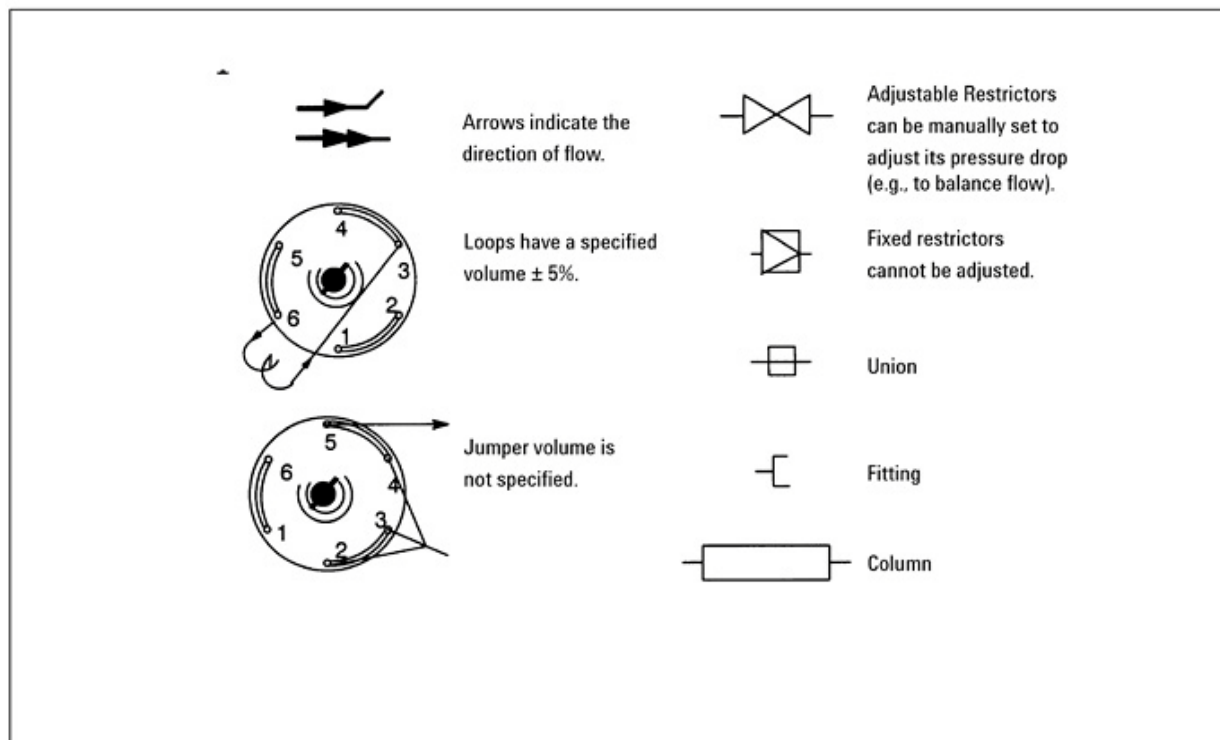


## Option numbers

Numbers of the form n0n are for standard temperatures (ambient to 200 °C).


Numbers of the form n3n are for high temperatures (100 to 350 °C).

## Custom Plumbing (diagram required), Option 700 or 730

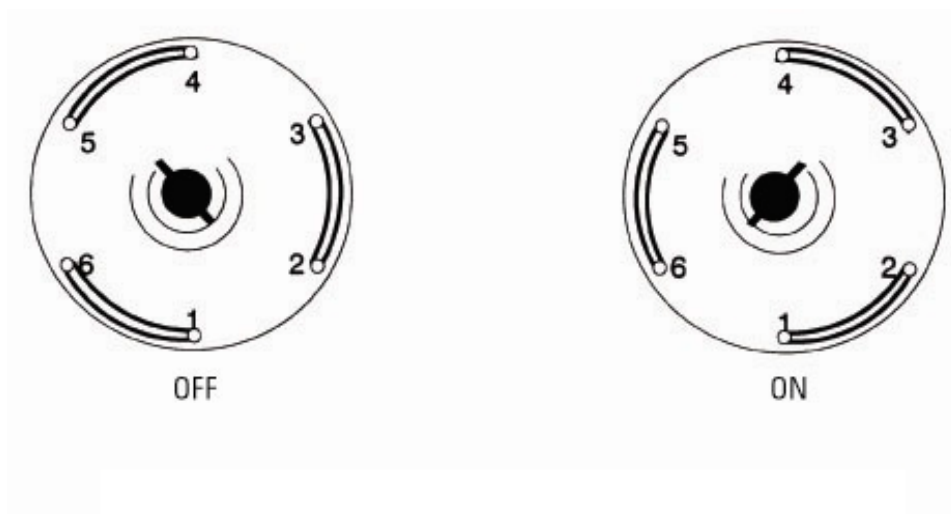

Custom Plumbing (Diagram Required), Option 700 or 730


## Gas Sampling Option, Option 701 or 731

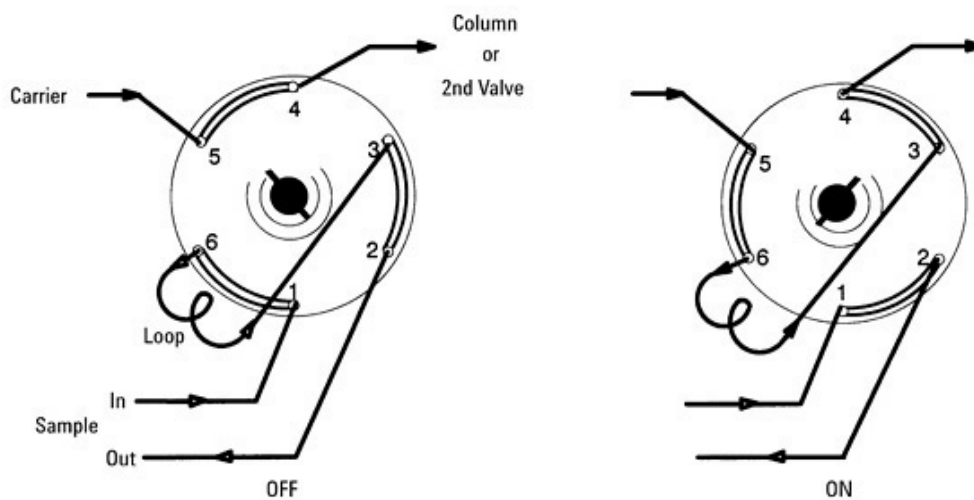

Gas Sampling Option, Option 701 or 731




### Column Isolation, Option 702 or 732

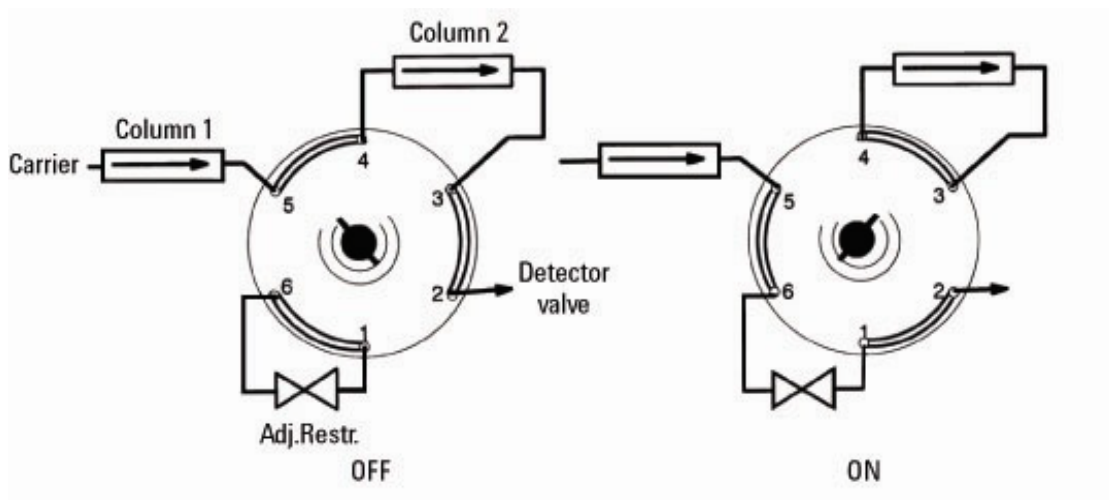

Column Isolation, Option 702 or 732

High temperature needle valve supplied with option 732


### Two Stream Selection (Requires Gas Sampling), Option 703 or 733

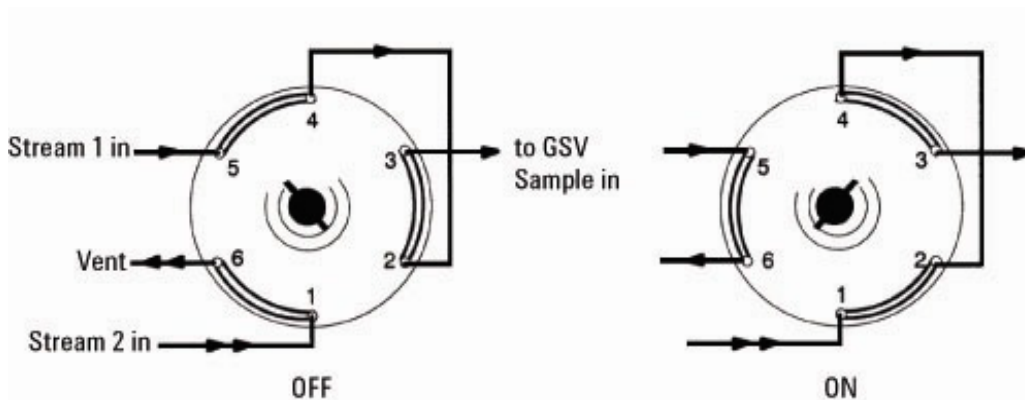

Two Stream Selection (Requires Gas Sampling), Option 703 or 753




### Backflush to Detector, Option 704 or 734

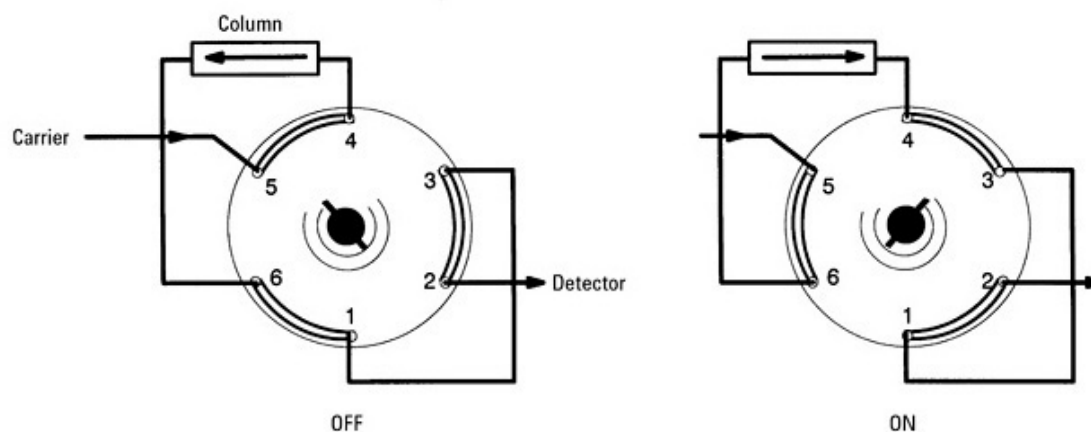


Backflush to Detector, Option 704 or 734

### Backflush a Precolumn to Vent, Option 705 or 735

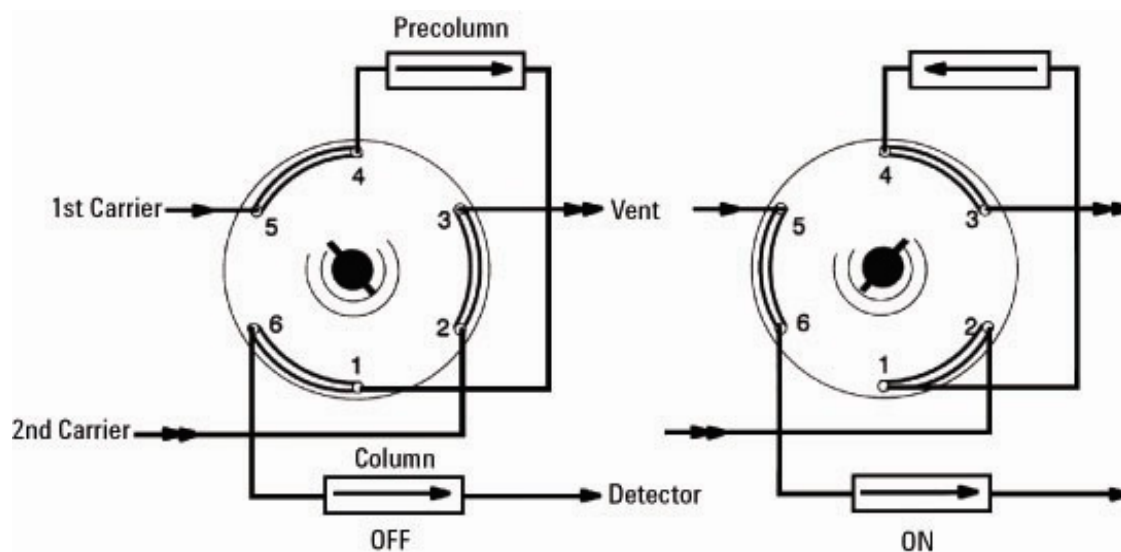


Backflush a Precolumn to Vent, Option 705 or 735

### Column Selection (Unused Column Isolated), Option 706 or 736

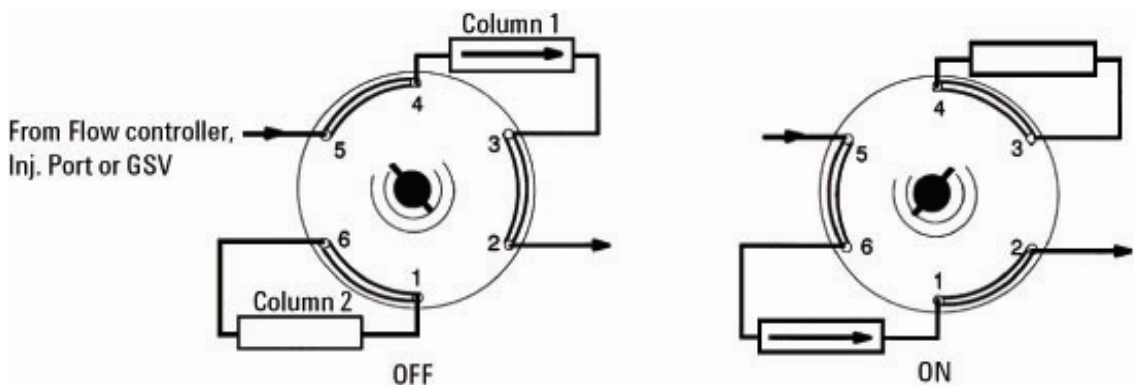


Column Selection (Unused Column Isolated), Option 706 or 736

### Sequence Reverse, Option 707 or 737

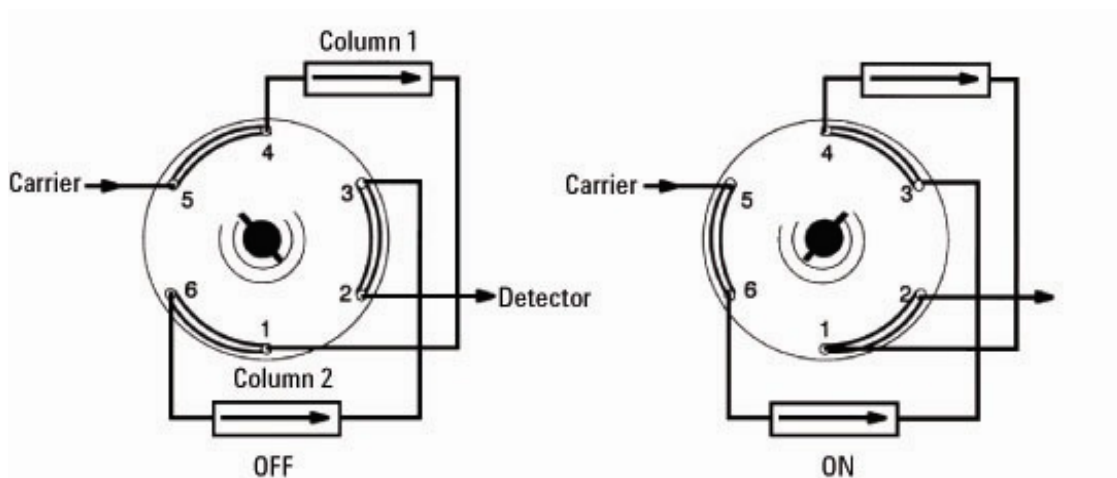


Sequence Reverse, Option 707 or 737

### Sequence Reverse with Backflush of Column 1, Option 708 or 738

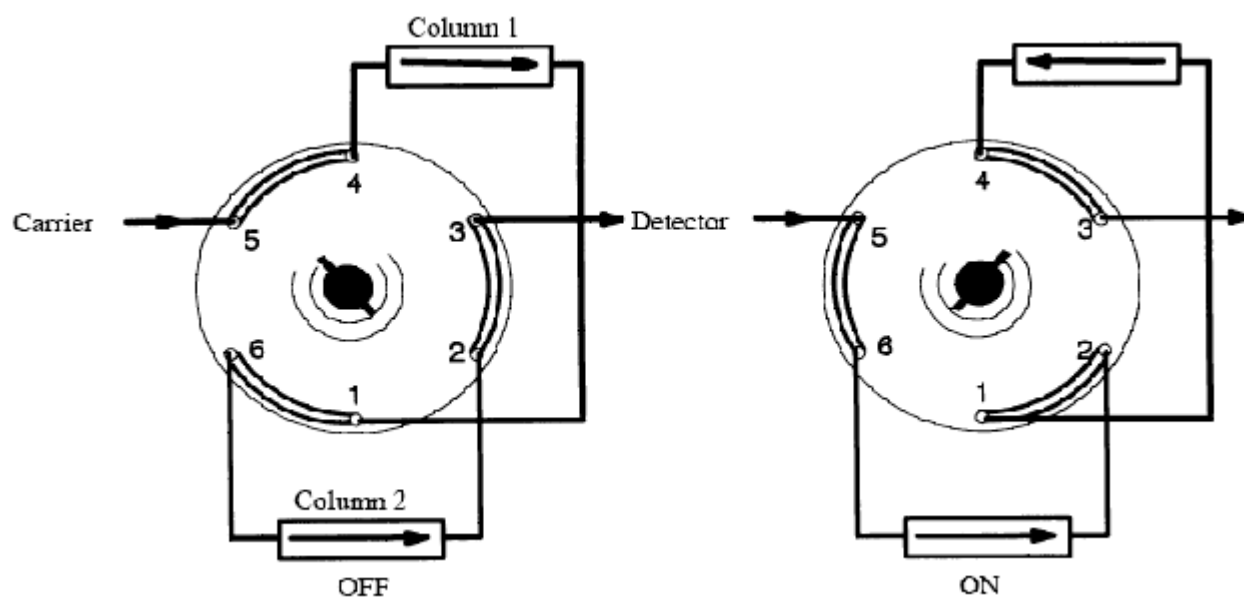


### Custom Plumbing (Diagram Required), Option 800 or 830

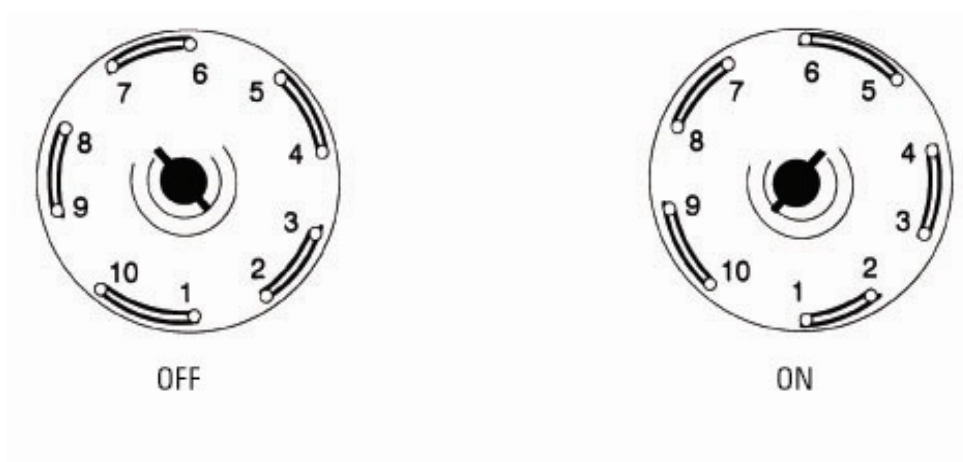


Custom Plumbing (Diagram Required)), Option 800 or 830

### Gas Sampling with Backflush of Precolumn to Vent, Option 801 or 831

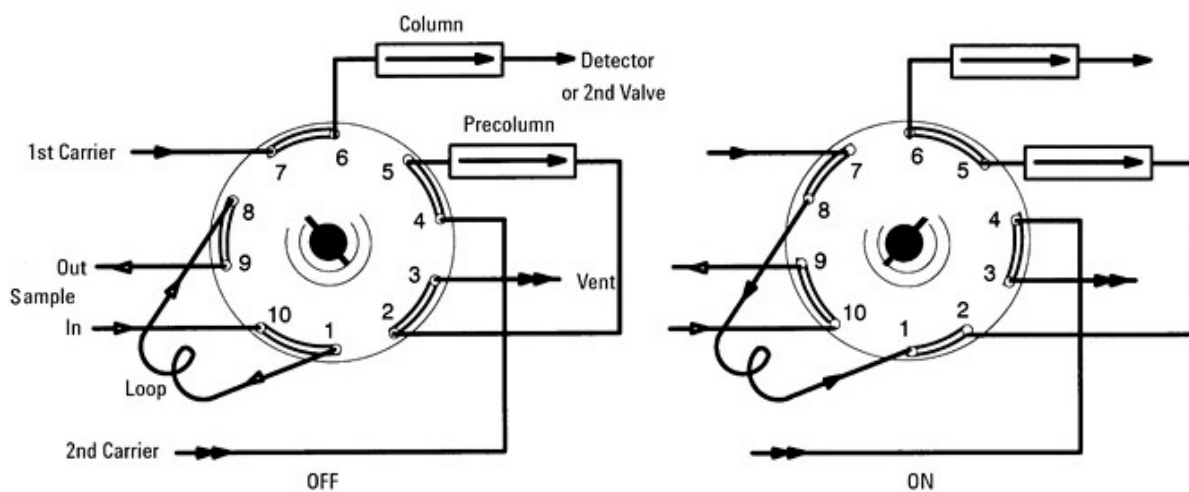


Gas Sampling with Backflush of Precolumn to Vent, Option 801 or 831

### Gas Sampling with Backflush to Detector, Option 802 or 832

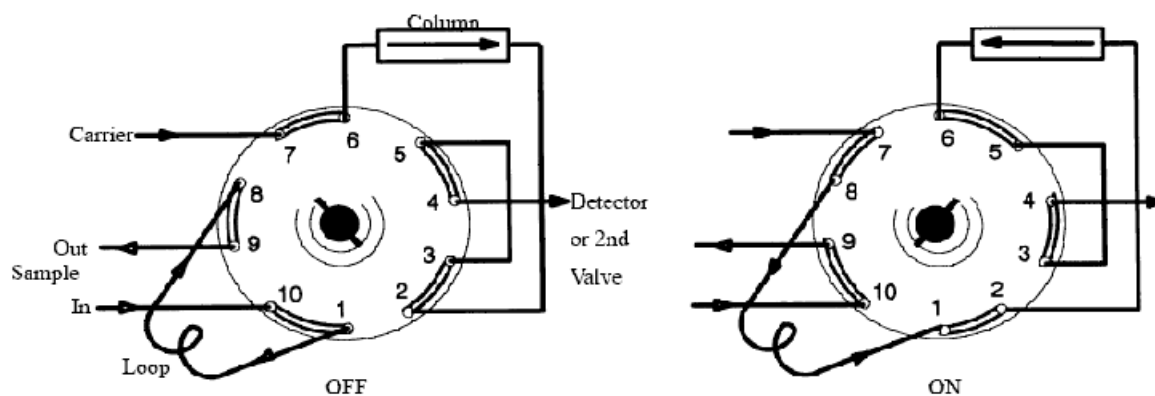


Gas Sampling with Backflush to Detector, Option 802 or 832


### Gas Sampling of Alternate Streams, Option 803 or 833

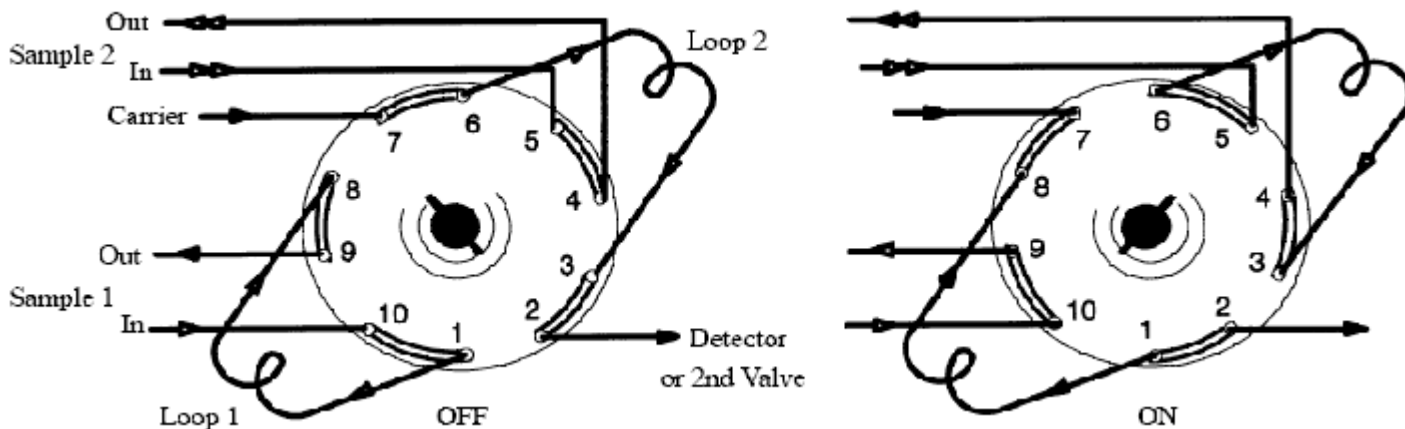

### Gas Sampling with Sequence Reverse, Option 804 or 834

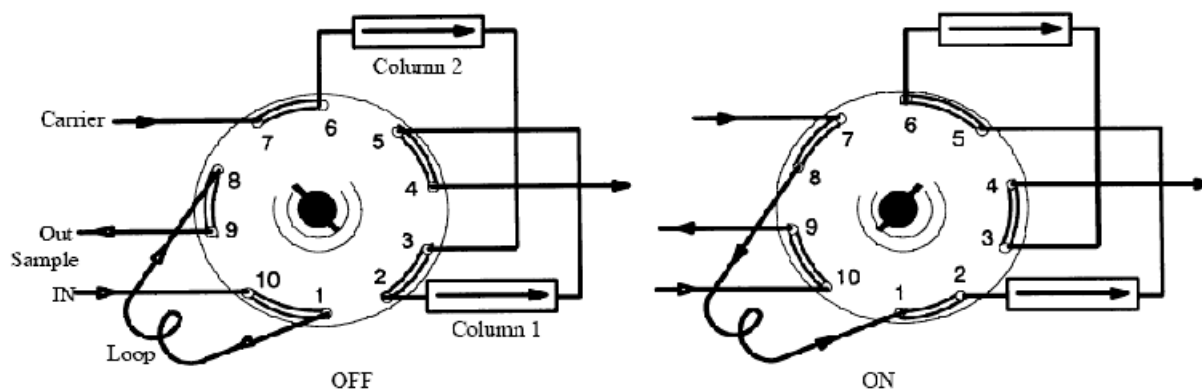



### Gas Sampling with Sequence Reverse and Backflush of Column 1, Option 805 or 835

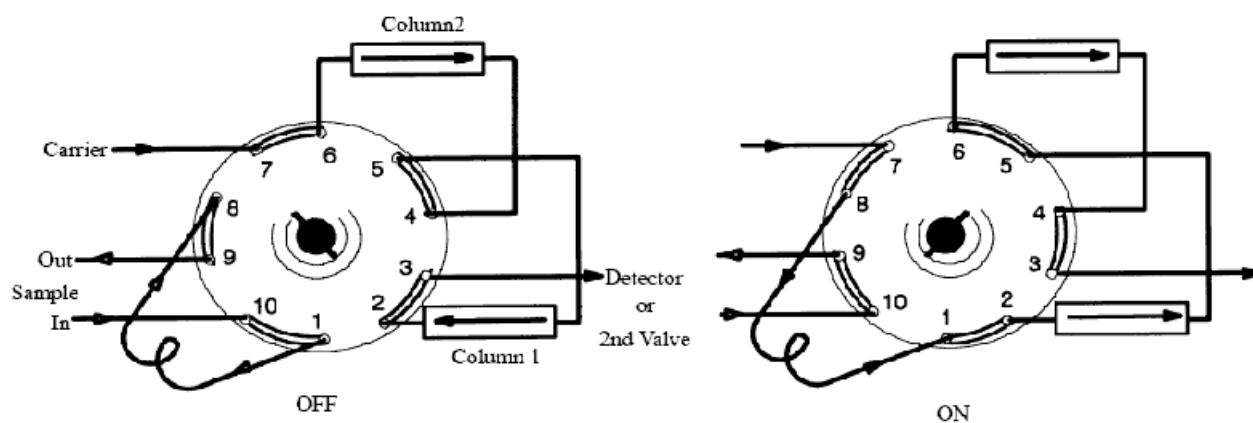

Gas Sampling with Sequence Reverse and Backflush of Column 1, Option 805 or 835

### Column Selection with Backflush to Vent, Option 806 or 836

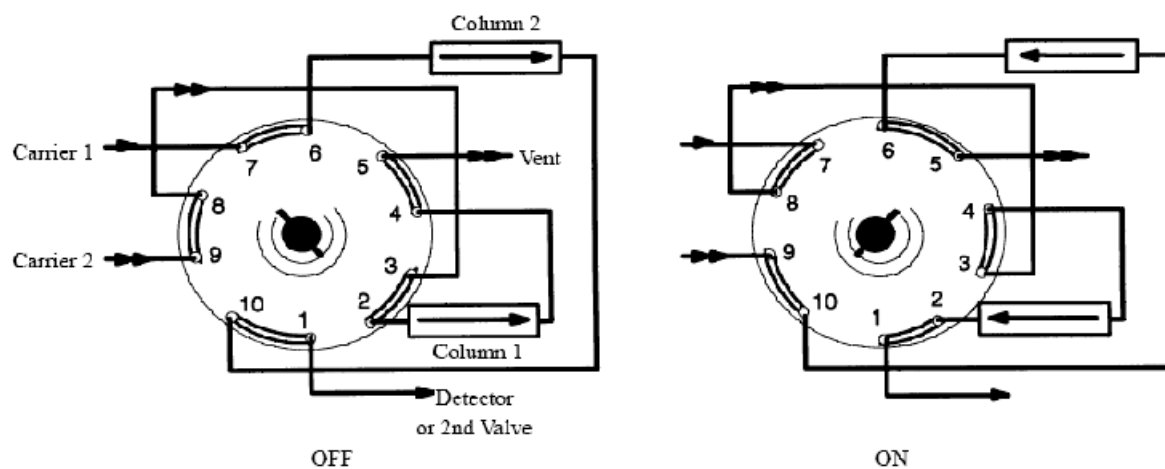

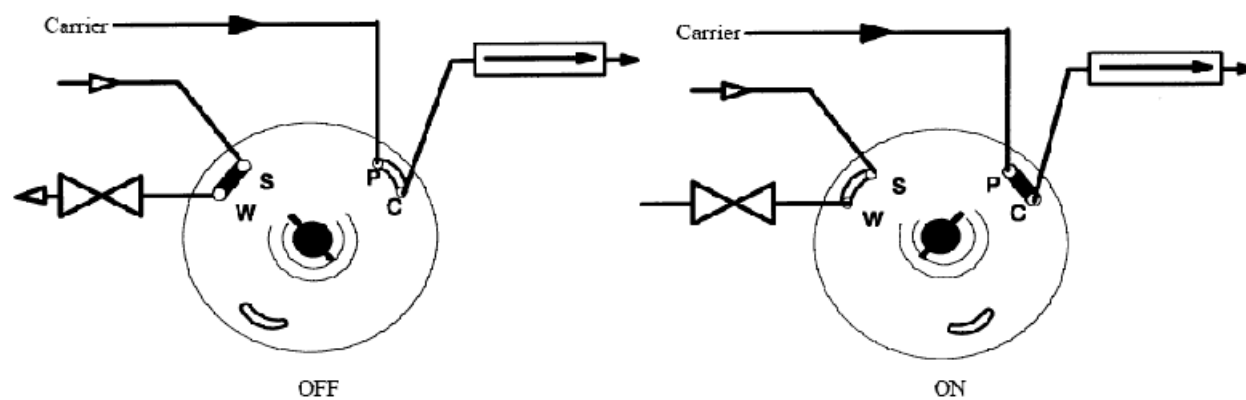**Liquid Sampling, Options 850, 852, 854, and 856**

## Troubleshooting

Most of the problems associated with sampling valves are related to peak broadening in transfer lines and inlets, sample adsorption by the valve or transfer lines, leaks, and perturbations in the baseline.

### Chromatographic symptoms

Troubleshooting valves and their related plumbing is primarily a matter of systematic checking and verification of unimpaired mechanical operation of any moving part. This requires an understanding of how the valve functions internally and how the plumbing is configured. A plumbing diagram is essential for effective troubleshooting.

The following "symptom-cause" list gives the most commonly encountered problems found with valves and their solution.

Symptom	Possible cause	Solution
Lost peaks (degradation)	Valve or transfer lines too hot	Reduce temperature 50 °C, reevaluate
	Transfer line activity	Use nickel or Hastelloy tubing
Lost or tailing peaks	Valve or transfer line too cold	Increase temperatures 50 °C, reevaluate
Baseline perturbation	Slow valve rotation	Increase actuator pressure
	Rotor distorted	Replace rotor
	Sample/column pressure too different	Add back-pressure regulator to sample drain
Peak tailing broad peaks	Column overload	Use smaller sample loop Increase split flow
	Flow too slow	Increase column flow Increase split flow
	System voids	Check connections Reduce volume of connecting tubing



## Loss of sensitivity or excessive drift

Several possible causes exist for overall deterioration of the chromatogram.

- Contamination in the valve requires a thorough cleaning.
- Internal leakage requires a complete disassembly and inspection of the mating surfaces.
- Poor temperature control may require a full check of electronic and thermal components.
- Lack of proper conditioning techniques, columns, etc.
- Failure or deterioration of other components (columns, detectors, etc.).

## Loss of peaks in specific areas of the chromatogram

Entire sections of chromatographic data can be lost due to a valve that does not rotate or one that rotates improperly. Other than obvious component failures (solenoid, actuator, etc.), improper adjustments and misalignments cause most problems.

- Check that adequate air (about 482 kPa or 70 psi) is supplied.
- Check the valve. Is it rotating?
- If the valve rotates, check for proper alignment of the actuator, mechanical binding or slippage of connecting parts.
- Check for blocked flow paths with valve in both positions.

## Extraneous peaks

Air peaks are sometimes seen in a chromatogram when leakage occurs because the valve rotor does not seal properly. These leaks may not be detectable using the soap-bubble method.

If a leak is suspected but cannot be located with soap bubbles, a pressure check will determine definitely if a leak exists. Extraneous peaks can occur due to contamination or improper conditioning of the valve. If leaks are not apparent, clean or condition the valve.

Other causes, totally unrelated to the valve, may produce similar symptoms. Impure carrier gas (i.e., containing water) can cause extraneous peaks.

## Peak broadening and tailing

Voids in the flow system (valve and connecting tubing) cause tailing and peak broadening. Use inlets and liners with small internal diameters and connect the valve to the inlet or column with short lengths of connecting tubing of narrow inner diameter.

If early-eluting peaks are too broad, stationary phase or thermal focusing effects should be used with packed-column ports or increased split flows when capillary split inlets are used. Inlets should be equipped with narrow inner diameter liners, and narrow-bore connecting tubing should be used between the valve and inlet.

## Baseline shifts

Baseline perturbations are caused by changes in column flow as the valve is rotated and as the sample loop equilibrates to system pressure. Slow valve rotation momentarily stops carrier gas flow and, when the valve stops rotating, a sudden increase in flow occurs which slowly returns to the set point. Check actuator pressure (usually 40 to 75 psi), valve rotor tension, and valve temperature to ensure that the valve rotates as quickly as possible. A restrictor or backpressure regulator can be added to the sample vent line to maintain the sample loop at system pressure. This will reduce the time it takes for the flow to stabilize after the valve is switched.

## Baseline upsets

Frequently, baseline upsets are seen on chromatograms when valves are switched. These upsets are caused by pressure changes within the system, injections of large volume samples, or by changing the amount of restriction in the flow path. These upsets will become more of a problem when high sensitivity is required. Addition of a fixed restriction downstream from the valve may help minimize the upset. Changes in column length may also help reduce the upsets.

Fixed restrictors are used immediately before flame detectors to prevent flameout and are used in some instances to prevent pressure surges from damaging TCD filaments. An adjustable restrictor (needle valve) can also be used where a matched restriction is desired but not for preventing pressure or flow surges.

Often confused with baseline upsets, an offset is a shift in the baseline that does not return quickly to the original level.

Baseline offsets may be caused by air leaks but more commonly are due to a change in gas purity or flow rate in the detector. Poor carrier gas or improperly conditioned filters and traps should be suspected whenever offsets occur.

## Variation in peak area and retention time

The amount of sample contained in the loop and, therefore, the amount injected onto the column is affected by loop pressure and temperature. Variations in pressure and temperature lead to variability in peak areas. Flow restrictors or back-pressure regulators help to maintain constant loop pressure, and valve boxes help maintain temperature.

Leaks can occur in the valve itself or at any of the connecting points with transfer lines. Leaks usually cause area irreproducibility, retention times changes, and increases in the area of air peaks (with thermal conductivity detectors). Leaks in rotors can sometimes be fixed by tightening the nuts holding the rotor in the valve body. Leaks in connections are usually found with an electronic leak detector or with a liquid leak detection fluid (e.g., Snoop).

## Pressure check


Leak checking the plumbing involved in a valve system must be done carefully and methodically. The pressure check method below will indicate, but sometimes not isolate, a leak in the flow path. Since this method does not necessarily isolate the leak, other leak check methods may be needed to locate the leak specifically.

### NOTE

Each valve in a system has two flow paths, ON and OFF. A leak sometimes occurs in only one of these two positions. Check both.

- 1 Disconnect the detector from the valve system.
- 2 Cap the valve system at its outlet and pressurize to 689 kPa (100 psi). Allow 2 to 5 minutes for pressure to equilibrate. If your instrument has flow control, it should read zero flow.
- 3 Turn off the gas supply at the source.
- 4 Generally, the pressure will drop quickly for approximately 30 to 60 seconds, then stabilize. After this initial pressure drop, the gauge should not indicate more than a 7 to 14 kPa (1 to 2 psi) drop during a 10 minute period.

- 5 If no leak is indicated, actuate all valves and repeat steps 2 to 4.
- 6 If a leak does show up, try to pinpoint the source using a soap bubble meter. Do not assume that the leak exists only at the valve. Often plumbing connections such as unions or bulkhead fittings are at fault. See “Valve Box” should it become necessary to expose the valve system.
- 7 If the leak cannot be found easily, divide the system in half and repeat the pressure check. Continue dividing in halves, and pressure check until the leak is isolated.

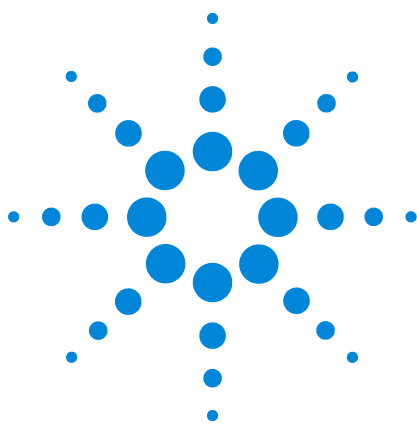

## 11 Capillary flow technology (CFT) devices

CFT overview	458
Ultimate union	460



## CFT overview

CFT devices consist of a plate and tubing. Ultimate unions are often used with CFT devices and are described here.

A CFT plate consists of two thin metal plates that are bonded together. Gas passages are left between the plates, and tubing fittings are brazed to the surface. Finally, the gas passages are deactivated. This results in a sturdy and inert unit with a very low thermal mass.

Accessories contain the plate, plus installation/use instructions and some related parts (ferrules, nuts, tubing, etc.) needed for use.

The following accessories and plates are available at the time of writing.

Description	Accessory	Plate only
Deans switch	G2855B	G2855-61500
Splitter, 2-way with make-up gas	G3180B	G3180-61500
Splitter, 2-way without make-up gas	G3181B	G3181-60500
Splitter, 3-way with make-up gas	G3183B	G3183-60500
Quick Swap for MSD	G3185B	G3185-60065

Tubing connects plates with one another and with other parts of the analytical path.

Description	Part Number
CPM-1/16-inch x 0.010 -inch id tube, 50cm	G1580-80060
CPM-1/16-inch x 0.20-inch id tube, 50cm	G1580-80062
Fused silica, deactivated, 0.20 mm x 5 m	160-2205-5
Fused silica, deactivated, 0.25 mm x 5 m	160-2255-5
Fused silica, deactivated, 0.18 mm x 10 m	160-2615-10
Fused silica, deactivated, 0.10 mm x 5 m	160-2635-5

Miscellaneous parts are used in making column connections.

Description	Part number
Oven bracket kit	G2855-60140
Reducing union, 1/8-inch to 1/16-inch	0100-0241
Restrictor, SS, 1/16-inch od x 0.01-inch id x 1 m	0100-2354
Union, SS, 1/16-inch	0100-0124
Screw, T-20, M4 x 8 mm	0515-2755
Column storage fitting	G2855-20590
Wire, SS, 0.015-inch x 40 mm	G2855-60593
Plug for CFT fittings or unions	G2855-60570
QuickSwap transfer line locking nut	G3185-20501
Deans switch supplies kit	G2855-60150

## Ultimate union

This union combines a stainless steel body, inert metal ferrules, and very low dead volume to provide durable connections of capillary columns, retention gaps, pre-columns, and other capillary devices.



### Tools and kits

The kits do not contain ferrules. The appropriate size must be ordered separately.

Part number	Description
G3182-61580	Union kit, deactivated
G3182-61581	Union kit, non-deactivated
G2855-20530	Internal nut
G2855-20555	Swaging nut
8710-0510	5/16 x 1/4 open end wrench (need 2)
5181-8836	Column cutter




## Metal ferrules and nuts

Identify the internal diameter of your column and select the part number for the appropriate sized nut and ferrule. Each packet contains an internal nut.


Part number	Description
5184-3569	Metal ferrules, 0.1- to 0.25-mm columns, pkg of 10
5184-3570	Metal ferrules, 0.32-mm columns, pkg of 10
5188-2789	Metal ferrules, 0.53-mm columns, pkg of 10

## Assembling nut, ferrule and swaging nut


- 1 Thread the column end through the internal nut and the ferrule leaving approximately 3 cm of column protruding beyond the ferrule. Thread the swaging nut onto the internal nut with the column protruding.

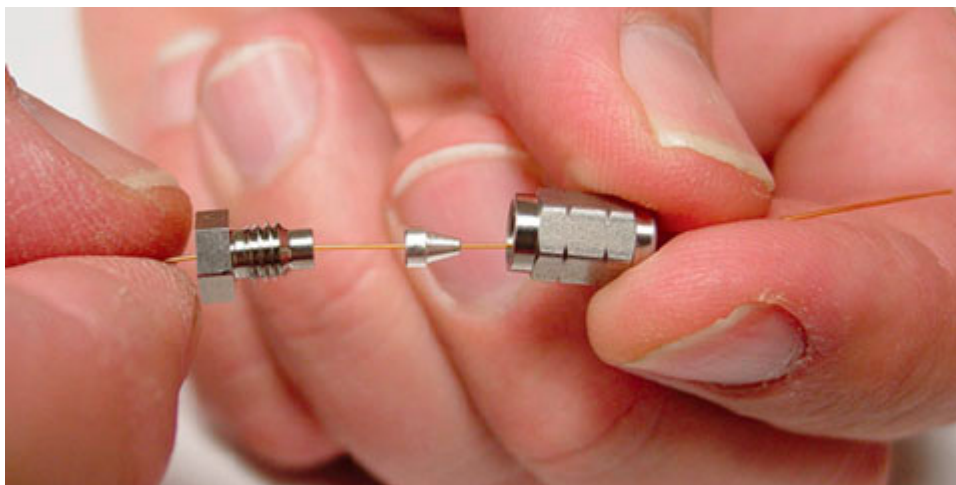


- 2 Tighten the two nuts together a little at a time, occasionally checking to see if the ferrule is grabbing the column. When the ferrule just starts to grab, notice the position of the nut and then tighten the nuts by turning 60 degrees (one flat).

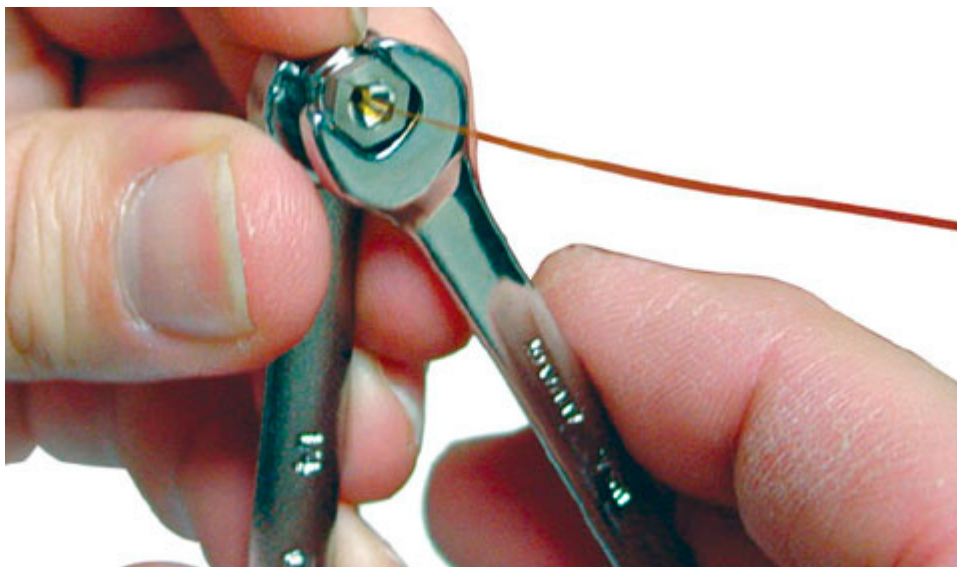


3 Remove the swaging nut.

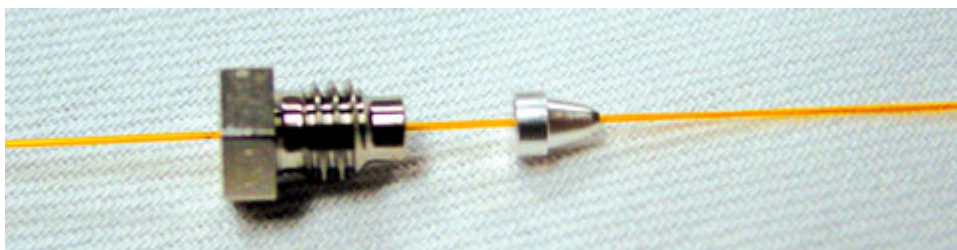

4 Rest the ceramic column cutter on the small end face of the ferrule at a 45-degree angle.

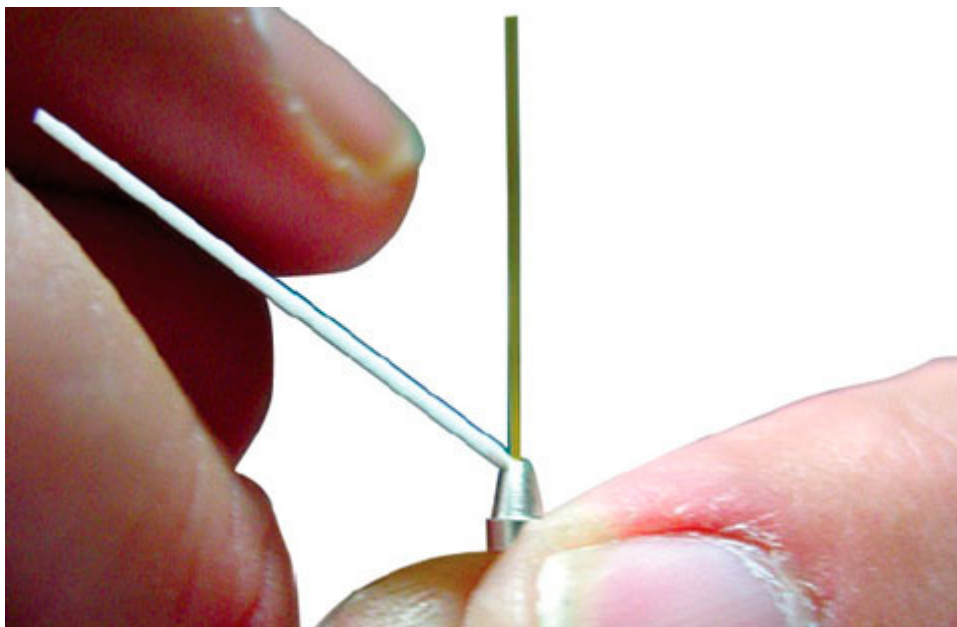


- 5 Scribe the column, then break off the end. It should break off with approximately 0.3 mm of column extending beyond the ferrule. It is important that you cut the column to a maximum of 0.5 mm beyond the swaged ferrule end. The end of the column need not be perfectly square, but should not have cracks that extend under the ferrule.

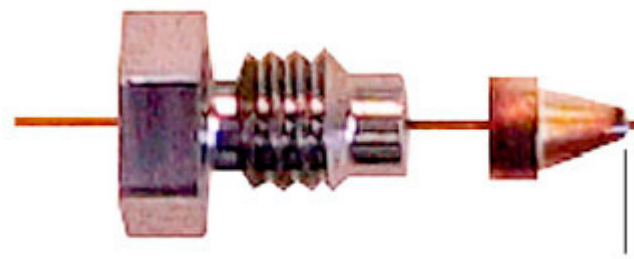

## Connecting and disconnecting columns

### Connecting columns

Connect columns to the union with internal nuts and preswaged SilTite ferrules. Finger-tighten the nuts. Further tighten with a wrench only 15 to 20 degrees.

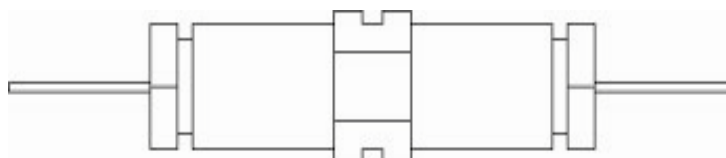



### Disconnecting columns

Loosen and remove the internal nut. If the column and ferrule do not come free, insert a pointed object (pen, paper clip) into the ferrule release hole and press firmly. You will hear a click as the ferrule releases.

## Mounting the union

The ultimate union is probably too heavy to be supported only by the capillary columns. Use the oven wall clip to support the weight of the union by attaching it to an opening in the oven wall.

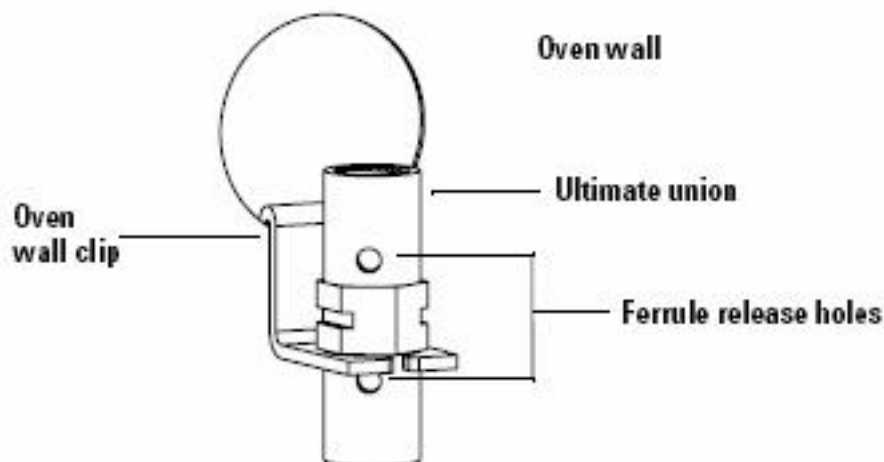

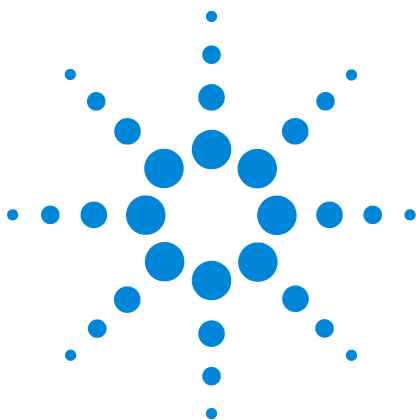

## 12 Electrical

Power options	466
Converting the power option	469
Remote start/stop connection	470
Remote control	471
Measuring inlet and detector heater and PRT resistance	472
Temperature sensor resistance	475
Oven heater coil resistance	476
Setting the instrument power configuration	477
Line voltage configuration plug	478
Ceramic and glass fuses	479
AC power board schematic	480
Oven shroud	481



## Power options

There are 8 power option kits. These contain the configuration plug, oven shroud, label, and power connector for each option. If you are converting from one power option to another you can order the kit or order the unique individual parts.

The power cord must be ordered separately. See “Power cords”.

### Regular oven

There are 4 power options for the regular oven:

- G1530-63800 Regular power kit 120V
- G1530-63810 Regular power kit 220V
- G1530-63820 Regular power kit 230V
- G1530-63830 Regular power kit 240V

Kits	120VAC G1530-63800	220VAC G1530-63810	230VAC G1530-63820	240VAC G1530-63830
Configuration plug	G1530-60690	G1530-60710	G1530-60720	G1530-60730
Oven shroud	G1530-61610	G1530-61230	G1530-61670	G1530-61240
Label	G1530-90928	G1530-90921	G1530-90922	G1530-90923
Power receptacle	G1530-61550	G1530-61560	G1530-61560	G1530-61560
Clamp	1400-1663			

Common to all regular power options:

- 0515-1084 Screw, M3x0.5x12mm T20, quantity 2 for power input assembly
- 0515-2711 Screw, M4X0.7x12mm T20, coated, quantity 4 for shroud
- 0535-0043 Nut with lock washer, M4X0.7x3.2mm 7mm for grounding wire
- 05890-00460 Insulation plug for heater/sensor tube
- G3430-00104 Heater Wire Insulation Sleeving

## Fast oven

There are 4 power options for the fast oven:

- G1530-63840 Fast power kit 200V (Japan)
- G1530-63850 Fast power kit 208 or 220V
- G1530-63860 Fast power kit 230V
- G1530-63870 Fast power kit 240V

<b>Kits</b>	<b>200VAC G1530-63840</b>	<b>208 220VAC G1530-63850</b>	<b>230VAC G1530-63860</b>	<b>240VAC G1530-63870</b>
<b>Configuration plug</b>	G1530-60700	G1530-60710	G1530-60720	G1530-60730
<b>Oven shroud</b>	G1530-61620	G1530-61630	G1530-61650	G1530-61640
<b>Label</b>	G1530-90924	G1530-90925	G1530-90926	G1530-90927

Common to all fast power options:

- 0515-1084 Screw, M3x0.5x12mm T20, quantity 2 for power input assembly
- 0515-2711 Screw, M4X0.7x12mm T20, coated, quantity 4 for shroud
- G1530-61550 Power input assembly
- 1400-1663 Clamp, 1.25x1.0 inch
- 0535-0043 Nut with lock washer, M4X0.7x3.2mm 7mm for grounding wire
- 05890-00460 Insulation plug for heater/sensor tube
- G3430-00104 Heater Wire Insulation Sleeving

Table of GC voltage by country

Country	Kit	Voltage	Oven type
Australia, 10 amp	G1530-63830	240 V	Regular
Australia, South Africa	G1530-63870	240 V	Fast
China	G1530-63810	220 V	Regular
China, Hong Kong	G1530-63850	220 V	Fast
Continental Europe, dual phase	G1530-63860	230 V	Fast
Continental Europe, single phase: Russia	G1530-63850	220 V	Fast
Denmark, Switzerland, 10 amp	G1530-63820	230 V	Regular
India, Denmark, Switzerland, 16 amp	G1530-63860	230 V	Fast
Israel	G1530-63850	220 V	Fast
Japan	G1530-63840	200 V	Fast
United Kingdom, Ireland	G1530-63870	240 V	Fast
USA	G1530-63800	120 V	Regular
USA	G1530-63850	208 V	Fast
USA	G1530-63870	240 V	Fast



## Converting the power option

This procedure outlines the steps you must complete to prepare a GC for use in another country or to configure for a different oven heating speed.

The basic requirements for use of the fast heating oven are:

- The electric service must be capable of providing  $\geq 200$  V at  $\geq 15$  amperes.
- In the United States, the electric service must be 208 V or 240 V.

### Procedure

- 1 Select the appropriate power option kit. See [“Power options”](#). There are 8 power option kits. These contain the configuration plug, oven shroud, label, and power connector for each option. If you are converting from one power option to another you can order the kit or order the unique individual parts. Compare the parts listings for the present and new power option to see what must be replaced.
- 2 Select the appropriate power cord. See [“Power cords”](#).
- 3 Cool down the oven.
- 4 Turn off the GC main power switch and disconnect its power cord.

#### **WARNING**


**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

#### **CAUTION**

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.

- 5 If necessary, replace the oven shroud. See [“Replacing the oven shroud, oven fan, and oven fan motor”](#).
- 6 Replace the oven configuration plug. See [“Line voltage configuration plug”](#).
- 7 If necessary, replace the power cord connector.
- 8 Attach the new label.
- 9 Power on the unit. Watch for any errors after self-test.

## Remote start/stop connection

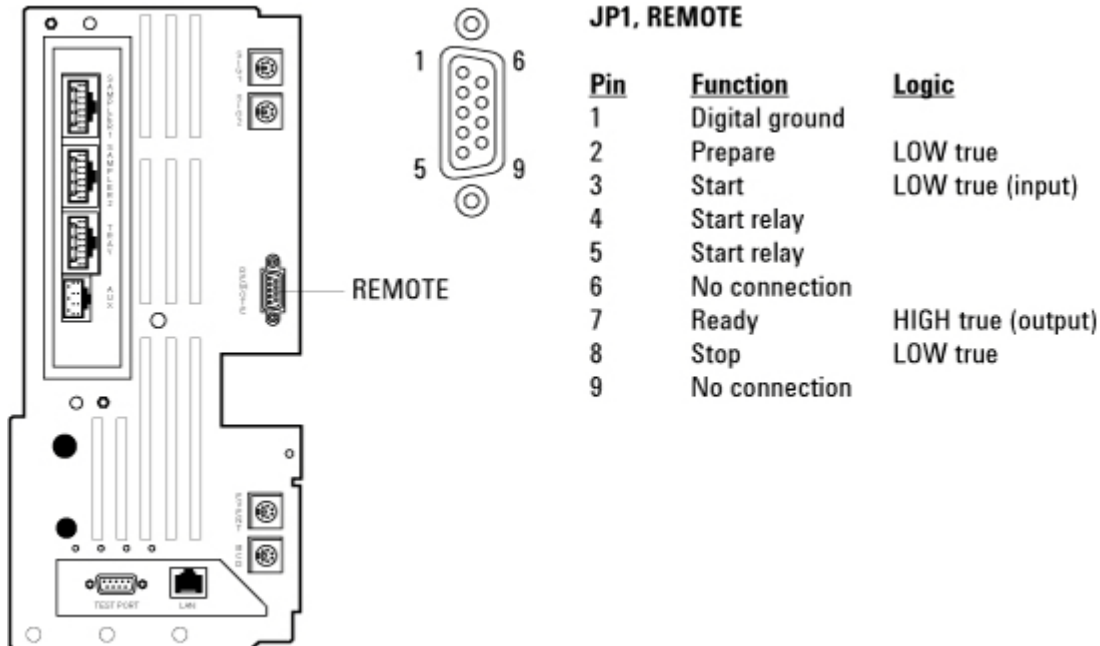

**Prepare (low)**—Request to prepare for analysis. Receiver is any module performing pre-analysis activities.

**Start (low)**—Request to start run/timetable. Receiver is any module performing runtime-controlled activities.

**Ready (high)**—System is ready for next analysis. Receiver is any sequence controller.

**Stop (low)**—Request to reach system ready state as soon as possible (for example, stop run, abort or finish, and stop injection). Receiver is any module performing runtime-controlled activities.

## Remote control

Remote control allows easy connection between single instruments or systems to ensure coordinated analysis with simple coupling requirements.

**Start Relay**—A 120 millisecond contact closure


Control of analysis is maintained by signal readiness READY for next analysis, followed by START of run and optional STOP of run triggered on the respective lines. In addition, PREPARE and START REQUEST may be issued.

The signal levels are defined as standard TTL levels (0 V is logic true, +5 V is logic false).

Input Load  $\geq 2.2$  kOhm against +5 V

Output type is open collector.

### Remote control timing

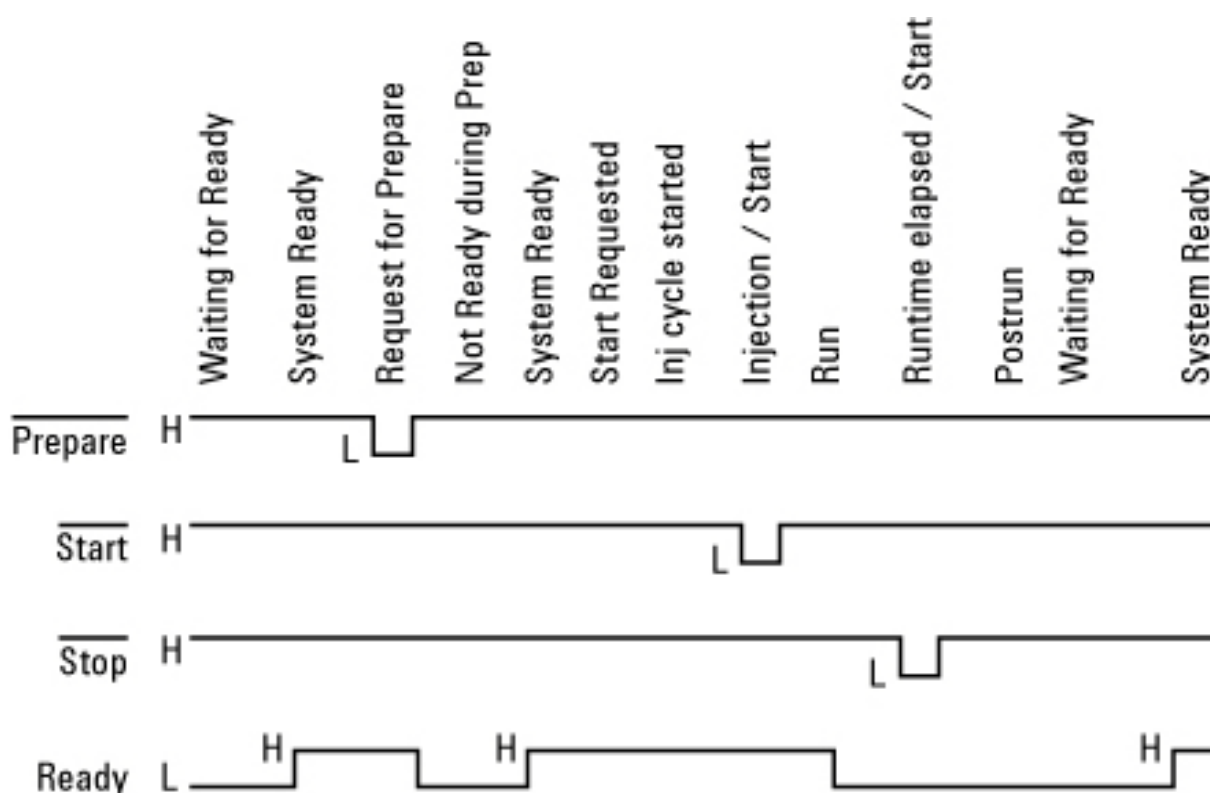

## Measuring inlet and detector heater and PRT resistance

If you suspect that an inlet or detector heater or PRT is defective, you can check its resistance and compare against the expected nominal values.

To check the resistance of an inlet or detector heater, or to check the resistance of an inlet or detector PRT:

**To measure the inlet or detector heater or PRT resistance:**

- 1 Turn off the GC main power switch and disconnect its power cord.

---

**WARNING**

**Hazardous voltages are present in the mainframe when the GC power cord is connected. Avoid a potentially dangerous shock hazard by disconnecting the power cord before removing any GC panels.**

---


**CAUTION**

Components can be damaged by static electricity: be sure to wear an ESD strap grounded to the GC chassis while performing this procedure.


---

- 2 Remove the GC panels required to access the zone's heater/sensor connector to be tested.
  - For inlets, remove the inlet cover and GC left side panel
  - For detectors, remove the electronics cover and GC right side panel


### 3 Remove the connector from the GC chassis.

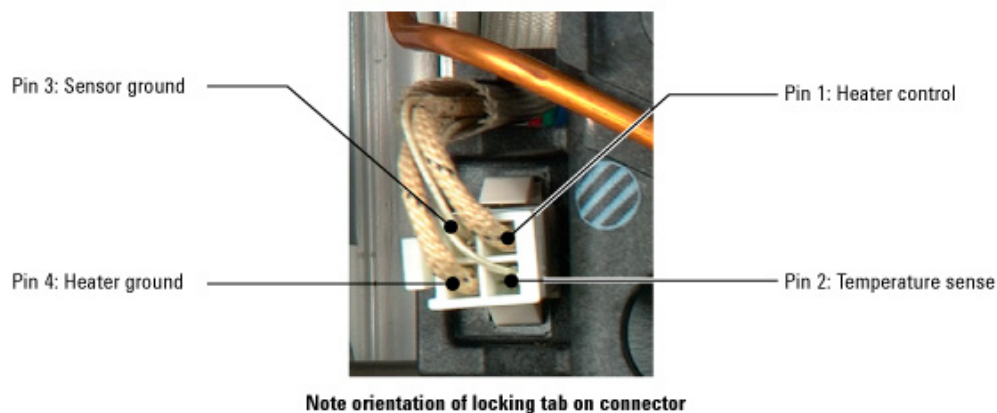

Inlet heater/sensor connection

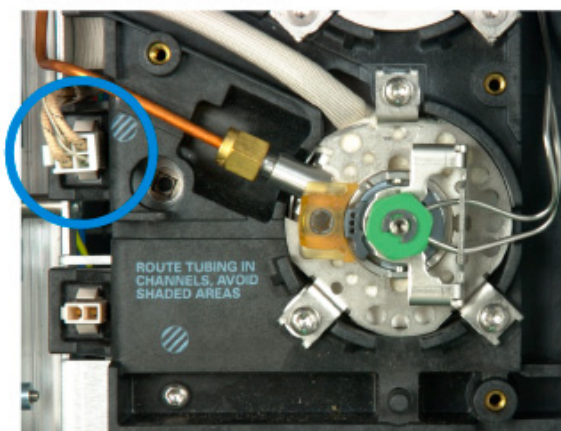

Detector heater/sensor connections

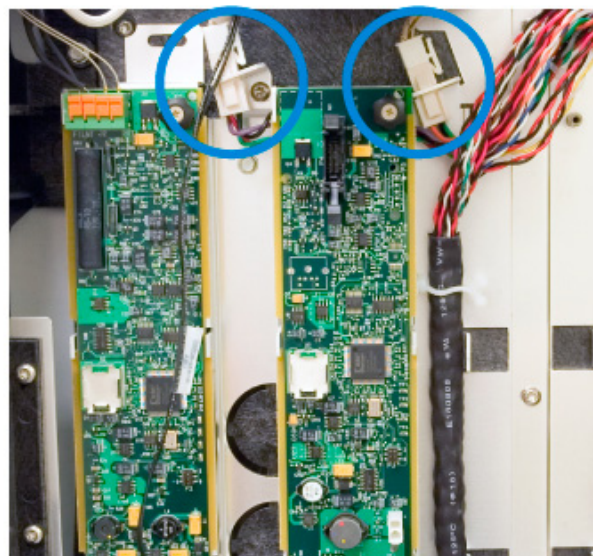

### 4 Use an Ohmmeter to measure the resistance across the connector pins of the heater/sensor cable.

- For heated zones the PRT is connected to the white wires on pins 2 and 3 of the connector.
- With the heater/sensor connector removed from the GC, there should be no connection to ground on any of the four pins. A reading of  $\sim 0$  indicates a possible shorted sensor or heater.

### 5 Compare the results against the resistance values in the tables below.

## Sensor resistance

The approximate resistance of a temperature sensor is:

$$R = 100 + (.35 \times t)$$

where R is resistance in ohms and t is sensor temperature in °C.

**Table 29** Temperature Sensor Resistance by Heater Temperature

°C	+0°	+10°	+20°	+30°	+40°	+50°	+60°	+70°	+80°	+90°	+100°
0°	100.00	103.90	107.79	111.67	115.54	119.40	123.24	127.07	130.89	134.70	138.50
100°	138.50	142.28	146.06	149.82	153.57	157.32	161.04	164.76	168.47	172.16	175.84
200°	175.84	179.51	183.17	186.82	190.46	194.08	197.70	201.30	204.88	208.46	212.03
300°	212.03	215.58	219.13	222.66	226.18	229.69	233.19	236.67	240.15	243.61	247.06
400°	247.06	250.50	253.34	257.34	260.75	264.14	267.52	270.89	274.25	277.60	280.93

Nominal inlet and detector heater/sensor resistances

Inlet or Detector	Heater type	Nominal heater resistance	Sensor type/comments
Split/splitless	70 W	23 Ω	100 Ω PRT
Purged packed	70 W	23 Ω	100 Ω PRT
Cool on-column	60 W	26-27 Ω	100 Ω PRT
Multimode	150 W	~ 9 Ω	Type K Thermocouple Resistance, cold: ~ 13 Ω
PTV	150 W	~ 9 Ω	Type K Thermocouple Resistance, cold: ~ 13 Ω
Volatiles interface	60 W	26-27 Ω	100 Ω PRT
FID	70 W	23 Ω	100 Ω PRT
TCD	70 W	23 Ω	100 Ω PRT
NPD	70 W	23 Ω	100 Ω PRT
uECD	70 W	23 Ω	100 Ω PRT
FPD	60 W (2 each)	26-27 ?	100 ? PRT (2 each)

## Temperature sensor resistance

The approximate resistance of a temperature sensor is:

$$R = 100 + (.35 \times t)$$

where R is resistance in ohms and t is sensor temperature in °C.

**Table 30** Temperature Sensor Resistance by Heater Temperature

°C	+0°	+10°	+20°	+30°	+40°	+50°	+60°	+70°	+80°	+90°	+100°
0°	100.00	103.90	107.79	111.67	115.54	119.40	123.24	127.07	130.89	134.70	138.50
100°	138.50	142.28	146.06	149.82	153.57	157.32	161.04	164.76	168.47	172.16	175.84
200°	175.84	179.51	183.17	186.82	190.46	194.08	197.70	201.30	204.88	208.46	212.03
300°	212.03	215.58	219.13	222.66	226.18	229.69	233.19	236.67	240.15	243.61	247.06
400°	247.06	250.50	253.34	257.34	260.75	264.14	267.52	270.89	274.25	277.60	280.93

## Oven heater coil resistance

If you believe that your heater coil is cracked or otherwise damaged and has caused an open circuit, you can check it by measuring its resistance.

### To measure the resistance:

- 1 Turn the instrument power off.
- 2 Disconnect the oven heater leads (P3, P4) from the AC power board.
- 3 Use an ohmmeter to measure resistance at the **connectors**.

**Acceptable resistance ranges** (in ohms) are given below. Acceptable resistances range from the nominal value for a new, cold heater to +5% from the nominal value.

### NOTE

Resistance goes up approximately +3% after heating the coil.

**Table 31** Resistances of the Heater Coil

Nominal cold heater resistances, ohms		
	Standard oven (1600 VA)	Fast-ramp oven (2250 VA)
120 V	9.07 – 9.52	n/a
200 V	n/a	17.78 – 18.7
220 V	n/a	21.51 – 22.6
230 V	33.06 – 34.71	23.51 – 24.7
240 V	n/a	25.60 – 26.9



## Setting the instrument power configuration

There are seven possible line voltage power configurations for the GC.

To change the power configuration for the instrument, you must install the appropriate types of the following components:

- Line voltage configuration plug
- Ceramic fuses on the AC power board
- Oven shroud assembly

These three components are explained in the topics that follow.

**Table 32** Voltage configuration information

Voltage (–10%, +5%)	Frequency (Hz)	Maximum power consumption (VA)	Power line requirement	Oven type
120 V	48-66	2,250	20-amp dedicated receptacle	Regular
200 V	48-66	2,950	15-amp dedicated receptacle	Fast-heating
208 V	48-66	2,950	15-amp dedicated receptacle	Fast-heating
220 V	48-66	2,950	15-amp dedicated receptacle	Fast-heating
230 V	48-66	2,250	10-amp dedicated receptacle	Regular
230 V	48-66	2,950	16-amp dedicated receptacle	Fast-heating
240 V	48-66	2,950	13- or 16-amp dedicated receptacle	Fast-heating

Line voltage configuration plug

There is a different line voltage configuration plug on the AC power board for each power configuration. Each configuration uses three or five jumper wires, each connecting to two different pins on the plug. The 120 VAC configuration uses five jumper wires and all other configurations use three jumper wires.

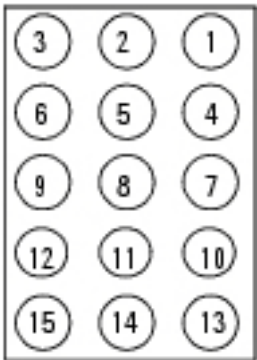



Figure 14 Top view

Jumper locations for power configuration plug					
	120 VAC	200 VAC	208/220 VAC	230 VAC	240 VAC
Transformer	3 « 13	2 « 13	3 « 13	5 « 13	3 « 13
	6 « 15	6 « 12	6 « 12	6 « 9	6 « 9
	9 « 14				
Oven fan	1 « 10	1 « 4	1 « 4	1 « 4	1 « 4
	4 « 8				

## Ceramic and glass fuses

The two sets of fuses on the G3440A AC power board:

- Two ceramic fuses to protect the board from the oven heating circuit (F1, F2)
- Two glass fuses to protect the board from the main supply (F3, F4)

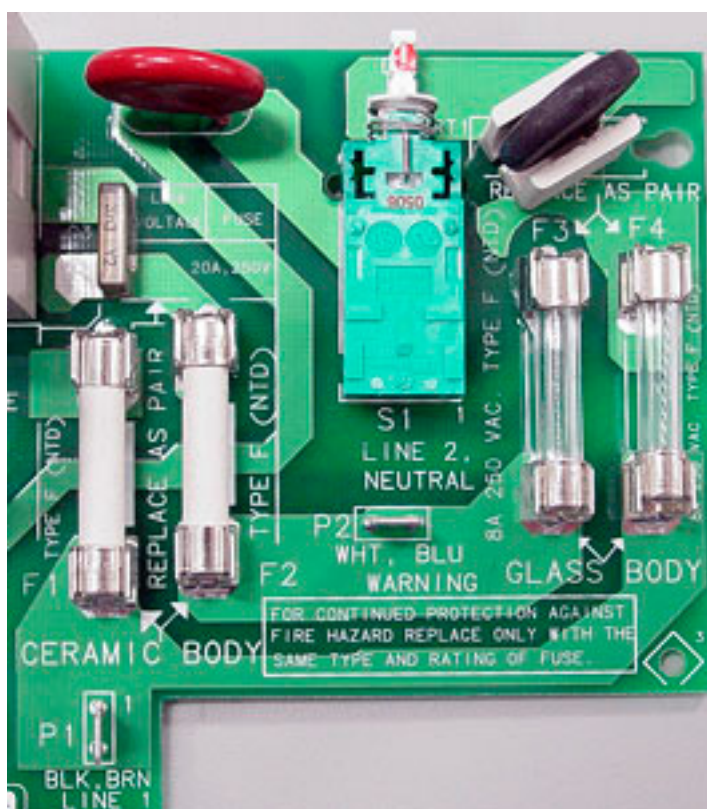


See the table below for the part numbers. When replacing both the glass and ceramic fuse types, always replace them in pairs.

I.D.	Description	System	Part number
F1	Ceramic Type F 20A/250 V	Oven heater	2110-0098
F2	Ceramic Type F 20A/250 V	Oven heater	2110-0098
F3	Glass Type F 8A/250 V	Main supply	2110-0036
F4	Glass Type F 8A/250 V	Main supply	2110-0036

## AC power board schematic

The AC power board supplies AC power to the oven heater and to the main transformer.

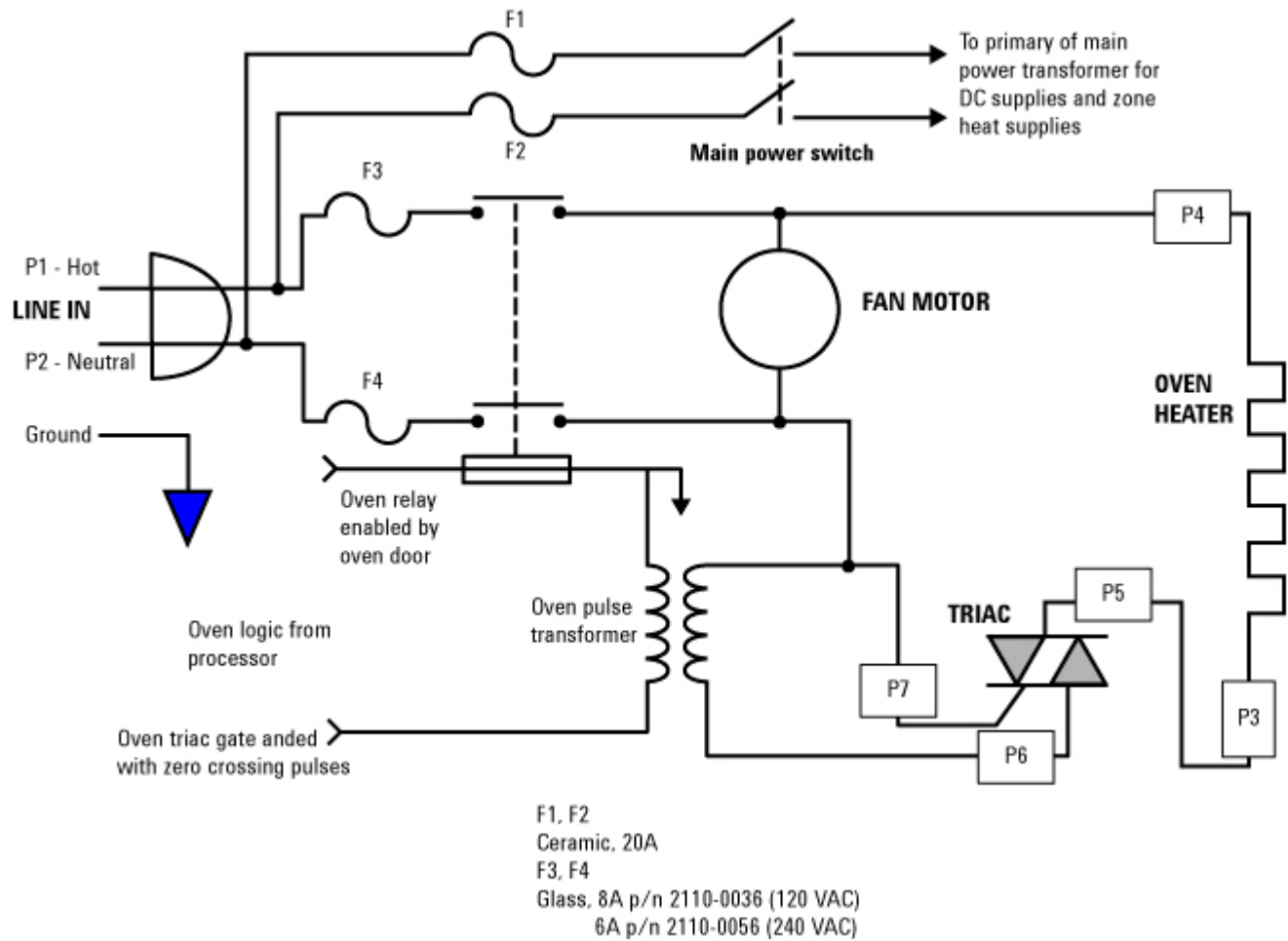


## Oven shroud

There are two different oven shrouds depending on the power option used. The oven shroud contains the oven heater and sensor as part of the assembly. If you need to replace the heater or sensor, you should replace the entire shroud assembly. See [“Replacing the oven shroud, oven fan, and oven fan motor”](#) for more details.

**Table 33** Part numbers for Oven Shrouds and Configuration Plugs

Voltage	Shroud part no.	Configuration plug part no.
<b>Regular oven shrouds</b>		
120 V	G1530-61610	G1530-60690
230 V	G1530-61670	G1530-60720
<b>Fast ramping oven shrouds</b>		
200 V	G1530-61620	G1530-60700
208 V	G1530-61630	G1530-60710
220 V	G1530-61630	G1530-60710
230 V	G1530-61650	G1530-60720
240 V	G1530-61640	G1530-60730

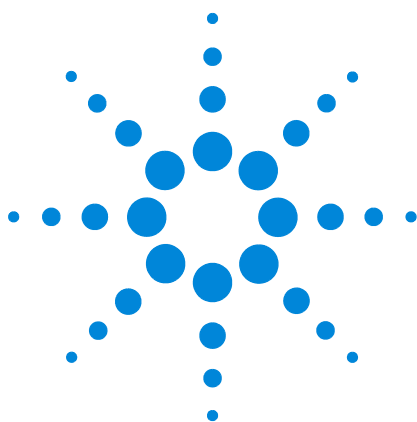

## 13 Firmware

Firmware overview	483
To Update GC Firmware	484
Firmware history	487
Confirm firmware update	495
Problems and solutions	496
Clearing battery-backed PF-RAM (7890A only)	497
Set the default method	500
Default method listing	501



## Firmware overview

The firmware is the programming that controls the entire GC. Periodic updates to correct errors and/or add features are available from the Agilent Web site and may be downloaded to the GC using the Agilent G4600AA Instrument Utilities software that shipped with the 7890A GC.

To successfully update GC firmware, you can use any version of the Instrument Utilities software. **Do not use any other firmware update utility.** To access the firmware update utility in Instrument Utilities, open the software and look for the Firmware Update link. You do not need to be monitoring the instrument.

About 160 7890A GCs were shipped to channel partners, distributors, and regional sales offices with very early versions of firmware. Contact the CAS factory problem manager (FPM) for instructions on how to update these units.

The latest version of Instrument Utilities software includes the ability to update the board serial number in the event that this was not done at the GC after replacing the logic board. (It is still best to update the serial number at the GC front panel after changing the logic board.) This version can also successfully update GCs without a manufacture date.

## To Update GC Firmware

To update the GC firmware, use Agilent Instrument Utilities, version B.01.06 or greater. When possible always use the latest version of Instrument Utilities.

The software maintains configuration information (including serial number) when updating a GC. However, any stored methods are lost. Click **Restore** to restore the active method after updating firmware.

**NOTE**

The Firmware Update utility requires full control of the instrument. Be sure to disconnect any Agilent data systems before using the Firmware Update utility.

- 1 Open Instrument Utilities.

**NOTE**

Once you connect to an instrument in the **Firmware Update** screen, no other tasks can be performed on the instrument until you **Disconnect**.

- 2 Make sure the instrument/component is turned on and not in use. Make sure any sampler or other device is also connected and turned on.
- 3 Record any information about the device stored in the device itself, for example, any local methods or service information (Early Maintenance Feedback counters, runtime details, and so forth).
- 4 Select **Firmware Update** from the explorer pane.
- 5 Select the desired instrument from the **Instrument** list and click **Connect**.
  - If it is not available, click **Add New Instrument** to open **Configure Instruments**. Add the new instrument, then return to firmware update.
- 6 Highlight the row of the component you wish to update.
- 7 Select the desired firmware version from the **Available** drop-down list. If the desired firmware version is not available, add it as follows:
  - a Go to the Agilent web site at <http://www.chem.agilent.com/en-US/Support/Downloads/firmware/Pages/default.aspx>.
  - b Download the updated firmware for your instrument type.



- c Unzip the firmware file(s).
  - d In the Firmware Update utility, click **Add**. Browse to and select the files unzipped in the previous step. The new firmware version will now appear in the **Available** drop-down list.
- 8 Click **Update**.
- 9 When prompted, save the active method and configuration data. You will restore the data from this file in a later step.

**WARNING**

**Do not turn off power to the device during firmware update.**

---

- 10 When the process completes, the instrument reboots. Reconnect to the instrument, then click **Restore**. Select the previously-saved data file to restore the active method.
- 11 Click **Disconnect** to disconnect the firmware update utility from the instrument.

## **Troubleshooting problems**

If the Instrument Utilities fails to connect to the GC, reboot the PC and try again.

If the firmware does not completely download to the GC, reboot the GC and try again. (Note that firmware download takes time, so be patient.)

If the firmware update fails, try again. Reboot the GC and PC if needed.

## Firmware history

### NOTE

The firmware history here is up to date as of the release of this manual. Check service notes and the firmware download site (Agilent internal: [http://cas.service.agilent.com/gc/Firmware/firmware\\_file\\_downloads.htm](http://cas.service.agilent.com/gc/Firmware/firmware_file_downloads.htm)) for more recent firmware versions and information.

Serial number  $\leq$  number in table - the unit could be updated to the latest revision.

Serial number  $>$  number in the table - the unit should already have this revision or greater.

Revision	CN serial number break	US serial number break	Service note
B.01.01			
A.01.10.2	TBD	TBD	G3440A-45
A.01.10.1	TBD	TBD	G3440A-41, -39
A.01.10	TBD	TBD	G3440A-37
A.01.09.2	TBD	TBD	G3440-31
A.01.09.1	CN10839078	TBD	G3440-30
A.01.09	CN10831XXX	US10831015	G3440A-27
A.01.08	TBD	TBD	G3440A-20
A.01.07	CN10751086	US10751031	G3440A-18
A.01.06	CN10738001	US10737020	
A.01.05	CN10731026	US10731014	
A.01.04	CN10722043	US10722003	G3440A-05
A.01.03	CN10717036	US10717009	G3440A-04
A.01.02.581	CN10716014	US10715013	G3440A-03
A.01.02	CN10708001	US10710002	G3440A-

Externally, firmware downloads are available from the Agilent web site.

Internally, go to:

[http://cas.service.agilent.com/gc/Firmware/firmware\\_file\\_downloads.htm](http://cas.service.agilent.com/gc/Firmware/firmware_file_downloads.htm) (link active as of October 2009)

**7890 A.01.10.2**

Described in Service note G3440A-45

Fixes: 7890 GC keypad now allows entry of serial numbers for injectors and tray. In certain circumstances the oven will not shut down if it cannot reach temperature for reasons such as an oven flap failure - This has now been fixed. Lab advisor can now connect to the GC whilst it is in MMON mode resolving firmware download connection problems. 7890 GC can now store barcode labels that contain more than 34 characters  
Injector light" is inconsistent with 7890 "Tower LED".

The error pertaining to the oven shut down may cause very hot air to be continuously released through the oven vent. It is advisable therefore that this firmware be installed on the GC at next visit. After installation of this firmware it is not necessary to requalify either GC nor the ALS. Firmware files can be obtained externally at [http://www.chem.agilent.com/\\_layouts/agilent/downloadFirmware.aspx?whid=50307](http://www.chem.agilent.com/_layouts/agilent/downloadFirmware.aspx?whid=50307) or internally at [http://cas.service.agilent.com/gc/Firmware/firmware\\_file\\_downloads.htm](http://cas.service.agilent.com/gc/Firmware/firmware_file_downloads.htm)

## 7890 A.01.10.1

Described in Service note G3440A-41, G3440-39

The 7890A GC firmware, A.01.10.1 has recently been released. Included in this firmware are not only enhancement and bug fixes but diagnostic "watchdogs" to help resolve and track possible future problems. Recent instrument failures have occurred that could have been avoided or troubleshooted more efficiently with the latest firmware installed.

GC Front/Back detector hardware fault (Team Track 02500) – Spurious error message. Error introduced in A.01.10. Evident at power on. NPD Analog signal always 0.0V (Team Track 02501). Error introduced in A.01.10. Digital Signal works OK. At times, GC keyboard, UI & LAN response is sluggish/stalled (Team Track 02502) PTV cryogenic cooling operating as per A.01.07. Cryo did not control inlet during run. (Team Track 00233) Fixes for External ready for 7820A. Fix method file upload and download of inlet temperature ramps. Fix power up procedure to not energize the small zone contactor when there is a shorted driver. Fix detector thermal shutoff status messages. Fix download of column calibration. Clean up small zone readiness calculation. Fix small zone shorted sensor detection. Add delay to valve state upload during a run. PVCS2476 Fix PTV Cryo for SOLVENT VENT Time. Fix method download during a run for non programmed zones. New PCOC PID's and control. MMI Temperature minimum for MMI multimode inlet set to -180 degC (Team Track 00234) LVI with a 500ul syringe set incorrect draw speed. (Team Track 00056) Fix Solvent Vent pressure calculation, use gage pressure. Give MMI inlets 10 ramps. 7693 ALS 7890 System Calibration should fail if G4514A internal calibration not completed (Team Track 00049) Tray LED does not flash when G4514A does not find a bottle during return all. (Team Track 00029). GC display updates BCR firmware revision correctly after a firmware update. Tray chiller info and text modified – chiller temperature upper limit is now 60C. Give MMI inlets 10 ramps Fix error in FPD single heater configuration. Tray Chiller configuration items added under Tray Config Tray Chiller temperature added to GC readiness Sampler devices reset after a firmware update/SN change causes 7890 to reconfigure QQQ Fix collision cell fault messages. Add QQQ to system readiness, Add QQQ shutdown test. Add QQQ diagnostic panels. Add Argon control mode to QQQ based on ECD Nitrogen and Argon frit calibrations. Add 'QQQ' – Collision cell gas controller Fix cannot pick collision cell from aux epc list

**7890 A.01.10**

Described in Service note G3440A-37

Enhancements: Support for the 7693A ALS and associated modules. Support for the 7000A Triple Quad GCMS. Support for the MMI. Fixes: Fix PTV Cryo for SOLVENT VENT Time. Fix method download during a run for non programmed zones. New PCOC PID's and control. Fix missing small zone fault messages and events. Fix download of column calibration. Clean up small zone readiness calculations. Fix small zone shorted sensor detection. Add delay to valve state upload during a run. Fix power up procedure to not energize the small zone contactor when there is a shorted driver. Fix detector thermal shutoff status messages. Fixed readiness for external devices Fix method file upload and download of inlet temperature ramps.

**7890 A.01.09.2**

Described in service note G3440A-31

Customers have been experiencing problems with oscillating temperatures with the PCOC inlets of the 7890. When the inlet was set to 50-80 degC, the temperature was seen to fluctuated +/- 5 degC causing a not ready status. In some cases the temperature would stabilize but only after a considerable length of time. Replacing the heater or even the whole weldment may or may not fix the problem.

Solution : The heating control coefficients (PIDs) have been modified in the firmware to better control and resolve these fluctuating temperatures. The new firmware revision is A.01.09.2 If you have have a customer with a PCOC inlet, regardless of whether they are experiencing temperature problems or not, it is advisable to update the firmware to A.01.09.2

**Additional Information**

In addition to the PCOC resolution one further problem has been fixed with this version. It has been noted that the DC210 led on the logic board intermittently flashes when the normal state is solid on. The status of this led does not affect the normal operation of the GC. To avoid all confusion firmware has now forced this led to lit all of the time.

**7890 A.01.09.1**

Described in service note G3440A-30

A list of defects fixed are as follows: Fix cryo trap for JAS Fix oven cryo shutdown. Fix download during run for oven thermal, and pressure zones. Fix download during a run for column flow programs. PCM AuxC reports as PCM AuxB Not Ready (p\_aux.c) Fix detector makeup mode for download. Fix status message for bad serial #. Fix entering zero value for analog out. Fix analog out for detector OFF.

**7890 A.01.09**

Described in service note G3440A-027

Changes applicable to Customers and Support Providers are as follows.

- Able to specify Cryo Type for User Aux Heaters
- Fixed analog out noise problem.
- Turn on PMT voltage when FPD flame is turned ON.
- Turn off FID air and H2 on ignition fail and gas shutdown
- Enable test signals for valve box heaters
- Add current limit to NPD SetBead voltage. If current exceeds 100pA during bead turn on, detector will turn off.
- Add ColComp setpoints to method.
- Add zone\_limit to AuxHeater config reporting
- Add several aux temp models for ICPMS
- Fix reporting of AuxDet Heater On/Off
- Remove Cryo Focus Valve configuration.
- Add Configure Collision Gas Control Module.
- Fix status display to not allow flow on COC.
- Add Column Compensation output in analog and digital modes (via config detector).
- Network settings (IP address, DHCP etc) can now be modified from the keyboard
- while in MMON mode.
- ALS appeared to go off line after 3 weeks of continuous use. Power reset resolved issue. Now fixed.
- (ALS firmware A.01.01)
- Gaps in Bar Code Reader can now be read. (ALS firmware A.01.01)

**7890 A.01.08**

Described in service note G3440A-020

**7890 A.01.07**

Described in service note G3440A-018

**7890A.01.04**

Described in service note.G3440A-005.

The following defects or perceived defects were fixed:

- Fixed PTV inlet configuration defect where Install Inlet (Cryo) was not a selection in the configuration list after a power cycle.
- Fixed the PCM Autoflow zero feature. Symptom would include PCM issues when in flow control mode. Measured flow is greater than flow displayed on screen. To confirm problem, navigate to Options, Calibration, PCM and scroll to Flow zero. Press Off / No to reset to factory defaults after updating firmware.
- Fixed symptom where pneumatic shutdown of PCM cannot be cleared without fixing the leak and power cycling the GC. Symptom may appear as PCM does not accept On / Yes or entering a new flow or pressure.
- Fixed symptoms where text alignment is incorrect and the cursor disappears when Options | Keyboard & Display | Language is set to Chinese.
- Fixed symptom where GC method is not downloaded from control software.
- Fixed the symptom where GC triggers a thermal shutdown with the flaps wide open. Chromatographic results indicate that the flaps were partially open for several runs.
- Fixed the symptom where the PP inlet behaves as if there is an oven track temperature mode.
- Fixed the symptom where no exception or error message is displayed when back detector heater, logic board, or EPC module is "broken."
- Fixed the symptom where you cannot program User Key 2 until after programming User Key 1. If you go through the steps of programming User Key 2 before User Key 1, the key will not work.



- Default configuration for the Auxiliary 2 Heater changed to MSD transfer line.

The following modifications were added to enable channel partners to configure and control their applications:

- Increased the maximum temperature of a small heated zone from 400 to 450°C. Using this maximum on our standard inlet or detector may cause damage to the instrument.
- Added a lock to the configuration settings so that someone cannot modify the configuration without turning the lock off. The lock is located under Options | Keyboard & display. At power On the lock is set to On.
- Added the capability to specify the wattage of the heater for a small heated zone: inlet, detector or auxiliary heater.
- Added the capability to configure 3 custom heaters with specific PID coefficients and store those coefficients.
- Added the capability for our channel partners to configure a variety of heaters.
- The Valve Box heater can be configured as an independent isothermal heater.
- Added the capability to configure an Auxiliary heated zone to have a maximum temperature between 40 and 450 °C.
- Adds the capability to configure a 140 watt inlet or auxiliary heater and program it to heat between 1 to 5 °C/min.

Added the capability to configure the PTV without a heater.

**7890A.01.03**

Described in service note G3440A-004.

The modifications in this firmware revision include:

- Adjusts control of oven flapper. Eliminates flapper not closing completely.
- Improves communication using DHCP with lower baud rate hubs. Lower baud rate hubs include the generation of early twisted-pair hubs. An example hub is the Hewlett-Packard J3128A AdvanceStack Hub-8e 10 base-T hub with the speed of 10 Mb/sec.
- Fixes this symptom: Aux EPC channel 6 always reads 0 pressure in the software Edit Parameters screen.

**7890A.01.02.581**

Described in service note G3440A-003.

The modifications in this revision include:

- Prevents interruption of firmware update process when using the standard tool.
- Fixes sequencing with multi-position and gas sampling valves.
- Fixes the displayed flow when column inlet is unspecified for SSL and PTV inlets. The flow was displayed as zero, even though the flow was accurate.
- Prevents inappropriate thermal shutdown for inlet, detector, and AUX heated zones when they are heating slowly.
- Fixes pressure setpoint ranges for Pneumatic Control Module (PCM).
- Eliminates cycling or multiple GC reboots after a power cycle.

## Confirm firmware update

Confirm the firmware update by pressing [**Status**][**Clear**]. The display should look like this.

Agilent 7890A GC  
A.01.09

Power on successful

The firmware update will become part of the GC's Event Log.

When upgrading 7890A GC firmware from:

- A.01.0x to A.01.07, or from
- A.01.03, A.01.04, and A.01.05 to A.01.08/A.01.09

configuration information can be lost. See document G3430-90035. The information is especially important if also upgrading OLCDS or GC ChemStation or MSD Productivity ChemStation.

Please check the following setpoints on the GC:

- Inlet and detector gas types - Config | <Front or Back><Inlet or Det>.
- Heater configuration. For example, the MSD transfer line may need reconfiguration.
- LAN configuration information - Options | Communications.
- Serial number and manufacturing date - Service Mode | Diagnostics | Instrument Status.

### CAUTION

If the GC does not have a valid serial number, you may have issues connecting with the ChemStation software. If the GC does not have a valid manufacturing date, the configuration settings will not be saved and will be lost if you power cycle the GC.

## Problems and solutions

If the update is interrupted or there is a defect in either of the boards, the update may fail. Click on the link to see more information about each of these symptoms.

- GC firmware update unsuccessful and GC reboots into MMON mode
- ALS firmware update incomplete

### Possible solutions to unusual behaviors

At installation and after a firmware update, we have had some communication issues between the GC and the software. The symptoms include unexpected changes to configuration or to method setpoints. Here are some suggested actions for eliminating these symptoms.

- Load the default method. Press Method.
- Replace the logic board.
  - a From the keyboard and display, check the configuration and method in the hardware.
    - 1 From the software, upload the GC method.
    - 2 From the software, select Instrument | GC Configuration. Review the GC configuration after clicking on Get GC configuration.
    - 3 From the software, select Instrument | Edit GC parameters. Review your method setpoints.
- Save the method under a new name.

## Clearing battery-backed PF-RAM (7890A only)

After updating to a firmware revision before A.01.06, clearing all or part of the PF-RAM memory has eliminated some symptoms, such as, heater faults and corrupted OLCDS or ChemStation methods.

Clearing some or all of the memory forces the GC to reload defaults. See [“Default method listing”](#) for a listing of the values.

### Partial clear

This procedure resets both configuration and current setpoints. Here is a listing of what is reset to defaults.

- Configuration items: Oven, Display, Inlets, Detectors, Aux EPC modules, PCM EPC modules, Valves, Thermal Aux, and Valve box.
- Sequence items: All items stored under the sequence including: Sample numbers, Overlap, Auto-injector, and Method number.
- Current setpoints (Workfiles): GC and ALS setpoints.
- Instrument state flags: Run states, Blank run, Power fail, Method directory dirty, Status table, Analog column compensation values.

#### Procedure

- 1 Power OFF the GC.
- 2 Press and hold down **[Clear]**.
- 3 Power ON the GC.
- 4 When the GC Status screen appears on the display, release **[Clear]**.


### Total clear

This procedure clears all of the battery-backed PF-RAM. In addition to what is reset by holding down the Clear key, these items are reset to defaults:


- ID items: Serial number, IP networking setpoints, Uptime clock, Column names
- Other: Service counters, Event log

### Procedure

- 1 Power off the GC.
- 2 Remove the left-side panel to access the logic board. Find the small DIP switch (G3430-60100) or DIP switches (G3430-60120 and G3430-60101) between the battery and the large connector.
- 3 For the G3430-60120 and G3430-60101 boards, move the switch to OFF. For the G63430-60100 (shown below), move both switches away from you toward the GC oven. This is the Off position.

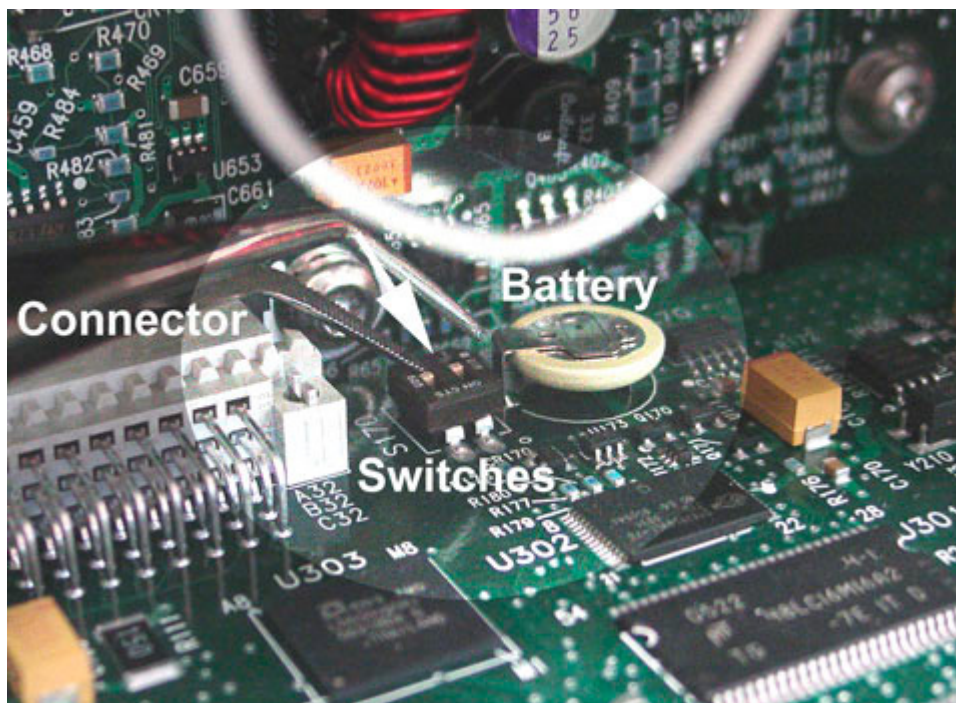

- 4 Wait 5 minutes to allow the battery-backed RAM to clear.
- 5 Reset the battery switches to On (G3430-60100 shown).

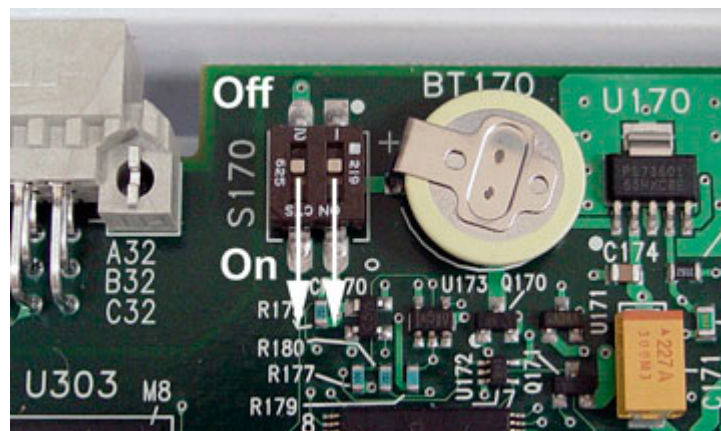

- 6 Power on the GC. The GC will post some messages about loading default configuration and method files.
- 7 Re-enter any lost configuration information, such as, gas type and method setpoints.

## Set the default method

To set the default method conditions, do the following:

- 1 Press **[Method]**.
- 2 Scroll to the last line of the table: **Set default method** and press **[Enter]**.
- 3 Press **[On/Yes]** to confirm the request.



## Default method listing

This page lists the default configuration and method setpoints for the GC.

### Instrument

- Post run hold time = 0
- Auto Prep Run = Unknown
- Keyboard lock = False
- Type of connection = LAN
- GCName = blank
- GCNotes = blank

**Table 34** Inlets

	COC	PP	PTV	SSL	VI
Determine readiness	True	True	True	True	True
COC oven track		False	False	False	
Septum purge flow		3	3	3	
Septum purge state		On	On	On	
Septum purge switching		False	False	False	
Purge switching time		0	0	0	
Total flow		3	18	3	
Total flow state		Off	Off	Off	
Control mode		Pressure	Pressure	Pressure	
Temperature state		Off	Off	Off	
Initial hold time		0	0	0	
Initial value		250	250	250	
Post run value		250	0	250	
Pressure state		Off	Off	Off	
Initial hold time		0	0		
Initial value		689475.7	689475.7	689475.7	
Post run value		689475.7	689475.7	689475.7	
Cryo		Off	Off	Off	
Cryo temperature		25	25	25	
Time out detection		Off	Off	Off	
Time out time		30	30	30	
Fault detection		Off	Off	Off	
Quick cool		Off	Off	Off	

	COC	PP	PTV	SSL	VI
Gas saver time		2	2	2	
Gas saver flow		20	20	20	
Gas saver state		Off	Off	Off	
Operation mode		Splitless	Solvent vent	Split	
Split ratio		0	100	100	
Purge flow		15	15	15	
Purge time		0.75	0.35	0.75	
Pressure pulse		0	689475	689475	
Pulse time		0	0.75	0.75	
Vi Inject time		0	0	0	
Vent flow		0	100	0	
Vent pressure		0	344737	0	
Vent time		0	0.25	0	

#### Columns 1 to 6

- Ignore ready = True; except for Columns 3 through 6 = False
- Mode = Constant pressure
- Flow state = Unknown
- Initial hold time = 0
- Initial value = 1; except for Column 1 where Initial value = 0
- Post run value = 1

#### Oven

- Ignore ready = True
- Cyro = Off
- Use Temperature = 0
- Time out detection = Off
- Time out = 0
- Fault detection = Off
- Quick cool = Off
- Equilibration time = 3
- Temperature state = Off
- Initial hold time = 0
- Initial value = 50

- Post run value = 50

#### Auxiliary Thermal Zones

- Thermal Aux 1, Determine readiness = True
- Thermal Aux 2, Determine readiness = True
- Thermal Aux 3, Determine readiness = False

**Table 35** Detectors

	FID	FPD	NPD	TCD	μECD
Determine readiness	True	True	True	True	
Temperature	250	200	250	250	
Heater	Off	Off	Off	Off	Off
Fuel flow	40	75	2	0	
Fuel	Off	Off	Off	Off	Off
Utility flow	450	100	60	10	
Utility	Off	Off	Off	Off	Off
Makeup flow	50	60	30	5	
Makeup	Off	Off	Off	Off	Off
Constant total flow	False	False	False	False	False
Total flow	0	0	0	0	0
Electronics	On	On	On	Off	
Flame	Off	Off	Off	Off	Off
State	Unknown	Unknown	Unknown	Unknown	Unknown

#### Auxiliary flow and pressure control modules

For each of the 9 possible Aux EPC channels (AUX).

- Determine readiness = False

For each of the PCM channels that function like a purged inlet.

- Determine readiness = True
- Mode = Forward pressure control
- Total flow = 0
- Pressure = 689475.7

- Pressure State = Off
- Initial hold time = 0
- Initial value = 689475.7
- Post run value = 689475.7

For each of the PCM channels that function like an Aux EPC channel.

- Determine readiness = False
- Pressure state = Off
- Initial hold time = 0
- Initial value = 689475.7
- Post run value = 689475.7

### Valves

For GC x GC Valve.

- GCxGC Valve sample time = 250
- Modulation time = 60000
- Valve idle state = Off
- Delay time = 60000

For each of the 8 possible valves.

- Valve initial state = Off

For Cyro-focusing valve.

- Pre-run time = 1
- Hold time = 1

**Table 36** Automatic liquid sampler

	Front injector	Back injector
Fan	On	On
<b>Pre-wash A Reps</b>	0	0
Volume	80	80
Draw speed	200	200
Dispense speed	10	10
Viscosity delay	0	0
<b>Pre-wash B Reps</b>	0	0
Volume	80	80
Draw speed	200	200
Dispense speed	10	10
Viscosity delay	0	0
<b>Sample wash Reps</b>	0	0
Volume	80	80
Draw speed	200	200
Dispense speed	10	10
Viscosity delay	0	0
Depth offset	0	0
<b>Sample pumps Reps</b>	6	6
Draw speed	200	200
Dispense speed	10	10
Viscosity delay	0	0
Depth offset	0	0
<b>Sample volume</b>	10	10
Draw speed	200	200
Viscosity delay	0	0
Depth offset	0	0
Injection dispense speed	10	10
Pre-injection dwell time	0	0
Post-injection dwell time	0	0
Large volume reps	1	1
Large volume delay	0	0
<b>Post-wash A Reps</b>	0	0
Volume	80	80
Draw speed	200	200
Dispense speed	10	10
Viscosity delay	0	0
<b>Post-wash B Reps</b>	0	0
Volume	80	80
Draw speed	200	200
Dispense speed	10	10
Viscosity delay	0	0

### Tray

- Temperature = 0
- Barcode reader = false
- Mixing intensity = Unknown
- Time = 0
- First sample preference = 2
- Enable ALS tray = False
- Overlapped injection mode = No overlap
- Overlapped injection delay = 0


### Service counters

- Twenty (20) counter IDs set to a default number.
- Enabled = False.

### Signals

Each of the 4 digital signals are set to

- Test Plot at 50 Hertz
- Signal is saved = False.

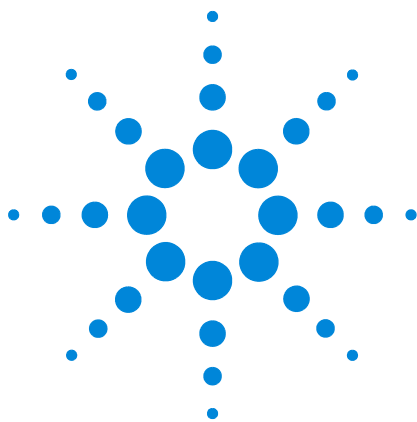

## 14 Firmware Instructions for Channel Partner Devices

Firmware overview for channel partners	508
Tuning a small heated zone	509
Thermal Aux Configuration	510
Valve Box Configuration	521
Inlet configuration	525
Cryo focus valve configuration	540
Cryo trap configuration	541
Detector configuration	543



## **Firmware overview for channel partners**

This section contains instructions for installing and configuring channel partner inlets, detectors, thermal aux zones, and cryo focus valves, on an Agilent 7890 Series GC. Each section includes information on:

- How to wire your hardware
- How to configure the GC to recognize your hardware
- How to access the method parameter screen for your hardware

For configuration of standard Agilent devices for the GC, see the Accessory Installation sheets and the Agilent 7890A GC Advanced User Guide or 7890B Operation Manual.



## Tuning a small heated zone

- 1 Disable the I and D terms. The D term is set to zero to disable it, but the I term has an inverse effect, i.e. the larger the I term the smaller its effect, so set the I term to its maximum value.
- 2 Increase the P term as needed to achieve oscillation, then decrease it until the oscillation is barely maintained. Note the value of the P term, which we call  $K_u$ , and the period of the oscillation in seconds, which we call  $P_u$ .
- 3 Use these values in one of the standard PID tuning methods, such as Ziegler-Nichols, to calculate the PID terms.

For example: using Ziegler-Nichols tuning the PID values would be:

$$P = K_u / 1.7$$

$$I = P_u / 2$$

$$D = P_u / 8$$

Ziegler-Nichols is probably a good tuning to use for any zone that needs to be temperature programmed. It has a small amount of overshoot but it settles fairly quickly. For zones that don't need to be temperature programmed there are other tunings that have less overshoot but are slower to settle, such as Tyreus-Luyben tuning. The PID values using Tyreus-Luyben tuning are:

$$P = K_u / 2.2$$

$$I = P_u * 2.2$$

$$D = P_u / 6.3$$

Reference: <http://www.chem.mtu.edu/~tbco/cm416/zn.html>

## Thermal Aux Configuration

There are several control models to select for the standard Agilent thermal aux zones. These include the AED transfer line, the MSD transfer line, and the Nickel catalyst. This section only covers the custom configuration of a thermal zone that is available by selecting **User configurable heater** as the auxiliary heater type during configuration. See the *Agilent 7890A GC Advanced User Guide* or *7890B Operation Manual* for configuration of standard thermal aux zones.

Up to 3 thermal auxiliary zones can be configured. All 6 of the GC's small heated zones are available for use. If cooling and heating are required in the thermal zone's control, a cryo valve is selected during configuration of a thermal zone.

### Heater and cryo connection locations

Each heater connector has a 40 to 180 watt heater power supply and a temperature sensor for the control loop. Labels on the wiring harness near the connector identify the connection. The labels match the GC display except for the FI or BI wiring harness labels displayed as FINLET and BINLET.

Connector label / GC display identifier	Connector location on GC
A1	Right side electrical compartment position A1 (far left) on valve bracket electrical harness
A2	Right side electrical compartment position A2 (far right) on valve bracket electrical harness
FI/FINLET	Left top side of GC frame adjacent to front inlet
BI/BINLET	Left top side of GC frame adjacent to back inlet
F-DET	Right side electrical compartment at top position adjacent to front detector
B-DET	Right side electrical compartment at top position adjacent to back detector
FV	Left top side of GC frame adjacent to front inlet
BV	Left top side of GC frame adjacent to back inlet

### Connections assigned to a thermal zone

Each Aux thermal zone can support one of two fixed heating connections containing a heater and temperature sensor. When one of the two allowed heater connections is configured to an Aux Zone the other connection is not available for thermal zone assignment. Additionally each thermal zone can support a cryo valve.

Thermal Zone	Thermal Aux 1	Thermal Aux 2	Thermal Aux 3
Heater Connection	Aux1 or BI/BINLET	Aux2 or B-DET	FI/FINLET or F-DET
Cryo Valve Connection	BV	BV	FV

The *User Configurable Heater* allows heated zones ranging from 40 watts to 180 watts.

## Auxiliary Type Selection

After assigning a heater or heater and cryo valve to an Aux thermal zone and rebooting the GC, you configure the type of device requiring thermal control. This is done by selecting the **Auxiliary type:** field from the configuration menu and pressing the [Mode/type] key. Several Agilent standard devices are listed along with a selection for a **User configurable heater**.

When configuring a device for thermal control that is not listed you can select the **User configurable heater** and create a custom thermal control model. All but the **User configurable heater** control model have a fixed maximum device setpoint of 400 deg C. The only auxiliary type to support a cryo valve is **User configurable heater**.

Auxiliary type	Use
AED transfer line	For the AED transfer line
MSD transfer line	For the Agilent GC/MSD transfer line
Nickel Catalyst	For the Agilent nickel catalyst
Unknown	
User configurable heater	For creation of a custom thermal device

## User Configurable Heater

The control algorithm used to maintain a stable temperature setpoint requires several constants in order to model the thermal response of a physical device. These constants can be modified to model a custom control device when the **Auxiliary type:** is set to **User configurable heater**.

### Custom heater parameters

**The PID and Mass fields** These fields allow entry of custom constants for the Proportional, Integral, Derivative, and Mass associated with a thermal device's control algorithm.

**The Power field** This field allows for entry of the maximum power that the heater requires. This value must be within 10% of the actual maximum or an error is generated.

**The Maximum Programming Rate field** Enter the maximum rate that the temperature of the device can change from one temperature to another in a fixed period of time. If this field is left at the default value of zero, the device cannot be temperature controlled.

**The Maximum Setpoint field** Enter the maximum setpoint allowed for the custom thermal device. The GC will not allow the user to enter a value higher than this number.

**Sensor** Pressing [Mode/Type] on this field brings up the sensor type for selection. Select from RTD and TC.

**The zone control model** Pressing [Mode/Type] on this field brings up the available zone control models summarized in the table.

Control Model	Use
SS	Isothermal zone with cryo
PTV	High power requirement with cryo and the ability to program temperature ramps during a run
COC	Low power requirement with cryo and the ability to program temperature ramps during a run

### Programmable temperature zones

The PTV and COC zone control models for the **User configurable heater** Auxiliary type and the **Unknown** Auxiliary type, allow a temperature program to be entered for the device similar to that available on an Agilent standard PTV and COC inlet.

## Cryo control model

This field is only displayed if a cryo valve is configured in the hardware.

Pressing [Mode/Type] on this field brings up the available cryo control models summarized in the table.

Cryo control model	Use
PTV	A control model that is similar to the Agilent PTV inlet.
CryoTrap	A control model that lowers a device's temperature anywhere during a run and then increase the device's temperature at a later time during the run.

### Configuring the thermal Aux zone (Firmware A.01.07)

- 1 With the GC disconnected from its power source, install the device containing the heater and temperature sensor and optionally install a cryo valve.
- 2 Plug the heater and sensor cable from this device into one of the six small heated zone connections.
- 3 Optionally plug the cryo valve 24-volt power supply into the FV or BV connector.
- 4 With all covers installed, attach the power supply to the GC and power the GC on.
- 5 Unlock the GC configuration. Press the [Options] key, select **Keyboard & Display** and press the [Enter] key. Scroll down to **Hard Configuration Lock** and press the [off] button.
- 6 From the section above titled "Connections assigned to a thermal zone" determine the Aux Zone (1, 2 or 3) that you are using based on the heater connection.
- 7 On the GC keypad press [Config] [Aux Temp #] press the number of the thermal zone 1, 2, or 3.
- 8 With **Unconfigured:** selected press the [Mode/Type] key on the GC keypad and if necessary scroll to see all selections. You must only make a selection that contains the GC identifier for the heater zone that your device is attached to. If a cryo valve is attached to the FV or BV connector, the GC also displays a selection for your heater zone and a cryo valve configuration. Press [Enter] to complete your selection.
- 9 When the GC prompts to power cycle the GC, reboot the GC by pressing the [Options] key and select **Communications**. From this option scroll down and select the reboot option and press [Enter] on the GC keypad. This soft reboot keeps the GC configuration unlocked for subsequent steps. Power cycling the GC forces you to unlock the GC configuration on power-up.
- 10 On the GC keypad, press [Config] [Aux Temp #] press the number of the thermal zone 1, 2, or 3. The display now shows your configured devices.
- 11 Scroll to **Auxiliary type:** and press the [Mode/Type] key on this field to bring up the available Auxiliary types. Select **User Configurable Heater** to create a custom thermal device and press [Enter].
- 12 If a cryo valve is included in this Aux Temperature zone, scroll to **Cryo type:** press the [Mode/Type] key and select the coolant supplied to the cryo valve. Next scroll to **Cryo**



**control model:** and press the [Mode/Type] key and select a suitable model from the list.

When the GC prompts to power cycle the GC, reboot the GC by pressing the [Options] key and select **Communications**. From this option scroll down and select the reboot option and press [Enter] on the GC keypad. This soft reboot keeps the GC configuration unlocked for subsequent steps.

On the GC keypad, press [Config] [Aux Temp #] press the number of the thermal zone 1, 2, or 3. The display now shows your configured devices.

- 13 Next scroll to and enter the initial values for **Proportional gain, Integral time, Derivative time** and **Mass**. These values are usually later modified based on empirical results until stable thermal performance is achieved.
- 14 Scroll to **Power (watts)** and enter the maximum power used to heat the device.
- 15 Scroll to **Maximum programming rate** and enter a rate of temperature rise per minute that your device is capable of achieving.
- 16 Scroll to **Maximum setpoint** and enter a value that is within the safe operating temperature of the device.
- 17 Scroll to the **Sensor:** field and make sure the sensor type agrees with your hardware. To change the sensor type, press the [Mode/Type] key and select the RTD or TC (thermocouple with conversion board) type.

When the GC prompts to power cycle the GC, reboot the GC by pressing the [Options] key and select **Communications**. From this option scroll down and select the reboot option and press [Enter] on the GC keypad. This soft reboot keeps the GC configuration unlocked for subsequent steps.

On the GC keypad, press [Config] [Aux Temp #] press the number of the thermal zone 1, 2, or 3. The display now shows your configured devices.

- 18 Scroll to **Zone control model:** and press the [Mode/Type] key on this field to bring up the available zone control models. Select a suitable control model and press [Enter].

When the GC prompts to power cycle the GC, reboot the GC by pressing the [Options] key and select **Communications**. From this option scroll down and select the reboot option and press [Enter] on the GC keypad. This soft reboot keeps the GC configuration unlocked for subsequent steps.

On the GC keypad, press [Config] [Aux Temp #] press the number of the thermal zone 1, 2, or 3. The display now shows your configured devices.

- 19** Power cycle the GC to lock this configuration.
- 20** On the GC keypad press [Aux Temp #]. If other Thermal Aux zones exist, press the keypad number (1, 2, or 3) for this zone. This display the Thermal Zone's parameter list for entering operating values for a method.

## Configuring the thermal Aux zone (Firmware A.01.08 and A.01.09)

- 1 With the GC disconnected from its power source, install the device containing the heater and temperature sensor and optionally install a cryo valve.
- 2 Plug the heater and sensor cable from this device into one of the six small heated zone connections.
- 3 Optionally plug the cryo valve 24 volt power supply into the FV or BV connector.
- 4 With all covers installed, attach the power supply to the GC and power the GC on.
- 5 Unlock the GC configuration. Press the [Options] key, select Keyboard & Display and press the [Enter] key. Scroll down to Hard Configuration Lock and press the [off] button.
- 6 From the section above titled "Connections assigned to a thermal zone" determine the Aux Zone (1, 2 or 3) that you are using based on the heater connection.
- 7 On the GC keypad press [Config] [Aux Temp #] press the number of the thermal zone 1, 2, or 3.
- 8 With Unconfigured: selected press the [Mode/Type] key on the GC keypad and if necessary scroll to see all selections. You must only make a selection that contains the GC identifier for the heater zone that your device is attached to. If a cryo valve is attached to the FV or BV connector, the GC also displays a selection for your heater zone and a cryo valve configuration. Press [Enter] to complete your selection.
- 9 When the GC prompts to power cycle the GC, reboot the GC by pressing the [Options] key and select Communications. From this option scroll down and select the reboot option and press [Enter] on the GC keypad. This soft reboot keeps the GC configuration unlocked for subsequent steps. Power cycling the GC forces you to unlock the GC configuration on power-up.
- 10 On the GC keypad, press [Config] [Aux Temp #] press the number of the thermal zone 1, 2, or 3. The display now shows your configured devices.
- 11 Scroll to Auxiliary type: and press the [Mode/Type] key on this field to bring up the available Auxiliary types. Select User Configurable Heater to create a custom thermal device and press [Enter].
- 12 If a cryo valve is included in this Aux Temperature zone, scroll to Cryo type: press the [Mode\Type] key and select the coolant supplied to the cryo valve. Next scroll to Cryo control

model: and press the [Mode/Type] key and select a suitable model from the list.

- 13 To display all the configuration parameters select View Custom Heater Conf and press [On/Yes].
- 14 Next scroll to and enter the initial values for Proportional gain, Integral time, Derivative time and Mass. These values are usually later modified based on empirical results until stable thermal performance is achieved.
- 15 Scroll to Zone control model: and press the [Mode/Type] key on this field to bring up the available zone control models. Select a suitable control model and press [Enter].
- 16 Scroll to the Sensor: field and make sure the sensor type agrees with your hardware. To change the sensor type, press the [Mode\Type] key and select the RTD or TC (thermocouple with conversion board) type.
- 17 Scroll to Power (watts) and enter the maximum power used to heat the device.
- 18 Scroll to Maximum setpoint and enter a value that is within the safe operating temperature of the device.
- 19 Scroll to Maximum programming rate and enter a rate of temperature rise per minute that your device is capable of achieving.
- 20 Power cycle the GC to lock this configuration.
- 21 On the GC keypad press [Aux Temp #]. If other Thermal Aux zones exist, press the keypad number (1, 2, or 3) for this zone. This display the Thermal Zone's parameter list for entering operating values for a method.

## Valve Box Configuration

The GC allows for a valve box hardware configuration consisting of one or two heaters in the valve box. An Agilent standard configuration uses one or two 70 watt heaters in the valve box. If you set the **Custom Heater Conf. Parameter** to [Off/ No] the PID, Mass, and Power fields use values optimized for these 70 watt heaters and any custom setting is ignored.

You may also design a custom heated valve box with a need for different heater parameters. To do this, set the **Custom Heater Conf. Parameter** to [On /Yes] and enter custom values for the PID, Mass, and Power fields.

## Heater connection locations

Each heater connector has a 40 to 180 watt heater power supply and a temperature sensor for the control loop. Labels on the wiring harness near the connector identify the connection.

Connector label / GC display identifier	Connector location on GC
A1	Right side electrical compartment position A1 (far left) on valve bracket electrical harness
A2	Right side electrical compartment position A2 (far right) on valve bracket electrical harness

## Custom heater parameters

### Custom Heater Conf.

Set the value to [On/Yes] to enable your custom heater parameters entered below. Setting the value to [Off/No] defaults to the Agilent standard valve box values for these numbers.

### View Custom Heater Conf.

Set the value to [On/yes] to display the custom heater field values of PID, Mass, and Power.

### The PID and Mass fields

These fields allow entry of custom constants for the Proportional, Integral, Derivative, and Mass associated with a thermal device's control algorithm.

### The Power field

This field allows for entry of the maximum power that the heater requires. This value must be within 10% of the actual maximum or an error is generated.

## Configuring the Valve Box Zone

- 1 With the GC disconnected from its power source, install the custom heater(s) and temperature sensor(s). You can install up to 2 heaters and temperatures for this zone.
- 2 Plug the heater and sensor cable(s) from the heater(s) into the A1 and/or A2 small heated zone connection.
- 3 With all covers installed, attach the power supply to the GC and power the GC on.
- 4 Unlock the GC configuration. Press the [Options] key, select **Keyboard & Display** and press the [Enter] key. Scroll down to **Hard Configuration Lock** and press the [off] button.
- 5 On the GC keypad press [Config] [Aux Temp #] and select **Valve Box**.
- 6 With **Unconfigured:** selected press the [Mode/Type] key on the GC keypad and select one of the following: Install Heater A1; Install Heater A2; Install 2 htr A1 & A2.. Press [Enter] to complete your selection.
- 7 When the GC prompts to power cycle the GC, reboot the GC by pressing the [Options] key and select **Communications**. From this option scroll down and select the reboot option and press [Enter] on the GC keypad. This soft reboot keeps the GC configuration unlocked for subsequent steps. Power cycling the GC forces you to unlock the GC configuration on power-up.
- 8 On the GC keypad, press [Config] [Aux Temp #] and select valve box. The display now shows your configured heater(s).
- 9 To allow use of the custom heater parameters, scroll to **Custom Heater Conf.** And press the [On/Yes] key.
- 10 To display all the configuration parameters select **View Custom Heater Conf** and press [On/Yes].
- 11 Next scroll to and enter the initial values for **Proportional gain**, **Integral time**, **Derivative time**, and **Mass**. These values are usually later modified based on empirical results until stable thermal performance is achieved.
- 12 Scroll to **Power (watts)** and enter the maximum power used to heat the device.
- 13 Power cycle the GC to lock this configuration.
- 14 On the GC keypad press [Aux Temp #] and select Valve Box. This displays the Valve box heating setpoint value for a method.



## Inlet configuration

A custom inlet can be configured using a channel partner PTV with or without a heater or cryo valve. Additionally a channel partner PCM configuration allows an Agilent PCM flow control to be used in a heated inlet control loop.

This section does not cover the configuration of the standard Agilent inlet EPC modules. See the *Agilent 7890A GC Advanced User Guide* or *7890B Operation Manual* for the configuration of Agilent standard inlets.

## Hardware configuration

Use this table to determine the allowable hardware configurations, where the hardware electrical connection must be made, and the configuration to select for this hardware as displayed by the GC. All hardware must be installed before attempting to configure the inlet.

GC configuration display selections	Hardware Configuration	EPC wiring connection	Heater wiring connection	Cryo valve connection
<b>Front inlet PTV</b>				
Install Front Inlet	PTV + HTR	EPC1	FI	None
Install Front Inlet (CRYO)	PTV+HTR+CRYO	EPC1	FI	FV
Install Front Inlet (NO HTR)	PTV	EPC1	None	None
<b>Front inlet PCM</b>				
Install Front Inlet	PCM + HTR	EPC1	FI	None
Install Front Inlet (CRYO)	PCM+HTR+CRYO	EPC1	FI	FV
<b>Back inlet PTV</b>				
Install Back Inlet	PTV + HTR	EPC2	BI	None
Install Back Inlet (CRYO)	PTV+HTR+CRYO	EPC2	BI	BV
Install Back Inlet (NO HTR)	PTV	EPC2	None	None
<b>Back inlet PCM</b>				
Install Back Inlet	PCM + HTR	EPC2	BI	None
Install Back Inlet (CRYO)	PCM+HTR+CRYO	EPC2	BI	BV

## EPC module installation

A PCM or PTV EPC installed in the standard EPC module location for a front or back inlet can be custom configured as a channel partner (CP) inlet by the GC. This requires removal of the top EPC module cover.

EPC wiring harness label	Location
EPC1	Left front EPC module slot
EPC2	Left back EPC module slot

## Heater and Cryo installation

Optionally install and attach your inlet device's heater and inlet controlled cryo valve to the GC's wiring harness. This requires the removal of the left side GC cover.

Wiring harness label	Location
FI	Left top side of GC frame adjacent to front inlet
FV	Left top side of GC frame adjacent to front inlet
BI	Left top side of GC frame adjacent to back inlet
BV	Left top side of GC frame adjacent to back inlet

## The Inlet configuration display

### Gas type

Pressing the [Mode/Type] key on this field brings up the choices for carrier gases connected to an inlet. Select the carrier gas going to this inlet and press the [Enter] key.

### Cryo Parameters

These parameters are displayed when a cryo valve is included in the hardware configuration.

### Cryo type

Pressing the [Mode/Type] key on this field brings up the choices for cryo gas connected to the inlet's cryo valve. Select the cryo gas connected to the inlet's cryo valve and press the [Enter] key. Changing this type requires rebooting the GC. This field displays when a cryo valve is configured to the inlet.

### Cryo

Use the [On/Yes] key to enable cryo cooling on this inlet. The [Off/No] key disables cryo cooling. This field displays when a cryo valve is configured to the inlet.

### Use Cryo temperature

Enter a maximum temperature for cryo cooling activation. This field displays when a cryo valve is configured to the inlet and the cooling type is not air.

### Cryo timeout

Enter a time in minutes for the cryo valve to close if no run has started before this time elapses. This field displays when a cryo valve is configured to the inlet and the cooling type is not air.

### Cryo fault

The GC can shut down temperature control of the inlet if setpoint is not reached in 16 minutes. Use the [On/Yes] key to enable temperature shutdown on this inlet. The [Off/No] key disables temperature shutdown. This field displays when a cryo valve is configured to the inlet.

### Custom heater parameters

The control algorithm used to maintain a stable temperature setpoint requires several constants in order to model the thermal response of a physical device. To access these constants, from the inlet's configuration display select **View Custom Heater Conf** and press [On/yes].

**The PID and Mass fields** These fields allow entry of custom constants for the Proportional, Integral, Derivative, and Mass associated with a thermal device's control algorithm.

**The Power field** This field allows for entry of the maximum power that the heater requires. This value must be within 10% of the actual maximum or an error is generated.

**The Maximum Setpoint field** Enter the maximum setpoint allowed for the custom thermal device. The GC will not allow the user to enter a value higher than this number.

**Sensor** Pressing the [Mode/Type] key on this field brings up the sensor type for selection. Select from RTD and TC.

**The Maximum Programming Rate field** Enter the maximum rate that the temperature of the device can change from one temperature to another in a fixed period of time. If this field is left at the default value of zero, the device cannot be temperature controlled. An inlet with a zone control model of PTV or COC requires a non-zero value to allow temperature programming.

**The zone control model** Pressing the [Mode/Type] key on this field brings up the available zone control models summarized in the table below.

Control Model	Use
SS	Isothermal zone with cryo
PTV	High power requirement with cryo and the ability to program temperature ramps during a run
COC	Low power requirement with cryo and the ability to program temperature ramps during a run

### Programmable temperature inlets

The PTV and COC zone control models allow a temperature program to be entered for the device similar to that available on an Agilent standard PTV and COC inlet. The SS zone control

model is for isothermal control and therefore can't be programmed.

### **Configured hardware**

The last few lines of the inlet's configuration screen lists the hardware controlled by the inlet. One line shows the EPC modules slot location and type. Another line indicates the heater used and its wattage, and the final line the cryo valve use.

### Configuring a custom inlet (Firmware A.01.05)

- 1 With the GC disconnected from its power source, install a PCM or PTV EPC module in the front or back inlet EPC module location. Plug the EPC module into the correct wiring harness connector. See EPC module location above.
- 2 Optionally install the inlet device containing the heater and temperature sensor.
- 3 Plug the heater and sensor cable from this device into one of the inlet's heated zone connections. See the "Heater and Cryo Installation" section above for the wiring harness connector to use.
- 4 Optionally install a cryo valve and plug it into the 24 volt power supply into the FV or BV connector.
- 5 With all covers installed, attach the power supply to the GC and power the GC on.
- 6 Unlock the GC configuration. Press the [Options] key, select **Keyboard & Display** and press the [Enter] key. Scroll down to **Hard Configuration Lock** and press the [off] button.
- 7 On the GC keypad press [Config] [Front Inlet] for the configuration of a front inlet device or [Config] [Back Inlet] for the configuration of a back inlet device.
- 8 With **Unconfigured:** selected press the [Mode/Type] key on the GC keypad and if necessary scroll to see all selections. See the table in the Hardware Configuration section for the correct display line to select. Press [Enter] to complete your selection.
- 9 When the GC prompts to power cycle the GC, reboot the GC by pressing the [Options] key and select **Communications**. From this option scroll down and select the reboot option and press [Enter] on the GC keypad. This soft reboot keeps the GC configuration unlocked for subsequent steps. Power cycling the GC forces you to unlock the GC configuration on power-up.
- 10 Select the configuration display for this inlet. The display now shows your configured hardware devices.
- 11 Scroll to **Gas type:** and press the [Mode/Type] key on this field to bring up the available carrier gas types. Select the gas type attached to the inlet and press [Enter].
- 12 If a cryo valve is included in this inlet, scroll to **Cryo type:** press the [Mode\Type] key and select the coolant supplied to the cryo valve. Set the **Cryo**, **Cryo timeout**, and **Cryo fault** parameters as required. See above.



- 13 To display all the configuration parameters select **View Custom Heater Conf** and press [On/Yes]. This displays all the heater parameters that can be set.
- 14 Scroll to **Zone control model:** and press the [Mode/Type] key on this field to bring up the available zone control models. Select a suitable control model and press [Enter].
- 15 Next scroll to and enter the initial values for **Proportional gain, Integral time, Derivative time** and **Mass**. These values are usually later modified based on empirical results until stable thermal performance is achieved.
- 16 Scroll to the **Sensor:** field and make sure the sensor type agrees with your hardware. To change the sensor type, press the [Mode\Type] key and select the RTD or TC (thermocouple with conversion board) type.
- 17 Scroll to **Power (watts)** and enter the power used to heat the device.
- 18 Scroll to **Maximum setpoint** and enter a value that is within the safe operating temperature of the device.
- 19 Scroll to **Maximum programming rate** and enter a rate of temperature rise per minute that your device is capable of achieving. This setting is not used with the SS zone control model which does not allow temperature programming.
- 20 Power cycle the GC to lock this configuration.
- 21 On the GC keypad press [Front Inlet] or [Back Inlet] to display the inlet's parameter list for entering operating values for a method.

## Configuring a custom inlet (Firmware A.01.07)

- 1 With the GC disconnected from its power source, install a PCM or PTV EPC module in the front or back inlet EPC module location. Plug the EPC module into the correct wiring harness connector. See EPC module location above.
- 2 Optionally install the inlet device containing the heater and temperature sensor.
- 3 Plug the heater and sensor cable from this device into one of the inlet's heated zone connections. See the "Heater and Cryo Installation" section above for the wiring harness connector to use.
- 4 Optionally install a cryo valve and plug it into the 24 volt power supply into the FV or BV connector.
- 5 With all covers installed, attach the power supply to the GC and power the GC on.
- 6 Unlock the GC configuration. Press the [Options] key, select **Keyboard & Display** and press the [Enter] key. Scroll down to **Hard Configuration Lock** and press the [off] button.
- 7 On the GC keypad press [Config] [Front Inlet] for the configuration of a front inlet device or [Config] [Back Inlet] for the configuration of a back inlet device.
- 8 With **Unconfigured:** selected press the [Mode/Type] key on the GC keypad and if necessary scroll to see all selections. See the table in the Hardware Configuration section for the correct display line to select. Press [Enter] to complete your selection.
- 9 When the GC prompts to power cycle the GC, reboot the GC by pressing the [Options] key and select **Communications**. From this option scroll down and select the reboot option and press [Enter] on the GC keypad. This soft reboot keeps the GC configuration unlocked for subsequent steps. Power cycling the GC forces you to unlock the GC configuration on power-up.
- 10 Select the configuration display for this inlet. The display now shows your configured hardware devices.
- 11 Scroll to **Gas type:** and press the [Mode/Type] key on this field to bring up the available carrier gas types. Select the gas type attached to the inlet and press [Enter].
- 12 If a cryo valve is included in this inlet, scroll to **Cryo type:** press the [Mode\Type] key and select the coolant supplied to the cryo valve.

When the GC prompts to power cycle the GC, reboot the GC by pressing the [Options] key and select **Communications**. From this option scroll down and select the reboot option and press [Enter] on the GC keypad. This soft reboot keeps the GC configuration unlocked for subsequent steps. Power cycling the GC forces you to unlock the GC configuration on power-up.

Select the configuration display for this inlet.

- 13 Set the **Cryo**, **Cryo timeout**, and **Cryo fault** parameters as required. See above.
- 14 Scroll to **Zone control model:** and press the [Mode/Type] key on this field to bring up the available zone control models. Select a suitable control model and press [Enter].

When the GC prompts to power cycle the GC, reboot the GC by pressing the [Options] key and select **Communications**. From this option scroll down and select the reboot option and press [Enter] on the GC keypad. This soft reboot keeps the GC configuration unlocked for subsequent steps. Power cycling the GC forces you to unlock the GC configuration on power-up.

Select the configuration display for this inlet.

- 15 To display all the configuration parameters select **View Custom Heater Conf** and press [On/Yes]. This displays all the heater parameters that can be set.
- 16 Next scroll to and enter the initial values for **Proportional gain**, **Integral time**, **Derivative time** and **Mass**. These values are usually later modified based on empirical results until stable thermal performance is achieved.
- 17 Scroll to **Power (watts)** and enter the power used to heat the device.
- 18 Scroll to **Maximum setpoint** and enter a value that is within the safe operating temperature of the device.
- 19 Scroll to **Maximum programming rate** and enter a rate of temperature rise per minute that your device is capable of achieving. This setting is not used with the SS zone control model which does not allow temperature programming.
- 20 Scroll to the **Sensor:** field and make sure the sensor type agrees with your hardware. To change the sensor type, press the [Mode\Type] key and select the RTD or TC (thermocouple with conversion board) type.
- 21 Power cycle the GC to lock this configuration.

- 22 On the GC keypad press [Front Inlet] or [Back Inlet] to display the inlet's parameter list for entering operating values for a method.

## Configuration procedure

- 1 With the GC disconnected from its power source, install the device containing the cryo focus valve.
- 2 Plug the 24 volt actuator cable from this device into the FV or BV connection.
- 3 With all covers installed, attach the power supply to the GC and power the GC on.
- 4 Unlock the GC configuration.
- 5 Press the [Config] key and select **cryo focus valve** from the list by scrolling to this entry and pressing [enter].
- 6 **For Firmware A.01.05:** Select the Cryo Focus valve attached to FV or BV.

**For Firmware A.01.07 through A.01.09:** With **Unconfigured:** selected, press [Mode/Type] and select the Cryo Focus valve attached to FV or BV.

- 7 After rebooting the GC, selecting the configuration for this cryo valve will show the valve properly associated to the FV or BV connector.
- 8 To display the method parameters for the cryofocus valve press the [Valve #] key. If other valves exist, scroll down to set the **Cryo focus equib time** and **Cryofocus delay time**.

### Configuring a custom inlet (Firmware A.01.08 and A.01.09)

- 1 With the GC disconnected from its power source, install a PCM or PTV EPC module in the front or back inlet EPC module location. Plug the EPC module into the correct wiring harness connector. See EPC module location above.
- 2 Optionally install the inlet device containing the heater and temperature sensor.
- 3 Plug the heater and sensor cable from this device into one of the inlet's heated zone connections. See the "Heater and Cryo Installation" section above for the wiring harness connector to use.
- 4 Optionally install a cryo valve and plug it into the 24 volt power supply into the FV or BV connector.
- 5 With all covers installed, attach the power supply to the GC and power the GC on.
- 6 Unlock the GC configuration. Press the [Options] key, select Keyboard & Display and press the [Enter] key. Scroll down to Hard Configuration Lock and press the [off] button.
- 7 On the GC keypad press [Config] [Front Inlet] for the configuration of a front inlet device or [Config] [Back Inlet] for the configuration of a back inlet device.
- 8 With Unconfigured: selected press the [Mode/Type] key on the GC keypad and if necessary scroll to see all selections. See the table in the Hardware Configuration section for the correct display line to select. Press [Enter] to complete your selection.
- 9 When the GC prompts to power cycle the GC, reboot the GC by pressing the [Options] key and select Communications. From this option scroll down and select the reboot option and press [Enter] on the GC keypad. This soft reboot keeps the GC configuration unlocked for subsequent steps. Power cycling the GC forces you to unlock the GC configuration on power-up.
- 10 Select the configuration display for this inlet. The display now shows your configured hardware devices.
- 11 Scroll to Gas type: and press the [Mode/Type] key on this field to bring up the available carrier gas types. Select the gas type attached to the inlet and press [Enter].
- 12 If a cryo valve is included in this inlet, scroll to Cryo type: press the [Mode\Type] key and select the coolant supplied to the cryo valve. Set the Cryo, Cryo timeout, and Cryo fault parameters as required. See above.

- 13 To display all the configuration parameters select View Custom Heater Conf and press [On/Yes]. This displays all the heater parameters that can be set.
- 14 Scroll to Zone control model: and press the [Mode/Type] key on this field to bring up the available zone control models. Select a suitable control model and press [Enter].
- 15 Next scroll to and enter the initial values for Proportional gain, Integral time, Derivative time and Mass. These values are usually later modified based on empirical results until stable thermal performance is achieved.
- 16 Scroll to the Sensor: field and make sure the sensor type agrees with your hardware. To change the sensor type, press the [Mode\Type] key and select the RTD or TC (thermocouple with conversion board) type.
- 17 Scroll to Power (watts) and enter the power used to heat the device.
- 18 Scroll to Maximum setpoint and enter a value that is within the safe operating temperature of the device.
- 19 Scroll to Maximum programming rate and enter a rate of temperature rise per minute that your device is capable of achieving. This setting is not used with the SS zone control model which does not allow temperature programming.
- 20 Power cycle the GC to lock this configuration.
- 21 On the GC keypad press [Front Inlet] or [Back Inlet] to display the inlet's parameter list for entering operating values for a method.

## Cryo focus valve configuration

Use the cryo focus valve when temperature control is not required and the run must start out with a relatively low temperature in a device which must be held for a fixed period of time after the run begins. This valve is programmed to remain open for a minimum time period prior to the run (**equib time**) and a fixed period of time after the run begins (**delay time**).

### Operation

The cryo focus valve control energizes the cryo valve actuator for a minimum programmable time before a run begins and a fixed amount of time after the run begins. These two time setpoints are entered by the user pressing the [Valve] key and then selecting the cryo focus valve from the valve list.

To start a run using manual injection, the user presses the [prep run] key to actuate the valve. The run cannot start until the **equib time** programmed elapses. The valve remains continually energized from the time [prep run] is pressed until the **delay time** programmed after the start run completes. The run starts when the user presses the [Start] key. Using an automatic liquid sampler or injection valve automates this process and it is only necessary for the user to press the [Start] key.

### Assigned connections

The cryo focus valve actuator is connected to the FV or BV 24 volt socket adjacent to the front and back inlet connections.



## Cryo trap configuration

This procedure applies to GCs with firmware version A.01.09. Firmware A.01.09 or greater is required.

This procedure assumes that the cryo trap is installed such that:

- the heater connects to the back inlet BI connector.
- the valve actuator connects to back inlet BV connector.

To configure the cryo trap:

- 1 Press **[Config]**, select **Aux Temp # > Thermal Aux 1**, press **[Enter]**.
- 2 Press **[Mode/Type]**, select **Install BINLET with BV Cryo**, press **[Enter]**.
- 3 Soft-reboot the GC. Press **[Options]**, select **Communications > Enter > Reboot GC > [On/Yes] > [On/Yes]**.
- 4 Press **[Config]**, select **Aux Temp # > Thermal Aux 1 > [Enter]**. Select **Auxiliary Type: Unknown**, press **[Mode/Type]**, select **User Configurable Heater**, press **[Enter]**.
- 5 Soft-reboot the GC. Press **[Options]**, select **Communications > Enter > Reboot GC > [On/Yes] > [On/Yes]**.
- 6 Press **[Config]**, select **Aux Temp # > Thermal Aux 1 > Cryo Type (Valve BV): N2**.

If **Cryo Type (Valve BV)** does not equal N2, then select **Cryo Type (Valve BV)**, press **[Mode/Type]**, select **N2 Cryo**, press **[Enter]** then **[Clear]**.

- 7 Configure Aux 1 by scrolling to each entry listed below and making the appropriate changes. After each step, press **[Clear]** to ignore the reboot request. You will reboot later.
  - a Proportional Gain: **5.30**, press **[Enter]** then **[Clear]**.
  - b Integral Time: **10**, press **[Enter]** then **[Clear]**.
  - c Derivative Time: **1.00**, press **[Enter]** then **[Clear]**.
  - d Mass (Watt-sec/deg): **18**, press **[Enter]** then **[Clear]**.
  - e Power (Watts):

To find the watts to set here, scroll down approximately 6 positions to Back Inlet Status (BINLET). As an example, it may read: **BINLET (OK): 125 Watts**. Use this reading (in this example 125 Watts) as the parameter for Power (Watts). Enter 125, press **[Enter]** then **[Clear]**.

**f** Cryo Control Mode > [Mode/Type]:

Cryo Control Type: PTV (ignore the first line, should already be PTV)

Cryo Control Type: Cryo Trap, press **[Enter]** then **[Clear]**.

**g** Zone Control Mode > [Mode/Type] > PTV, press **[Enter]** then **[Clear]**.

**h** Sensor > [Mode/Type] > Thermocouple, press **[Enter]** then **[Clear]**.

**i** Maximum Setpoint: 400, press **[Enter]** then **[Clear]**.

**j** Maximum Programming Rate: 720, press **[Enter]** then **[Clear]**.

**8** Soft-reboot the GC. Press **[Options]**, select **Communications** > **Enter** > **Reboot GC** > **[On/Yes]** > **[On/Yes]**.

## Detector configuration

Use the following 3 tables to determine the component wiring connections for the hardware configuration you are designing. See the *Agilent Advanced User guide* for configuring Agilent standard detector types.

## FID Signal 1 Detector Hardware Connections

The following table lists the hardware connections for an Agilent standard FID. Additionally listed is the hardware connections required for a channel partner primary FID signal board with or without a secondary FID or AIB signal board.

GC configuration display selection	EPC wire harness connection	Heater wire harness connection	Signal 1 Detector Board wire harness connection	Signal 2 Detector Board wire harness connection
Front detector				
FID	EPC3	F-DET	F-DET	None
FID + AIB	None	None	F-DET	AUX DET1
CPDET FID, No htr, No EPC	None	None	F-DET	None
CPDET FID & Htr, No EPC	None	F-DET	F-DET	None
CPDET FID & EPC, No Htr	EPC3	None	F-DET	None
CPDET FID & EPC & Htr	EPC3	F-DET	F-DET	None
CPDET FID & AIB & Htr, No EPC	None	F-DET	F-DET	AUX DET1
CPDET FID & AIB & EPC, No htr	EPC3	None	F-DET	AUX DET1
CPDET FID & AIB & Htr & EPC	EPC3	F-DET	F-DET	AUX DET1
CPDET FID & FID & Htr, No EPC	None	F-DET	F-DET	AUX DET1
CPDET FID & FID & EPC, No htr	EPC3	None	F-DET	AUX DET1
CPDET FID & FID & Htr & EPC	EPC3	F-DET	F-DET	AUX DET1
Back detector				
FID	EPC4	B-DET	B-DET	None
FID + AIB	None	None	B-DET	AUX DET2
CPDET FID, No htr, No EPC	None	None	B-DET	None

<b>GC configuration display selection</b>	<b>EPC wire harness connection</b>	<b>Heater wire harness connection</b>	<b>Signal 1 Detector Board wire harness connection</b>	<b>Signal 2 Detector Board wire harness connection</b>
CPDET FID & Htr, No EPC	None	B-DET	B-DET	None
CPDET FID & EPC, No Htr	EPC4	None	B-DET	None
CPDET FID & EPC & Htr	EPC4	B-DET	B-DET	None
CPDET FID & AIB & Htr, No EPC	None	B-DET	B-DET	AUX DET2
CPDET FID & AIB & EPC, No htr	EPC4	None	B-DET	AUX DET2
CPDET FID & AIB & Htr & EPC	EPC4	B-DET	B-DET	AUX DET2
CPDET FID & FID & Htr, No EPC	None	B-DET	B-DET	AUX DET2
CPDET FID & FID & EPC, No htr	EPC4	None	B-DET	AUX DET2
CPDET FID & FID & Htr & EPC	EPC4	B-DET	B-DET	AUX DET2

In the above table, FID indicates a detector signal board with igniter control support. CPDET FID indicates a channel partner FID detector signal board without support for an igniter.

## AIB Signal 1 Detector Hardware Connections

The following table lists the hardware connections required for a channel partner primary AIB signal board with or without a secondary FID or AIB signal board. The first two entries under Front Detector and Back Detector entries are for the standard legacy Agilent AIB installations.

GC configuration display selection	EPC wire harness connection	Heater wire harness connection	Signal 1 Detector Board wire harness connection	Signal 2 Detector Board wire harness connection
<b>Front detector</b>				
Install AIB with No Heater	None	None	F-DET	None
Install AIB with Heater	None	F-DET	F-DET	None
CPDET AIB, No htr, No EPC	None	None	F-DET	None
CPDET AIB & Htr, No EPC	None	F-DET	F-DET	None
CPDET AIB & EPC, No Htr	EPC3	None	F-DET	None
CPDET AIB & EPC & Htr	EPC3	F-DET	F-DET	None
CPDET AIB & AIB & Htr, No EPC	None	F-DET	F-DET	AUX DET1
CPDET AIB & AIB & EPC, No htr	EPC3	None	F-DET	AUX DET1
CPDET AIB & AIB & Htr & EPC	EPC3	F-DET	F-DET	AUX DET1
CPDET AIB & FID & Htr, No EPC	None	F-DET	F-DET	AUX DET1
CPDET AIB & FID & EPC, No htr	EPC3	None	F-DET	AUX DET1
CPDET AIB & FID & Htr & EPC	EPC3	F-DET	F-DET	AUX DET1

GC configuration display selection	EPC wire harness connection	Heater wire harness connection	Signal 1 Detector Board wire harness connection	Signal 2 Detector Board wire harness connection
<b>Back detector</b>				
Install AIB with No Heater	None	None	B-DET	None
Install AIB with Heater	None	B-DET	B-DET	None
CPDET AIB, No htr, No EPC	None	None	B-DET	None
CPDET AIB & Htr, No EPC	None	B-DET	B-DET	None
CPDET AIB & EPC, No Htr	EPC4	None	B-DET	None
CPDET AIB & EPC & Htr	EPC4	B-DET	B-DET	None
CPDET AIB & AIB & Htr, No EPC	None	B-DET	B-DET	AUX DET2
CPDET AIB & AIB & EPC, No htr	EPC4	None	B-DET	AUX DET2
CPDET AIB & AIB & Htr & EPC	EPC4	B-DET	B-DET	AUX DET2
CPDET AIB & FID & Htr, No EPC	None	B-DET	B-DET	AUX DET2
CPDET AIB & FID & EPC, No htr	EPC4	None	B-DET	AUX DET2
CPDET AIB & FID & Htr & EPC	EPC4	B-DET	B-DET	AUX DET2

## Aux Detector Hardware Connections

GC configuration display selection	EPC wire harness connection	Heater wire harness connection	Signal 1 Detector Board wire harness connection	Signal 2 Detector Board wire harness connection
<b>Aux Detector 1</b>				
TCD	EPC5	A1 or A2 or BI/BINLET	AUX DET1	None
AIB	None	None	AUX DET1	None
<b>Aux Detector 2</b>				
AIB	None	None	AUX DET2	None

The Aux Detector 1 TCD installation requires an extension cable to reach the EPC module and detector board located in a side mounted enclosure. The heater and sensor for the Aux TCD also requires an extension cable.

## EPC module installation

Optionally install and attach your detector device's EPC module to the GC's wiring harness.

EPC wiring harness label	Use
EPC3	Front detector EPC module
EPC4	Back detector EPC module
EPC5	3 <sup>rd</sup> detector TCD EPC module



## Detector Board Installation

Wiring harness label	Use
F-DET	Front detector board
B-DET	Back detector board
AUX DET1	3rd detector TCD board or AIB
AUX DET2	Aux detector board AIB

## Heater installation


Optionally install and attach your detector device's heater to the GC's wiring harness.

Wiring harness label	Use
F-DET	Front detector heater
B-DET	Back detector heater
A1	3 <sup>rd</sup> detector heater
A2	3 <sup>rd</sup> detector heater

## Configuring a custom detector

- 1 With the GC disconnected from its power source, remove the electronics side panel and install the detector boards in the detector slots. See the section on configuration and plug the correct communications cable into the board.
- 2 Install the detector optionally containing a heater and temperature sensor. Plug the heater and sensor cable from this device into one of the six small heated zone connections. See the "Heater Installation" section above for the wiring harness connector to use.
- 3 Optionally install an FID EPC module in the front, back, or aux detector EPC module location. Plug the EPC module into the correct wiring harness connector. See EPC module location above.
- 4 With all covers installed, attach the power supply to the GC and power the GC on.
- 5 Unlock the GC configuration. Press the [Options] key, select **Keyboard & Display** and press the [Enter] key. Scroll down to **Hard Configuration Lock** and press the [off] button.
- 6 On the GC keypad press [Config] [Front Det] for the configuration of a front detector device or [Config] [Back Det] for the configuration of a back detector device or [Config] [Aux Det #] for the configuration of a Aux detector device.
- 7 With **Unconfigured:** selected press the [Mode/Type] key on the GC keypad and if necessary scroll to see all selections. Select your hardware configuration from the list and press [Enter] to complete your selection.
- 8 When the GC prompts to power cycle the GC, reboot the GC by pressing the [Options] key and select **Communications**. From this option scroll down and select the reboot option and press [Enter] on the GC keypad. This soft reboot keeps the GC configuration unlocked for subsequent steps. Power cycling the GC forces you to unlock the GC configuration on power-up.
- 9 Select the configuration display for this detector. The display now shows your configured devices.
- 10 For an FID installation, scroll to **Gas type:** and press the [Mode/Type] key on this field to bring up the available makeup gas types. Select the gas type and press [Enter].
- 11 Power cycle the GC to lock this configuration.

- 12 On the GC keypad press [Front Det] or [Back Det] or [Aux Det #] to display the detector's parameter list for entering operating values for a method.

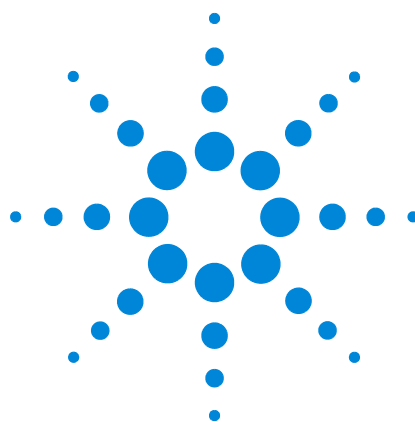

## 15 Service Mode Diagnostics

Service mode overview	553
Detector diagnostics	554
Signal board diagnostics	555
Multiplexed ADC	556
Pneumatics	557
Power diagnostics	558



## Service mode overview

Service mode generates a series of tables that include information that can help with diagnosing problems. For example, the detector signal board diagnostic includes information to verify a short in an FID ignitor. It also contains a measurement of the Noise in the signal output.

To access the service mode diagnostics:

- 1 Press [**Service Mode**].
- 2 Scroll to Diagnostics and press [**Enter**].

The list of diagnostics that are available vary with your configuration. Here is an example of a 7890A GC with 2 inlets, 2 detectors, an AUX module, and a PCM module.

```
DIAGNOSTICS
Instrument status
Front inlet
Back inlet
Front detector
Back detector
Auxiliary pressure (1,2,3)
Pcm B
Thermal PID control
Electronics
```

## Detector diagnostics

For example, the detector signal board diagnostic includes information to verify a short in an FID ignitor. It also contains a measurement of the Noise in the signal output.

To access the detector diagnostics:

- 1 Press [**Service Mode**].
- 2 Scroll to **Diagnostics** and press [**Enter**].
- 3 Scroll to **<Front, Back, or AUX detector>** and press [**Enter**].

The list of diagnostics that are available vary with your configuration. Here is an example of an FID detector.

### **FRONT FID DIAGNOSTICS**

Detector Signal  
Pneumatics module  
Hydrogen gas  
Air gas  
Makeup gas

## Signal board diagnostics

Service mode is a series of tables that include information that can help with diagnosing problems. For example, the detector signal board diagnostic includes information to verify a short in an FID ignitor. It also contains a measurement of the Noise in the signal output.

To access the detector signal board diagnostics:

- 1 Press [**Service Mode**].
- 2 Scroll to **Diagnostics** and press [**Enter**].
- 3 Scroll to **<Front, Back, or AUX detector>** and press [**Enter**].
- 4 Scroll to **Detector Signal** and press [**Enter**].

The list of diagnostics that are available vary with your detector. Here is an example of an FID detector.

### FRONT SIGNAL (FID)

Signal Value		0.0
Noise		0
High Voltage		12.9V
Ignitor:	0.00V	0.0A
ADC Offset:	60.1mV	Ref: 2.48V
Input disable		Off
Ignitor voltage		Off
Mode:		Normal

## Multiplexed ADC

This information is based on Service Mode | Diagnostics | Electronics | Multiplexed ADC.

### MULTIPLEXED ADC

Definition	Oven rtd
Reading noise	159.8
Reading value	1801821
Offset	54907
Gain counts/volt	4173672
Gain counts/ohm	14607
Reference	95.446
Temperature	
Reading value	1801887
Reading value	-95
Reading value	-180
Reading value	198
Reading value	64
Reading value	-45

There is a multiplexed analog-to-digital circuit (ADC) for each of the following resistance temperature detectors (rtd):

- Oven rtd reading
- Front detector rtd reading
- Back detector rtd reading
- Front inlet rtd reading
- Back inlet rtd reading
- Auxiliary 1 rtd reading
- Auxiliary 2 rtd reading
- Line Sense reading
- Temperature ref reading
- Atmospheric pressure reading
- Minus 24V power supply rdg
- Logic current reading
- 24 volt poser supply reading
- Valve current reading

To view the diagnostics for a particular reading,

- 1 With the cursor on the Definition line, press **[Mode/Type]**.
- 2 Scroll to the reading that you want and press **[Enter]**.



# Pneumatics

This information is based on Service Mode | Diagnostics | Electronics | Pneumatics.

<b>PNEUMATICS</b>	
Atmospheric pressure	14.715
24 Volt valve drive	26.71
Zero all pressure sensors	

To zero all of the pressure sensors, do the following.


- 1 Set all flows and pressures to zero.
- 2 Disconnect the supply gases.
- 3 Wait 30 minutes.
- 4 Press [**Enter**] and confirm.

## Power diagnostics

This information is based on Service Mode > Diagnostics > Electronics > Power Info.

### POWER DIAGNOSTICS

Front inlet	166.5 watts
Front detector	70.1 watts
Auxiliary 1	0.9 watts
Back inlet	73.9 watts
Back detector	63.6 watts
Auxiliary 3	65.4 watts
Logic current	1.370 amps
Valve current	0.326 amps
24 Volts	26.66 volts
-24 Volts	-27.03 volts
Zone supply	42.92 volts

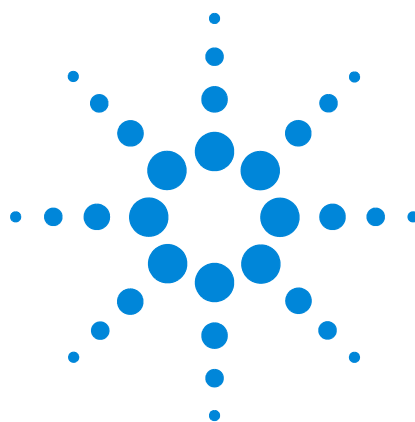

## 16 Accessories and options

For convenience, the table below lists the accessories and options for the 7890 Series GCs. **Please refer to ELSA for the latest and most accurate selections.**

**Table 37** 7890 Series GC accessories and options

Accessory	Option	Description	Notes
<b>Inlets</b>			
G3451A	102	Packed inlet with EPC, septum purge. 0–100 psi operation. Configured for 0.53-mm columns.	Incompatible with MSD software G1701CA or later. Includes liners and adapters/connectors for 0.53-mm, 1/8-inch, and 1/4-inch columns.
G3452A	112	0–100 psi Split/splitless with EPC.	
G3460A	113	0–150 psi Split/splitless with EPC.	Select for use with 0.100 id columns. Not optimal for 0.530 id columns.
G3454A	122	On-column capillary inlet with EPC. 0–100 psi temperature-programmed, on-column capillary inlet system with EPC.	Order a cryo accessory kit for use with subambient oven temperatures or for faster cycle time to near ambient. Order 19245-60590 for manual injection with a fused silica needle.
G3510A		<b>7890A only</b> Multimode inlet with liquid CO <sub>2</sub> cooling	Only one type of cryogenic source (LN <sub>2</sub> or LC02) is allowed per instrument regardless of number of cryogenic accessories. Not compatible with PTV (or 2nd MMI on 7890A). Incompatible with 2nd FPD, #651, #887 or #221.
G3511A		<b>7890A only</b> Multimode inlet with LN <sub>2</sub> or air cooling.	



**Table 37** 7890 Series GC accessories and options

Accessory	Option	Description	Notes
G3510B		<b>7890B only</b> Multimode inlet with liquid CO <sub>2</sub> cooling	Only one type of cryogenic source (LN <sub>2</sub> or LC02) is allowed per instrument regardless of number of cryogenic accessories. Not compatible with PTV (or 2nd MMI on 7890A). Incompatible with 2nd FPD, #651, #887 or #221.
G3511B		<b>7890B only</b> Multimode inlet with LN <sub>2</sub> or air cooling.	
G3500A	130	Cap PTV Inlet, with EPC, Septumless, LC02 0–100 psi, septumless head, LC02 cryo.	Temperature range – 60 °C to 450 °C with 3 ramps and 3 operating modes. Split, Splitless, solvent venting. Only one PTV inlet allowed per GC. If using oven cooling, it must also be LC02 cooled. CTC and septumless PTV requires a 23" syringe (not included).
G3501A	131	Capillary PTV Inlet with EPC, Septum Head, LC02 0–100 psi, w/septum head, LC02 cryo.	Temperature range –60 to 450 C with 3 ramps and 3 operating modes. Split, Splitless, solvent venting. Only one PTV inlet allowed per GC. If using oven cooling, it must also be LC02 cooled.
G3502A	132	Capillary PTV Inlet with EPC, Septumless, LN <sub>2</sub> 0–100 psi, septumless head, LN <sub>2</sub> cryo.	Temperature range –160 C to 450 C with 3 ramps and 3 operating modes. Split, Splitless, solvent venting. Only one PTV inlet allowed per GC. If using oven cooling, it must also be LN <sub>2</sub> cooled. CTC and septumless PTV requires a 23" syringe (not included).
G3503A	133	Capillary PTV Inlet with EPC, Septum head, LN <sub>2</sub> 0–100 psi, septum head, LN <sub>2</sub> cryo.	Temperature range –160 C to 450 C with 3 ramps and 3 operating modes. Split, Splitless, solvent venting. Only one PTV inlet allowed per GC. If using oven cooling, it must also be LN <sub>2</sub> cooled.
G3504A	142	Volatiles interface with EPC. 0–100 psi. Split, splitless, or direct injection. Deactivated.	Calibrated for use with H <sub>2</sub> or He carrier gas only.
G2399A		Solvent Vapor Exit Accessory for use on Cool on column inlet	Includes retention gap, pre-column, analytical column, valve and 50 uL autosampler syringe.

**Table 37** 7890 Series GC accessories and options

Accessory	Option	Description	Notes
<b>Detectors</b>			
G3461A	210	FID with EPC. Flame ionization detector with EPC for packed and capillary columns.	Includes adapters and connectors for capillary, 1/8-inch and 1/4-inch columns.
G3462A	211	FID with EPC, optimized for capillary columns only.	Includes digital control of all detector gases. Includes connectors for capillary columns.
G3432A	220	TCD with EPC. Thermal conductivity detector with EPC.	Includes adapters/connectors for capillary, 1/8-inch and 1/4-inch columns.
G3437A	221	TCD with EPC mounted on left side of GC.	Includes adapters/connectors for capillary, 1/8-inch and 1/4-inch columns. Incompatible with MSD or CTC
G3437B	221	<b>7890B Only.</b> Side-mounted TCD with EPC mounted on left side of GC.	Includes adapters/connectors for capillary, 1/8-inch and 1/4-inch columns. Incompatible with MSD.
G2397AD	231	Micro Electron Capture Detector with EPC	Includes adapters/connectors for 1/4-in. and capillary columns.
G2398AD		Micro Electron Capture Detector with EPC for Japan.	Includes adapters/connectors for 1/4-in. and capillary columns.
G3435A	240	Single flame photometric detector (FPD) with EPC.	For use with packed and capillary columns (1/4-in. glass columns require adapter p/n G1532-20710, not included). Requires two heated zones.
G3435B	240	FPD Plus with EPC. Single flame photometric detector with EPC.	<b>7890A+ and 7890B only.</b> For use with packed and capillary columns (1/4-in. glass columns require adapter p/n G1532-20710, not included). Requires: <ul style="list-style-type: none"> <li>• Two heated zones.</li> <li>• G3444A PCB enhancement kit.</li> <li>• A.01.05 OLCDS, B.07.00 MassHunter or higher.</li> </ul>

**Table 37** 7890 Series GC accessories and options

Accessory	Option	Description	Notes
G3436A	241	Dual flame photometric detector (FPD) with EPC with make-up gas.	For use with packed and capillary columns. For use with packed and capillary columns (1/4" glass columns require adapter p/n G1532-20710, not included). Not Compatible with any heated valve options. Requires two heated zones and two detector circuit board slots. Not compatible with rear ALS.
G3436B	241	Dual flame photometric detector (FPD) with EPC with make-up gas.	<b>7890A+ and 7890B only.</b> Not compatible with any heated valve box options. For use with packed and capillary columns (1/4-in. glass columns require adapter p/n G1532-20710, not included). Requires: <ul style="list-style-type: none"> <li>• Two heated zones. and two detector board slots.</li> <li>• G3444A PCB enhancement kit.</li> <li>• A.01.05 OLCDS, B.07.00 MassHunter or higher.</li> </ul>
G3495A	250	NPD (white ceramic bead) for packed and capillary columns.	For use with packed and capillary columns. Includes adapters and connectors for capillary and 1/4-inch columns.
G3496A	251	NPD (white ceramic bead) optimized for capillary columns only.	Includes connectors for capillary columns.
G3463A	252	NPD (Blos bead) for packed and capillary columns.	For use with packed and capillary columns. Includes adapters and connectors for capillary and 1/4-inch columns.
G3464A	253	NPD (Blos bead) optimized for capillary columns only (not for use with packed columns).	Includes connectors for capillary columns.
G3490A	201	MSD interface.	For use with Agilent 5975 Series MSD or Agilent 7000A Triple Quad GC/MS or 7200 Q-TOF GC/MS). Uses one heated zone.
G3456A	206	Analog input board. Single-channel analog input board to accommodate an external detector.	Mounts in detector board position. Cable included.
19232C	925	TCD to FID Series Adapter	Requires FID for packed columns (#210 or G3461A).

**Table 37** 7890 Series GC accessories and options

Accessory	Option	Description	Notes
<b>Pneumatics</b>			
G3470A	301	Aux EPC.	Provides three channels of auxiliary 0–100 psi electronic pneumatics control.
G3471A	309	Pneumatics Control Module. Provides two channels of EPC. One channel is either flow or pressure, the second is pressure only.	Occupies one of two aux slots and/or one of two inlet slots. Maximum of three.
G3469A	205	EPC headspace interface.	Required for the EPC-ready Agilent 7694 headspace sampler (G1289B). Optional for the G1290B and G1888A headspace sampler.
<b>CFTs</b>			
G2855B	888	Capillary Flow Technologies Dean Switch Requires two detectors, two columns and either PCM (preferred) or aux EPC.	Enables "heartcutting". Switch includes high temperature ferrules and calculation SW CD.
G3180B	889	Capillary Flow Technologies Splitter with Makeup Gas	Requires either aux epc or PCM. Typical application is to split column effluent to two detectors. Includes hw and Calculator SW. Requires MS Excel.
G3183B	890	Capillary Flow Technologies 3 way splitter with Makeup Gas.	Typical application is to split column effluent to three detectors. Includes hw and Calculator SW. Requires either aux epc or PCM. Requires MS Excel.
G3181B		Capillary Flow Technologies Splitter, unpurged	Enables leak free splitting of column effluent to two detectors. Does not include installation.
G3184B		CFT Unpurged Tee for MSD Transfer Line.	
G3486A	887	CFT Flow Modulator Accessory requires 7890A GC w/SSL inlet, FID , PCM, ALS injector, G2070BA SW (Ver B.03.02).	2D Data analysis SW is neither provided nor supported by Agilent. Also requires G3487A Checkout Kit.
G3487A		CFT Flow Modulator Checkout Kit, Includes various parts for checkout.	Required use with 7890A CFT Flow Modulator ( G3440A#887 or G3486A).

**Table 37** 7890 Series GC accessories and options

Accessory	Option	Description	Notes
G1472A		Purged Ultimate Union kit for GC/MS system with mounting and instructions. Installation and Familiarization are included. Optional installation service is not required.	
G3182-61580		CPM Union kit	High temperature, inert, leak free column connector with oven mounting bracket
G1580-61060		CPM-1/16" Tube Assembly Kit	This kit provides leak-free, inert connection from GC valve to CPM system
<b>Cryo</b>			
G3465A	311	Cryogenic oven and inlet cooling w/ CO2	Provides cooling for oven and up to two cool-on-column inlets if installed. Not compatible with GC QQQ (7000 MS).
G3467A	311	Cryogenic oven cooling w/CO2	Provides cooling for oven only. If PTV inlet is present, it must also be LC02 cooled. Not compatible with GC QQQ (7000 MS).
G3468A	312	Cryogenic oven and inlet cooling w/ LN2	Provides cooling for oven and up to two cool-on-column inlets if installed. Not compatible with GC QQQ (7000 MS).
G3466A	312	Cryogenic oven cooling w/ LN2	Provides cooling for oven only. If PTV inlet is present, it must also be LN2 cooled. Not compatible with GC QQQ (7000 MS).
<b>Miscellaneous</b>			
G3478A	307	Adds a nickel catalyst (methanizer)	Requires availability of rear injection port position. If two injection ports are required, contact your sales assist engineer for alternatives. Compatible with PCM in rear position. Not compatible with capillary-only FID.
G1530-80650	306	Exhaust deflector	
G3485A		Manual Pneumatics Carrier	<b>7890A and 7890A+Only</b> Fits on left side of 7890A for housing unheated valves. Incompatible with MSD.



**Table 37** 7890 Series GC accessories and options

Accessory	Option	Description	Notes
G3485B		Manual Pneumatics Carrier	<b>7890B Only</b> Fits on left side of 7890A for housing unheated valves. Incompatible with MSD.
G2887BA		Simulated Distillation software for ChemStation (included with #651). Requires ChemStation. See G3506A below.	<b>7890B Only</b>
G6586AA		Agilent LTM Control Software. Provides software control for up to four LTM modules on one Agilent 7890/6890 GC.	<b>7890B Only</b>
G2646-60500		Oven Insert Accessory for fast oven temperature ramping. Cuts the oven volume in half.	Only compatible when using rear inlet/detector. If using MSD, inlet and MSD must be in rear position.
G3444A		7890A Logic Board Enhancements Kits.	<b>7890A Only</b> Makes 7890A ready for FPD+, Integrated Intelligence features, serial port enhancements. Requires software upgrades for A.01.05 OLCDS and MassHunter B.07.00 or higher.
G3494A		Bar Code Reader with USB Connection for 7890B	7890A requires: <ul style="list-style-type: none"> <li>• G3444A enhancement kit</li> <li>• A.01.05 OpenLAB or later, or B07.00 MassHunter or later</li> </ul>
G3494B		Bar Code Reader with Serial Connection for 7890B	7890A requires: <ul style="list-style-type: none"> <li>• G3444A enhancement kit</li> </ul>
<b>Valves</b>			
G1580A	751	Heated, automated valve box for one valve	
G1581A	752	Heated, automated valve box for two valves	
G3480-67585	860	Kit to interface an upstream gas or liquid sampling valve to a capillary or Multi-Mode inlet for sample splitting and introduction into a capillary or megabore column.	
G2739A		Liquid sampling valve plumbing kit	
G2740A		Gas valve plumbing kit	
G2741A		Adds an adjustable restrictor to control gas flows.	Max temp of 225 °C.

**Table 37** 7890 Series GC accessories and options

Accessory	Option	Description	Notes
G2742A	874	Kit to interface a capillary inlet or multimode inlet to a valve located downstream	
G2743A		7890A actuator kit. Adds automation to a 7890 valve. It consists of an actuator and solenoid valve.	
G2743B		7890 Actuator Kit. Adds automation to 7890 valve system for 7890 with S/N = CN10201096 or later, or for 7890 S/N = US10281001 or later.	
G2744A		Adds a heated zone to a 7890 valve box	
G2748A	861	Hardware kit to mount an unheated gas sampling or liquid sampling valve on the left side of the oven.	Incompatible with MSD. Requires side carrier.
G3478A	307	Adds a nickel catalyst (methanizer) to a 7890 Gas Chromatograph	Requires availability of rear injection port position. If two injection ports are required, contact your sales assist engineer for alternatives. Compatible with PCM in rear position.
G3505A		High Pressure Injection Device Accessory. Requires SSL inlet and one flow source from Aux EPC or PCM module. User must provide appropriate hardware to interface sample to HPID.	Only one can be installed, front position only.
G4338A	872	0.53 mm Column to Valve Interface	Requires mounting plate 05890-80660. Each plate can hold 5 interface adaptors.
G4339A	873	0.32 mm Column to Valve Interface	Requires mounting plate 05890-80660. Each plate can hold 5 interface adaptors.



## 17 Printed circuit boards

Summary of PCB Part Numbers	568
AC power board	569
Analog and power board	570
Logic and communications board (7890A)	572
Logic and communications board (7890B)	573
FID signal board	574
TCD signal board	575
uECD signal board	576
FPD signal board	577
NPD signal board	578
FID electrometer board	579
uECD electrometer board	580
NPD electrometer board	581
Analog input board	582

This section describes the following circuit boards and their connections.

Some connectors are marked MTC, which stands for Manufacturing Test Connector. There may or may not be a connector loaded on the board. MTCs are used at the factory to connect to test equipment so that the board can be exercised without the need for a GC. These connectors are not used when the board is installed in an instrument.




## Summary of PCB Part Numbers

Description	Part number
AC power board	G3430-60050
Analog and power board	G3430-60151 (replaces G3430-60150)
Logic and communications board	G3430-60120 (replaces G3430-60101 and G3430-60100)
Keyboard/display connection board	G3430-60003
FID signal board	G3431-60020
TCD signal board	G3432-60010
uECD signal board	G3433-61020 (replaces G3433-60020)
FPD signal board	G3435-60010
NPD signal board	G3434-60021 (replaces G3434-60020)
FID electrometer board	G3431-60010
uECD electrometer board	G3433-60010
NPD electrometer board	G3434-60010
Analog input board	G3456-60010 (obsolete) shown, G3433-61020 has similar connectors)

# AC power board


Part number G3430-60050

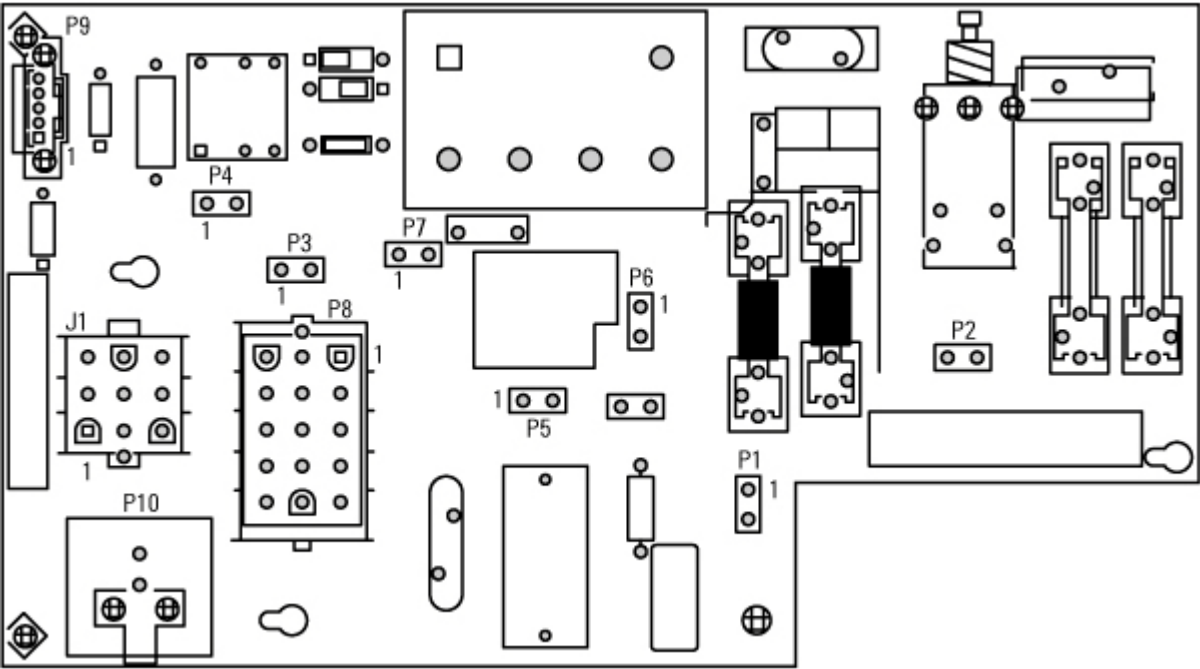

Connector	Function	See also
J1	Supplies power to transformer primary windings.	<a href="#">“Main transformer”</a>
P1 and P2	Supply line power to the AC board.	<a href="#">“AC power board schematic”</a>
P3 and P4	Supply power to the oven heater.	
P5, P6 and P7	These three connectors take jumpers to the oven triac. Each jumper connects to the nearest connector on the triac, which must be properly oriented.	<a href="#">“Replacing the oven triac”</a>
P8	This takes the power configuration plug, which must have the correct jumpers for the voltage supplied.	<a href="#">“Line voltage configuration plug”</a>
P9	Oven control from logic board.	
P10	MTC	

## Analog and power board


Part number G3430-60150 (shown). Connectors for the G3430-60151 are in the same location.

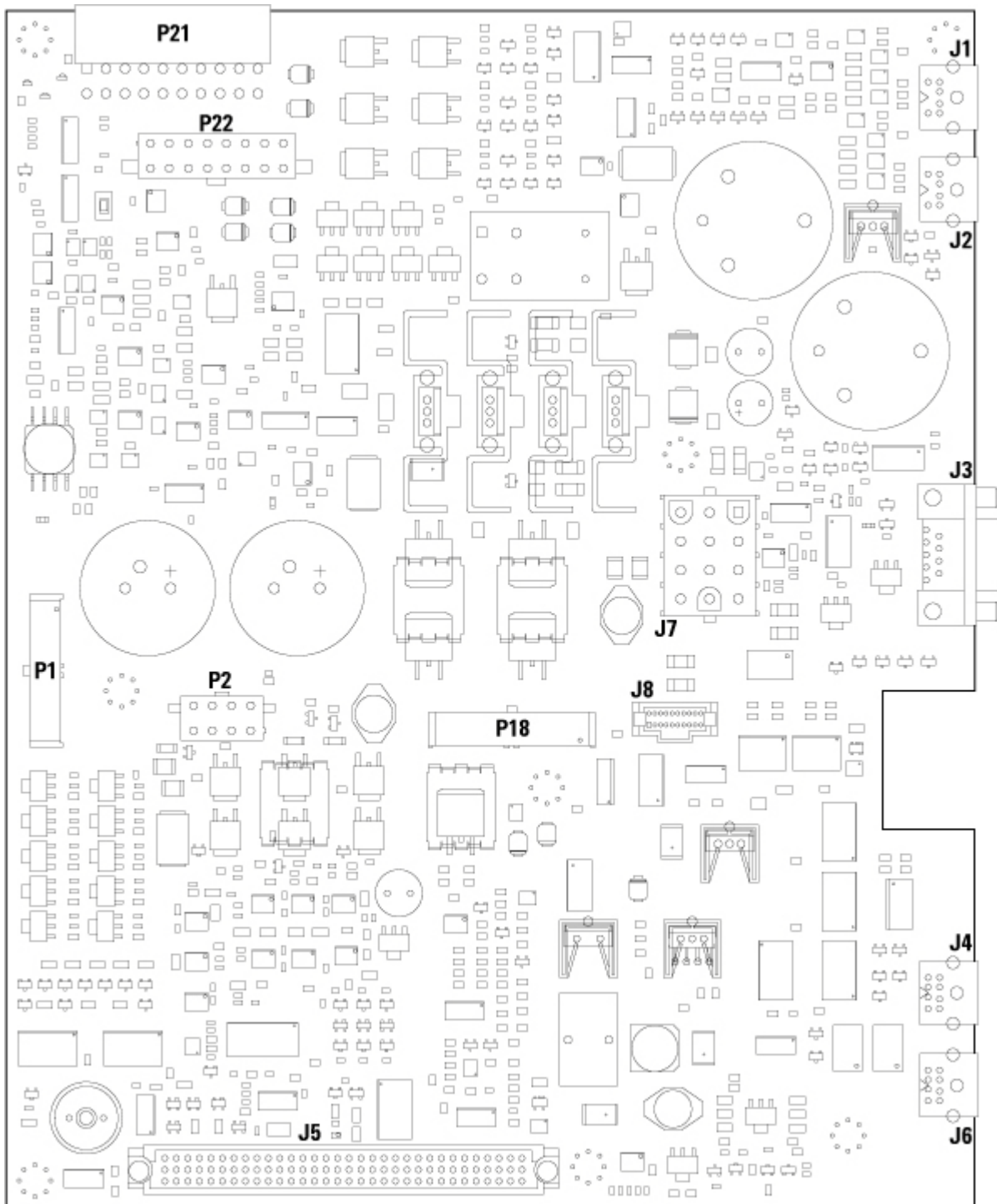

Connector	Function	See also
J1	SIG 1. Analog output for signal 1.	<a href="#">“Analog signal outputs”</a>
J2	SIG 2. Analog output for signal 2.	<a href="#">“Analog signal outputs”</a>
J3	REMOTE. Synchronize start/stop functions for up to 10 devices.	<a href="#">“Remote Start/Stop”</a>
J4	EVENT. Control external devices, either active (24 VDC) or passive )contact closures).	<a href="#">“External event”</a>
J5	96-pin data connection to analog and power board.	
J6	BCD. Accept BCD-coded input from selection valves and other devices.	
J7	Receive power from transformer secondary.	<a href="#">“Main transformer”</a>
J8	Ribbon cable to automatic liquid sampler board.	<a href="#">“RS-232 cable, ALS controller”</a>
P1	Supply power to the Communication harness.	<a href="#">“Communication harness”</a>
P2	Supply power to the oven fan.	
P18	Connects to the Motor Harness and AC Control cable.	<a href="#">“Motor Harness &amp; AC Control Assembly”</a>
P21	Connects to Inlet/Detector harness.	<a href="#">“Inlet/Detector harness”</a>
P22	Connects to the Auxiliary Zone/Valve Box Harness.	<a href="#">“Auxiliary zone/Valve box harness”</a>

Logic and communications board (7890A)

Part number G3430-60100 shown. Connectors for the G3430-60120 and G3430-60101 boards are in the same locations.

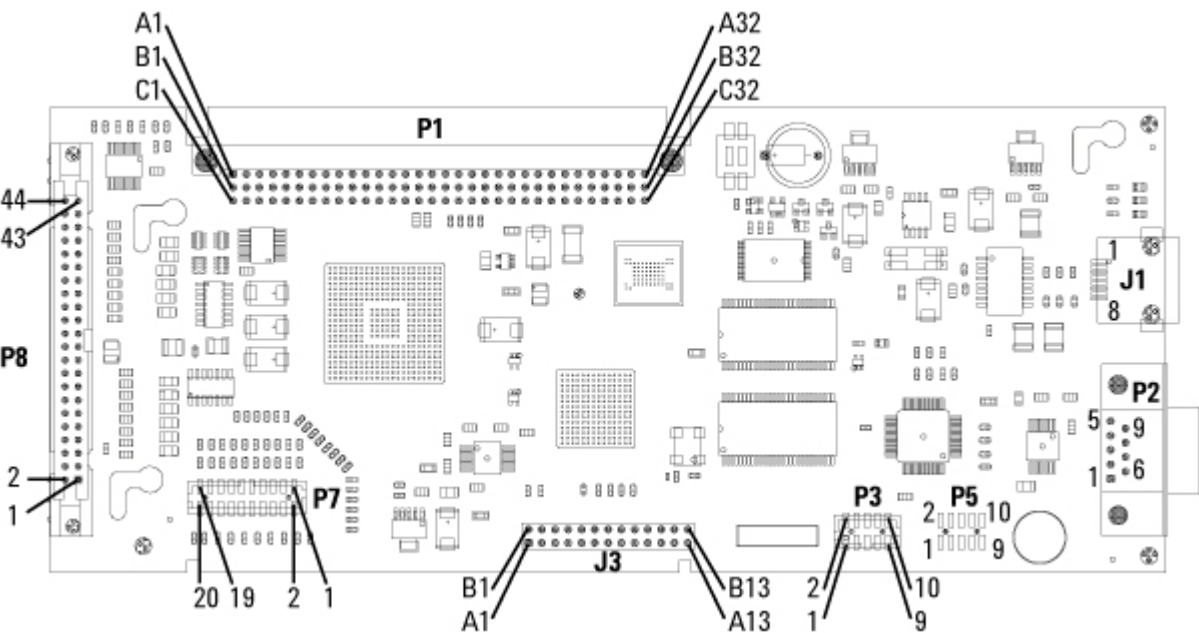


Connector	Function	See also
J1	External LAN port.	
J3	MTC	
P1	96-pin buss to analog and power board.	
P3+	MTC	
P5	MTC	
P7	Connects to Comm Buss on Communications Harness.	"Communication harness"
P8	44-pin ribbon cable to keyboard and display.	"Keyboard/Display harness", "Keyboard interconnect board"



Logic and communications board (7890B)

Part number G3430-61010 shown.


NEED GRAPHIC FOR 7890B

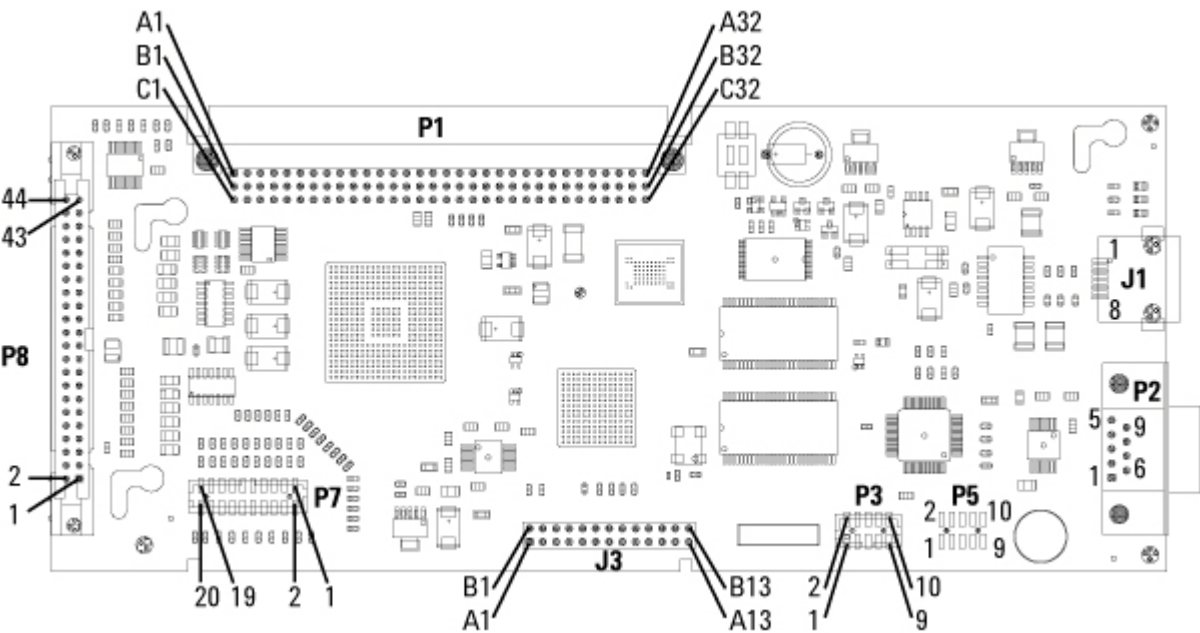

Connector	Function	See also
J1	External LAN port.	
J3	MTC	
P1	96-pin buss to analog and power board.	
P2	96-pin buss to analog and power board.	
P3+	MTC	
P4	96-pin buss to analog and power board.	
P6		
P7	Connects to Comm Buss on Communications Harness.	<a href="#">“Communication harness”</a>
P8	44-pin ribbon cable to keyboard and display.	<a href="#">“Keyboard/Display harness”</a> , <a href="#">“Keyboard interconnect board”</a>

FID signal board


Part number G3431-60020

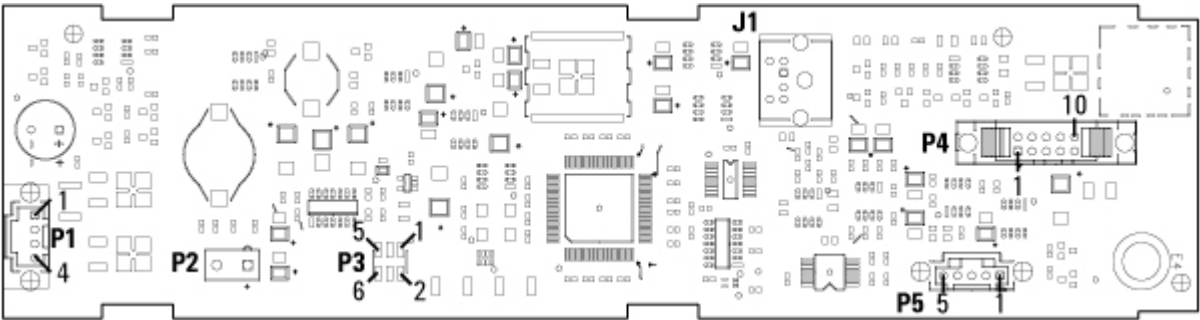

Connector	Function	See also
J1	MTC	
P1	Communications buss.	<a href="#">“Communication harness”</a>
P2	Ignitor.	
P3	MTC	
P4	Signal from electrometer.	
P5	MTC	

TCD signal board


Part number G3432-60015

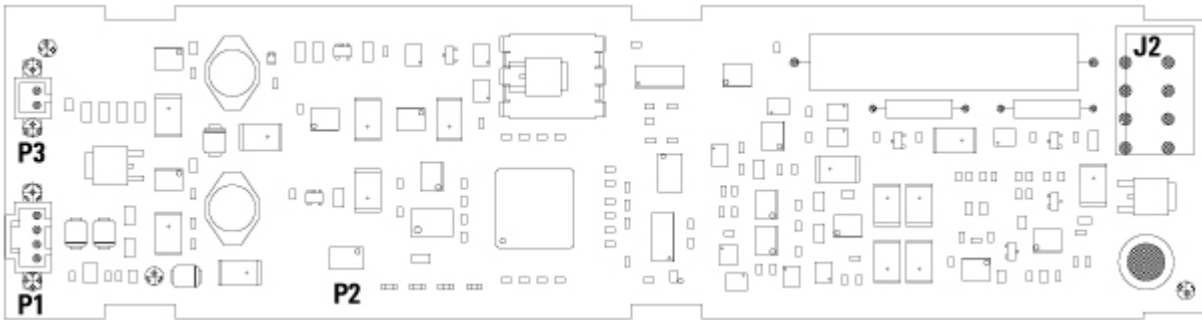

Connector	Function	See also
J2	Filament (upper 2) and ?PRT (lower 2) connections.	
P1	Communications buss.	“Communication harness”
P2	MTC	
P3	Switching valve.	“Replacing the TCD switching valve”

uECD signal board


Part number G3433-60020 shown below. (Replaced by G3433-61020.)

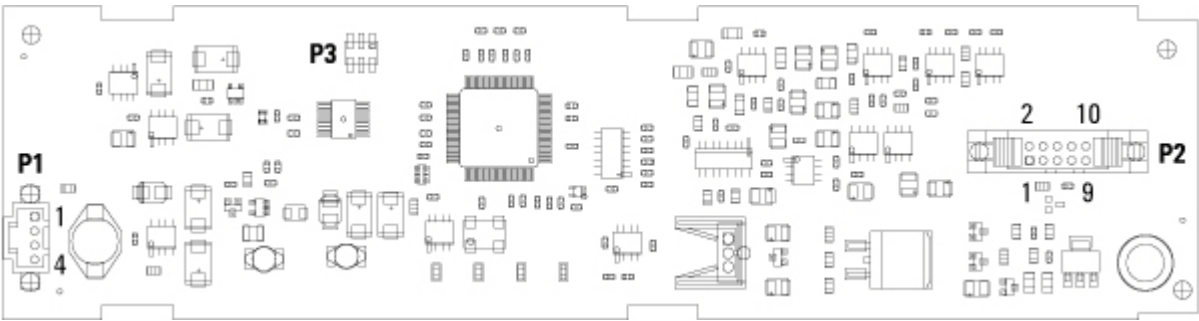

Connector	Function	See also
P1	Communications buss.	<a href="#">“Communication harness”</a>
P2	Signal from electrometer.	
P3	MTC	

FPD signal board


Part number G3435-60010

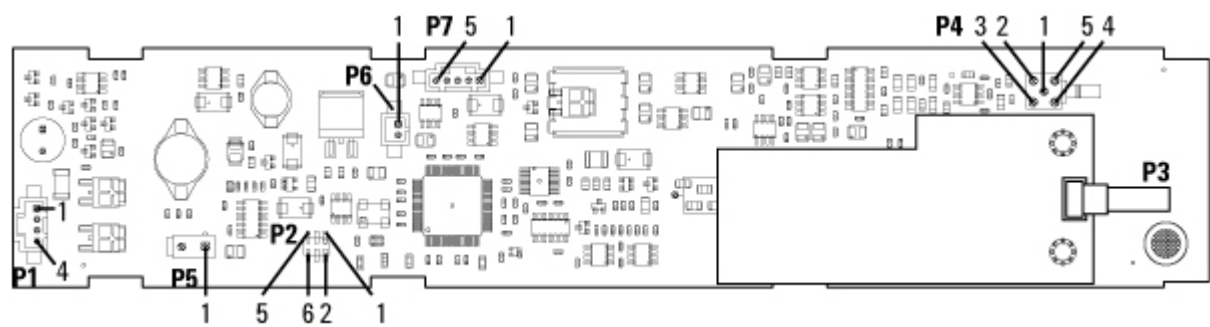

Connector	Function	See also
P1	Communications buss.	<a href="#">“Communication harness”</a>
P2	MTC	
P3	High voltage to photomultiplier.	
P4	Signal from electrometer.	
P5	Ignitor.	
P^	MTC	
P7	MTC	

NPD signal board


Part number G3434-60021 (replaces G3434-60020)

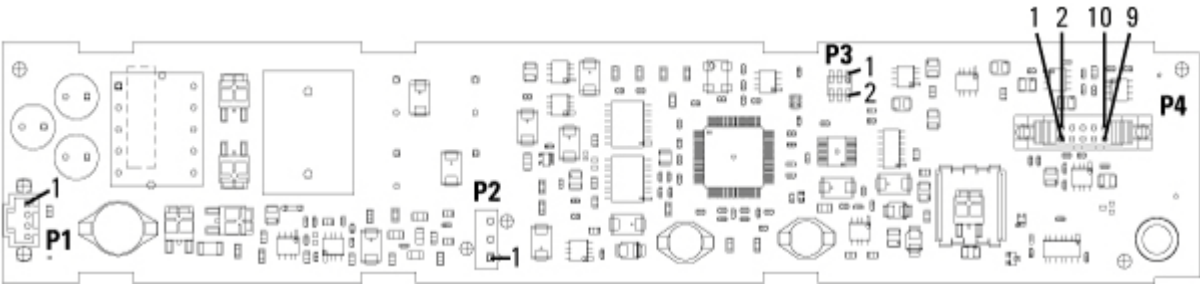

Connector	Function	See also
P1	Communications buss.	<a href="#">“Communication harness”</a>
P2	Bead power.	<a href="#">“NPD power cable”</a>
P3	MTC	
P4	Signal from electrometer.	

FID electrometer board


Part number G3431-60010

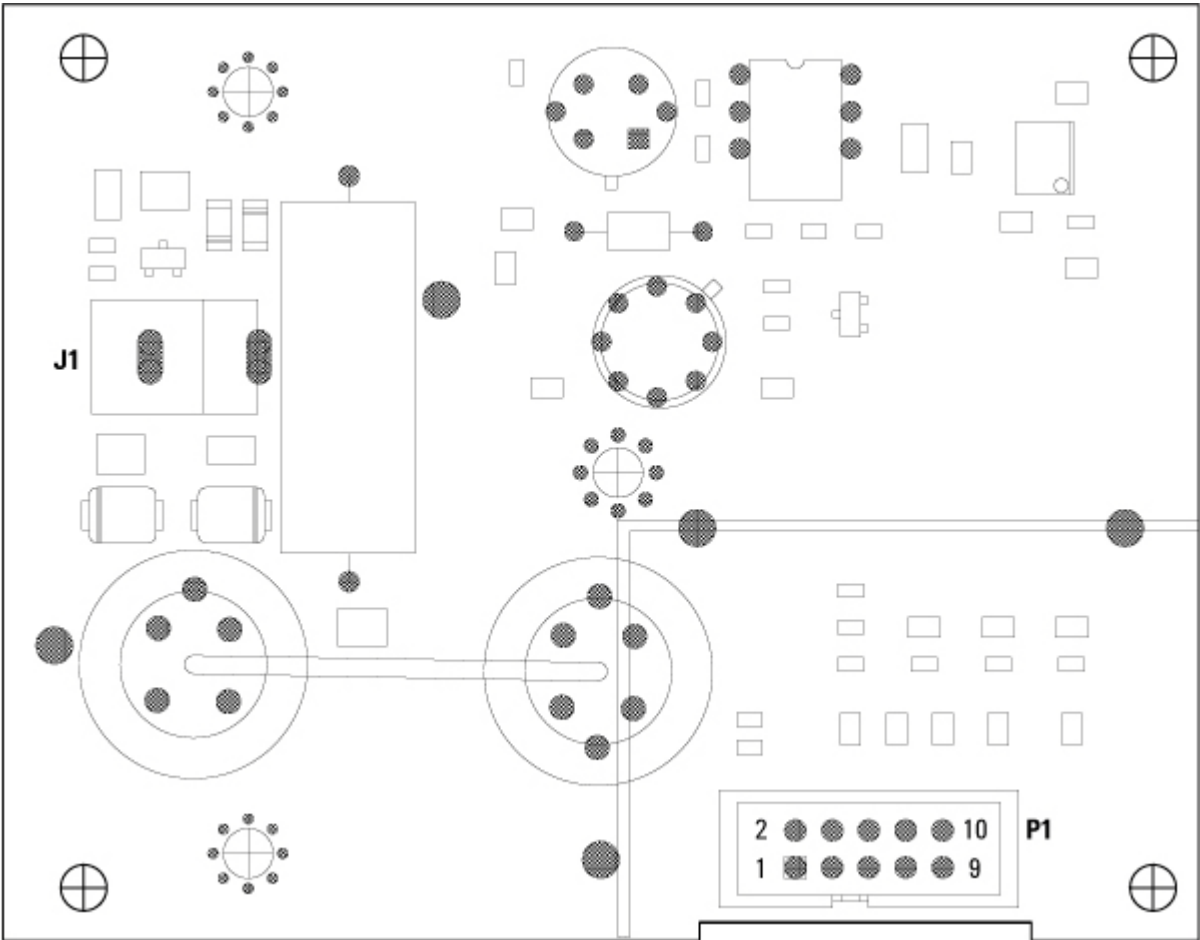

Connector	Function
J1	Input from detector.
P1	Output to signal board.

## uECD electrometer board

Part number G3433-60010

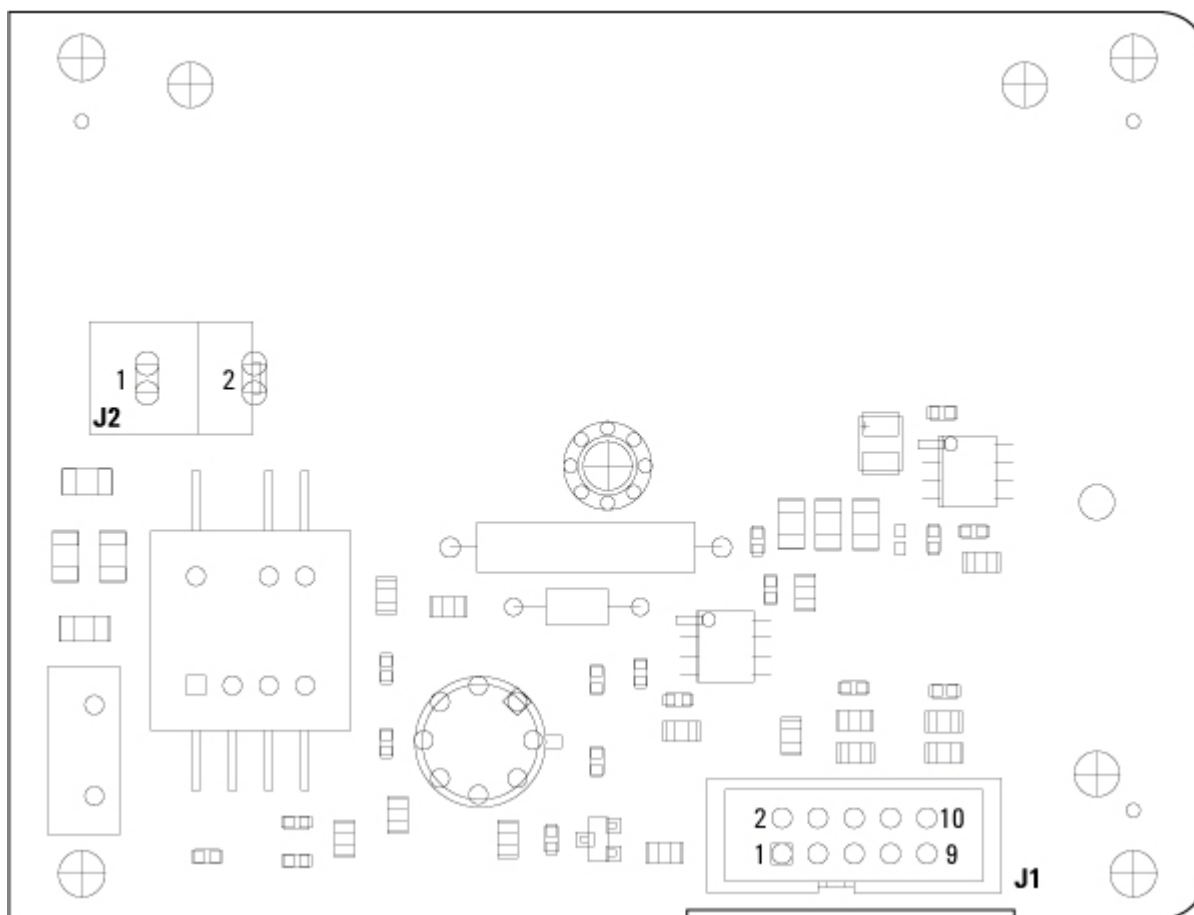


Connector	Function
J1	Output to signal board.
J2	Input from detector.



NPD electrometer board


Part number G3434-60010

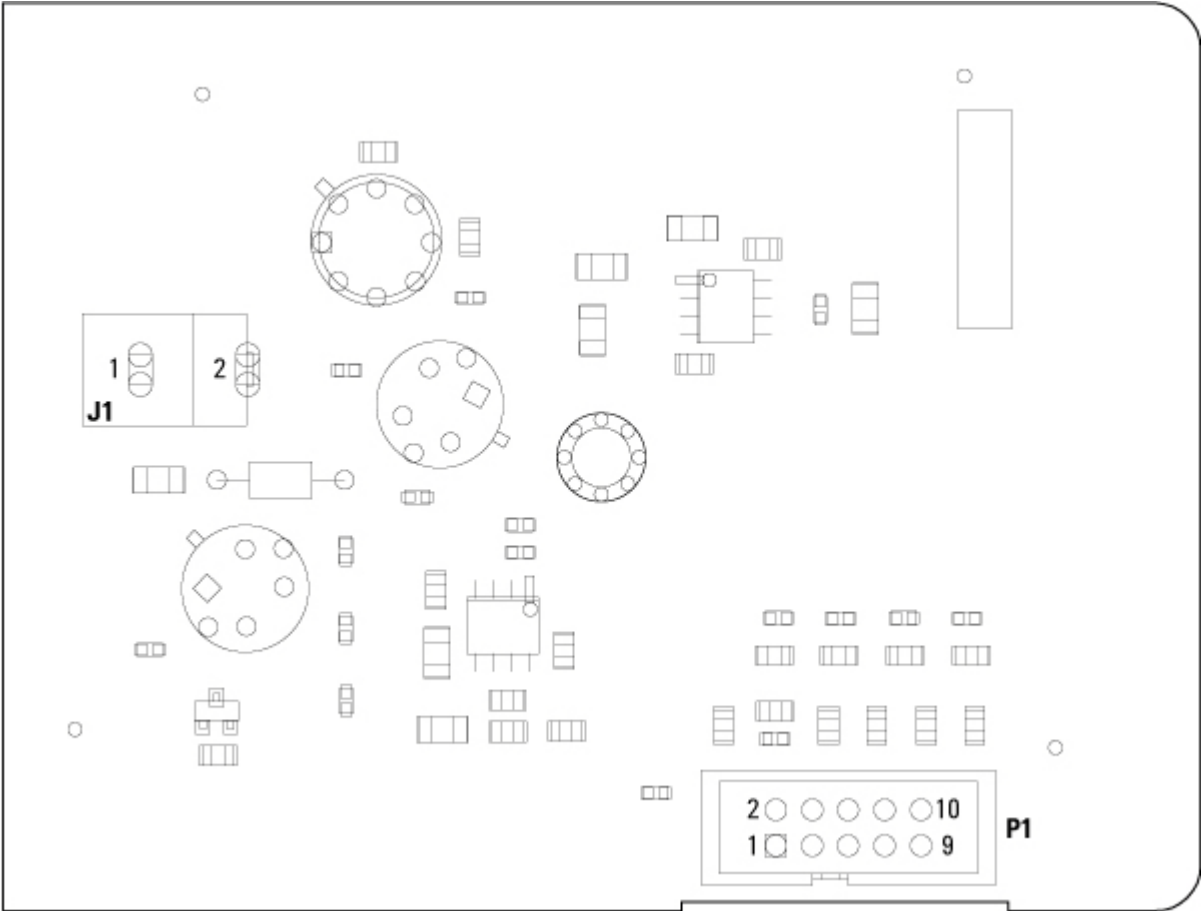

Connector	Function
J1	Input from detector.
P1	Output to signal board.

# Analog input board

Part number G3456-60010

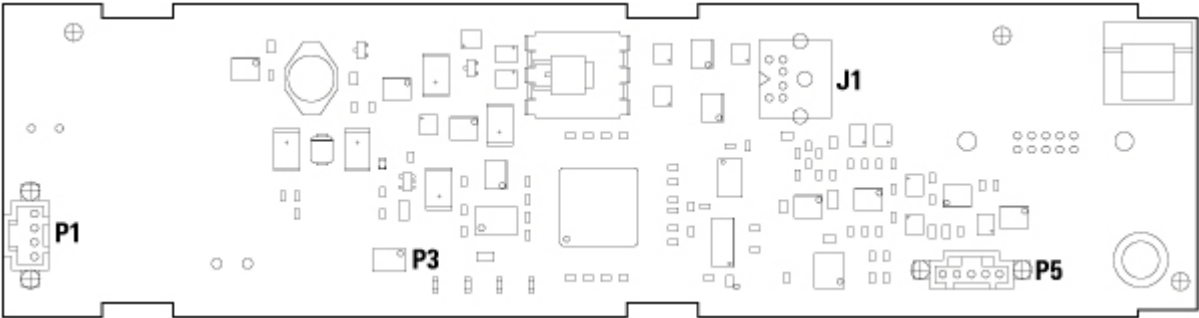

Connector	Function	See also
J1	External input cable.	<a href="#">“Analog signal outputs”</a>
P1	Communications buss.	<a href="#">“Communication harness”</a>
P3	MTC	
P5	MTC	



## 18 Illustrated Parts Breakdown

Inlets	584
Detectors	627
EPC Module Parts	669
Nickel Catalyst Accessory	674
Oven	676
Valves	685
Electronics and Fans	700
Covers	720
Tools	729

This section consists of parts lists and exploded views of the instrument modules. It is a major reference for part identification and part numbers.

For consumables and other customer-orderable parts, see the latest Agilent Parts Finder software version and IPB models for the 7890 Series GC.

In any table listing parts, a part with ID “NS” means the part is “not shown” in the accompanying graphics.

Quantities listed are the quantities required in the assembly, not necessarily the quantity to order (the quantity to order is not the same as the quantity used when the parts are sold in packages of more than one).



## **Inlets**

This section contains illustrated parts breakdowns for GC inlets and related components.


- Split/Splitless Inlet (SSL)
- Multimode Inlet (MMI)
- Purged Packed Inlet (PP)
- Cool on-column inlet (COC)
- Programmed Temperature Vaporization Inlet (PTV)
- Volatiles Interface (VI)

Split/Splitless Inlet (SSL)

Top level subassemblies for SSL inlet accessories:


- G3452A 100 psi split/splitless with EPC, Kit
- G3460A 150 psi split/splitless inlet with EPC, Kit
- G3453B Inert flow path, 0–100 psi

Description	Part number
SSL Inlet Accessory 0-100 PSIG without ship kit (includes inlet weldment, EPC module, and split vent filter assy)	G3452-60500
SSL Inlet Accessory 0-150 PSIG without ship kit (includes inlet weldment, EPC module, and split vent filter assy)	G3452-60509
SSL Weldment Assembly without EPC module (includes inlet weldment and split vent filter assembly)	G3452-60512
SSL inlet body with heater/sensor and insulation (includes inlet weldment only)	G3452-67000
SSL EPC module, 100 PSI	G3452-60552
SSL EPC module, 150 PSI	G3452-60510

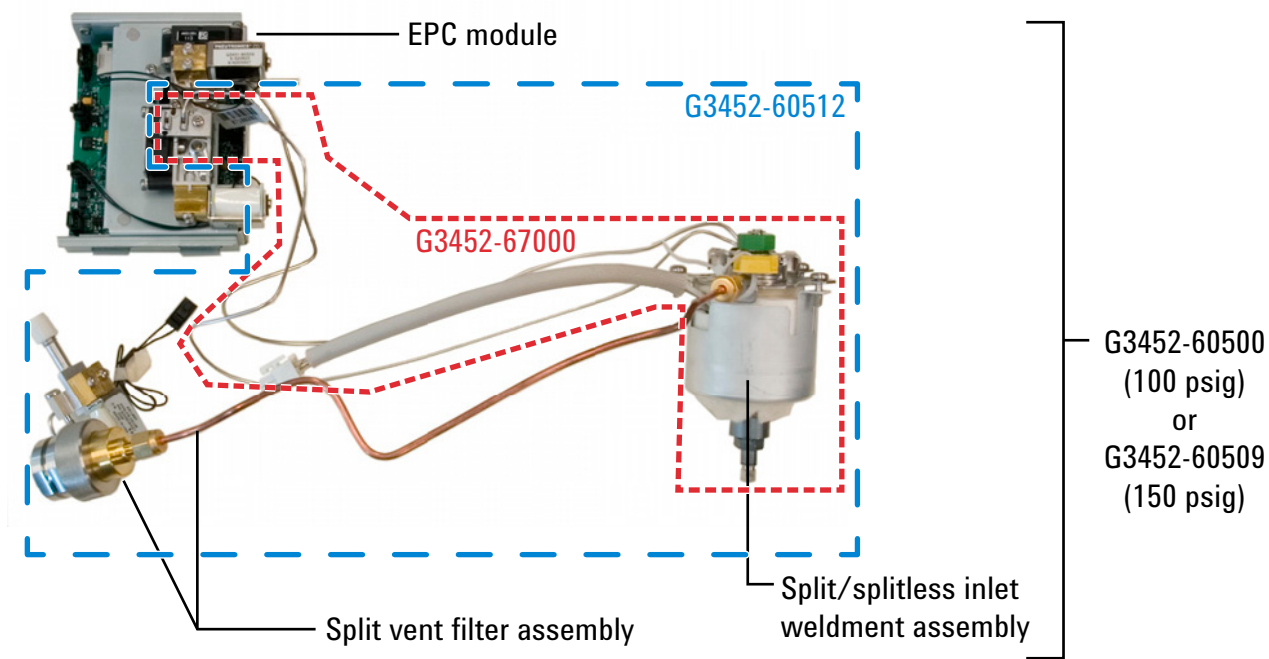


**Split/Splitless Inlet, Top**

Item	Description	Part number	Qty
1	Cap, Merlin Microseal, low pressure	See Parts Finder	1
	Cap, Merlin Microseal, general purpose	See Parts Finder	1
2	Merlin Microseal	See Parts Finder	
3	Top insert assembly, standard	G3452-60730	1
NS	Insert weldment, inert, 7890	G3452-60586	
NS	Top SSL insert weldment assembly, headspace	G3452-60100	
NS	Top insert, AC gang fitting weldment	G3430-60011	
NS	Top insert assembly, valve	G3480-67585	
4	Septum nut, standard	18740-60835	1
5	Septum, standard	See Parts Finder	1
6	O-ring, 12/pk (2 required for EPC module seal, 2 required for valve seal)	5180-4181	4
7	Screw, captive, M3	G1946-20168	2

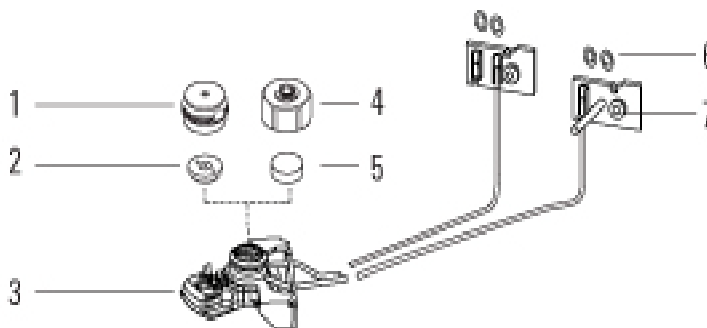


Item	Description	Part number	Qty
1	Tubing assembly	G3430-60011	1
2	Split ring	0510-1306	1
3	Turn-top	G3430-40035	1

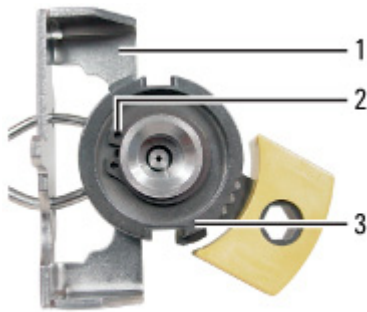

Bottom view

## SSL Split Vent Trap


Item	Description	Part number	Qty
1	Vent line tube assembly	G3430-60526	1
2	Screw, captive, M4 x 20 mm, T-20	1390-1024	2
3	Trap shell, back assembly (does not include item 9 or 10)	G3452-60506	1
4	Trap nut	G1544-20650	1
NS	Split vent trap	See Parts Finder	1
5	Threaded trap shell, front	G1544-20675	1
6	Proportional valve, split vent	G3430-60527	1
7	Screw, M3 x 16 mm T-10 to attach valve	0515-1141	2
NS	O-ring, 12/pk (2 required for EPC module seal, 2 required for valve seal)	5180-4181	4
8	1/8-inch copper tubing with length, 50 ft. Trim to 325 mm or 12.8 in.	5180-4196	
9	Screw, M4 x 12 mm, T-20	0515-2496	1
10	Trap retainer (key)	G3452-00005	1

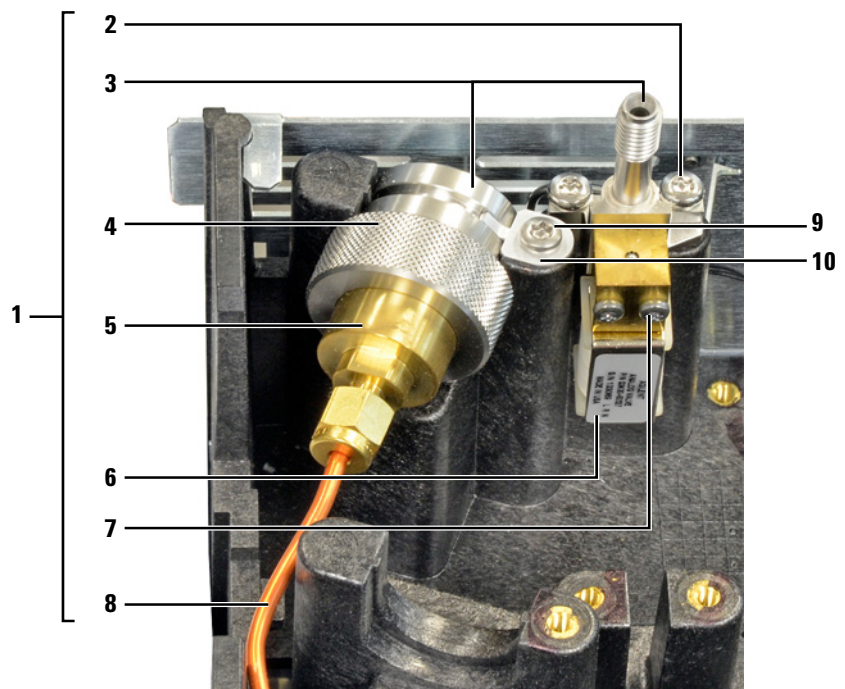




## Split Splitless Inlet Body

Item	Description	Part number	Qty
1	Screw, M4 x 12 mm, Torx T-20 (captive) C-clip	1390-1023	3
		0510-1799	3
2	Inlet weldment	G3452-80570	1
	Cap Inlet Shell Weldment Assembly, Inert	G3452-60570	
3	Top insulation	G1544-00010	1
4	Heater/Sensor assembly	G1530-67950	1
5	Screw, M4 x 10 mm	0515-4788	1
6	Washer, 12/pk	5061-5869	1
7	Heater block	G3452-20500	1
8	Insulation block	G3452-00002	1
9	Insulation bottom	G3432-00003	1
10	SSL insulation cover	G3431-00005	1
11	Retaining nut for heater block	G1544-20590	1
12	Gold seal kit, includes washer	5188-5367	1
	Gold seal, Ultra Inert, with washer	5190-6144	1
	Gold plated seal with cross, split only	5182-9652	1
13	Reducing nut	18740-20800	1
14	Ferrule	See Parts Finder	
15	Column nut	See Parts Finder	

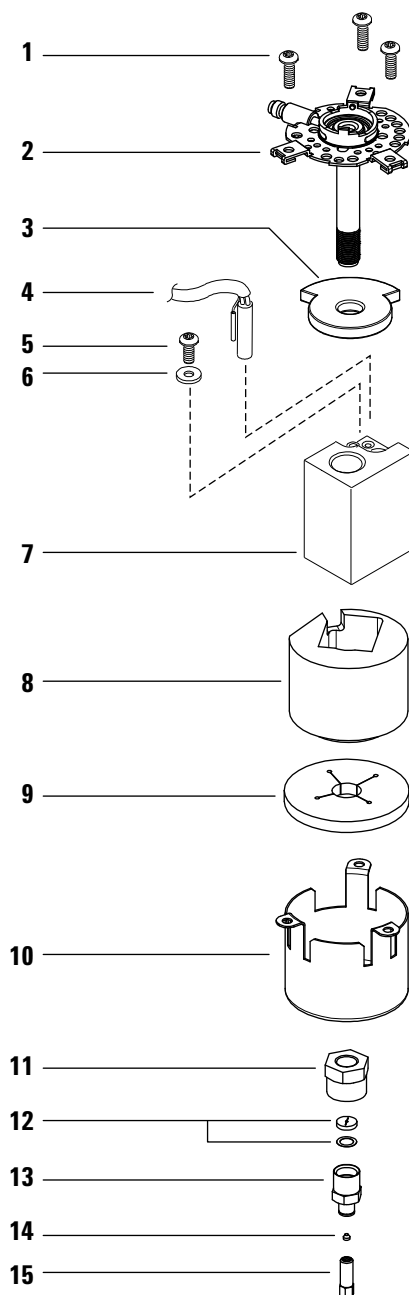

Split Splitless Inlet Warmer

Item	Description	Part number	Qty
1	SSL insulation kit, 3 pieces	5188-5241	1
2	Bottom insulation cover	19243-00070	1
3	Screw, M4 x 12 mm, Torx T-20, chrome-plated	0515-2711	2

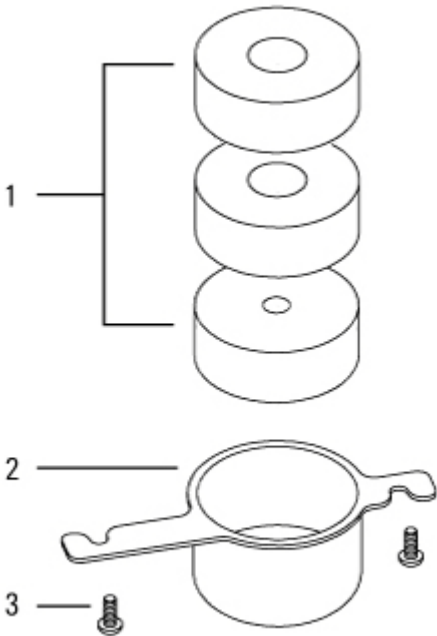

## Multimode Inlet (MMI)

Top level subassemblies for the 7890B:


- G3510B Multimode inlet, liquid carbon dioxide (CO<sub>2</sub>), kit
- G3511B Multimode inlet, liquid nitrogen (N<sub>2</sub>), kit



Top level subassemblies for the 7890A:





- G3510A Multimode inlet, liquid carbon dioxide (CO<sub>2</sub>), kit
- G3511A Multimode inlet, liquid nitrogen (N<sub>2</sub>), kit

See the figures below.

Item	Description	Part number
<b>7890A and 7890B</b>		
1	Inlet chassis	G3430-40011
2	MMI EPC module	G3510-60500
<b>7890A Only</b>		
3	Multimode inlet, liquid carbon dioxide (CO <sub>2</sub> ), kit	G3510-64000
4	Multimode inlet, liquid nitrogen (N <sub>2</sub> ), kit	G3511-64000
<b>7890B only</b>		
5	MMI LCO <sub>2</sub> Module 7890B (includes EPC module, inlet body, cables, cryo components, and thermocouple PCA)	G3510-65200
6	MMI LN <sub>2</sub> Module 7890B (includes EPC module, inlet body, cables, cryo components, and thermocouple PCA)	G3510-65300

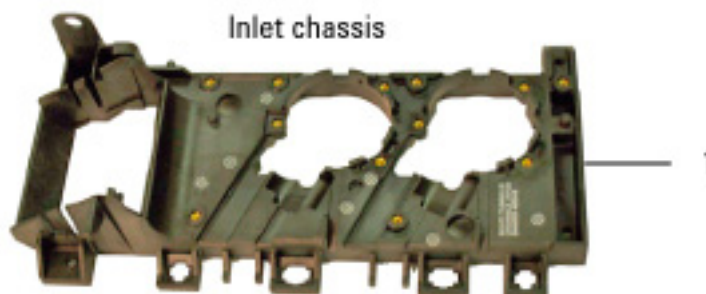

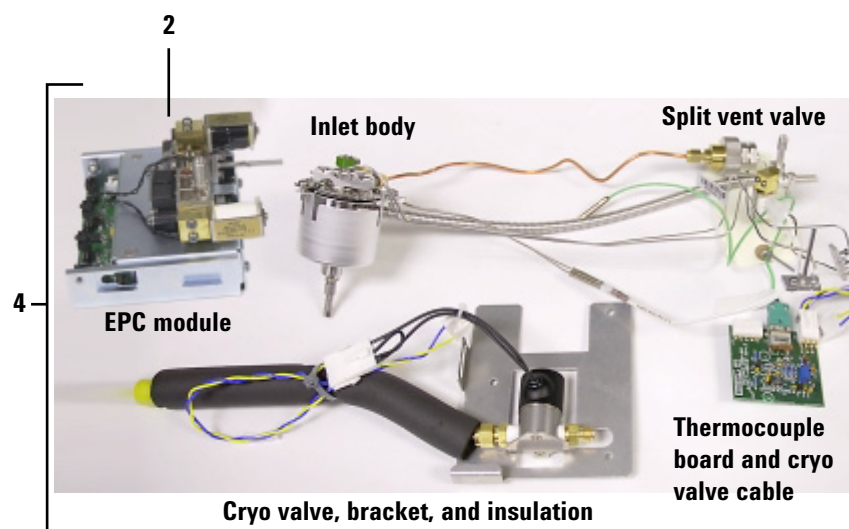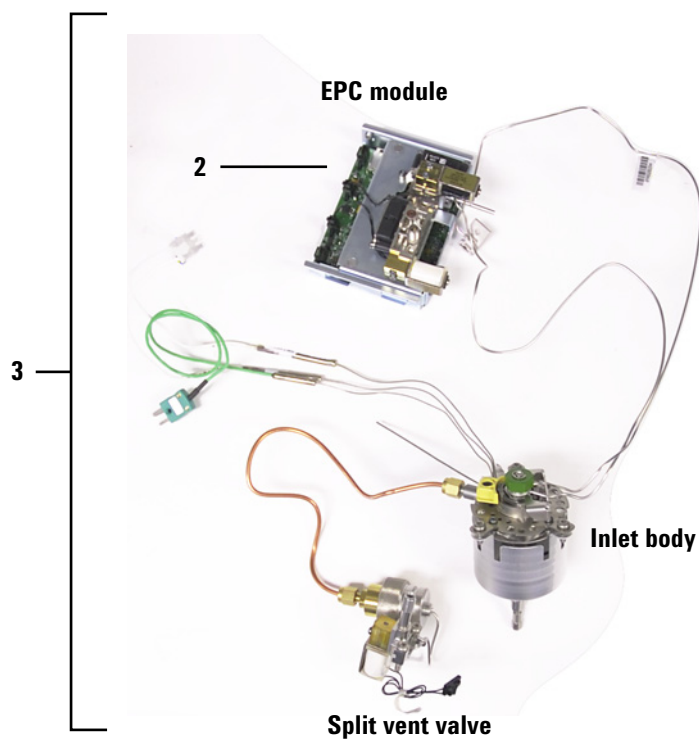


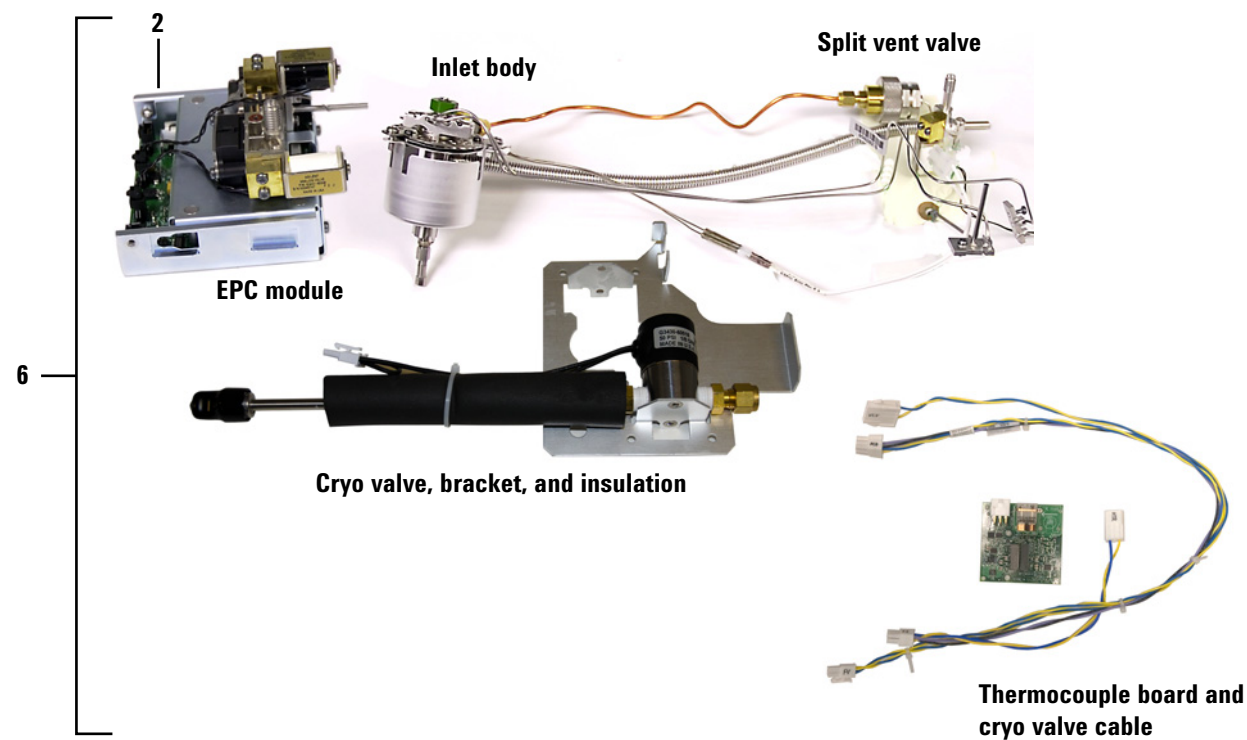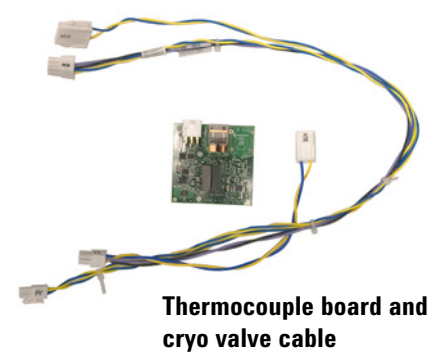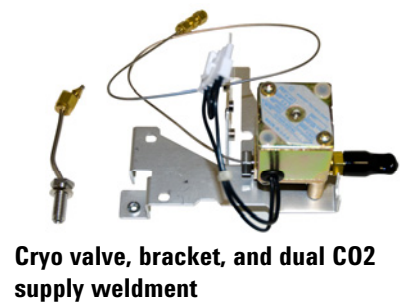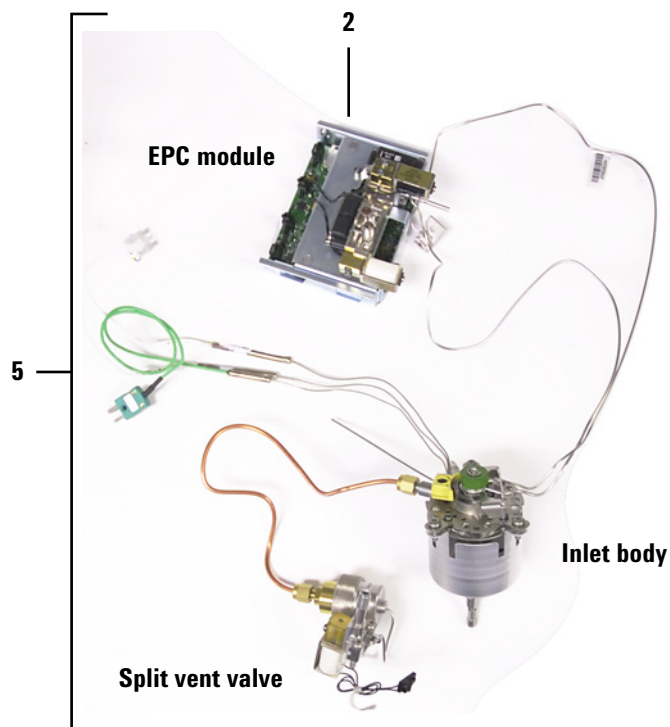

### Multimode Inlet Body

Item	Description	Part number	Qty
NS	MMI EPC manifold	G3510-60500	1
1	MMI Full Weldment* (includes cryo jacket)	G3510-64100	1
	<b>MMI cryo jackets</b>		1
2	MMI LN <sub>2</sub> jacket, <b>7890A only</b>	G3510-81310	
2	MMI LN <sub>2</sub> jacket, <b>7890B only</b>	G3510-81011	
3	MMI CO <sub>2</sub> jacket, 7890A and 7890B	G3510-81210	
4	Column nut adapter	G3510-20018	1
5	Foil (LN <sub>2</sub> only)		
6	MMI Insulation (LN <sub>2</sub> or LCO <sub>2</sub> )	19243-00069	1
7	MMI Metal Insulation Cup (LN <sub>2</sub> )	G3510-00001	1
	MMI Metal Insulation Cup (LCO <sub>2</sub> )	G3510-00002	
8	Tube, precise, 0.010-in. id x 1/16-in. od x 15 cm long	G2855-20700	1
NS	MMI CO <sub>2</sub> kit, <b>7890A only</b>	G3510-67520	1
NS	MMI LN <sub>2</sub> kit, <b>7890A only</b>	G3510-67530	1
NS	Brazed shell	G3510-81000	1
* UltraClean part			


See the Agilent Parts Finder software and the Agilent chemical analysis consumables and accessories catalog for Merlin assemblies.

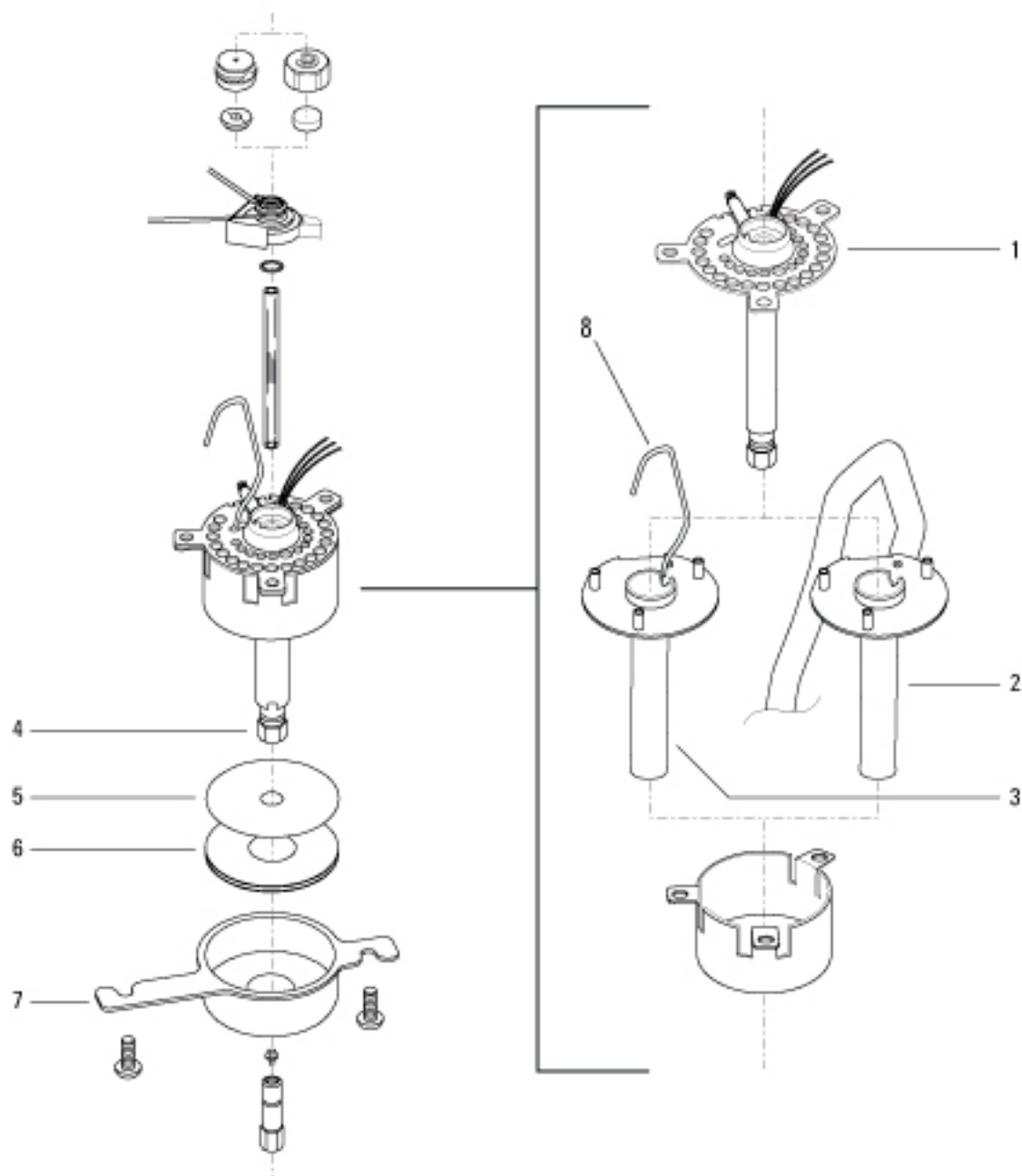




### MMI Top Insert Assembly

Item	Description	Part number	Qty
1	Cap, Merlin Microseal, low pressure	See Parts Finder	1
	Cap, Merlin Microseal, general purpose	See Parts Finder	1
2	Merlin Microseal	See Parts Finder	
3	Top insert assembly, standard	G3452-60730	1
NS	Top SSL insert weldment assembly, headspace	G3452-60100	
NS	Top insert, AC gang fitting weldment	G3430-60011	
NS	Top insert assembly, valve	G3480-67585	
4	Septum nut, standard	18740-60835	1
5	Septum, standard	See Parts Finder	1
6	O-ring, 12/pk (2 required for EPC module seal, 2 required for valve seal)	5180-4181	4
7	Screw, captive, M3	G1946-20168	2

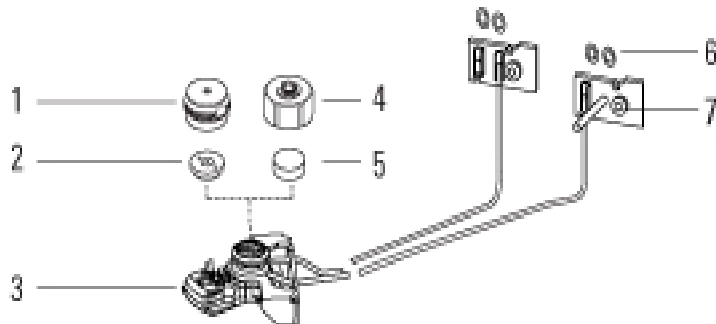


Item	Description	Part number	Qty
1	Tubing assembly (top insert AC gang fitting weldment)	G3430-60011	1
2	Retainer ring	0510-1306	1
3	Turn-top (helix plate with molded handle)	G3430-40035	1

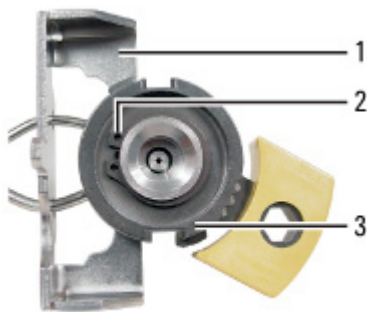

Bottom view


## MMI Split Vent Trap

Item	Description	Part number	Qty
1	Vent line tube assembly	G3430-60526	1
2	Screw, captive, M4 x 20 mm, T-20	1390-1024	2
3	Trap shell, back assembly (does not include item 9 or 10)	G3452-60506	1
4	Trap nut	G1544-20650	1
NS	Split vent trap	See Parts Finder	1
5	Threaded trap shell, front	G1544-20675	1
6	Proportional valve, split vent	G3430-60527	1
7	Screw, M3 x 16 mm T-10 to attach valve	0515-1141	2
NS	O-ring, 12/pk (2 required for EPC module seal, 2 required for valve seal)	5180-4181	4
8	1/8-inch copper tubing with length, 50 ft. Trim to 325 mm or 12.8 in.	5180-4196	
9	Screw, M4 x 12 mm, T-20	0515-2496	1
10	Trap retainer (key)	G3452-00005	1

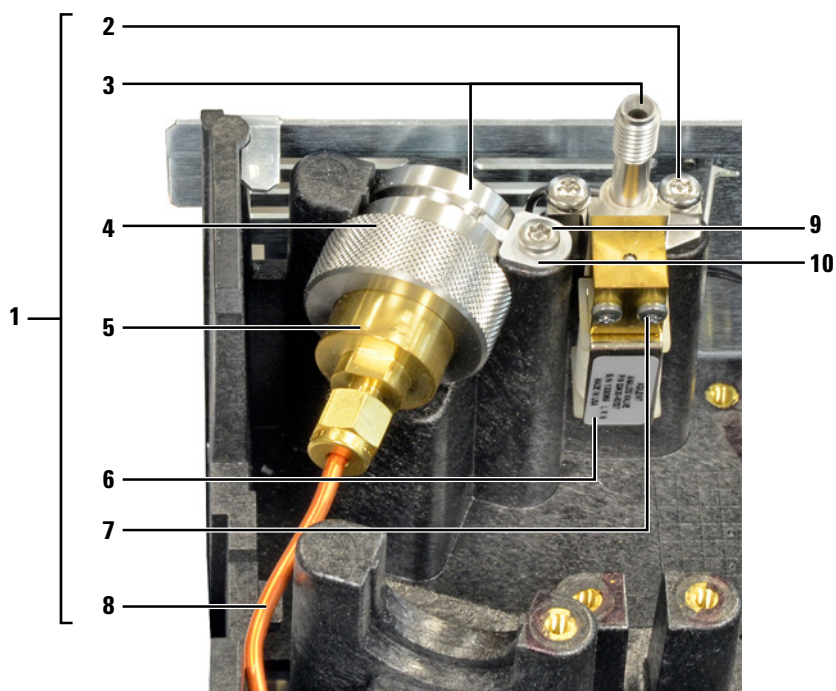


**7890A MMI Temperature controller**

Item	Description	Part number	Qty
1	Wiring harness cable	G3510-60110	1
2	Thermocouple board	G3430-60095	1
3	Screw, M4 x 6 mm, T-20 (attach board to bracket)	0515-2832	2
4	Screw, M4 x 12 mm, T-20 (attach bracket to carrier)	0515-2496	2
5	Bracket, PTV thermocouple board and LC02 valve (single CO2 bracket)	G3430-00058	1
6	Kit, containing thermocouple board and wiring harness cable	G3510-67000	-

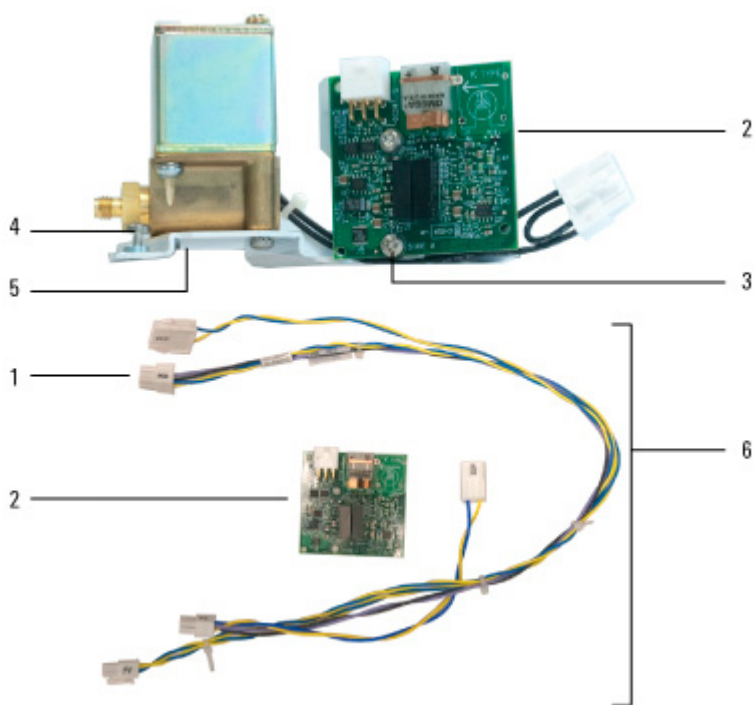**Figure 15** 7890A Thermocouple board, bracket, and cables


**7890B MMI Temperature controller**

Item	Description	Part number	Qty
1	Thermocouple board	G3430-60095	1
	Bracket, cryo		1
2	Bracket, PTV thermocouple board, LC02 cryo	G3430-00092	
4	Bracket, PTV thermocouple board, LN2 cryo (Single CO2 bracket)	G3430-00058	
3	Screw, M4 x 12 mm, T-20 (attach board to bracket)	0515-2496	2
5	Screw, M4 x 12 mm, T-20 (attach bracket to GC carrier)	1390-1023	2
6	Wiring harness cable	G3510-60110	1
7	Kit, containing thermocouple board and wiring harness cable	G3510-67000	—

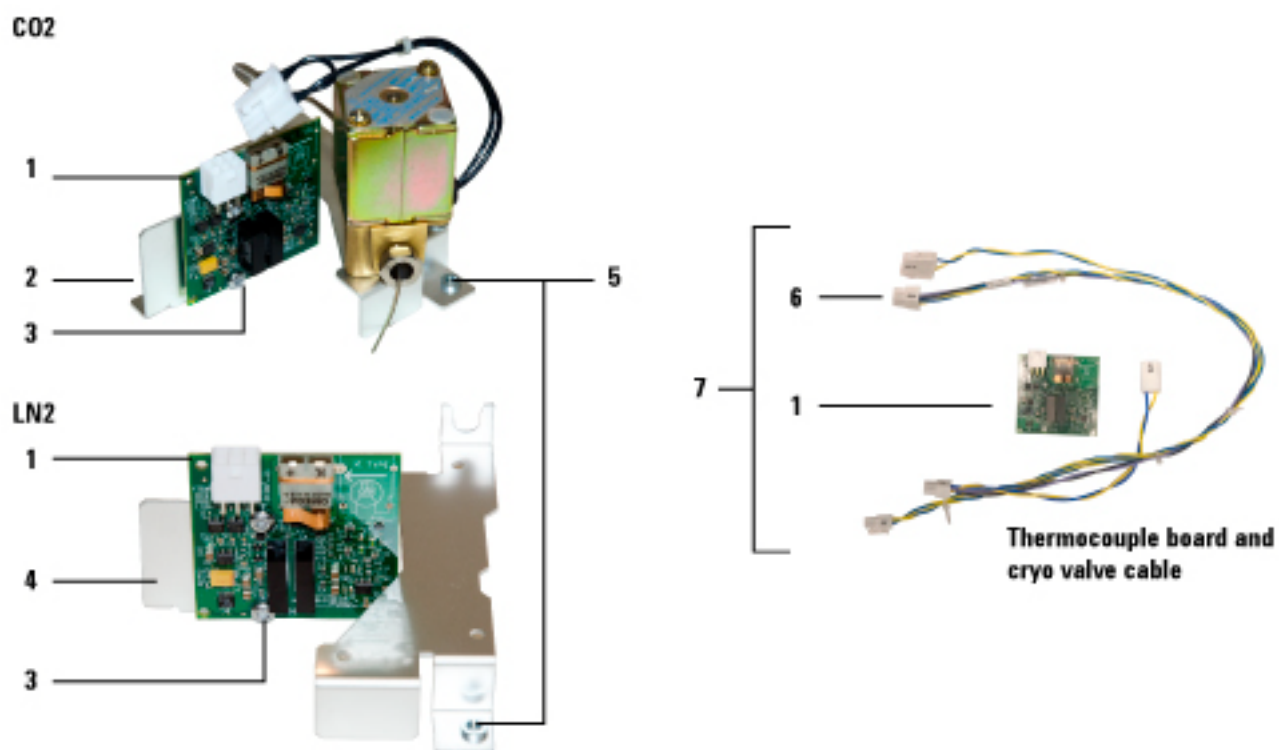**Figure 16** 7890B Thermocouple board, bracket, and cable

**7890A MMI Cryo assemblies**

Item	Description	Part number	Qty
NS	7890 MMI LC02 kit (includes CO2 cryo valve assembly, cryo bracket, insulation cup)	G3510-67520	1
NS	7890 MMI LN2 kit (includes insulation cover, cryo bracket, insulation cup, and LN2 cryo valve assemblies)	G3510-67530	1
1	CO2 brazement (valve to IP tube)	G2617-80510	1
NS	Cap, protective without flange	1401-0044	1
NS	1/8-inch male connector without O-ring	G1543-80025	1
NS	Gland seal CO <sub>2</sub>	G1565-20590	1
NS	PTFE seal	G1565-20840	1
2	LN2 Noise reduction Restrictor	G3510-27000	

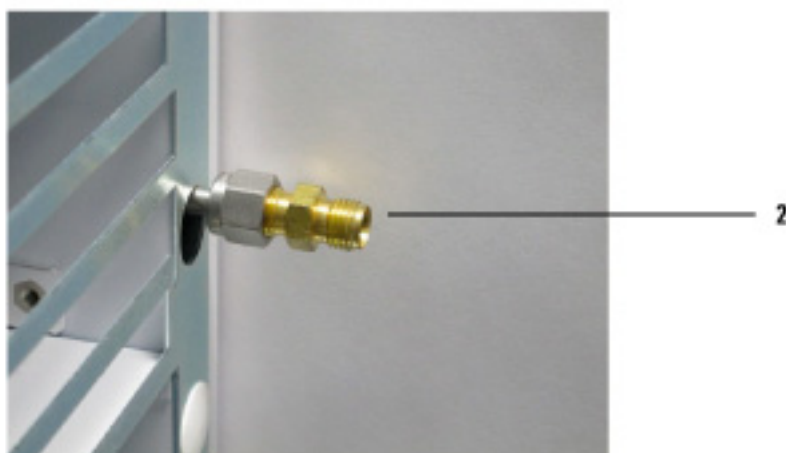**CO2****LN2**

### 7890B MMI Cryo assemblies


Item	Description	Part number	Qty
	MMI LCO2 Module (includes MMI Full Weldment - CO2, EPC module, grommet, TC board, and cable)	G3510-65200	
	MMI Full Weldment - CO2 (includes CO2 cryo assembly, MMI body and split vent valve, insulation and insulation cup)	G3510-60520	
	MMI LN2 Module 7890B (includes EPC module, grommet, cryo valve bracket, TC board, cable, and MMI Full Weldment LN2 assembly)	G3510-67300	
	MMI Full Weldment LN2 7890B (includes LN2 cryo assembly, MMI body and split vent valve, insulation and insulation cup)	G3510-60535	

## Purged Packed Inlet (PP)

Top level subassemblies for Purged Packed inlet accessories:

- G3451A Purged-packed with EPC, Kit


Description	Part number
Purged-packed Inlet Accessory without ship kit	G3451-65520
Purged packed inlet weldment assembly without EPC module	G3451-80506
Purged packed inlet EPC module	G3451-60551

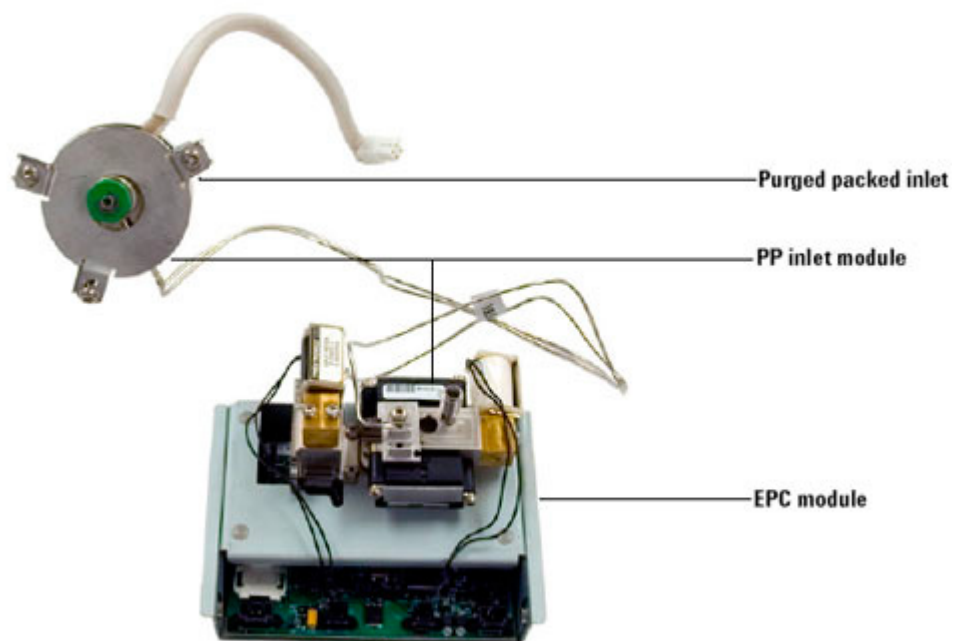




## Purged Packed Inlet Upper Body

Item	Description	Part number	Qty
1	Screw, Torx T-20, M4 x 12 mm	1390-1023	3
2	Top cover plate	G1543-00085	1
3	Top insulation	G1543-00100	1
4	Nut, Merlin Microseal	See Parts Finder	
5	Septum, Merlin Microseal	See Parts Finder	
6	Top insert weldment	19243-80570	1
7	Viton O-ring (12/pk)	5080-8898	1
8	Screw, Torx T-10, M3 x 16 mm	G1946-20168	2
9	Inlet weldment	G3451-80501	1
10	Thermal block	G1543-20765	1
11	Bottom nut	G1543-20580	1
12	Septum nut	18740-60835	1
13	Septum, 11 mm, low bleed (50/pk)	See Parts Finder	-
14	Heater/Sensor assembly	G1543-61540	1
NS	Top gasket insulation	G1543-00155	1

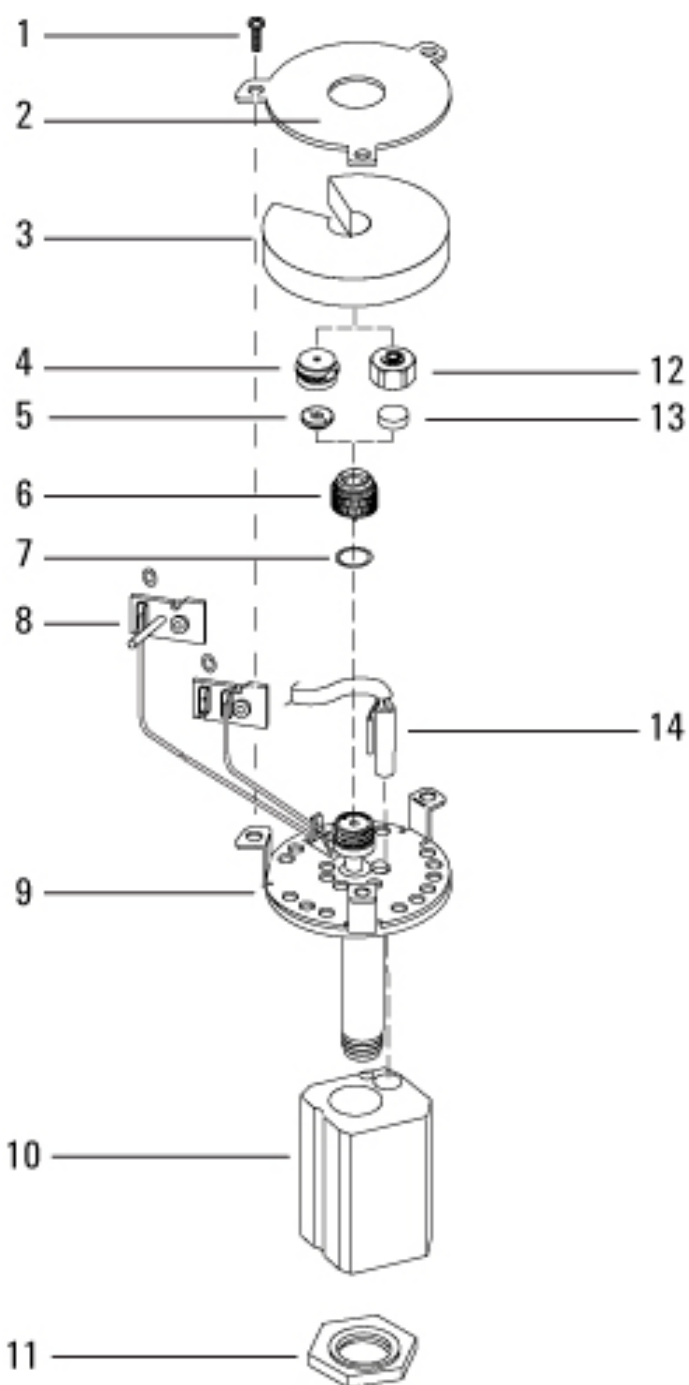

## Purged Packed Inlet Lower Parts

Item	Description	Part number	Qty
15	PP insulation	G3432-00004	1
16	PP bottom insulation	G3432-00003	1
17	PP insulation cup	G3431-00005	1
18	C-clip	0510-1799	3
19	Glass liner (25/pk)	See Parts Finder	1
20	Vespel ferrule (10/pk)	5080-8774	1
21	Tubing nut (1/4-inch brass), 10/pk	5180-4105	1
22	Column adapter	See Parts Finder	1
23	Bottom plate for oven	G1543-00060	1
24	Screw, T-20, M4 x 12mm	0515-2711	2
25	Adapter insulation, 1/pk	19234-60715	3
26	Insulating cup	19234-60700	1
27	Ferrules	See Parts Finder	
28	Column nut (capillary, universal, 2/pk)	See Parts Finder	1

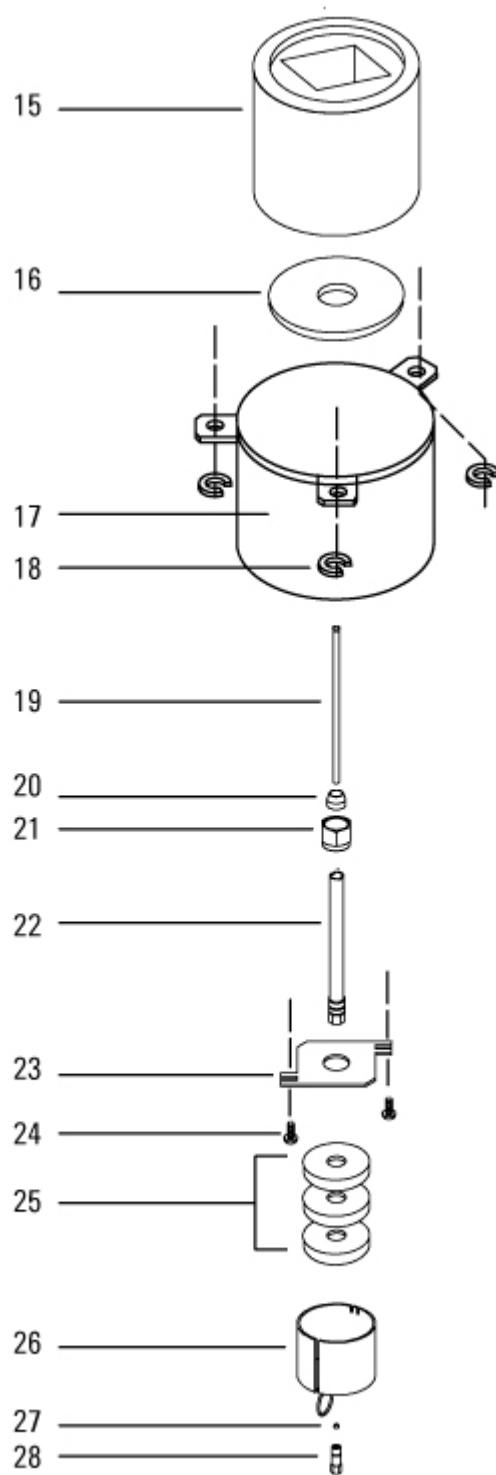


## Cool on-column inlet (COC)

Top level subassemblies for Cool on-column inlet accessories:


- G3454A Cool on-column inlet with EPC, Kit
- For cryoblast parts for G3467A and G3468A see “CO2 cryogenic cooling” on page 682 or “Liquid nitrogen (LN2) cryogenic cooling” on page 684.
- 19245-60590 for manual injection with a fused silica needle

Description	Part number
Cool On-column Inlet Accessory without ship kit	G3454-60501
Cool On-column Inlet weldment assembly without EPC module	G3454-60502
Cool On-column Inlet EPC module	G3454-60554

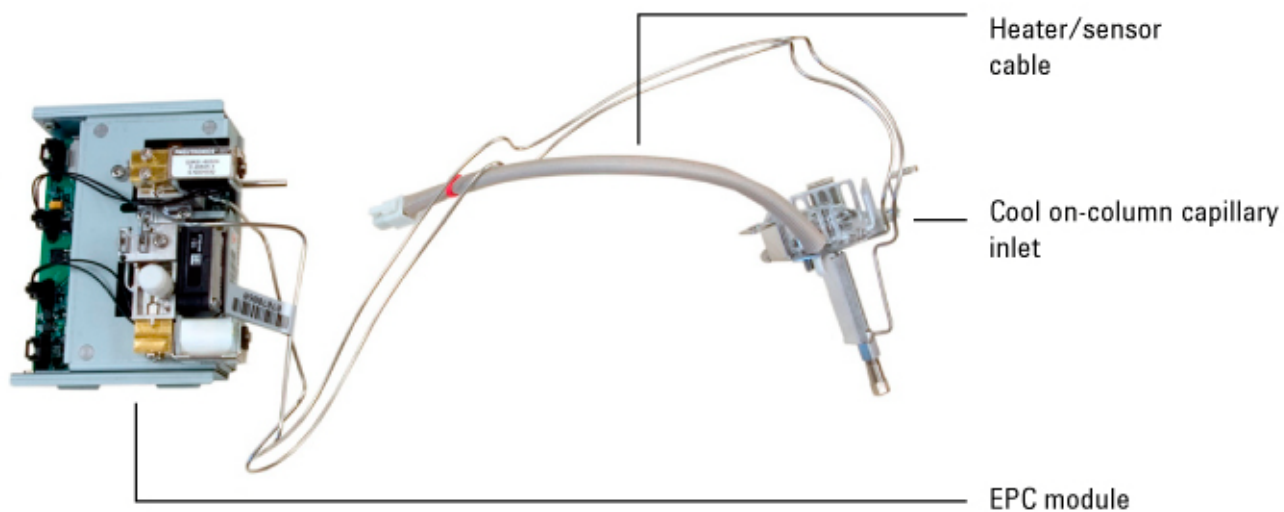

**Cool on-column inlet body**

<b>Item</b>	<b>Description</b>	<b>Part number</b>	<b>Qty</b>
2	Septum nut base assembly, for .250/.320-mm columns	19245-80521	1
3	Septum nut base, for .530 mm columns	G1545-80520	1
4	Manual cooling tower assembly (optional), for .200 mm columns	19320-80625	1
5	Duckbill for use with cooling tower, 10/pk	19245-40050	
6	Screw, M4 × 12 mm long	0515-2711	3
7	Heatsink fin	G1545-00010	1
8	Heater/Sensor assembly	G1545-60520	1
9	Screw, M3 × 16 mm, T-10	1390-1022	1
10	Septa	See Parts Finder	
11	Insert spring	19245-60760	1
12	Inserts (identify by number of rings on insert)	See Parts Finder	
13	Screw, M4 × 12 mm, T-20	0515-2496	3
14	Inlet weldment	G3454-80500	1
15	Inlet weldment insulation	G1545-20630	1
16	Cavity sleeve	19245-00060	1
17	Ferrules (identify by internal diameter)	See Parts Finder	
18	Column nut	See Parts Finder	

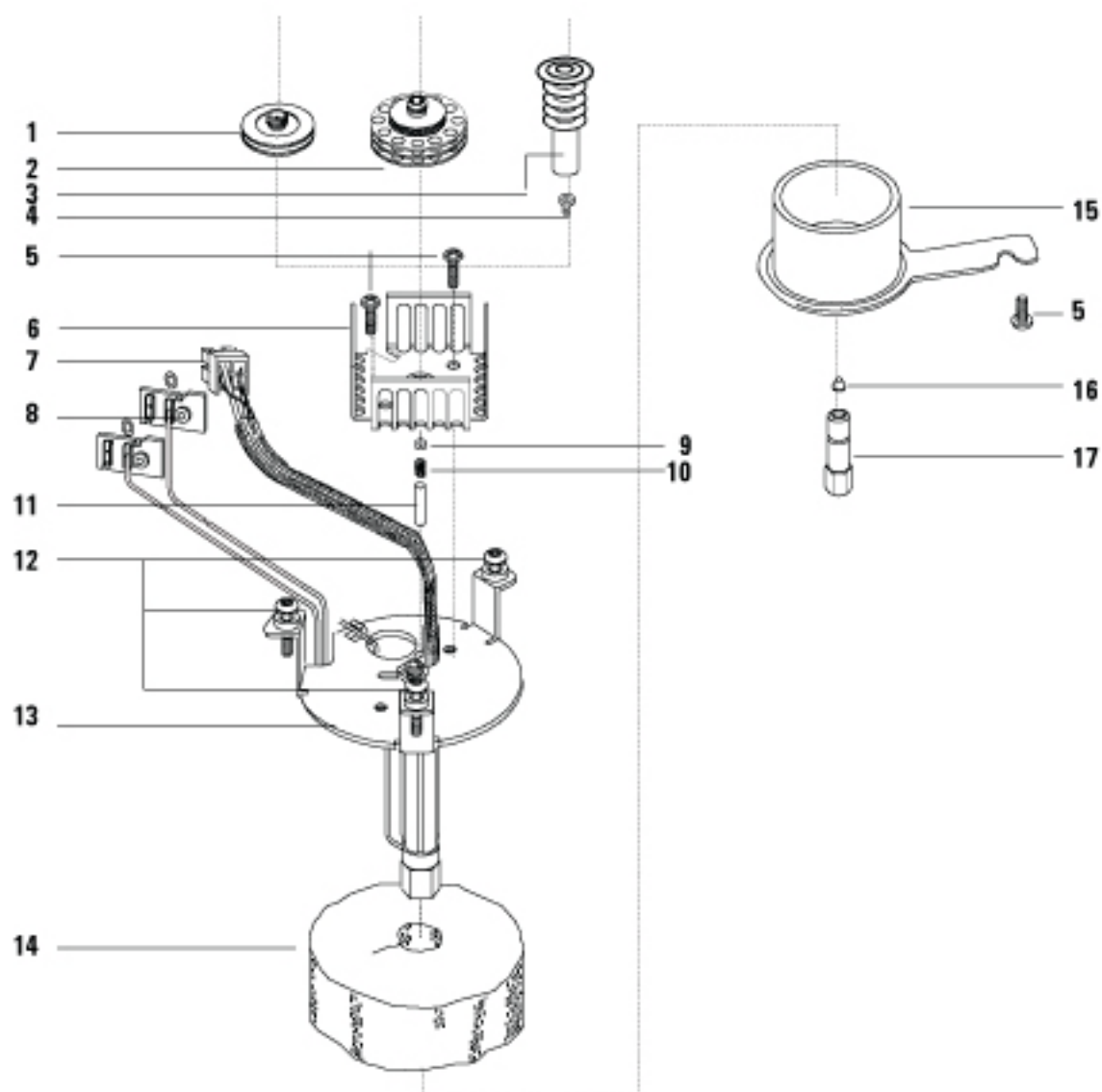

## Programmed Temperature Vaporization Inlet (PTV)


Top level subassemblies for PTV inlet accessories:

- G3500A Septumless head PTV, LCO2 Accessory with EPC, Kit
- G3501A Septum head PTV, LCO2 Accessory with EPC, Kit
- G3502A Septumless head PTV, LN2 Accessory with EPC, Kit, for 7890A only
- G3502B Septumless head PTV, LN2 Accessory with EPC, Kit, for 7890B only
- G3503A Septum head PTV, LN2 Accessory with EPC, Kit, for 7890A only
- G3503B Septum head PTV, LN2 Accessory with EPC, Kit, for 7890B only


Description	Part number
Septumless head PTV inlet, LCO2 accessory without ship kit	G3500-65500
Septum head PTV inlet, LCO2 accessory without ship kit	G3501-65500
Septumless head PTV inlet, LN2 accessory without ship kit	G3502-65500
Septum head PTV inlet, LN2 accessory without ship kit	G3503-60500
PTV EPC assembly, 100 PSI	G3500-60500
PTV Inlet chassis	G3430-40011

### Septumless head



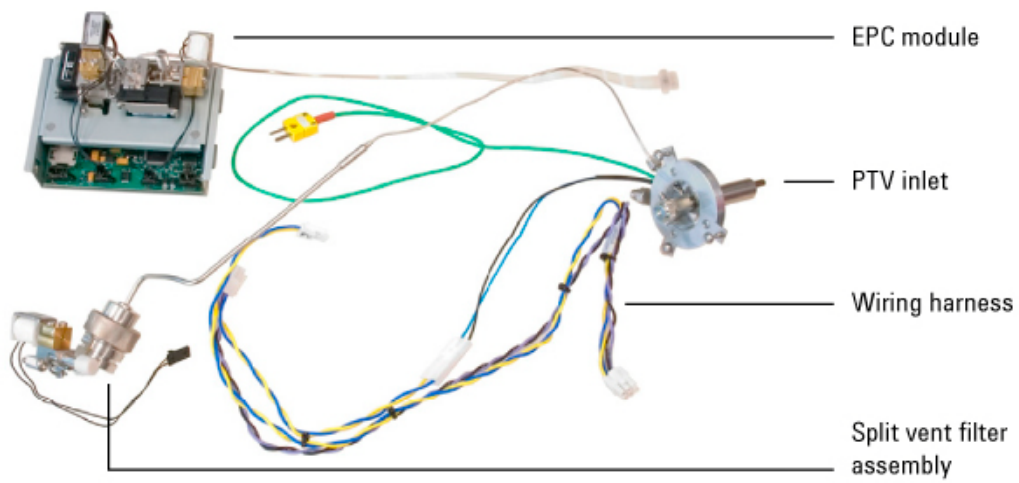


**Septum head**

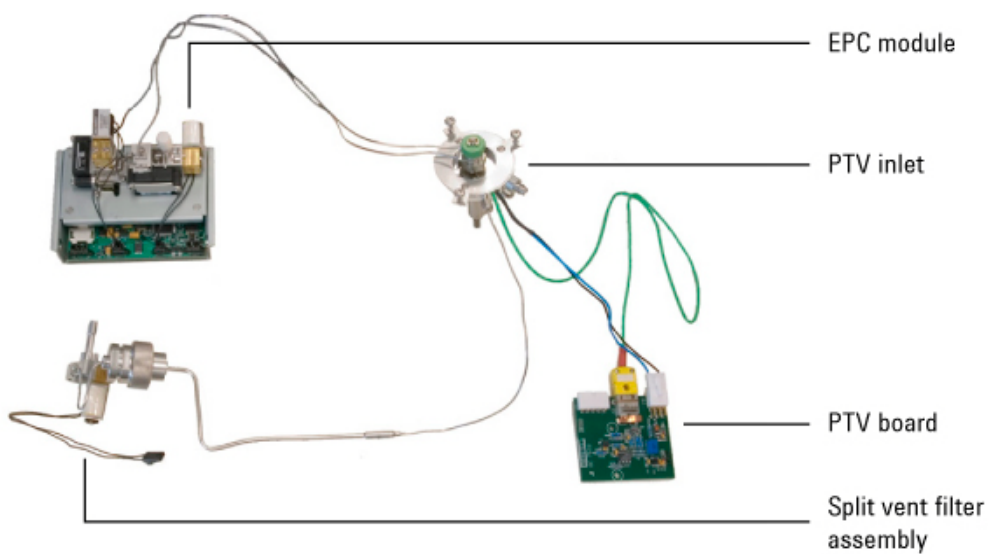

**PTV Inlet Body**


<b>Item</b>	<b>Description</b>	<b>Part number</b>	<b>Qty</b>
1	Septumless head assembly	G2617-60507	1
2	Septumless head weldment	G3500-80000	1
3	Septum nut (see also Parts Finder.)	18740-60835	1
4	Septum head (see also Parts Finder.)	G3500-80001	1
5	PTV inlet assembly	G2617-60506	1
6	Cooling jacket		
	PTV, CO2 cooling jacket	G2617-60508	1
	PTV N2 cooling jacket	G2619-60501	1
NS	PTV inlet assembly with cooling jacket		
	PTV inlet with CO2 cooling jacket	G2617-60518	
	PTV inlet with LN2 cooling jacket	G2617-60517	
7	Silver seal (5/pk)	5182-9763	1
8	Graphpak inlet adapter (includes silver seal and ferrule)	See Parts Finder	1
9	Ferrules for graphpak inlet	See Parts Finder	1
10	Split nut for inlet adapters	5062-3525	1
11	PTV liner	See Parts Finder	1
12	PTV insulation block	G2617-20510	1
NS	PTV cryo insulator	G2617-60510	1
NS	PTV ferrule assembly tool	G2617-80540	—

\* See the Parts Finder and the Agilent chemical analysis consumables and accessories catalog for Merlin assemblies and other consumables, such as the Graphpak adapter and ferrules.

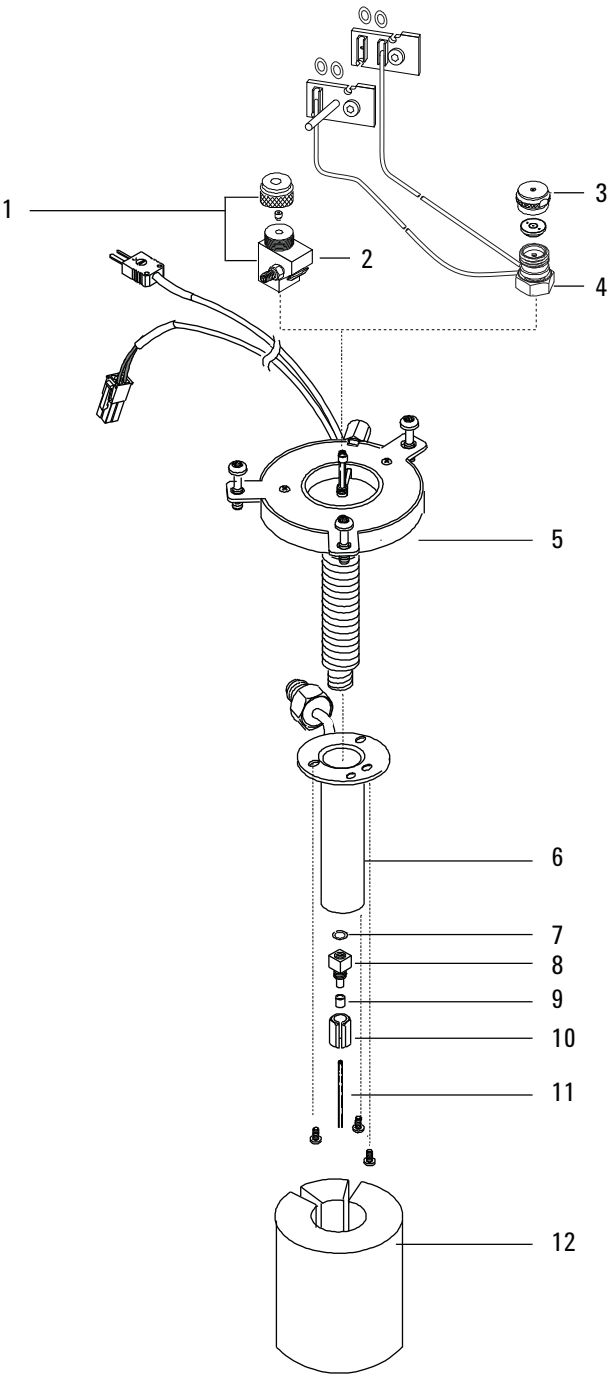


**PTV Septumless head**

	Description	Part number	Qty
	Septumless head	G2617-60507	
1	PTFE ferrule (needle seal)	5182-9748	1
2	Service kit for SLH (Septumless head repair kit: includes Kalrez seal, spring, and Viton seal)	5182-9747	1
(3)	Kalrez seal (normally order item 2 instead)	5182-9759	1
(4)	Pressure spring (normally order item 2 instead)	5182-9758	1
(5)	Viton seal (normally order item 2 instead)	5182-9775	1
6	Valve body	5182-9757	1
7	Sealing element	5182-9760	1
8	Ferrule, 1/16-inch PTFE, for septumless head carrier gas tube	0100-1375	1
9	PTV column adapter tube (includes 1/6-inch nut)	G2617-80550	1
10	PTV septumless head weldment (includes jumper tube to septum purge channel)	G3500-80000	1

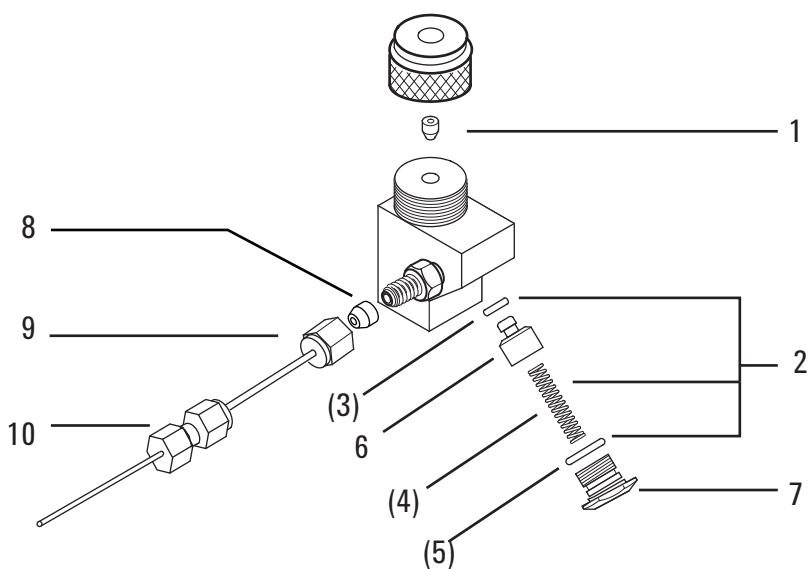

**PTV Split vent**


	Description	Part number	Qty
NS	Replacement filter kit	G1544-80530	1
1	Proportional valve, split vent	G3430-60527	1
2	Screw, M3 x 16 mm, T-10	0515-1141	2
NS	O-ring, 12/pk (2 for EPC module seal, 2 for valve seal)	5180-4181	2
3	Flow block, split vent trap	G3452-60506	1
4	Trap nut	G1544-20650	1
5	PTV front trap tubing assembly	G2617-80520	1
6	Screw, M4 x 12 mm, T-20	0515-2496	1
7	Nut, 1/16-inch SS	0100-0053	1
8	Ferrule, 1/16-inch SS	0100-1490	1
NS	Split vent trap	See Parts Finder	
NS	Screw, M4 x 12 mm, T-20 (for trap retainer)	0515-2496	1
NS	Trap retainer	G3452-00005	1

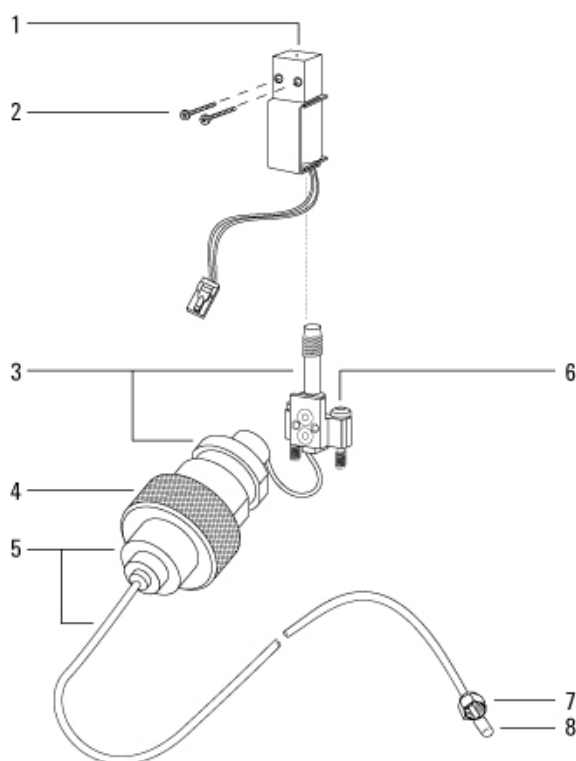

**PTV Temperature controller**


Item	Description	Part number	Qty
1	Cable, thermocouple	G2617-60505	1
3	Thermocouple interface board	G2617-60010	1
4	Screw, M4 x 12 mm long, T-20 (attach board to bracket)*	0515-2496	2
5	Screw, M4 x 12 mm long, T-20 (attach bracket to carrier)	0515-2496	2
6	Bracket, PTV thermocouple board and LC02 valve	G3430-00058	1

\*. Some 7890A assemblies use Screw, M4 x 6 mm long, T-20 (attach board to bracket), p/n 0515-2832.

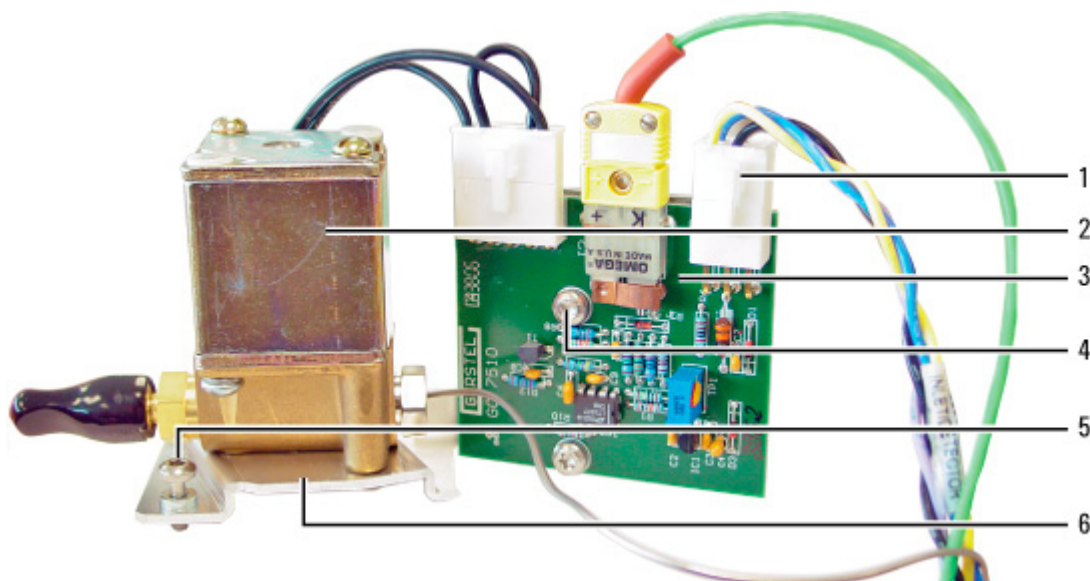


## PTV CO2 cryo assembly

Item	Description	Part number	Qty
1	7890 PTV LC02 Cyro Assembly (recommended repair level)	G3500-60515	1
2	Screw, #8-32 x .375 inch long (to attach valve to bracket)	2510-0102	2
3	Bracket, PTV thermocouple board and LC02 valve	G3430-00058	1
4	Valve, PTV LC02 cryo (included in G3500-60515)		1
5	Cap, protective without flange	1401-0044	1
6	PTV CO2 brazement (valve to IP tube)	G2617-80510	1
NS	1/8-inch male connector without O-ring	G1543-80025	1
NS	Gland seal CO <sub>2</sub>	G1565-20590	1
NS	PTFE seal	G1565-20840	1

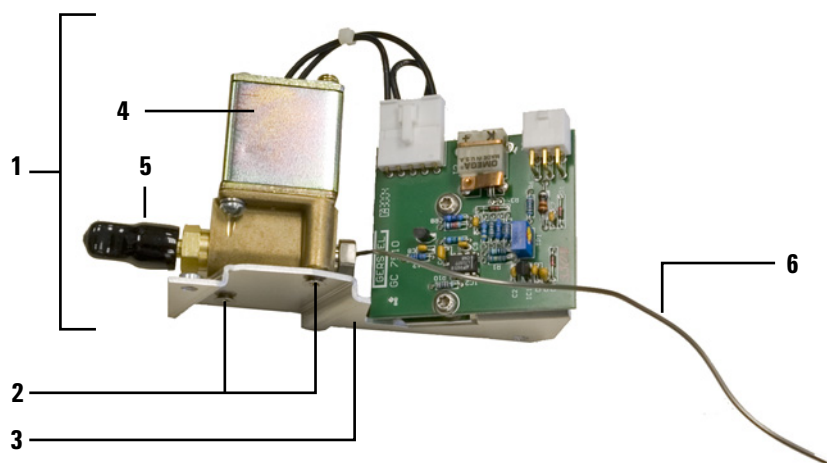

PTV LN<sub>2</sub> cryo assembly, 7890A

Item	Description	Part number	Qty
1	PTV LN <sub>2</sub> nut plate	G2619-00040	1
2	PTV LN <sub>2</sub> supply tube (included in 8)	G2619-20520	
3	Screw, M4 x 12 mm, T-20 (attach nut plate)	0515-2496	1
4	PTV LN <sub>2</sub> insulation cover	G2619-00030	1
5	Screw, M4 x 12 mm, T-20 (attach cover)	1390-1024	2
6	PTV LN <sub>2</sub> cell inlet tube, for front inlet PTV LN <sub>2</sub> cell inlet tube, for back inlet	G2619-20510 G3500-20006	1
7	Screw, 8-32 x .312-inch, pan head, pozi (included in 8)	2510-0043	2
8	7890 PTV LN <sub>2</sub> Cryo Assembly (recommended valve repair level). Includes cryo valve, mounting bracket, PTV jumper cable, LN <sub>2</sub> supply tube, and mounting screws.	G3502-60510	1
9	PTV inlet tube insulation (included in 8)	G2619-00010	1
10	Mounting bracket, LN <sub>2</sub> cryo valve (included in 8)	G3502-00001	1
11	Screw, M4 x 12 mm long, T-20	0515-2496	3
NS	PTV jumper cable (included in 8)	G2619-60502	1
NS	Connector, male (included in 8)	0100-0208	1

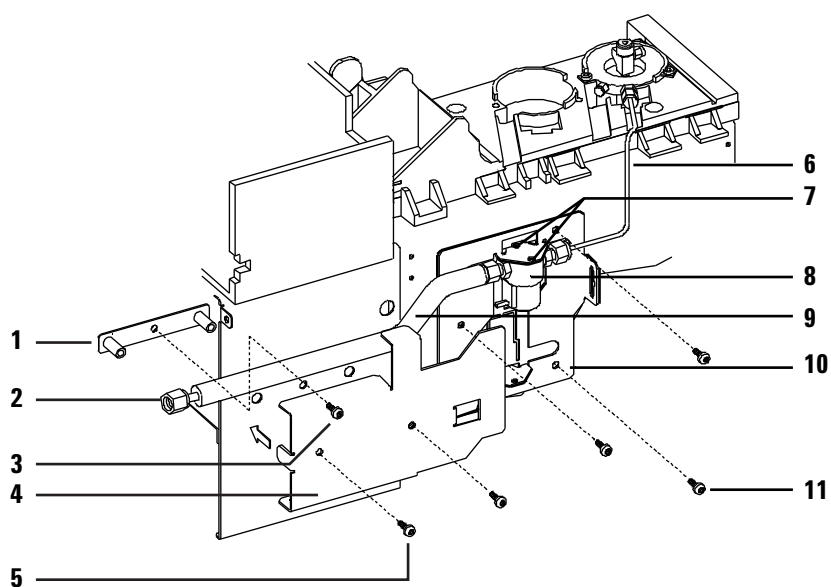




# PTV LN2 cryo assembly, 7890B

Item	Description	Part number	Qty
1	PTV LN <sub>2</sub> supply tube assembly, back	G3510-60008	1
2	PTV supply tube insulation	G3510-00021	1
3	PTV LN <sub>2</sub> cell inlet tube, for front inlet PTV LN <sub>2</sub> cell inlet tube, for back inlet	G3502-20001 G3502-20002	1
4	Screw, 8-32 x .312-inch pan head, pozi (included in 5)	2510-0043	2
5	7890 PTV LN2 Cryo Assembly (Recommended valve repair level. Includes: cryo valve mounted on cryo bracket, with LN2 supply tubes and connector between the LN2 supply tube and the cryo valve.)	G3502-60515	1
	Mounting bracket, LN2 cryo valve (included in 5)	G3430-00021	1
	PTV jumper cable (included in 5)	G2619-60502	1
6	Connector, male, 1/4-in. brass	0100-0208	1
8	Screw, M4 x 12 mm, T-20	0515-2496	4

## 18 Illustrated Parts Breakdown

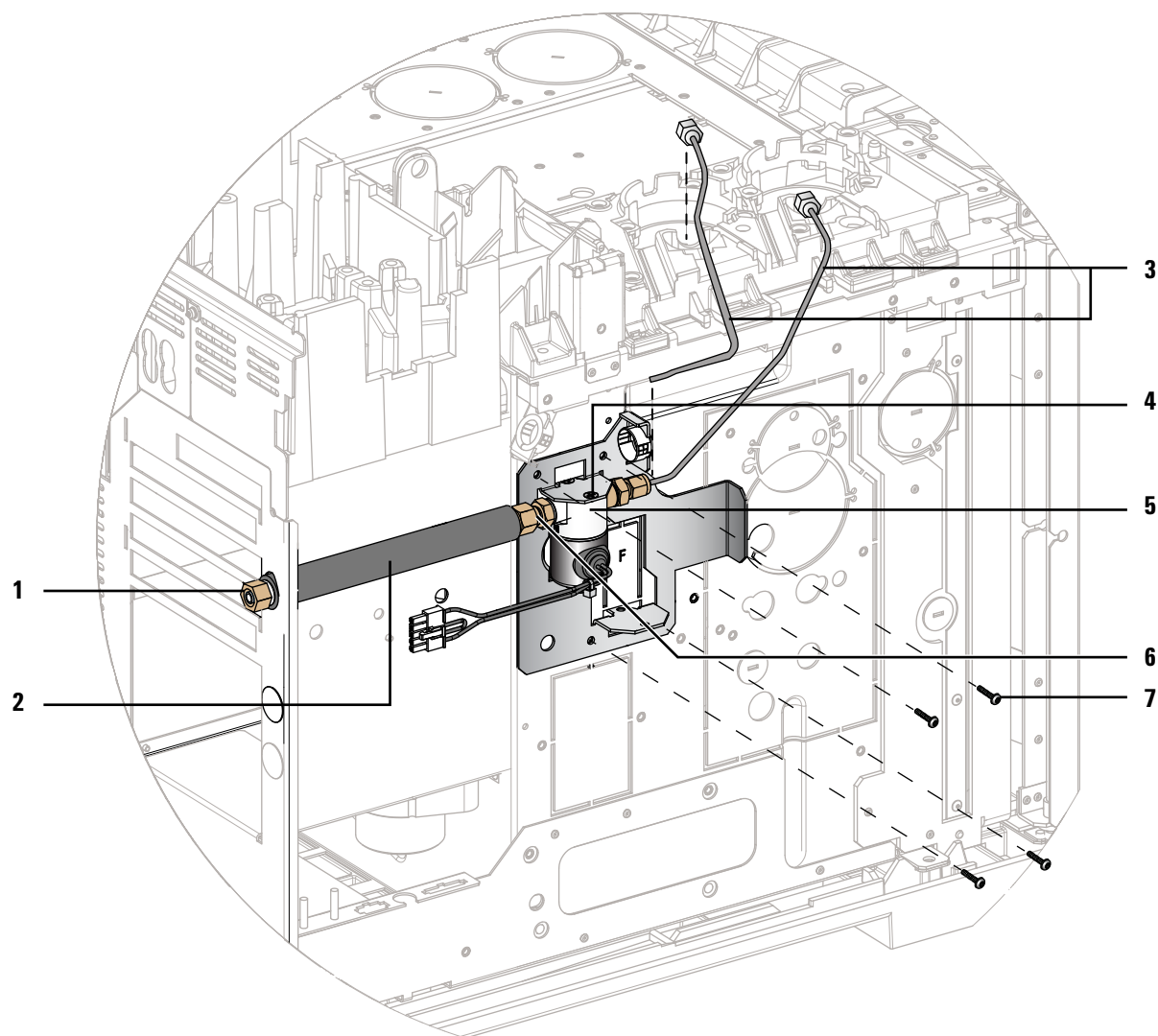


## Volatiles Interface (VI)

Top level subassemblies for VI inlet accessories:


- G3504A 100 psi Volatiles Interface with EPC, Kit

Description	Part no.
VI Inlet Accessory 0-100 PSIG without ship kit	G3504-65500
VI Weldment Assembly without EPC module	G3504-60506
VI EPC module, 100 PSI	G3504-60501

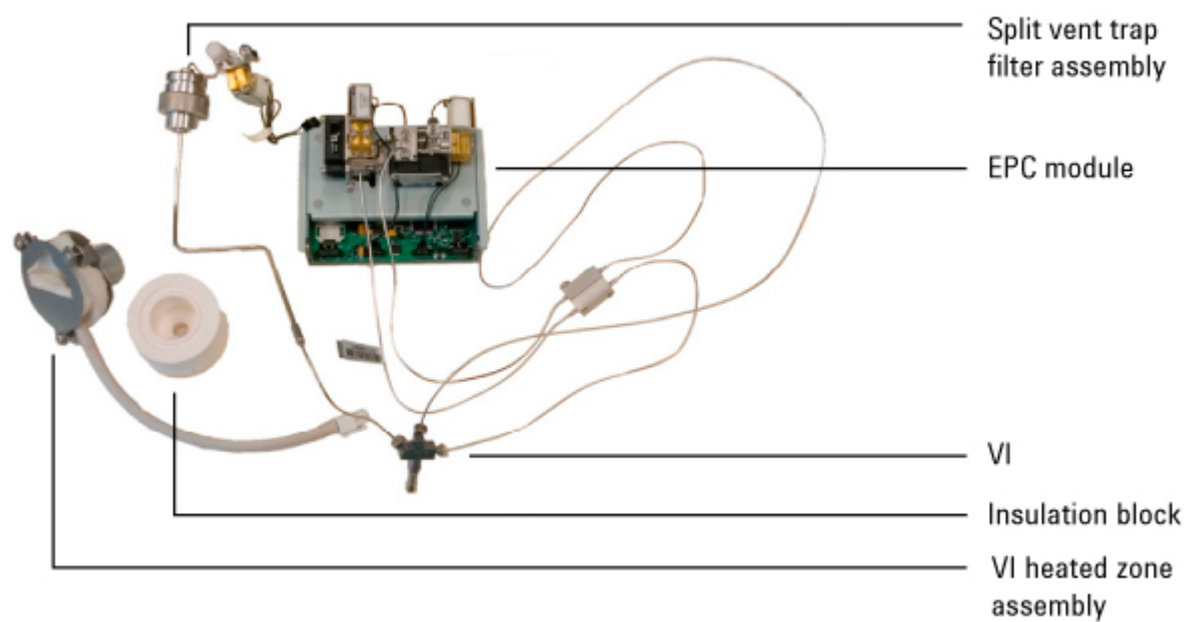


**Volatiles Interface (VI)**

<b>Item</b>	<b>Description</b>	<b>Part number</b>	<b>Qty</b>
1	VI top cover	G2319-00020	1
2	Top insulation	G2319-20580	1
3	Screw, M4 x 12 mm	1390-1023	4
4	Clamping plate	G2319-20540	1
5	Volatiles interface	G2319-60505	1
6	VI heater/sensor assembly	G2319-60503	1
7	Heater block assembly	G2319-60507	1
8	VI insulation, heater block	G2319-20530	1
9	Pneumatic fitting assembly	G3504-60500	1
10	Nut, 1/16-inch	19258-20830	1
11	SilTite Metal Ferrules, 1/16-in. x 0.5 mm ID (10/pk)	5184-3570	1

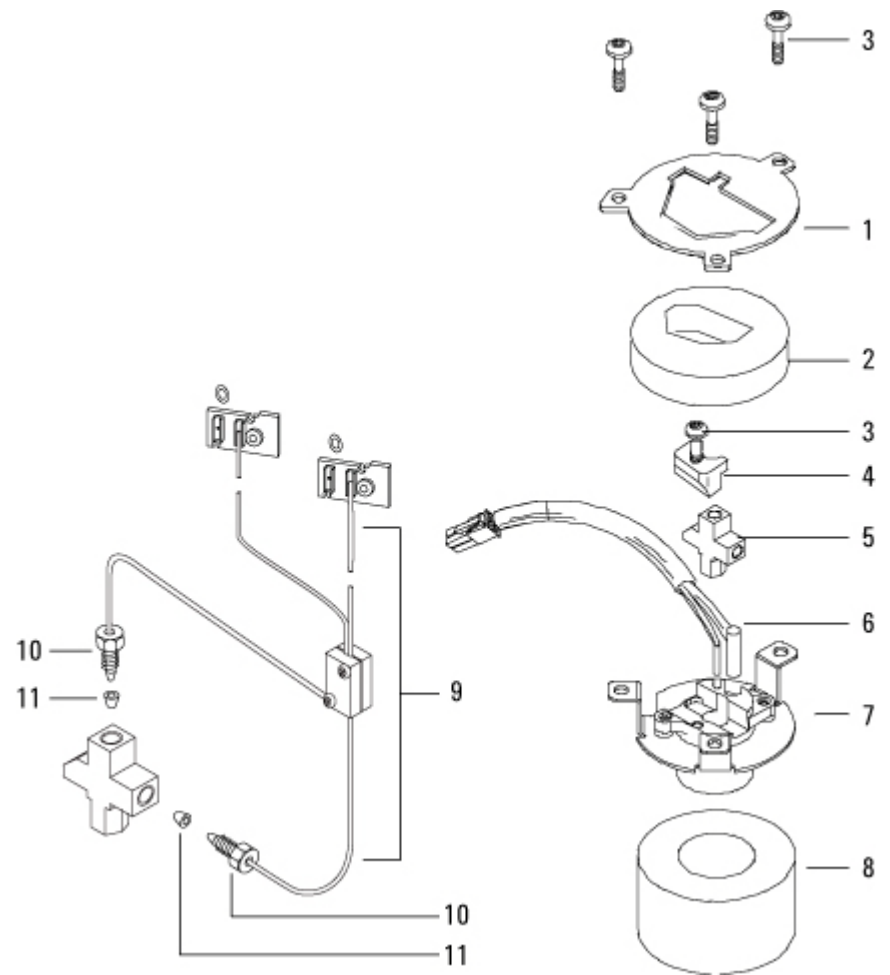

**Volatiles Interface Trap**

<b>Item</b>	<b>Description</b>	<b>Part number</b>	<b>Qty</b>
1	T08PC Valve Chemraz 0.01	G3430-60527	1
2	Screw, M3 x 16 mm, T-10, to attach valve	0515-1141	2
NS	O-ring, 12/pk for valve seal	5180-4181	2
3	Flow block, split vent trap	G3452-60506	1
NS	Screw, M4 x 12 mm, T-20	0515-2496	1
NS	Trap retainer	G3452-00005	1
4	Trap nut	G1544-20650	1
5	VI front trap assembly	G2319-80530	1
NS	Split Vent Trap	See Parts Finder	
6	Screw, captive, M4 x 20 mm, T-20	1390-1024	2
7	Nut, 1/16-inch	19258-20830	1
8	SilTite Metal Ferrules, 1/16" x 0.5mm ID (10/pk)	5184-3570	1

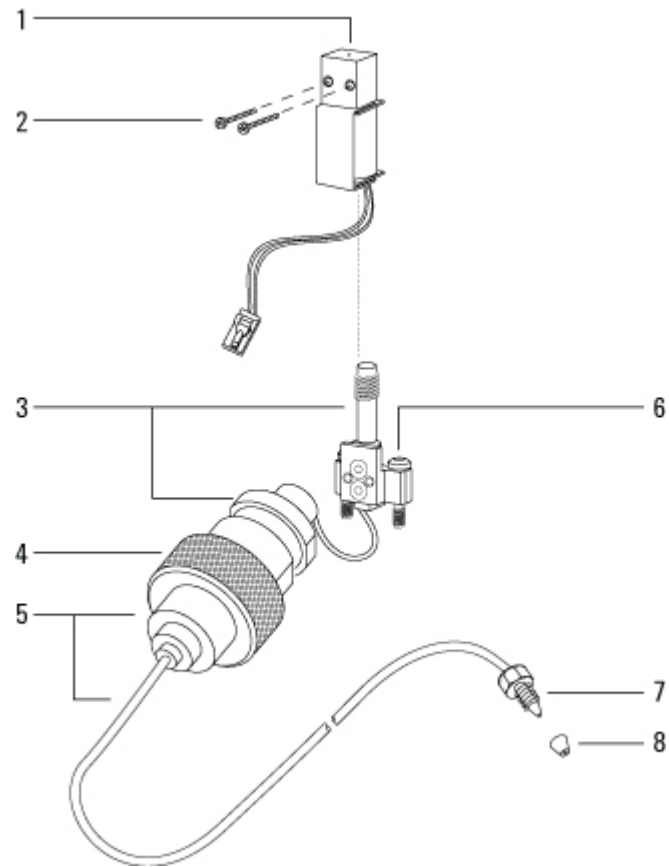

## **Detectors**

This section contains an illustrated parts breakdown for each of the following 7890 GC detectors and related components.

- Flame Ionization Detector (FID)
- Flame Photometric Detector (FPD)
- FPD+ (High Temperature FPD)
- Microcell Electron Capture Detector (uECD)
- Nitrogen Phosphorus Detector (NPD)
- Thermal Conductivity Detector (TCD)




## Flame Ionization Detector (FID)

Top level subassemblies for FID detector accessories:


- G3461A FID with EPC, Kit adaptable to either packed or capillary columns.
- G3462A Capillary FID with EPC

Description	Part number
FID Detector Module, adaptable without ship kit	G3431-60506
FID Detector Module, capillary without ship kit	G3431-60500
FID Packed Weldment Assembly with pallet, electrometer, and base	G3431-80508
FID Capillary Weldment Assembly with pallet, electrometer, and base	G3431-80506
FID EPC Module	G3431-60531
FID Signal Board	G3431-60020

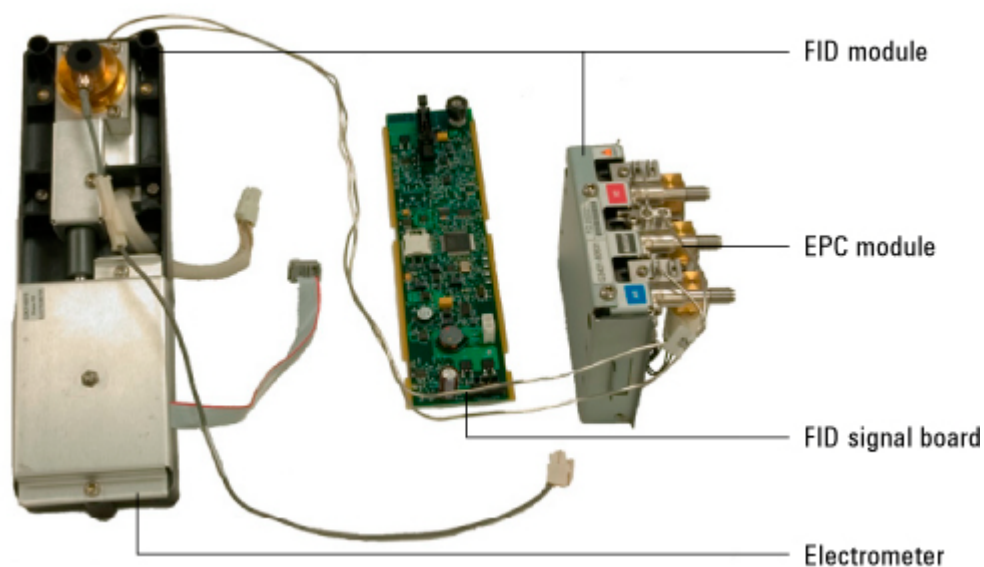

**FID detector**

Item	Description	Part number	Qty
1	Screw, M4 × 25 mm, Torx T-20	0515-2712	4
2	FID interconnect cover	G1531-00220	1
3	PTFE chimney (optional)	19231-21050	1
4	Collector assembly For individual collector parts see Parts Finder.	G1531-60690	1
5	Screw, M4 × 10 mm, Torx T-20	0515-2495	3
6	Interconnect J-clamp	19231-00040	1
7	Mounting pallet, FID	G1531-40020	1
8	<b>FID Base Assembly</b>		
9	Cylindrical EMI suppressor, ferrite	9170-1730	1
10	Spring, electrometer interconnect	1460-2142	1
11	FID interconnect assembly	G1531-60715	1
12	FID electrometer	G3431-60015	1

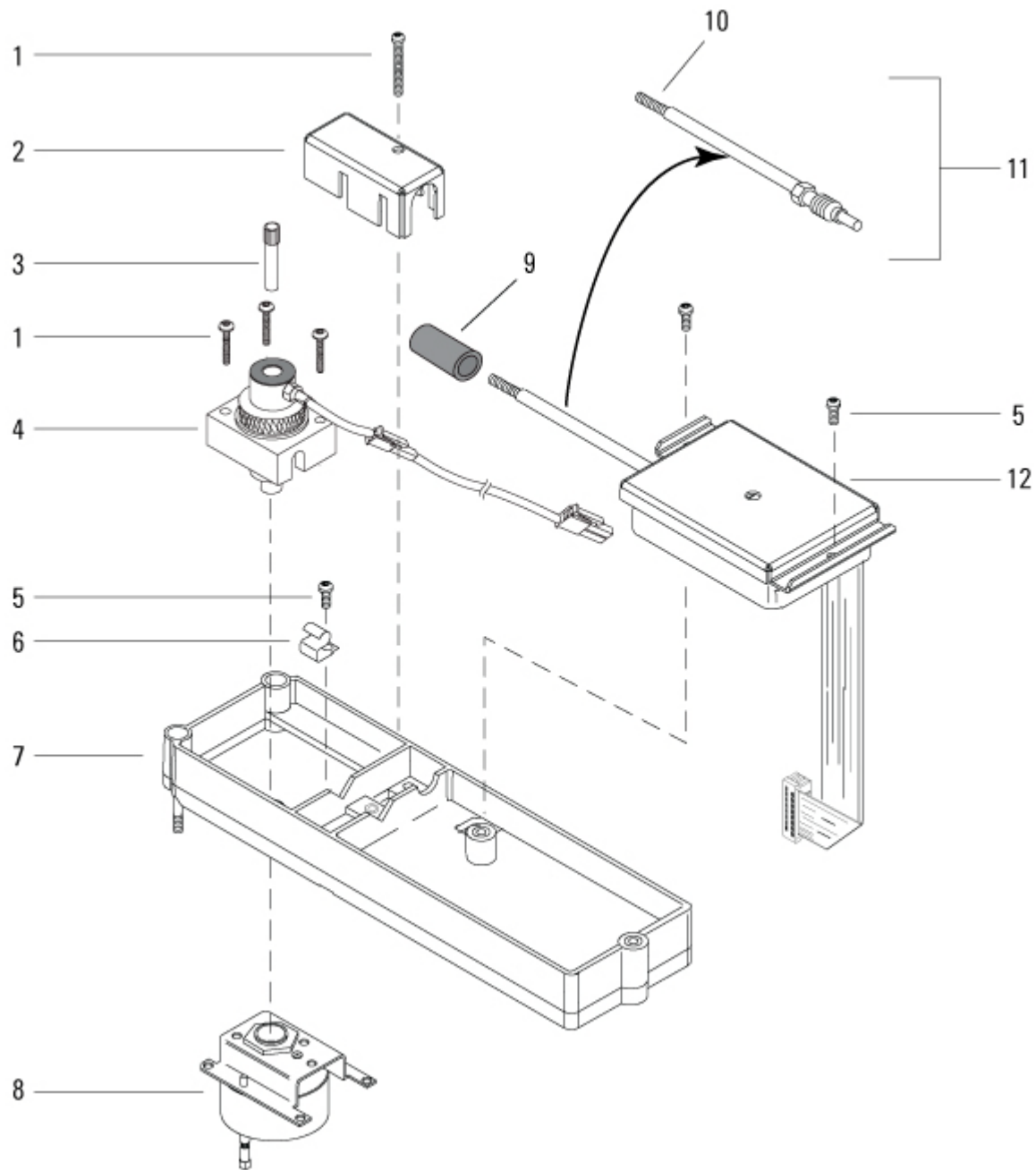


**FID Base Assembly**

<b>Item</b>	<b>Description</b>	<b>Part number</b>	<b>Qty</b>
1	Base spanner nut	19231-20990	1
2	Thermal strap	G1531-00105	1
3	Screw, captive, M3, T-10	G1946-20168	3
	O-ring, size 2-006, fluorocarbon, 12/pk	5180-4181	3
4	Jets	See Parts Finder	1
5	Adaptable FID, base weldment with jet	G3431-80509	—
	Capillary column FID, base weldment with jet	G3431-80507	
6	Screw, M4 × 25 mm, Torx T-20	0515-2712	2
7	FID block insulation	G1531-60700	1
8	Adaptable FID column adapters:		1
	— FID/NPD capillary column	19244-80610	—
	— FID/NPD 1/8-inch packed column	19231-80520	—
	— FID/NPD 1/4-inch packed column	19231-80530	—
9	Nut warmer insulation and cup assembly.	19234-60700	1
10	Heater/sensor assembly	G1530-67950	1
NS	Insulation, nutwarmer (1/pk)	19234-60715	3

G3431-80509 includes jet 19244-80560, capillary adaptable, 61.5 mm

G3431-80507 includes jet G1531-80560, capillary optimized, 43 mm, .29 mm id

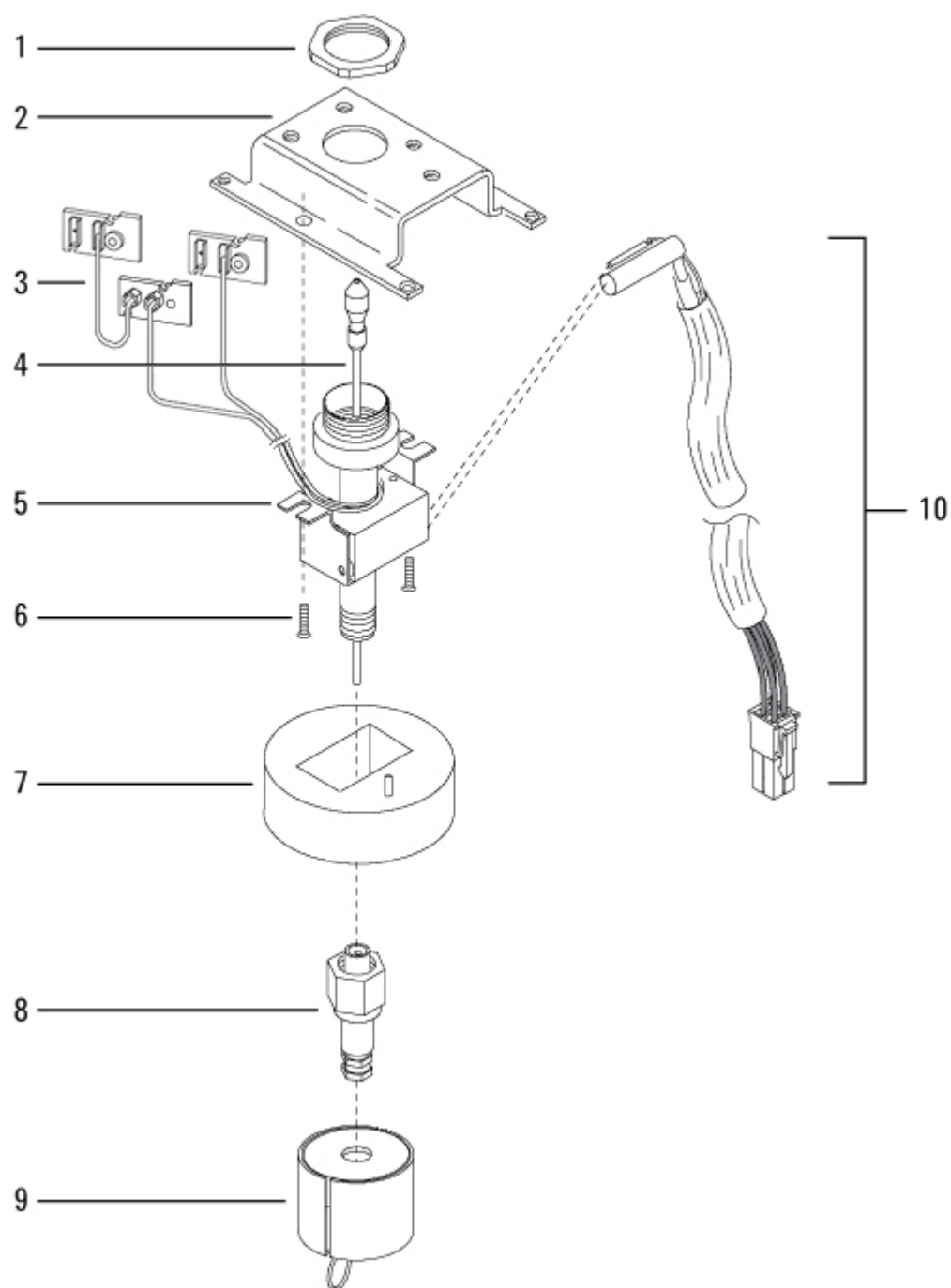

**PM kits for the FID**

See Parts Finder and the Agilent catalog for consumables and supplies for consumable parts, or visit the Agilent Web site for the latest information ([www.agilent.com/chem/supplies](http://www.agilent.com/chem/supplies)).

**FID PM kits**


Description	Part no.
Maintenance kit, FID rebuilding	G1531-67001
Maintenance kit, FID cleaning	G1531-67000

## Flame Photometric Detector (FPD)

Top level subassemblies for FPD detector accessories:

- G3435A Single FPD with EPC, Kit
- G3436A Dual FPD with EPC, Kit

Description	Part number
FPD Detector Module, Single, without ship kit	G3435-60500
FPD Detector Module, Dual, without ship kit	G3436-60500
FPD Body Assembly (without EPC, signal board, or PMT)	G3435-60575
Dual FPD Body Assembly (without EPC, signal board, or PMT)	G3435-60585
FPD EPC module	G3435-60535
FPD Signal board (1 for single FPD, 2 for dual FPD)	G3435-60010

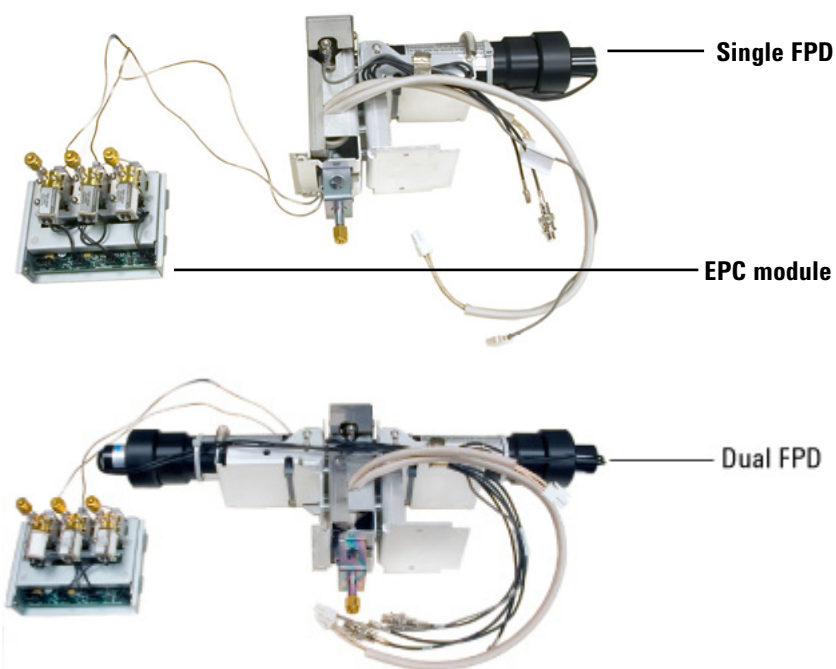




FPD Signal board



### FPD inert transfer line parts


Item	Description	Part number	Qty.
1	O-ring, FPD jet seal, white	0905-1608	1
2	Screw, M3 x 12 mm, T-10	0515-1084	1
3	Retaining clip	19256-00340	1
4	Inert transfer line	G3435-60555	1
5	Heater/Sensor assembly	G1535-60620	1
6	Capillary adapter seat	19256-21140	1
7	Capillary adapter nut	19256-21150	1
8	1/8-inch ferrule (10/pk)	0100-1332	1
9	1/8-inch nut (10/pk)	5180-4103	1
10	1/4-inch packed adapter	G1532-20710	1

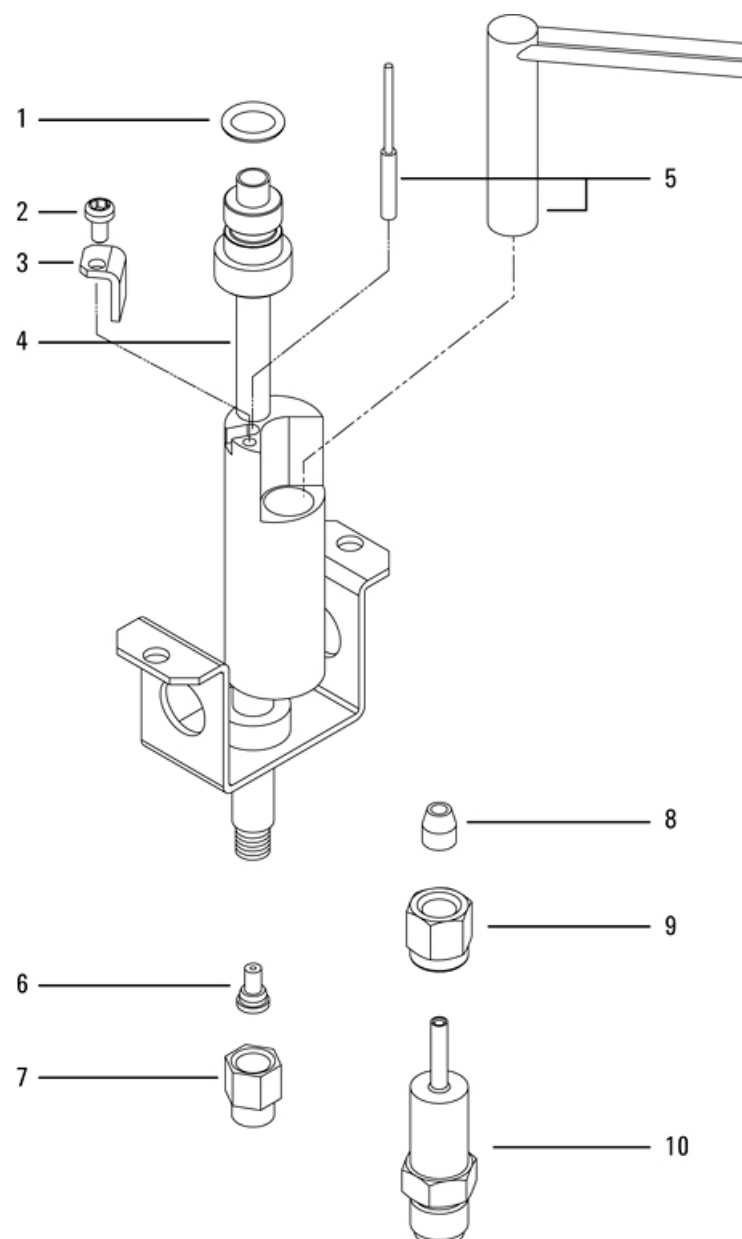

**FPD ignitor and heat shield assembly**

<b>Item</b>	<b>Description</b>	<b>Part number</b>	<b>Qty.</b>
1	FPD tube assembly: Aluminum Stainless Steel	19256-60700 19256-20705	1
2	Vespel ferrule, 1/4-inch id, 20/pk	5080-8774	1
2	Front ferrule, 1/4-inch id, brass, 10/pk	5180-4111	1
2	Back ferrule, 1/4-inch id, brass, 10/pk	5180-4117	1
3	Emissions block assembly FPD, single FPD, dual (not shown)	19256-80560 19256-80600	1
4	Ignitor replacement kit A. O-ring B. Spacer C. Glow plug	19256-60800	1
13	Screw, M3 x 66 mm, T-10	0515-0680	1
12	Collet, ignitor	19256-20690	1
5	Ignitor cable assembly	G3435-60600	1
6	Heat shield gasket, white	19256-80045	1
7	First heat shield window	19256-80030	1
8	Heat shield disk	19256-20580	1
9	Stainless steel coupling	19256-20550	1
10	Lock washer (4 required)	2190-0584	4
11	Screw, M3 x 12, T-10 (4 required)	0515-1084	4


## 18 Illustrated Parts Breakdown

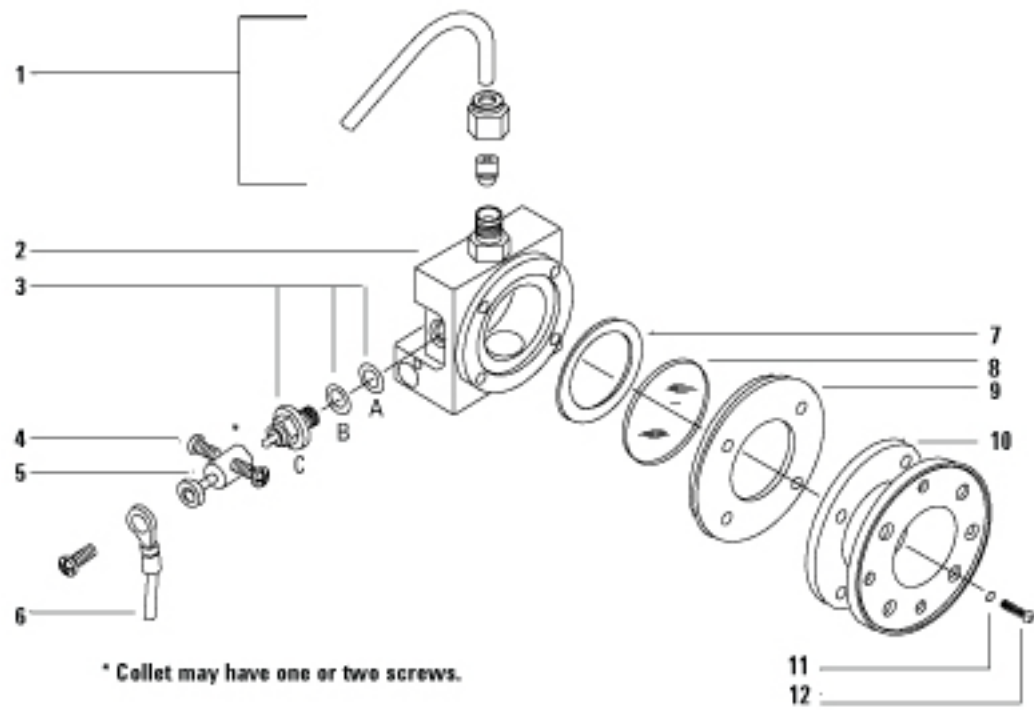

**PMT and bracket assemblies****Single FPD**


Item	Description	Part number	Qty.
1	Chimney back cover	G1535-80520	1
2	Heater/Sensor assembly	G1535-60610	1
3	Transfer line support bracket, single	19256-00320	1
4	Bracket/Support	G1535-00010	1
5	Emissions block assembly, single	19256-60650	1
6	Filters:		1
	Sulfur	1000-1437	
	Phosphorus	19256-80010	
7	Filter spacer (used only with sulfur filter)	19256-20910	1
8	Sping, PMT	1460-1160	1
9	PMT housing assembly	19256-60510	1

## 18 Illustrated Parts Breakdown

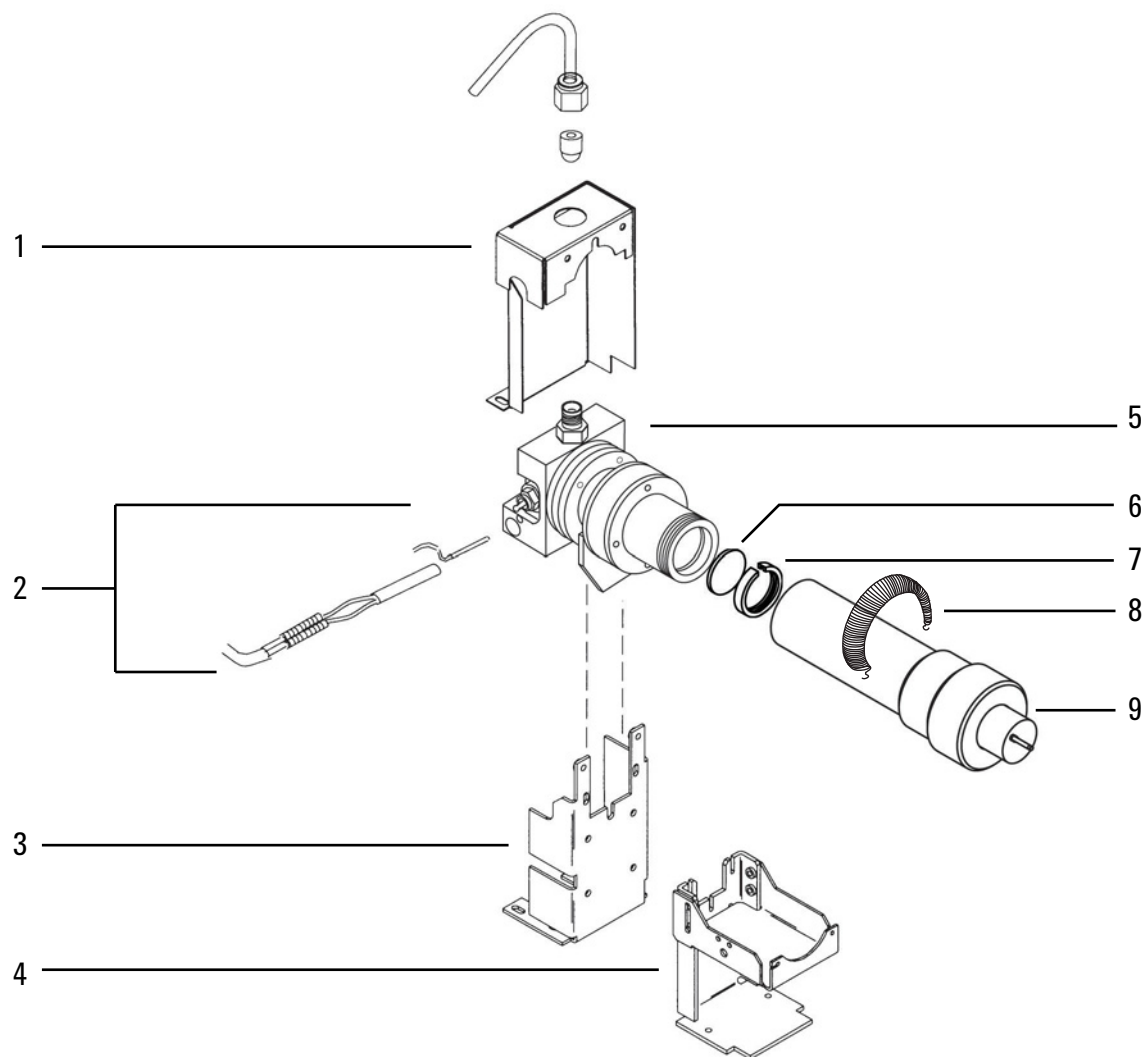


## Dual FPD

Item	Description	Part number	Qty.
11	Dual cover assembly	G1535-80530	1
12	Emissions block assembly, dual	19256-60690	1
13	Dual FPD chimney back	19256-00330	1
14	Dual main bracket	G1535-00040	1
15	Dual FPD chimney front	G1535-00030	1

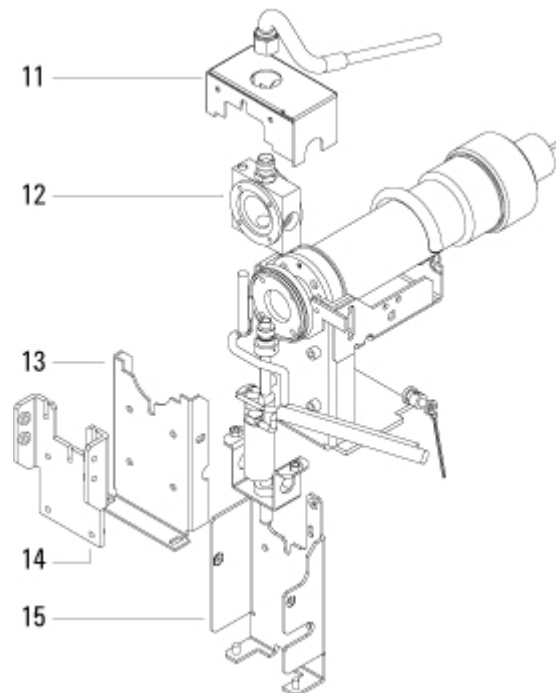

## FPD lens assembly


Item	Description	Part number	Qty.
1	Clamp	19256-00090	1
2	Screw, M3 × 25 (4 required)	0515-0683	4
3	Silicone O-ring, 0.926-inch id (orange)	0905-1099	1
4	Convex lens	1000-1438	1
5	Lens housing	19256-20900	1
6	Flange ring	19256-00200	1
7	O-ring 12/pk, Fluoroelastomer, 1.239-in. ID	5061-5890	1

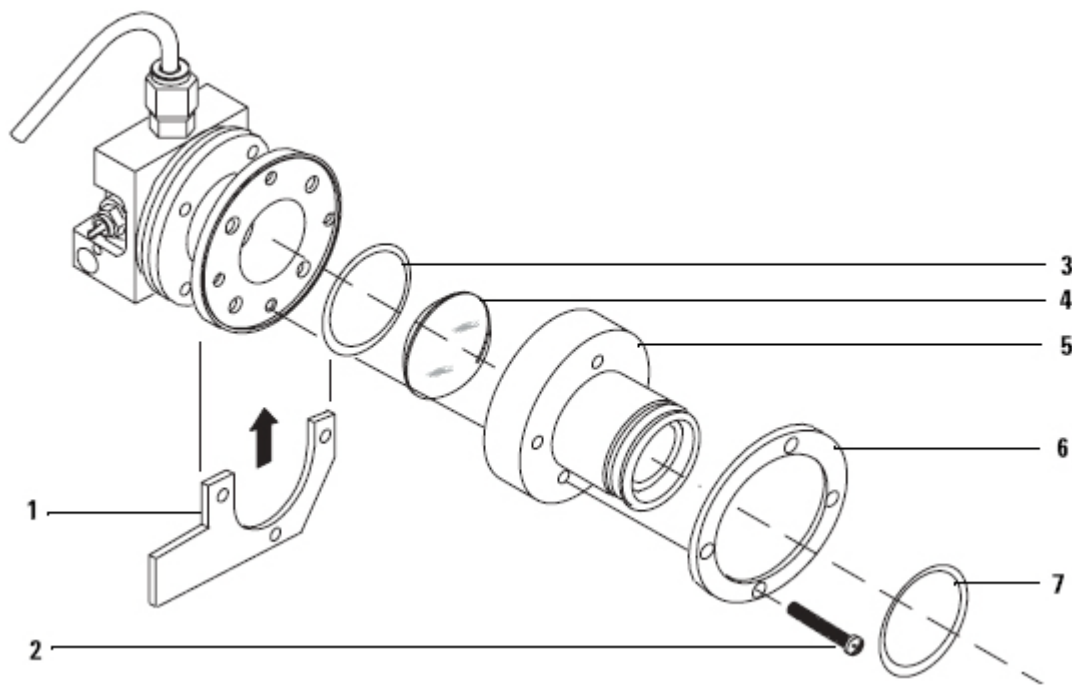




## FPD Covers

Item	Description	Part number
1	Single FPD top cover	G3435-80500
2	Dual FPD top cover	G3435-80501
3	Electronics top cover, with ground strap Ground strap	G3435-60503 19256-60730
4	Auxiliary heater and valve control assembly Aux zone/valve box harness	 G1530-60660

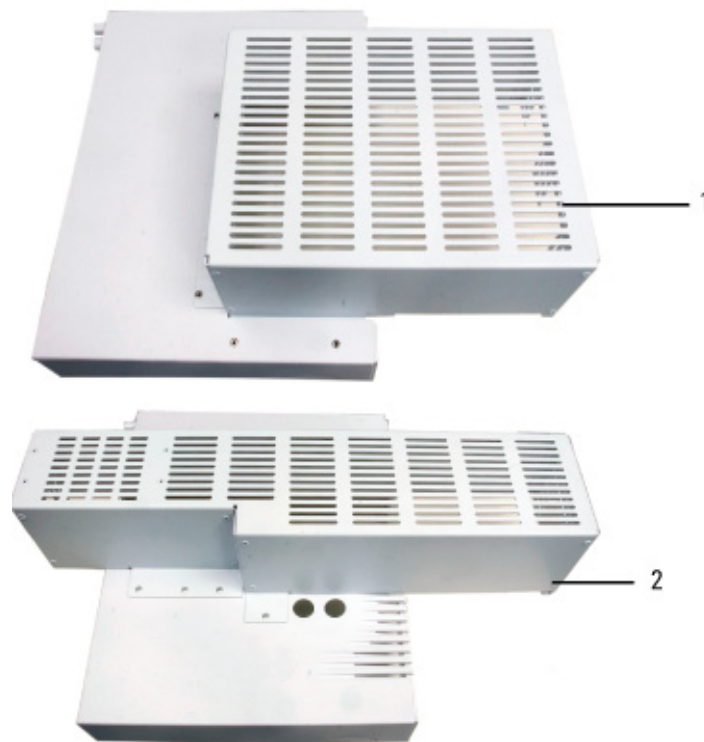

## 18 Illustrated Parts Breakdown

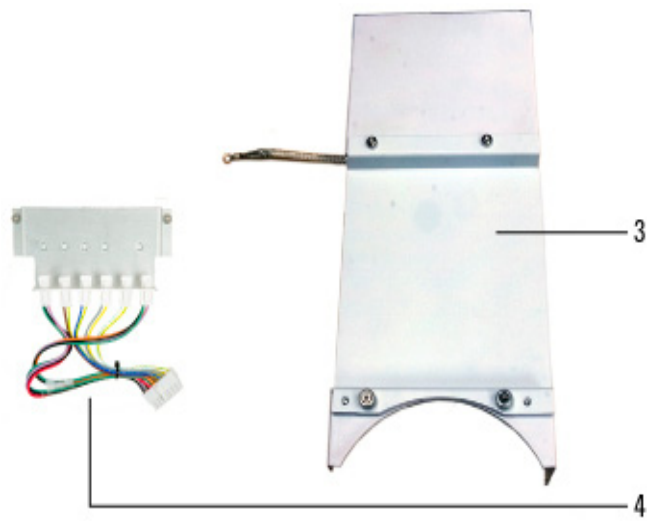

FPD<sup>+</sup> (High Temperature FPD)

Top level subassemblies for FPD<sup>+</sup> detector accessories:

- G3435B Single FPD<sup>+</sup> with EPC, Kit
- G3436B Dual FPD<sup>+</sup> with EPC, Kit

Description	Part number
FPD High Temp, Single	G3435-60300
FPD High Temp, Dual	G3436-60300
FPD <sup>+</sup> EPC module	G3435-60335
FPD Signal board (1 for single FPD, 2 for dual FPD)	G3435-60010



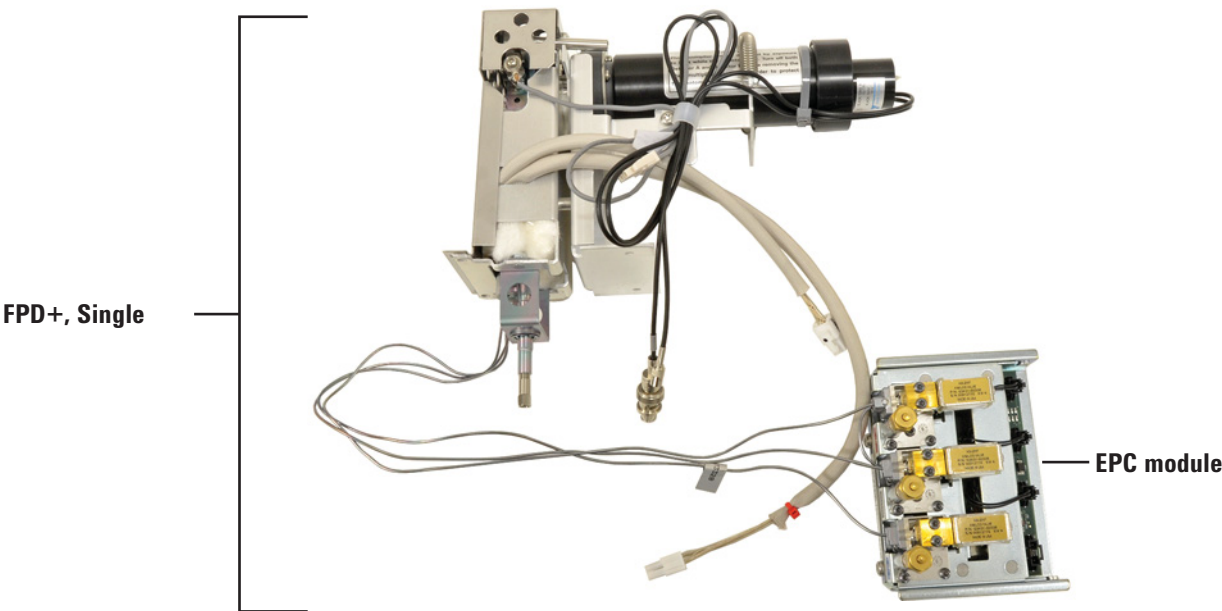


Figure 17 FPD<sup>+</sup>, Single

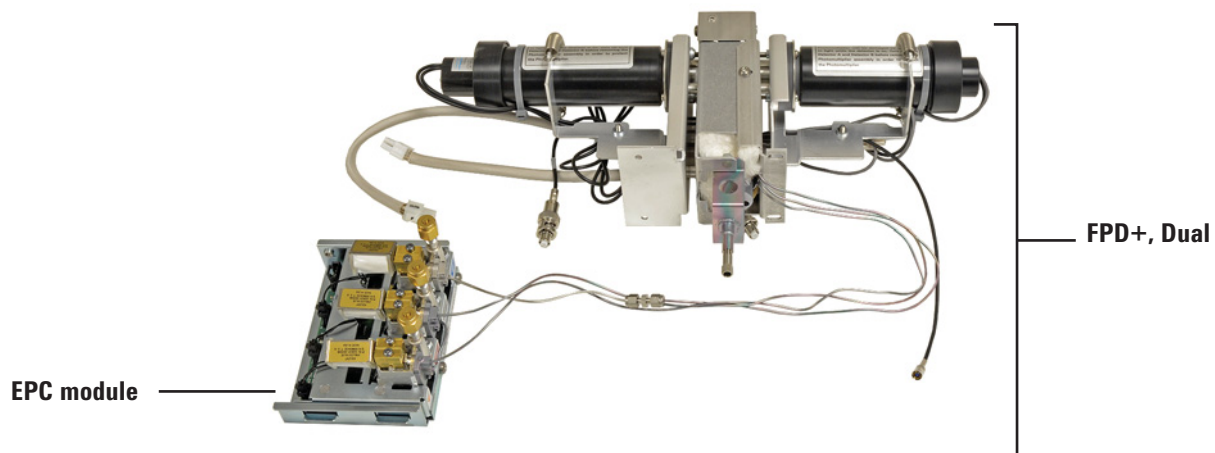


**Figure 18** FPD+, Dual

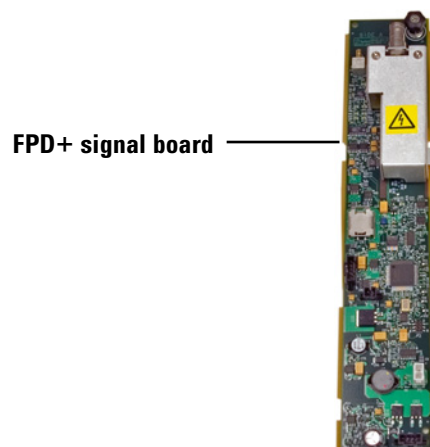

## FPD<sup>+</sup> Assemblies


**Table 38** Coupling, ignitor, and exit tube

ID	Description	Part number	Qty
1	O-Ring, 2-010, Fluorocarbon, brown (exit tube)	0905-1014	1
2	Exit tube assembly, FPD <sup>+</sup>	G3435-60330	1
3	Screw, M3 X 6 mm long, pan head	0515-0680	2
4	FPD Ignitor cable assembly	G3435-60600	1
5	Collet, glow plug	19256-20690	1
6	Ignitor, cleaned w/copper gasket (for FPD <sup>+</sup> , discard copper gasket)	19256-60750	1
7	Heat shield gasket, perfluoroelastomer	19256-80045	1
8	Window, first	19256-80030	1
9	Shield, heat	19256-20580	1
10	Lens housing coupling	G3435-20345	1
11	Screw, machine, M3 X 12 mm long, pan head	0515-1084	4
12	Washer, lock, HLCL 3.0 3.1 MM-ID 6.2 MM-OD	2190-0584	4
NS	Washer-Flat metallic 3.2 mm-ID 9 mm-OD 1mm	3050-2326	
NS	Spacer-Round 3.2 mm-ID 6 mm-OD 6mm-LG Br Z	0380-4876	

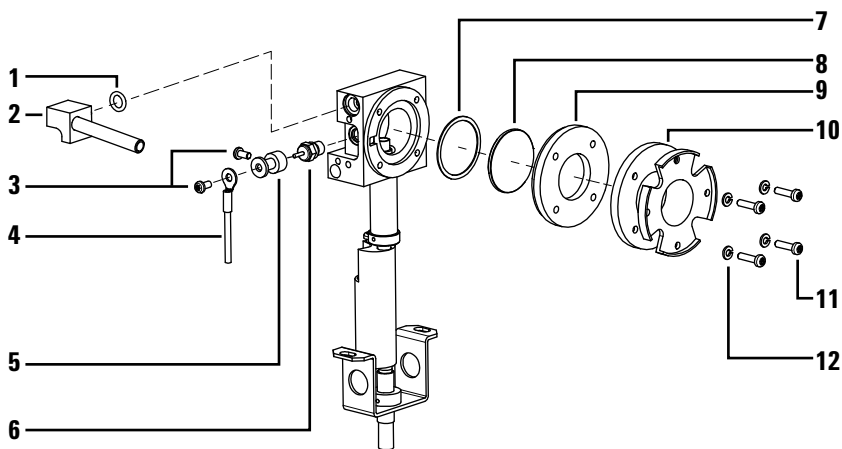


**Table 39** Lens housing

ID	Description	Part number	Qty
1	BK7 Convex glass lens, 31.2mm f.l.	1000-1438	
2	O-Ring 2-021 Silicone Rust Color	0905-0955	1
3	Clamp	19256-00090	1
4	Lens housing	G3435-20391	
5	RING, FLANGE	19256-00200	
6	SCREW, 3.0 25 PN TX	0515-0683	4
7	O-Ring, Size 2-123, Silicone, Rust Color	0905-1159	1

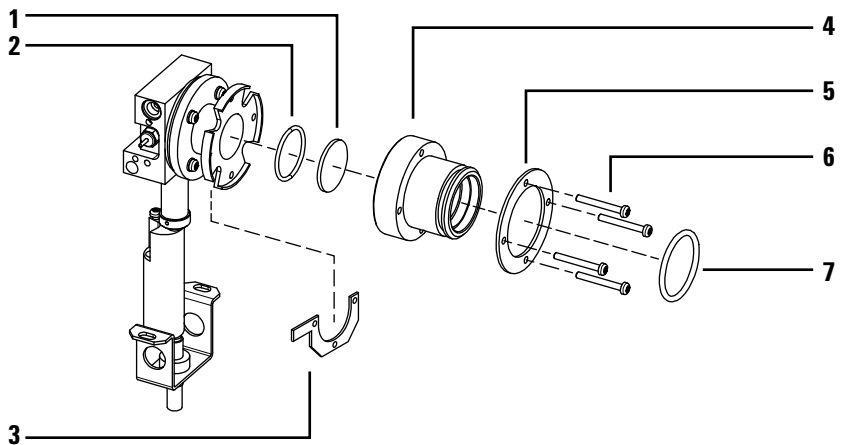
**Table 40** Brazements and heaters

ID	Description	Part number	Qty
1	WS2 Coated M3x12 Screw	G3435-20357	1
2	Retaining Clip, FPD Heater-Sensor	19256-00340	
3	FPD UP HEAT/SENS ASSY	G1535-60610	1
4	FPD Cell Brazement - Inert	G3435-60558	1
5	FPD Emission Block Heater Assembly	G3435-60610	1

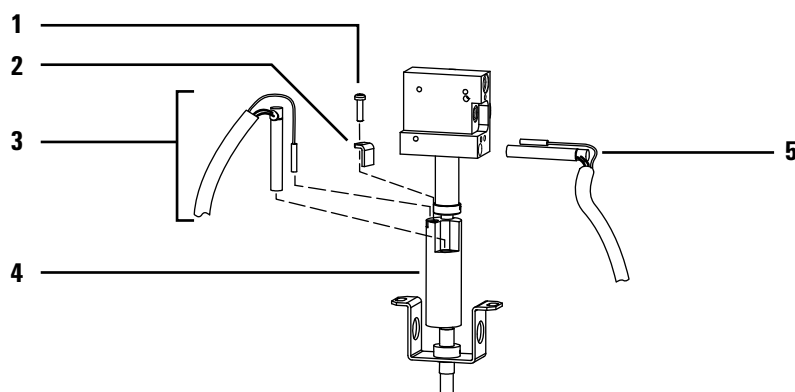

**Table 41** Brackets, insulation, PMT, and filters

ID	Description	Part number	Qty
1	Front detector cover assembly	G3435-81330	1
2	Chimney insulation	G3435-00340	1
3	Nut, hex, with lockwasher, M4	0535-0043	2
4	Chimney Insulation, Middle	G3435-00342	1
5	Screw, machine, M4 x 10-mm long	0515-2495	5
6	FPD Chimney back	19256-00320	1
7	WHC-375 Wire Clip	1400-3196	1
8	Filters		1
	FPD Sulphur filter 394nm	1000-1437	
	Filter, P 5890 FPD	19256-80010	
9	FPD Filter spacer (for sulfur filter only)	19256-20910	1
10	Spring, compression	1460-1160	1
11	PMT Assembly (single FPD+ and dual FPD+ right side)	G3435-60360	1
12	FPD Main bracket	G3435-00370	1

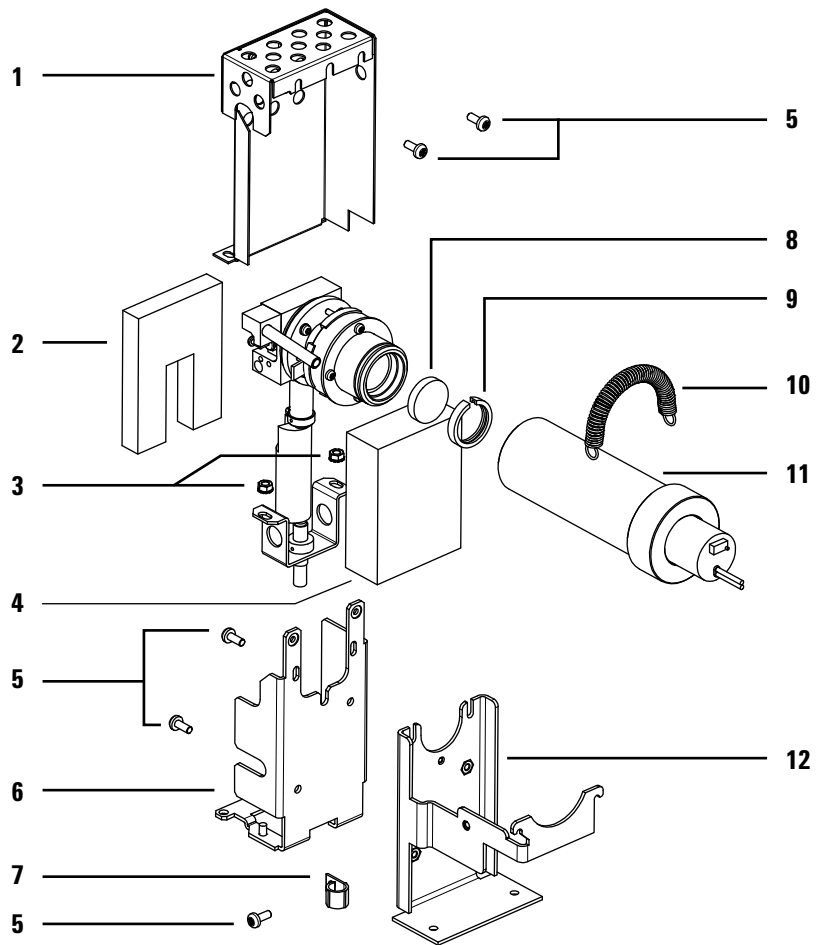

**Table 42** Covers (external)

ID	Description	Part number	Qty
	M3 X 8mm Captive WS2 coated screw	G3430-20028	
1	Single FPD+ top cover	G3435-81350	1
3	Ground strap	19256-60730	1
4	Electronics top cover, with ground strap	G3435-60503	1






See [Table 43](#) below for additional parts used only on the dual FPD<sup>+</sup>.

**Table 43** Dual FPD<sup>+</sup> parts

ID	Description	Part number	Qty
1	Dual FPD cover assembly	G3435-81360	1
2	FPD Cell brazement, dual - inert	G3435-60559	1
3	Dual wavelength main bracket	G3435-00375	1
4	PMT Assembly (Dual FPD+ only, left PMT, short)	G3435-60370	1
5	DW FPD Chimney front	G1535-00030	1
6	PMT Assembly (right PMT, long)	G3435-60360	1
7	WHC-375 Wire Clip	1400-3196	2
8	Screw, machine, M4 x 10-mm long	0515-2495	5

## 18 Illustrated Parts Breakdown

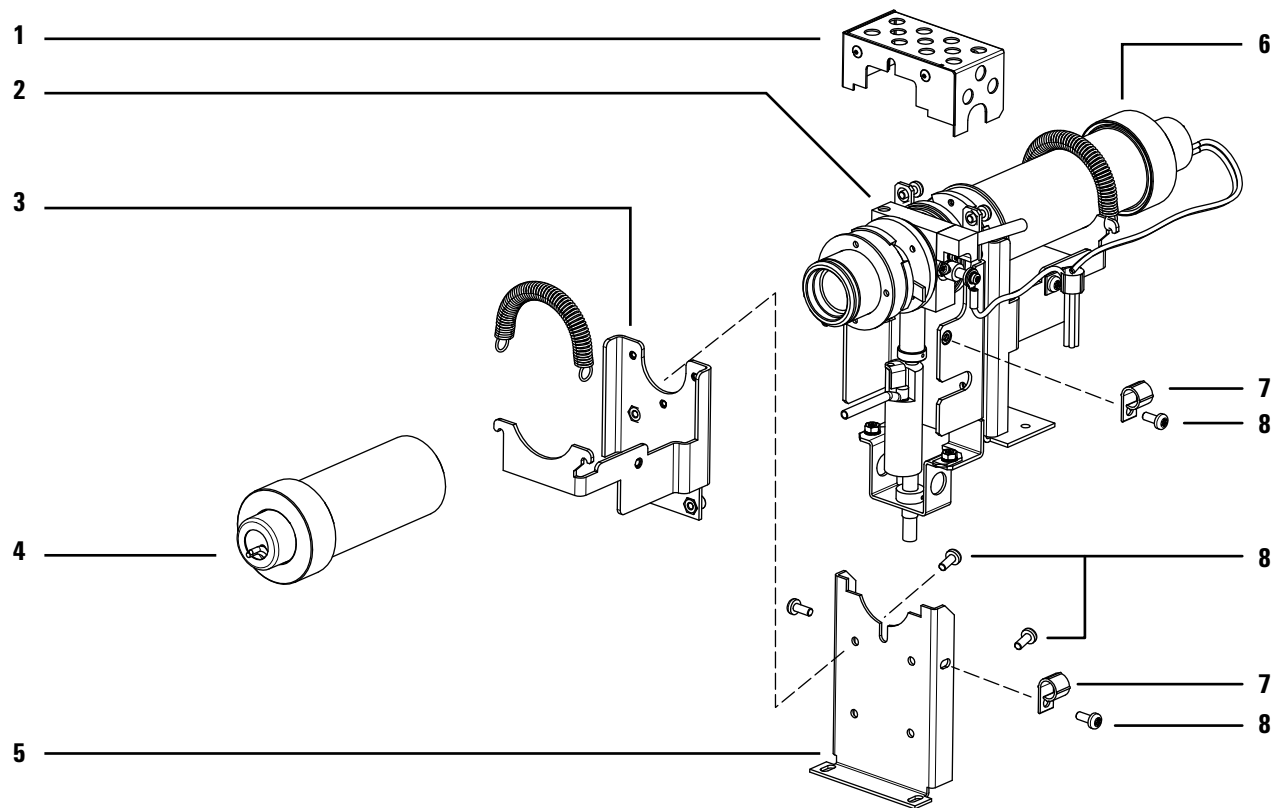


Microcell Electron Capture Detector (uECD)

Top level subassemblies for  $\mu$ ECD detector accessories:


- G2397AD  $\mu$ ECD with EPC Accessory
- G2398AD  $\mu$ ECD with EPC Accessory for Japan

Description	Part number
$\mu$ ECD EPC module	G3433-60533
$\mu$ ECD Signal board	G3433-61020*
* Replaces $\mu$ ECD Signal board G3433-60020.	

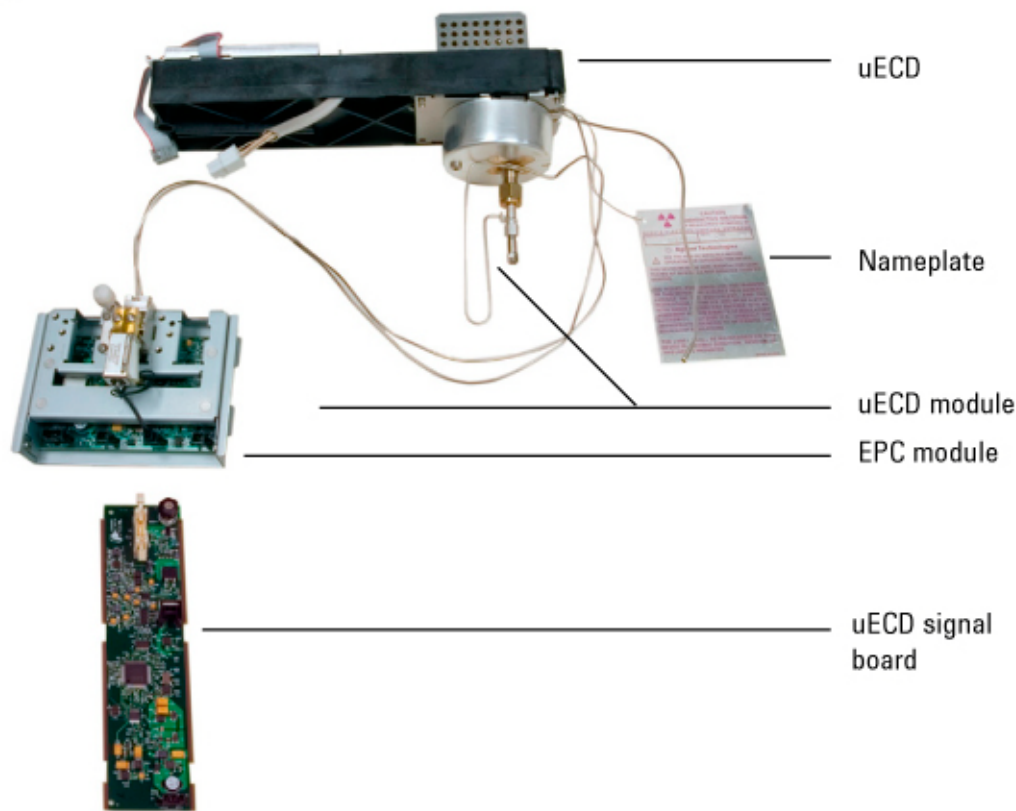

## Microcell Electron Capture Detector (uECD) parts

Item	Description	Part number	Qty
1	uECD top cover	G2397-00011	1
2	uECD signal wire assembly	19233-60635	1
3	Anode/Ferrule/Nut assembly (specific license required)	—	1
4	uECD thermal cover	G1533-00030	1
5	Thermal cover clip	19233-00095	1
6	uECD top insulation	G1533-00020	1
7	Screw, M4 × 10 mm, Torx T-20	0515-2495	1
8	Clamp, interconnect	19231-00040	1
9	Captive screw, detector pallet, M4 × 20 mm	1390-1024	4
10	FID pallet	G1531-40020	1
13	Tubing, Tygon, 30 inch	0890-0934	1
14	Screw, M4 × 45 mm, Torx T-20	0515-2484	2
15	Screw, M4 × 12 mm, Torx T-20	0515-2496	
16	ECD Insulation cup kit	G3433-00004	1
17	Screw, M4 × 40mm T-20	0515-4793	1
18	uECD make-up gas adapter weldment assembly, includes:		1
	— Fused silica indented mixing liner	G3433-63000	
	— End cap	G2397-20540	
		19233-20755	
19	Ferrules, capillary columns	See Parts Finder	
20	Capillary column nut	See Parts Finder	
21	Replacement ECD Cell (general license)	G2397-60610	1
	Replacement ECD Cell (JP Labelled)	G2397-60615	
	Replacement ECD Cell (JP General)	G2398-60615	
22	uECD interconnect assembly	G1533-60510	1
NS	Screw, M4 × 10 mm, Torx T-20 to attach electrometer	0515-2495	2
23	uECD electrometer without interconnect	G3433-60010	1
24	Upper heated block	G1533-20525	1
25	Heater/Sensor assembly	G1533-60625	1
26	Lower heated block	19233-20515	1
27	uECD mounting plate	G3433-60501	1
28	Nut warmer insulation and cup assembly. Includes:	19234-60720	1
	Nutwarmer insulation	19234-60715	
	Nutwarmer cup assembly	19234-60700	
29	Ferrule, 1/4 in. Vespel (10/pkg)	5080-8774	1
30	1/4 in. nut, brass (10/pkg)	5180-4105	1

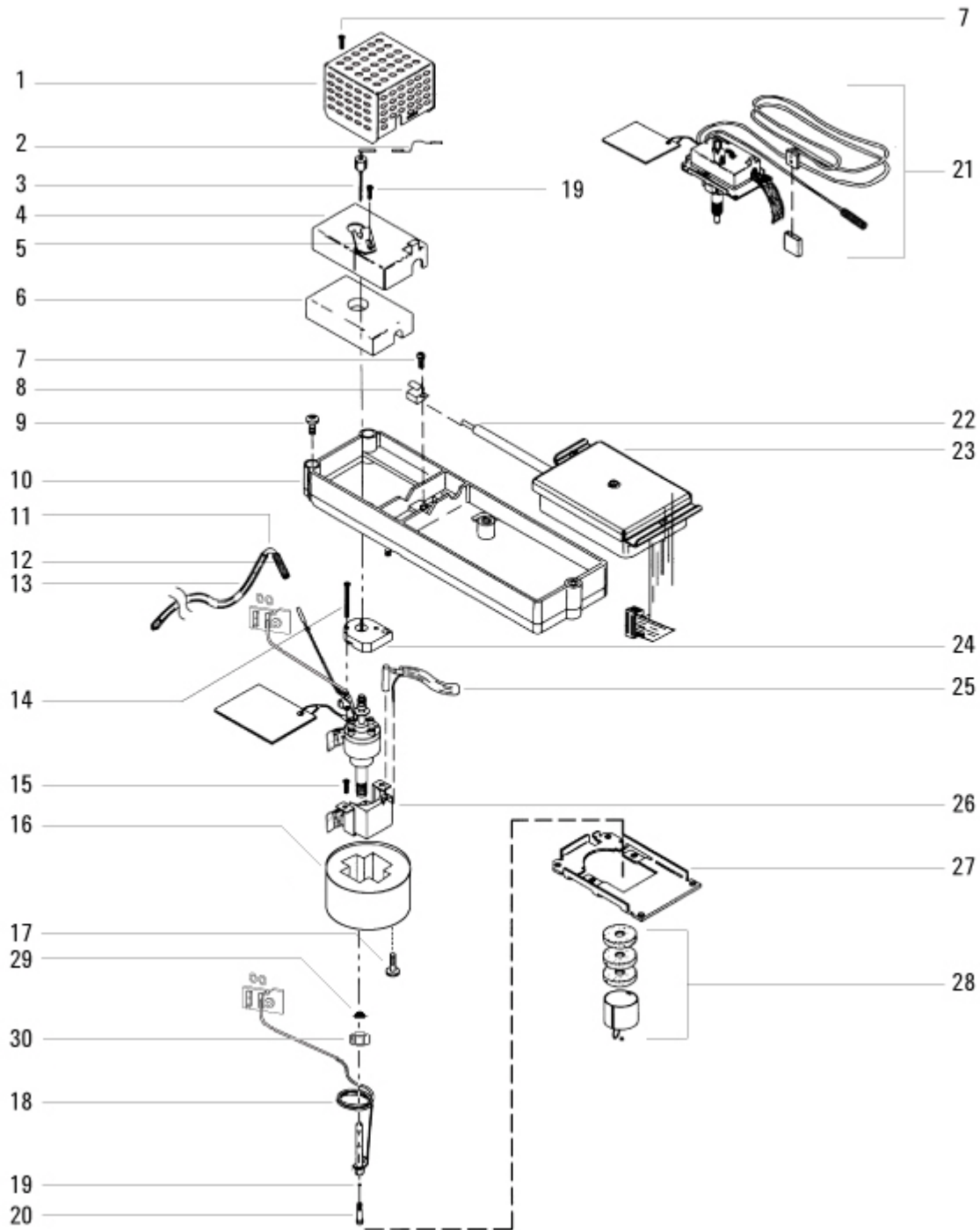


**uECD Side Mount (3rd Detector) Parts**

This section lists parts specific to the  $\mu$ ECD mounted in the side carrier as a third detector. For other service parts, see “[Microcell Electron Capture Detector \(uECD\) parts](#)”.

- G2397AD Micro ECD wit EPC
- G2398AD Micro ECD with EPC (**Japan only**)
- G3485B Manual Pneumatics Carrier for 7890B
- G3440-68011 Mounting bracket kit

Item	Description	Part number	Qty
1	Third Detector Heated Zone Cable	G3432-60506	1
2	Third Detector EPC Communication Cable	G3432-60505	1
3	1/8-inch Bulkhead union, 316 SST	0100-0133	1
4	aT-Manifold Assy	G1176-60569	1
5	Screw, machine, M4 x 12 mm long	0515-2496	7
6	Screw, machine, M4 x 25 mm long	0515-2712	2
7	Bulkhead fitting retainer, 316 SST	0100-0137	1
8	Nutwarmer cup assembly	19234-60700	1
9	Label, ECD Caution	19233-90730	1
NS	Nut Warmer Insulation	19234-60715	1

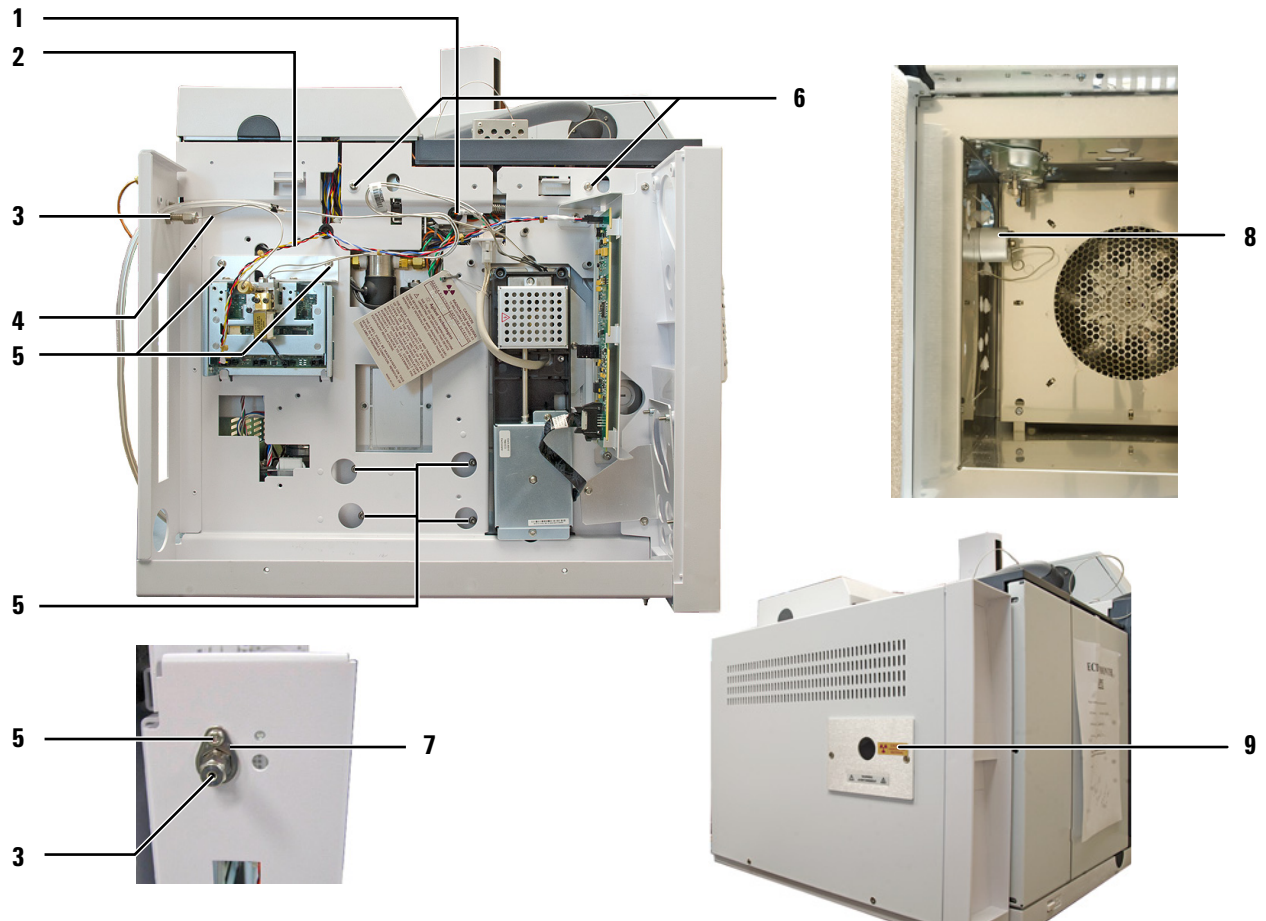


## Nitrogen Phosphorus Detector (NPD)

Top level subassemblies for NPD detector accessories:

- G3495A NPD (white ceramic bead) for packed and capillary columns
- G3496A NPD (white ceramic bead) optimized for capillary columns only
- G3463A NPD (Blos bead) for packed and capillary columns
- G3464A NPD (Blos bead) optimized for capillary columns only

Description	Part number
NPD Detector Module, adaptable, without ship kit, includes bead G1534-60670 (white)	G3434-60506
NPD Detector Module, capillary, without ship kit	G3434-60501
NPD Adaptable Weldment Assembly with pallet, electrometer, and base	G3434-80508
NPD Capillary Weldment Assembly with pallet, electrometer, and base	G3434-80507
NPD EPC Module	G3434-60502
NPD Signal Board	G3434-60021

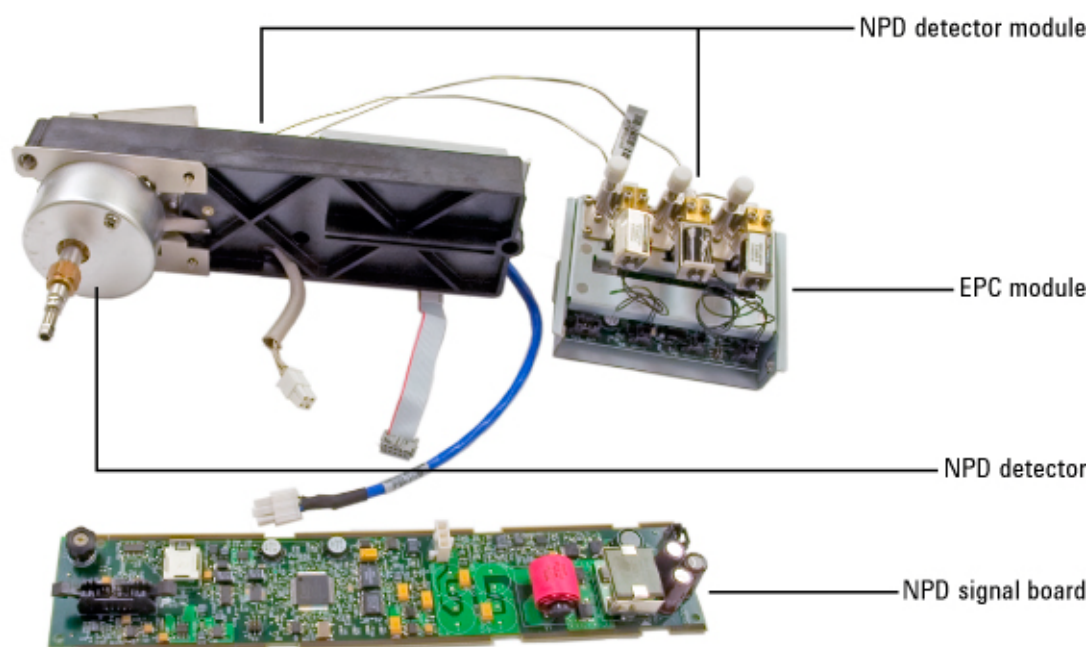




## Nitrogen Phosphorus Detector (NPD) parts

Item	Description	Part number	Qty
1	Screw, M4 × 10 mm, Torx T-20	0515-2495	1
2	Hinged cover assembly	G1534-80521	1
3	Screw, M3 × 8 mm, Torx T-10	0515-0655	3
4	NPD bead assembly	See Parts Finder	1
5	Screw, M4, Torx T-20, 10 mm long	0515-2495	3
6	Power cable assembly	G3434-60600	1
7	Lid weldment	G1534-80510	1
8	Collector funnel, standard	G1534-20530	1
9	NPD ceramic replacement kit Includes top and bottom C-rings and upper and lower insulators.	5182-9722	—
10	Screw, M4 × 10 mm long	0515-2495	1
11	J-clamp	1400-0015	1
12	Screw, M4 × 10 mm long	0515-2495	2
13	NPD interconnect assembly	G1534-60610	1
NS	Spring, interconnect	1460-2142	1
14	NPD electrometer	G3434-60010	1
15	Pallet captive screws, M4 x 20 mm long	1390-1024	4
16	O-ring, size 2-006, fluorocarbon, 12/pk	5180-4181	3
17	Screw, captive, M3, T-10	G1946-20168	1
18	Mounting pallet	G1531-40020	1
19	Jets	See Parts Finder	1
20	Base weldment with jet - Capillary column NPD - Packed column NPD	G3434-80505 G3434-80506	1
21	Lid stop/standoff	G1534-20590	3
NS	Insulation cap	G1531-00130	1
22	NPD Insulation cup kit. Includes:	G3434-00004 G3434-00001 G3434-00002	1
23	Heater/sensor assembly	G1530-67950	1
24	Column adapters for packed NPD	See Parts Finder	
25	Cup insulation kit. Includes: Insulation Nutwarmer cup assembly	19234-60720 19234-60715 19234-60700	1
NS	Flow measurement adapter	G1534-60640	—

## 18 Illustrated Parts Breakdown

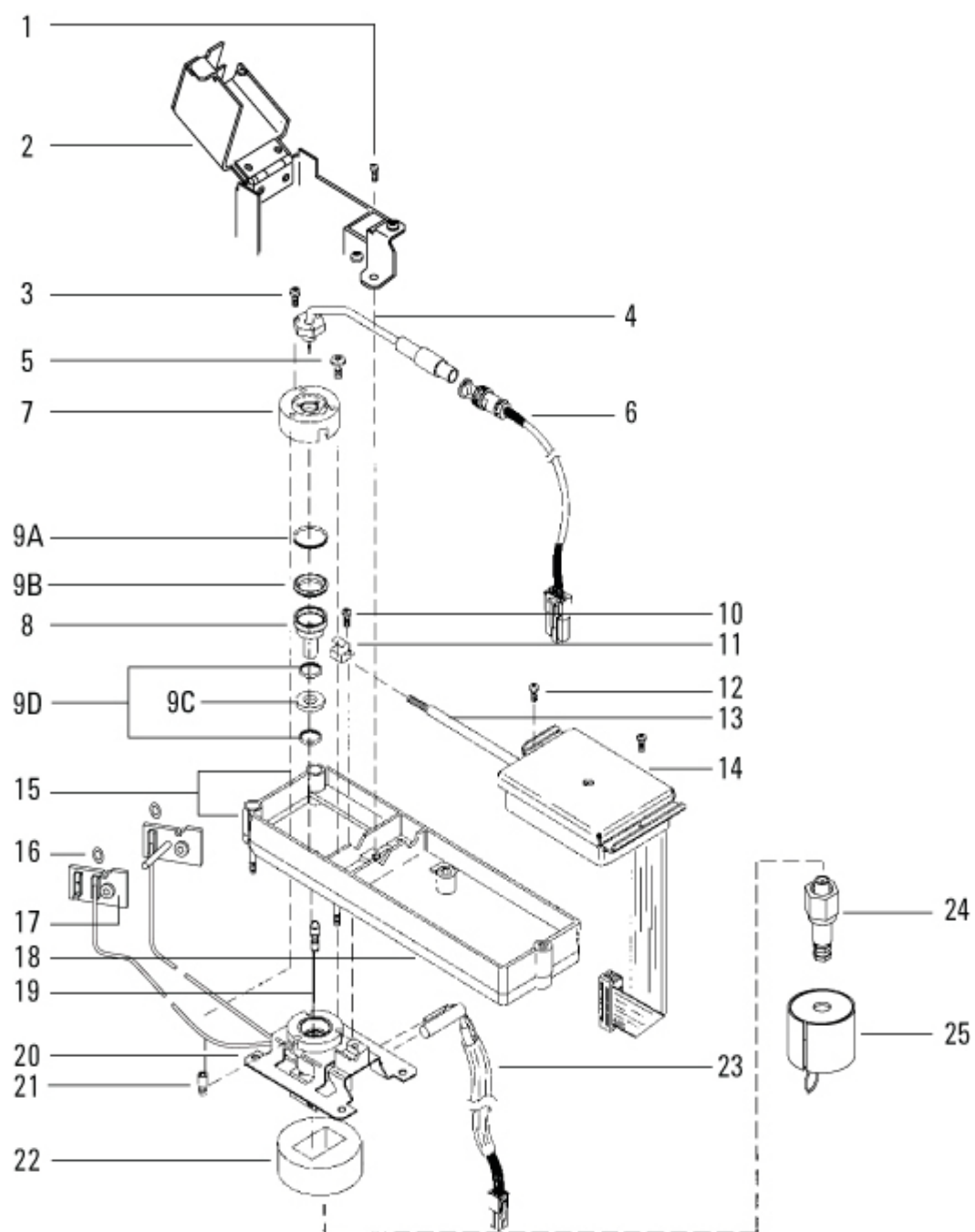

## Thermal Conductivity Detector (TCD)

Top level subassemblies for TCD detector accessories:

- G3432A TCD Accessory with EPC, Kit
- G3437A TCD mounted on left-hand side of GC, Kit, **for 7890A GC**
- G3437B TCD with EPC mounted on left side of **7890B GC**
- 19232C TCD to FID adapter kit


Description	Part no.
TCD Detector Module without ship kit (includes TCD cell, heater and sensors, EPC module without switching valve, signal board PCA, insulation cup assembly, and mounting pallet)	G3432-60500
TCD base assembly, tested (includes filament block assembly, delta PRT, heater/sensor, TCD switching valve, detector baseline noise results)	G3432-60220
Acid tested TCD (includes filament block assembly and TCD switching valve). <b>Use only for specific applications</b> , such as fatty acid analysis.	G3432-60509
TCD EPC Module (without switching valve)	G3432-60532
TCD Signal Board	G3432-60015 *

\*. Replaces G3432-60010, see service note G3440A-40.


### Selecting a replacement TCD cell assembly

When replacing a TCD cell, normally order the TCD base assembly (G3432-60220). This is the standard replacement part for most users.


Customers performing analyses of certain acidic compounds, for example, fatty acid analysis, may need the Acid Tested TCD (G3432-60509). The Acid Tested TCD costs about twice as much as a standard cell, and should only be ordered for the specific applications that require it.

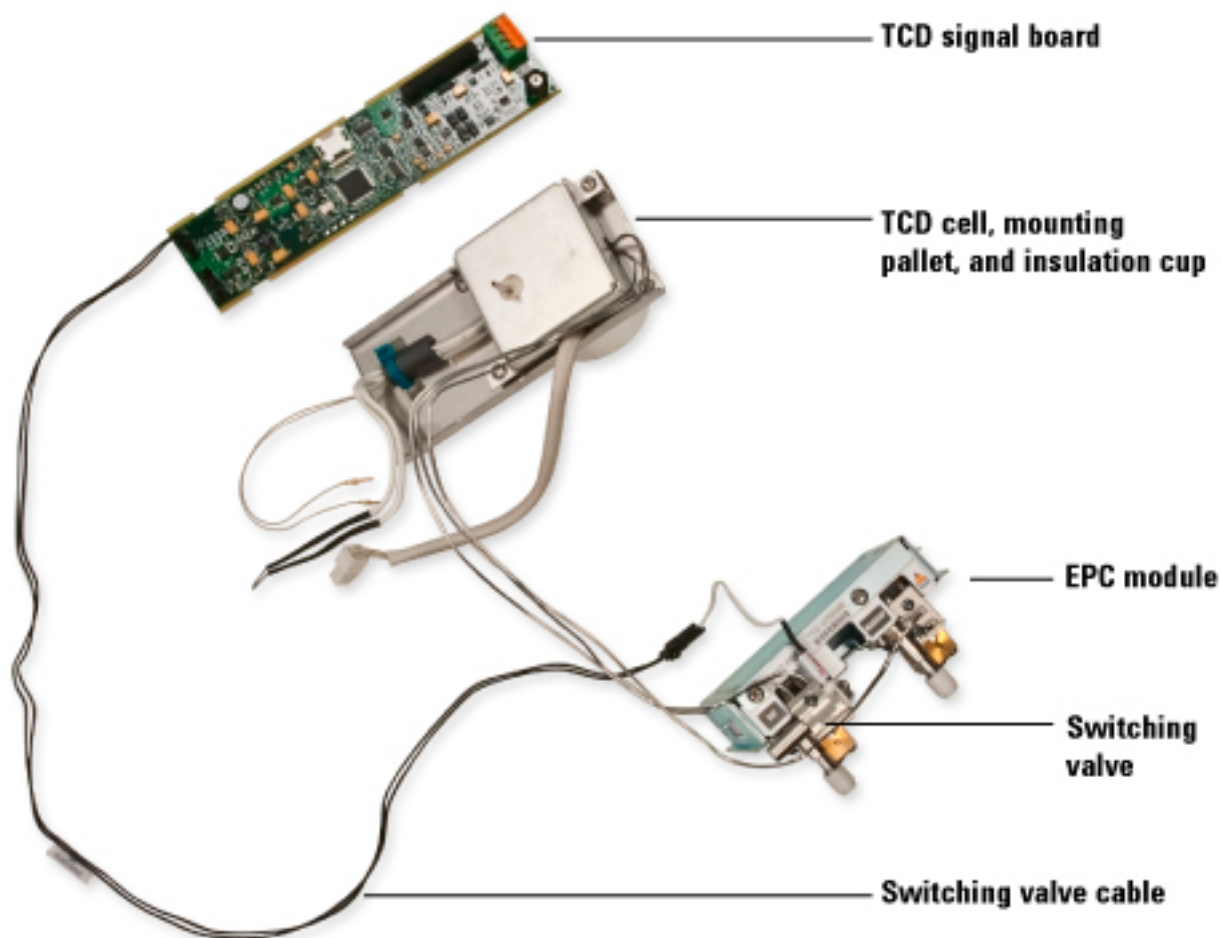

**Figure 19** TCD major component identification (top-mounted TCD, G3437B, shown)

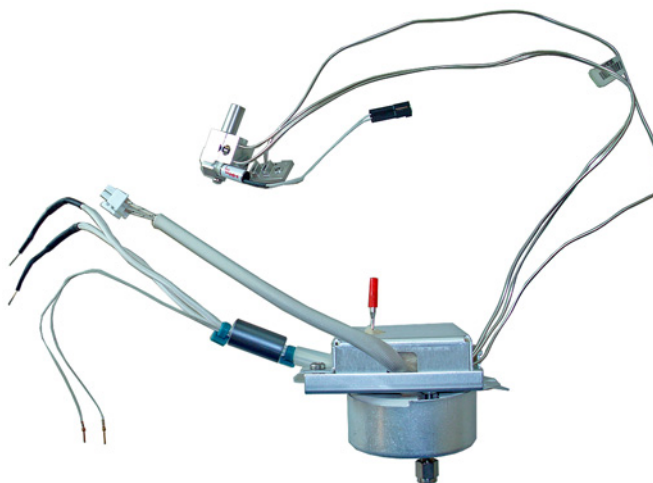


**Figure 20** G3432-60220, TCD base assembly

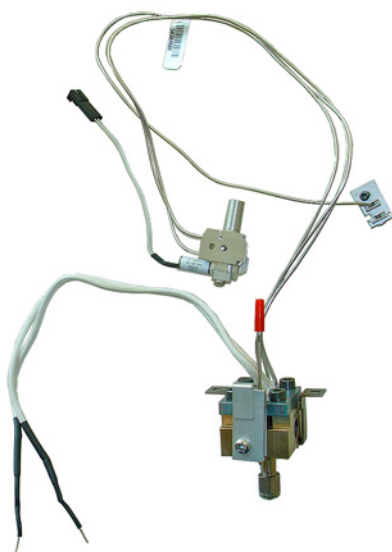

**Figure 21** G3432-60500, Acid tested TCD


## Thermal Conductivity Detector (TCD) parts

Item	Description	Part number	Qty
1	Screw, M4 × 12 mm, T-20	0515-2496	2
2	Thermal cover and insulation assembly	G1532-00027	1
3	O-ring, 12/pk (2 for EPC module seal, 4 for valve seals)	5180-4181	6
4	Screw, captive, M3, T-10	G1946-20168	2
5	Screw, M4 × 10 mm, Torx T-20	0515-2711	2
6	TCD mounting pallet (top mount)	G3432-00002	1
7	TCD Insulation Cup Kit	G3432-00004	1
8	Screw, M4 × 40 mm, T-20	0515-4793	1
9	TCD seal insulation	G1532-00080	1
10	Column adapters, ferrules, and nut	See Parts Finder	
11	TCD cell assembly	See “Thermal Conductivity Detector (TCD)” on page 662.	
12	PRT Sensor, TCD (Delta PRT)	G1532-60660	1
13	Heater/Sensor assembly	G1530-67950	1
14	Cylindrical EMI suppressor, ferrite	9170-1730	1
15	Screw, M4 × 12 mm, T-20 (attach detector to oven top)	0515-2496	2
NS	TCD vent restrictor kit (See <a href="#">Table 44</a> and service note G1530-25.)	G1532-60700	

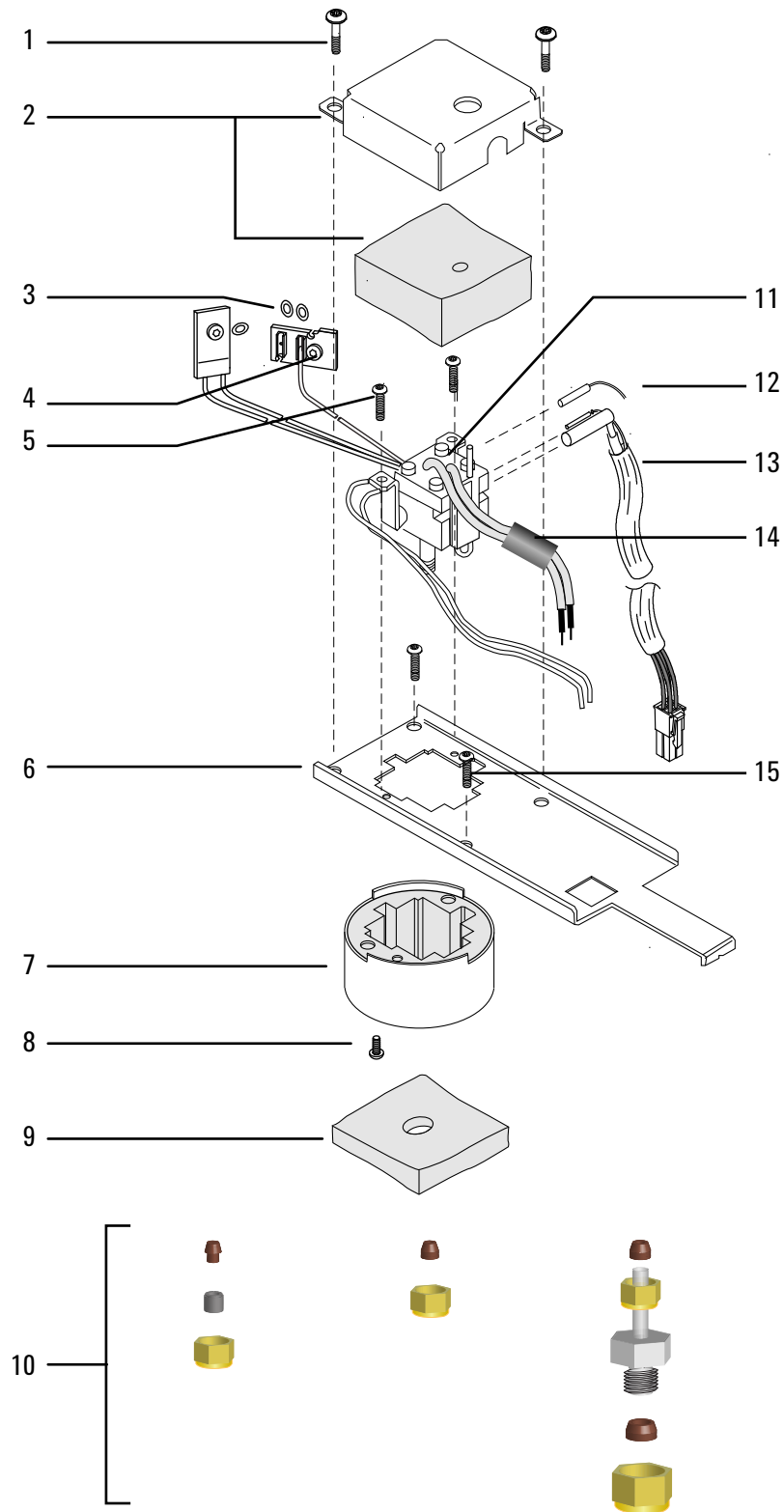

**Reference gas switching valve** This valve is located on the EPC module.

Item	Description	Part number	Qty.
1	Clamp, switching valve	G3432-20523	1
2	Valve, TCD Switching	G1532-60620	1
3	Screw, M3 x 12 mm T-10	0515-1084	1
4	Cable, switching valve	G1532-60550	1

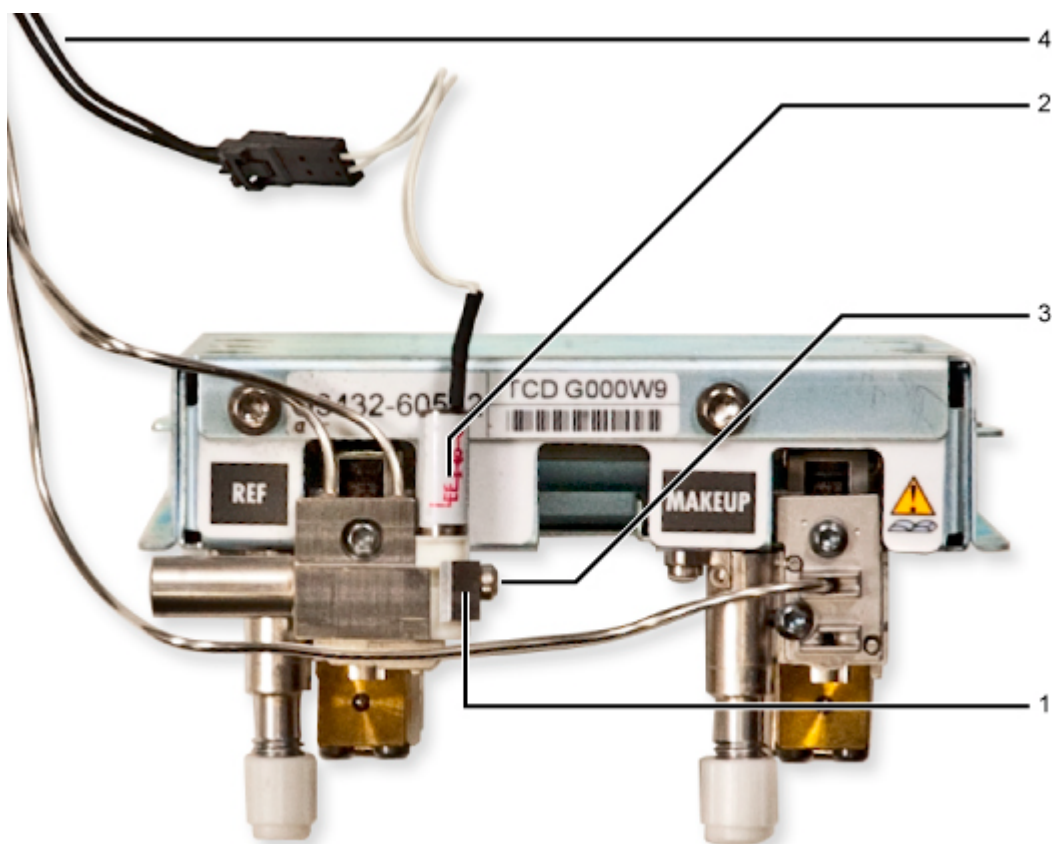



**TCD vent restrictor kit** Use this kit to eliminate baseline signal spikes due to changes in laboratory air pressure due to doors closing and similar events.

**Table 44** TCD outlet restrictor kit G1532-60700 parts

Item	Description	Part number	Qty.
1	Swagelok reducer	0100-0121	1
2	Restrictor	19232-60555	1
3	1/8-inch Graphite/Vespel ferrule (10/pk)	0100-1107	1

## EPC Module Parts

Lookup table [669](#)  
Inlet modules [670](#)  
Detector modules [671](#)  
PCM modules [672](#)  
Aux EPCs [673](#)

This section provides a quick reference to the EPC-related part numbers.

Each proportional valve requires 2 O-rings (5180-4181, 12/pk). The valve part numbers are for kits, which contain the valve, screws, and a package of O-rings.

### Lookup table

Each proportional valve has a part number on the valve. This table lists the reorder part number for each proportional valve.

Part number on valve	Reorder
G3431-60508	G3430-67013
G3431-60509	G3430-67014
G3430-60527	G3430-67015
G3430-60528	G3430-67016

**Inlet modules****EPC modules**

Description	Part number	O-rings	Slot
COC	G3454-60554	3	1 or 2
MMI	G3452-60552	3	1 or 2
PP	G3451-60551	3	1 or 2
PTV	G3500-60500	3	1 or 2
S/SL 100 PSI	G3452-60552	3	1 or 2
S/SL 150 PSI	G3452-60510	3	1 or 2
VI	G3504-60501	4	1 or 2

**Proportional valves**

Description	Carrier	Septum purge	Split vent
COC	G3430-67013	G3430-67014	
MMI	G3430-67013	G3430-67014	G3430-67015
PP	G3430-67013	G3430-67014	
PTV	G3430-67013	G3430-67014	G3430-67015
S/SL 100 PSI	G3430-67013	G3430-67014	G3430-67015
S/SL 150 PSI	G3430-67013	G3430-67014	G3430-67015
VI	G3430-67013	G3430-67014	G3430-67015

## Detector modules

### EPC modules

Description	Part number	O-rings	Slot
uECD	G3433-60533	1	3 or 4
FID	G3431-60531	3	3 or 4
FPD	G3435-60535	3	3 or 4
FPD <sup>+</sup>	G3435-60335	3	3 or 4
NPD	G3434-60502	3	3 or 4
TCD	G3432-60532	2	3 or 4 or AUX Det 2

### Proportional valves

Description	Air	Makeup	Hydrogen	Reference
uECD		G3430-67013		
FID	G3430-67013	G3430-67014	G3430-67014	
FPD	G3430-67013	G3430-67013	G3430-67013	
FPD <sup>+</sup>	G3430-67013	G3430-67013	G3430-67013	
NPD	G3430-67013	G3430-67014	G3430-67014	
TCD		G3430-67013		G3430-67013

## PCM modules

### EPC module

Description	Part number	O-rings	Slot
PCM	G3476-60501	3	1, 2, 5, or 6

### Proportional valves

Description	Carrier	AUX
PCM	G3430-67013	G3430-67016

## Aux EPCs

### EPC module

Description	Part number	Slot
AUX	G3470-60501	5 or 6 or AUX

### Tubing weldment

Description	Part number
Pneumatic tubing with block connector and captive screw	G3430-60550



**Figure 22** Tubing weldment

### Proportional valves


Description	Channel	Channel	Channel
AUX	G3430-67013	G3430-67013	G3430-67013

## Nickel Catalyst Accessory

Top level subassemblies for FID detector accessories:

- G3478A Nickel Catalyst Kit
- G3440-60510 Nickel Catalyst Kit without ship kit

Item	Description	Part number	Qty
	Nickel catalyst assembly (items 2 through 8)	G3440-63002	1
1	Nickel catalyst hydrogen mix weldment	G3440-81000	1
NS	Screw,		
2	Screw, M4 x 12 mm, captive	0515-2711	3
3	Top cover plate	G1543-00085	1
4	Heater/sensor assembly	G1580-61160	1
5	Top insulation	G1543-00100	1
6	Nut, hex, with lockwasher	0535-0043	2
7	Heater block	18900-20835	1
8	Screw, socket M4 x 20 mm	0515-0038	2
9	PP base insulation	G1543-00030	1
10	PP bottom insulation	G1543-00070	1
11	Adapter insulation	19234-60715	3
12	Insulating cup	19234-60700	1

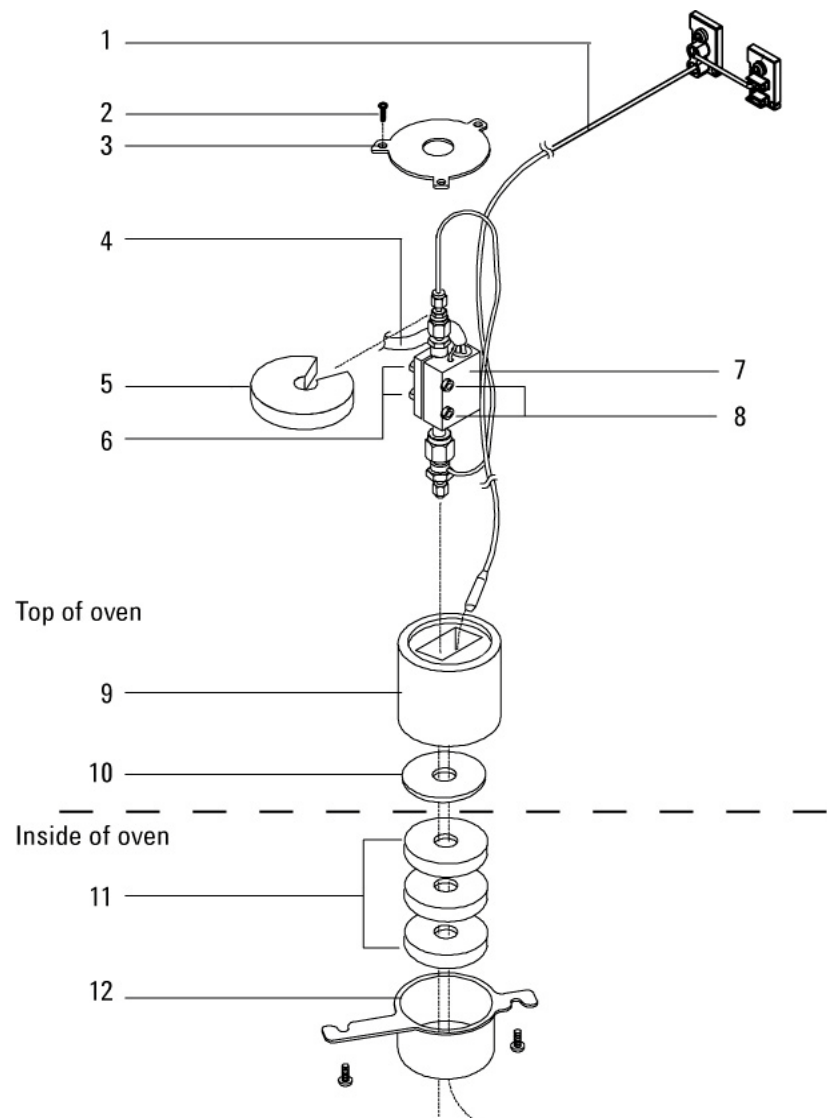






## Oven

This section contains illustrated parts breakdowns for the following 7890A GC components:

- Oven assembly
- Oven flapper assembly
- CO2 cryogenic cooling
- Liquid nitrogen (LN2) cryogenic cooling


## Oven assembly

Item	Description	Part number	Qty
1	Door hinge rod	G3430-20518	1
2	Oven door (non-removable)	G3430-60502	1
3	Oven heater shroud assemblies:		1
	–120V USA power	G1530-61610	
	–200V Japan power	G1530-61620	
	–220V Single phase power, Europe, and –208V USA power, single phase	G1530-61630	
	–220V Single phase power, Hong Kong	G1530-61630	
	–220V Single phase power, China	G1530-61630	
	–220V China power, slow ramp	G1530-61230	
	–220V Israel power	G1530-61630	
	–230V Single phase power, Switzerland, 10 amp	G1530-61670	
	–230V Denmark power, 10 amp	G1530-61670	
	–230V Denmark/Switzerland power, 16 amp	G1530-61650	
	–230V Continental Europe power	G1530-61650	
	–240V Australia power	G1530-61640	
	–240V South Africa/India power	G1530-61640	
	–240V Great Britain/Ireland power	G1530-61640	
	–240V USA power	G1530-61640	
	–240V Australia power, slow ramp	G1530-61640	
4	Screw, M4, Torx T-20, chrome plated	0515-2711	4
5	On/Off switch rod	G3430-40023	1
6	Hex nut, M4 x .7 3.2 mm	0535-0043	3
7	Motor, oven fan	G3430-60504	1
8	Washer	2190-0712	6
9	Door sensor assembly	G1530-80670	1
10	Oven door button	G3430-80513	1
11	Column clip kit (includes 4 spring clips and screws)	G1530-61580	1
12	Oven sensor assembly	G1530-61030	1
13	Column hanger	1460-1914	
14	Fan blade, stainless steel, and set screw	05890-80270	1

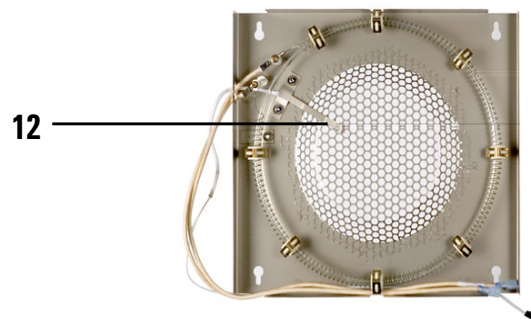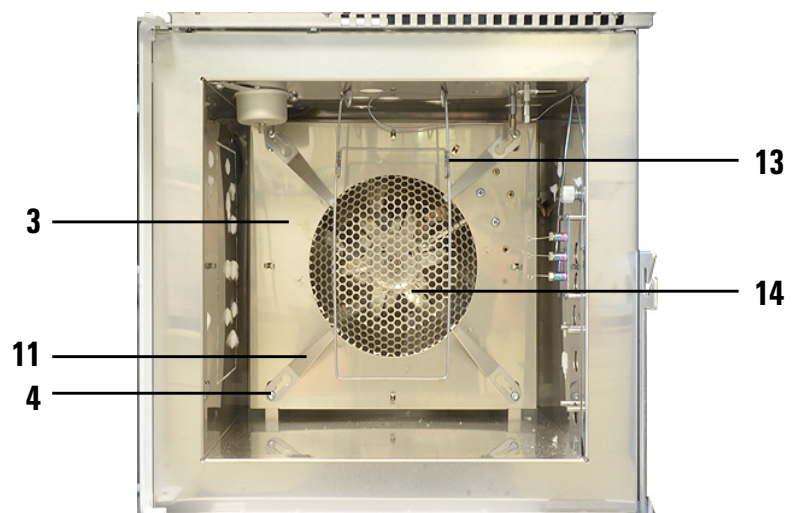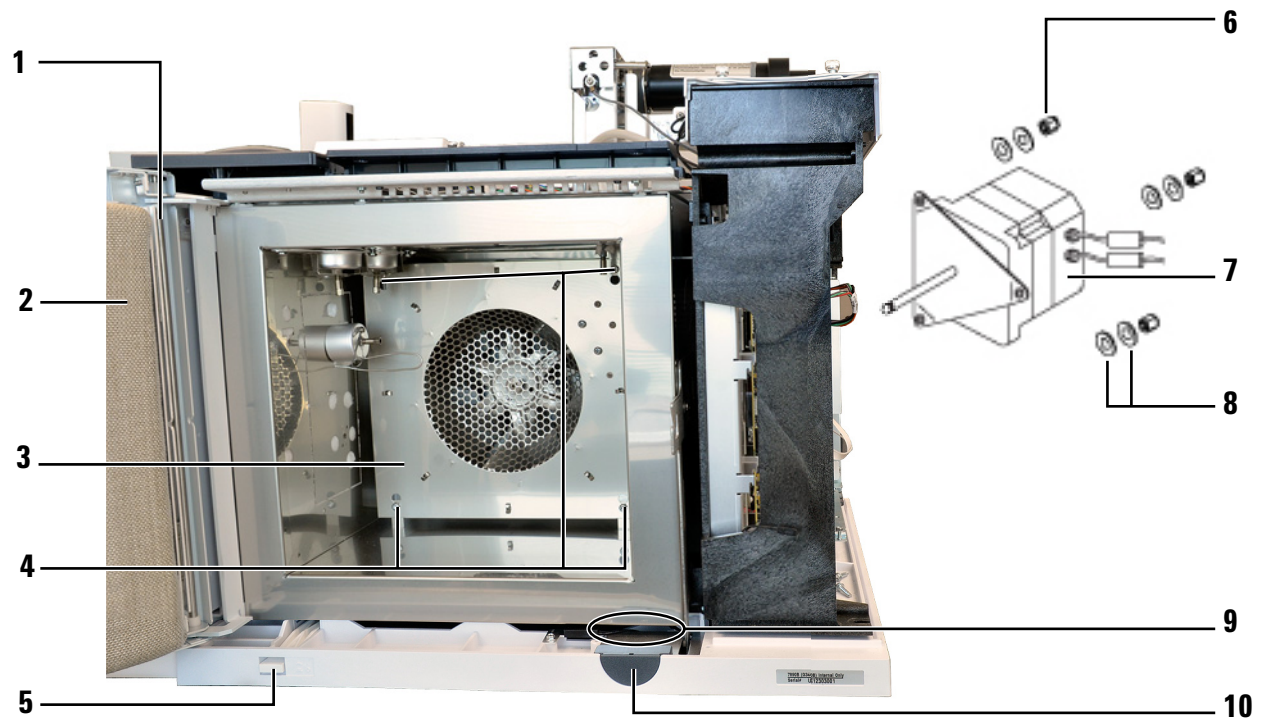

Removable oven door kit parts


Item	Description	Part number	Qty
NS	Removable oven door kit (includes screw, door pin, hook plate, base plate, door, and support spacer)	G3430-68640	1
1	Rod	G3430-20643	1
2	Screw: Mach M4 x 0.7, 20 mm, LG T-20 Pan HD SST	0515-0686	1
3	Lower Oven Flange Support Spacer	G3430-20645	1
4	Hook Plate	G3430-00647	1
5	Base Plate	G3430-00646	1
NS	Door, removable	G3430-60640	1

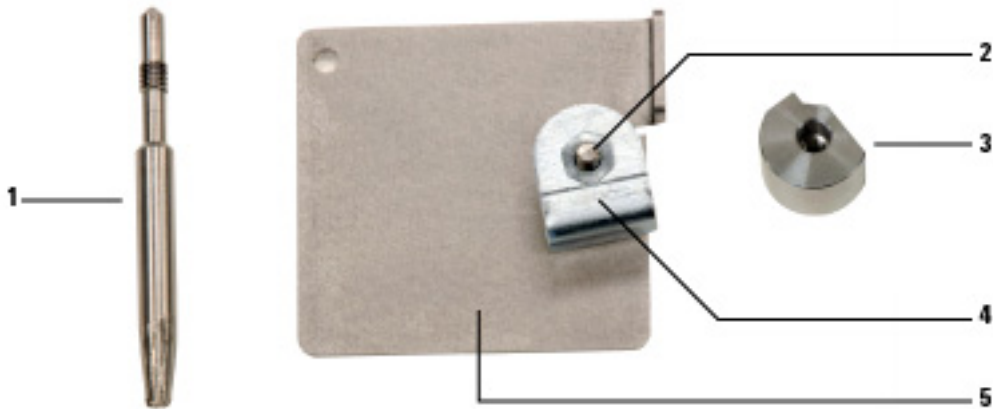

## Oven flapper assembly


Item	Description	Part number	Qty
NS	Oven exhaust duct and flapper assembly (includes 3 captured screws 1390-1023)	G3430-60008	1
1	Dual duct assembly	G3430-81001	1
2	Screw, Torx T-10, M3 × 12 mm (to attach duct to oven)	1390-1023	3
3	Screw, Torx T-20, M4 x 12 mm (to attach motor to shaft)	0515-2496	1
4	Screw, Torx T-10, M3 × 8 mm (to attach motor to duct)	0515-0655	2
5	Flapper bushing	G1530-20550	1
6	Flapper shaft assembly	G3430--81002	1
7	Stepper motor assembly	G1530-60945	1

## 18 Illustrated Parts Breakdown

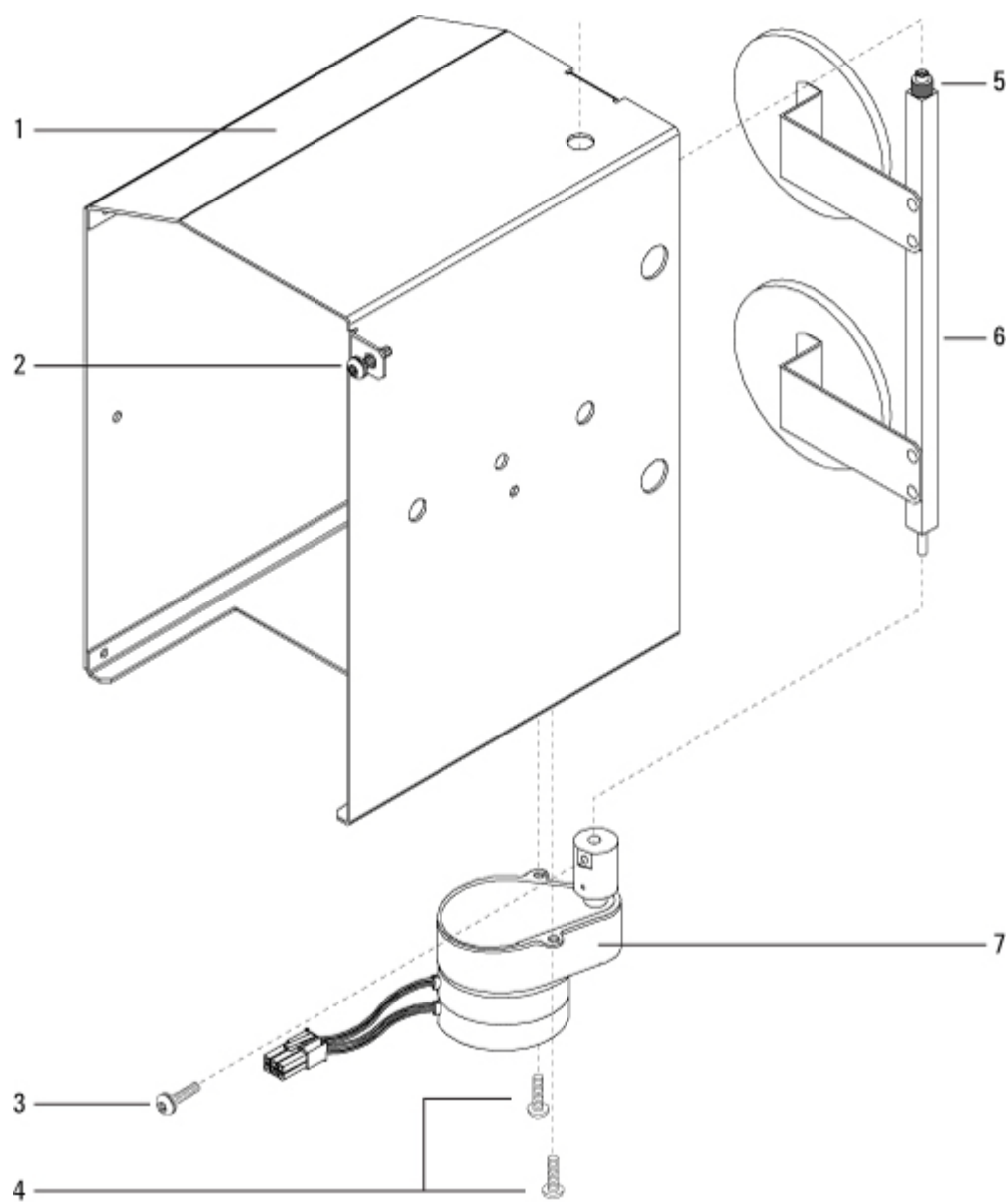

**CO2 cryogenic cooling**


Item	Description	Part number	Qty
1	Tee, brass, 1/8-inch	0100-0090	1
2	CO2 blast brazement (Cryoblast)	G1565-80550	1
3	Cryo bracket	G1565-00010	1
4	CO2cryo valve assembly (recommended level for valve)	G1565-65510	1
5	1/8-inch male connector without O-ring	G1543-80025	1
6	CO2 inlet tube	G1565-20600	1
7	CO2 cryo inline filter	3150-0602	1
8	Cryoblast restrictor weldment	G1565-80590	1
9	CO2 weldment (oven cryo)	G1565-80505	1
10	Clamp screw	G1565-20560	1
11	Gland seal, CO2	G1565-20590	1
12	Screw, 8-32, .312-inch	2510-0043	2
NS	Screw Torx T-20, M4 × 12 mm	0515-2496	2
NS	PTFE seal	G1565-20840	1

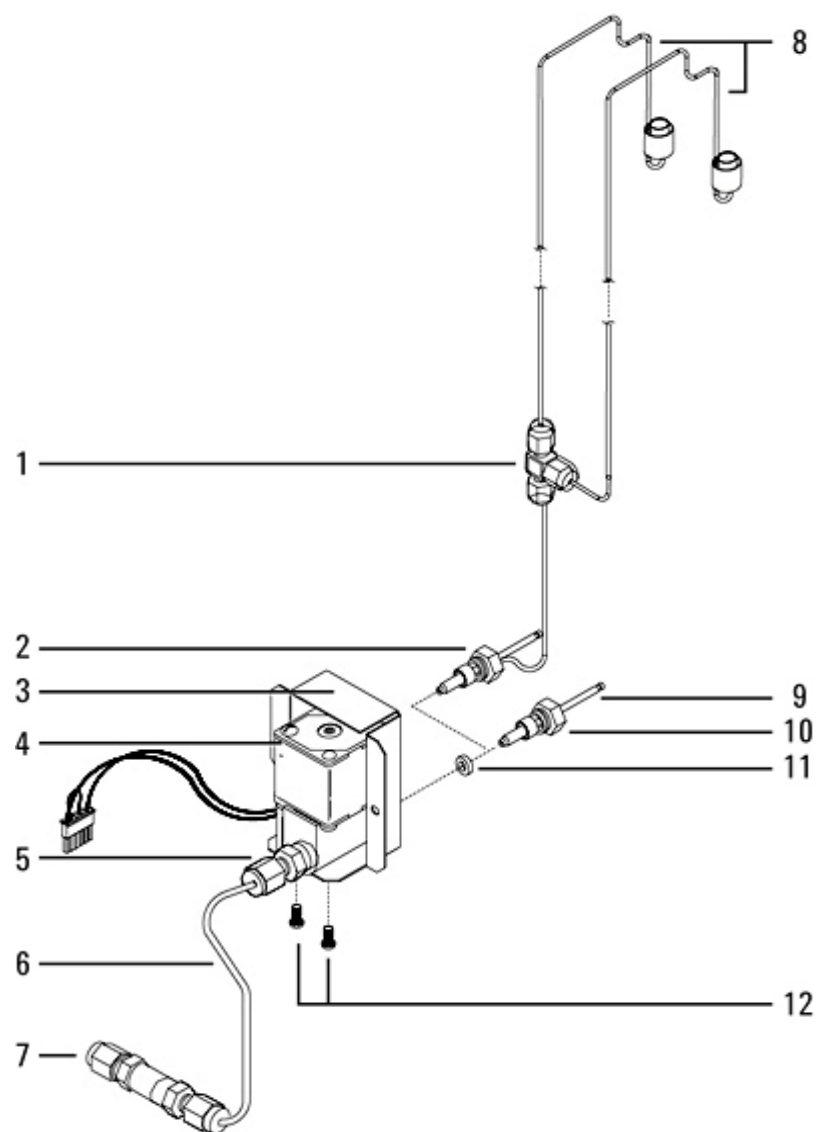



## Liquid nitrogen (LN<sub>2</sub>) cryogenic cooling

Item	Description	Part number	Qty
1	N2 blast brazement	G1566-80535	1
2	Tee, 1/8-inch brass union, 2/pk	5180-4160	1
3	Connector, male, 1/4-inch	0100-0208	1
4	Screw, 8-32, .312-inch	2510-0043	2
5	Cryo bracket	G1565-00010	1
6	N2 restrictor tube	G1566-20575	2
7	N2 cryo valve	G3430-60018	1
8	Liquid N2 nozzle	19310-20500	1
NS	Screw Torx T-20, M4 × 12 mm	0515-2496	2
NS	Pipe thread sealant tape, 1/4-inch wide	0460-1266	
NS	N2 cryo valve/dual blast assembly	G1566-65508	
NS	N2 cryogenic kit	G1566-65517	

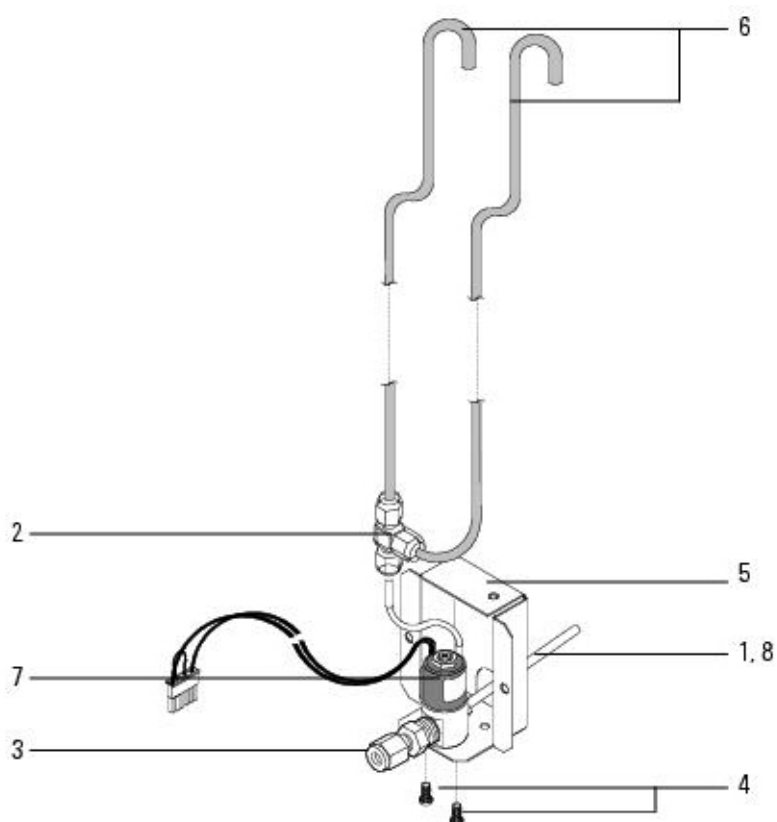

## Valves

Top level subassemblies for Valve accessories:


- G1580A Valve box for one valve
- G1581A Valve box for two valves

This section contains illustrated parts breakdowns for the following 7890 Series GC valves and related components. In general, the valve parts are the same for both 7890 Series and 6890 Series GCs. The solenoid valve stacks have changed, however. Be sure to order the correct version of parts or upgrade the GC to the newer version if needed.

- [Valve box assembly](#)
- [Valve driver assembly](#)
- [Valve actuator assembly \(1 of 2\)](#)
- [Valve actuator assembly \(2 of 2\)](#)
- Valco W-Series minivalve (See Parts Finder.)
- [Miscellaneous Valve Parts](#)

## Valve box assembly

Item	Description	Part number	Qty
1	Screw, M4 × 12 mm, Torx T-20, chromeplated	0515-2711	6
2	Valve box cover	G1580-00030	1
3	Valve box insulation, top	G1580-00050	1
4	Insulation retainer plate	G1580-00040	1
5	Hex nut, insulation plate		2
6	Screw, M3 × 30 mm, Torx T-10, chromeplated	0515-2525	4
7	Heater block	G1580-20520	1 or 2
8	Standoff, valve box	G1580-20500	2
9	Valve box bottom plate	G1580-00010	1
10	Screw, M3 × 8 mm, Torx T-10, chromeplated	0515-2726	2
11	Heater/Sensor harness cable clamp	1400-0015	2
12	Heater/Sensor assembly	G1580-61140	2
NS	Cable-tie strap	1400-0249	4
NS	Aluminum tube (split/splitless inlet only)	18900-20320	1
NS	In-line filter (sample in line)	0101-0532	1

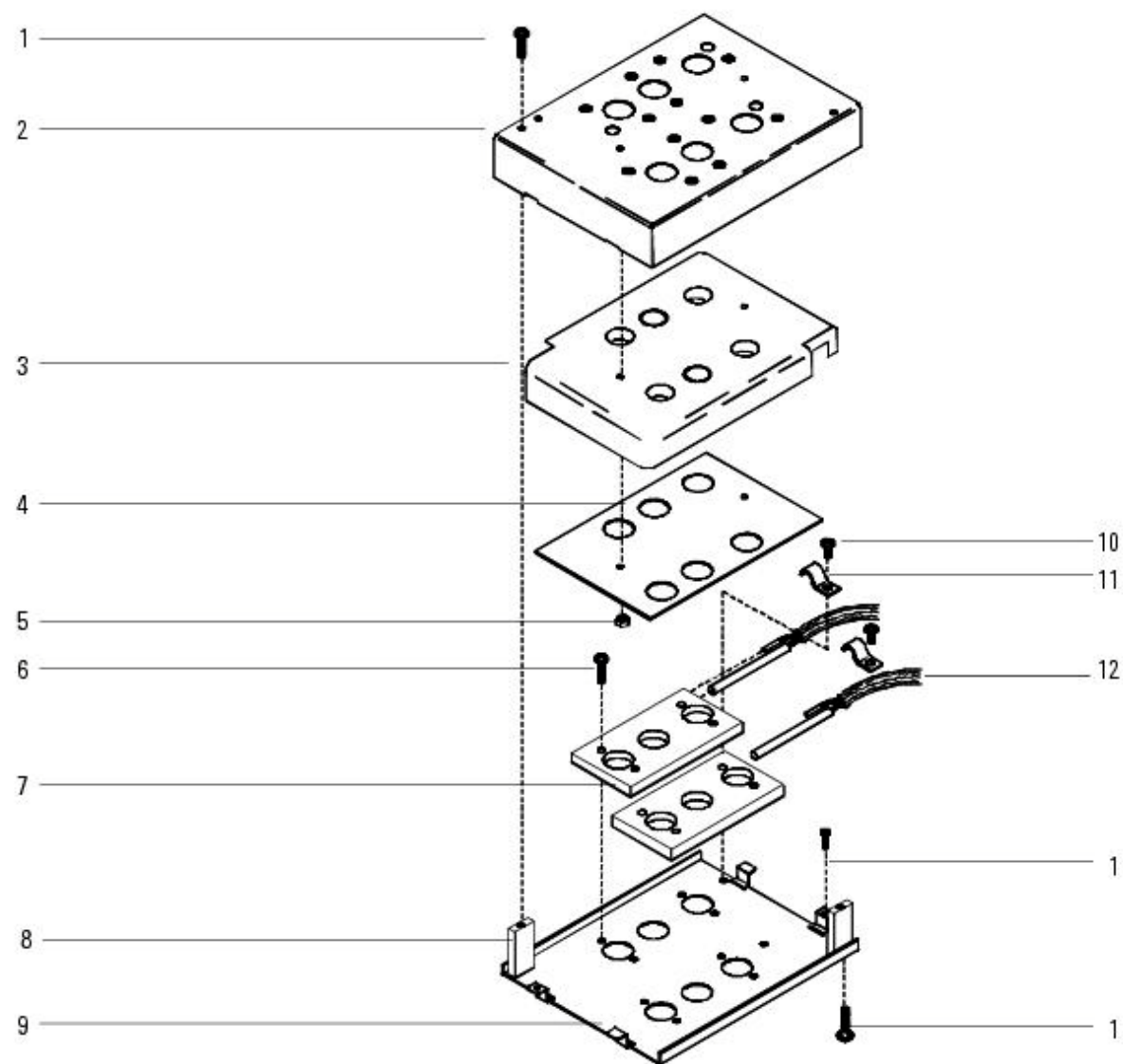

## Valve driver assembly


**For 7890A:** If your GC has serial number CN10201096 or greater, or serial number US10281011 or greater, use the new-style valve driver assembly parts.

If your GC has a serial number less than listed above, use the old-style valve driver assembly parts (if available).

**For 7890B:** All 7890B GCs use the new style valve driver assembly parts.

**Valve driver assembly (original style)**

Item	Description	Part number	Qty
1	Valve driver bracket	G1580-00070	1
2	Captured screw, M4 × 12 mm, Torx T-20	1390-1023	2
3	Valve driver wiring harness	G1530-60660	1
4	PTFE tubing, 5 m, 1.5 id, 3 mm od	5062-2483	8
5	Adapter fitting	5185-8386	8
6	Standoff (included with valves)		8
7	4-way Solenoid valves	G1580-61095	4
8	O-ring (with valves and end plate kit)		10
9	Solenoid valve end plate kit	05890-61097	1
10	Screw, socket head, M5 x 16 mm long	0515-1214	2
11	Elbow fitting, 1/4-inch, male	0100-1632	2
12	Exhaust tubing, 1/4-inch od, 120 inch	0890-1489	1
NS	Heater cable assembly	G1530-60790	1


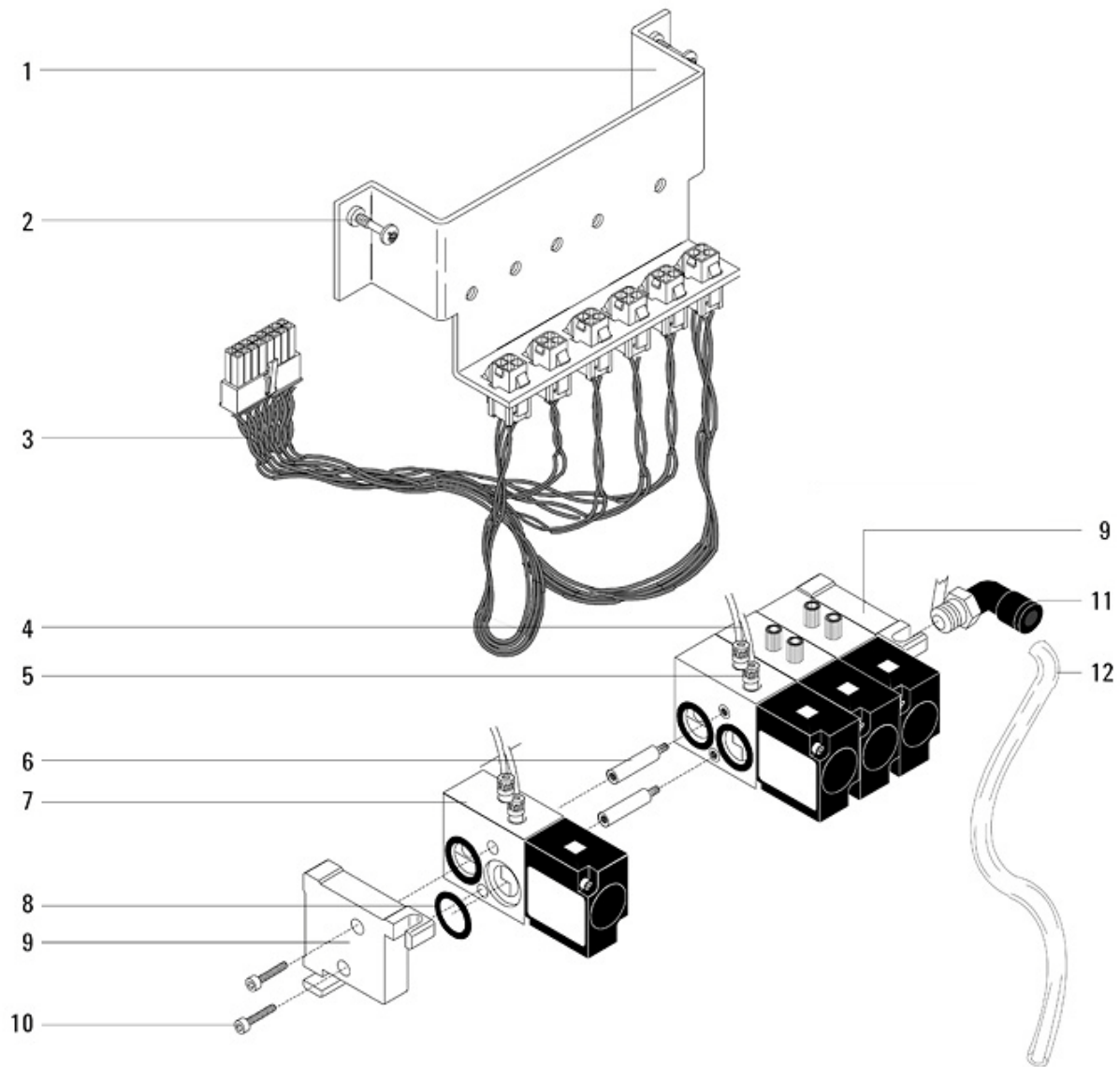


**Valve driver assembly (new style)**


Item	Description	Part number	Qty
1	Valve driver bracket	G1580-00071	1
2	Captured screw, M4 × 12 mm, Torx T-20		2
3	Valve driver wiring harness	G1530-60660	1
4	PTFE tubing, 5 m, 1.5 mm id, 3 mm od	5062-2483	1
5	Exhaust tubing, 1/4-inch od, 120 inch	0890-1489	1
NS	Heater cable assembly	G1530-60790	1
6	Adapter fitting	5185-8386	2
7	Airtec Solenoid Valve w screws	G1580-68000	1
NS	Screws for One Valve Stack (M4 x 30)	5185-8381	2
8	Screws for Two Valve Stack (M4 x 50)	5185-8382	2
9	Screws for Three Valve Stack (M4 x 70)	5185-8383	2
10	Screws for Four Valve Stack (M4 x 85)	5185-8384	2
NS	One Valve Stack*	G1580-65111	
NS	Two Valve Stack*	G1580-65112	
NS	Three Valve Stack*	G1580-65113	
11	Four Valve Stack*	G1580-65114	
12	PT 1/8-inch male elbow fitting for tube 1/4-inch	5185-8385	1
13	Solenoid Valve End Plate Kit	G1580-61097	1

- Includes bracket, solenoid(s), fittings, end plates, and harness cable.



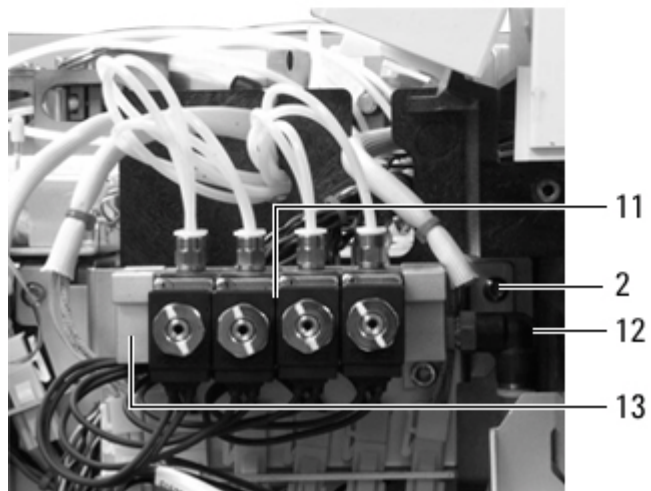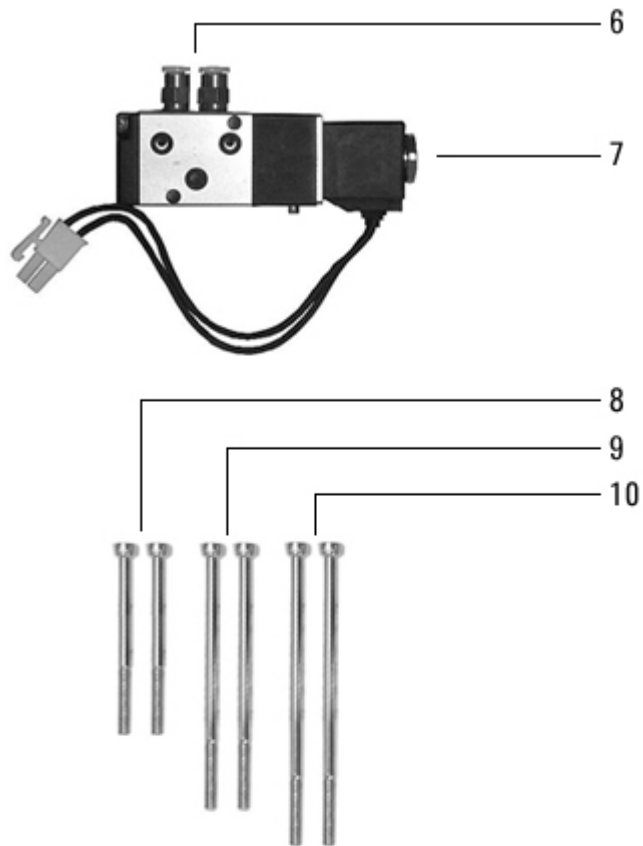


## 18 Illustrated Parts Breakdown

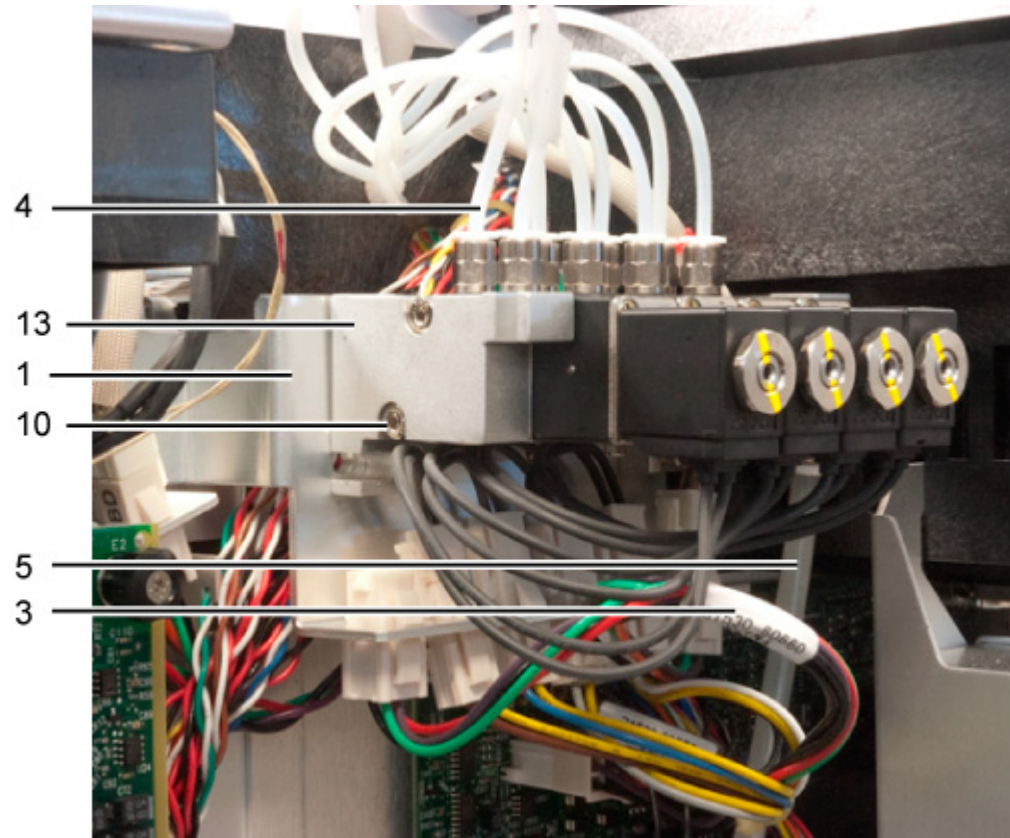

**Valve actuator assembly (1 of 2)**

Item	Description	Part number	Qty
1	Valve actuator assembly	19325-60660	1
2	Modified screw	19325-80030	2
3	Quick-release pin		1
4	Hose fitting	0100-1707	1
5	PTFE tubing 5m, 1.5mm ID, 3mm OD	5062-2483	3
6	36° Actuator limiter	18900-21000	1
7	End cap		1
8	Elbow fitting	0100-2264	1
9	Coupler/Shaft assembly		1
NS	Hex key, 3 mm	8710-0911	


## 18 Illustrated Parts Breakdown

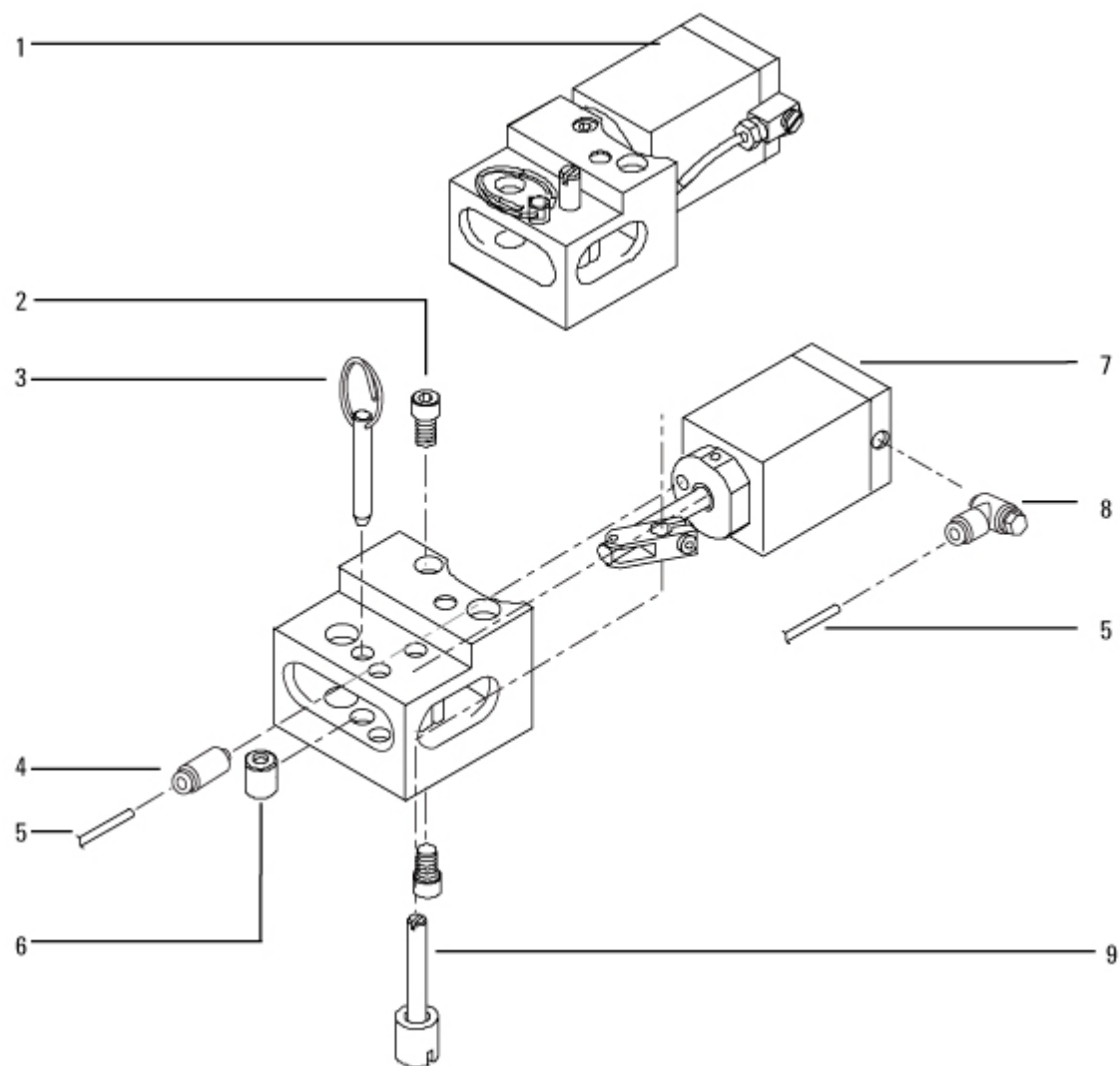

## Valve actuator assembly (2 of 2)

Item	Description	Part number	Qty
1	Cylinder end cap		1
2	O-ring, 1.176-inch id	0905-1405	1
3	O-ring, 1.046-inch id	0905-0463	1
4	Piston	19325-20640	2
5	Actuator cylinder		1
6	Hose fitting	0100-1707	1
7	PTFE tubing 5m, 1.5mm ID, 3mm OD	5062-2483	1
8	Piston rod		1
9	Dowel pin	1480-0017	1
10	Link	19325-80010	1
11	Screw, socket head, M4 x 8 mm		3
12	Elbow fitting	0100-1707	1

## 18 Illustrated Parts Breakdown

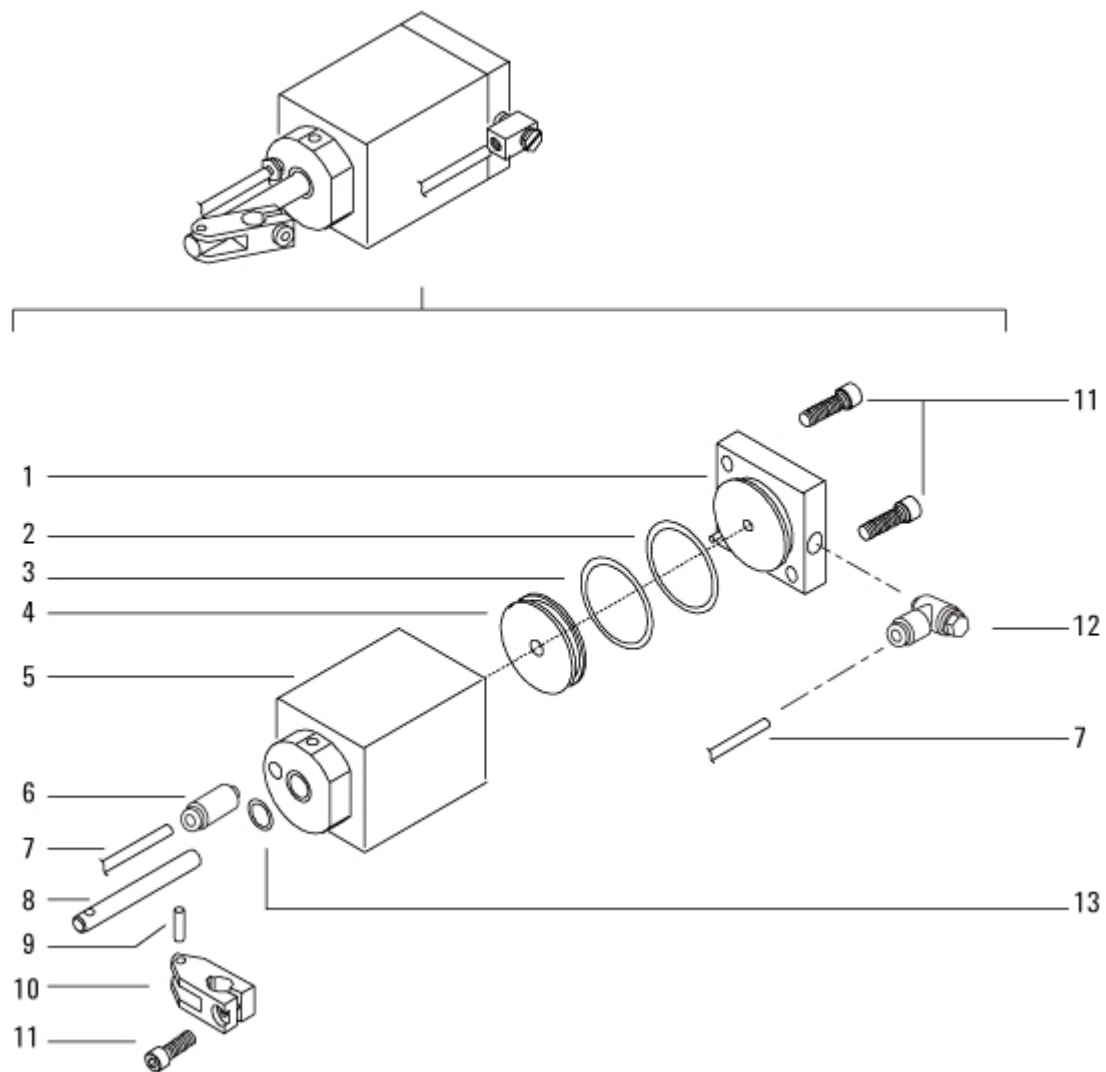

## Miscellaneous Valve Parts

The table below lists several valve accessory kits and parts that are useful when installing valves, valve boxes, and related parts.

Description	Comment	Part number
G1580A	Heated valve box for one valve	G1580A
G1581A	Heated valve box for two valves.	G1581A
Nut plate assembly		05890-80660
Capillary (SSL or MMI) inlet to upstream sampling valve kit	Kit to interface an upstream gas or liquid sampling valve to a capillary or Multi-Mode inlet for sample splitting and introduction into a capillary or megabore column.	G3480-67585
G2739A Liquid valve plumbing kit	This kit contains the hardware needed to install a liquid sample valve into a GC. It does not contain the valve or any automation.	G2739A
G2740A Gas valve plumbing kit	This kit contains the hardware needed to install any 6 port or 10 port valve into a GC and build any of the standard configurations. It does not contain the valve, automation, loops or needle valves.	G2740A
G2741A Needle valve kit	Adds an adjustable restrictor to control gas flows. Max temp of 225 °C. This kit is needed when the customer requires to have a column isolation configuration built into a GC.	G2741A
G2742A Downstream capillary port to valve kit	Kit to interface a split/splitless inlet or multi-mode inlet to a valve located downstream.	G2742A
G2743B 7890 Actuator kit	7890 Actuator Kit. Adds automation to 7890A/B valve system <ul style="list-style-type: none"> <li>• 7890A with S/N = CN10201096 or later and 7890A S/N = US10281001 or later</li> <li>• 7890B (all S/N's)</li> </ul>	G2743B
G2744A Heated zone kit	This kit is needed when adding a third heated valve to a GC that already has two heated valves installed. The heated zone will allow a 4th valve to be added if required.	G2744A
G2748A	Hardware kit to mount an unheated gas sampling or liquid sampling valve on the left side of the oven. Requires side carrier accessory (7890A: G3485A or 7890B: G3485B). Incompatible with MSD.	G2748A
G3478A Nickel Catalyst	Adds a nickel catalyst (methanizer) to a 7890 Gas Chromatograph. Requires back inlet position. Compatible with PCM in rear position.	G3478A
G3505A high pressure injection device accessory	Requires an SSL inlet and one flow source from an Aux EPC or PCM module. User must provide appropriate hardware to interface sample to HPID. One only. Front inlet only.	G3505A

## 18 Illustrated Parts Breakdown

Description	Comment	Part number
G4338A Inert 0.530 mm id capillary column to valve interface kit	Requires mounting plate 05890-80660. Each plate can hold 5 interface adaptors.	G4338A
G4339A Inert 0.320 mm id capillary column to valve interface kit	Requires mounting plate 05890-80660. Each plate can hold 5 interface adaptors.	G4339A
G2855-60200 Swaging wrench	Swaging Nut or Swaging Wrench is required to swage SilTite ferrule to column.	G2855-60200
G2855-20555 Swaging nut	Swaging Nut or Swaging Wrench is required to swage SilTite ferrule to column.	G2855-20555



## Electronics and Fans

This section contains illustrated parts breakdowns for the following 7890 GC electrical components.

- [Electronics carrier](#) (logic, A&P, and LAN boards)
- [Keyboard Assemblies](#)
- [AC circuit board components](#)
- [AC power](#)
- [Power cords](#)
- [Main transformer](#)
- [Analog input board \(G1556A accessory\)](#)
- [Chassis fans](#)

## Electronics carrier

The table below lists the 7890 Series GC electronics carrier parts. To identify the parts, see [Figure 23](#) for 7890B and [Figure 24](#) for 7890A.

Item	Description	Part number	Qty
NS	Fuse, ALS controller PCA, 7A 125V (ALS control PCA)	2110-0961	1
1	Analog & Power board	see <a href="#">Table 45</a>	1
2	Logic board	see <a href="#">Table 46</a>	1
3	Harness, 4-wire communication buss (LVDS communications harness)	G3430-60513	1
4	Harness, keyboard and display	G3430-60514	1
5	Backplane for Analog and power board	G3430-00114	1
6	Backplane for Detector signal boards	G3430-00085	1
7	Grounding screw, M4 × 25 mm, T-20 (for #1)	0515-2712	1
8	Screw, M4 x 6 mm	0515-2832	10
9	Screw, self-tapping	0515-4897	3
10	Nut, hex, with lock washer M4 x 3.2 mm thick	0535-0043	1
11	Grounding screw, M4 x 12 mm, T-20 (for #6)	0515-2496	1
12	ALS Controller Assembly	G3430-60529	1
13	Cable, RS-232,7683 ALS controller	G2612-60510	1
14	Backplane right logic board	G3430-00114	1
15	Screw, T-20, M4 × 12 mm (attach ALS board to carrier)	0515-2496	2


The tables below list the logic and A&P board changes for the 7890 series GCs. Entries in *italics* are obsolete.

**Table 45** Analog & Power board changes


Part number	Model	Min. firmware	Service Note	Features
<i>G3430-60150</i> (Obsolete)	<i>7890A</i>			<i>Original board. Power supply LEDs inside GC. LVDS FETs not overcurrent protected.</i>
<i>G3430-60151</i> (Obsolete)	<i>7890A</i>	<i>A.01.09</i>	<i>G3440A-025A</i>	<i>Power supply LEDs viewable on the back of the GC. LVDS FETs overcurrent protected.</i>
G3430-61050	7890B, 7890A			Modifications to improve the power line hold up time. If the lights blink in the lab (due to a few missed line cycles), the new board will continue to supply power to the system whereas the older board would have caused the GC to reboot.

**Table 46** Logic board changes

Part number	Model	Min. firmware	Service Note	Features
<i>G3430-60100</i> (Obsolete)	<i>7890A</i>			<i>Original board. LAN LEDs inside GC. Two battery switches.</i>
G3430-60101	7890A		G3440A-025A	LAN LEDs on back of GC. One battery switch. Used by OEMs that require old firmware.
G3430-61020	7890A	A.01.11	G3440A-058A	LAN LEDs on back of GC. One battery switch.
G3430-61110	7890B, 7890A+	B.02.00		Real Time Clock module can be replaced. Extra LVDS channels.

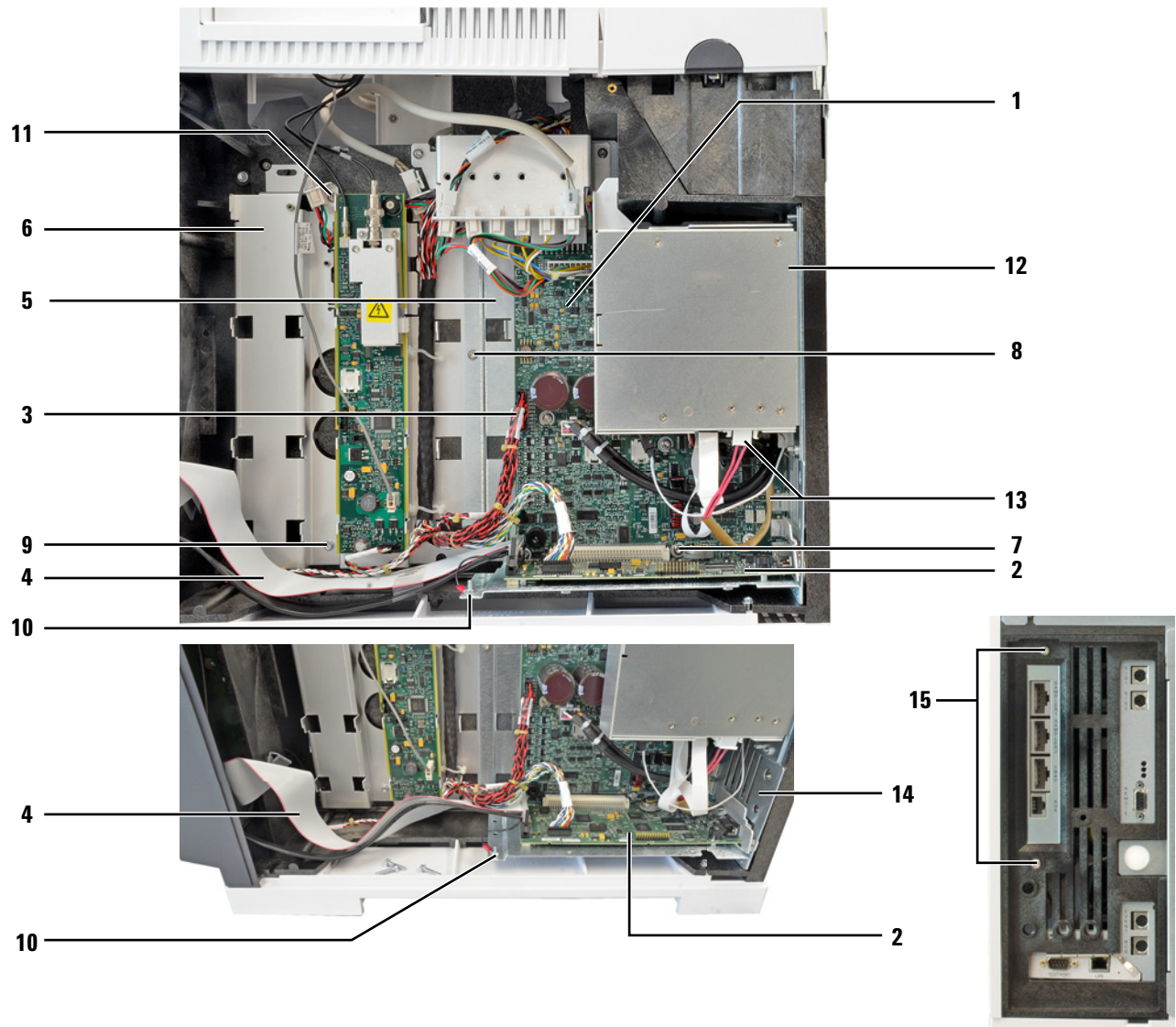

**Figure 23** 7890B Electronics carrier parts

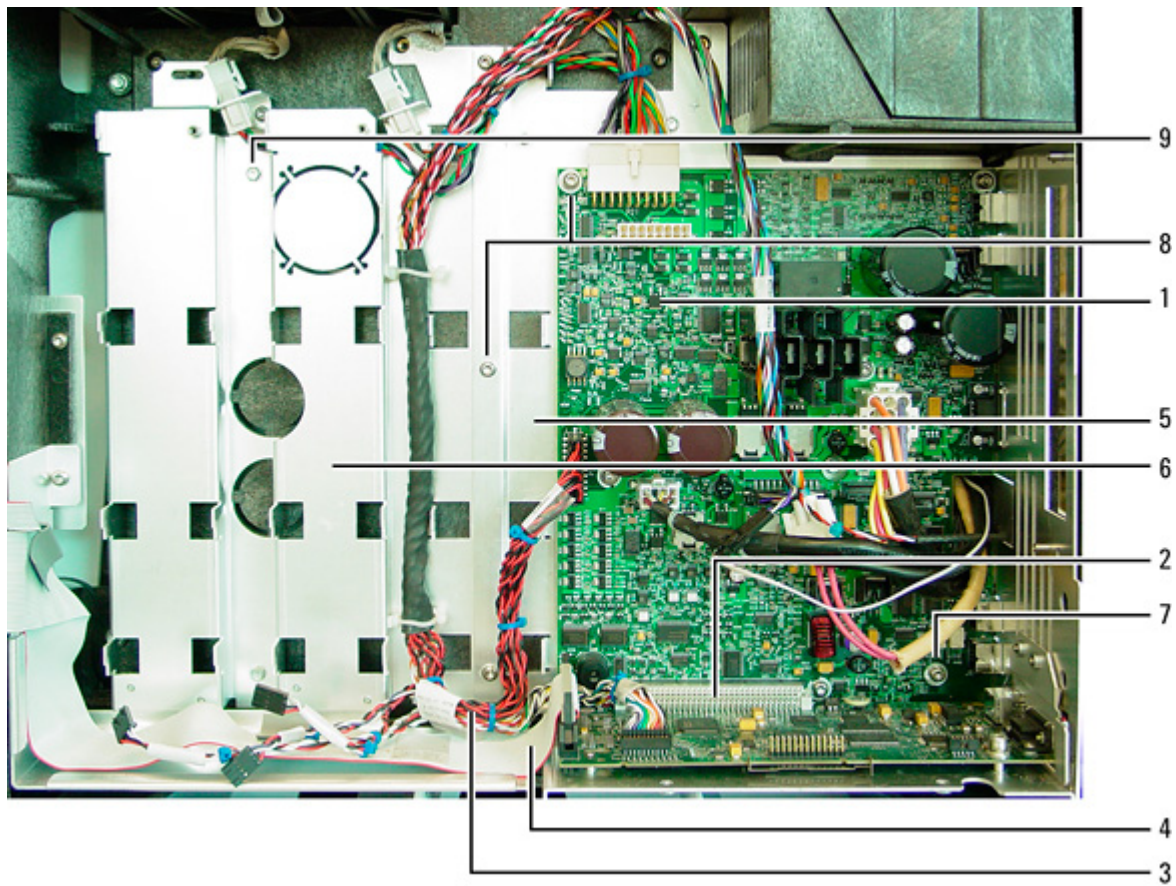


**Figure 24** 7890A Electronics carrier parts

## Keyboard Assemblies

### 7890A Keyboard and Display (original style)

Item	Description	Part number	Qty
1	Keyboard/display assembly with display and harness	G3430-60555	1
2	Replacement keyboard assembly without display	G3430-67500	1
3	Display	G3430-80003	1
4	Self-tapping screws	0624-1076	4
5	Keyboard/display harness	G3430-60514	1
6	Bracket	G3430-00098	1
7	Interconnect board	G3430-60003	1
8	Screw, M4 x 12 mm, T-20	1390-1023	1
NS	Screw, self-tapping, 1/4-inch (connect left and bottom sides of assembly to electronics carrier)	0515-4897	3
NS	Screw, M4 x 12 mm, T-20 (connects top right corner of assembly to electronics carrier)	0515-2496	1

Much of the keyboard/display is concealed by its plastic shell, so that a conventional exploded view is not very helpful. Use the schematic diagram below to identify parts.

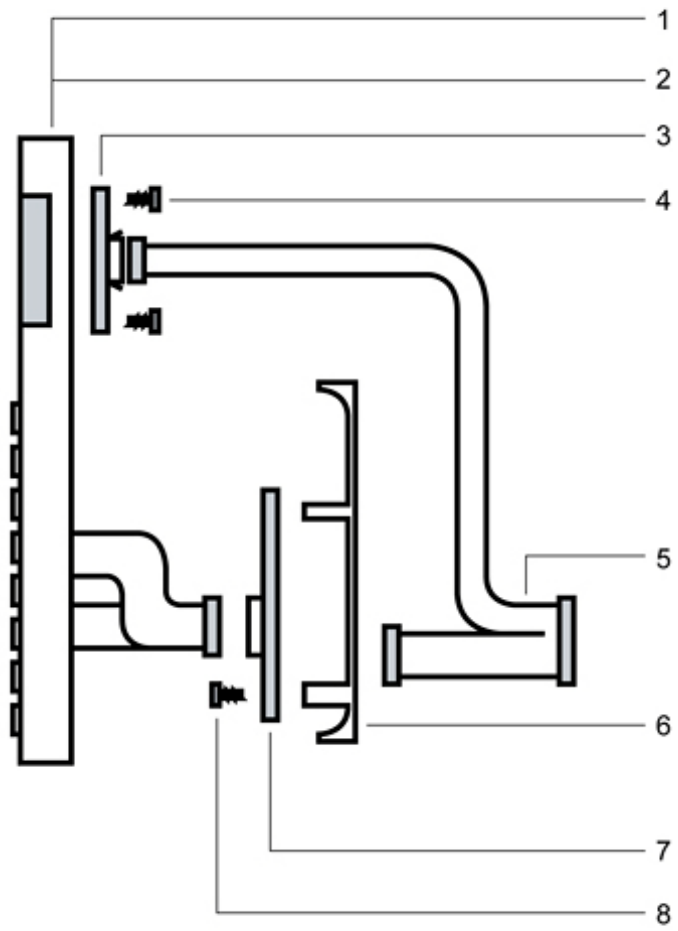


**7890A Keyboard and Display (new style)**

Item	Description	Part number	Qty
1	Keyboard/display assembly with display and harness (includes G3430-67500, G3430-60554, and G3430-80003 assemblies)	G3430-60555	1
2	Replacement keypad assembly without display, without keys	G3430-61080	1
NS	Replacement key and bezel assembly (no cable harness, no board, no display)	G3430-67500	1
3	Display	G3430-80003	1
4	Self-tapping screws	0624-1076	4
5	Keyboard/display harness	G3430-60554*	1
6	Door sensor interconnect cable	G3430-60561	1
7	Grounding nut, for display harness	0535-0043	1
8	Keypad support spacer	G3430-60022	8
9	Keypad support screw	G3430-20034	8
NS	Screw, self-tapping, 1/4-inch (connect left and bottom sides of assembly to electronics carrier)	0515-4897	3
NS	Screw, M4 x 12 mm, T-20 (connects top right corner of assembly to electronics carrier)	0515-2496	1

\*. Use for G3430-67500 and G3430-61080 keyboard assemblies.




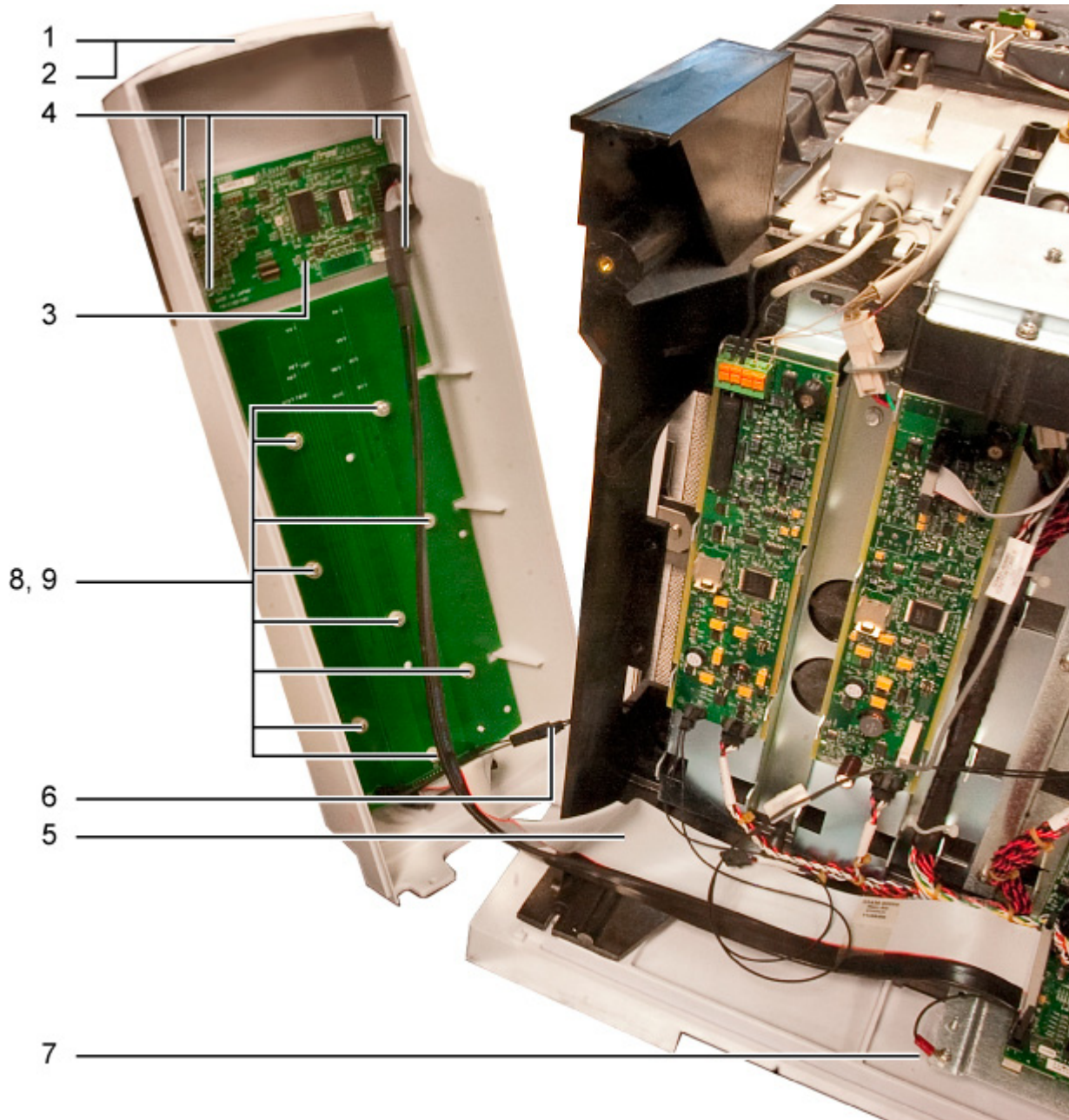

**7890B Keyboard and Display**


Item	Description	Part number	Qty
1	Keyboard Assembly 7890B	G3440-60555	1
2	Display	G3430-80003	1
3	Screw, M4 x 12 mm, T-20 (connects top right corner of assembly to electronics carrier)	0515-2496	1
4	Screw, self-tapping, 1/4-inch (connect left and bottom sides of assembly to electronics carrier)	0515-4897	3
5	Keyboard/display harness	G3430-60554*	1
6	Replacement key and bezel assembly (no cable harness, no board, no display)	G3430-67500	1
7	Self-tapping screws	0624-1076	4
8	Replacement keypad assembly without display, without keys	G3430-61080	1
9	Keypad support spacer	G3430-60022	8
9	Keypad support screw	G3430-20034	8
10	Door sensor interconnect cable	G3430-60561	1
NS	Grounding nut, for display harness	0535-0043	1

\*. Use for G3430-67500 and G3430-61080 keyboard assemblies.

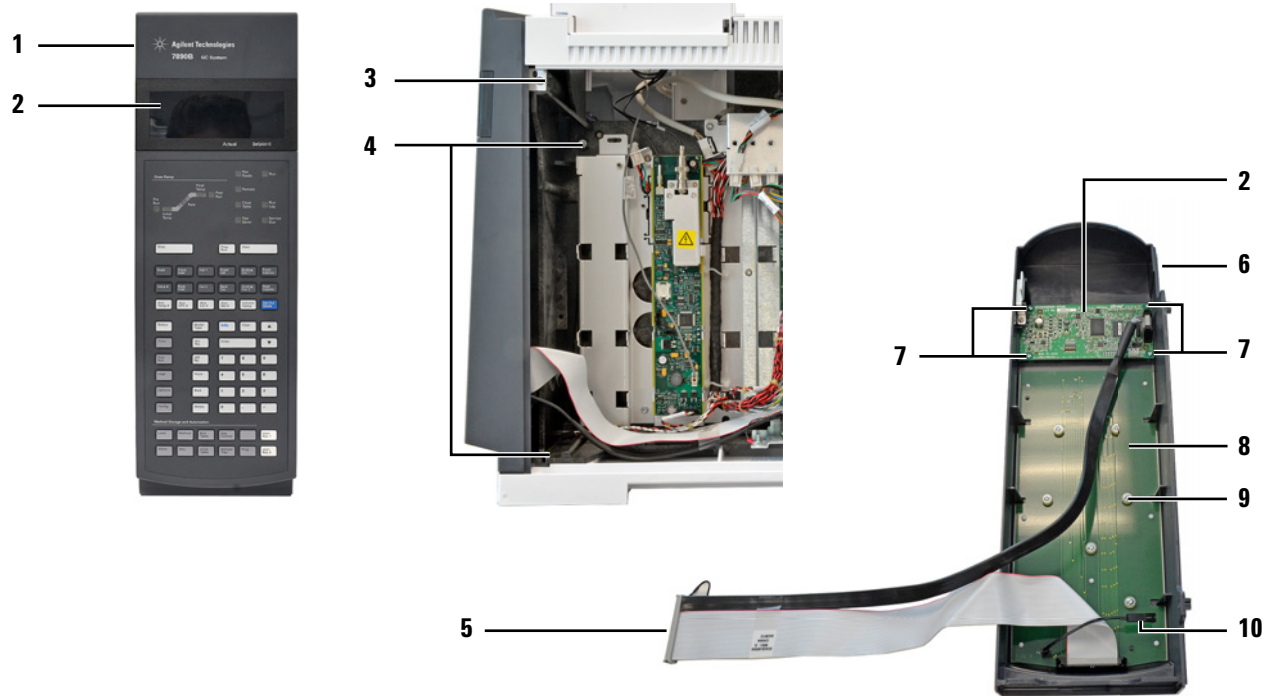


## AC circuit board components

Item	Description	Part number	Qty
1	Triac jumper (part of 5)		3
2	Screw, M4 x 20 mm long	1390-1024	1
3	Washer		1
4	AC configuration plug		1
	– USA and Canada 120V	G1530-60690	
	– Japan 200V	G1530-60700	
	– Continental Europe, 220V single phase power	G1530-60710	
	– Hong Kong, 220V, single phase	G1530-60710	
	– China 220V	G1530-60710	
	– China 220V, 10 amp slow ramp	G1530-60710	
	– Israel 220V	G1530-60710	
	– Switzerland/Denmark, 230V, 16A	G1530-60720	
	– Switzerland power cable, 230V, 10A	G1530-60720	
	– Denmark, 230V, 10A	G1530-60720	
	– Continental Europe, 230V	G1530-60720	
	– Australia, 240V	G1530-60730	
	– South Africa/India, 240V	G1530-60730	
	– United Kingdom/Ireland, 240V	G1530-60730	
	– USA and Canada, 240V	G1530-60730	
	– Australia 240V, slow ramp	G1530-60730	
5	AC board (without triac kit)	G3430-60050	1
6	Ceramic fuse, type x, 20A/250V	2110-0098	2
7	Glass fuse, type F, 8A/250V	2110-0036	2
8	Triac kit, includes triac (item 9) and thermal link pad (item 10)	G1530-61340	–
9	Triac	Order item 8	1
10	Pad (under triac)	1205-0934	1
11	Hex nut	0535-0031	2

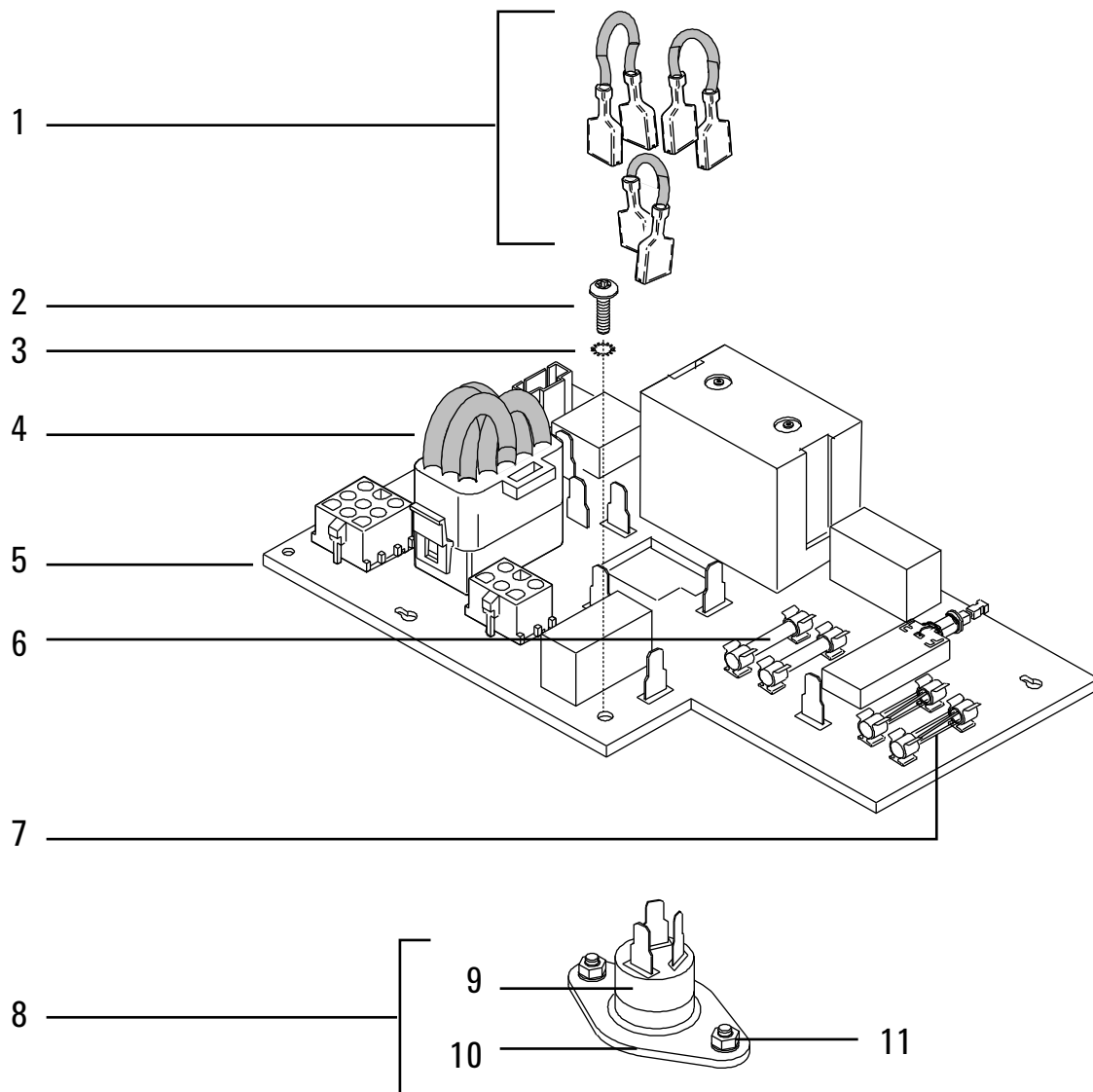

**AC power**

Item	Description	Part number	Qty
1	EPC pneumatics carrier assembly with fan	G3430-80503	1
NS	Screw, T-20, M4 × 12 mm (attach carrier to oven)	0515-2496	3
2	Fan, pneumatics carrier	G3430-60006	1
NS	Screw, M4 X 30 mm, T-20	0515-0669	4
3	Washer, steel to attach motor	2190-0712	6
3	Nut with lock washer, M4 to attach motor	0535-0043	3
4	Oven flapper assembly	G3430-60008	1
5	Motor, oven fan, DC	G3430-60504	1
6	Nut with lock washer, M4 X 3.2 mm for grounding straps (7 mm nut driver)	0535-0043	2
7	Transformer kit with large washers	G3430-60975	1
8	AC power board	G3430-60050	1
NS	Screw, M4, T-20	1390-1024	1
9	AC power inlet assembly (includes power receptacle)		1
<b>Regular oven</b>			
	120 VAC (requires a C19 power cord)	G1530-61550	
	220, 230, 240 VAC (requires a C13 power cord)	G1530-61560	
<b>Fast oven</b>			
	All (200/208/220/230/240 VAC) (requires a C19 power cord)	G1530-61550	

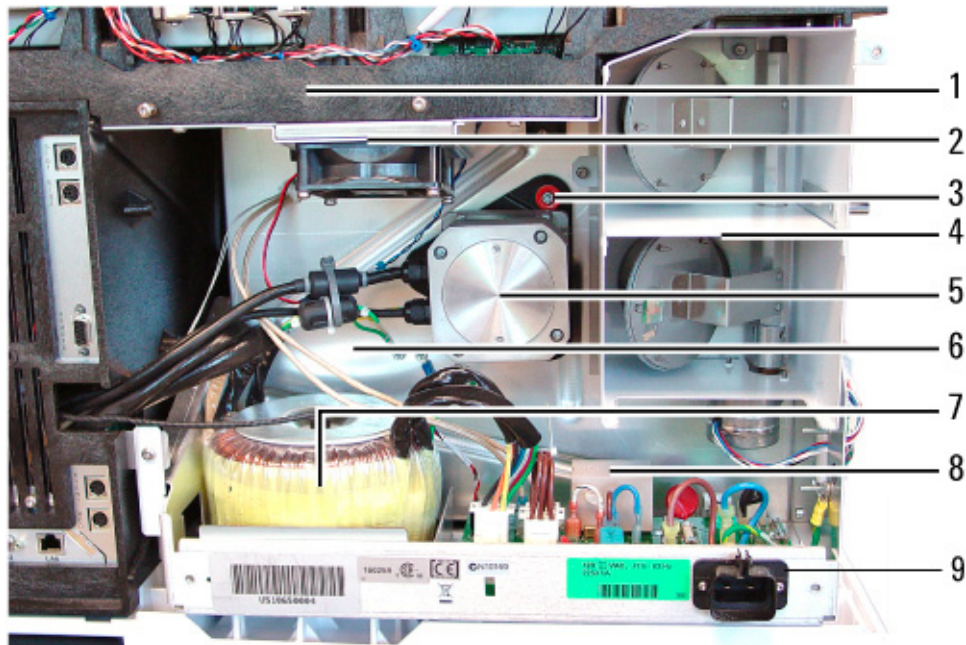

## Power cords



The 7890 Series GC uses detachable power cords.

Countries	VAC	Description	Part number
Australia, New Zealand		C13, 10 amp	8120-1369
Americas	120	10 amp	8120-1378
Europe		C13, 10 amp	8120-1689
US / Canada / Taiwan / Thailand		C13, 13 amp	8120-1992
Switzerland		C13, 10 amp	8120-2104
Denmark, Greenland		C13, 10 amp	8120-3997
India / South Africa		C13, 10 amp	8120-4211
Japan		10 amp	8120-4753
Israel		C13, 10 amp	8120-5182
Taiwan / South America		C19, 20A	8120-6360
Argentina			8120-6869
US	120, 208, 240	C19, 20 amp	8120-6894
Japan		C19, 20 amp	8120-6903
Chile		C13, 10 amp	8120-6978
Australia		C19, 16 amp	8120-8619
Great Britain / Hong Kong / Singapore / Myanmar		C19, 13 amp	8120-8620
Europe		C19, 16 amp	8120-8621
Switzerland, Denmark		C19, 16 amp	8120-8622
Great Britain / Hong Kong / Singapore / Myanmar		C13, 10 amp	8120-8705
China		C19, 15 amp, Fast	8121-0070
US	240	C19, 15 amp	8121-0075
Israel		C19, 16 amp	8121-0161
Argentina		C19, 20 amp	8121-0675
India / South Africa		C19, 15 Amp	8121-0710
China		C13, 10 amp	8121-0723
Korea		C19, 16 amp	8121-1222
Korea		C13, 10 amp	8121-1226
Thailand	220	15 A, 1.8M, C19	8121-1301



Main transformer

Part number G3430-60975

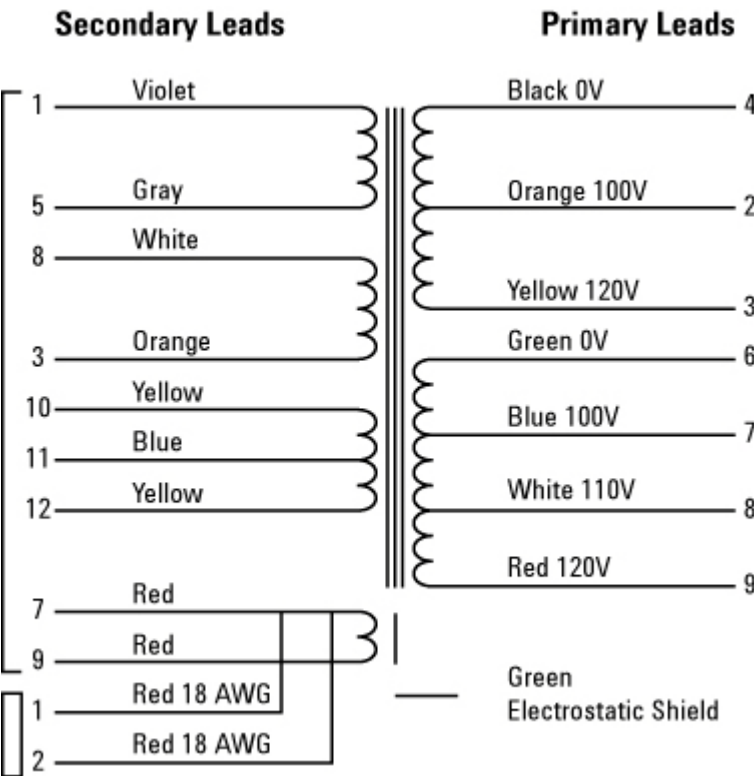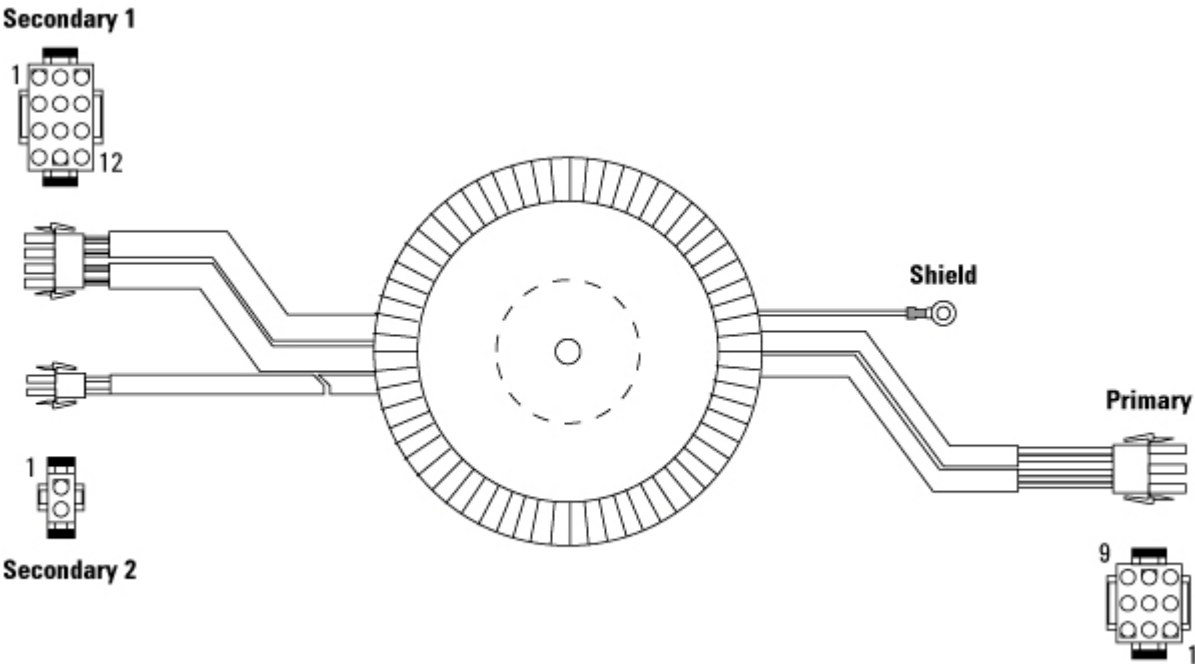


Note: Some transformers will have solid color wires as given; others may have black wires with color markers.

Connector (pin)	Wire color
Primary (1)	NC
Primary (2)	Orange
Primary (3)	Yellow
Primary (4)	Black
Primary (5)	NC
Primary (6)	Green
Primary (7)	Blue
Primary (8)	White
Primary (9)	Red
Secondary 1 (1)	Violet
Secondary 1 (2)	NC
Secondary 1 (3)	Orange
Secondary 1 (4)	NC
Secondary 1 (5)	Gray
Secondary 1 (6)	NC
Secondary 1 (7)	Red
Secondary 1 (8)	White
Secondary 1 (9)	Red
Secondary 1 (10)	Yellow
Secondary 1 (11)	Blue
Secondary 1 (12)	Yellow
Secondary 2 (1)	Red
Secondary 2 (2)	Red



Analog input board (G1556A accessory)

tem	Description	Part number	Qty
1	Analog input board	G3456-60010	1
2	General purpose analog output cable assembly	G1530-60560	1

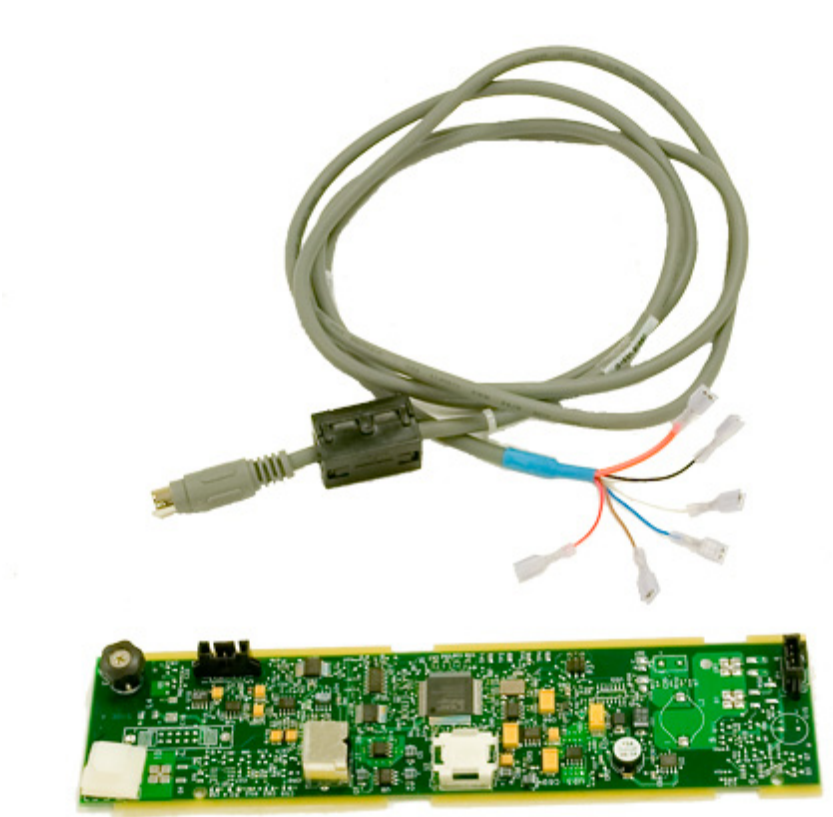

## Chassis fans

Item	Description	Part number	Qty
1	Inlet fan assembly	G3430-60590	1
2	Pneumatics area fan	G3430-60006	1

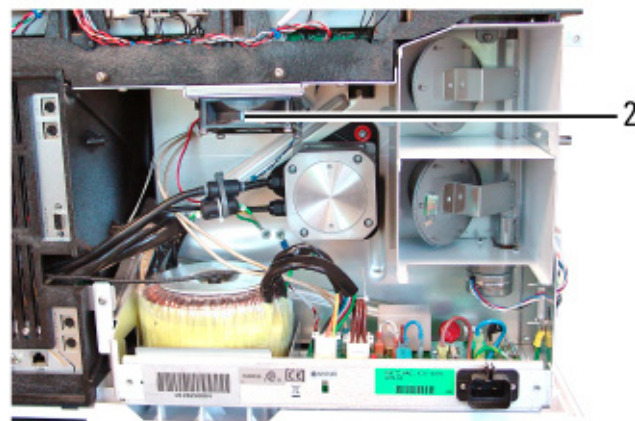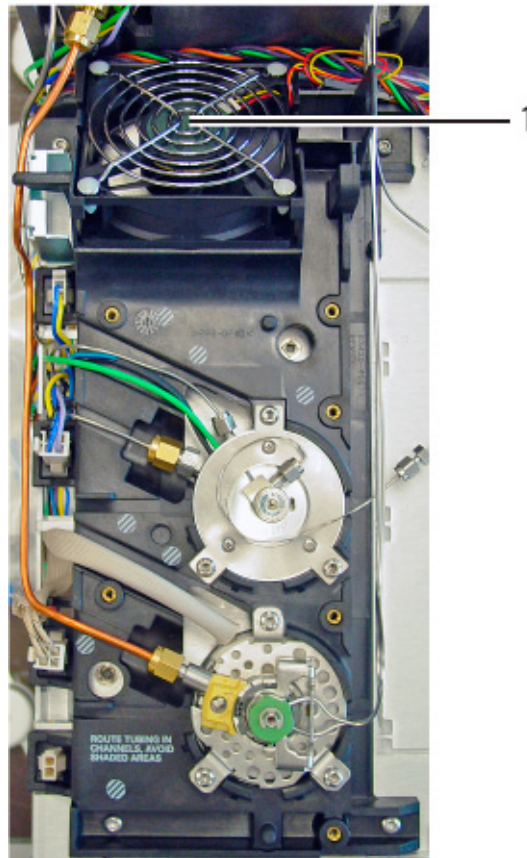

## Covers



Plastic covers	721
Back covers	723
Auto-injector mounting and parking posts	724
7890A Side Carrier	726
7890B Side Carrier	727
Mainframe	728

## Plastic covers

**NOTE**

7890A GC shown. Plastic covers listed are the same for 7890A and 7890B unless otherwise noted.


Item	Description	Part number	Qty
1	Pneumatics cover	G3430-60546	1
2	Detector cover	G3430-40007	1
NS	Detector cover for HS	G3430-60017	—
NS	Removable detector top cover for CTC	G3440-60013	—
3	Inlet cover (7 screws to attach to inlet carrier)	G3430-60540	1
4	Cover, left side with plugs	G3430-60541	1
5	Rotating insert	G3430-40024	1
6	1 1/2-inch hole plug	5040-4643	1
7	1-inch hole plug	5040-4642	2
8	Bezel, oven, black plastic, with 2 retained screws	G3430-60545	1
9	Electronics top cover (1 screw to attach)	G3430-60543	1
10	Cover, right side (1 screw to attach)	G3430-60542	1
NS	Screw, Captive M4 × 12 mm, T-20	1390-1023	5
NS	Screw, M4 × 12 mm, T-20	0515-2496	3
NS	Screw, M4, T-20	1390-1024	1
NS	Detector top cover repair kit	G3430-67648	—

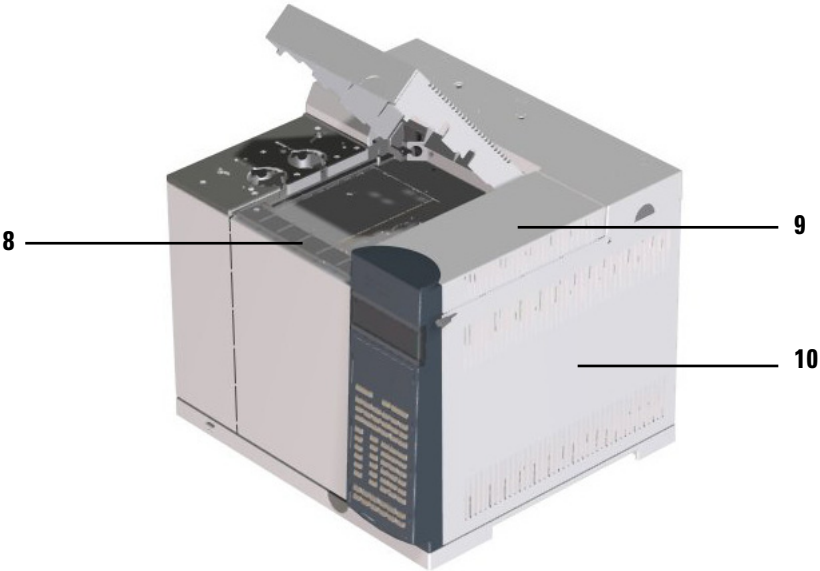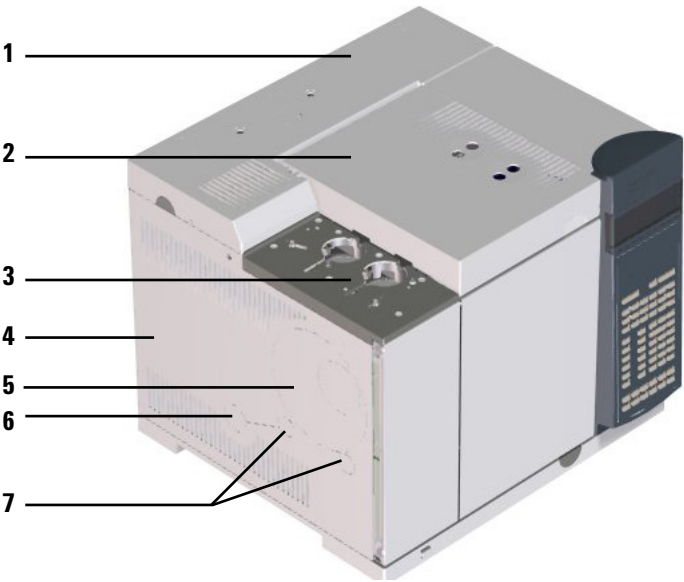

## Back covers

### Covers

Item	Description	Part number	Qty
1	Cover, back upper	G3440-01090	1
2	Oven exhaust deflector (optional)	G1530-80650	1
3	Cover, back lower	G3430-00023	1
4	Screw, M4 × 12 mm, T-20	0515-2496	6
5	Lock washer, No. 8, 0.168 in.		2

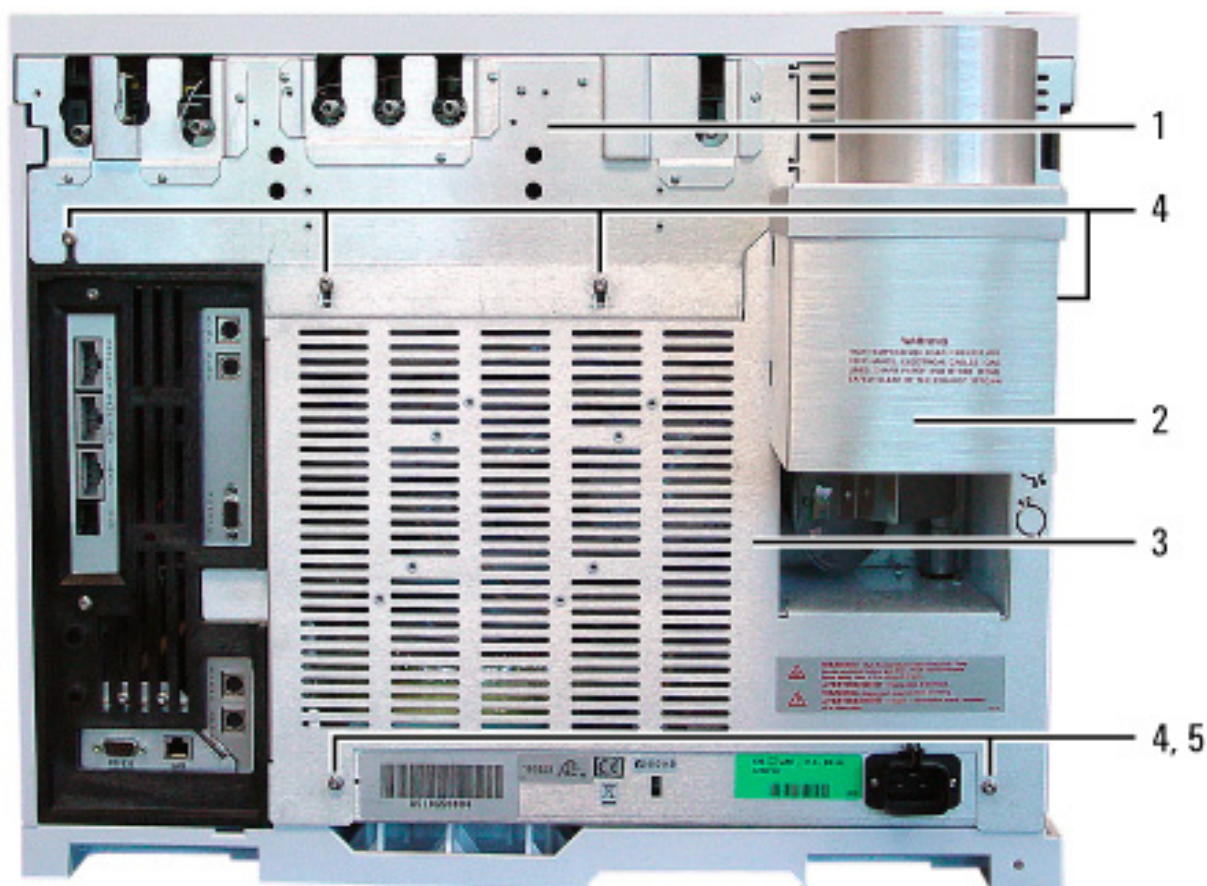




Auto-injector mounting and parking posts


7683B ALS

Item	Description	Part number	Qty
1	Injector mounting post	G2613-20500	2
2	Parking post, dual purpose	05890-61525	2

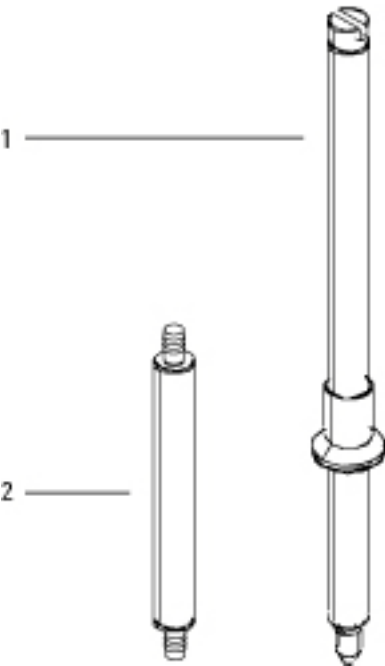

**7893A ALS, 7650 ALS**

Item	Description	Part number	Qty
1	Injector mounting post	G4513-20561	2
2	Parking post, dual purpose	05890-61525	2

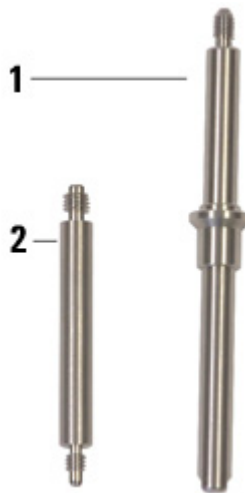**ALS Post Repair Kit**

If the threaded portion of the ALS post is broken off in the inlet chassis, order ALS Post repair Kit G3440-67001. The kit contains:


**Table 47** ALS Post Repair Kit contents

Description	Part number	Qty
Unique 7683 ALS mounting post	G2613-20505	2
Unique 7693 ALS mounting post	G4513-20574	2
Unique injection port cover	G3430-60022	1

The ALS posts in this kit are **not** compatible with the standard injection port cover. The unique injection port cover in this kit requires the unique ALS mounting posts.


7890A Side Carrier

	Description	Part number	Qty
1	Screw, M4 × 25 mm, Torx T-20	0515-2712	2
2	Pneumatic carrier assembly	G1530-60950	1
3	Screw, M4 × 12 mm, Torx T-20	0515-2496	4
4	Flow side cover	G1530-01240	1
5	Screw, M4 × 10 mm, Torx T-20, flathead	0515-2725	2

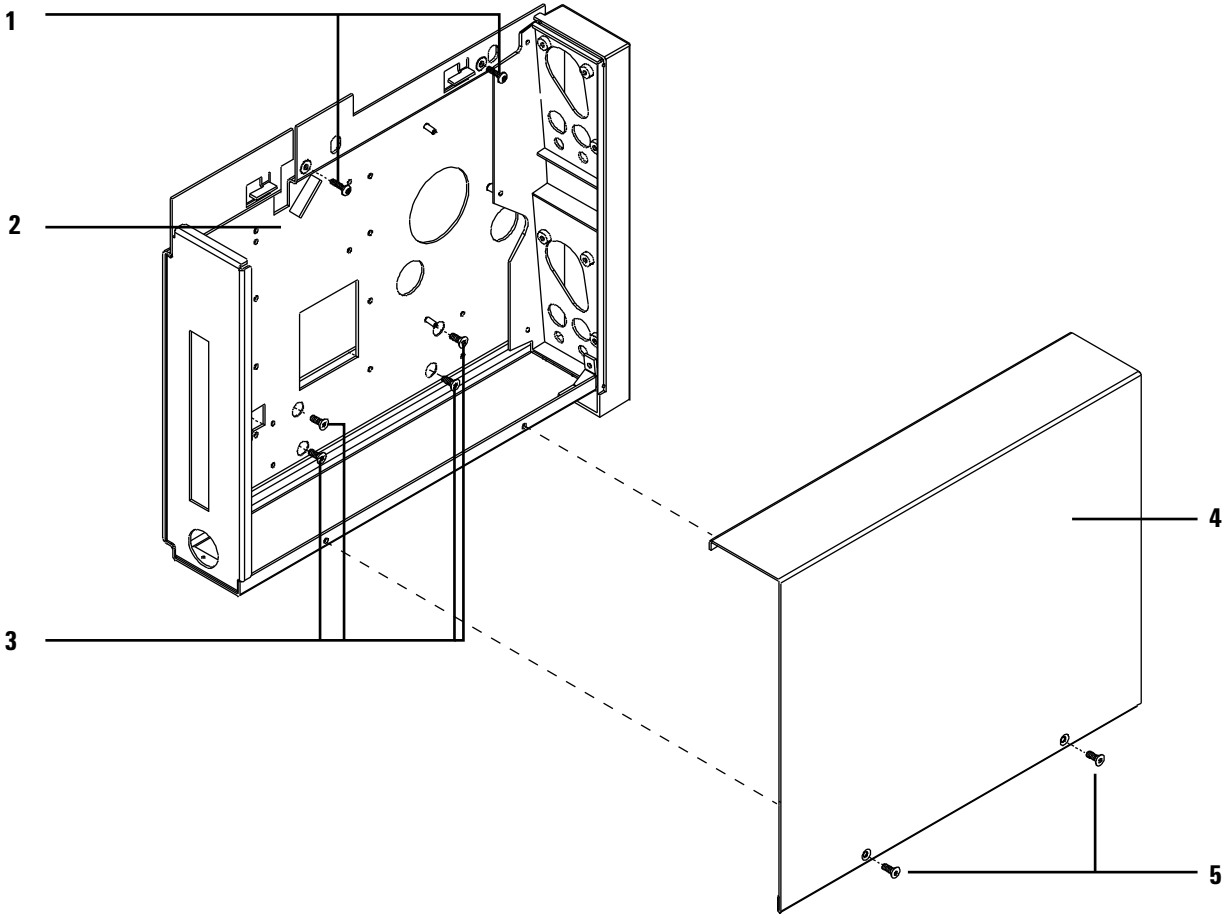

**7890B Side Carrier**


	Description	Part number	Qty
1	1/8-inch Bulkhead union, 316 SST	0100-0133	
2	aT-Manifold assembly	G1176-60569	
3	Screw, machine, M4 x 12 mm long	0515-2496	
4	Screw, machine, M4 x 25 mm long	0515-2712	
5	Bulkhead fitting retainer, 316 SST	0100-0137	
6	3rd Detector box cover assembly	G3432-60001	
7	Screw, machine, M4 x 10 mm long	0515-1269	
NS	Blank label plate	G1530-91705	

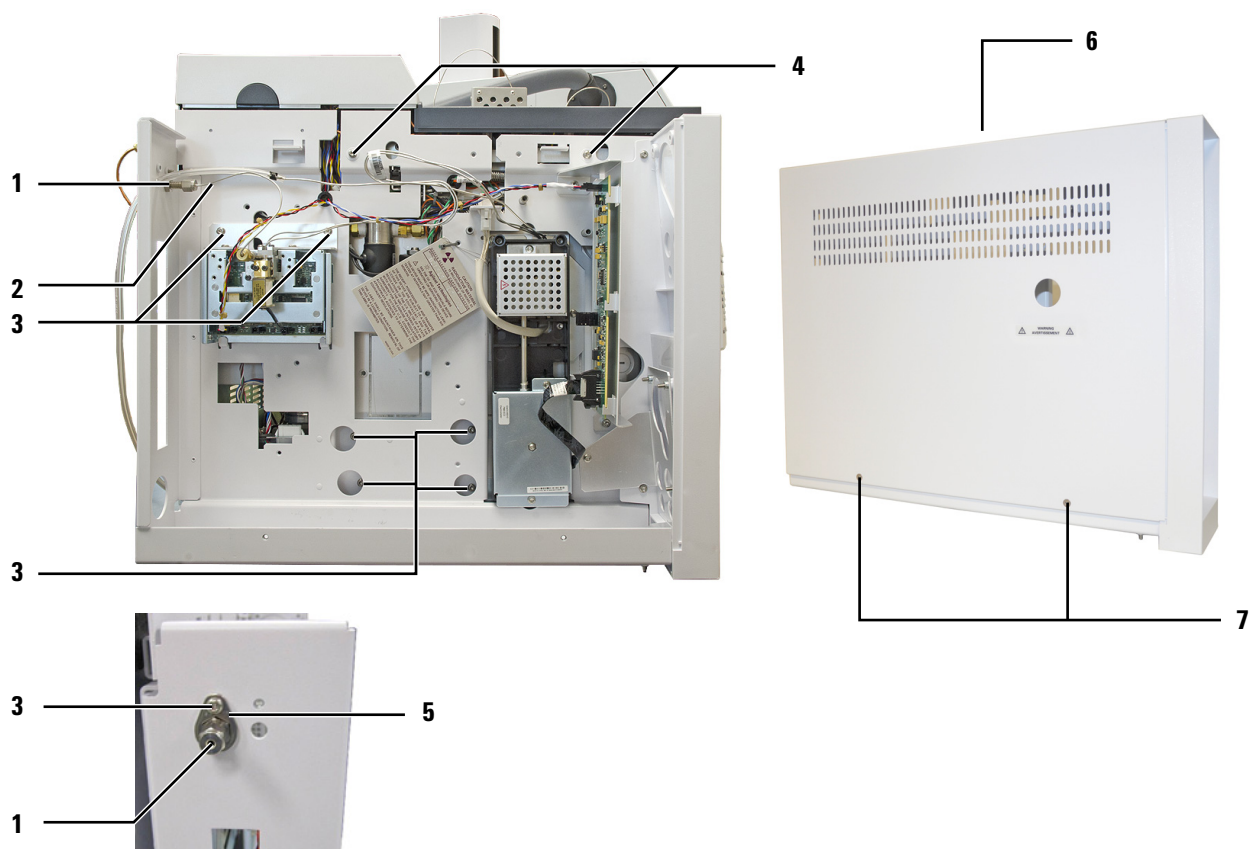

Mainframe

	Description	Part number	Qty
1	Inlet chassis, for SSL, PP, COC, and VI inlets	G1530-40020	1
2	PTV Inlet chassis, required for MMI and PTV. (Accepts all inlet types.)	G3430-40011	1

Standard inlet chassis

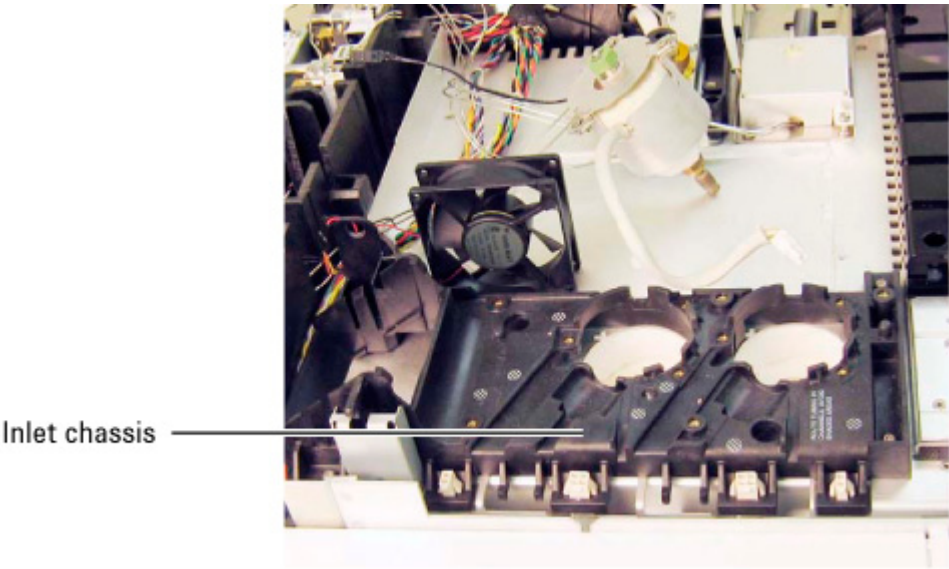

MMI and PTV inlet chassis



## Tools

This section contains parts not applicable to other sections.

### Shipping Crates

Description	Part number	Where used
MS shipping crate	G2589-63001	for 5973 MSD
MS shipping crate	G3170-80000	for 5975 MSD
MS shipping crate	G3870-80100	for 5977 MSD
GC shipping container, 7890	G3430-98000	for 7890 series GC
7683 A/B Tray package	G2614-60500	for 7683A/B Tray, G2614A
7683 A/B Injector package	G2913-60480	for 7683B injector

Note: MS shipping crates include the wooden pallet, cardboard carton, polyethylene bag, and protective foam. They do not include transfer line packing, straps, or packing for a rough pump.

## Hand tools

We assume you have a well-equipped toolbox, but some of the special tools listed here will be helpful.

**Table 48** 7890 Series GC tool bag, 5182-3456

Description	Part no.
Wrench, 1/4-inch to 5/16-inch	8710-0510
Wrench, open-end, 7/16-inch and 9/16-inch	8710-0803
Wrench, open-end, 1/2-inch and 7/16-inch	8710-0806
Wrench, open-end, 7/16-inch and 3/8-inch	8710-0972
Driver, nut, 7-mm	8710-1217
Cutter, tube	8710-1709
Screwdriver, Torx T-20	8710-1615
Screwdriver, Torx T-10	8710-1623
Tool bag	9222-3253
Tool, ferrule, .25–.32 mm, 2 x 3 x 3/4-inch	RFT-2500
Tool, ferrule, .53 mm	RFT-5300
Lens, magnifying	G2855-40001
Septum tool, knurled handle	450-1000



### Drivers

- T-20 Torx Driver — 5182-3465
- T-20 Torx Key (for close-quarters work) — 8710-1807
- T-10 Torx Driver — 5182-3466



- T-10 Torx Key (for close-quarters work) – 8710-2140
- #1 Pozidriv screwdriver – 8710-0899
- #2 Pozidriv screwdriver – 8710-0900
- Flat-blade screwdriver – 8730-0008
- 1/4-inch Nut Driver – 8710-1561 (FID Jets)
- 7-mm Nut Driver – 8710-1217

### **Wrenches**

- Wrench, 1/4-inch to 5/16-inch – 8710-0510
- Wrench, open-end, 7/16-inch and 9/16-inch – 8710-0803
- Wrench, open-end, 1/2-inch and 7/16-inch – 8710-0806
- Wrench, open-end, 7/16-inch and 3/8-inch – 8710-0972
- Wrench, adjustable, 12-inch – 8710-1712
- Wrench, open-end, 3/16-inch x 1/4-inch – 8710-2618
- Wrench, open-end, 9/16-inch x 5/8-inch – 8720-0010
- Wrench, 7-mm x 8-mm – 1340407010
- Wrench, 6-mm x 7-mm – 1340407011
- Wrench, 8-mm x 10-mm – 1340407012
- Wrench, 6-mm x 6-mm – 8710-2156
- Wrench, 5-mm x 5-mm – 8710-2157

### **Cutters/Crimpers/Strippers**

- Diagonal Cutters
- Crimper/Wire Strippers
- Tubing Cutters
  - Restek #20193 <http://www.restekcorp.com>
  - Agilent 8710-1709
  - Agilent precision tubing cutter for 1/16-in. SS tubing 5190-1442
  - Alternate for Europe:

HICHROM Ltd - Part No: HI-196

Description: Tubing Cutters - Price: £17.00

1 The Markham Centre Station Road – Theale, Berkshire

RG7 4PE, UK Telephone: +44 (0) 118 930 3660

[mailto: salesw@hichrom.co.uk](mailto:salesw@hichrom.co.uk)

**Pliers**

- Channel lock pliers
- Slip-joint pliers
- Needle nose pliers
- Agilent needle nose pliers – 8710-0004
- Agilent Truarc pliers – 8710-0018

## Specialized tools for GC support



- FID Flow Tool – 19301-60660
- FID Cleaning Kit – 9301-0985, (.010-inch Guitar String)
- Injection Port Cleaning Kit – 480-0003
- 12 Piece File Kit – RSF-1200

- (Ferrule Removal tool, not shown – 440-1000)
- Fiberglass Tape – 0460-0186
- PTFE Tape – 0460-0016
- White Gloves – 8650-0030



- NPT Adapter – 0100-0118
- 5/16-inch by #20 Thread Chaser (Die), to clean 1/8-inch Swagelok threads – Other Supplier
- On/Off Valve, 1/8-inch Ball Valve – 0100-2144
- Exacto Knife
- Metric (8710-0641) and English Allen Wrench – Hex Keys
- Pin Vise and small drill bit
- Tweezers – 8710-0007
- Micro Probes (Sharp Object) – RMP-5005 – Qty 5
- Inspection Mirror – 707-0027
- 1.25-inch socket and ratchet for removing FID brass retaining nut
- Column cutter wafer, 4/pk – 5181-8836
- 6890 inlet wrench – 19251-00100
- 7890 inlet wrench (for SSL and MMI) – G3452-20512




## Electronic tools for GC support




- Digital Multimeter – (Agilent U1231A True RMS shown)
- Power Outlet Test Tool (Radio Shack/Sears/Electrical Supply)
- Static Strap – 9300-0969 (small) or 9300-0970 (large), 5 ft cord/clip – 9300-0980
- Wrist strap, disposable, 4-LG 1-W – 9300-1408
- Various jumpers and clip leads – electronic supply
- ESD mat – 9300-1484
- Useful Cables – not shown:
  - Crossover LAN Cable – 5183-4649
  - 9 Pin RS-232 Null Modem Cable – G1530-60600
  - USB/RS-232 adapter cable – 8121-1013
- Sensor pin extractor – 8710-1542

## 18 Illustrated Parts Breakdown

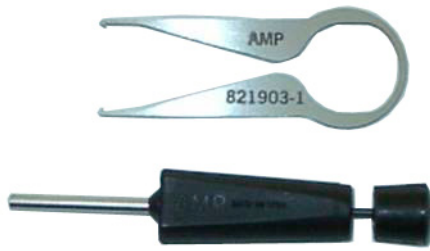

## Pneumatic tools for GC support

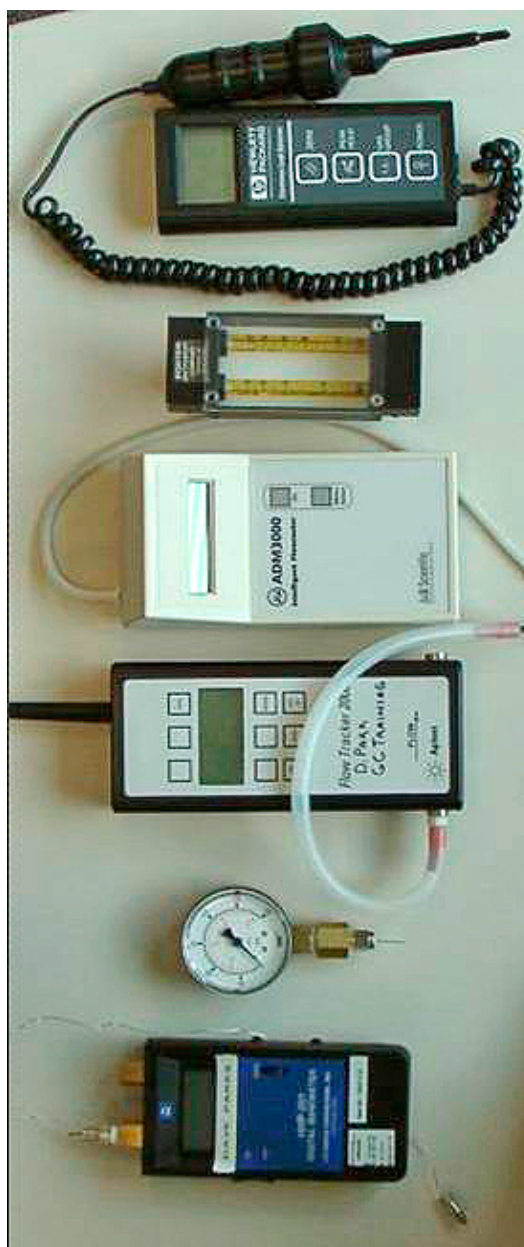



- Electronic Leak Detector
  - G3388B (or G3388A)
  - 5182-9646, 120 VAC
  - 5182-9648, 220 VAC (Agilent)
- Rotameter, Porter Model 65
 

Tube 1: 10-70 cc/min He, Tube 2: 85-850 cc/min Air

[http://www.porterinstrument.com/indust\\_21.php](http://www.porterinstrument.com/indust_21.php)
- Electronic Flow Meter
  - Preferred flow meter: Precision gas flow meter, flow range 5-500 mL/min – 5067-0223



- ADM 1000 – 220-1170
- ADM 2000 – 220-1171-U (Mass Flow Version)
- Electronic "Mass" Flow Meter – Flow Tracker
  - Model 1000 Flow Only – 5183-4779
  - Model 2000 with Leak Detector – 5183-4780
  - Universal Power Adapter – 5183-4781
- White Silicone Tubing:
  - 4 meters – 701-0016
  - 7 cm – 220-1179
- PTFE tubing, 1/8-inch, 1 m long – G1290-80050



### For Checking Inlet Pressure

- Quality Analog Pressure Gauge (0-100 or 0-60 PSIG) adapted with a Headspace Probe (301-013-HSP) – need NPT adapter – 0100-0118 and TCD Ferrules Set – 5180-4103, 5182-9673, 5182-3477
- Digital Pressure Gauge (0-100 PSIG) – Omega – Model HHP-201 Adapted with 1/8-inch NPT adapter 0100-0111, SS Capillary 5021-1831, Reducing Ferrule – 0100-1342 and 530  $\mu$ M Column Ferrule – 5080-8773



NPD Flow adapter - G1534-60640



G1530-20610 - FID/NPD Jet Plug



"No-Hole" Column Ferrules – 5181-7458



FID Flow adapter - 19301-60660



5060-9055 uECD/TCD Detector Plug



Plastic 1/8-inch Swagelok Cap – P/N 0100-2414 Qty 3

Replacement Septum Purge Fitting for EPC Module – P/N  
G2630-61720 (Not included in kit)



Blanking plug, finger tight, 5020-8294. (Can be used with any ferrule to duplicate using a no-hole ferrule.)

